



VistA Blood Establishment Computer Software (VBECS) Version 2.1.0

User Guide

April 2016

Department of Veterans Affairs
Product Development

This page intentionally left blank.

Revision History

Date	Revision	Description	Author
6/3/15	1.0	<p>Modified VistA Blood Establishment Computer Software (VBECS) User Guide for VBECS 2.0.0, Version 7.0 to create the VistA Blood Establishment Computer Software (VBECS) User Guide for VBECS 2.1.0, Version 1.0:</p> <ul style="list-style-type: none"> • Global: Replaced "2.0.0" with "2.1.0" • Global: Replaced "April 2015" with "June 2015". • Table 1: Removed last row "VBECS System Administrator (not a VBECS user)". • Testing: Revised seventh bullet. • Hardware and Infrastructure Architecture section: Changed "vSphere 5" to "vSphere 5.5". Changed "AlwaysOn cluster" to "Windows Server Failover Clustering with AlwaysOn cluster". • Commonly Used System Rules, Date and Time, first bullet: Changed "System Administrator" to "National Help Desk". • Configuring Site Parameters, Assumptions second bullet and Additional Information fifth bullet: Changed "VistA System Administrator" to "VBECS Application Administrator". • Abo/Rh Confirmation: Revised first sentence. • ABO/Rh Confirmation, Limitations and Restrictions: Revised second bullet to read "When building a worklist by invoice number, only invoiced units not previously tested and confirmed appear on the worklist. The user may deselect or add additional units." • ABO/Rh Confirmation, Step 4: Revised last note to read "Upon successful completion of ABO/Rh confirmation, the unit status is updated. Successful completion means that the ABO and Rh (if required) interpretation is consistent with recorded observed results and matches the ABO/Rh at unit login and any repeat ABO/Rh tests on the unit that are not invalidated." • Return Issued Units To Blood Bank, Step 7: Removed first item under VBECS column and revised second note. • Discard or Quarantine section, Step 6: Revised eighth note. • Discard or Quarantine by Invoice section, Step 4: Revised eighth note. • Edit Unit Information Section, Step 3: Revised third note. • Table 13: Added row for Serologic Crossmatch Test Required. • Table 19: Removed row "Unit out of controlled storage for more than 30 minutes found acceptable". • Appendix D: Revised first item under ABO-Rh Confirmation. Moved "Users may perform testing for no more than four patients in one VBECS session." to Patient Testing, Pending Task List. 	BBM team
11/30/15	2.0	<ul style="list-style-type: none"> • Global: Replaced "June 2015" with "December 2015". • Global: Removed notes in the reports which read "VBECS prints scheduled reports once an hour.", because VBECS checks every 5 seconds and prints when scheduled. • Table of Contents: Under Appendix B area, removed "Table" and "Figure" wording. • Customer Support section: Revised to mirror VBECS 2.1.0 Release Notes, Version 1.0. • Configuring Site Parameters, Additional Information: Added new seventh bullet to see VBECS Technical Manual-Security Guide for Active Directory tool installation. 	BBM team
4/4/16	3.0	<ul style="list-style-type: none"> • Global: Replaced "December 2015" with "April 2016". 	BBM team

This page intentionally left blank.

Table of Contents

REVISION HISTORY	3
INTRODUCTION.....	1
FREEWARE DISCLAIMER.....	1
RELATED MANUALS AND MATERIALS	2
ORIENTATION.....	3
HOW THIS USER GUIDE IS ORGANIZED.....	3
Terms.....	3
Background and Basic Knowledge	3
Figures and Tables	3
Options	3
Supplemental Information.....	4
Appendices	4
SECURITY.....	4
APPLICATION ARCHITECTURE.....	6
HARDWARE AND INFRASTRUCTURE ARCHITECTURE.....	6
Options That Require VistALink	7
CUSTOMER SUPPORT	10
Problems?	10
USING THE SOFTWARE.....	10
SOFTWARE BASICS	10
COMMONLY USED SYSTEM RULES.....	12
Date and Time	12
Locking and Time-Outs.....	12
Working with Data	12
Searching the Database	13
Testing.....	14
Creating and Viewing Reports	16
REPORT TOOLBAR	18
Navigation Arrows and Report Pages	18
Report Layout.....	18
Report Zoom-in/Zoom-out	19
Report Search	19
Report Export	20
Copying an Exported Report.....	23
USING VBECS	25
USING ONLINE HELP.....	26
Contents Tab	28
Index Tab.....	29
Search Tab.....	30
Glossary Tab	31
VIEWING THE VBECS VERSION NUMBER.....	32

VISTA RECORDS IN VBECS	34
ACCEESSING THE SYSTEM	35
Log Into VBECS and Vista.....	35
Exit VBECS	38
Screen Settings.....	39
Screen Display.....	42
Sound.....	43
Volume	45
Keyboard Shortcuts for Windows	46
Drop-Down Menus.....	47
Tool Tips	47
Icons and Buttons	48
Toolbars.....	48
PATIENT INFORMATION TOOLBAR	49
Patient Blood Availability.....	49
<i>Medication Profile</i>	50
<i>View Special Instructions and Transfusion Requirements</i>	51
View Recent Orders, Recent Transfusions/Issued Units, and Transfusion Reaction History	53
<i>Patient Information Toolbar: Recent Orders</i>	53
<i>Patient Information Toolbar: Recent Transfusions/Issued Units</i>	54
<i>Patient Information Toolbar: Transfusion Reaction History</i>	54
OTHER VBECS FUNCTIONS	55
Collecting Workload	55
CONFIGURING SITE PARAMETERS	57
UPDATE USER ROLES	57
COMPONENT CLASSES	59
CONFIGURE DAILY QC	62
CONFIGURE DIVISION	65
Configure Division.....	65
Configure Testing.....	67
Product Modifications	69
Order Alerts.....	71
Login Message	73
LOCAL FACILITIES	75
Blood Products	78
Antibodies	81
Canned Comments	83
SETTING TRANSFUSION PARAMETERS.....	85
Transfusion Complications.....	85
Transfusion Effectiveness	86
MSBOS	88
Workload Codes.....	90
REAGENTS AND SUPPLIES.....	92
Enter Daily QC Results	92
REAGENTS	99
Supplies	105
Equipment	109

COMPONENT PROCESSING	115
SHIPMENTS	115
Incoming Shipment	115
Edit Invoice Text	133
UNIT SEARCH SCREEN	134
TEST UNITS	136
ABO/Rh Confirmation	136
Unit Antigen Typing	139
MODIFY COMPONENTS	147
Modify Units	147
Modify Units: Split a Unit.....	157
Modify Units: Pool Units	160
Add and/or Remove Units from a Pool	164
MAINTAIN UNIT RECORDS	167
Edit Unit Information	167
Edit Financial Data.....	172
Free Directed Unit For Crossover	174
Discard or Quarantine	176
SELECT UNITS	182
ISSUE BLOOD COMPONENTS.....	192
RETURN ISSUED UNITS TO BLOOD BANK	204
RELEASE UNITS FROM PATIENT ASSIGNMENT.....	207
PRINT UNIT CAUTION TAG & TRANSFUSION RECORD FORM	209
INVALIDATE TEST RESULTS	217
REMOVE FINAL STATUS	221
PROCESSING ORDERS	225
ACCEPT ORDERS: ACCEPT AN ORDER	225
Accept Orders: Pending Order List.....	230
Display Patient Update Alerts	233
Display Patient Merge Alerts	234
MAINTAIN SPECIMEN.....	235
PATIENT TESTING: CANCEL AN ACTIVE ORDER	239
PATIENT TESTING: PENDING TASK LIST.....	241
Accept Orders: Cancel a Pending Order	243
ORDER REFLEX TESTS	245
PATIENT SERVICES.....	249
PATIENT SEARCH SCREEN	249
Patient Merge	251
Updated Patients and Deceased Patients	255
Patient Testing: General Instructions	256
Patient Testing: Record a Patient ABO/Rh	261
Patient Testing: Record a Patient Antibody Screen	264
Patient Testing: Record a Direct Antiglobulin Test	266
Patient Testing: Record a Crossmatch.....	269
Patient Testing: Record a Patient Antigen Typing	274
Patient Testing: Enter Antibody Identification Results.....	280

Patient Testing: Record a Transfusion Reaction Workup	283
Finalize/Print TRW	291
MAINTAIN PATIENT RECORDS.....	295
Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement.....	295
Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions.....	300
POST-TRANSFUSION INFORMATION	303
DOCUMENT ABO INCOMPATIBLE TRANSFUSIONS.....	308
Justify ABO/Rh Change.....	310
REPORTS.....	313
ADMINISTRATIVE DATA REPORT	313
AUDIT TRAIL	315
BLOOD AVAILABILITY	317
COST ACCOUNTING REPORT	319
C:T RATIO REPORT.....	321
DIVISION TRANSFUSION REPORT	323
DIVISION WORKLOAD REPORT.....	325
EXCEPTION REPORT.....	327
FINALIZE/PRINT TRW (SEE PATIENT SERVICES, FINALIZE/PRINT TRW).....	329
INAPPROPRIATE TRANSFUSION REQUEST REPORT.....	329
ISSUED/RETURNED UNIT REPORT	331
MEDICATION PROFILE (SEE PATIENT INFORMATION TOOLBAR)	333
ORDER HISTORY REPORT	333
PATIENT HISTORY REPORT.....	335
PROLONGED TRANSFUSION TIME REPORT	337
PATIENT TESTING WORKLIST AND TESTING WORKLIST REPORTS	339
TRANSFUSION COMPLICATIONS REPORT.....	342
TRANSFUSION EFFECTIVENESS REPORT	344
TRANSFUSION REACTION COUNT REPORT	346
TRANSFUSION REQUIREMENTS REPORT.....	348
UNIT HISTORY REPORT	350
REFERENCES.....	353
GLOSSARY.....	357
APPENDICES	367
APPENDIX A: DOWNTIME FORMS AND INSTRUCTIONS	367
Patient Testing Form	369
Unit ABO/Rh Confirmation Form	372
Unit Issue and Inspection Log.....	373
Unit Modification Form	374
APPENDIX B: DATABASE TABLE INFORMATION	377
Antibody and Antigen Tables.....	377
Antibody and Antigen Table: Irregular Antibodies.....	377
Antibody and Antigen Table: Database Conversion	380
Antibody Screen Test Interpretation	381
Canned Comment Category Types and Text.....	382
CPRS Orderable Blood Components (Component Classes) Mapped to ICCBBA Component Classes	390

Details in Audit Trail Report.....	391
Details in Exception Report	396
Enable Crossmatch Option for OTHER Product Types.....	398
Allowable Product Modifications by Original Product Type.....	399
National Treating Specialty Table.....	404
Order Status Flowchart.....	407
Rules for Electronic and Serologic Crossmatch.....	408
Unit Status Flowchart.....	409
VBECS Maximum Surgical Blood Order Schedule	411
BCE COTS Message Updates.....	412
Workload Process to Use Case Mapping Table	416
APPENDIX C: VBECS WORK FLOW	427
APPENDIX D: LIMITATIONS AND RESTRICTIONS	431
APPENDIX E: KNOWN DEFECTS AND ANOMALIES	445
APPENDIX F: SYSTEM VALIDATION.....	447
Introduction	447
Validation Plan.....	447
Prerequisites	448
Record System Settings.....	449
Configure Division.....	449
Configure Testing: Full Service	449
Product Modifications: Full Service.....	450
Order Alerts.....	451
Configure Daily QC: Full Service.....	452
Maintain Minimum Levels: Full Service	453
Log In Reagents: Full Service	454
Blood Products	455
Local Facilities	456
Component Classes	457
Transfusion Complications.....	458
Transfusion Effectiveness	459
Workload Codes.....	460
Canned Comments	460
MSBOS	460
Edit Invoice Text.....	460
Antibodies	461
Update User Roles.....	462
Testing.....	462
VBECS Forms and Reports.....	463
APPENDIX G: WORK PROCESS FLOWCHARTS	465
Configure Division Settings: Full-Service Divisions	465
Configure Division Settings: Transfusion-Only Divisions	466
Maintain Reagents and Daily QC: Full-Service Divisions.....	467
Process Incoming Shipment: Full-Service Divisions	468
Process Incoming Shipment: Transfusion-Only Divisions	469
Modify Units	470
Issue, Return, or Transfuse Units: Full-Service Divisions	471
Issue, Return, or Transfuse: Transfusion-Only Divisions	472

Process a Patient Order.....	473
Patient Testing.....	474
APPENDIX H: FREQUENTLY ASKED QUESTIONS	475
APPENDIX I: COMMON ERROR CORRECTIONS.....	477
Correct Patient Test Results	477
Invalidate Test Results	477
UN-Pooling a Pooled Unit (in the computer only).....	478
UN-Modifying a Unit.....	478
APPENDIX J: RECOMMENDED REPORT USAGE	479
APPENDIX K: BARCODE SCANNER CONFIGURATION AND TROUBLESHOOTING	481
APPENDIX L: UNITED BLOOD SERVICES (UBS) FACILITY BARCODES	483
APPENDIX M: SYSTEM RESPONSES TO ACTIVE TRANSFUSION REQUIREMENTS IN SELECT UNIT AND ISSUE BLOOD COMPONENT.....	485
APPENDIX N: WORKLOAD PROCESS MAPPING TO APPLICATION OPTION TABLE.....	495
INDEX.....	503

Introduction

The main purpose of the VistA Blood Establishment Computer Software (VBECS) is to automate the daily processing of blood inventory and patient transfusions in a hospital transfusion service.

VBECS is an improved Blood Bank application that facilitates ongoing compliance with Food and Drug Administration (FDA) standards for medical devices and enhances the Department of Veterans Affairs (VAs) ability to produce high-quality blood products and services to veterans. The system follows blood bank standards, standards of national accrediting agencies, FDA regulations, and VA policies.

VBECS supersedes VistA Blood Bank v5.2 for blood bank operations. VistA Blood Bank v5.2 blood unit records remaining after the transfer of patient information to VBECS are available for reference only and cannot be edited. VistA Blood Bank v5.2 validation records must be maintained for five years after the last of the blood unit records is transferred to VBECS.

Limitations and Restrictions

VBECS does not support:

- Blood donor collection activities (donor module)
- Direct data entry of antibody identification
- Direct data entry of transfusion reaction workups
- Patient testing, transfusions, and blood components recorded in Legacy VistA are not available in VBECS
- Pediatric blood product preparation
- Printing of blood product labels
- User defined customized queries or reports
- User editing blood product, antigen, antibody, or transfusion reaction tables
- Validation record storage

Limitations and Restrictions are listed in Appendix D: Limitations and Restrictions by option as well as throughout this document as part of the detailed information associated with the functionality.

 Unauthorized access or misuse of this system and/or its data is a federal crime. Use of all data, printed or electronic, must be in accordance with VA policy on security and privacy.

 The U.S. Food and Drug Administration classifies this software as a medical device. Unauthorized modifications will render this device an adulterated medical device under Section 501 of the Medical Device Amendments to the Federal Food, Drug, and Cosmetic Act. Acquiring and implementing this software through the Freedom of Information Act require the implementer to assume total responsibility for the software and become a registered manufacturer of a medical device, subject to FDA regulations.

Freeware Disclaimer

Requests for VBECS may be submitted through the Freedom of Information Act (FOIA). Effective on September 19, 2011, the Department of Veterans Affairs will accept FOIA requests electronically. Most VA Central Office (VACO) Offices are in the process of establishing an electronic mailbox. In the meantime, the FOIA Service will accept electronic FOIA requests for any VACO FOIA Office that does not yet have a mailbox in place, and will forward your request to the appropriate office for processing.

Submit requests to http://www.foia.va.gov/FOIA_Offices.asp

Related Manuals and Materials

- *VistA Blood Establishment Computer Software (VBECS) Technical Manual-Security Guide*
- *VistA Blood Establishment Computer Software (VBECS) Release Notes*

Copies of *VistA Blood Establishment Computer Software (VBECS) User Guide* may be obtained from the VDL: VistA Documentation Library (VDL), Clinical Section, Laboratory: VistA Blood Establishment Computer Software (VBECS) Web site.

Orientation

This section covers:

- How This User Guide Is Organized
- Security: who can use VBECS
- Application Architecture: the architectural requirements for VBECS
- Hardware and Infrastructure Architecture: the system environment for installation and operation of VBECS
- VistA and VBECS: how VBECS compares to VistA
- Customer Support: how to get help

How This User Guide Is Organized

Terms

For consistency and space considerations, the pronouns “he,” “him,” and “his” are used as pronouns of indeterminate gender equally applicable to males and females.

In many instances, a user may scan a barcode or enter data manually (by typing). The term “enter” is used throughout this guide to mean “enter manually.”

See the Glossary for definitions of other terms and acronyms used in this guide.

Background and Basic Knowledge

The Introduction, Orientation, and Using the Software sections provide the background and basic knowledge VBECS users need.

Figures and Tables

If you refer to figures and tables from the user guide in your local policy and procedure documents, you may wish to use their titles only, without figure or table numbers: as the user guide is updated, those numbers may change.

Options

The Configuring Site Parameters, Component Processing, Processing Orders, Patient Services, and Reports sections cover VBECS options.

Each option section starts with an introduction, a short paragraph that explains what a function allows a user to do, followed by:

- Assumptions: tasks that must be completed before performing a function.
- Outcome: results of performing a function.
- Limitations and Restrictions: features of the system that are restricted by and exist within this function.
- Additional Information: information that helps the user perform a function.
- User Roles with Access to This Option: which users may perform this option in whole or in part.
- A detailed description of what the function allows each user to do. Two or more arrows (►) indicate functions limited to a specific role.
- Table of steps for the basic path (user and VBECS actions).
- Tables of steps for alternate paths (user and VBECS actions).
- Figures (examples of reports, screen shots, etc.).

- Table of alerts that may occur during an option.

A  (bullhorn) indicates that VBECS emits an audible alert. Keep the sound on to hear the alerts (see Volume for instructions on adjusting the volume).

- Outlined, italicized text highlights tips, warnings, and cautions, for example:

 signals potential hazards throughout this guide as they occur.

Supplemental Information

The third group of sections (References, Glossary) provides knowledge that may be needed beyond the standard, daily use of VBECS.

Appendices

The appendices contain truth tables, diagrams, and other materials for reference.

Click a section name or page number in the table of contents to move to a section. The index does not incorporate this feature.

Security

User roles identify who can access which functions. Not every VBECS user has access to all VBECS options or functions within options. Two or more arrows () indicate functions limited to a specific role.

 It is critical that each site evaluate its staff and determine appropriate security levels based on training and qualifications.

User roles also govern which options are available and which warnings can be overridden. VBECS displays options that are not available to a user in gray text.

Each of the six user roles is associated with a security level (see Table 1). Privileges accumulate as the security level increases. For example, a Lead Technologist's privileges include those of a Blood Bank Technologist and Enhanced Technologist.

Table 1: User Roles and Functions

Security Levels and User Roles	Accessible Functions	Who Should Fill These Roles?
Level 1: Blood Bank Technologist (all users)	<ul style="list-style-type: none"> • Activate/edit some shipper information. • Activate/edit some blood product information. • Standard access (includes patient testing, accepting and canceling orders, modifying units, processing shipments, processing transfusion reaction workups, accessing report functions) 	Rotating technologists, new employees, and students trained to perform associated tasks and overrides.

Security Levels and User Roles	Accessible Functions	Who Should Fill These Roles?
»» Level 2: Enhanced Technologist	<ul style="list-style-type: none"> • Edit unit cost. • Modify/issue expired blood products. • Release blood products from quarantine. • Edit unit login prior to defining patient associations. • Edit verified unit confirmation testing. • Edit patient record/verified data and test results (except the ABO/Rh for historic record and the antibody field, which results in an antigen negative requirement). • Release directed units to the main blood supply. 	Experienced technologists trained to process associated overrides and options.
»»» Level 3: Lead Technologist	<ul style="list-style-type: none"> • Select/issue ABO incompatible blood. • Select/issue antigen positive or untyped red blood cells (for the clinically significant antibodies in the table). • Print internal reports. • Maintain login message. 	Experienced technologists trained to process associated overrides and options (e.g., evening- or night-shift supervisors).
»»»» Level 4: Traditional Supervisor	<ul style="list-style-type: none"> • Maintain minimum levels. • Edit permanent record of the patient's historic ABO/Rh. • Edit Transfusion Requirements (TRs) and Special Instructions (SIs). • Edit patient record red cell antibody permanent fields (match antigens in the antigen table). • Edit a blood unit's record (verified unit data when unit has any previous patient associations). • Edit a Patient's Transfusion Record. • Remove a Blood Unit's Final Status. 	Experienced blood bank supervisors trained to process associated overrides and options
»»»»» Level 5: Enhanced Supervisor	<ul style="list-style-type: none"> • Edit antigen frequency table to reflect local population. • Add new local shipper, not already on the table issued with VBECS. • Edit and create canned comments for the division. • Edit and create consultative reports, templates, MSBOS, reagents, and equipment maintenance logs for the division. 	Experienced blood bank supervisors or VBECS administrators trained to process associated overrides and options

Security Levels and User Roles	Accessible Functions	Who Should Fill These Roles?
►►►►► Level 6: Administrator/Supervisor	<ul style="list-style-type: none"> • Define VBECS users' security settings. • Define the division as full service or transfusion only. • Enable electronic crossmatch at the division. • Populate International Council for Commonality in Blood Banking Automation, Inc. (ICCBBA) number for the division. • Define valid unit modifications for the division. • Assign workload codes for the division. • Configure testing. • Configure division. 	VBECS administrators trained to process all overrides and options not designed to change frequently. This role may be assigned temporarily to a blood bank supervisor at setup to configure the division.

The VA Wide Area Network (WAN) is protected from unauthorized access. The hardware running the VBECS application and its database will be installed at the Austin Information Technology Center [AITC (primary site)] and the Philadelphia Information Technology Center [PITC (disaster recovery site)]. Both sites are secure data centers where they are protected against unauthorized modifications and accidental damage.

Application Architecture

The primary purpose of VBECS is to automate the daily activities of a hospital transfusion service: blood inventory and patient transfusion processing. The major goals of the application are:

- Improve the safety of blood/blood component transfusions by decreasing the risk of errors through effective use of current technology such as barcode scanning, retrieval of previous records to be used for comparison, and detection of inconsistencies in data input
- Improve the quality of patient care through evaluation of transfusion appropriateness and follow-up on the effectiveness of transfusions
- Provide comprehensive reporting capabilities for quality monitoring within the transfusion service and for clinicians

VBECS operates as a discrete application, in accordance with FDA regulations regarding change to a medical device. It interfaces with the current VistA applications, using Health Level 7 (HL7) and VistALink. Areas of operation that are part of the VistA Blood Bank v5.2 software medical device (except for the donor module) are contained in this encapsulated VBECS package. Patient orders originate and are received from the Computerized Patient Record System (CPRS) and Laboratory applications.

Hardware and Infrastructure Architecture

The VBECS application requires that hardware and system software serve five users in a standard configuration and up to twenty-five users in an integrated Veterans Integrated Service Network (VISN) environment. Figure 1: System Schematic describes the major system components: a Windows 2008 Server system (the execution environment for the VBECS application) and Windows 7 workstations for access to VBECS through Windows Remote Desktop Connection. The VBECS server communicates with and exchanges information with VistA applications through Extensible Markup Language (XML) and Health Level 7 (HL7) messages over Transmission Control Protocol/Internet Protocol (TCP/IP) networking. The VBECS server communicates with BCE through HL7 messages over TCP/IP networking.

VBECS 2.1.0 is installed in a virtualized network using vSphere 5.5 as the virtualization platform.

SQL Servers exist in Windows Server Failover Clustering with AlwaysOn cluster which consists of three nodes. The Primary and High Availability server reside at the primary site while a Disaster Recovery server resides at an alternate location:

- Primary server: This server fields all requests. Its data is replicated to the High Availability and Disaster Recovery servers.
- High Availability (HA) server: This server provides backup database services through synchronous replication. Its data are guaranteed to be consistent with the Primary.
- Disaster Recovery (DR) server: This server resides at a remote site and provides backup database services through asynchronous replication. It will become the Primary server should both the Primary and HA server fail.

Options That Require VistALink

These options require a connection through VistALink. When VistALink is not active and VBECS requires access to VistA data, VBECS prompts the user to activate VistALink. When VistALink is not available, VBECS warns the user and asks him to try again later.



Contact your local support staff immediately when VistALink is not available. VBECS, VistA and CPRS all depend on VistALink to exchange information.

Main Menu: Orders

- Accept Orders: Pending Order List

Main Menu: Blood Units

- Edit Unit Information (used when editing patient restrictions)
- Return Issued Units To Blood Bank

Main Menu: Patients

- Post-Transfusion Information
- Return Issued Units To Blood Bank
- Issue Blood Components (to include a hospital location)

Main Menu: Shipments

- Incoming Shipment (used when restricting units to patients)

Main Menu: Tools

- Blood Products
- Component Classes
- Transfusion Complications
- Transfusion Effectiveness
- Workload Codes

Main Menu: Supervisor

- Document ABO Incompatible Transfusion

Main Menu: Reports


- Transfusion Complications Report
- Transfusion Effectiveness Report

- Transfusion Requirements Report
- Medication Profile (Reports or from the Patient Tool Bar icon)

Main menu window, icon in the bottom of the window:

- Display Order Alerts

Figure 1: System Schematic

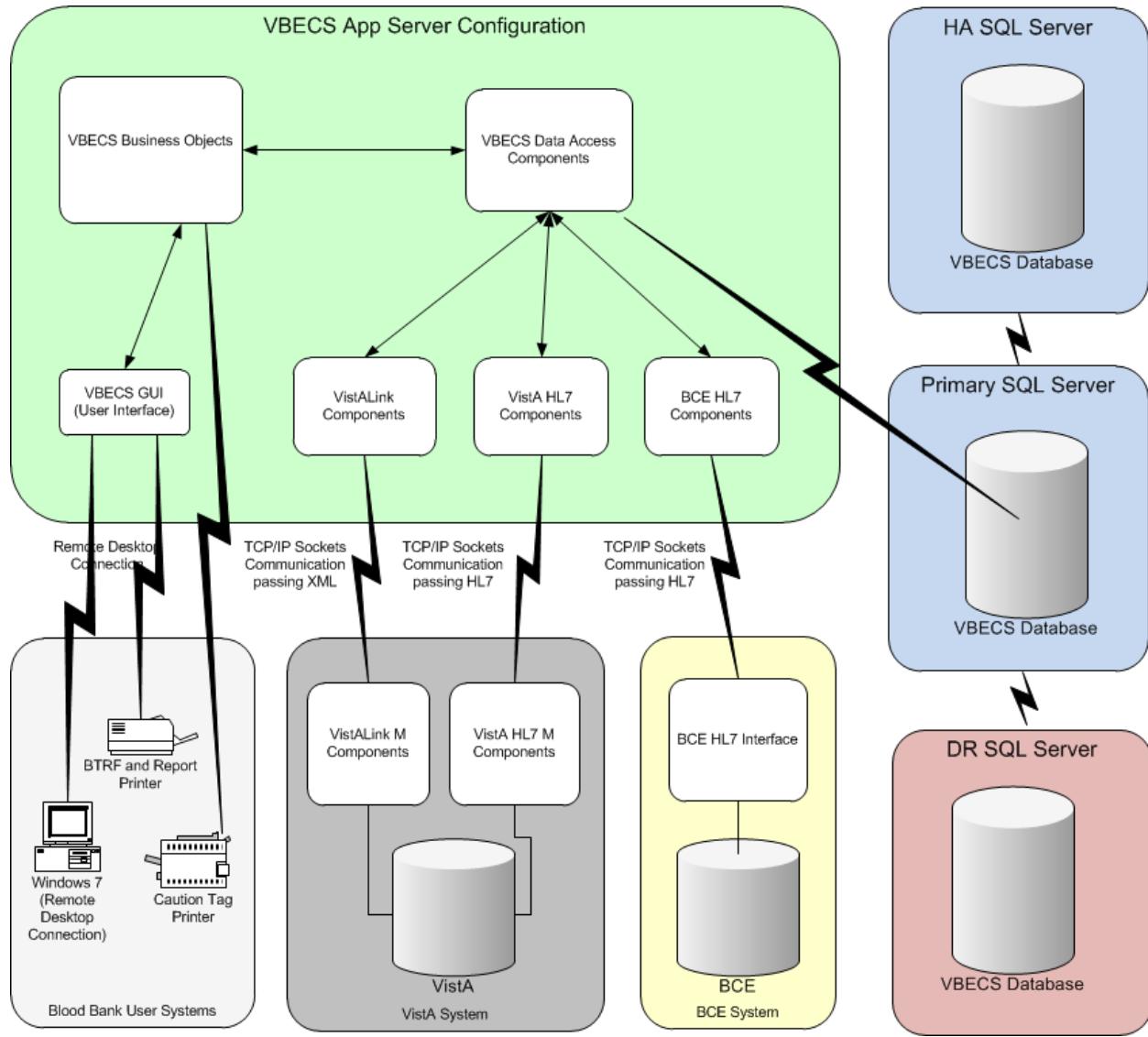

Users may not change or replace operating platforms.

Customer Support

Problems?

Contact your **Local** Information Resource Management (IRM) or Laboratory Automated Data Processing Application Coordinator (ADPAC) if you encounter VistA or CPRS connection problems and for training support before contacting the National Service Desk (NSD).

Please ensure local VistA Support contact information is available at all times. National Support will engage Enterprise Operations (EO) personnel as needed.

Problems with connectivity to VistA and CPRS may require personnel from EO with VBECS server administrator access and VistA IT support access.

If the problem remains unresolved after local VistA triage, contact the national Information Technology (IT) support service or request a region ticket to be reassigned to:

National VA Service Desk Contact

For Information Technology (IT) support, call the National Service Desk (NSD), 888-596-HELP (4357) (toll free), 24 hours per day, 7 days per week. [Users with access to the VASD-supported request tool (e.g., CA SDM) may file a ticket in lieu of calling the NSD.]

National VA Service Desk Alternate Contacts

- Web site: <http://vaww.aac.va.gov/servicedesk/> (National Service Desk Tuscaloosa).
- Email : NSDTuscaloosaUSD@va.gov

Using the Software

Using the Software explains:

- How to access VBECS and VistA
- What users will see on the screen
- What the icons and other symbols mean
- Which conventions are used in this guide
- How to exit VBECS

Software Basics

- When a message window requires a user response and VBECS displays it behind the active window, the user must minimize the active window or click and drag dialog boxes to expose the message window.
- Press the **F1 key** to view screen-specific online help.
- A user may click **Cancel**, **Close**, or **☒** in the upper right corner at any time to exit without saving.
- All operations depend on the user being logged into the appropriate division.
- A user may click a column heading to resort the items in that column (Figure 2 and Figure 3).
- Reaction results entries are allowed for only those reagents defined as part of the reagent rack.
- Select a date from the calendar in the drop-down list to edit a date, or click each part of the date to edit the date and time (Figure 4 and Figure 5).
- Labels shown in figures are visual aids only and do not represent actual labels.

- VBECS does not offer a spell-check feature. (The user may write text in Microsoft Word, spell check the text, and then copy and paste the text into VBECS.)



Although the copy-and-paste function is enabled, do not use it to enter patient identification information. It is highly recommended that the user scan this information when possible; otherwise, enter it.

Figure 2: Click to Sort by Patient Name

Patient Name	Patient ID	Lab Order No.	Urgency	Specimen	Request	Ordered By	Patient Location	Received
VBPATIENT, EIGHTY	000-00-0080	65014	Routine	Required	TAS	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYONE	000-00-0081	65022	Routine	Required	RED BLOOD CELLS x 2	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYONE	000-00-0081	65022	Routine	Required	TAS	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYTHREE	000-00-0083	65060	Routine	Not required	PLATELETS x 3	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYTHREE	000-00-0083	65060	Routine	Required	TAS	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYTWO	000-00-0082	65028	Routine	Not required	FRESH FROZEN PLASMA x 2	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/

Figure 3: Click to Sort by Request

Patient Name	Patient ID	Lab Order No.	Urgency	Specimen	Request	Ordered By	Patient Location	Received
VBPATIENT, EIGHTYTWO	000-00-0082	65028	Routine	Not required	FRESH FROZEN PLASMA x 2	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYTHREE	000-00-0083	65060	Routine	Not required	PLATELETS x 3	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYONE	000-00-0081	65022	Routine	Required	RED BLOOD CELLS x 2	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYTWO	000-00-0082	65028	Routine	Required	RED BLOOD CELLS x 2	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTY	000-00-0080	65014	Routine	Required	TAS	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/
VBPATIENT, EIGHTYONE	000-00-0081	65022	Routine	Required	TAS	USER, ONE	VBEC LEAVENWORTH OUTP	01/31/

Figure 4: Select a Date from the Calendar


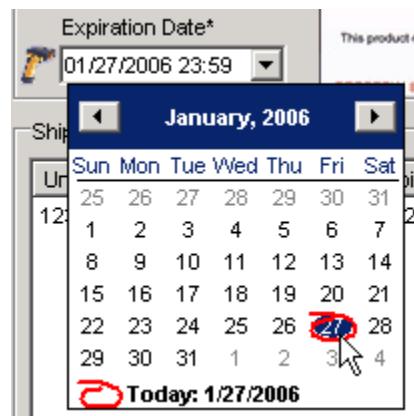


Figure 5: Edit the Date and Time

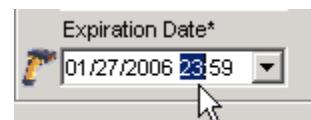

Commonly Used System Rules

This section includes system rules that apply to several or all options.

Date and Time

- Sites are not notified when an administrator changes the date or time on the server. Users must contact the National Help Desk when the date or time is incorrect.
- The date, time, and time zone in VBECS is defined where the blood bank is located. Any requests for data from a particular blood bank database will reflect the time zone where the data was created and not the time zone where the report is generated.
- Refer to *VistA Blood Establishment Computer Software (VBECS) Technical Manual-Security Guide* to change the time zone.

Locking and Time-Outs

- In many cases, VBECS's locking feature prevents more than one user from accessing data. When a user attempts to access locked data, VBECS displays a message indicating which user has the locked data.
- A local configurable setting (Lock Inactivity Timeout field) specifies how long a user can be idle and in control of data being edited. VBECS warns the user 60 seconds before the lock inactivity period expires that he will lose priority for the data. When he responds within 60 seconds, VBECS clears the warning and resets the lock activity timer. Otherwise, VBECS informs him that his lock was released and he must reenter his changes.

Working with Data

- VBECS will emit an audible alert when a user exceeds the maximum character limit for a field.
- The user may enter interpretations and reactions in full-service facilities. The user may enter only interpretations in transfusion-only facilities.
- Information displayed on an active window is updated with information from other options only when it is refreshed (i.e., closed and reopened) in the active window.
- When VBECS displays a scanner icon next to a data entry field, a user may enter or scan the data. The Caps Lock key must not be engaged.
- VBECS accommodates Codabar and ISBT 128-labeled units. VBECS uses special data identifiers embedded in an ISBT 128 product code to determine the label type scanned. The lack of data identifiers indicates that the barcode is for a Codabar unit (default: unknown).
- VBECS warns the user when a Codabar unit ID contains “A,” “B,” “I,” “O,” or “U” (invalid characters) and does not save the unit.
- VBECS captures all technologist IDs and the name of the technologist who displays and/or prints testing information on reports associated with entering test results on a worksheet.
- With each valid keystroke entry in the antibody screen and crossmatch test grids, VBECS moves the cursor from top to bottom in the far left column within the cells in which results are being entered. The cursor returns to the top cell in the next column to the right and so on across the worksheet grid. The cursor stays in the reaction results cells by default (the technologist may navigate independently) until the last entry is made in the reaction result portion of the grid, then the cursor moves to the first (top) interpretation cell.
- To clear an entry in a test grid cell, press the **Backspace** key.
- VBECS allows selection of deceased patients without a warning message presented when VBECS has not received a patient death update message from VistA.

- VBECS displays the data that a user enters in a session. The user may edit the data and save them. When a user cancels, VBECS warns that it will not save the data. VBECS closes the form and returns the user to the main menu screen that may include unrelated open windows.
- VBECS displays and enforces active Transfusion Requirements (TRs) to all divisions in a multidivisional database.
- The user may not edit the title column of the data grid that specifies the reagent used in the test with a patient specimen.
- Unless otherwise indicated, data entry and retrieval are restricted to the division logged into.
- To enter data after VBECS has been unavailable, a user may enter a past date, time, and the technologist who performed the work.
- VBECS displays changes in reverse chronological order.
- VBECS calculates and displays the days (hours, if less than one day) remaining until a unit's expiration based on the date and time of unit selection entered during modification.
- The format of an expiration date may be "mmddyy" or "mmddyyyy" (Codabar).
- Data are never deleted: they are saved indefinitely.
- Although a user may appear to select and delete data in gray fields, VBECS does not save such edits.
- A column's width may be adjusted or hidden for a session. This view returns to the default presentation when the view is opened again.

Searching the Database

- A user may search for a patient by:
 - Last name
 - ID (omit hyphens) ("Full")
 - Initial of the last name and last four numbers of the ID ("Last 5")
 - Last four numbers of the ID ("Last 4")
- A user must scan or enter a unit ID and a product code to select a unit from inventory.
- When VBECS does not find a unit record, it displays an error message stating that the unit was not found and instructing the user to reenter unit information. A user may close this message, click **Clear**, and enter another unit ID and product code.
- VBECS associates the technologist ID, date, time, and division with each process for retrieval by division.
- VBECS lists all or part of a division's blood component inventory, based on system rules. A user may enter a varied combination of unit information to create a customized list of selectable units, including:
 - Unit ID
 - Product type
 - ABO/Rh
 - Special testing
 - RBC antigens (available only for the RED BLOOD CELLS component class)
 - Received date range
 - Expiration date range

VBECS displays this unit information in the selection list whether full inventory, subset, or filtered subset:

- Unit ID
- Unit ABO/Rh
- Product short name
- Product code

- o Unit status
- o Unit expiration date
- o Donation type
- o Date and time received (standard format)

Testing

- ➤➤➤➤➤ Administrator/Supervisors may configure reagent racks to use a tube, solid phase, or gel method.
 - o VBECS allows a user to perform QC testing at AHG and later perform patient testing at all phases using a rack configured for AHG only.
 - o VBECS does not verify that the phase of testing matches the reagents used for testing.
 - o Users working with reagent racks containing different reagents must do so in accordance with local policy and procedures.
- Reaction result cells in the data grid contain only symbols. A legend of the symbols and their meanings, specific to the data grid and its truth tables, is located in the bottom left corner of the form and is available to the user throughout data entry.
- Entry of a symbol (not a number or letter) into a test grid reaction result cell or interpretation field will not elicit a system response. The entry of the symbol is not allowed and does not overwrite the previous entry, if any. See Table 2 for valid observed test result entries.
- After entering a canned or free-text comment in the testing comment field of the testing grid, the user must press **Enter** to store the comment.
- The user may enter canned comments in comment cells by selecting from a pick list, or enter free-text comments by selecting “Other.”

Table 2: Valid Observed Test Result Entries

Entry	VBECS Translation	
H	Hemolysis	Positive
W	Weak	Positive
F	Mixed Field	Positive
M	Microscopic	Positive
1	1+	Positive
2	2+	Positive
3	3+	Positive
4	4+	Positive
X	Not Tested	Not Tested
0	No Agglutination	Negative
R	Rouleaux, No Agglutination	Negative

- When recording test results:
 - o In the Test Details window, click a check box to select a rack.
 - o In the Date Tested field, select the default date and time (current date and time) for testing performed, which the user may edit to a past date and time.
 - o In the Tested By field, select the default identification (current user) for the testing technologist or another user from the drop-down menu.
 - o Click the appropriate radio button to determine the grid configuration for antibody screen and crossmatch grids [All phases or Anti-Human Globulin (AHG) only (pre-checked default varies with local configuration)].
 - o When an automated instrument is used, click the **Using Automated Instrument** check box, enter the instrument name in the Instrument Name field, and click the **Instrument QC'd?** check box, as appropriate.

- VBECS requires units in ICCBBA component classes Red Blood Cells, Whole Blood, Granulocytes and Leukocytes to have ABO confirmation testing (forward typing), and Rh confirmation testing for units labeled as Rh negative.
- The user may not edit VBECS-populated grid worksheet fields that contain unit or patient data.
- The user may sort worksheet entries only before beginning data entry, preferably during unit selection. When the user wishes to sort a worksheet after entering the first reaction result, he must exit without saving the work and begin again.
-  When a user enters an invalid symbol in a cell, VBECS emits an audible alert and prevents the cursor from moving to the next cell.
- VBECS moves the cursor to the comment cell automatically, when applicable; the user may also use the mouse to move to another cell.
- VBECS displays **E**(error) at the beginning of invalid rows and **V**(valid reagent) or **W** (reagent will expire within 24 hours) at the beginning of valid rows. These warning icons change as the user updates information. When the user holds the mouse over an icon, VBECS displays the message associated with the truth table interpretation (as with a tool tip).
- VBECS compares user-entered interpretations with truth table data.
- The user performs tests and enters data simultaneously.
- Throughout this option, the user may view the patient's VBECS demographics, including the information entered in the patient Special Instructions (SIs), Transfusion Requirements (TRs), and other clinical lab results related to the transfusion through CPRS.
- The historical ABO/Rh display includes the date and division of creation or indicates that the record is from the VistA database conversion.
- The test record captures the entering and testing technologists' ID and the display/print technologists' name, the division, date, and time of the save, the content (details of the verified testing), QC rack ID (if not by automated instrument), test method, automated test instrument information, if indicated, data entries, interpretations, and comments. (When the configuration changes in Configure Testing after the test data are entered, the record is retrievable as originally configured.)
-  When a site is defined as "full service" and daily QC was not performed on the selected rack or the Not in Use Today check box was checked for the current 24-hour period (beginning with the time set for the alert defined in Configure Daily QC), VBECS emits an audible alert, warns the user, and asks whether he wishes to continue. **Yes** requires a comment and captures details for inclusion in the Exception Report (exception type: QC not performed on rack used for testing). **No** allows the user to select a different rack or cancel and return to the main menu.
- Repeat tests are not billable. Repeat tests are not available to the patient chart unless they are part of a corrected report test.
- When reagents or antisera are within 24 hours of expiration (23:59 on the expiration date), VBECS warns the user when the reagent will expire and asks whether the user wishes to continue. **Yes** allows the user to continue testing. **No** moves the cursor to a cleared lot number cell, which is ready for the user to enter a new lot number.
- The unit record will capture the entering tech ID, and display and/or print the entering tech name; the selecting tech ID, and display and/or print the selecting tech name; division, date, and time of the save; and the content (details of the unit selection), if indicated, data entries, interpretations, and comments.
- VBECS allows a user to save partially completed tests and allows the same or different user to complete the tests later. A first or second user may enter only new results and may not edit prior entries. VBECS allows the user to save partially completed testing. The user may return to complete the pending tasks, but may not edit saved interpretations.

Creating and Viewing Reports

- A division's printer name cannot exceed 50 characters in length.
- Click in the Include Report Sections area to select all report sections (Figure 6).
- Click in the Include Report Sections area to clear all report sections (Figure 7).

Figure 6: Click to Select All Report Sections


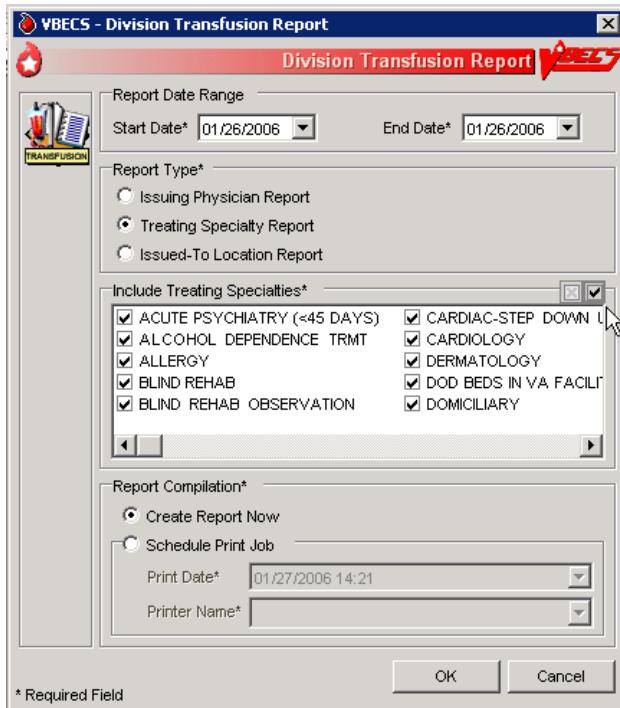


Figure 7: Click to Clear All Report Sections

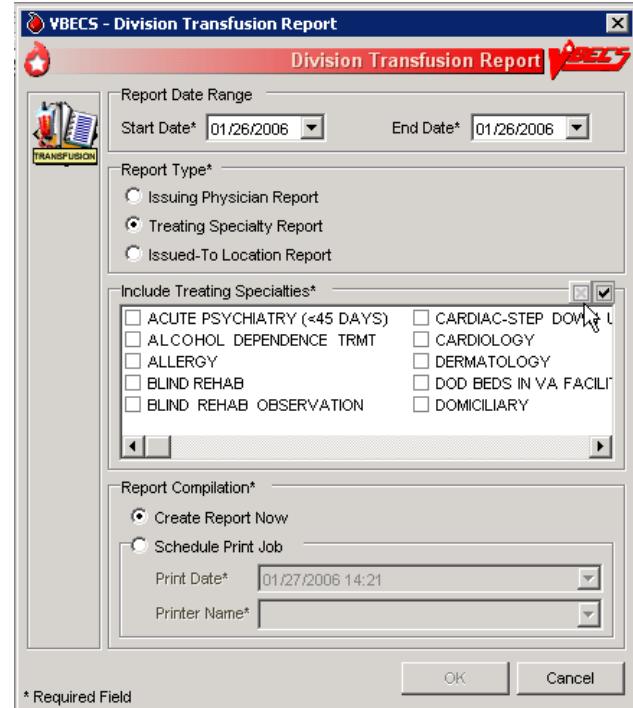

- In the Report Compilation area, click the **Create Report Now** radio button and **OK** to view the report, or click the **Schedule Print Job** radio button and **OK**, then enter or select a print date and printer name, then click **OK** to view and print the report.
- Reports are printed only in their entirety.
- A single-division report includes information from one division.
- When a site is a single division or a user does not have multidivisional access and a report has a multiple-division option, the user may not select a multidivisional report.
- VBECS maintains a cumulative history of all transactions, including testing and status changes, for every blood unit entered in the database (Unit History Report).
- VBECS captures changes to verified data for inclusion in an audit trail report (Audit Trail Report, Unit History Report, Patient History Report, and Order History Report).
- VBECS displays patient ABO/Rh results in reverse chronological order, regardless of division, then from the VistA database conversion. When no ABO/Rh records exist, VBECS displays "NR" (no records). When there are discrepant ABO/Rh results from the same universal ID (UID) or from two specimen UIDs, VBECS displays the results recorded just prior to the latest result that created the discrepancy and warns the user.


- When a search will return a large number of items on the report, VBECS notifies the user that the report will be longer than 50 pages and asks whether he wants to display and perhaps print the report. **Yes** notifies the user that the process will take a few minutes. **No** clears the screen and returns to the option's main menu.
- When no information that fulfills search criteria is available, VBECS notifies the user that no information is available, per entered search criteria.
- The user may change the default start and end dates (the previous day). The end date and time must be later than or the same as the start date and time. Report activity begins at and includes midnight on the start date, and ends at and includes midnight on the end date. When the date range includes the current date and time, VBECS displays "Preliminary" in the header on each page of the report.
- VBECS displays headers on every page of a report; it displays search and sort criteria and date ranges at the end of the report. VBECS identifies whether reports are generated from a test environment or from a production environment.
- A user may schedule a report to run and print in the future. When the activity date range includes a date that precedes the date and time of printing, VBECS processes the data on the date and time the report is scheduled to print. When a date range is greater than 30 days, VBECS warns the user and asks whether he wishes to continue.
- Local record retention policy instructs users to print, review, document, and save the Audit Trail, Exception, and Testing Worklist Reports. The reports contain only data entries generated in the user's login division. Each entry includes an area for supervisor acknowledgement (signature or initials and date of review) and comment.
- The search functionality (binoculars icon) may not work in a displayed report.
- When the processing technologist differs from the login technologist, the report includes both names.
- When the processing date differs from the data entry date, the report includes both dates.
- VBECS displays the patient historical ABO/Rh, antibodies identified, antigen negative requirements, and transfusion reaction history (date, division of testing, whether from the VistA database conversion), and whether the historical ABO/Rh was changed.
- When a patient has a history of an ABO/Rh change justified through Justify ABO/Rh Change, VBECS displays "History of justified ABO/Rh change" next to the historical ABO/Rh information.

Reports in VBECS are created using a Microsoft technology called SQL Server Reporting Services (SSRS). Reports are rendered on screen using ReportViewer control. Reports are not deployed on the reporting server but as a part of the main application. All reports can be exported. The user has the option to export reports in three different formats (PDF, Excel and/or Word).

Report Toolbar

Each report in VBECS will have a toolbar on top that allows the user to interact with the report (Figure 8). The following sub-sections elaborate on the controls found on the Report Toolbar. It is suggested that users create a local directory for exporting their reports (see the *VistA Blood Establishment Computer Software (VBECS) Technical Manual-Security Guide*, Configure the Report Share and the Configure a Shortcut to the Report Share sections for more information).

Figure 8: Report Toolbar

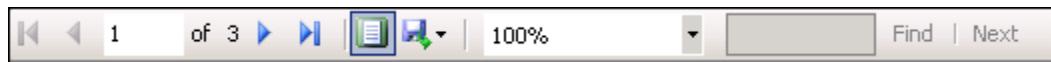

Navigation Arrows and Report Pages

Navigation arrows allow the user to move to the next or previous page of the report. If the user is on the first page of the report, the arrows pointing to the left will be disabled. If the user is on the last page of the report, the arrows pointing to the right will be disabled.

The first arrow button moves the report by one page forward (or backward for left facing arrow). The second moves the report to either the last or first page (in case of the left facing arrow). The number X of Y tells user which page is currently being viewed (X) and what is the total number of pages (Y).

Report Layout

There are two ways (layouts) in which reports can be displayed. By default, reports are displayed in Print Layout. It allows reports to appear on the screen exactly how they will look when printed. However when reports are displayed in Print Layout mode, the Report Search function (see Report Search) is disabled. In order to be able to use Report Search, the user has to turn off Print Layout mode by pressing this icon:




. When Print Layout is disabled, the icon will look like this: . It can be pressed again to return to Print Layout mode.

Page size of reports not displayed in Print Layout mode differs from the regular letter size (8.5in. x 11in.) and that is why the total number of pages for the same report displayed in different layouts may differ.

Report Zoom-in/Zoom-out

Clicking the drop-down menu showing 100% will give the user a choice of the report display size (Figure 9).

Figure 9: Report Zoom-in/Zoom-out Control

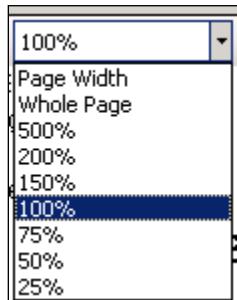

The default value is 100%. The option Page Width will resize the report to fit the screen horizontal size. The option Whole Page will display the entire page on the screen. Selecting an option less than 100% will zoom out the report making the font smaller and allowing more data to fit on the screen. Selecting an option higher than 100% will zoom in the report making the font bigger. If the report does not fit on the screen, the user will be able to use scroll bars to scroll up or down and left or right to view the entire report.

Report Search


This option allows the user to search for specific text in a report and is only enabled if Print Layout is disabled (see Report Layout). By default, the **Find** and **Next** buttons are disabled. Once the user enters text in the search box, the **Find** button becomes enabled (Figure 10).

Figure 10: Report Search Control



Clicking the **Find** button will search the entire report. If the searched text is found, the report will display the page where it was found. The **Next** button will then become enabled allowing the user to search for subsequent instances of the search string. If the searched text cannot be found, the user will be presented with the message shown in Figure 11.

Figure 11: Report Search Warning Message

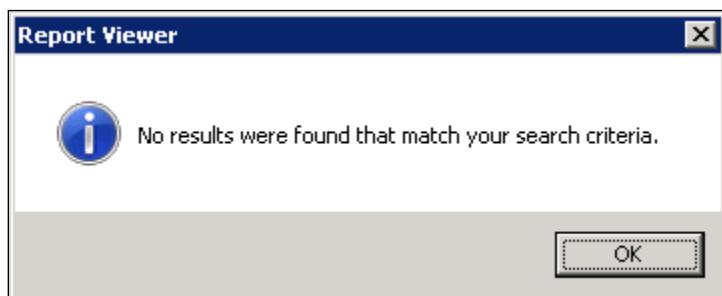

Report Export

The user has an option to export the report in three different formats: PDF, Excel and Word. The export to PDF creates a report that is Read-Only (may not be modified) and will match the report that users can see directly in VBECS in terms of data and format. The export to Excel and Word creates a report that may be modified. **VBECS reports exported to Excel or Word can be altered by a user. Only copies printed directly from VBECS should be shared with clinicians to ensure data integrity.**

The data output from VBECS is not controlled by the medical device and must be strictly regulated by local facility personnel. Many of the reports contain patient sensitive information and must be treated as such. Use of the exported reports, including alteration of data, is the responsibility of the user.

Export of the reports to Excel and Word comes with known problems. Reports exported to Excel or Word contain the same data as reports that can be seen directly in VBECS; however, their formatting may differ from the original report. Not all reports can be exported to all the identified formats. Table 3 lists the reports and their exportable formats.

Table 3: Exportable Formats for Reports

Report Name	PDF	Word	Excel
Maintain Reagent Inventory	Yes	Yes	No
Testing Worklist Report	Yes	Yes	No
Blood Availability Report	Yes	Yes	No
Division Transfusion Report	Yes	Yes	No
Blood Bank Supplies Report	Yes	Yes	No
Unit History Report	Yes	Yes	No
Patient History Report	Yes	Yes	No
Other VBECS reports	Yes	Yes	Yes

There is also a bug in Microsoft Word that causes the message “Word cannot start the converter MSWRD632.wpc” to display when opening files in Microsoft Word that were not created using Microsoft Word (e.g., a text file that was renamed as a .doc file). This message also shows for reports exported to Microsoft Word from VBECS. Microsoft published an article that describes the bug and steps on how to fix it (<http://support.microsoft.com/kb/973904>). The fix needs to be applied to the machine in which the exported reports are being opened and not to the VBECS application server. The PC used to view reports will get this error. Users can either click on the error to dismiss it each time an exported report is opened or apply the Microsoft hotfix to correct the issue.



By clicking on the report export icon, the user has the option to select and export to an Excel, PDF or Word file (Figure 12). Choosing any of these options will open the dialog shown in Figure 13.

Figure 12: Report Export Toolbar


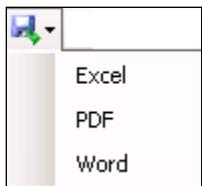


Figure 13: Example of Report Export Dialog with Assumption that Report Export Directory is D:\VBECSReports

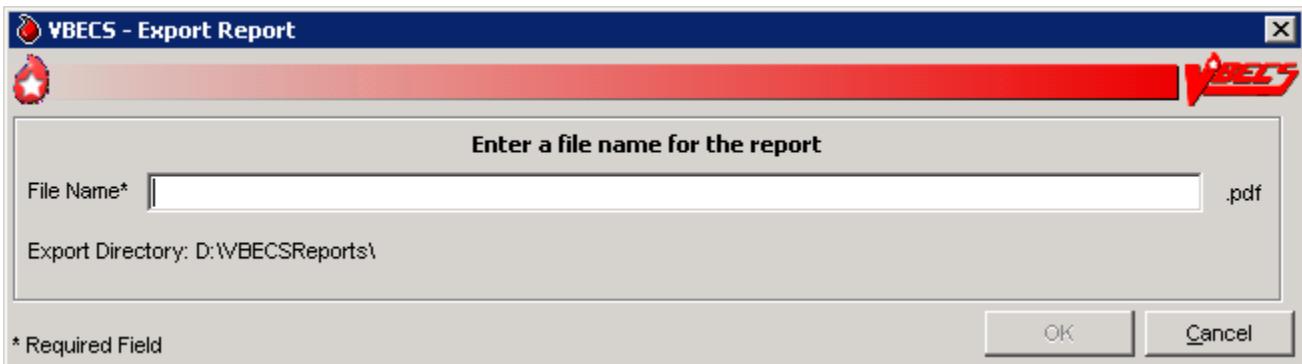


This dialog allows the user to specify the name of the file to which the report will be exported. The extension of the file is pre-filled for the user and depends on which export option was selected (.xls for Excel, .pdf for PDF or .doc for Word). The file name length is limited to 200 characters and will not accept special characters as part of the name (e.g., \ / ? % * | : ' " < >).

If the user attempts to use a special character in the file name, it will be ignored and not appear on the screen.


Once the user enters a proper file name, the **OK** button will become enabled. Clicking the **OK** button will present the user with the message shown in Figure 14.

Figure 14: Report Exporting in Progress

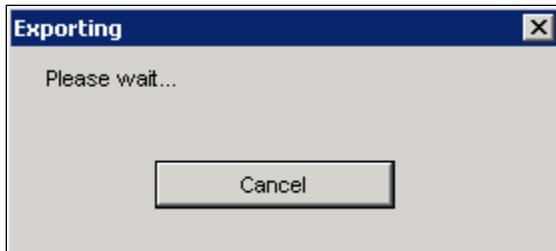

Clicking the **Cancel** button will stop the export and the user will be reverted back to viewing the report. If the export is successful, the user will be presented with the message shown in Figure 15.


Figure 15: Example of Report Exporting Success

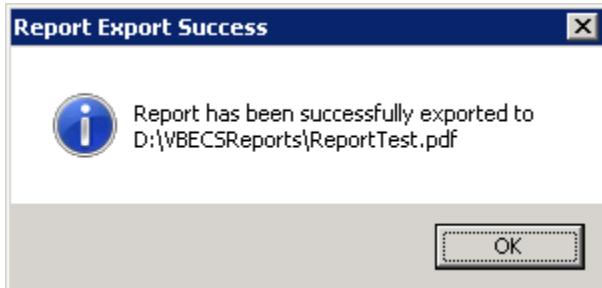

This dialog shows the full path to the report and its name allowing the user to find it quickly using Windows Explorer.

If the user enters the name of a report that already exists, the system will present the message shown in Figure 16.

Figure 16: Report File Already Exists

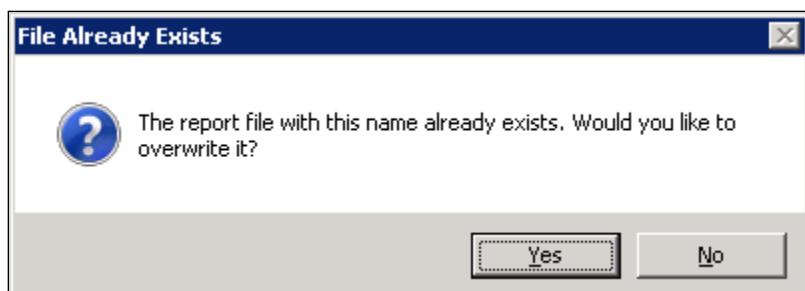

Clicking **Yes** will overwrite the existing file. Clicking **No** will allow the user to specify a new name.

Copying an Exported Report


Here are some tips for copying a file or folder of files from one location to another.

Before you begin:

- Verify that you or the person copying the file has permissions to access both the folder you are copying files from and the folder you are copying the files to.
- Determine the type of file (i.e., PDF, Word or Excel) you want to use for your local document before you export. Do not change the file type during this process.

You may have a shortcut defined on your desktop that opens the report export folder on the server (Figure 17). This allows you to open files on the server on your local PC. It is the server's D: drive, which is (Y) in this example.

Figure 17: Example of D: Drive

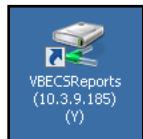

Option One: Drag and drop a copy of the file in "MY LOCAL FOLDER".

This works well when you do not need to change the file name when placing it in your local folder.

- 1) When in 'Windows Explorer' (not Internet Explorer), left-click on the file and while holding down the mouse button, drag the file to the destination folder. The dark-blue "Admin Data Report 071514" is selected (Figure 18).

Figure 18: Example of Selection of "Admin Data Report 071514"



- 2) Release the left mouse button when hovering over "MY LOCAL FOLDER". The file may copy without any notification or may have a copying window appear if it is a large file (Figure 19).

Figure 19: Example of a Copying Window




- 3) Skip to Step 3 of Option Two.

Option Two: Open the document and save a copy to another file.

- 1) Open the document from the server folder.
- 2) Select the **Save As** option and then the local folder, which in this example is called MY LOCAL FOLDER.


You may edit the file name at this point. The name you use for the file in your local folder should be something that allows you to manage those files as you use them. You can name the exported file accordingly or edit the file name in your local folder. A generic file name allows you to save and refresh the information. It may be appropriate for something like the Cumulative Transfusion Requirements report (Figure 20).

Figure 20: Example of File Name

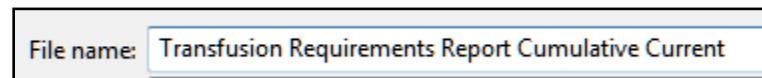


You may also opt to save with a date in the file name to allow a snapshot report for historic purposes (Figure 21).

Figure 21: Example of Dated File Name

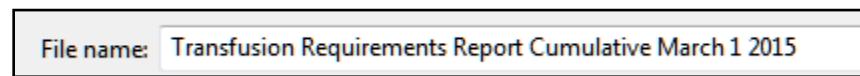

- 3) Verify that “MY LOCAL FOLDER” contains the saved file. You can edit and save to this folder with those changes (Figure 22). Unlike the server file, you may also delete this file.

Figure 22: Example of Saved File

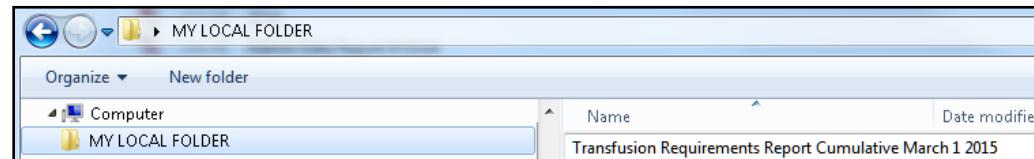

Note: The copy of the file you save to your local folder is separate and does not update the file in the server folder.

Using VBECS

VBECS is a Windows-based application.

One VBECS system is implemented for each VistA database.

Subject matter experts (SMEs) at user sites will help implement VBECS and train blood bank staff and other authorized users.

The VA Software Document Library (VDL) site contains VBECS user documents (<http://www.va.gov/vdl/application.asp?appid=182>).

The VBECS Sharepoint site contains training, Frequently Asked Questions (FAQs), Technical Bulletins, ListServ messages and sign up, etc. to support the implementation and use of VBECS (<http://vaww.oed.portal.va.gov/projects/vbecs/default.aspx>). Users may access links to VBECS training modules and frequently asked questions, and join a ListServ (an email distribution list to manage discussion groups) for announcements and updates on the VBECS Sharepoint site. Contact information is provided for comments and questions.

Documentation

Users may retrieve VistA documentation from the VistA Documentation Library (VDL) Web page.

Using Online Help

Select **VBECS Help** from the Help menu (Figure 23) to view a directory of help contents (Figure 24).

Figure 23: VBECS Help



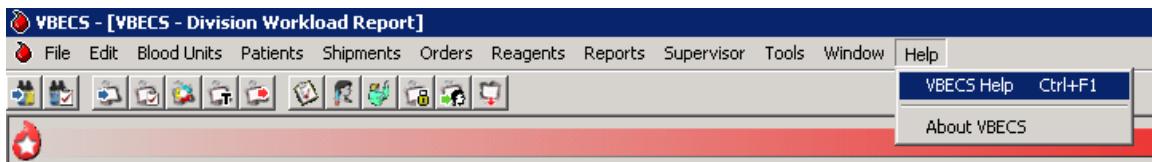


Figure 24: Example of VBECS Online Help Contents

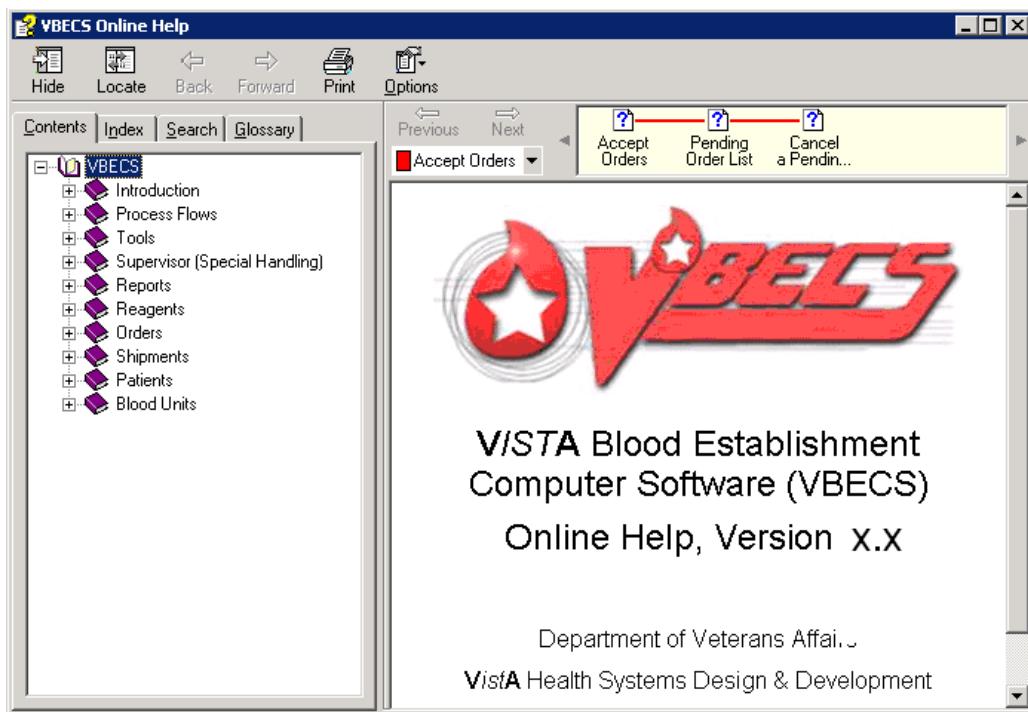

Press the **F1** key to view screen-specific help (Figure 25). Click text in blue to expose more information.


Figure 25: Example of Screen-Specific Online Help

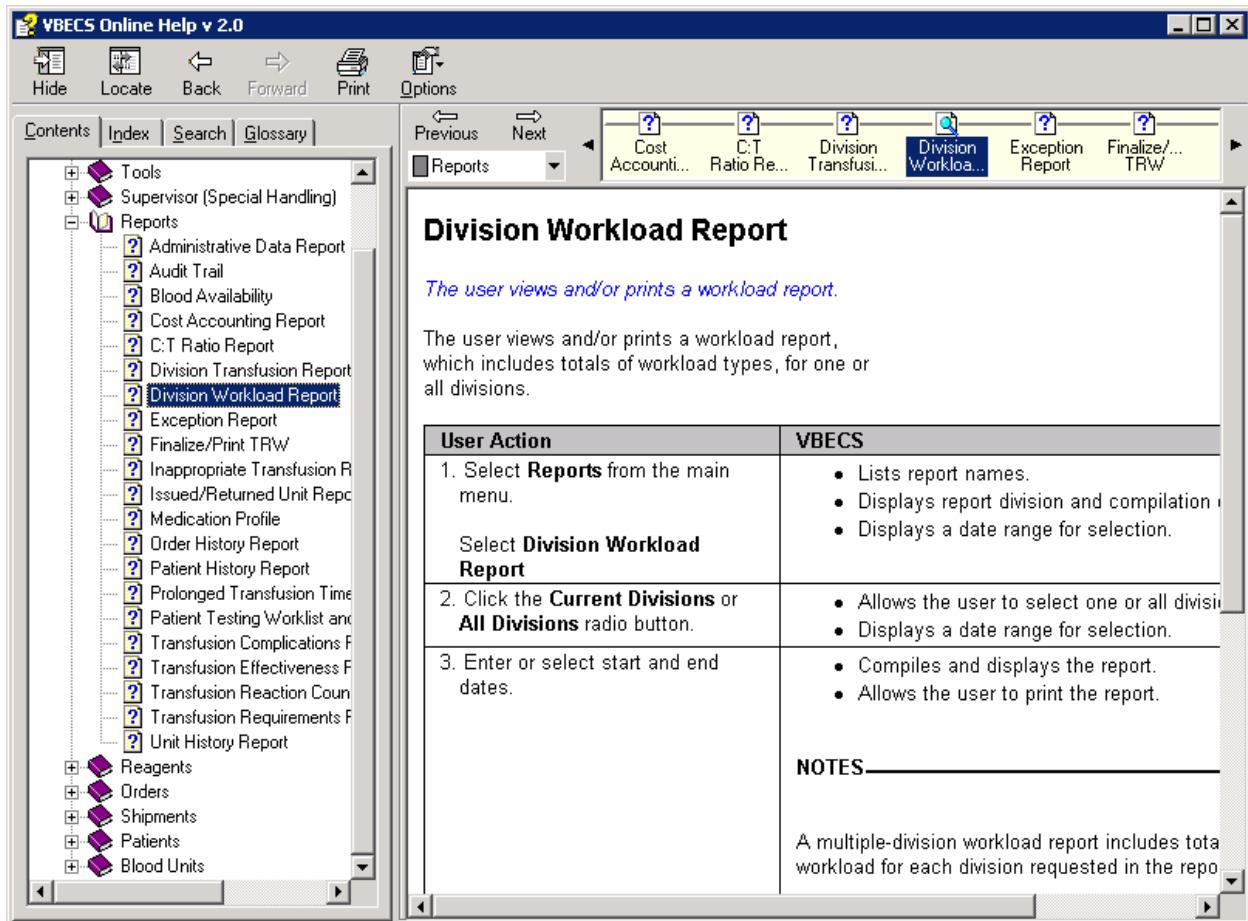

Contents Tab

Click the **Contents** tab in the navigation panel on the left to display categories of VBECS functions. Click a plus sign (+) to view the topics in each category. Click a topic preceded by a question mark to view information on that topic on the right (Figure 26). Click a minus sign (–) to condense a list.


Figure 26: Using the Contents Tab

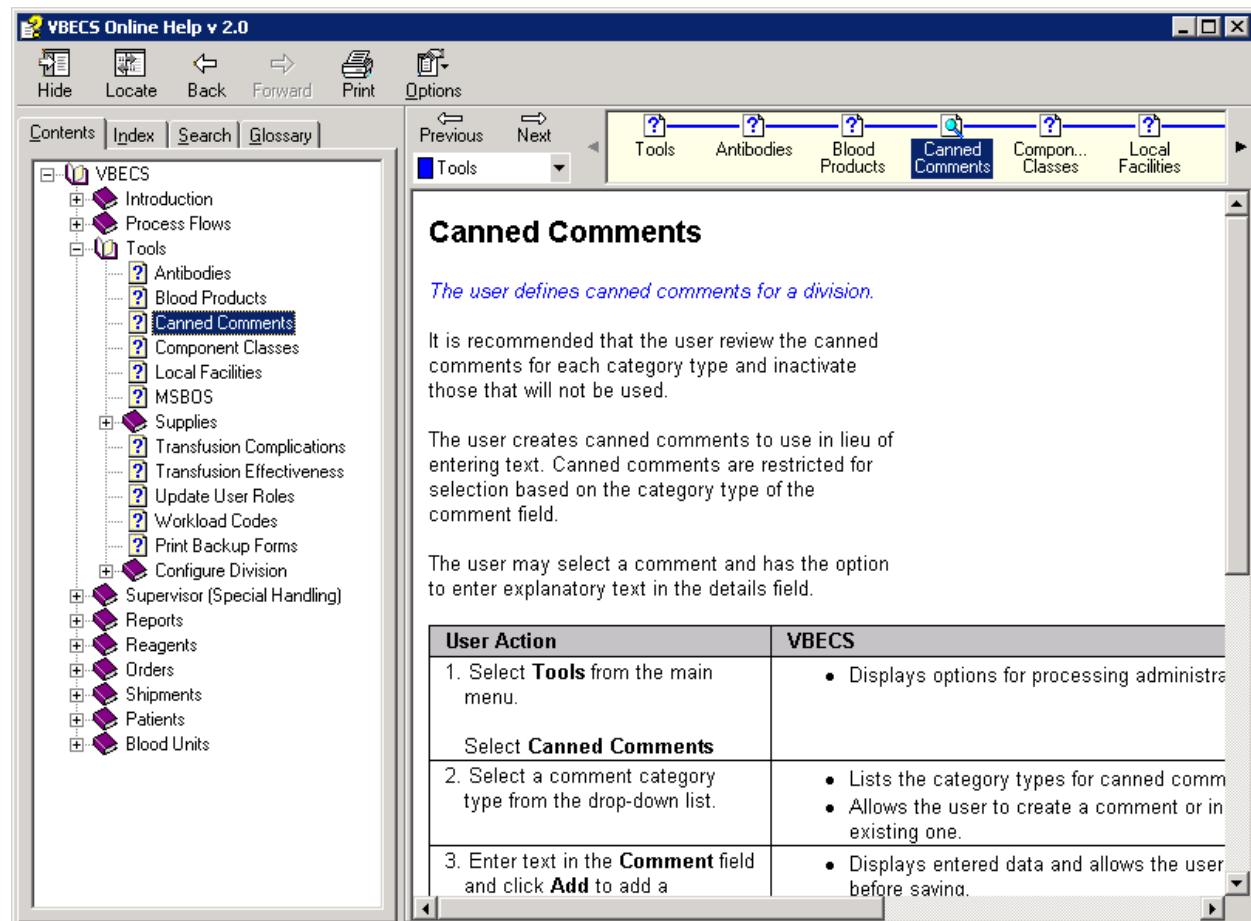

Index Tab

Enter a keyword (or one or more letters of a keyword) in the text box to find a topic. Double-click the topic to view the topic on the right (Figure 27).


Figure 27: Using the Index Tab

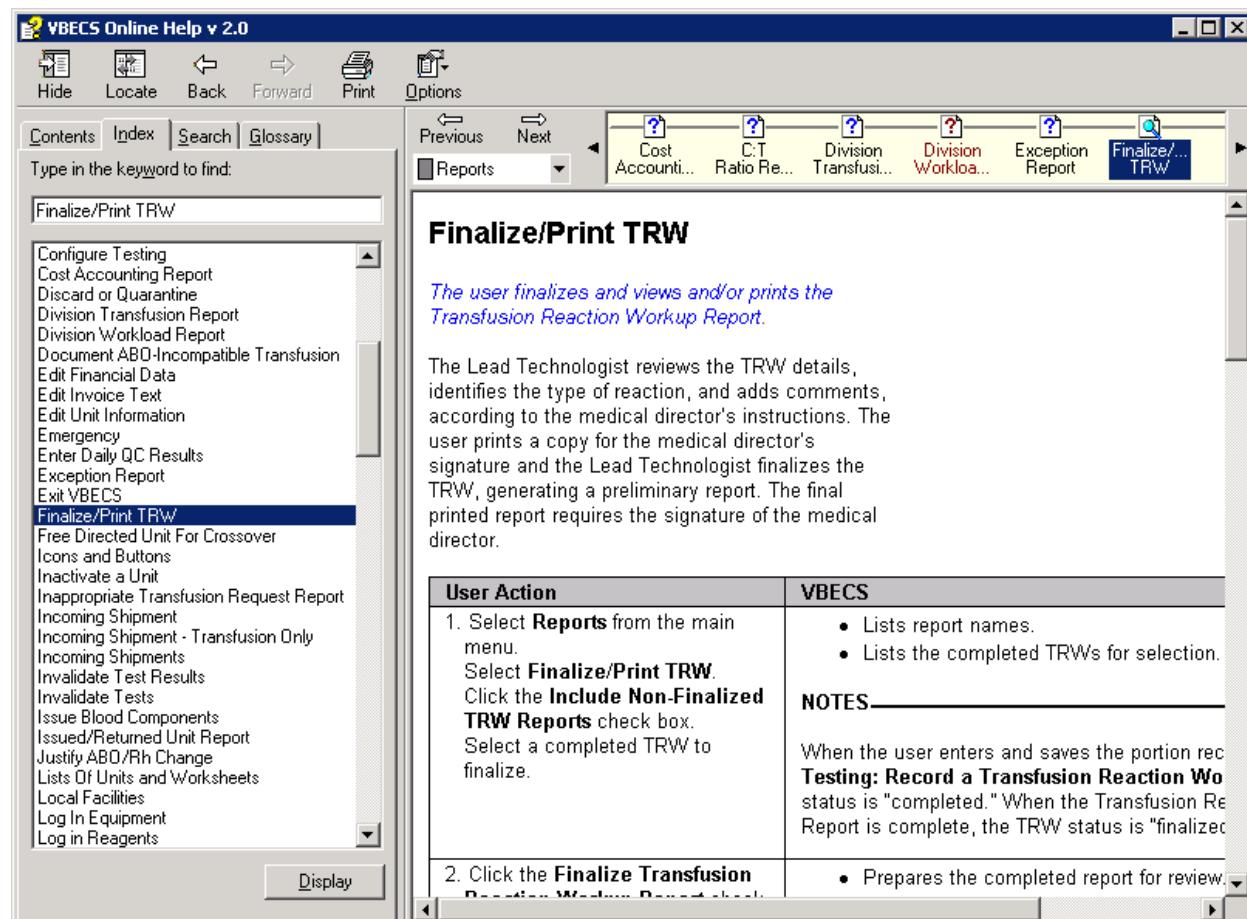

Search Tab

Enter a keyword or phrase in the text box and click **List Topics**. Select a topic and click **Display**, or double-click the topic. The topic appears on the right; each instance of the keyword entered is highlighted (Figure 28). (To turn the highlight feature on or off, click **Options** at the top of the screen. Click **Search Highlight On** or **Search Highlight Off**.)


Figure 28: Using the Search Tab

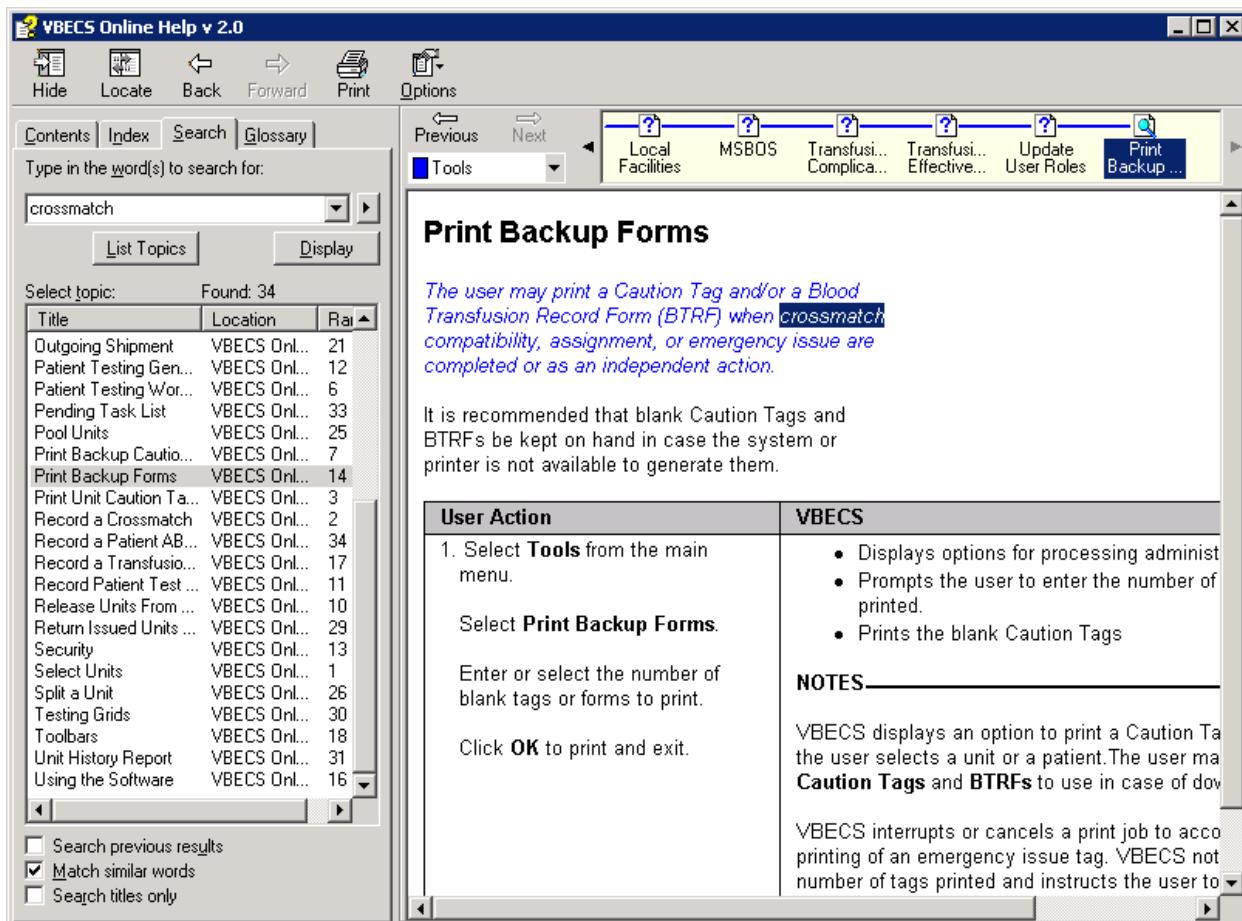

Glossary Tab

Click and drag the scroll bar up or down to find a word or phrase. Click the word or phrase to display a definition at the bottom (Figure 29).

Figure 29: Using the Glossary Tab

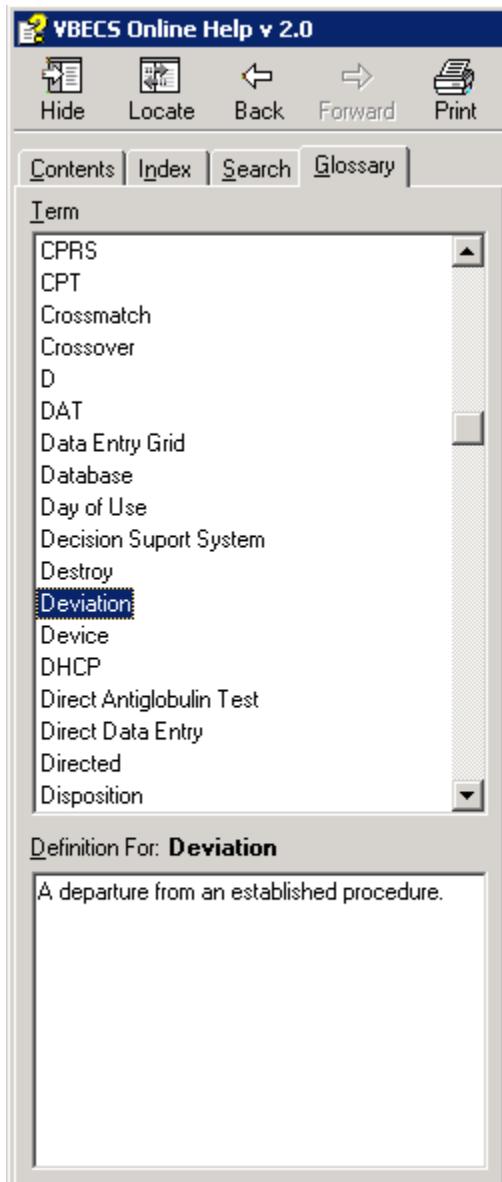

Viewing the VBECS Version Number

In earlier VBECS patch releases, the user documentation referred to the VBECS version in a 4-digit format (e.g., 2.0.0.2 – where 2.0.0 represents the patch version and the last 2 is the patch build number). The build number is used by VBECS Product Support for diagnostic and troubleshooting purposes.

The VBECS version will be represented with only the first three digits (e.g., 2.0.0) and will appear that way in all user documentation to simplify readability. The full 4-digit version can still be found under the **Help, About VBECS** window in VBECS (Figure 30 and Figure 31) and will appear in patch installation guides where build specific files are referenced. The revision letter tracks database-only updates (e.g., blood product table updates, canned comments updates). The revision letter is normally a single alpha character (e.g., C), but can be two characters (e.g., AA, AB, AC) in the unlikely event that more than 25 database updates are made before a code change is implemented. The revision letter starts at A with each new code release which means that when the first database-only update is made, it is incremented to B. The revision letter is then updated by one character in the alphabet for every successive database-only update until a new code release is implemented at which time the revision letter reverts back to A.

The VBECS **Help, About VBECS** window displays the release number (e.g., 1.5.0), the build number (e.g., 1.5.0.4) and the revision letter (e.g., Revision A).

When a version number is incremented, each number set is considered independent of the others. For example, “1.10” follows “1.9,” even though 1.10 is arithmetically smaller than 1.9. The exception to this rule is when a new release is started. The four-digit set is separated by decimal points. The first number is the Major release number. The second number is the Minor release number. The third number is the Patch release number. The fourth number is the Build number.

- At a major release start, the major number is incremented by one, the minor and patch numbers return to zero, and the build number returns to one (e.g., 1.7.0.4 → 2.0.0.1).
- At a minor release start, the major number remains the same, the minor number is incremented by one, the patch number returns to zero, and the build number returns to one (e.g., 1.5.1.2 → 1.6.0.1).
- At a patch release start, the major and minor numbers remain the same, the patch number is incremented by one, and the build number returns to one (e.g., 1.6.0.3 → 1.6.1.1).

Figure 30: About VBECS


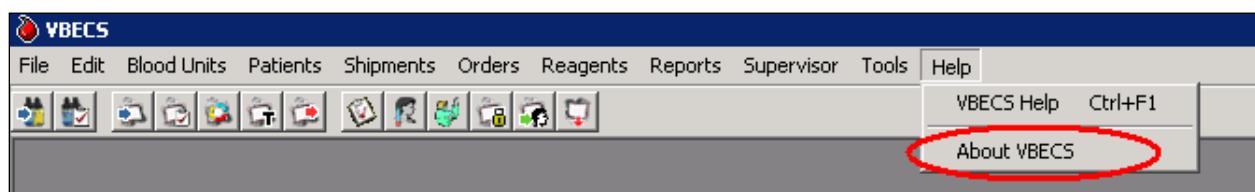


Figure 31: Example of VBECS Version and Build Numbers



VistA Records in VBECS

Historic VistA blood bank records were transferred to VBECS and validated during the database conversion. VBECS handles these records in accordance with system rules:

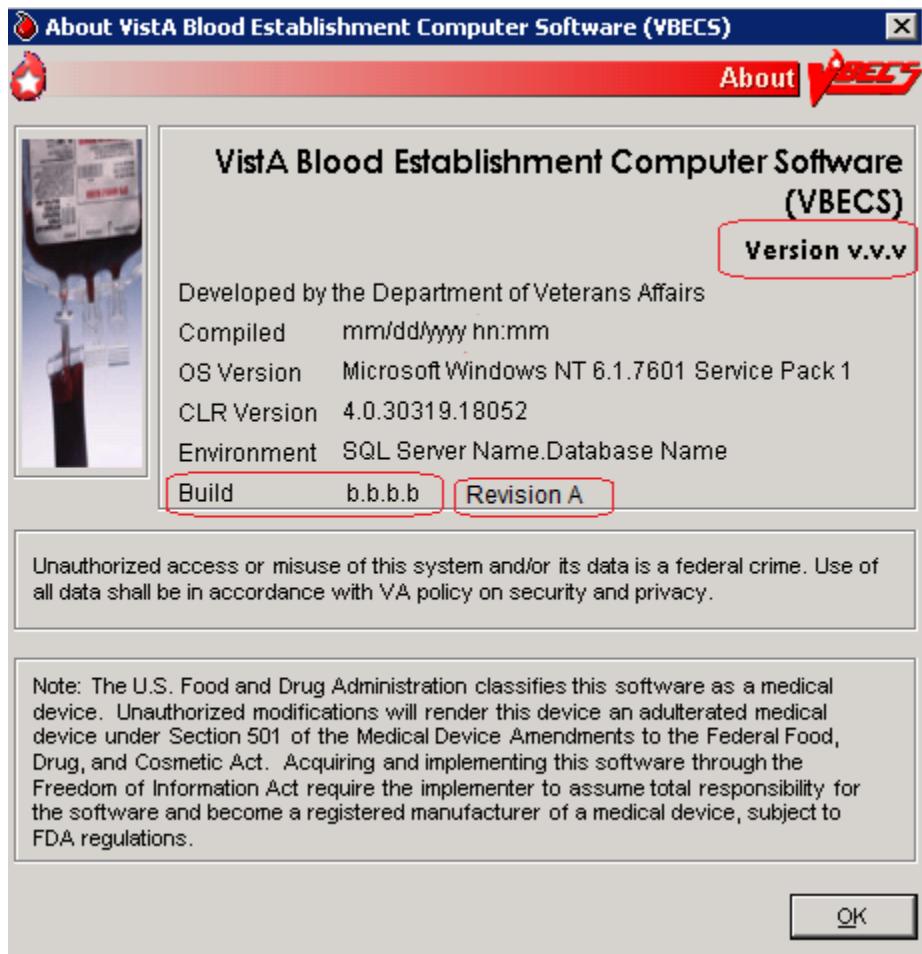

VBECS does not store or display patient name prefixes or suffixes from VistA.



Patient names in VistA that exceed the database standard for File #2 (30 characters) or do not have a first or last name do not convert to VBECS.

- VBECS may display patient information differently from VistA. For example, VBECS displays patient names as “LAST NAME, FIRST NAME, MIDDLE NAME or INITIAL.”
- VBECS will truncate Special Instructions (SIs) text from the VistA database conversion.
- VBECS displays historic ABO/Rh types from the database conversion, but doesn’t use them for issuing blood or for comparing test results.
- SIs from database conversion entries are labeled as such.
- VistA SIs from the database conversion may include component or antigen negative requirements. The VistA SIs remain SIs and are not enforced unless they are entered in the Transfusion Requirements tab in Special Instructions & Transfusion Requirements.
- Patient antibody identifications from VistA are included in the database conversion.
- Antigen negative requirements properly set in VistA will be transferred to antigen negative requirements in VBECS.
- VBECS informs the user when the displayed historical ABO/Rh, transfusion reaction history, and/or antigen typing result from the database conversion. Transfusion reaction types without their associated units are transferred to VBECS.

Accessing the System

Log Into VBECS and VistA

A user may work in only one division at a time. After login¹, a user with access privileges in several divisions may change divisions. There is no option to change divisions in a single-division installation. A user must have access privileges in at least two divisions to change divisions.

Assumptions

- The user's role is configured in both VistA and VBECS. (See Update User Roles.)
- A user must be associated with at least one division to gain access to VBECS.
- VBECS allows access only to valid divisions associated with the user.

Outcome

- A user is logged into the selected division.

Limitations and Restrictions

- VBECS automatically logs a user with access to only one division into that division.

Additional Information

- It is recommended that local policy limit a user to only one VBECS session per terminal.
- The user may log into multiple terminals simultaneously.
- A user with multidivisional access privileges may change divisions during a VBECS session.
- VBECS provides an interface (VistALink) with VistA to allow the user to log onto VistA (not required for most options; see the list of options that use VistALink in Hardware and Infrastructure Architecture). VBECS checks the characters in the Access and Verify Codes and warns when there are fewer than six characters or when a semicolon (;), colon (:), or caret (^) is entered. VBECS sends the entered Access and Verify Codes to VistA. VistA checks the codes and responds to VBECS with a message stating whether the entered Access and Verify Codes are invalid or valid for the user attempting to log on. VBECS warns when VistA does not accept the Access and Verify Codes entered.
- After a number of unsuccessful attempts to log onto VistA (defined in VistA server settings), the VistA server terminates the logon session and displays a message prompting the user to log onto VistA again.
- The user initially logs onto VistA by entering VistA Access and Verify Codes. The system reuses the security credentials (VistA Access and Verify Codes) initially entered in the user's session to perform VistA logoff and logon without user interaction.
- After 15 minutes of inactivity, VBECS disconnects the user from the server. When the user logs back in, he will see the last screen he worked on. After 30 minutes of inactivity, work done after the last save will be lost. The user must log in and start over.

User Roles with Access to This Option

All users

¹ There is a slight difference in terminology between VistA and VBECS: VistA uses “log on” and “logon,” and VBECS uses “log in” and “login.” Therefore, both terms are used throughout this guide. “Log in” and “login” are used generically when referring to both systems at one time.

Log Into VBECS

User Action	VBECS
1. To log into VBECS, double-click the Remote Desktop Connection icon on the Windows desktop.	<ul style="list-style-type: none"> Displays the Remote Desktop Connection screen.
2. Enter or select the name of the VBECS server on the Remote Desktop Connection (Figure 32) screen and click Connect .	<ul style="list-style-type: none"> Displays Security Warning screen.
3. Click OK on the Security Warning screen.	<ul style="list-style-type: none"> Displays Log On to Windows screen.
4. Enter your user name, password and domain (log on to field) on the Log On to Windows screen. Click OK .	<ul style="list-style-type: none"> Displays VBECS server desktop with VBECS.exe, VBECS Test.exe, VBECS Administrator.exe, and VBECS Administrator Test.exe icons. <p>NOTES</p> <p>Users with VBECS Administrator access can select to start VBECS administrator or VBECS.</p>
5. Click on the VBECS.exe or VBECS Test.exe icon to access VBECS.	<ul style="list-style-type: none"> Displays the Select Division screen when the user has access privileges to multiple divisions in a single VBECS installation. <p>NOTES</p> <p>A user may exit without logging into a division.</p>
6. Multidivision users will be prompted to select a division from the Select Division list and click OK , or double-click the division.	<ul style="list-style-type: none"> Displays an error message to the user when logging into VBECS, or switching divisions, if the default printer for the division was removed. <p>NOTES</p> <p>When VBECS is started, the user may have access to more than one division.</p>
7. Continue to the VistA logon screen (Figure 33).	<ul style="list-style-type: none"> Opens the VistA Logon – Authorization screen. The user may log onto VistA or continue and log on as needed. Displays a division login message. <p>NOTES</p> <p>VBECS allows the user to log in without establishing a VistALink connection.</p> <p>The division login message allows the display of data as each user logs into his division. When no login message text is entered, VBECS does not display a login message.</p>
8. Log onto VistA at VBECS startup or at the invocation of any option that uses VistALink when VistALink is not connected.	<ul style="list-style-type: none"> Allows a user to log on by entering VistA Access and Verify Codes, separated by a semicolon (;), in the Access Code data entry field. When a user accesses an option that requires a VistALink connection and the connection becomes unavailable, allows the user to restore the connection. When a reconnection attempt is successful, VBECS closes the connection status window and returns to the desktop. The VistALink Connected icon in the status bar indicates a successful connection. When a reconnection attempt is unsuccessful, VBECS attempts to reconnect to VistALink until the user cancels. <p>NOTES</p> <p>When a user logs into VBECS, the connection to VistA is</p>

User Action	VBECS
	<p>established, and the data in the cache tables are out of date, VBECS retrieves the current data from VistA and saves them. This update process may delay user login.</p> <p>When the VistALink connection is not restorable, VBECS displays a message stating that the requested use cannot be executed because VistALink is unavailable.</p>
9. Enter the VistA Access and Verify Codes.	<ul style="list-style-type: none"> Verifies that the user credentials for the VBECS and VistA Access and Verify Codes belong to the same user. VBECS will attempt to resend queued DSS information when VistALink is available after any user logs in. Emits a tone to verify that the workstation's sound works.
10. Continue working in VBECS in the selected division.	<ul style="list-style-type: none"> Displays the main menu.
11. Users with multidivisional access: to change division, select File . Select Change Division . Continue at Step 2.	<ul style="list-style-type: none"> Ensures that there are no incomplete activities for the user's session. Displays a division login message. Opens the selected division. <p>NOTES</p> <p>A user may log into or change divisions without being logged onto VistA.</p> <p>VBECS displays a message advising the user to close dialogs and transactions with pending data changes, and closes open dialogs with no pending data changes before the user changes division.</p> <p>VBECS automatically logs onto VistA for the selected division.</p>

Figure 32: Remote Desktop Connection Options


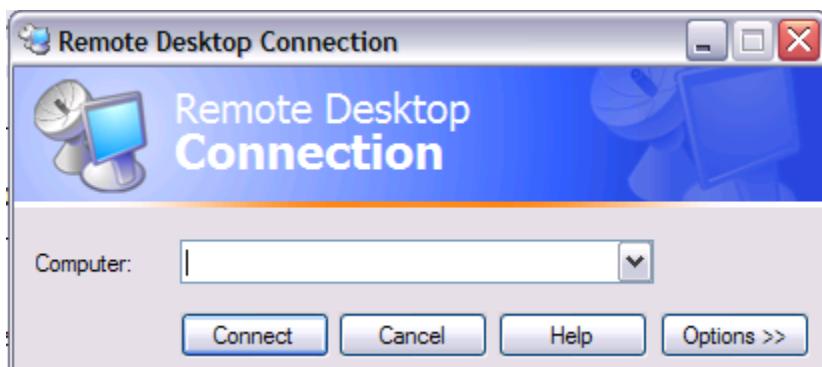


Figure 33: VistA Logon

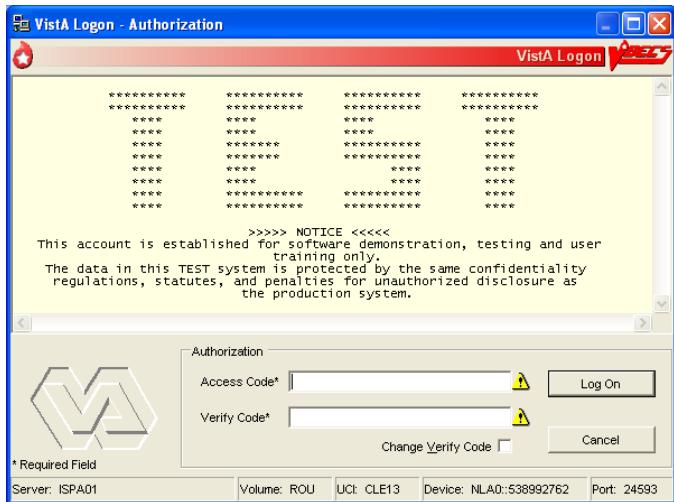

Restore VistA Connection

When a user is logged into a VBECS division and a VistA connection is not in place, VBECS attempts to reconnect.

Restore VistA Connection

User Action	VBECS
1. Attempt to use a VBECS option that requires a connection to VistA for data retrieval.	<ul style="list-style-type: none"> Displays a message that the VistA connection is not available and asks whether the user wishes to attempt to reconnect. Allows the user to enter his VistA logon information.
2. Enter personal logon information or cancel.	<ul style="list-style-type: none"> Attempts to connect using the user's information when logon information is entered. When a connection is achieved, closes the pop-up and takes the user to the selected option. When a successful connection cannot be achieved, does not allow the user to continue in the option. Closes the pop-up when the user clicks Cancel.
3. Continue in the selected option.	

Exit VBECS

A user may shut down the VBECS session.

Shut Down VBECS

User Action	VBECS
1. To shut down VBECS, select File , Exit VBECS . Click Yes at the prompt.	<ul style="list-style-type: none"> Exposes the server desktop.
2. Select Start , Log Off . At the prompt, select Log Off or Cancel .	<ul style="list-style-type: none"> Shuts down the session. Terminates open applications.

Screen Settings



VBECS was designed to operate with the settings described in this section to display complete screens. Failure to maintain these settings may result in partial display of safety-critical information.

Local Machine Screen Resolution




Set the local machine's desktop resolution to 1024 by 768 pixels or higher for optimal display.

To adjust the resolution of the screens that you view:

- 1) Click **Start**, **Settings**, **Control Panel**.
- 2) Click **Appearance and Themes**.
- 3) Click **Change the screen resolution**.
- 4) Select the **Settings** tab (Figure 34).
- 5) Click, hold, and slide the pointer to a screen resolution of 1024 by 768 pixels.
- 6) Click **Apply**.
- 7) Click **OK** to save the setting.

Figure 34: Local Machine Screen Resolution

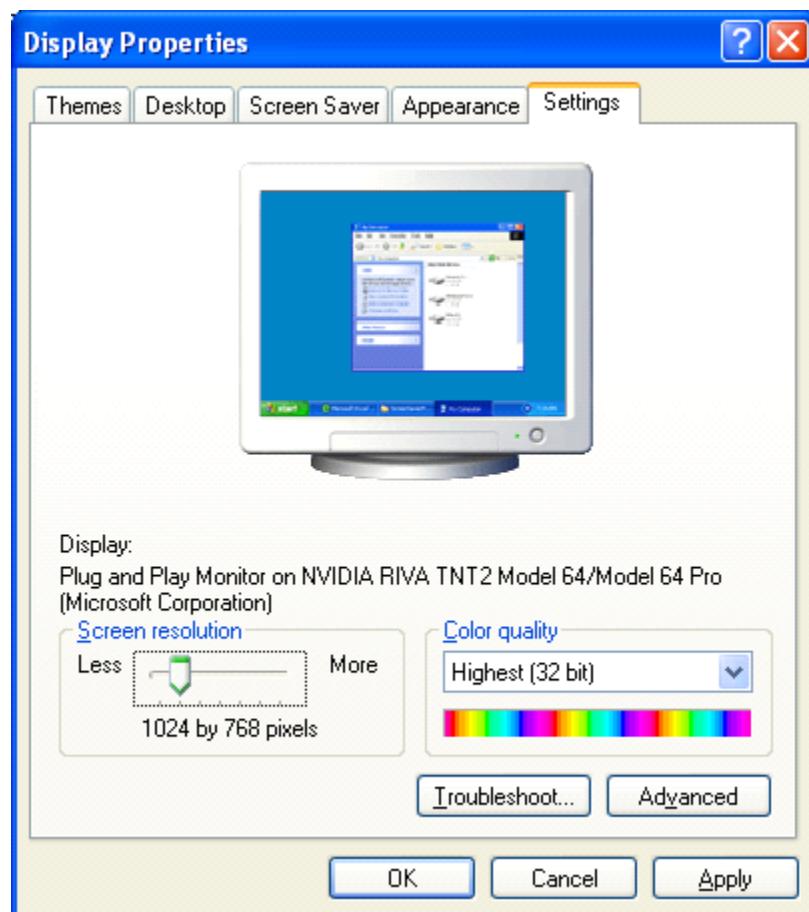

Remote Desktop Connection Screen Size



The Remote Desktop Connection desktop size must be set to “Full Screen.”

To adjust the dimension of the screens that you view:


- 1) Double-click the VBECS (**Remote Desktop Connection icon**) located on the local workstation desktop.
- 2) Click **Options** (Figure 35).

Figure 35: Remote Desktop Connection Options



- 3) Click the **Display** tab (Figure 36).

Figure 36: Remote Desktop Connection Remote Desktop Size

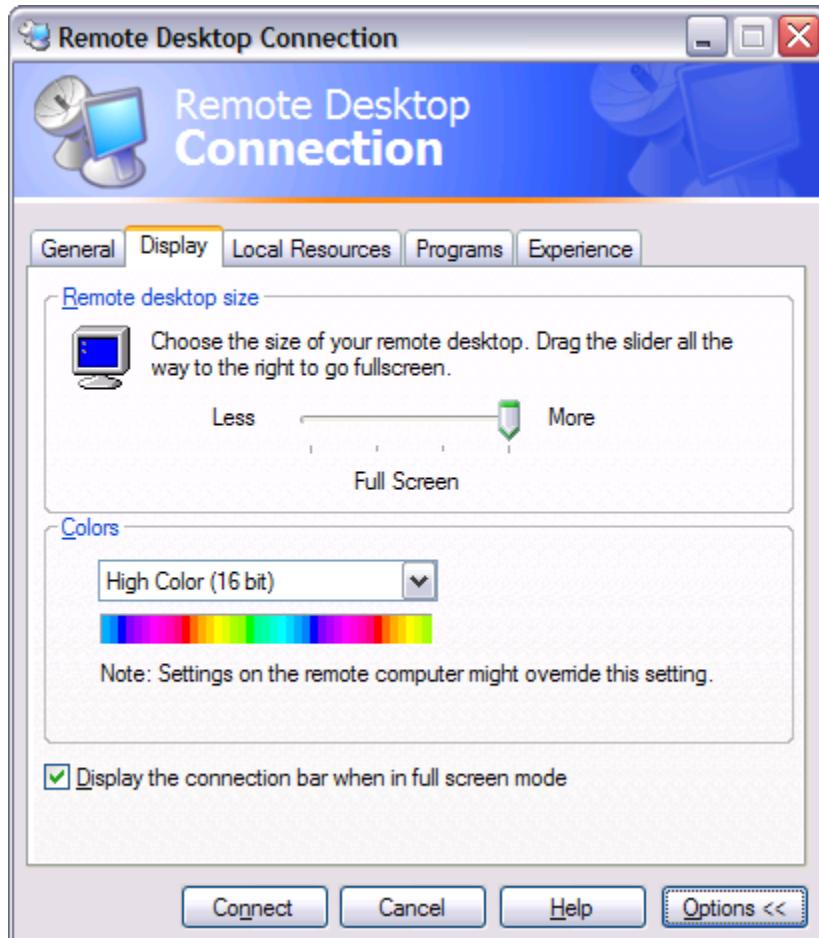

- 4) Click, hold, and slide the pointer to the right until “Full Screen” appears.
- 5) Click on the **General** tab.
- 6) Click **SAVE** to save the setting.

Screen Display

A *window* is the section of the screen that displays the active application or file.

Click the control buttons in the upper-right corner (Figure 37) of an open file or application window to:

- Minimize: reduce the window to appear as a button (rectangle) at the bottom of the screen. Click the button to enlarge the window to its original size.
- Maximize: enlarge a window to fill the entire screen
- Restore: return a window to a smaller size that can be adjusted (Figure 38)
- Close: shut down a document or application quickly




Maximize option windows to display safety-critical information in full.

Figure 37: Resize or Close a Window



Figure 38: Restore a Window

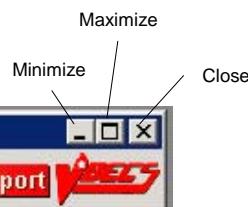

When a window is partially displayed or freezes; or when you move from one window to another, enter information, return to the original window, and the entered information is not available, minimize then maximize the window to refresh it.

To resize an open window:

- Make sure the Maximize symbol is visible.
- Change the width: move the mouse cursor to the left or right window border until the pointer changes into a horizontal double-headed arrow, then drag the border to the right or left.
- Change the height: move the mouse cursor to the top or bottom window border until the pointer changes into a vertical double-headed arrow, then drag the border up or down.
- Change the height and width together: move the mouse cursor to any window corner until the pointer changes into a diagonal double-headed arrow, and then drag the border in any direction.

A maximized window cannot be resized or moved.

To move a window, click the dark blue bar at the top of the window. While depressing the mouse button, drag the window to the desired position.

Select **Help**, **VBECS Help**, or press F1 at any time for screen-specific online help.

The first time a user opens a form [such as the Blood Transfusion Record Form (BTRF)], VBECS displays it in cascade style. When a user opens the form again, VBECS displays it in the size and position last selected.


Forms can be resized and moved in the same way that windows are resized and moved.

Sound

To enable sound on Remote Desktop Connection:

- 1) Double-click  (the **Remote Desktop Connection icon**).
- 2) Click **Options** (Figure 39).
- 3) Click the **Local Resources** tab.
- 4) Select **Bring to this computer** from the Remote computer sound drop-down list.


Figure 39: Remote Desktop Connection Remote Computer Sound

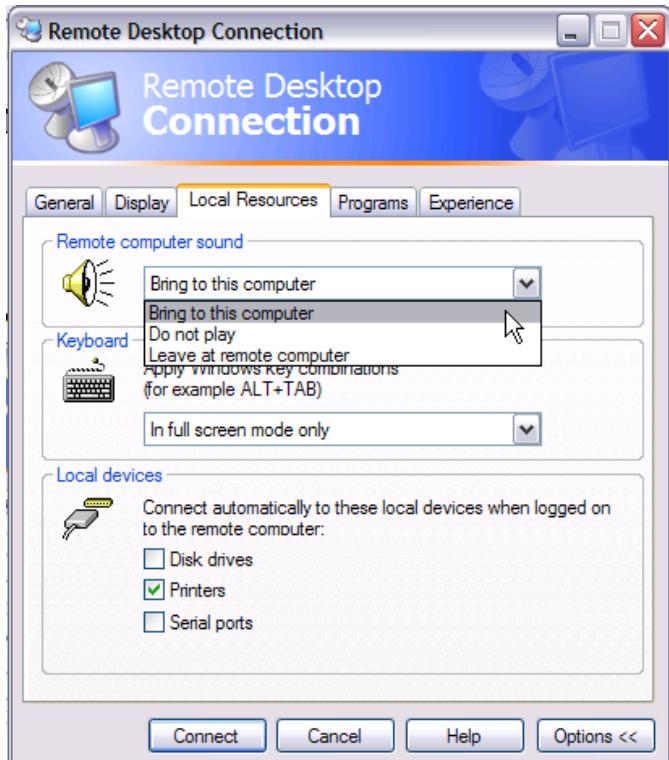

- 5) Click the **Experience** tab.

- 6) Select **LAN** from the **Choose your connection speed to optimize performance** drop-down list (Figure 40).

Figure 40: Remote Desktop Connection Choose Connection Speed

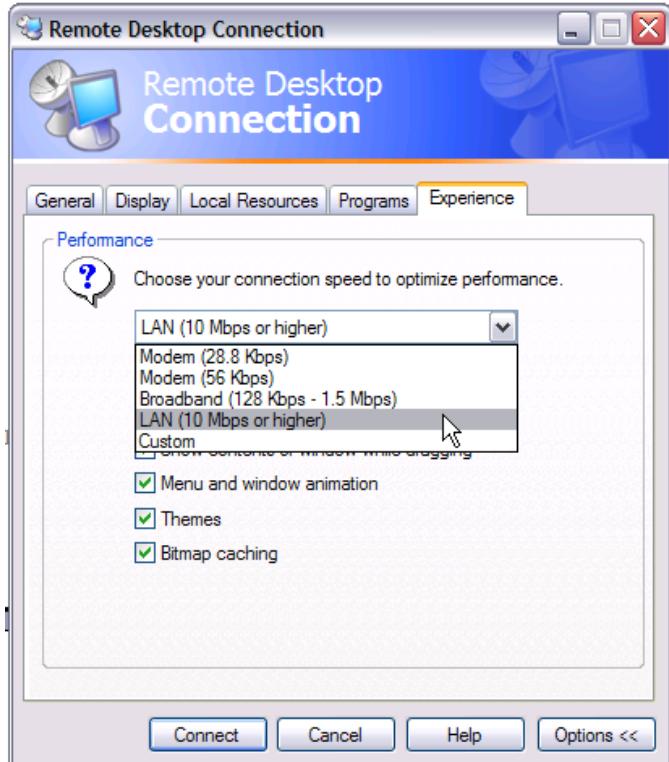

- 7) Click **Connect**.

Volume




Do not change volume settings unless alerts are inaudible (muted) or low.

VBECS uses sound as one way to alert users. To restore the volume setting on your workstation so that alerts are audible:

- 1) Click **Start, Settings, Control Panel**.
- 2) Click **Sounds, Speech, and Audio Devices**.
- 3) Click **Sounds and Audio Devices**.
- 4) Select the **Volume** tab (Figure 41).
- 5) Click, hold, and slide the pointer to adjust the volume.
- 6) Click **OK** to save the setting.

Figure 41: Sounds and Audio Devices: Volume

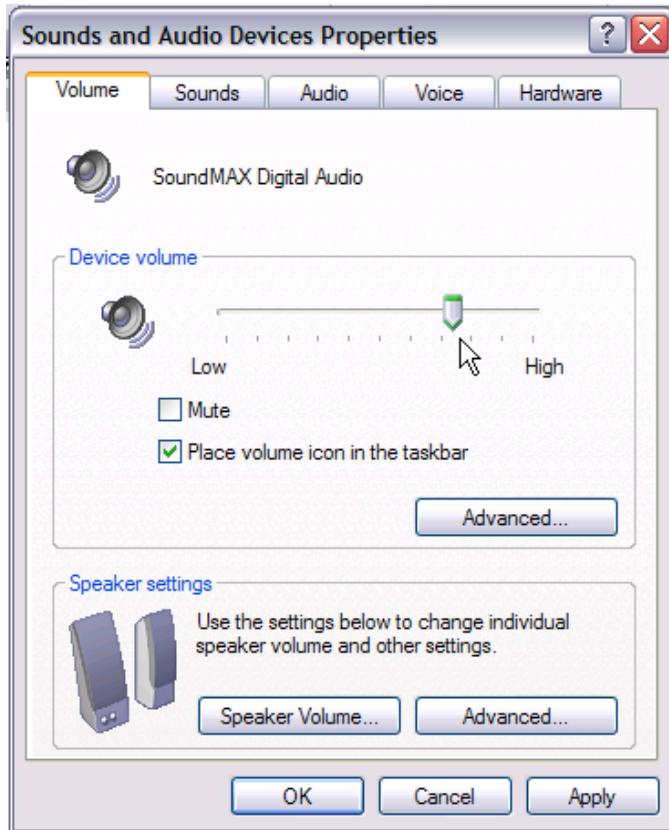

Keyboard Shortcuts for Windows

Using a mouse to navigate VBECS requires fewer steps than using the keyboard. A user who prefers using the keyboard may refer to these keyboard shortcuts. Keys preceded and/or followed by “+” indicate that they must be pressed simultaneously:

Keyboard navigation of VBECS does not accommodate the presentation of error providers or tool tips. A mouse is required to reveal this information as needed.

CPRS and VBECS

- Toggle between **VBECS** and **CPRS** or another application: **Alt + Tab**

Menus

- Display the **Start** menu (use the arrow keys to select an item): **Alt + Home**
- Start VBECS online help: **F1**
- Activate a menu bar option: **F10**
- Open a menu: **Alt + underscored letter**
- Move to and select a menu item: Arrow key + **Enter**
- Open a submenu item: **Shift + letter** (case insensitive)
- Open a drop-down menu: **Alt + down arrow key**
- Close a drop-down menu: **Alt + up arrow key**
- Display a shortcut menu for the selected item: Application key (

Buttons, Check Boxes, and Fields

- Move to and select fields and buttons: **Tab + Enter**
- Select or clear a check box: **Tab**, arrow key (or **Alt + underscored letter**), spacebar
- Select or clear a radio button: **Tab**, arrow key (or **Alt + underscored letter**)

Windows

- Display the **Control** window to restore, move, resize, minimize, maximize, or close the window: **Alt + Delete**
- Switch between a window and a full screen: **Ctrl + Alt + Break**
- Close the current window: **Alt + F4**
- Display the **Windows Security** dialog box to lock the computer, log off, shut down, change a password, and access **Task Manager**: **Ctrl + Alt + End**

Screen Prints

- Place a snapshot of the desktop on the local clipboard and paste it in another application: **Ctrl + Alt + Plus (+ on the numeric keypad; Ctrl + V to paste)**
- Place a snapshot of the active window on the local clipboard and paste it in another application: **Ctrl + Alt + Minus (- on the numeric keypad; Ctrl + V to paste)**

Drop-Down Menus

A *drop-down menu* may appear when a button is clicked or text is selected. Click or hold the mouse pointer over an item in a drop-down menu to reveal other items and submenus. Press an underscored letter (mnemonic) to select a menu item. To view underscored letters:

- 1) Click **Start** at the bottom left of the screen.
- 2) Click **Settings**.
- 3) Click **Control Panel**.
- 4) Double-click **Display**.
- 5) Click the **Appearance** tab.
- 6) Click **Effects**.
- 7) Clear the **Hide underlined letters for keyboard navigation until I press the Alt key** check box.
- 8) Click **OK** to save the setting.
- 9) Click **OK** to close the Display Properties window.
- 10) Close the Control Panel.

A ▶ (submenu arrow) next to a menu item indicates that a submenu is available. An option that is not available to a user is visible but grayed out.

Tool Tips

A *tool tip* is a description or explanation that appears when the mouse pointer is held over (not selected or highlighted) a tool, button, or other object (without clicking a mouse button).

When a tool tip does not appear, wait a few moments, then reposition the mouse pointer over the tool, button, or other object. If necessary, exit and reenter the window.

Tool tips appear with:

- Some control buttons
- Icons in the Patient Information Toolbar

Text enclosed in angle brackets (< >) in tool tips and messages will change to reflect the current activity.

The user may click column headings in some data displays, such as the Pending Task List, Pending Order List, and Select Units, to sort the information in those columns.

Icons and Buttons

VBECS uses several types of images to help you move through the system:

- *Icons* represent objects (information used in an application). The rectangular area in the lower right corner of the screen (system tray) contains icons of programs that always run in the background, such as the clock.
- *Buttons* are icons that represent tools or commands.
- *Radio buttons* are small circles preceding text. Related radio buttons are grouped together. Only one button in a group may be selected: clicking a button turns the other radio buttons off.

Toolbars

A *tool* is a shortcut button, usually shown on a bar near the top of a window. Tools provide quick access to commonly used functions. A *toolbar* displays groups of tools. Toolbars may contain buttons, menus, or combinations of both. (See also Patient Information Toolbar.)

A toolbar with tools for several options (Main Toolbar) appears near the top of the screen.

Click a tool to open a function quickly.

A function may also be opened by pressing Alt and the key for the underscored letter in the tool name.

Main Toolbar

The *Main Toolbar* appears at the top of the screen. Click the icons in the Main Toolbar to access²:

-  Log In Reagents
-  Enter Daily QC Results
-  Incoming Shipment
-  ABO/Rh Confirmation
-  Modify Units
-  Unit Antigen Typing
-  Outgoing Shipment
-  Accept Orders: Accept an Order
-  Maintain Specimen
-  Patient Testing
-  Select Units
-  Issue Blood Components
-  Post-Transfusion Information

² Larger versions of these icons appear throughout VBECS in the corresponding option windows. They cannot be used to access options.

Clicking an icon in the Main Toolbar saves a step: it's the same as clicking the main menu item, and then on the option.

VBECS uses tree views to displays list items. Click a plus sign (+) to expand a list. Click a minus sign (-) to condense a list.

Patient Information Toolbar

The *Patient Information Toolbar* appears in the active window when a patient record is displayed. VBECS displays only the icons appropriate to the option. Click the icons in the Patient Information Toolbar to display:

-  a Blood Availability Report (Blood Availability)
-  a patient's current medication profile (Medication Profile)
-  a patient's recent orders (Recent Orders)
-  a patient's Special Instructions and/or Transfusion Requirements (View Special Instructions and Transfusion Requirements)
-  a patient's recent transfusion history (Recent Transfusions/Issued Units)
-  a patient's transfusion reaction history (Transfusion Reaction History)

Patient Blood Availability

The user views blood products available to a patient. The Blood Availability Report does not include unassigned, restricted blood products. To view unassigned, restricted blood products and other custom reports, see Blood Availability.

Medication Profile

The user views a patient's current medication profile.

Assumptions

- The connection to VistA is active.

Outcome

- VBECS displays a patient's current medications.

Limitations and Restrictions

- None

Additional Information

- The user may access CPRS for additional information on a patient medication order.
- These data are accessible to the user through the Patient Information Toolbar on all patient-specific screens.

User Roles with Access to This Option

All users

Medication Profile

At the user's request, VBECS queries VistA to search for active prescriptions within a user-specified time frame or a default time frame of six months, regardless of division, in a multidivisional database.

User Action	VBECS
1. Select Reports from the main menu. Select Medication Profile . Select a patient, or Select an option that displays the Patient Information Toolbar. Click  to view a patient's current medication list.	<ul style="list-style-type: none">• Lists report names.• Displays the report date range, sections, and compilation criteria.• Displays options for processing patient-related functions.• Compiles the active patient identification data. <p>NOTES —————</p> <p>The user may edit the default start date (six months before the current date) for the patient medication search.</p>
2. Accept or edit the default start and end dates.	<ul style="list-style-type: none">• Displays each medication name status, order number, and issue/start date returned from the VistA query.• When a query does not return any medications, notifies the user that he may repeat the query with a different date range.
3. View the display. When desired, go to Step 1 to repeat the query,	<ul style="list-style-type: none">• Allows the user to print the displayed patient medication profile.
4. Print the profile or click Cancel to exit.	<ul style="list-style-type: none">• Displays and prints the medication information profile.

View Special Instructions and Transfusion Requirements

The user views patient-specific Special Instructions (SIs) and Transfusion Requirements (TRs).

Assumptions

- SI entries were created through Special Instructions & Transfusion Requirements: Enter a Special Instruction or from database conversion.
- TR entries were entered through Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement or Patient Testing: Enter Antibody Identification Results.

Outcome

- VBECS displays a patient's SIs and TRs.

Limitations and Restrictions

- Changing a patient's antigen negative requirement or antibody identified after a unit is crossmatched does not provide an alert at Issue Unit regarding the type of crossmatch (electronic or serologic) performed on a unit available for issue.

Additional Information

- VBECS displays a patient's known antibodies in all divisions in a multidivisional database. VBECS does not restrict viewing SI and TR entries by division; they are linked to the patient's record and not to an order, a specimen, or a unit.
- VBECS does not enforce SIs; VBECS enforces TRs and persistent antigen negative requirements throughout the database.

User Roles with Access to This Option

All users

View Special Instructions and Transfusion Requirements

The user views a patient record. When SIs and TRs are available for the patient, VBECS activates the Patient Information Toolbar icons. VBECS alerts the user when a patient's record includes SIs and TRs.

 It is highly recommended that the user click every active icon when reviewing patient data.	
User Action 1. Select Patients from the main menu. Select Special Instructions & Transfusion Requirements , or Click  in the Patient Information Toolbar, when it is active.	VBECS <ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the option to select a patient, when one is not already selected. <p>NOTES _____</p> <p> When a patient's VBECS record includes an active TR, SI, or antigen negative requirement, VBECS emits an audible alert and notifies the user by enabling the icon.</p>
2. Select a patient and click OK . View the data.	<ul style="list-style-type: none">• Displays the data. <p>NOTES _____</p>

User Action	VBECS
Click OK to continue.	<p>VBECS displays active TRs to all divisions within a multidivisional database.</p> <p>VBECS alerts the user when an active TR, SI, or antigen negative requirement is in a patient's record.</p> <p>VBECS displays component requirements to all divisions in a multidivisional database. No division allows issue until requirements are met.</p> <p>When VBECS displays the Patient Information Toolbar, it calculates the overall compatibility percentage of units to be screened to find compatible units.</p> <p>The compatibility percentage displayed is based on the following calculation:</p> <p>Type Specific ABO % x Antigen compatibility % as defined in Antibodies equals % compatible.</p> <p>For example, if a patient is A Positive and has both Anti-E (compatibility percentage = 75%) and Anti-K (compatibility percentage = 91%) antigen negative requirements, the percentage of the population considered compatible for this patient would be .45 x .75 x .91 = .30 (rounded). Rh type of the patient is not used in the calculation when the patient does not have an Anti-D as Rh Pos or Rh Neg units may be screened.</p> <p>The blood type percentages in VBECS are as follows and are general percentages for the American population: A is 45%, B is 15%, AB is 6% and O is 45%.</p>
3. Select another patient and return to Step 2, if desired.	
4. Click OK to exit.	

View Recent Orders, Recent Transfusions/Issued Units, and Transfusion Reaction History

The user views portions of a patient's data previously recorded in VBECS in three formats: recent orders, recent transfusions/issued units, and transfusion reaction history.

Assumptions

- A VBECS record must be established for a patient in this database in at least one division, if multidivisional.

Outcome

- The user views:
 - A patient's active order data (recent orders).
 - A patient's transfusions and currently issued units (recent transfusion history).
 - A snapshot of the patient's VBECS clinical record (transfusion reaction history).

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option

All users

Patient Information Toolbar: Recent Orders

The user views current (pending and active) order data (including reflex tests), which include clinical data from all divisions within a multidivisional database.

User Action	VBECS
1. In the Patient Information Toolbar, click  to view the patient's recent orders.	<ul style="list-style-type: none">• Allows the user to click an icon to view a patient's recent orders. <p>NOTES _____</p> <p>When the user selects , VBECS displays all orders associated with the selected patient, including open orders and orders placed and processed within the last 14 days.</p> <p>Diagnostic Test orders sort and display by the date the order was placed in CPRS. Component orders sort and display by the date the order was accepted in VBECS.</p> <p>Orders with a status of Canceled, Completed, Expired or Filled are not displayed in this view. See Order History Report.</p>
2. View the data and click Close to exit.	

Patient Information Toolbar: Recent Transfusions/Issued Units

The user views transfused and currently issued units for a patient.

User Action	VBECS
1. In the Patient Information Toolbar, click  to view the recent transfusion history.	<ul style="list-style-type: none">Allows the user to click an icon to view a patient's recent transfusion history. <p>NOTES _____</p> <p>When the user selects , VBECS displays blood units issued and transfused within the last seven days in reverse chronological order. The date and time issued or the location issued to does not display.</p>
2. View the data and exit.	

Patient Information Toolbar: Transfusion Reaction History

The user views a snapshot of the patient's transfusion reaction record.

User Action	VBECS
1. In the Patient Information Toolbar, click  to view the transfusion reaction history.	<ul style="list-style-type: none">Allows the user to click an icon to view a patient's transfusion reaction history. <p>NOTES _____</p> <p>When the user selects , VBECS displays all transfusion reaction entries in the database.</p>
2. View the data and exit.	

Main Status Bar

The following icons are displayed in the main status bar:

- ⌚ The Order Alert icon blinks three times and remains visible when orders are added to the system. Hold the mouse pointer over the icon to see order information. Click the associated button to open the Pending Order List.
- ❗ The Patient Alert icon flashes three times and remains displayed when patient information is updated. Hold the mouse pointer over the icon to see what type of patient change occurred. Click the associated button to view updated patient information.
- 📶 The VistALink icon indicates that the connection to VistA is active.
- 🚫 The VistALink Failure icon indicates that the connection to VistA is severed.

Throughout VBECS

- ✖ The Delete/Invalidate icon appears on buttons that cause entries to be deleted or invalidated in the database.
- Click to select all items in a list.
- Click to clear all items in list.
-  Click the **find** button to display a list of items for selection.
-  Click the **Print** icon's associated button to print a document.
-  Click the **Save** icon's associated button to update the database.
-  The Scanner icon indicates which screen prompts permit scanned entries.
-  Click the **Search** icon's associated button to search for entries meeting search criteria entered.
-  Click the **Order Details** icon's associated button to display the full details of a selected order.
-  The Input Error icon indicates that a user made an invalid entry at a screen prompt. The icon appears and flashes next to the invalid entry. Hold the mouse pointer over the icon to see a message explaining the problem. Enter valid information to make the icon disappear.

Other VBECS Functions

Collecting Workload In some options, VBECS collects workload (statistical data used to calculate productivity) when the user saves data:

- ABO/Rh Confirmation
- Discard or Quarantine
- Document ABO Incompatible Transfusion
- Edit Unit Information
- Enter Daily QC Results
- Finalize/Print TRW
- Free Directed Unit For Crossover
- Invalidate Test Results

- Issue Blood Components
- Modify Units
- Modify Units: Pool Units
- Modify Units: Split a Unit
- Outgoing Shipment
- Patient Testing: Enter Antibody Identification Results
- Patient Testing: Record a Patient Antigen Typing
- Patient Testing: Record a Transfusion Reaction Workup
- Patient Testing: Record Patient Test Results
- Post-Transfusion Information
- Processing Incoming Shipment
- Release Units From Patient Assignment
- Remove Final Status
- Return Issued Units To Blood Bank
- Select Units
- Unit Antigen Typing

Configuring Site Parameters

Configuration must be done as part of the VBECS installation.

Update User Roles

The user reviews and changes the user roles and associated security levels that determine which options and functions are available to a user in the division.

Assumptions

- The user is logged into VBECS in the division of the users to be edited.
- A VBECS Application Administrator added the users to the VBECS database and assigned a division-specific role to each user.

Outcome

- A user's role code changed.

Limitations and Restrictions

- One role at a time may be assigned to a user at a division. The role may be reset.
- VBECS does not allow users to edit or delete the sole Administrator/Supervisor.

Additional Information

- An inactivated user remains in the VBECS database.
- There must be at least one Administrator/Supervisor in the division to maintain users.
- A VBECS user ID is the user's NT Logon ID.
- One or more users may be selected for batch update.
- A VBECS Application Administrator activates and inactivates VBECS users.
- A user may be defined at multiple divisions within a multidivisional database. The settings for a user may differ from division to division and are mutually exclusive.
- See *VistA Blood Establishment Computer Software (VBECS) Technical Manual-Security Guide*, VBECS Maintenance Operations, Configure Users section for Active Directory tool installation.
- See *VistA Blood Establishment Computer Software (VBECS) Technical Manual-Security Guide* to add a user to a mail group that receives messages when VistALink errors occur.

User Roles with Access to This Option

- Administrator/Supervisor
➤➤➤ Lead Technologist

Update User Roles

➤➤➤➤➤ An Administrator/Supervisor maintains assigned user roles, including associated levels of system access and the ability to process system overrides. These roles determine the available options and the ability to override certain warnings.

It is recommended that only one staff member be assigned the role of Administrator/Supervisor: when more than one Administrator/Supervisor updates roles simultaneously, changes made by one user are not visible to the other.

One of six user roles may be assigned to users, by division:

- Level 1: Technologist
- Level 2: Enhanced Technologist
- Level 3: Lead Technologist
- Level 4: Traditional Supervisor
- Level 5: Enhanced Supervisor
- Level 6: Administrator/Supervisor

The functions associated with each user role are described in Table 1: User Roles and Functions.

User Action	VBECS
1. Select Tools from the main menu. Select Update User Roles .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Lists active users and associated information, including user name, ID, role, initials and email address.• Displays an option to print a VBECS User Report. <p>NOTES —————</p> <p>Click Print at any time in this process to print a VBECS User Report.</p>
2. Click one or more check boxes to select users.	<ul style="list-style-type: none">• Allows the assignment of a security level to one or more users at a time. <p>NOTES —————</p> <p>The role that a user selects is assigned to all selected users.</p>
Select a user role from the drop-down list in the Use Role field.	<ul style="list-style-type: none">• Requests confirmation of the names of users to be changed and the role to be assigned.
3. Click OK .	
4. Click Yes to confirm the save or No to change the users selected.	<ul style="list-style-type: none">• Creates an Audit Trail Report entry. <p>NOTES —————</p> <p>See the sample VBECS User Report.</p>
5. Repeat Steps 2–4 to change the roles of additional users.	
6. Click Cancel to exit.	

VBECS User Report

The user may generate a division report that includes the names, VBECS user IDs, and roles of active users. ➤➤➤ A Lead Technologist may generate a report of current VBECS users, by division, that includes the users' name, user ID, division, and role.

Component Classes

The user configures the component class parameters for the selected division.

Assumptions

- The connection to VistA is active.
- The WHOLE BLOOD and RED BLOOD CELLS component classes are deployed as requiring current specimen ABO/Rh results before issue and may not be edited.
- The FFP, PLT, CRYO and OTHER component classes do not have a specimen requirement defined at deployment.

Outcome

- Divisions have local settings for all component classes.

Limitations and Restrictions

- The user may not change the specimen requirement for a component class when there are incomplete orders for that component class.
- There is no retrospective data entry or user selection for this option.

Additional Information

- None

User Roles with Access to This Option

►►►►► Administrator/Supervisor


Component Classes

The user configures the parameters for the division's component classes. These parameters are applied to all blood units that descend from this class. Component classes may be reset when business practices affecting the use of component classes change at a division.

User Action	VBECS
1. Select Tools from the main menu. Select Component Classes (Figure 42).	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Lists component classes.• Displays options for configuring division component classes:<ul style="list-style-type: none">◦ Require Current Specimen ABO/Rh Results?◦ Maximum Transfusion Time.• Displays options for defining inappropriate request indicators. <p>NOTES —————</p> <p>See Appendix B: Table 17: CPRS Orderable Blood Components (Component Classes) Mapped to ICCBBA Component Classes.</p>
2. Select a class from the Select Component Class list. Click OK to accept the list, or click the Require Current Specimen ABO/Rh Results? check box, when appropriate.	<ul style="list-style-type: none">• Requires current specimen ABO/Rh results for issue.• Defines the maximum transfusion times.• Defines inappropriate request indicators.• Allows the user to edit this parameter for component classes that do not contain only red blood cells.• Saves the component to the division list. <p>NOTES —————</p>

User Action	VBECS
	<p>Component class configuration is required. Component classes: RED BLOOD CELLS, FFP, CRYO, PLT, WHOLE BLOOD, and OTHER.</p>
3. Edit the maximum transfusion time.	<p>NOTES</p> <p>VBECS records the maximum transfusion time from the transfusion start to the transfusion end time in minutes with a default of 270 minutes. The user may edit the maximum transfusion time to a minimum time of 5 minutes and a maximum time of 999 minutes.</p> <p>VBECS uses the maximum transfusion times in minutes to determine whether the time used to transfuse a unit of blood is inappropriate for a class.</p>
4. Define the inappropriate request indicators. Click the find button to search for a new laboratory test name. Enter a partial laboratory test name and click Search . Select a laboratory test name and click OK to add it. Click the Show Inactive Inappropriate Request Indicators check box, if desired.	<ul style="list-style-type: none"> • Facilitates data entry of one or more inappropriate request indicators: <ul style="list-style-type: none"> ◦ Laboratory Test Name ◦ Specimen Type ◦ Threshold Result • Uses VistALink to retrieve the laboratory test name and specimen type indicators. The user enters the threshold result indicator. • Displays the entry and allows the user to edit it. • Displays inactive tests and allows the user to reactivate them. <p>NOTES</p> <p>VBECS uses the laboratory test name as the only input parameter to search for all possible matches and return the Laboratory Test Name, Specimen Type(s), and VistA Lab Test Internal Entry Number (IEN).</p> <p>A numeric threshold result must include "<" or ">" and a numerical result, e.g., ">8.0."</p> <p>The settings here indicate predetermined levels used to determine the appropriateness of a component order. Settings are different for each component class. The example of a ">8.0" hemoglobin level is set according to local policy.</p> <p>An example of an inappropriate request indicator:</p> <ul style="list-style-type: none"> • The test name is "HGB." • The specimen is entered as "Blood." • The threshold result is set to ">8.0."
5. Click Add to add the VBECS component class site parameters to the database.	<ul style="list-style-type: none"> • Displays the component class parameters and allows the user to edit them before saving.
6. Repeat Steps 2–6 for each component class to be edited.	
7. Click OK to save. Click Yes to confirm changes and exit.	<ul style="list-style-type: none"> • Updates the database and captures changes to previously defined parameters for inclusion in an Audit Trail Report.

Figure 42: Component Classes

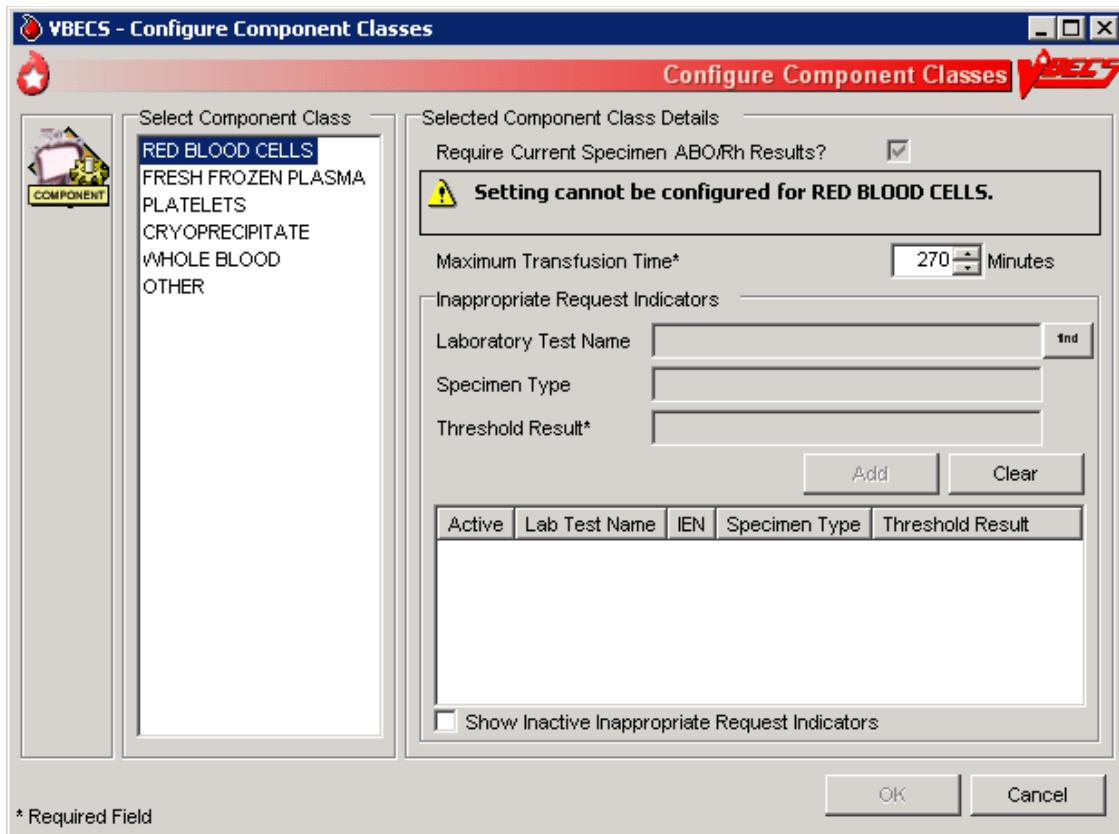

Configure Daily QC

The user sets up the routine reagent racks and QC template for the division before performing daily reagent QC.

Assumptions

- The division is “full service.”
- In the first use of this option, there are no defaults or data in place until the user selects the type of QC template (commercial or non-commercial).

Outcome

- A testing template is available for Enter Daily QC Results; reagent racks and identifiers were created.
- Minimum reaction strengths for each reagent type are defined and used by Enter Daily QC Results.
- A historical record for each rack will be maintained.
- Rack names created are available for use by other options.

Limitations and Restrictions

- This option pertains only to routine reagents.
- When there are partially or currently QC'd racks in the system, the user may not change the QC template.

Additional Information

- This option is disabled for transfusion-only facilities.
- The naming convention for the active template is not editable. When fewer racks are indicated, the system will inactivate existing rack names to keep the list in order without skipping a number or letter. Rack names are reusable.
- The daily use field for QC Rack template entries exported as REQUIRED for Daily Use are not editable by the user. The daily use field for entries exported with the daily use as OPTIONAL may be changed to be REQUIRED and remain editable by the user at any time.
- If Screening Cell 3 (SC3) and/or (SC4) are set to “required” as part of the configuration of the antibody screening cell set they are enabled on the antibody screen data grid.
- The user must identify an enhancement media reagent type as the primary. In addition, the user must identify the Anti-Human Globulin serum used during daily quality control testing of the screening cells and primary enhancement media.
- VBECS does not evaluate maximum strength of reaction for daily QC, as required by CAP. Local policy must define maximum strength of reaction.

User Roles with Access to This Option


►►►► Enhanced Supervisor

Configure Daily QC

The user sets up the daily reagent QC template for reagent rack testing for the division. In the first use of this option, all fields default to blank with no template preferred. Additional tests are available and may be included in the daily QC testing template in accordance with local practice. When a user changes the active template, there is no carryover of field settings to the new template. All are set to the default and configuration must be completed.

User Action	VBECS
1. Select Reagents from the main menu. Select Configure Daily QC .	<ul style="list-style-type: none"> • Displays options for processing reagents. • Displays an option to configure online Daily Reagent Rack Quality Control.
2. Click the Commercial or Non-Commercial radio button (Figure 43).	<ul style="list-style-type: none"> • Lists the two system-supplied Daily QC templates. <ul style="list-style-type: none"> ◦ Commercial Daily QC ◦ Non-commercial Daily QC • When there are partially or currently QC'd racks in the system, warns the user to wait until QC expires (passes the daily alert time) and then make the changes.
3. Set the daily alert time.	<ul style="list-style-type: none"> • Asks the user to set the time the daily alert to all users at the site will begin. <p>NOTES</p> <p>The user indicates the local time that defines the start of the 24-hour period (the QC day) during which a rack is considered QC'd. VBECS alerts users when a rack was not performed after this local time passed.</p> <p>The user may test a rack no earlier than the alert time for the next 24 hours.</p>
4. Enter the number of racks at the site. Select a naming convention.	<ul style="list-style-type: none"> • Asks the user to set the number of racks to be QC'd daily at the division. • Asks the user to set the naming convention for racks. <p>NOTES</p> <p>A user may indicate 1–26 racks per site and names them according to the convention selected (A–Z or 1–26).</p>
5. Accept or edit the default settings of the template. Select primary and secondary enhancement medium options from the drop-down menus.	<ul style="list-style-type: none"> • Lists the default reagent types for inclusion in the daily reagent rack quality control and expected minimum reactions for each test type. • Displays changes and edits for review before saving. <p>NOTES</p> <p>One setting covers all racks at the site.</p> <p>Reagent types selected display on the rack lot number template and the reaction result entry worksheet for each rack in Enter Daily QC Results with some exceptions: gel kit, selected AHG, and automated testing lots numbers.</p> <p>When previously entered parameters for daily QC rack definitions are changed, VBECS requires a comment before the user saves the update.</p>
6. Click OK to save.	

Figure 43: Commercial and Non-Commercial Radio Buttons

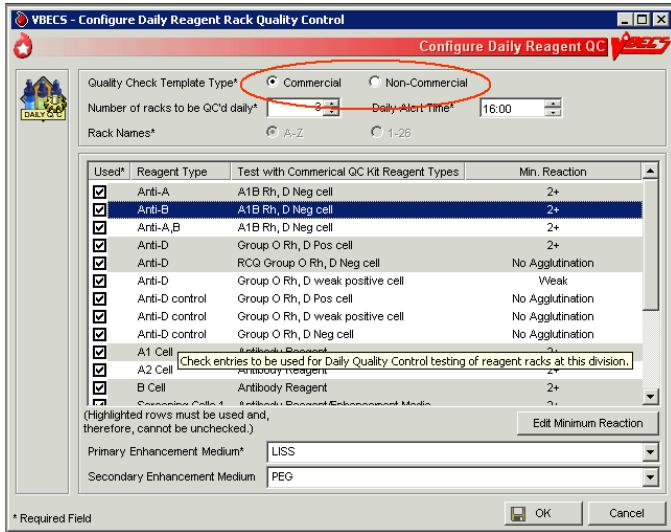

Configure Division

Configure Division

A user defines site parameters that affect how VBECS behaves at a specific division. VBECS may be configured to function differently at each division within a consolidated database. Site parameters may be redefined.

Assumptions

- The site administrator defines the VBECS configuration, adds users, and assigns privileges for each division so that they mirror the associated VistA database.

Outcome

- Functional parameters are set for divisions that use the VistA database.

Limitations and Restrictions

- Changes made to the division setting are not applied to a user until they re-connect to VBECS after the change is made. It is highly recommended that all users log off VBECS prior to any configuration setting change.

*It is recommended that users configure a division before configuring testing (click the **Configure Testing** tab).*

Additional Information

- To incorporate changes in division codes and names from VistA to VBECS, and to change time zones, see *VistA Blood Establishment Computer Software (VBECS) Technical Manual-Security Guide*.
- VBECS records configuration changes (such as time zone) on the Audit Trail Report.

User Roles with Access to This Option

►►►►► Administrator/Supervisor


Configure Division

Once VBECS is configured to identify all of the divisions, the user sets up critical fields required for routine use of VBECS.

User Action	VBECS
1. Select Tools from the main menu. Select Configure Division (Figure 44).	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays the Configure Division tab with fields for entering division-specific blood product information.
2. Enter the ICCBBA Registration Number, the medical director's name (or select it from the drop-down list), and the maximum number of days until the specimen expires in their respective fields.	<p>NOTES —————</p> <p>An ICCBBA Registration Number contains an alphabetical character followed by four numeric characters.</p> <p>A user may enter the name of the blood bank medical director or acting medical director to appear on reports. The user may edit the name and title but may select only one name at a time for</p>

User Action	VBECS
	<p>use. VBECS saves these edits for retrieval.</p> <p>The user may edit the default setting for the maximum number of specimen expiration days from "3" to any value from "3" to "93." This setting allows the user to extend the expiration date of any one specimen to the maximum date indicated. This does not change the normal specimen expiration date of all specimens received.</p>
3. Click OK to save and click OK again to confirm the save and exit, or click another tab.	<ul style="list-style-type: none"> • Displays the information entered. • Requests confirmation to update the database. • Updates the database.

Figure 44: Configure Division Tabs

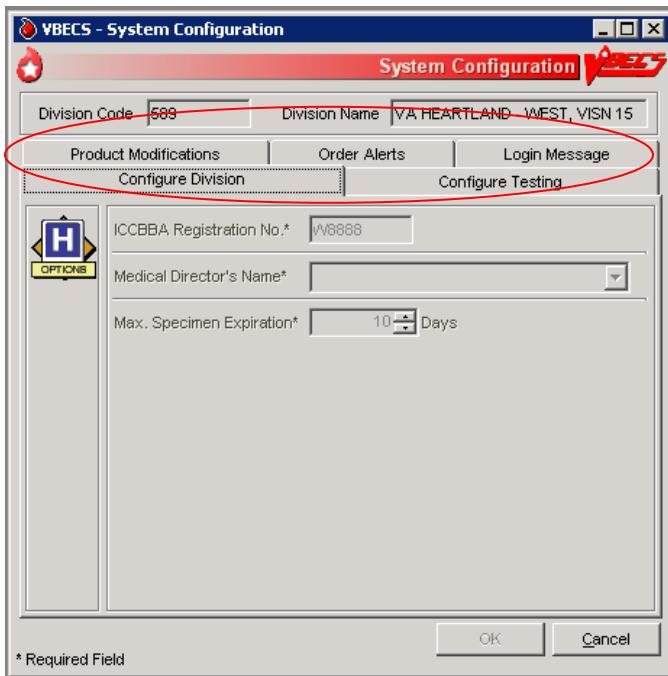

Configure Testing

The user sets testing preferences according to a division's medical policies.

Assumptions

- A division's worksheets affected by the option are not open or incomplete at the time of configuration.

Outcome

- Division testing preferences are configured.

Limitations and Restrictions

- None

*It is recommended that users configure testing before modifying products (click the **Product Modifications** tab).*

Additional Information

- None

User Roles with Access to This Option


>>>> Administrator/Supervisor

Configure Testing

The user sets the system to display testing grids to the users in their division. The user determines whether an autocontrol is used and whether a Direct Antiglobulin Test (DAT) is added as part of the routine Antibody Screen Test (ABS) described in Patient Testing: Record a Patient Antibody Screen.

User Action	VBECS
1. Select Tools from the main menu. Select Configure Division .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays the Configure Division tab with fields for entering division-specific testing information.
2. Click the Configure Testing tab (Figure 45).	<ul style="list-style-type: none">• Displays the Configure Testing tab with fields for entering division-specific blood product information.• Displays options to set standard operating procedures for a division's laboratory testing.
3. Click the appropriate radio buttons to: Enable Electronic Crossmatch Determine antibody screen (display antibody screen only, or with a DAT or Autocontrol) Determine the default test grid displayed for ABS and XM Testing Phases	<ul style="list-style-type: none">• Displays default configuration settings of:<ul style="list-style-type: none">○ Electronic Crossmatch: Disabled○ Antibody Screen: Antibody Screen Only○ ABS and XM Testing Phases: All Phases• Changes to the configuration settings are captured for the Audit Trail Report.
4. Click OK to save and click OK again to confirm the save and exit, or click another tab.	<ul style="list-style-type: none">• Displays the appropriate testing grids for subsequent patient orders.

Figure 45: Configure Testing Tab

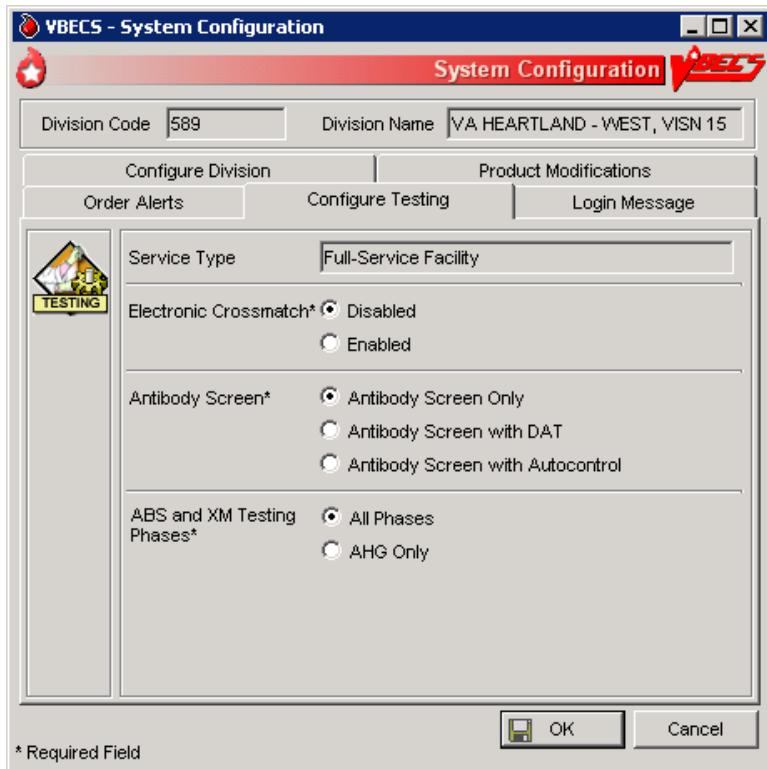

Product Modifications

The user enables the types of modifications made at a division.

Assumptions

- An ICCBBA number was entered for the division.
- VBECS locks a user out when the user attempts to change the configuration and there are open or partially completed worksheets or processes in the division.

Outcome

- This option dictates the available modification types displayed in Modify Units: Pool Units, Modify Units: Split a Unit, and Modify Units.

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option

►►►►► Administrator/Supervisor


Product Modifications

The user enables specific modification type procedures performed in the blood bank in the hospital division. The user may enable one, all, or some of the modification processes. Once the types of modifications are enabled, individual blood unit modification is controlled by the rules of modification truth tables. VBECS is exported with all modification types disabled. VBECS allows for full customization upon deployment as well as changes to the enabled modifications when a division's capabilities and procedures change.

User Action	VBECS
1. Select Tools from the main menu. Select Configure Division .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays the Configure Division tab with fields for entering division-specific testing information.
2. Click the Product Modifications tab (Figure 46).	<ul style="list-style-type: none">• Lists valid modification procedures. <p>NOTES _____</p> <p>A user may enable one, more, or all modification types.</p>
3. Click one or more check boxes to enable modification types. Enter the processing cost associated with each modification procedure.	<ul style="list-style-type: none">• Lists modification procedures. <p>NOTES _____</p> <p>Modification procedures:</p> <ul style="list-style-type: none">• Thaw• Pool• Thaw/Pool Cryo• Split/Divide• Irradiate• Leukoreduce

User Action	VBECS
	<ul style="list-style-type: none"> • Volume reduce • Wash • Rejuvenate • Freeze • Deglycerolize <p>The modification process fee is added to the base unit cost of the first target unit and special testing costs, creating a patient bill for a transfused blood product, not including on-site antigen typing performed. The modification process fee may be U.S. \$0.00 to \$9,999.99. VBECS will generate an audit trail report entry when this field is filled or changed.</p>
4. Click OK to save and click OK again to confirm the save and exit, or click another tab.	Updates the database and captures changes to any previously defined parameters for inclusion in the Audit Trail Report.

Figure 46: Product Modifications Tab

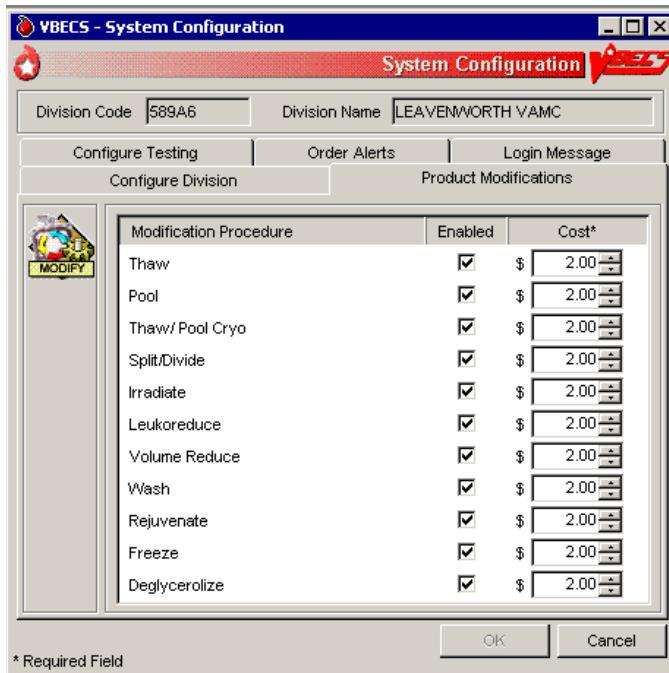

Order Alerts

The user configures VistA alerts that appear at login.

Assumptions

- None

Outcome

- User alerts of physician orders for blood bank diagnostic tests and blood components are configured.

Limitations and Restrictions

- This option does not create a mechanism to accept or reject orders or specimens or define reports for viewing orders.
- The user may select the option to display printer alerts only when the designated printer is configured and available.
- The printer must be configured by the System Administrator.

Additional Information

- A user may review current user alert configurations.
- VBECS uses other alerts (patient update alerts, VistALink status alerts) that do not require configuration.
- The user may set VBECS to automatically print an order alert on the designated printer when VBECS acknowledges the order. This is an alert only and contains limited information from the CPRS order. The CPRS can be configured to print a copy of the order on a printer in the blood bank.

User Roles with Access to This Option

►►►►► Administrator/Supervisor


Order Alerts

The user configures diagnostic test or component order alerts, and turns them on or off. The alerts notify users when new orders are entered, updated, or canceled in VistA.

User Action	VBECS
1. Select Tools from the main menu. Select Configure Division.	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays configuration tabs.
2. Click the Order Alerts tab (Figure 47). Click the Printer Alerts On and/or the Icon Alerts On check box. Select the order alert type from the drop-down menus in the Order Alert Type fields.	<ul style="list-style-type: none">• Displays the existing order alert configuration and allows the user to edit it.• Requires the user to select a printer name and to define which orders to print.• Displays the selected printer, the selected order alert type, and the number of hours (default: 48 hours) entered for review. <p>NOTES</p> <p>VBECS lists local and network printers.</p> <p>The user must select an order alert type:</p> <ul style="list-style-type: none">• All orders (orders received by VBECS)• Hours limit to fill order (orders received by VBECS needed in

User Action	VBECS
	<p>the next <i>n</i> hours)</p> <ul style="list-style-type: none"> Specimen OK or N/A (orders where no specimen is required)
3. Enter or select a refresh rate (the interval between searches for new orders) for the order alert in the Refresh Rate field.	<ul style="list-style-type: none"> Requires the user to enter a refresh rate (default: 60 seconds) (Icon Alerts On). Displays the updated refresh rate. Displays the selected order alert type.
4. Click OK to save	<ul style="list-style-type: none"> Saves the user alert configuration settings.
5. Click Yes to confirm the changes and exit.	

Figure 47: Order Alerts Tab

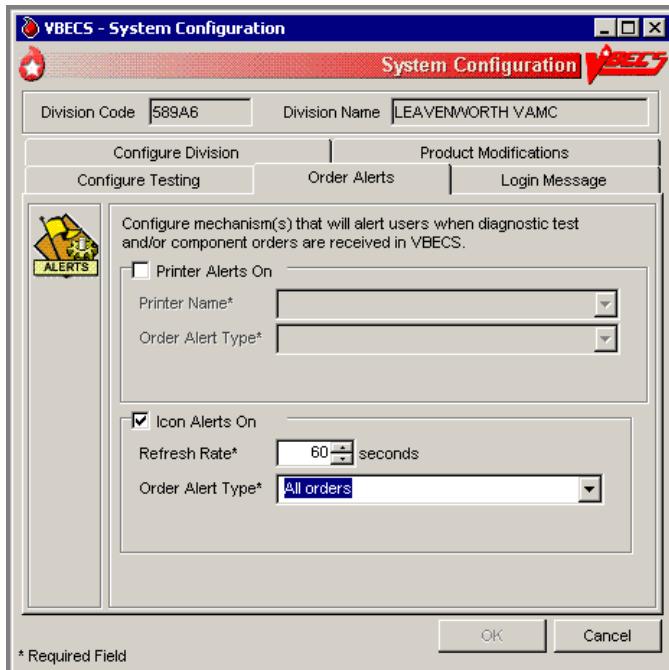

Login Message

The user creates a division-wide VBECS message that appears at login.

Assumptions

- None

Outcome

- A message was created for display when users log into a VBECS division.

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option


►►► Lead Technologist

Login Message

The user may add, update, or delete a message that VBECS displays to users logging into the division.

User Action	VBECS
1. Select Tools from the main menu. Select Configure Division .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays configuration tabs.
2. Click the Login Message tab (Figure 48). Enter a message in the message field to display at login.	<ul style="list-style-type: none">• Displays the existing login message and allows the user to edit or delete it.• Displays the updated login message. <p>NOTES —————</p> <p>VBECS does not display a blank login message.</p> <p>The user may click Clear to clear the message field and enter a new login message.</p>
3. Click OK to save.	<ul style="list-style-type: none">• Saves the login message changes.
4. Click Yes to confirm the changes and exit.	

Figure 48: Login Message Tab

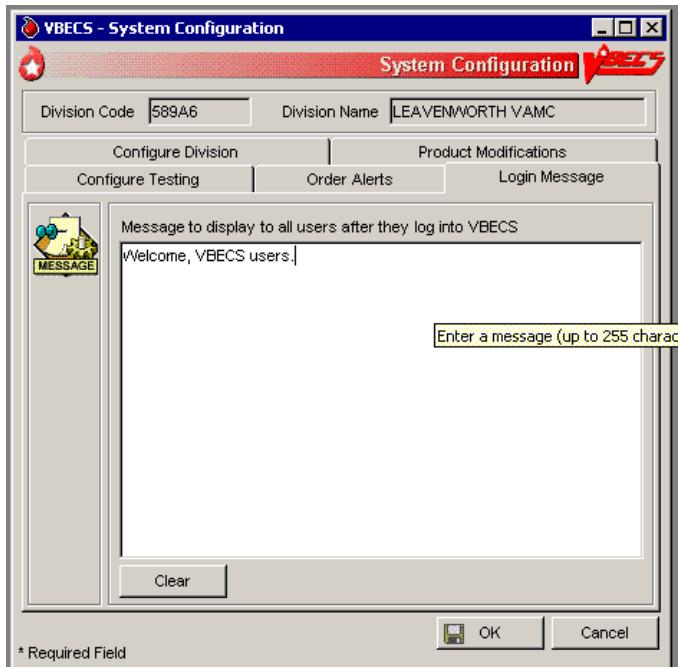

Local Facilities

The user customizes collection facility information. (A “collection facility” is not a facility for blood collection: it is for blood distribution only.)

Assumptions

- A table of known FDA-registered facilities accompanies VBECS.

Outcome

- The supplier is configured for a division.

Limitations and Restrictions

- A user may:
 - Activate and edit collection facility records distributed with VBECS.
 - Add and locally edit active collection facilities.
 - Activate and/or deactivate local collection facilities.
 - Edit a division name, which may result in a mismatch between VistA and VBECS division names.
 - The Local Facility Address maximum length is 45 characters including spaces to allow the address to display on VBECS reports. The third line of the facility address is not displayed on VBECS reports.

Additional Information

- A facility must be activated and have at least one blood product defined to be a shipper.
- A user may identify his site as a collection facility through this option.
- An ICCBBA Registration Number has one letter, indicating the country of origin (“W” for the U.S.) and four numbers. The field is five characters long.
- The user may change facility addresses that appear on reports.

User Roles with Access to This Option

All users

Local Facilities

This option allows a user to activate an entry from the national collection facility table, thereby adding it to a division’s local collection facility table, and add or update the facility’s information (except for the FDA Registration Number). A user may select an active collection facility to associate blood products (which makes the facility a shipper).

This option also allows a user to create a collection facility.

Verify that site information entered during installation is correct prior to adding local facilities.

User Action	VBECS
1. Select Tools from the main menu. Select Local Facilities .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays fields for entering facility information.• Allows a user to access an existing facility record.
2. To find a collection facility, scan or enter an FDA Registration	<ul style="list-style-type: none">• Displays the collection facilities that match the search criteria.

User Action	VBECS
Number in the FDA Reg. No. field and go to Step 4 (or click the find button to open the Facility Search window).	<ul style="list-style-type: none"> Displays properties of the selected collection facility.
3. In the Facility Search window, enter a full or partial facility name in the Partial Facility Name field (or enter a full or partial FDA Registration Number in the FDA Reg. No. field), and click Search to view a list of active facilities. Select a facility name from the list and click OK .	<p>NOTES</p> <p>The collection facility is tied to the user's current login division. A user may activate a facility. By default, VBECS marks active facilities as collection facilities.</p> <p>►►►► The Enhanced Supervisor may enter search criteria that do not exist in the national or the local collection facility tables to add a collection facility to the division's local list. By default, VBECS marks such entries as a collection facility.</p> <p>When a user clears the Active Facility check box for a collection facility in a division, warns the user, and asks whether he wishes to inactivate this facility and continue.</p> <p>Yes allows the action to continue. No leaves the collection facility marked as active for the division.</p> <p>VBECS does not allow inactivation of the product type for a shipper with in-date blood component units. VBECS warns that it is unable to deactivate the product type or shipper as there are in-date products for the shipper.</p>
4. Edit or enter information in the collection facility information fields (for example, address, telephone number, whether the facility is a testing and/or collection facility) (Figure 49) and click OK to save the updated information.	<ul style="list-style-type: none"> Displays entries for review and acceptance. Prompts to save. <p>NOTES</p> <p>OK remains disabled until the user edits or enters information or simply clicks in a field or check box.</p> <p>The Codabar eye-readable prefix (when used) is always two characters.</p>
5. Click Yes to confirm the save.	<ul style="list-style-type: none"> Saves the active collection facility record and prompts to add or edit another facility. <p>NOTES</p> <p>VBECS saves an active collection facility only when it is identified as a collection facility.</p> <p>A facility must be activated and have at least one blood product defined to be a shipper. Refer to Blood Products section for more information.</p>
6. Repeat these steps to process another facility.	
7. Click Cancel in the Maintain Facility screen to exit.	

Figure 49: Maintain Facility

VBECS - Maintain Facility

Maintain Facility

 **FACILITY**

FDA Reg. No.*	1474346	Find	Clear
Facility Name*	VAMC Hines, IL		
ICCBBA Reg. No.	W0384		
Facility Address	5th Avenue and Roosevelt Road		
Address Line 2			
Address Line 3			
City*	Hines		
State*	IL	Phone	(<input type="text"/>) <input type="text"/> - <input type="text"/>
Zip Code	60141-0000		
Codabar Unit Format			
Unit ID contains alpha characters?	<input checked="" type="checkbox"/>		
2-digit eye-readable prefix (if applicable)	<input type="text"/>		
Facility Type*			
Collection Facility?	<input checked="" type="checkbox"/>		
Testing Facility?	<input checked="" type="checkbox"/>		
Active Facility?	<input checked="" type="checkbox"/>		

*** Required Field**

OK Cancel

Blood Products

The user activates, deactivates, and edits a limited number of variables for an existing blood product site parameter. VBECS provides a Codabar and an ISBT 128 Blood Product reference table to facilitate these and other actions.

Assumptions

- The connection to VistA is active.
- The user understands Healthcare Common Procedure Coding System (HCPCS) codes. (See Additional Information.)
- A shipper must be marked as an active collection facility for a division before a user may select it as a valid shipper for a blood product for that division. This includes marking the division facility itself as an active collection facility to accommodate the creation of new blood products through modification.

Outcome

- Blood product codes are activated in the user's division.
- Each product's HCPCS code and its activation and deactivation dates are maintained in a permanent historical record.

Limitations and Restrictions

- Blood product records distributed through VBECS are national standards.
- The division must have activated the valid target blood products for modification.

Additional Information

- HCPCS codes are required for medical and surgical services and supplies for the outpatient setting under the final code set standards for the Health Insurance Portability and Accountability Act of 1996 (HIPAA). See your local Current Procedural Terminology (CPT) coordinator for more information.
- At least one valid shipper and product cost must be defined for each active product code in use in the division; multiple shippers may be associated with a blood product code.

User Roles with Access to This Option

All users

Blood Products

Before processing a blood unit in an incoming shipment, or activating or creating a blood unit in-house, the user identifies the shipper(s), cost, and return credit percentage for each shipper. The user also tags each product as active for use within the division and assigns an HCPCS code to each. These are the only fields the user may edit when maintaining an existing blood product code.

A "shipper" is the collection facility that delivers and bills for the blood product, regardless of the facility named on the unit's label. The collection facility that delivers the product to the shipper is not considered a shipper nor does it have a product cost associated with it, but it needs to be activated in the division for VBECS to interpret the unit ID. When components are modified within the division, the shipper is defined as the division. The cost may be \$0.00 or other U.S. dollar value assigned by the shipper.

User Action	VBECS
1. Select Tools from the main menu. Select Blood Products.	<ul style="list-style-type: none"> • Displays options for processing administrative functions. • Displays fields for entering blood product information.
2. To identify the product, scan or enter a Codabar or ISBT product code. Click the find button to search for a product. Select a product type from the drop-down Product Type list or enter a partial product name in the Partial Product Name field. Click Search .	<ul style="list-style-type: none"> • Displays current blood product information. <p>NOTES</p> <p>When a user enters a product code that VBECS does not recognize, VBECS warns that the blood product was not found and that it does not support entry of new blood products.</p> <p>OK returns the user to a blank screen to enter a different product. VBECS clears the screen. The user may enter a different product code.</p>
3. Select a product from the Search Results list and click OK (Figure 50). Select an HCPCS code from the drop-down HCPCS Code list or enter an HCPCS code.	<ul style="list-style-type: none"> • Allows a user to enter the HCPCS code. <p>NOTES</p> <p>A blood product may have only one HCPCS code at a time.</p>
4. Select or enter product shipper information: shipper name, FDA Registration Number, cost, and return credit percentage.	<p>NOTES</p> <p>Site parameters determine which products are active for a division. Products marked as active can be received in an incoming shipment or created in-house through modification.</p> <p>A user with no increased security can maintain an existing blood product for use only within the user's current division.</p> <p>The return credit default is "100%."</p> <p>VBECS does not allow inactivation of the product type for a shipper with in-date blood component units. VBECS warns that it is unable to deactivate the product type or shipper as there are in-date products for the shipper.</p>
5. Repeat Steps 3 and 4 for each selected product's shipper and for other products.	<ul style="list-style-type: none"> • Displays the information entered. • Requests confirmation to update the database. • Updates the database or warns the user with no override for discrepancies.
6. Click OK to exit.	

Figure 50: Maintain Blood Products

VBECS - Maintain Blood Products

Maintain Blood Products

Blood Product

Product Code*	04210	1nd	Clear
Short Product Name	RBC CPD>AS1		
Full Product Name	AS-1 RED BLOOD CELLS		
Product Type	RED BLOOD CELLS		
Volume	350	mL	
Maximum Storage	42	Days	0 Hours

HPCPS

HPCPS Code*	P9023
HPCPS Text*	FROZEN PLASMA, POOLED, SD

Product Shippers

Shipper Name*			
FDA Reg. #*			
Cost (\$)*	0.00	Up	Down
Return Credit %*	100	Up	Down
Add			

Active	FDA Reg #	Shipper	Cost	R
<input checked="" type="checkbox"/>	1474396	VAMC North Chicago, IL	\$1.00	1
<input checked="" type="checkbox"/>	1472204	LifeSource - Glenview	\$1.00	1

* Required Field

OK Cancel

Antibodies

The user edits certain fields in the supplied antibody table.

Assumptions

- VBECS is configured to reflect antibody reference parameters. Only some of these parameters may be changed locally.
- The default settings are enforced without user interaction.

Outcome

- Each division has a populated antibody table whether or not the user makes any edits.

Limitations and Restrictions

- None

Additional Information

- The antibody table defines an antigen negative requirement for all instances of the antibody for all patients in the division. When an antibody specificity does not generally require antigen negative blood but a patient has a unique appearance of that antibody that requires antigen negative blood, the user may enter the antigen negative requirement for that patient in Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Antibodies

The user chooses to edit the system-supplied blood antibody table. The Compatibility Percentage field is fully editable (except for ABO antibodies and anti-D) to allow facilities to reflect local populations. The user may edit only the Clinically Significant entries initially exported as “No.”

User Action	VBECS
1. Select Tools from the main menu. Select Antibodies .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Lists available antibodies based on the system-supplied table. <p>NOTES _____</p> <p>See Appendix B: Table 13: Antibody and Antigen Table: Irregular Antibodies.</p> <p>The table's presentation order is predefined at implementation and may not be sorted.</p>
2. Select the antibody, edit each field, and click Update (Figure 51).	<ul style="list-style-type: none">• Displays data and allows the user to edit them.• Accommodates data entry of whole numbers (0–100) and edit of the Antigen Negative Compatibility Percentage and Higher-Level Override fields. The user may not edit fields displayed on a gray background. <p>NOTES _____</p> <p>When Table 13: Antibody and Antigen Table: Irregular Antibodies “antibody sets an antigen negative requirement” field is defined as REQ (required) at installation, it is not editable.</p>

User Action	VBECS
	When the default settings is OPT (optional), the user may change the name to REQ at any time and it will remain editable.
3. Repeat for each antibody to be edited.	
4. Click OK and Yes to confirm changes and exit.	

Figure 51: Maintain Antibodies

The screenshot shows the 'VBECS - Maintain Antibodies' window. At the top, there's a toolbar with icons for back, forward, search, and other functions. Below the toolbar is a title bar with the window name and a close button. The main area is divided into two sections: 'Antibody List' and 'Selected Antibody'.

Antibody List: A grid table showing a list of antibodies with columns: Antibody Name, Clinically Significant?, AG Neg %, and Higher-Level Override?. The rows include: Anti-C, Anti-E, Anti-c, Anti-e, Anti-f, Anti-G, Anti-K, Anti-k, Anti-kp(a), Anti-kp(b), Anti-Jk(a), Anti-Jk(b), Anti-Jk3, and Anti-EvRn. The 'Anti-C' row is highlighted.

Antibody Name	Clinically Significant?	AG Neg %	Higher-Level Override?
Anti-C	Yes	53	Yes
Anti-E	Yes	75	Yes
Anti-c	Yes	12	Yes
Anti-e	Yes	2	Yes
Anti-f	Yes	N/A	Yes
Anti-G	Yes	N/A	Yes
Anti-K	Yes	91	Yes
Anti-k	Yes	1	Yes
Anti-kp(a)	Yes	98	Yes
Anti-kp(b)	Yes	2	Yes
Anti-Jk(a)	Yes	15	Yes
Anti-Jk(b)	Yes	42	Yes
Anti-Jk3	Yes	N/A	Yes
Anti-EvRn	Yes	63	Yes

Selected Antibody: A panel on the right showing details for the selected antibody (Anti-C). It includes fields for Antibody Name (Auto-filled), Clinically Significant? (checked), Antigen Negative Compatibility Percentage (53%), and a checkbox for 'Require high-level user to override antigen negative requirement?' (checked). There are 'Update' and 'Cancel' buttons at the bottom.

Canned Comments

The user defines canned comments for a division.

Assumptions

- None

Outcome

- Canned comments for a division are available for selection.

Limitations and Restrictions

- The user may not edit comments provided with VBECS but may create additional comments.
- Site personnel must enter comments in accordance with local policy and procedures.

Additional Information

- Canned comments are specific to each division. See Appendix B: Table 16: Canned Comment Category Types and Text defines the category types and comments.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Canned Comments

It is recommended that the user review the canned comments provided with VBECS for each category type and inactivate those that will not be used.


The user creates canned comments to use in lieu of entering text. Canned comments are restricted for selection based on the category type of the comment field.

The user may select a comment and has the option to enter explanatory text in the details field.

User Action	VBECS
1. Select Tools from the main menu. Select Canned Comments .	<ul style="list-style-type: none">• Displays options for processing administrative functions.
2. Select a comment category type from the drop-down list (Figure 52).	<ul style="list-style-type: none">• Lists the category types for canned comments.• Allows the user to create a comment or inactivate an existing one.
3. Enter text in the Comment field and click Add to add a comment. Click one or more check boxes to inactivate existing comments.	<ul style="list-style-type: none">• Displays entered data and allows the user to edit them before saving. <p>NOTES —————</p> <p>A user may:</p> <ul style="list-style-type: none">• Add canned comments.• Inactivate canned comments.• Reactivate inactive comments.• View inactive comments and dates of inactivation.• Control the order in which VBECS displays comments in a context.
4. Click the Show inactive Comments check box to view inactive comments.	<ul style="list-style-type: none">• Displays inactive comments and allows the user to reactivate them.

User Action	VBECS
Click one or more check boxes to reactivate inactive comments.	
5. Repeat Steps 2–4, as needed. Click OK and Yes to save and exit, or Cancel to exit without saving.	<ul style="list-style-type: none"> Updates the database.

Figure 52: Maintain Comments

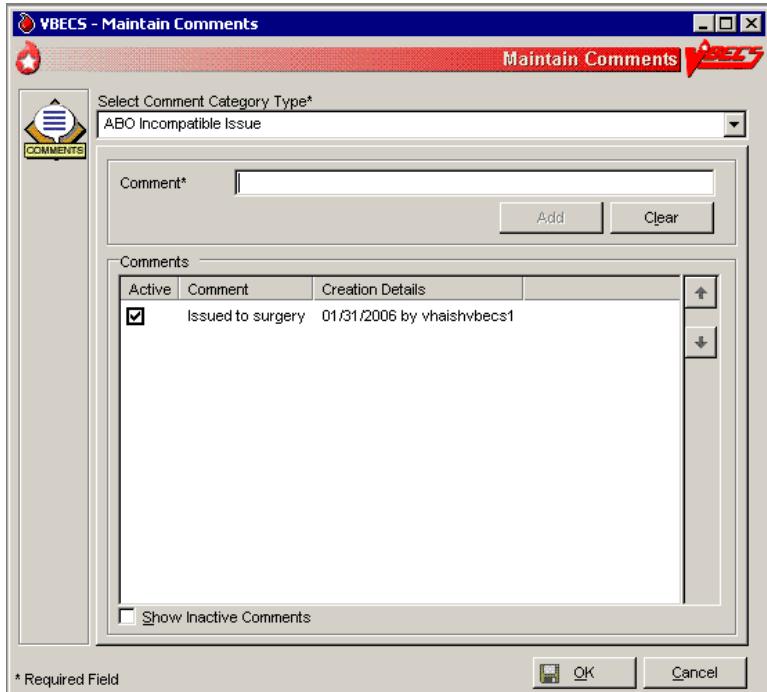

Setting Transfusion Parameters

Transfusion Complications

The user defines Laboratory tests and threshold values used to generate a report listing patients who were transfused within a specified period and who may have developed a transfusion-associated disease.

Assumptions

- VistA tests are defined.
- The connection to VistA is active.

Outcome

- Parameters are defined for future Transfusion Complication Reports.

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option

►►►► Enhanced Supervisor


Transfusion Complications

The user defines Laboratory tests for use in the transfusion complication table.

User Action	VBECS
1. Select Tools from the main menu. Select Transfusion Complications .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays fields for entering transfusion complication threshold data.
2. Click the find button to search for a Laboratory test. Select a test. Click OK .	<ul style="list-style-type: none">• Allows the user to edit or add previously defined tests for use in a Transfusion Complication Report. <p>NOTES —————</p> <p>VBECS displays the entered laboratory test names in all uppercase letters. The retrieval of lab tests occurs in VistA and results for the Transfusion Complications Report are not case sensitive. There is a limit of 25 active entries. (If the user enters all of the available results for the test, he can get results that are not filtered by VBECS.)</p>
3. Click the Numeric Threshold radio button and enter or edit test threshold information (Figure 53). Click Add to add the test to the list, or Click Update to save changes to a previously saved test.	<ul style="list-style-type: none">• Allows the user to enter one or more tests.• Uses VistALink to retrieve test data. <p>NOTES —————</p> <p>When defining a test reported as a text value, VBECS accommodates the input of specific text (1–255 characters) that will be used to search for matching values.</p>

User Action	VBECS
Click OK to save and Yes to confirm the save.	
4. Click OK to save and Yes to confirm the save and exit.	<ul style="list-style-type: none"> Updates the database.
5. Repeat Steps 2–4 for each transfusion complication test to be defined or click Cancel to exit.	

Figure 53: Transfusion Complication Thresholds

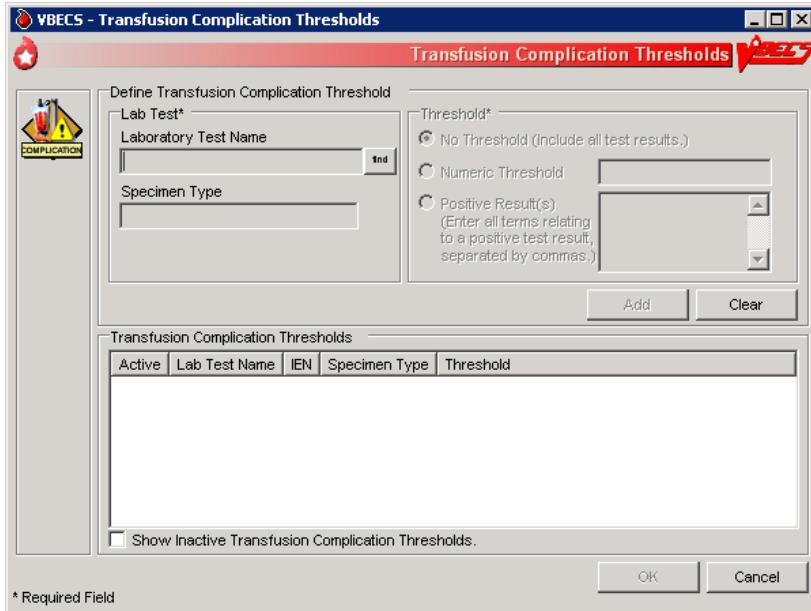

Transfusion Effectiveness

The user defines Laboratory tests to monitor the effectiveness of transfusions.

Assumptions

- VistA lab tests are defined.
- The connection to VistA is active.

Outcome

- Parameters are defined for future reports.

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option


►►►► Enhanced Supervisor

Transfusion Effectiveness

The user enters the tests to be included in the Transfusion Effectiveness Report.

User Action	VBECS
1. Select Tools from the main menu. Select Transfusion Effectiveness .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Allows the user to edit a previously defined Transfusion Effectiveness Tests or add a test.
2. Click the find button to search for a Transfusion Effectiveness Test. Enter part or all of a laboratory test name. Click Search . Click one to four check boxes to select tests. Click OK . Click Add to add the tests (Figure 54), or Click the Show Inactive Transfusion Effectiveness Report Tests check box, if desired. Click Update to edit the tests.	<ul style="list-style-type: none">• Allows the user to enter one or more Transfusion Effectiveness Tests.• Uses VistALink to retrieve Laboratory test data for each of these indicators:<ul style="list-style-type: none">o Laboratory test nameo VistA Laboratory test Internal Entry Number (IEN)o Specimen type
3. Repeat Step 2 for each Transfusion Effectiveness Test to be defined, or click Cancel to exit.	
4. Click OK to save and Yes to confirm the save.	<ul style="list-style-type: none">• Updates the database, including the date and time and the technologist ID. Saves changes to previously defined tests for inclusion in an Audit Trail Report.

Figure 54: Transfusion Effectiveness Report Tests

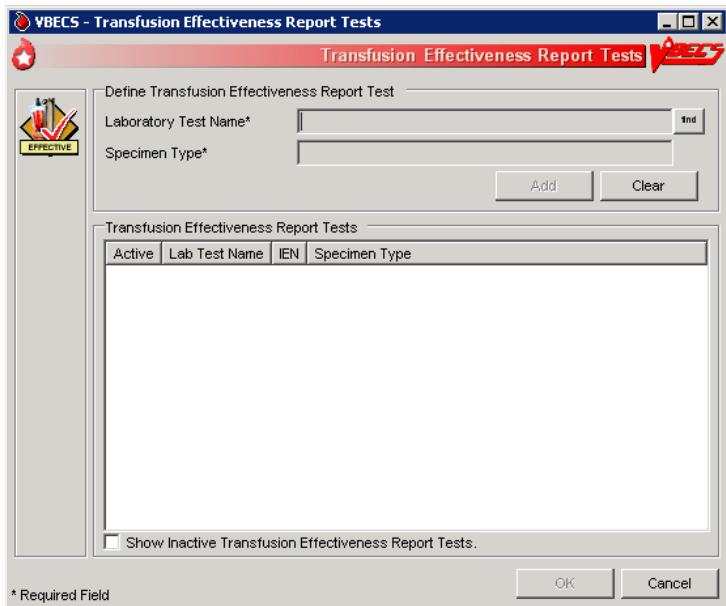

MSBOS

The user defines parameters for comparing blood component orders placed by a clinician with the Maximum Surgical Blood Ordering Schedule (MSBOS).

Assumptions

- A table of generic surgery names is available for MSBOS activation.

Outcome

- CPRS users will see only active surgery names when placing a blood bank order.
- MSBOS will be available to VBECS users to evaluate the appropriateness of pre-op orders.

Limitations and Restrictions

- None

Additional Information

- The maximum number of blood products serves as a guideline. Orders in excess of the MSBOS are not prohibited.
- The VBECS user may view active and inactive surgery names.

User Roles with Access to This Option


►►►► Enhanced Supervisor

MSBOS

The user selects a surgery name from the generic surgery list or creates a new surgery name, then defines the associated MSBOS recommendation. MSBOS is used when VBECS acknowledges orders from VistA that are written for preoperative patients.

User Action	VBECS
1. Select Tools from the main menu. Select MSBOS .	<ul style="list-style-type: none"> • Displays options for processing administrative functions. • Displays the schedule and component recommendation options.
2. To activate an existing surgery name, click the Show Inactive MSBOS entries check box (Figure 55). Select a surgery name. Click the Active MSBOS Entry? check box, or To add and activate a surgery name, click New Surgery and enter a surgery name in the Surgery Name field. Click Add , or To edit and/or deactivate an existing surgery name, click the appropriate options in the MSBOS Component Recommendation area of the screen edit. Click Update .	<ul style="list-style-type: none"> • Lists active surgery names or allows the user to create an entry. • Activate, add, or edit a MSBOS entry and its associated data. • Displays surgery names available for activation. • Displays entries for review and acceptance. <p>NOTES</p> <p>MSBOS is specific to the division.</p> <p>See Appendix B: Table 24: VBECS Maximum Surgical Blood Order Schedule (MSBOS).</p> <p>Each surgery name has a Type & Screen (TAS) (no blood order recommended) or one or more component classes defined as recommended orders.</p> <p>The user may select one or more component classes for MSBOS recommendation. Each defined component class must have an associated maximum number recommended. The user may deactivate previously activated component classes for an MSBOS recommendation. The OTHER component class is not available for selection.</p>
3. Repeat Step 2 until all new entries are created or existing entries are edited. Click OK to save the updated data.	<p>NOTES</p> <p>VBECS does not save changes to the database when the user clicks Cancel without clicking OK.</p>
4. Click Cancel to exit.	

Figure 55: Maintain MSBOS

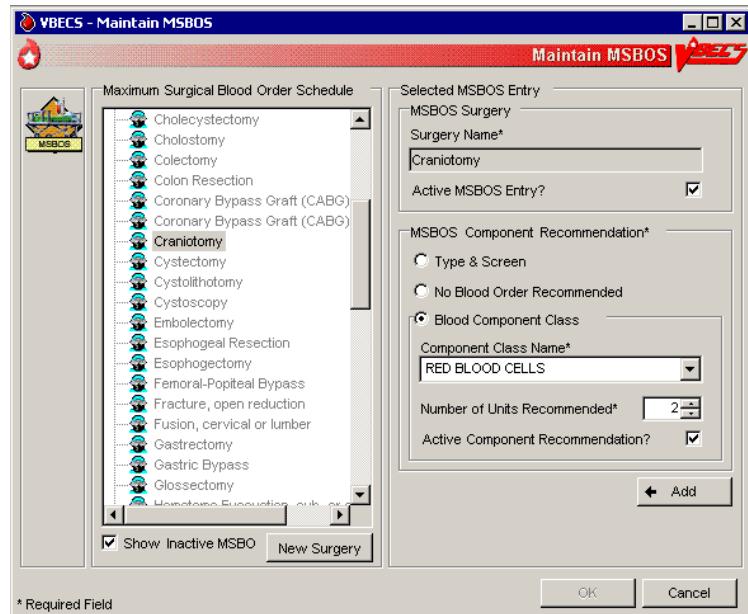

Workload Codes

The user assigns workload Laboratory Management Index Program/National Laboratory Test (LMIP/NLT) tests to processes performed through VBECS. VBECS test names are associated with VistA workload codes, which VBECS uses when compiling and storing data to provide workload-related reports to VistA.

Assumptions

- The VistA database was updated to include all required workload codes for the VBECS process.
- The VistA database was updated to associate CPT codes and costs with the workload codes to be used by VBECS.
- The connection to VistA is active.

Outcome

- Data are available to the VistA workload package.

Limitations and Restrictions

- Transfusion-only sites must not map workload codes to processes, such as ABO/Rh testing, not performed on site.
- Retrieved LMIP/NLT codes from VistA may have multiple CPT codes. Inactivating any of the associated CPT codes in VBECS will deactivate the LMIP/NLT code. Editing of the LMIP/NLT code associations must be done in VistA.

Additional Information

- VistALink is required to query VistA files 64 (NLT/LMIP) and 81 (CPT) for appropriate entries.
- The division provides CPT codes and costs associated with the workload codes when they need updating.
- Each site is responsible for notifying users when incorrect CPT codes are used.
- Workload events are placed in a queue for transmission via VistALink.

User Roles with Access to This Option

>>>> Administrator/Supervisor


Workload Codes

As part of record-keeping activity, VBECS allows the user to account for workload units associated with each activity and to assign costs incurred by the blood bank due to these activities. Data gathered because of this workload collection are available to VistA, as needed, for various reports.

User Action	VBECS
1. Select Tools from the main menu. Select Workload Codes .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays fields for entering workload-related data.• Lists and retrieves processes that may involve workload.
2. Select a process for which to define workload codes.	<p>NOTES —————</p> <p>The user may assign multiple LMIP and CPT workload codes to a process.</p> <p>VBECS checks for active and inactivated workload associations</p>

User Action	VBECS
	<p>and does not allow the user to enter a duplicate association. VBECS allows the user to reactivate previously inactivate entries and add workload associations.</p>
3. Select an LMIP/NLT code. Click Add to add the code, or click Update to select a current association to edit.	<ul style="list-style-type: none"> Lists previously assigned LMIP/NLT codes for the selected process. Allows the user to select an LMIP/NLT code. <p>NOTES</p> <p>VBECS uses VistALink to query the VistA WKLD CODE File (#64) for valid LMIP/NLT codes.</p> <p>VistA uses the LMIP/NLT code to create workload reports as well as reports to patient care encounter (PCE) and the decision support system (DSS) database.</p>
4. Enter the start or stop date (whichever is appropriate) of the workload code association with the VBECS process (Figure 56).	<ul style="list-style-type: none"> Allows the user to enter a date. <p>NOTES</p> <p>The user may enter an inactivation date only when an activation date was previously entered (default start date and time: current date and time).</p> <p>The inactivation date must be later than or the same as the activation date. When the user attempts to enter an inactivation date that is earlier than the activation date, VBECS warns the user. There is no override.</p>
5. Repeat Steps 4–5 for each LMIP/NLT code to be associated with the selected process.	
6. Add the selected parameters to the selected process.	<ul style="list-style-type: none"> Allows the user to add the selected parameters.
7. Repeat Steps 1–8 for to define workload parameters for another process, or exit.	

Figure 56: Assign Workload Codes

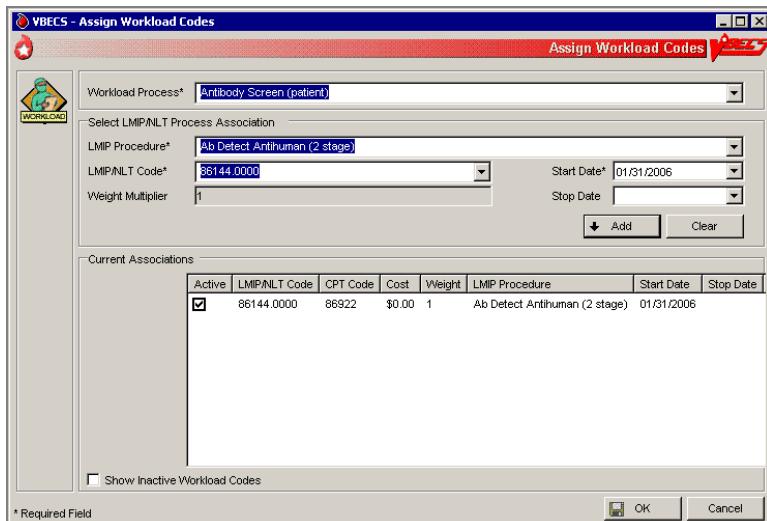

Remapping Workload Codes

When workload codes or a workload code multiplier is changed VBECS must be updated. The active workload process to LMIP/NLT Code association must be inactivated and remapped for VBECS to track workload correctly.

User Action	VBECS
1. Select Tools from the main menu. Select Workload Codes .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays fields for entering workload-related data.• Lists and retrieves processes that may involve workload.
2. Select a process to edit workload codes.	<ul style="list-style-type: none">• Allows the user to edit existing data or add data.
3. Click one or more check boxes to deactivate LMIP/NLT Code. Click OK and Yes to save and exit.	<ul style="list-style-type: none">• Removes the code from the Current Associations window. <p>NOTES —</p> <p>Inactive codes can be viewed by clicking the Show Inactive Workload Codes checkbox.</p>
4. Click OK and Yes to save and exit.	<ul style="list-style-type: none">• Saves the changes to the workload code mapping.
5. Select Tools from the main menu. Select Workload Codes	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays fields for entering workload-related data.• Lists and retrieves processes that may involve workload.
6. Select a process to edit workload codes.	<ul style="list-style-type: none">• Allows the user to edit existing data or add data.
7. Select an LMIP/NLT code. Click Add to add the code.	<ul style="list-style-type: none">• Allows the user to select an LMIP/NLT code. <p>NOTES —</p> <p>LMIP/NLT weight multipliers values are set in VistA and that value is imported into VBECS.</p> <p>Changes made to VistA files 64 (NLT/LMIP) and 81 (CPT) are updated in VBECS on a rolling 24 hour clock. The first user to log in 24 hours after the last update will trigger the next update.</p>
8. Repeat Steps 1–7 for to edit workload parameters for another process.	

Reagents and Supplies

Enter Daily QC Results

The user enters observed serological reactions for all phases of testing and interpretations of the various routine reagents.

Assumptions

- The division is “full service.”
- Site parameters regarding reagent quality control are set up.
- Minimum inventory levels for reagent types are defined.
- Inventory control functionality is in use and reagents are in the inventory of the user’s division.
- Daily QC testing does not include reagents used in special antigen typing.

Outcome

- Routine reagents are available for use, shutting down the alert for the day associated with the rack.
- The user may enter and save reaction results below the minimum expected reaction; an Exception Report may be generated for supervisory review; the rack may or may not be used by the site for routine testing throughout that day.

Limitations and Restrictions

- Additional user-defined reagent types not available in Configure Daily QC or Enter Daily QC Results are not part of a reagent rack's daily QC testing.

Additional Information

- VBECS allows the user to save a partially filled worksheet for a specific QC rack ID and allows the same or a different user to redisplay the form for the entry of remaining test results at a later time. The user, first or second, may enter new results but may not edit prior entries. VBECS records the ID of each user who enters test observations on a worksheet.
- VBECS issues warnings when they are defined in Maintain Minimum Levels.
- VistALink does not need to be active when a user enters daily quality control results.
- A testing template is generated for daily quality control testing of routine reagents for each rack in a division. (Multiple templates may be used on a given day.)
- System rules determine the validity of the test interpretations and determine if all phases of testing were recorded for both scenarios.
- A user must enter and may view one day's QC data at a time. The supervisor or administrator must review these entries and the reagent information daily.
- VBECS compares the current rack QC results with the prior day's QC results.
- Multiple technologists may use the same rack.
- Rack number and reagents used to test a specimen are linked to the lot number, manufacturer, and expiration date of the reagent type as well as the rack identifier.
- A user may use more than one testing method in accordance with local policies and procedures.
- If a rack has never been used with reagents prior to using the QC Documentation Offline option and the rack is not so checked when the user attempts to use this untested rack, the message appears that this rack has never been tested and cannot be used. The user may not proceed. To use this rack on this date, the QC Offline check needs to be indicated or testing performed.

User Roles with Access to This Option

All users

Enter Daily QC Results

The American Association of Blood Banks' *Standards for Blood Banks and Transfusion Services* requires that commercial reagents in the blood bank—antisera and routinely-used red cell reagents—must be tested for performance “on the day of use.” As routine reagents are always used, they are tested daily. Records of these reagents, their lot numbers, and their testing results must be stored for reference. VBECS verifies and indicates whether the reagent is satisfactory for use.

This option facilitates the direct entry of reaction results and interpretations, constituting daily reagent quality control testing for a rack of reagents.

VBECS allows the user to change individual reagents without repeating the QC for the full rack during the day.

User Action	VBECS
<p>1. Select Reagents from the main menu.</p> <p>Select Enter Daily QC Results.</p> <p>2. Verify the default date tested and tester, and edit them, if needed.</p>	<ul style="list-style-type: none"> Displays options to perform daily reagent rack QC. <p>NOTES</p> <p>Daily QC testing does not include reagents used in special antigen typing.</p>
<p>3. Select the Test Rack(s), Partially Tested Rack(s), or Retest Rack(s) tab (Figure 57).</p> <p>Click one or more check boxes to select rack numbers.</p> <p>Click OK.</p>	<ul style="list-style-type: none"> Displays the lot number verification template for each rack selected. <p>NOTES</p> <p>A reagent rack is displayed on one of three tabs based on its testing status and the daily alert time setting:</p> <ul style="list-style-type: none"> Test Rack(s): the current alert time has passed and is/are considered "untested" Partially tested: racks that have incomplete test results for the current 24 hour time clock and are considered "untested" Retest Racks: racks that are fully tested within the current alert time and are available for new lot number selection and its testing. <p>The user may test a rack no earlier than the alert time for the next 24 hours.</p> <p>On occasion, a user may save a partially completed QC worksheet and finish it later. The same user or a different user may complete the testing.</p> <p>The user may retest a rack when a reagent lot number changes after initial testing. VBECS does not allow the user to retest saved valid results.</p> <p>VBECS requires the user to indicate the rack identifier to provide a testing template for each daily reagent rack QC.</p> <p>Site settings determine reagent types and lot numbers present. VBECS displays only reagent types previously entered for the division using the reagent inventory function. VBECS allows entry of lot numbers only. When VBECS displays lot numbers for selection, it checks the inventory system for the reagent name and lot numbers in use in a rack. VBECS does not offer a reagent with an unsatisfactory value in the "Inspection" field or with a quantity of zero as a selectable reagent. The user cannot force the entry of a lot number that VBECS disallowed on the selection list. The user must choose a different reagent.</p> <p> VBECS compares the expiration date of each reagent on the selected rack to the testing date, emits an audible alert, and warns when the reagent is expired.</p>

User Action	VBECS
	<p>Yes requires a comment and VBECS captures details for inclusion in an Exception Report (exception type: expired reagent QC'd).</p> <p>No removes the unsatisfactory reagent from the daily QC rack. The user must select a new one. VBECS does not display unsatisfactory lot numbers.</p> <p>The user may enter the lot numbers of the reagents used for the daily QC rack and the lot numbers of the enhancement media reagents defined in Configure Daily QC. VBECS does not display these reagents on the reaction results grid for a rack.</p>
<p>4. Select a worksheet to complete the testing (Figure 58), or</p> <p>Click the Rack not in use today check box to indicate that testing of a rack will not be performed.</p> <p>Click the QC Testing Documentation Offline check box to acknowledge QC is not performed in VBECS but has been properly documented. Click Yes to save changes. Close the window.</p>	<ul style="list-style-type: none"> Displays the selected worksheet. Displays the selected partially completed worksheet with the option to cancel the worksheet or continue to enter more testing information. <p>NOTES</p> <p>When the QC Testing Documentation Offline check box is checked for a reagent rack, VBECS allows the rack to be used for patient testing until the next daily alert time.</p> <p>The primary lot number can be tested and changed and the results for both will display on the Testing Worklist Report.</p> <p> A user may cancel or invalidate a partially completed QC rack worksheet when it is available. VBECS emits an audible alert and warns that invalidating this worksheet will necessitate starting over and asks whether the user wishes to continue.</p> <p>No returns the user to the partially partially-completed worksheet. Yes invalidates the worksheet. VBECS does not delete reaction results from the database.</p> <p>An Exception Report comment documents the invalidation of the rack's worksheet. VBECS does not record workload calculation for invalidated testing. The invalidated reaction results and associated testing information remain in the database for future retrieval (exception type: deletion of partially completed QC).</p> <p> VBECS compares the results entered to those of the previous day for the same rack number and alerts the user when there is a change in reactivity of two or more or when the result is changed from positive to negative for each reagent type by lot number. VBECS emits an audible alert and warns that there is a decrease in reagent reactivity of two or more and asks whether the user wishes to continue to use this vial.</p> <p>Yes allows the user to continue and requires a comment. VBECS captures details for inclusion in an Exception Report and allows this reagent to be used with the change (within acceptable limits) (exception type: QC decrease reagent reactivity ≥ 2).</p> <p>No allows the user to return to the lot number entry for this reagent type and approve this lot number (different vial) or enter a new lot number. VBECS returns to the testing worksheet where the user left</p>

User Action	VBECS
5. Verify or select from the drop-down menu the reagent type lot number and test method for each rack chosen for processing. Click OK .	<p>off. (An individual vial of the same lot number may replace the one losing its reactivity.)</p> <p>The screening cell and reverse typing reagents are delivered in sets and share the same lot number; therefore, VBECS displays only one reagent type on the lot number template.</p> <ul style="list-style-type: none"> • Displays the lot numbers and expiration dates for verification, if any. • Displays each rack in its own tab. <p>NOTES</p> <p>VBECS allows for one lot number entry per reagent type per rack template. When an entry was made, VBECS fills in the reagent type lot number, the expiration date of that lot number, and the manufacturer. The user may edit the lot number. The user must verify the lot numbers of racks in the selected batch before VBECS allows the user to move to the reaction results grids.</p> <p>Each reagent rack has a separate reaction results grid to collect the user's direct observations.</p> <p>VBECS displays reagent types (except for the gel and the instrument lot numbers) listed as REQUIRED (by VBECS or by the user) in the daily use field on the rack lot number template and the reaction result entry worksheet for each rack in Enter Daily QC Results.</p> <p>When a rack is tested, VBECS displays a retest option at the daily reagent QC prompt with the untested rack numbers. When a user requests to retest a rack on a given date, VBECS displays the option to reenter a lot number for reagent types for all or part of the daily QC with a comment. All other reagent type results previously entered are locked and may not be edited. VBECS saves the new information and the original data, then generates a report entry.</p> <p>When a user changes a reagent lot number, VBECS clears the testing template for that reagent type and allows the user to enter new results.</p>
6. Enter the observed serologic reactions (Figure 59). Click OK and Yes to save and exit. Repeat for each rack selected in Step 3.	<ul style="list-style-type: none"> • Displays the testing template reaction results grid for each rack selected, when no lot number changes were made. • Determines whether the entered results are acceptable and indicates acceptability of the lot number for the reagent type. • Verifies (by type and lot number) that the reagent is satisfactory or unsatisfactory for use. • Saves the partially or fully completed data. <p>NOTES</p> <p>The user may enter results on multiple reagent racks. VBECS displays each rack in its own tab.</p> <p>VBECS verifies satisfactory and unsatisfactory QC testing based on site parameters and user data entry. If the user-entered test results indicate a satisfactory interpretation, VBECS displays "Satisfactory." VBECS saves results that are satisfactory only when the testing for a row is complete and the worksheet is saved. VBECS clears results associated with an unsatisfactory interpretation from the</p>

User Action	VBECS
	<p>screen/worksheet and the testing for the affected QC rack is considered incomplete until the rack is retested with satisfactory results for the affected row. (Warnings generated by other options when QC testing is incomplete apply).</p> <p>The daily reagent QC for a date requires a valid entry for each reagent displayed on the worksheet. VBECS considers a rack ready and available for patient testing on that date when all reagents on a rack have valid test results.</p>

Figure 57: Enter Daily QC Results

VBECS - Enter Daily QC Results

Date Tested* 02/02/2006 11:58

Tested By* vhaishvbechs1

Select Rack(s)*

Test Rack(s) | Partially Tested Rack(s) | Retest Rack(s)

Rack
<input checked="" type="checkbox"/> A

* Required Field

OK Cancel

Figure 58: Record Lot Numbers for Daily QC Testing

The screenshot shows a Windows application window titled "VBECS - Enter Daily QC Results". The main area contains a table titled "Record Lot Numbers for Daily QC Testing". The table has columns for Reagent Type, Lot Number*, Manufacturer, Expiration Date, and Test Method*. There are 10 rows in the table, each with a green checkmark icon in the first column. The rows list various reagents: Anti-A, Anti-B, Anti-D, Screening Cells, Reverse ABO Typing Cells, QC kit, LISS, AHO, PS, and CC. All entries have the same manufacturer (Immucor), expiration date (07/25/2012), and test method (Tube). A key at the bottom left defines the icons: a green checkmark for Valid Reagent, a yellow warning sign for Reagent will expire within 24 hours, and a red exclamation mark for Expired Reagent. A note at the bottom left indicates "* Required Field". Buttons for "OK" and "Cancel" are at the bottom right.

Record Lot Numbers for Daily QC Testing		QC Testing Documentation Offline		<input type="checkbox"/>	<input type="checkbox"/> Rack not in use today	
Reagent Type	Lot Number*	Manufacturer	Expiration Date	Test Method*		
Anti-A	RA000001	Immucor	07/25/2012 11:15	Tube		
Anti-B	RA000003	Immucor	07/25/2012 11:16	Tube		
Anti-D	RA000007	Immucor	07/25/2012 11:17	Tube		
Screening Cells	RA000015	Immucor	07/25/2012 11:20	Tube		
Reverse ABO Typing Cells	RA000007	Immucor	07/25/2012 11:14	Tube		
QC kit	RA000099	Immucor	07/25/2012 11:22	Tube		
LISS	RA000028	Immucor	07/25/2012 11:21	Tube		
AHO, PS	RA000101	Immucor	07/25/2012 11:23	Tube		
CC	RA000019	Immucor	07/25/2012 11:20	Tube		

Key
V Valid Reagent
W Reagent will expire within 24 hours.
E Expired Reagent

* Required Field

Figure 59: Enter Daily QC Results with Serologic Reactions

The screenshot shows a Windows application window titled "Daily QC - Enter Results". The main area contains a table with columns: Reagent Type, Test With, Min Result, Phase, IS*, IAT*, and CC*. Below the table is a key and a note about valid entries. A note at the bottom left indicates "* Required Field". Buttons for "Invalidate Test Results", "OK", and "Cancel" are at the bottom right.

Reagent Type	Test With	Min Result	Phase	IS*	IAT*	CC*
S Anti-A	A1B Rh, D Neg cell	2	IS	2		
S Anti-B	A1B Rh, D Neg cell	2	IS	2		
S Anti-D	Group O Rh, D Pos cell	2	IS	2		
S Anti-D	RCQ Group O Rh, D Neg cell	0	IS/MAT	0	0	2
S Screening Cells 1	Antibody Reagent/Enhancement Media	2	IAT		2	
S Screening Cells 2	Antibody Reagent/Enhancement Media	2	IAT			2
S A1 Cell 1	Antibody Reagent	2	IS	2		
S B Cell 2	Antibody Reagent	2	IS	2		

Key
S Satisfactory
W Warning
U Unsatisfactory

Valid Entries
0 = No Agglutination M = Microscopic W = Weak 1 = 1+ 2 = 2+ 3 = 3+ 4 = 4+

* Required Field

Reagents

The options related to reagent inventory control record the receipt of specific reagents by reagent type for each division.

Log In Reagents, Update Inventory, View/Print Inventory, Maintain Minimum Levels

The reagent log is maintained by individual division even in a multidivisional database. Update Inventory allows the user to update reagent inventory by changing the number of vials available. The user may view inventory online using various search criteria and print a copy of the search results.

Assumptions

- The division is “full service.”
- Reagents are in the inventory of the user’s division.

Outcome

- A user may select previously logged-in reagents through Enter Daily QC Results, Patient Testing: Record a Patient Antigen Typing, and Unit Antigen Typing.
- Reagents entered are accessible only within the division in which they were entered.

Limitations and Restrictions

- VBECS does not display a reagent report by invoice number.
- The unsatisfactory portion of a shipment must be recorded individually.
- There is no retrospective data entry or user selection for these options.
- Each site must set a standard for recording quantity: it may use the number of cases or the number of units within the cases.
- When a user enters a lot number in the Lot # field and clicks the magnifying glass, VBECS displays only the first instance of a duplicate lot number.
- Free text details can be added only when “Other” is selected as the canned comment.

Additional Information

- This option is disabled for transfusion-only facilities.
- VBECS displays alerts, as appropriate, during execution of Enter Daily QC Results to notify the user when reagent inventory levels fall below minimum inventory levels, as defined in Maintain Minimum Levels.
- Local policies and procedures must define the use of patient antisera as reagents.

User Roles with Access to This Option

All users

Log in Reagents

A user records the receipt of specific reagents by type. The reagent record includes the date the shipment was received, reagent name, lot number, expiration date, manufacturer, and number of vials received in a shipment.

User Action	VBECS
1. Select Reagents from the main menu.	<ul style="list-style-type: none">• Displays options for Quality Control functions.• Displays options for processing reagents.
2. Select Log In Reagents .	<ul style="list-style-type: none">• Displays fields for entering invoice and reagent information.

User Action	VBECS
3. Invoice Information: Select or edit the date and time in the Date Received field. Enter the invoice number in the Invoice Number field. Enter the name of the manufacturer in the Manufacturer field, or select a name from the drop-down list.	NOTES More than one technologist at a time may process a large shipment of reagents; only one technologist at a time may process individual reagents in the shipment. The default date received is the current date and time. The user may edit it to a past date and time only. When the reagent type and lot number were already entered for this invoice, VBECS warns that the user is trying to enter a reagent that was already logged into the system for this invoice. Yes takes the user to the “Update Reagent Inventory” screen where VBECS displays the reagent type and lot number for update. No clears all fields on the screen, allowing the user to log in another reagent.
4. Reagent Information: Enter the reagent name and reagent type in the Manufacturer Reagent Name and Reagent Type fields, or select from the drop-down lists. Enter the lot number in the Lot Number field. Select or edit the date and time in the Lot Expiration Date field. Enter (or select) a quantity (1–999) in the Vials Received per Lot Number field.	NOTES The default lot expiration date is one minute before midnight (23:59) of the current date.
5. Visual Inspection Information: Click the Satisfactory radio button, or click the Unsatisfactory radio button and select a comment from the drop-down list in the Comment field. Click Add to add the reagents to the inventory list.	NOTES Satisfactory reagents are available for later selection, unsatisfactory reagents are not. When a user selects Other from the drop-down list in the Comment field, entry of additional text in the Details field is required.
6. Repeat Steps 3–5 to log in more supplies, as desired.	NOTES Click Clear to clear all but the Invoice Number, Date Received, Manufacturer, and Lot Expiration Date fields.
7. Click OK to save. Click Yes to confirm the save and exit.	<ul style="list-style-type: none"> Requests confirmation to update the database.

Update Reagent Inventory

The user changes the number of available items in the reagent inventory.

User Action	VBECS
1. Select Reagents from the main menu.	<ul style="list-style-type: none"> • Displays options for Quality Control functions. • Displays options for processing reagents.
2. Select Update Inventory .	<ul style="list-style-type: none"> • Lists existing satisfactory reagents by default.
3. Select a reagent to edit. To view reagents with a specific lot number, enter the lot number in the Lot # field and click the magnifying glass .	<p>NOTES</p> <p>Click a reagent type to reveal manufacturer names. Click a manufacturer name to reveal lot numbers.</p>
4. Click a lot number to enable the Number of Vials Available field. Enter (or select) a quantity that is less than or equal to the current entry in the Quantity Available field. Select or edit the date and time in the Date Changed field (Figure 60).	<p>NOTES</p> <p>The default change date is the current date and time. The user may edit it to a past date and time only.</p> <p>When a user selects Other from the drop-down list in the Reason for Change field, he may enter additional text in the Comments field.</p>
5. Repeat Steps 3 and 4 to edit more reagents, as desired.	
6. Click OK to save. Click Yes to confirm the save and view the updated list.	<ul style="list-style-type: none"> • Prompts the user to confirm the save. • Saves the record and redisplays current inventory, as in Step 3. • Displays an option to exit.
7. Click Close to exit.	

Figure 60: Update Reagent Inventory

VBECS - Update Reagent Inventory

Update Reagent Inventory

Inventory

Selected Reagent

REAGENTS

Lot #

Anti-A

Anti-B

Anti-A,B

Anti-D

Anti-D control

LISS

Reverse ABO Typing Cells

Screening Cells

CC

QC kit

Anti-Complement

AHG, PS

Anti-IgG

Complement CC

Anti-K

Other, Set

Reagent Type: Anti-A

Minimum Vial Stock Level: 0

Manufacturer Reagent Name: A

Manufacturer: AAA

Lot Number: 1

Expiration Date: 02/28/2006 23:59

Inspection: Satisfactory Unsatisfactory

Number of Vials Available*: 4

Date/Time Changed*: 02/01/2006 10:17

Reason for Change*

Comments

Update

Close


* Required Field

View/Print Reagent Inventory

The user views inventory online based on various search criteria and prints the search results.

User Action	VBECS
1. Select Reagents from the main menu.	<ul style="list-style-type: none"> • Displays options for Quality Control functions. • Displays options for processing reagents.
2. Select View/Print Inventory .	<ul style="list-style-type: none"> • Lists reagent types.
3. Click one or more check boxes in the Reagent Type list to select reagent types.	<ul style="list-style-type: none"> • Displays or allows the user to print a report of all current in-date, unsatisfactory, expired reagent inventory and minimum stock levels as a read-only display based on the selection criteria entered: <ul style="list-style-type: none"> o Reagent type o In-date (active) reagents o Unsatisfactory reagents o Expired (inactive) reagents o Expiration date range (the user enters a date to begin and end the search) o Date received (the user enters a date to begin and end the search)
NOTES —————	
VBECS display only reagents selected based on search criteria.	
4. Click OK to accept the displayed selection criteria (Figure 61).	<ul style="list-style-type: none"> • Compiles the report and displays an option to print or exit the report.
5. Click Print to print the list of reagents.	<ul style="list-style-type: none"> • Sends the report to print on the printer selected by the user.
6. Print another report or click Close to exit.	

Figure 61: Reagent Search

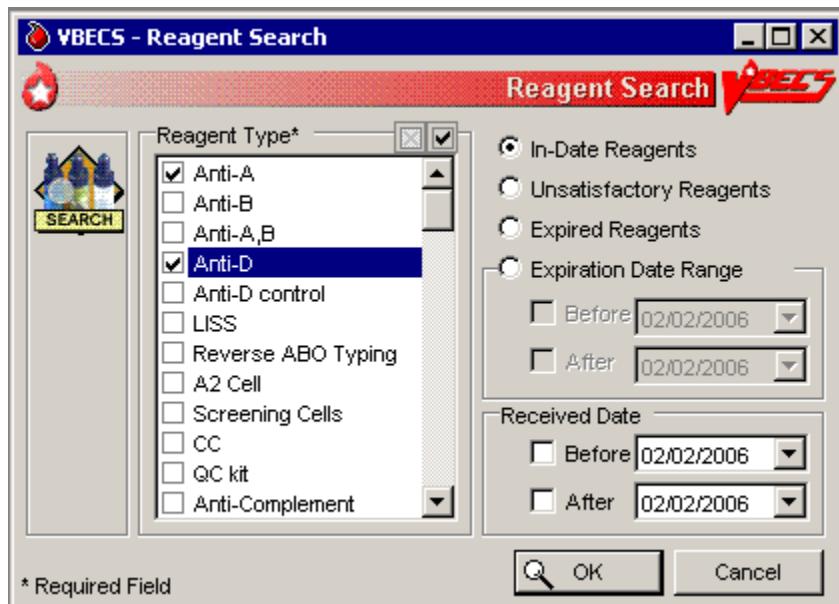


Maintain Minimum Levels

The user sets minimum stock levels for reagent types used in the division.

User Action	VBECS
1. Select Reagents from the main menu.	<ul style="list-style-type: none">• Displays options for processing reagents.
2. Select Maintain Minimum Levels .	<ul style="list-style-type: none">• Lists available reagent types and their associated minimum stock levels.
3. Select a Reagent Type to edit the minimum level threshold.	<ul style="list-style-type: none">• Displays the reagent type selected and allows the user to edit the Minimum Stock Level field or select another reagent. <p>NOTES</p> <p>A minimum inventory level (a whole number from 0 to 999) is set for each reagent type. This is related to the reagent type, independent of manufacturer name and lot numbers.</p>
4. Enter (or select) a minimum number of vials in the Minimum Stock Level field.	<ul style="list-style-type: none">• Displays the changes and allows the user to select another reagent to edit. <p>NOTES</p> <p>When the number of vials in stock is at or is less than the stock minimum, VBECS notifies the user and instructs him to check the inventory.</p> <p>Inventory levels are tied to the reagent type, not to the manufacturer reagent name or the lot number. A particular type may have multiple manufacturer reagent names and multiple lot numbers associated with that manufacturer reagent name, all of which are summed to create the inventory of that reagent type. VBECS tracks the total number of vials in the system for a particular reagent type.</p> <p>Each reagent type has one minimum inventory level established regardless of manufacturer name and lot numbers. Valid entries are whole numbers from 0 to 999.</p>
5. Review the displayed information and click Update (Figure 62).	<ul style="list-style-type: none">• Displays the changes and allows the user to select another reagent to edit.• Prompts the user to confirm the save.
6. Repeat Steps 3–5 until all minimum reagent levels are set.	
7. Click OK to save, and Yes to confirm the save.	<ul style="list-style-type: none">• Saves the record and exits.

Figure 62: Maintain Minimum Reagent Levels

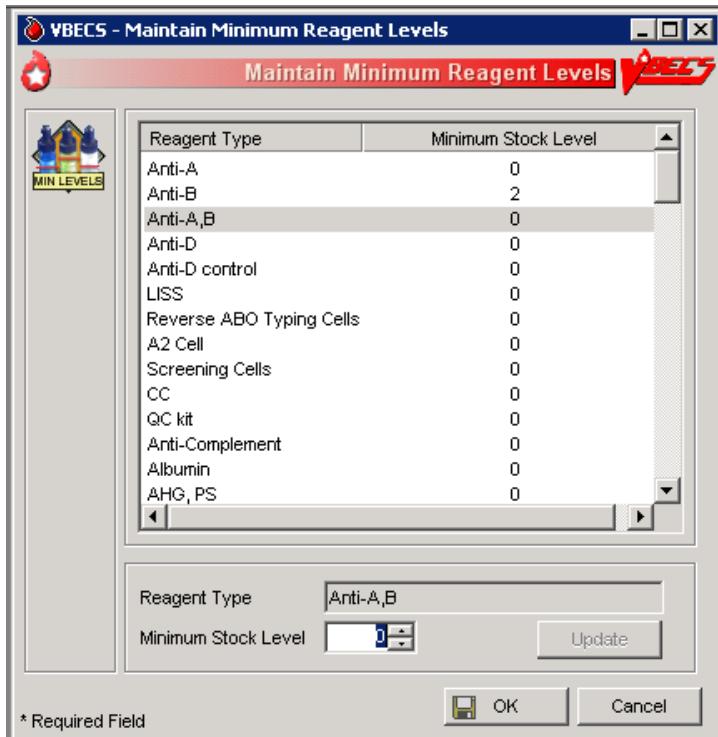

Supplies

The user documents the receipt of supplies, and updates and views the supply inventory.

Assumptions

- A user may maintain supply inventory only in the division logged into.

Outcome

- A supply inventory that is accessible only within the division of entry was established.

Limitations and Restrictions

- The division's modification processes may allow or require the user to select previously entered supplies for association with a blood unit's modification. An inventory of containers (transfer bags) and sterile connection device (SCD) wafers is *required* for various modification processes.
- There is no retrospective data entry or user selection for this option.
- The unsatisfactory portion of a shipment must be recorded individually.
- Each site must set a standard for recording quantity: it may use the number of cases or the number of units within the cases.

Additional Information

- Maintaining Blood Bank Equipment (Equipment: Log In Equipment, Equipment: Maintain Equipment) and Reagents (Log In Reagents, Update Inventory, View/Print Inventory, Maintain Minimum Levels) accommodate the maintenance of reagents and equipment, respectively.
- VBECS saves the lot number, manufacturer name, and expiration date of supplies used in a unit's modification process as part of a unit's modification record when associated during the modification transaction.
- VBECS creates a list of manufacturers as they are entered and makes the list available for selection. Users may not edit entries.

User Roles with Access to This Option

All users


Log in Supplies

The user records the receipt of specific supplies by type. The supply record includes the date the shipment was received, supply name, lot number, expiration date, manufacturer, and number of items received in a shipment.

User Action	VBECS
1. Select Tools from the main menu. Select Supplies .	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Displays options for processing supplies.
2. Select Log In Supplies .	<ul style="list-style-type: none">• Displays fields for entering invoice and supply information.
3. Invoice Information: Enter the invoice number in the Invoice Number field. Select or edit the date and time in the Date Received field.	<p>NOTES</p> <p>More than one technologist at a time may process a large shipment of supplies; only one technologist at a time may process individual supplies in the shipment.</p>

User Action	VBECS
Enter the name of the manufacturer in the Manufacturer field, or select a name from the drop-down list.	<p>The default date received is the current date and time. The user may edit it to a past date and time only.</p> <p>When the supply type and lot number were already entered for this invoice, VBECS warns that the user is trying to enter a supply that was already logged into the system for this invoice. Yes takes the user to the “Update Supply Inventory” screen where VBECS displays the supply type and lot number for update. No clears all fields on the screen, allowing the user to log in another supply.</p>
4. Supply Information:	<p>NOTES</p> <p>The default lot expiration date is one minute before midnight (23:59) of the current date.</p>
Enter the name of the manufacturer and supply type in the Manufacturer and Supply Type fields, or select from the drop-down lists. Enter the lot number in the Lot Number field. Select or edit the date and time in the Lot Expiration Date field. Enter (or select) a quantity (1–999) in the Quantity Received field.	
5. Visual Inspection Information:	<p>NOTES</p> <p>Satisfactory inventory is available for later selection, unsatisfactory inventory is not.</p> <p>When a user selects Other from the drop-down list in the Comment field, entry of additional text in the Details field is required.</p>
6. Repeat Steps 3–5 to log in more supplies, as desired.	<p>NOTES</p> <p>Click Clear to clear all but the Invoice Number, Date Received, Manufacturer, and Lot Expiration Date fields.</p>
7. Click OK to save. Click Yes to confirm the save and exit.	<ul style="list-style-type: none"> • Requests confirmation to update the database.

Figure 63: Log In Supplies

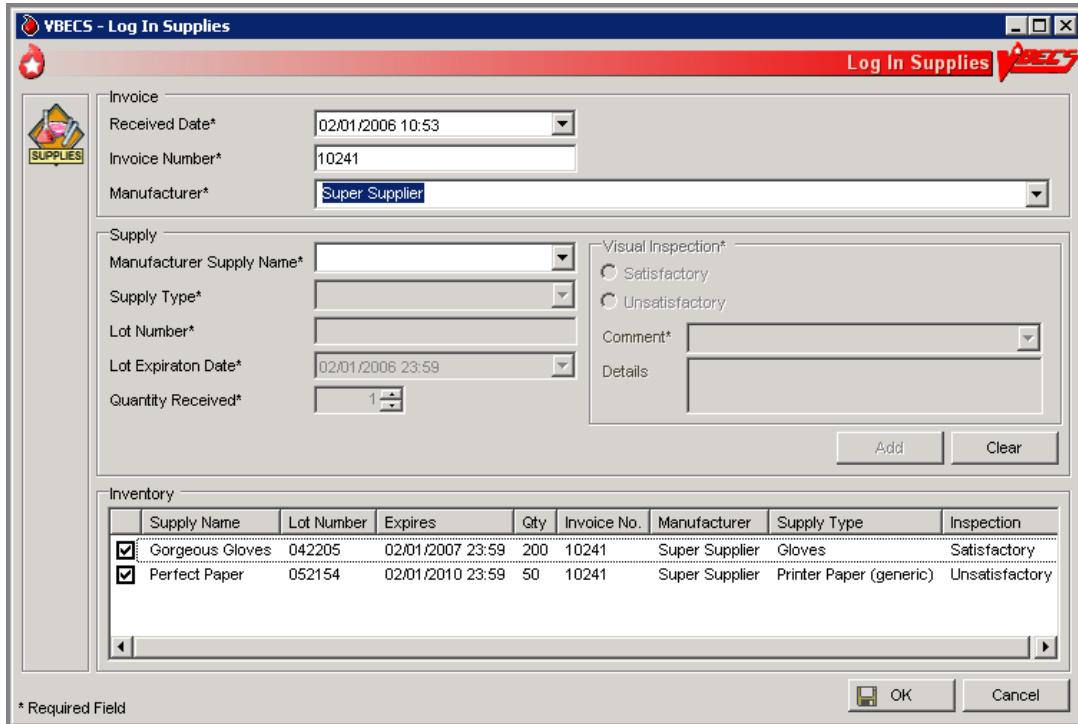


Update Supply Inventory

The user changes the number of available items in the supply inventory.

User Action	VBECS
1. Select Tools from the main menu. Select Supplies .	<ul style="list-style-type: none"> Displays options for processing administrative functions. Displays options for processing supplies.
2. Select Update Supply Inventory . Click the Include Expired Supplies check box to include expired supplies in the list.	<ul style="list-style-type: none"> Lists existing satisfactory supplies by default.
3. Select a supply to edit. To view supplies with a specific lot number, enter the lot number in the Lot # field and click the magnifying glass .	NOTES — Click a supply type to reveal manufacturer supply names. Click a manufacturer supply name to reveal lot numbers.
4. Click a lot number to enable the Quantity Available field. Enter (or select) a quantity that is less than or equal to the current entry in the Quantity Available field. Select or edit the date and time in the Date Changed field (Figure 64).	NOTES — The default change date is the current date and time. The user may edit it to a past date and time only. When a user selects Other from the drop-down list in the Reason for Change field, he may enter additional text in the Comments field.

User Action	VBECS
5. Repeat Steps 3 and 4 to edit more supplies, as desired.	
6. Click OK to save. Click Yes to confirm the save and view the updated list.	
7. Click Print Supply Inventory to print the list of supplies.	<ul style="list-style-type: none"> • Displays the Supply Inventory window (see View/Print Supply Inventory).
8. Click Cancel to exit.	

Figure 64: Update Supplies

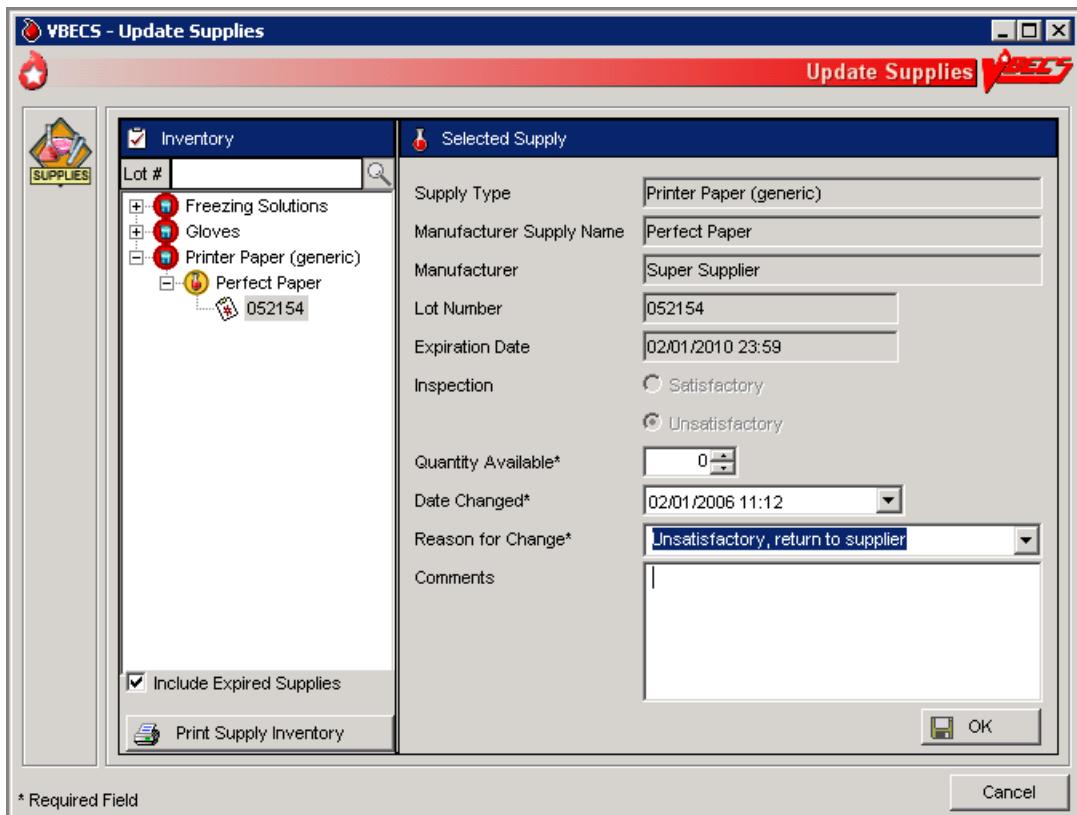

View/Print Supply Inventory

The user views inventory online based on various search criteria and prints the search results.

User Action	VBECS
1. Select Tools from the main menu. Select Supplies .	<ul style="list-style-type: none"> • Displays options for processing administrative functions. • Displays options for processing supplies.
2. Select View/Print Supply Inventory .	<ul style="list-style-type: none"> • Lists supply types.
3. Click one or more check boxes in the Supply Type list to select supply types.	<ul style="list-style-type: none"> • Displays or allows the user to print a report of all current in-date, expired, and unsatisfactory supply inventory and minimum stock levels as a read-only display based on the selection criteria entered: <ul style="list-style-type: none"> o Date range (the user enters a date to begin and end the search) o Supply type o In-date (active) Supplies

User Action	VBECS
	<ul style="list-style-type: none"> o Expired (inactive) Supplies o Expiration date of the supply o Date received <p>NOTES —————</p> <p>A user may view the supply inventory by one, some, or all of the selection criteria.</p>
4. Click OK to accept the displayed selection criteria.	<ul style="list-style-type: none"> • Compiles the report and displays an option to print or exit the report.
5. Click Print Supply Inventory to print the list of supplies.	<ul style="list-style-type: none"> • Sends the report to print on the printer selected by the user. <p>NOTES —————</p> <p>VBECS displays and prints the report in the supply type order : in-date by lot number in reverse chronological order of expiration date, followed by expired and unsatisfactory supplies. Details included in the report:</p> <ul style="list-style-type: none"> • Supply Type (Unit of Measure) • Manufacturer • Manufacturer's Supply Name • Lot Number • Expiration Date • Number Available • Inspection upon receipt • Invoice Number • Received Date/Time • Received by Tech Name • Comments
6. Print another report or click Cancel to exit.	

Equipment

The user registers blood bank instruments and other equipment (refrigerators, freezers, thawing bath, incubators, centrifuges, etc.) and records associated maintenance activities.

Assumptions

- A user may maintain equipment inventory only in the division currently logged into.

Outcome

- An equipment maintenance record, accessible within the division of entry, was established.

Limitations and Restrictions

- VBECS does not check for duplicate entries of equipment names and identifiers.
- There is no retrospective data entry or user selection for this option.

Additional Information

- VBECS does not calculate inventory levels for supplies or equipment.

User Roles with Access to This Option

All users

Equipment: Log In Equipment

A user adds new equipment or maintains existing equipment. After equipment is registered, the user must record periodic maintenance. The equipment records are updated and retained according to local procedure and policy. Once an equipment record is established, it may be selected as equipment used during a blood product modification.

User Action	VBECS
1. To add new equipment, select Tools from the main menu. Select Equipment .	<ul style="list-style-type: none"> Displays options to maintain equipment maintenance.
2. Select Log In Equipment to create an equipment record or Maintain Equipment to record maintenance (see Equipment: Maintain Equipment).	<ul style="list-style-type: none"> Displays an empty equipment property template. Displays all in-service equipment items for the blood bank associated with the user's division and offers the option to display out-of-service equipment.
3. Equipment Info area: Enter the equipment information in the appropriate fields. Enter or select an equipment type from the drop-down menu. Installation Info area: Select or edit the date and time in the Date Received/Installed field. Click the Expiration/Renewal Date check box when there's an expiration or renewal date, and edit the date as above. Manufacturer Info area: Enter the manufacturer name or select a manufacturer from the drop-down list. Enter the manufacturer's contact information in the appropriate fields or click the find button to select the manufacturer's address. Click Add to add the equipment to inventory (Figure 65).	<ul style="list-style-type: none"> Displays fields for entering equipment name, ID, serial number, and type in the Equipment Info area. Displays fields for entering the date the equipment was received or installed, warranty or service contract number, and expiration or renewal service contract date in the Installation Info area. Displays fields for entering the manufacturer name, address, and phone in the Manufacturer Info area.
4. Repeat Step 3 to add new items to the equipment list. Click OK to save.	<ul style="list-style-type: none"> VBECS maintains a historical record for each piece of equipment entered.
5. Click Yes to confirm and exit.	

Figure 65: Log In Equipment

The screenshot shows the 'VBECS - Log In Equipment' window. On the left, there's a sidebar with a magnifying glass icon and the word 'EQUIPMENT'. The main area has two main sections: 'Equipment Info' and 'Manufacturer Info'. 'Equipment Info' includes fields for 'Equipment Name*', 'Local Equipment ID', 'Serial Number*', 'Equipment Type*', 'Date Received/Installed*' (set to 08/11/2011 14:36), 'Warranty/Service Contract No.', and an optional 'Expiration/Renewal Date'. 'Manufacturer Info' includes fields for 'Manufacturer Name*', 'Address Line 1' (with a 'Find' button), 'Address Line 2', 'City', 'State', 'Zip', and 'Phone'. Below these are 'Add' and 'Clear' buttons. A large grid titled 'Equipment Logged In' displays columns for Equipment Name, Serial No., Type, Manufacturer, Received, and Contract No. At the bottom, there are 'OK' and 'Cancel' buttons, and a note: '* Required Field'.


Equipment: Maintain Equipment

The user maintains an online equipment record to record periodic maintenance. Equipment records are updated and retained according to local procedure and policy.

User Action	VBECS
1. To maintain equipment, select Tools from the main menu. Select Equipment.	<ul style="list-style-type: none"> Displays options to maintain equipment maintenance.
2. Select Maintain Equipment to record maintenance or Log In Equipment to create an equipment record (see Equipment: Log In Equipment).	<ul style="list-style-type: none"> Displays entered in-service equipment items for the blood bank associated with the user's division and offers the option to display out-of-service equipment.
3. In the Equipment tab, enter the equipment serial number (or enter part of the serial number and click the magnifying glass to view serial numbers), or click [+] in front of each equipment type and name to view equipment items. Select the serial number of the item whose record will be maintained.	<ul style="list-style-type: none"> Displays fields for entering the type of maintenance event performed, the date the maintenance was performed, the results of the event, the next scheduled event, the next scheduled maintenance date, and pending maintenance issues in the Maintenance tab, and allows the user to select a maintenance event type and event details. Displays entered information.

User Action	VBECS
Click the Include Out of Service Equipment check box, if desired.	
4. Review and edit the equipment data in the Equipment tab.	<ul style="list-style-type: none"> Displays entered information.
5. Select the Maintenance tab . To add a maintenance event, enter the maintenance type or select it from the drop-down list. Select or edit the date and time in the Maintenance Date field. Enter details in the Maintenance Results and Pending Issues fields. Click Add to add the maintenance event to the Maintenance Events list (Figure 66). Click Print to print the record.	<ul style="list-style-type: none"> Displays recorded maintenance information and allows the user to update information. Displays an option to print an Equipment Report.
6. Click the Out of Service Equipment check box to remove equipment.	<ul style="list-style-type: none"> Removes the equipment from the list of currently maintained equipment.
7. Repeat Steps 3–6, as desired.	
8. Click OK and Yes to save and continue, or No to return to the screen and continue. Click Cancel to exit without saving.	<ul style="list-style-type: none"> Prompts the user to confirm the save. Saves the updates to the equipment record.

Figure 66: Maintain Equipment

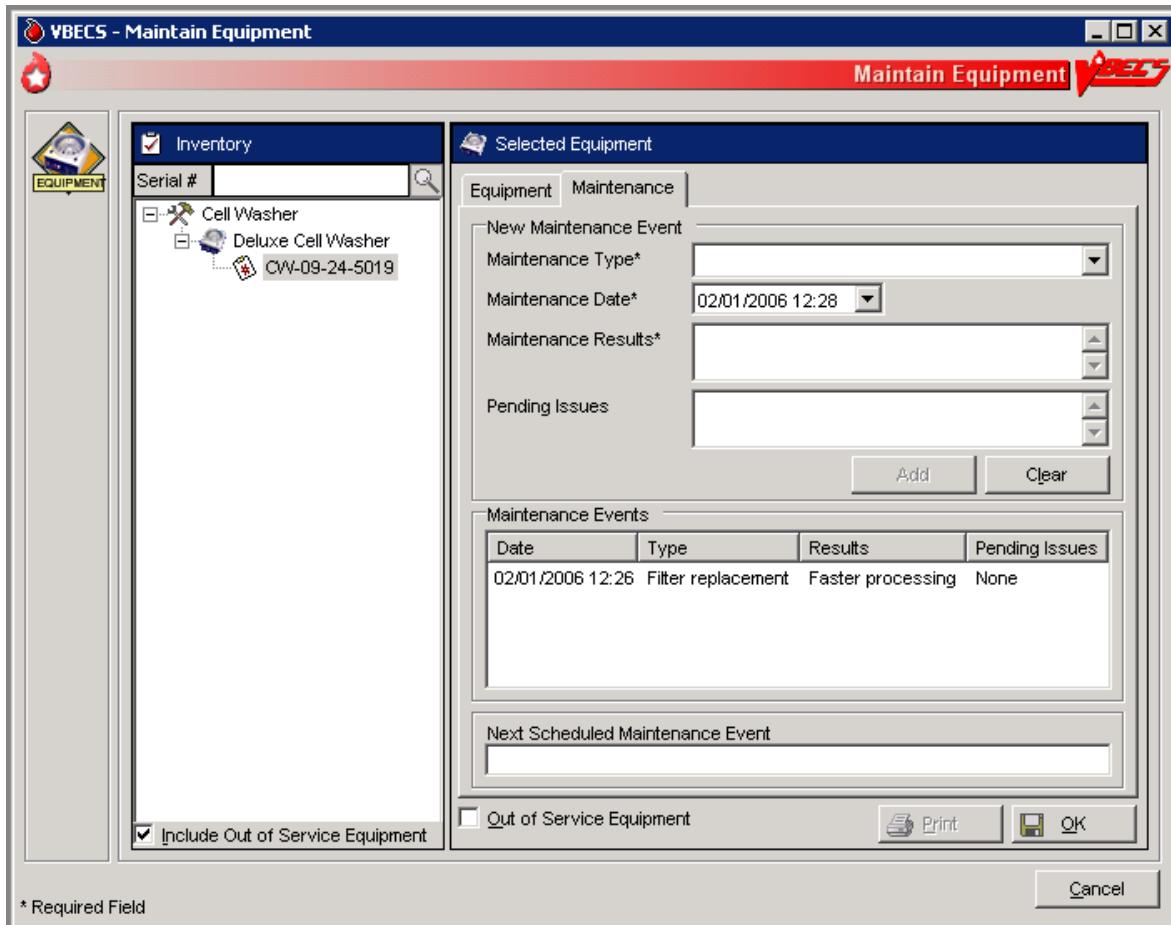

This page intentionally left blank.

Component Processing

Shipments

Incoming Shipment

A shipment handled in this option contains one or more containers of human blood products.

The user may use a barcode scanner to read information from a unit blood bag (recommended) or enter this information. The user reviews and confirms the information before adding a unit to the database.

Assumptions

- The connection to VistA is active.
- The division is configured.
- The shipper is active.
- Units are inspected and deemed satisfactory for receipt before they are entered in VBECS.

Outcome

- A shipment or partial shipment of blood units was processed and added to the blood bank inventory.
- The initial status for each unit was set in accordance with component class.

Limitations and Restrictions

- VBECS does not compare the ABO/Rh of an autologous unit or display a message when the ABO/Rh is incompatible with the restricted for patient.

When a product type, code, or local code is not in VBECS or a blood supplier creates a product code, contact the Help Desk. See Customer Support for contact information.

- The user must enter the donation type for Codabar-labeled units.
- The unit expiration [Maximum Storage Time (MST)] is set to the longest possible storage time based on product type, core condition, attributes, and available references as a guide. The MST is set at a higher level than an individual collection product code. VBECS does not know the specific collection details. It is the responsibility of each user to accurately enter and process any one blood unit based on its expiration date/time as labeled by the collection facility.
- When a user inadvertently enters a unit ID that is longer than seven characters for a Codabar unit, or 13 characters for an ISBT 128 unit, VBECS accepts the entry but will not associate it with the scanned unit ID.
- When blood products intended for quarantine are scanned, VBECS does not automatically quarantine them: the user must do so.
- The display of Special Testing for a unit is limited to 253 characters in Incoming Shipment. All of the special testing is saved in VBECS and is available on reports.
- The display label that is built when logging a blood product into VBECS does not change from volunteer donor to autologous when the donation type is changed.
- VBECS does not accept scanned information for some modified ISBT 128 blood units. VBECS accepts entered information in these instances.
- VBECS does not recognize scanned labels for para-Bombay or Bombay ISBT units.

- VBECS does not check whether blood products with the same unit ID and different product codes are the same blood type.
- VBECS does not allow duplicate IDs for the same product code.
- VBECS reads Codabar prefixes and translates them into digits. When this results in identical Codabar unit IDs, keep the original unit and return the other unit to the supplier.
- When a user scans a product code or expiration date barcode in the Unit ID field, VBECS does not display the entered information or an error icon.
- When a user scans an ISBT product ID with “0” as the donation type (sixth character), VBECS indicates that the donation type is not specified. When the user selects a donation type from the drop-down menu, VBECS replaces “0” with a recognizable character in the product code text box and stores it as such. When the user rescans the product ID, he must manually replace “0” or search using the unit ID alone, or the unit ID and the first five characters of the product ID.
- Blood product code entries in the blood product table include those from ICCBBA for ISBT 128 and from the AABB Codabar table. Product codes created after the release of VBECS cannot be entered until the updated tables are incorporated through patches to VBECS. Refuse acceptance of the unit, if possible. Users may not edit the blood product table.
- ISBT 128 tables include blood product codes labeled as not for manufacturing or transfusion. VBECS does not restrict such products from being entered, issued, or transferred (local policy controls their use), but does not accommodate blood modification for these products.
- A user cannot bring in pooled ABO or pooled Rh types for plasma or FFP. Cryoprecipitate is **not** affected.

Additional Information

- A unique blood unit record is defined by the following data elements:
 - VBECS Division
 - Unit ID
 - Product type code (Codabar Only)
 - For example: Product Code: 12000 will present a choice of Product type code:
 - Platelets
 - Washed Platelets
 - Product code
 - Codabar (5 digits)
 - ISBT 128 (8 alphanumeric characters).
- VBECS stores scanned ISBT 128 product codes and displays them as eight characters:
 - Characters 1–5: product type
 - Character 6: donation type
 - Characters 7 and 8: division code (split)
- VBECS may display a scanned Codabar product code as seven characters.
- When ISBT 128 units received during downtime must be issued, it is recommended that the user record the check character. When a blood unit is to be logged in after downtime and the unit is unavailable for scanning, the user must enter the ISBT unit ID and the case-sensitive check character (the character in the square following the unit ID).
- To calculate a check character, see the Quick K Calculator on the ICCBBA Web page under ISBT 128, Other Tools. Follow the instructions to download the application.
- Each division may have only one blood unit record of the same product type, product code, and unit ID. A unit shipped out (transferred) and subsequently shipped to (received by) a division in the same database will maintain individual records of the unit ID and product code for each division.

- Several data elements, including the unit ID and product code, define a blood unit record as unique. The unit record also includes expiration date and time, blood type, and donation type. These data elements are not editable.
- Some suppliers relabel blood products and do not affix their own FDA Registration Number to those units. It is recommended that the user obtain FDA Registration Number barcodes from their suppliers and use them for scanning the shipper.

User Roles with Access to This Option

All users

Incoming Shipment

In a “transfusion-only” facility, units must be restricted for use by a specific patient while an incoming shipment is processed. VBECS identifies the original collection facility to determine the presence of an eye-readable prefix and/or alphabetical characters in a Codabar unit ID. A user may enter unique unit information regarding special testing and antigen phenotyping. The user confirms that all entered units from the shipment were added to the database in the blood bank.

User Action	VBECS
1. Select Shipments from the main menu. Select Incoming Shipment .	<ul style="list-style-type: none"> • Displays options for processing blood product shipments. • Displays invoice information fields for entering shipment information.
2. Scan the invoice number barcode or enter the invoice number. To search for an existing invoice, click the find button. Enter an invoice number and/or select a shipper from the drop-down list, then click Search . Click an invoice number to select it and click OK .	<ul style="list-style-type: none"> • When an identified invoice does not exist, asks the user whether he wants to create an invoice and add units to it. • Displays the option to search for existing invoices by invoice number and/or by active shipper. When no match is found, notifies the user. <p>NOTES</p> <p>When a user chooses to select an existing invoice to add units, VBECS ensures that the shipper is marked as active. VBECS displays and allows the user to edit existing invoices associated with active shippers.</p> <p>To correct an invoice associated with an inactive shipper, the user must reactivate the shipper.</p> <p>Multiple users may process different units from a shipment.</p> <p>When more than one user intends to work on an invoice, one user must create the invoice. The other users must select that invoice and add units to it.</p>
3. Select a source from the Source (Shipper) drop-down list. If needed, edit the date received in the Date Received field.	<p>NOTES</p> <p>The date and time a shipment is received cannot be in the future.</p> <p>After a user enters or selects an invoice number, the user may enter information in the shipper, date received, and processing technologist fields.</p> <p>During retrospective entry, the user must enter the shipment date and time before entering any other unit transaction.</p>
4. Click the ISBT or Codabar radio button to select a labeling type, or	<ul style="list-style-type: none"> • When the barcode is scanned, determines the labeling type (autodiscriminates).

User Action	VBECS
scan the unit ID barcode. Click OK again to continue.	
5. Click Yes to add units to the invoice or click No to select another invoice.	<ul style="list-style-type: none"> When a user selects an existing invoice, displays the shipment information and blood unit details associated with the invoice.
6. Select a unit to add. For a Codabar-labeled unit, scan the FDA Registration Number barcode or enter the FDA Registration Number. For an ISBT 128-labeled unit, continue at Step 7.	<ul style="list-style-type: none"> Displays shipment data and prompts the user to scan individual units. Displays information corresponding to the FDA registration barcode scanned or allows a user to search for the FDA Registration Number and facility information (all users may activate a facility). When a facility receives a blood product from a shipper for the first time, prompts the user to activate the blood product for that shipper. <p>NOTES</p> <p>VBECS may display a scanned FDA Registration Number with two additional characters.</p> <p>For United Blood Services (UBS) facilities other than Scottsdale, see Appendix L: United Blood Services (UBS) Facility Barcodes.</p> <p>When VBECS does not recognize an FDA Registration Number, it warns the user. The user may edit or reenter the number. If VBECS still does not recognize the number, the user must add the unknown facility to continue.</p>
7. For an ISBT 128-labeled unit, scan the label, or enter the unit ID and check character (in that order).	<ul style="list-style-type: none"> Disables the FDA Registration Number field. Displays information corresponding to the unit ID barcode scanned. When a user enters an ISBT 128-labeled unit ID, VBECS prompts the user to enter the check character from the blood bag. <ul style="list-style-type: none"> When a user scans ISBT 128 labels units and the cursor is in an inappropriate field, VBECS warns that the incorrect barcode was scanned. <p>NOTES</p> <p>For Codabar-labeled units:</p> <p>VBECS may display a scanned FDA Registration Number with two additional characters.</p> <p>When VBECS does not recognize an FDA Registration Number, it warns the user. The user may edit or reenter the number. If VBECS still does not recognize the number, the user must add the unknown facility to continue.</p> <p>When no FDA Registration Number barcode was scanned, VBECS displays an error message to remind the user to scan it so that VBECS can properly interpret the Codabar unit ID.</p> <p>When a user scans a unit ID barcode but enters other unit information, VBECS marks the unit record as scanned and does not indicate that other data field information was entered during incoming shipment of that unit.</p> <p>When a unit ID or its prefix does not match the format defined for the collection facility, VBECS warns the user.</p> <p>The ID of a unit associated with a local supplier with a facility prefix and alphabetic characters has two numbers for the facility prefix, and one or two alphabetic characters. When there is one</p>

User Action	VBECS
	<p>alphabetic character, at least six numbers must follow. When there are two alphabetic characters, at least five numbers must follow. Codabar unit IDs cannot contain "A," "B," "D," "I," "O," or "U."</p> <p>For ISBT 128-labeled units: When the user-entered check character from the label does not match the VBECS-generated check character, VBECS warns the user.</p> <p>When VBECS does not find a matching entry for the ICCBBA Registration Number on the local supplier table, the user must activate or add a local supplier before continuing.</p>
8. Scan the product code barcode or enter the product code.	<ul style="list-style-type: none"> • When a user scans a product code barcode that is not activated in the user's division, prompts the user to activate the code. (See Blood Products.) • Displays the corresponding product code, product name, and (for ISBT) donation type (see Step 11). • When a unit ID of the same product type and code exists in a division, the expiration date is in the future, and there is no record of inactivation or transfer of that unit ID and product code in the division's record: <ul style="list-style-type: none"> ◦ Warns that the blood unit already exists in the division. ◦ Clears the data entry fields. The user may enter a different unit ID and product code. • When a unit ID of the same product type and code exists in a division and the unit in the division has an expiration date in the past, warns that the blood unit already exists in the division. • When a unit ID of the same product type and code exists in a division, the unit status is "Transferred," and the ABO/Rh and unit expiration date are identical to the existing unit of record: <ul style="list-style-type: none"> ◦ Makes the existing unit record available. ◦ Warns that the user is attempting to reenter a unit in inventory. No clears the unit. Yes and saves the unit. ◦ Updates the unit status to "Limited" or "Available." ◦ Updates the unit record to reflect the reentry of the unit in the division's inventory. The user may add new information about this unit through Edit Unit Information.

NOTES

Scanned Product Codes

When a unit is transferred between facilities in a consolidated database or is reentered after transfer to an unrelated facility, VBECS records all unit activity.

Codabar units may not be added to an ISBT 128 shipment in the same session and ISBT 128 units may not be added to a Codabar shipment in the same session. However, a completed incoming shipment invoice may contain both types of units.

When a blood product code in the RED BLOOD CELLS, WHOLE BLOOD, or APHERESIS RED BLOOD CELLS ICCBBA component class is entered, VBECS enables data entry of sickle cell negative and antigen typing. The cytomegalovirus (CMV) negative data entry option is always enabled.

User Action	VBECS
	<p>When a user enters a product code that VBECS does not recognize, VBECS warns that the blood product was not found and that it does not support entry of new blood products.</p> <p>OK returns the user to a blank screen to enter a different product. The user may enter a different product code.</p> <p>VBECS requires the full 8-digit ISBT 128 product code. The user should enter the 8-digit ISBT 128 product code as it appears on the blood unit label of the actual blood unit being entered into VBECS. When a user does not enter the sixth through eighth characters, VBECS automatically completes the required information with defaults, V and 00, respectively. This may be changed by clearing the unit information and entering the complete 8-digit ISBT 128 unit product code.</p> <p>VBECS assigns a volume to the unit based on the average volume field from the Blood Product Table entry corresponding to the unit's blood product code.</p> <p>For each product type in a shipment, site parameters must indicate that the product is marked as active for the division and that a cost is associated with the shipper. If not, VBECS prompts the user to mark the component as active.</p> <p>For Codabar-labeled units: VBECS displays a scanned Codabar product code as seven characters.</p> <p>In some instances, Codabar product codes are not unique when they are modified.</p> <p>For ISBT 128-labeled units: VBECS stores scanned ISBT 128 product codes and displays them as eight characters:</p> <ul style="list-style-type: none"> • Characters 1–5 are used to search the ISBT Blood Product Table to identify the product type. • Character 6 indicates donation type. • Character 7 indicates a division at the time of collection. • Character 8 indicates an additional aliquot.
9. For Codabar- and ISBT 128-labeled units, scan the ABO/Rh number barcode, enter the ABO/Rh, or select the ABO/Rh from the drop-down list.	<ul style="list-style-type: none"> • Displays the corresponding ABO/Rh. • Allows the selection of only ABO/Rh types acceptable for ICCBBA component classes. • Enables the ABO/Rh field when a valid product code is entered.
10. Scan the expiration date barcode or enter the expiration date and time when an expiration time is required.	<ul style="list-style-type: none"> • Displays the corresponding expiration date and time. <p>NOTES ——————</p> <p>The default unit expiration time is 23:59.</p> <p>VBECS compares the entered blood product's expiration date to the maximum storage date and time limit for the product type. The date and time of shipment receipt are used as the collection date and time. When the expiration date entered exceeds the defined maximum storage time, VBECS displays an error message and clears the field. The user may reenter an expiration date and time for the product.</p>

User Action	VBECS
	<p>Products with a maximum expiration date and time of \leq24 hours default to an expiration time \leq24 hours from the time of the incoming shipment. Product types with a maximum expiration date and time of $>$ 24 hours have a default expiration time of 23:59.</p> <p> When a unit expiration date precedes the incoming shipment date and time, VBECS emits an audible alert, warns that the user is attempting to enter an expired unit, and asks whether the user wishes to continue.</p> <p>No clears the unit expiration date field. The user must enter a new, acceptable expiration date. Yes requires the user to enter a comment.</p> <p>VBECS captures details of this override in an Exception Report.</p>
<p>11. For Codabar-labeled units, accept “voluntary allogeneic” or select another donation type.</p> <p>For ISBT 128-labeled units, accept the donation type based on character 6 of the product code (cannot be edited) or, when character 6 is “0,” indicate a donation type.</p>	<ul style="list-style-type: none"> Assigns the donation type to the unit or requires user input, as needed. For Codabar, assigns the “voluntary allogeneic” donation type to the unit. Requires data entry of the “reserved for” patient name and/or patient ID and the disease market testing status for “autologous” and “directed” donation types. <p>NOTES</p> <p>For ISBT 128, donation types other than “Special Testing Barcode” and “Not Specified” cannot be edited. When the donation type is “Special Testing Barcode” or “Not Specified,” the user must select a different donation type.</p> <p>Codabar barcodes do not include donation type codes, so a user must select a valid donation type. When a donation type requires restriction (directed or autologous), the user must associate the unit with a patient name.</p> <p>When a site is a transfusion-only facility, each unit processed must be restricted to a specific patient regardless of donation type.</p>
<p>12. Click the find button next to the Restricted For Patient field to indicate the patient for whom the unit is reserved, as needed.</p>	<ul style="list-style-type: none"> Allows the user to search for active VistA patients to associate with the unit. Requires data entry of a patient name and/or patient ID to associate the unit with that patient when a unit has a “directed” or “autologous” donation type or the receiving facility is a “transfusion-only” facility. When no match exists, warns the user. When the unit is autologous, requires a user to indicate the disease marker testing status: <ul style="list-style-type: none"> NEG: unit fully tested and negative for all disease markers POS: unit tested and positive for one or more disease markers (biohazard) NFT: unit not fully tested for one or more disease markers (biohazard) <p>NOTES</p> <p>Displays a message when the user selects a patient with a blank</p>

User Action	VBECS
	<p>first or last name, or the name is too long. The user is instructed to contact the VBECS Administrator for assistance.</p> <p>The Patient Blood Availability Report does not include restricted units not assigned to a patient. To track such units, the user may enter restricted units' IDs in the Special Instruction field in Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions. When a user restricts a unit to a patient and VBECS does not find that patient in its database, VBECS searches VistA. When it finds a valid patient, VBECS adds the patient to the database.</p>
13. Do one or more:	<ul style="list-style-type: none"> Displays an option to save or clear unit information. <p>NOTES</p> <p>VBECS accommodates data entry of:</p> <ul style="list-style-type: none"> Special testing results (antigen phenotype results, CMV negative, or sickle cell negative). A single special typing cost (the sum of charges from \$0.00 to \$999.00, other than the base unit cost) to associate with a unit. <p>An antigen type may be positive or negative for any antigen listed as a VBECS corresponding antigen in Appendix B: Table 13: Antibody and Antigen Table: Irregular Antibodies. An antigen may not be positive and negative for the same unit.</p>
14. Review and click Save to accept the unit.	<ul style="list-style-type: none"> Asks the user whether he has inspected the units and that they are satisfactory. When the user clicks yes, the unit record is added to the division's inventory. When the user clicks no, the unit information is cleared and the unit is not entered into the division's inventory. Adds the unit to the database. Lists summary information for the accepted unit. Allows the user to begin processing a new unit. <p>NOTES</p> <p>If the unit is not satisfactory per the visual inspection, the user must click no. The unit is not entered into VBECS inventory and should be processed according to local procedure and policy.</p> <p>VBECS assigns unit statuses based on a site's configuration in Configure Division. For a:</p> <ul style="list-style-type: none"> "Full-service" facility, VBECS assigns "Limited" status to blood products that contain red cells and "Available" status to all other components. "Transfusion-only" facility, VBECS assigns "Available" status to the units. ABO/Rh confirmation testing is not required. <p>Units are available only to users in the division in which the units were processed (separate inventories are maintained in a multidivisional site).</p> <p>To edit incorrect saved information, deactivate the unit through Edit Unit Information, then reenter unit information in Incoming Shipment.</p>
15. Repeat until all units are entered.	<ul style="list-style-type: none"> Displays all unit data entries.

User Action	VBECS
16. Click Close .	

Figure 67: Process Incoming Shipment

Outgoing Shipment

The user processes and ships available blood units from his division to an outside facility, blood center, or hospital.

The user may edit return credits (additional fees or credits may be issued from one site to another) before processing the shipment.

VBECS generates a shipping invoice.

Assumptions

- Users are authorized to release units from inventory.
- Blood units being shipped:
 - Were entered in the division.
 - Are available for scanning.
 - Are located in the user's division.

Outcome

- An invoice for a completed shipment is closed and may not be modified.
- An invoice for a shipment that is not confirmed is open and may be edited.
- The units are not available in the division's inventory.

- When a unit is transferred to another division in the database, VBECS displays the history of the unit at both divisions. The unit is active in one division at a time.
- A credit for a returned product may be assigned.

Limitations and Restrictions

- A user may edit only unconfirmed invoices.
- The number of tests displayed on an invoice is limited to 8,000 characters per unit.
- When the user creates an invoice and cancels before saving it, VBECS excludes the invoice number from use.

Additional Information

- None

User Roles with Access to This Option

All users

Outgoing Shipment

A user ships blood units from his division's inventory to a facility that is not a transfusion location of his division. These units were removed from inventory and are available for scanning and packaging.

The user identifies the shipment destination, scans the unit ID and product code for the first unit in the shipment, and repeats this process for each unit added to the shipment. The user completes the data entry and then reviews and accepts the shipment information. When the shipment is confirmed, the user may print a shipping invoice.

VBECS generates an invoice that includes information for the blood units on that invoice.

A user may edit return credits (additional fees or credits may be issued from one site to another) before processing the shipment.

User Action	VBECS
<p>1. Select Shipments from the main menu.</p> <p>Select Outgoing Shipment.</p> <p>2. To process a new invoice, select New Invoice. (To process an existing invoice, continue at Step 3.)</p> <p>Enter the invoice number provided by the outside facility or leave the Invoice Number field blank.</p> <p>Select a shipment destination from the drop-down list.</p> <p>Click the Return to Shipper or Send to Unrelated Shipper radio button to select a shipment type (Figure 68).</p>	<ul style="list-style-type: none"> • Displays options for processing shipments. • Displays options for selecting a new invoice or an existing (recent) invoice. <p>NOTES</p> <p>To select an existing invoice while in the new invoice window, click Use Existing Invoice.</p>

User Action	VBECS
Continue at Step 4.	<ul style="list-style-type: none"> Displays a list of open (unconfirmed) invoices. <p>NOTES</p> <p>To select a new invoice while in the existing invoice window, click Create New Invoice.</p> <p>A user may edit only saved, not confirmed, outgoing shipment invoices.</p> <p>A user may select shipping destination (ship-to) facilities based on site parameters that define a shipper for a division.</p> <p>A shipper must be marked as an active collection facility for a division before a user may select it as a valid shipper for a blood product for that division. This includes marking the division facility itself as an active collection facility to accommodate the creation of new blood products through modification.</p>
3. To process an existing invoice, select Recent Invoice from the Outgoing Shipment menu. Select an invoice.	<ul style="list-style-type: none"> Displays fields for entering unit information. <p>NOTES</p> <p>When a unit selected for inclusion in a shipment is quarantined, VBECS warns the user and asks the user whether he wishes to continue.</p> <p>No clears the unit information and allows entry of a new unit or shipment process. Yes adds the unit to the shipment.</p> <p>VBECS allows the shipment or transfer of units located only in the same division as the user and in the blood bank (for “Limited,” “Available,” “Assigned,” and “Crossmatched” unit statuses only).</p> <p>When a unit is expired, VBECS warns the user.</p> <p>When units selected for shipment are assigned or crossmatched to a patient, VBECS warns the user and asks the user whether to release the unit before shipping.</p> <p>Yes releases all patient assignments from the unit prior to adding the unit to the invoice (no assignment information is displayed on the invoice). No indicates that the outgoing shipment invoice must include the unit’s patient association details.</p>
4. Click OK to save the information entered (Figure 68). 5. For each unit to be added to the shipping invoice, scan the unit ID and product code barcodes or enter them in the Unit ID and Product fields (Figure 69 and Figure 70). Click OK .	<ul style="list-style-type: none"> Displays corresponding unit data. <p>NOTES</p>
6. Respond to warnings, when indicated. Click Select to add the selected unit to the shipping invoice or click Clear to delete a unit (Figure 71).	<ul style="list-style-type: none"> Adds unit to or removes unit from a temporary list of selected units.
7. Repeat Steps 5 and 6 until all units for outgoing shipment are entered. Click Save to save the information	<ul style="list-style-type: none"> Displays shipment information. When a user clicks Cancel Shipment, asks whether the user wants to remove all units from the invoice. Yes deletes the invoice. No returns to the invoice.

User Action	VBECS
entered (Figure 72) and Yes to confirm the save. To cancel the shipment, click Cancel Shipment , then Yes to delete the invoice and return to Step 1 to begin again, or click No to return to the invoice (Figure 73).	NOTES <ul style="list-style-type: none"> A user may add or delete units from a temporary list before accepting the information entered.
8. To change the return credit amount for a unit in a “Return to Shipper” unconfirmed shipment, select a unit and right-click the mouse. Select Return Credit Update (Figure 74). Enter the return credit amount. Click OK to save the information entered or Cancel to return to the invoice.	NOTES <ul style="list-style-type: none"> Prompts the user to accept the new amount. Returns to the invoice and displays the return credit amount on the invoice. Displays the unit's financial information and allows the user to edit the return credit amount. NOTES <p>VBECS applies a return credit (equal to the return credit percentage of the original cost of the unit) to units processed in an outgoing shipment to the original shipper. Users may edit return credits for individual units.</p> <p>Instead of right-clicking the mouse, a user may press the application key () on the keyboard.</p>
9. To confirm the shipment, click Save , then Yes to save the invoice, or Cancel to return to the invoice (Figure 75). Click Print Invoice to generate and view an invoice. Click Printer to print an invoice to send with the shipment or click Close to return to the invoice. Click Cancel to exit or move to Step 9 to change a return credit amount.	NOTES <p>When units selected for shipment are currently restricted to a patient and the donation type of the unit is Voluntary Allogeneic, in addition to releasing the patient assignment, VBECS also removes the unit restriction in a Transfusion Only division configuration.</p>
10. Cancel or confirm the shipment as in Step 9.	<ul style="list-style-type: none"> Displays a “Confirm Shipment” option and stores the open invoice until the user chooses to complete the shipment. NOTES <p>When a shipment is saved but not confirmed, VBECS displays “NOT CONFIRMED” as the shipment date on the print version of the invoice.</p> <p>The shipping invoice includes patient assignment, restriction, quarantine indicator, biohazard information, antigen typing (specificity and interpretation only), special testing information, and (when it is in a unit record) the autologous unit testing status.</p>
11. When the shipment is confirmed, select or edit the date and time in the Shipment Date/Time field.	<ul style="list-style-type: none"> Displays shipment information, including the shipment date and time. When a user accepts the entered date and time, displays the invoice to be printed. Updates the status of each unit in the shipment to show that each unit

User Action	VBECS
Click OK to save the information entered or Cancel to return to the invoice.	<p>was transferred to the ship-to facility and assigns the shipment date and time as the date and time of the transaction.</p> <p>NOTES</p> <p>When VBECS confirms an outgoing shipment, it:</p> <ul style="list-style-type: none"> • Adds the shipment information to the unit record, changes the status of each unit in the shipment to "Transferred," and assigns the date and time of the transaction to the shipment. • Saves the technologist ID, the name of the printing technologist, the division, and the date and time the shipment was processed. It also saves the confirmation date and time as the date and time of shipment. The user may accept or edit the current date and time displayed (no future date or time may be entered). The unit status is changed to transferred upon invoice confirmation.
12. Click Print Invoice to generate and view an invoice as in Step 9, or click Cancel to exit.	<ul style="list-style-type: none"> • Displays the information to be printed on the invoice. • Prints the shipping invoice. <p>NOTES</p> <p>When a user prints an invoice and VBECS does not find information to include in the invoice, VBECS notifies the user and asks whether he wants to continue to print. The user may cancel or print the report, which will include "No pending orders found."</p>

Figure 68: Process Outgoing Shipment

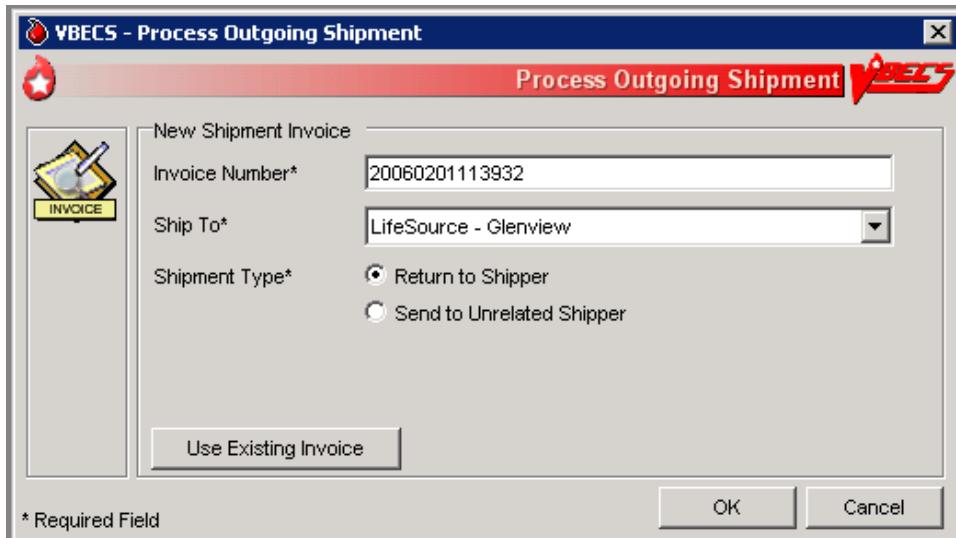


Figure 69: Select Units 1

The screenshot shows the 'VBECS - Process Outgoing Shipment' window. At the top, it displays 'Invoice 959', 'Ship To IVAMC North Chicago, IL', and 'Date Processed 08/11/2011 14:46'. Below this is a 'Unit Search' section with fields for 'Unit ID*' and 'Product*'. To the right is a preview area showing a barcode and the text 'Unit Number' and 'EXPIRES'. Underneath are fields for 'Details', including 'Donation Type', 'Date Received', 'Unit Status', and 'Volume (mL)'. A note states 'See circular of information for indications, contraindications, cautions and methods of infusion.' and 'This product may transmit infectious agents. Rx Only'. A red warning at the bottom right says 'PROPERLY IDENTIFY INTENDED RECIPIENT'. At the bottom left is a 'Selected Units' table with columns: Unit ID, Product, ABO/Rh, Expires, Donation, Status, Date/Time Received, Patient, Biohazardous, Quarantined, and Return Credit. The bottom of the window has buttons for 'Select', 'Clear', 'Cancel Shipment', 'Confirm Shipment', 'Print Invoice', 'Save', and 'Cancel'.

Figure 70: Select Units 2

The screenshot shows the 'VBECS - Select Units for Shipment' window. It features a 'Search Criteria' section with dropdown menus for 'Product Type' (listing options like Deglycerized RED BLOOD CELLS, FRESH FROZEN PLASMA, Thawed FRESH FROZEN PLASMA, Washed RED BLOOD CELLS, RED BLOOD CELLS, and WHOLE BLOOD), 'ABO/Rh' (listing O Pos, A Pos, B Pos, AB Pos, O Neg, A Neg, and B Neg, with O Pos checked), 'Special Testing' (listing CMV Negative, Sickle Cell Negative, and RBC Antigens), and date filters for 'Received Date' (Before and After) and 'Expiration Date' (Before and After). Below these is a table titled 'Units within Criteria : 4' with columns: Unit ID, Product, Code, ABO/Rh, Expires, Donation, Status, and ID. The table lists four units: 0100007A, 0100007B, 0100009, and 0100010, all marked with a checkmark in the first column. The bottom of the window has 'OK' and 'Cancel' buttons.

Figure 71: Select Units 3

The screenshot shows the VBECS - Process Outgoing Shipment window. At the top, it displays 'Invoice 959', 'Ship To VAMC North Chicago, IL', and 'Date Processed 08/11/2011 14:46'. In the center, there's a 'Unit Search' section with fields for 'Unit ID*' (W000000000009) and 'Product*' (E0291V00). Below this is a 'Details' tab showing 'Donation Type: Volunteer Allogeneic Donor', 'Date Received: 07/27/2011 12:33', 'Unit Status: Available', and 'Volume (mL): 350'. To the right, a large barcode is displayed with the identifier 'W000000000009' and a label indicating 'A Rh NEGATIVE VOLUNTEER DONOR'. At the bottom, there are buttons for 'Select' and 'Clear'. A 'Selected Units' table is shown below, which is currently empty. At the very bottom, there are buttons for 'Cancel Shipment', 'Confirm Shipment', 'Print Invoice', 'Save', and 'Cancel'.

Figure 72: Select Units 4

This screenshot is identical to Figure 71, showing the VBECS - Process Outgoing Shipment window. The main difference is that the 'Selected Units' table now contains one row with a checked checkbox. The row details are: Unit ID W000000000009, Product RBC CPD>AS1, ABO/Rh A Neg, Expires 08/06/2011, Donation Volunteer Allogeneic Donor, Status Available, and Date/Time Received 07/27/2011 12:33.

Figure 73: Select Units 5

Unit ID	Product	ABO/Rh	Expires	Donation	Status	Date/Time Received
W000000000009	RBC CPD>AS1	A Neg	08/06/2011	Volunteer Allogeneic Donor	Available	07/27/2011 12:33

Figure 74: Return Credit Update

Unit ID	Product	ABO/Rh	Expires	Donation	Status	Date/Time Received
W000000000009	RBC CPD>AS1	A Neg	08/06/2011	Volunteer Allogeneic Donor	Available	07/27/2011 12:33

Figure 75: Confirm Invoice

The screenshot shows the 'Process Outgoing Shipment' window with the title bar 'VBECS - Process Outgoing Shipment'. The main area contains fields for 'Invoice' (959), 'Ship To' (VAMC North Chicago, IL), and 'Date Processed' (06/11/2011 14:46). A 'Unit Search' section includes 'Unit ID*' and 'Product*' search boxes. Below these are 'Details' fields for 'Donation Type', 'Date Received', 'Unit Status', and 'Volume (mL)'. To the right, there's a barcode area with instructions: 'Property Identify Intended Recipient', 'Do Not Dispense by patients, non-pharmacy technicians and', and 'Caution: Federal Law prohibits dispensing without a prescription'. A 'Selected Units' table lists one unit: W000000000009, RBC CPD>AS1, A Neg, 08/06/2011, Volunteer Allogeneic Donor, Available, 07/27/2011 12:33. At the bottom are buttons for 'Cancel Shipment', 'Confirm Shipment' (with a checkmark), 'Print Invoice', 'Save', and 'Cancel'.

Figure 76: Modify Return Credit

The screenshot shows the 'Return Credit' window with the title bar 'VBECS - Return Credit'. The main area is titled 'Modify Return Credit'. It features a 'Unit Financials' section with fields for 'Unit ID' (W038403000041), 'Unit Cost' (\$1.00), 'Special Testing Cost' (\$), and 'Return Credit Amount*' (\$2.00). A note at the bottom left says '* Required Field'. At the bottom are 'OK' and 'Cancel' buttons.

Figure 77: Confirm Shipment

The screenshot shows the 'Process Outgoing Shipment' window. At the top, there are fields for 'Invoice' (959), 'Ship To' (VAMC North Chicago, IL), and 'Date Processed' (08/11/2011 14:46). Below these are sections for 'Unit Search' (with 'Unit ID*' and 'Product*' fields) and 'Details' (with fields for 'Donation Type', 'Date Received', 'Unit Status', and 'Volume (mL)'). To the right, there is a barcode area with instructions: 'Property of Intended Recipient', 'Do Not Dispense by Patients, Interpatients, Nurses and Other Health Care Personnel', and 'Caution: Federal Law prohibits dispensing without a prescription'. Below the search section is a table titled 'Selected Units' containing one row:


	Unit ID	Product	ABO/Rh	Expires	Donation	Status	Date/Time Received
<input checked="" type="checkbox"/>	W000000000009	RBC CPD>AS1	A Neg	08/06/2011	Volunteer Allogeneic Donor	Available	07/27/2011 12:33

At the bottom, there are buttons for 'Select', 'Clear', 'Cancel Shipment' (with a red X icon), 'Confirm Shipment' (with a checkmark icon), 'Print Invoice' (with a printer icon), 'Save', and 'Cancel'.

Figure 78: Confirm Shipment Date

The screenshot shows the 'Process Outgoing Shipment' window with a 'Shipment Date Entry' dialog box overlaid. The dialog box has a title 'Shipment Date Entry' and contains fields for 'Invoice Number' (959), 'Ship-To Facility' (VAMC North Chicago, IL), 'City, State' (North Chicago, IL), and 'Shipment Date*' (08/11/2011 14:54). The background window shows the same 'Selected Units' table as Figure 77, with the same single row selected. The bottom buttons are identical to Figure 77: 'Cancel Shipment' (red X), 'Confirm Shipment' (checkmark), 'Print Invoice' (printer), 'Save', and 'Cancel'.

Figure 79: Confirmation of Invoice

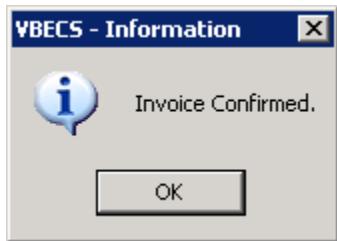

Edit Invoice Text

VBECS includes a shipping invoice. The user may edit the explanatory text above the signature on the last page of a shipping invoice.

Assumptions

- Information on the outgoing shipment invoice is generated during Outgoing Shipment.

Outcome

- A division customized the text displayed on an outgoing shipment invoice.

Limitations and Restrictions

- None

Additional Information

- The default explanatory text reads: "I certify that the blood products listed were properly maintained, in accordance with the Code of Federal Regulations, while in storage at this institution. Components were inspected when packed for shipment and found to be satisfactory in color and appearance."

User Roles with Access to This Option

►►►► Enhanced Supervisor

Edit Invoice Text

The blood bank user needs to edit the shipping invoice text box to conform to changes in regulations or division requirements. When the edit is complete, the user saves the altered text for future shipping invoices. An audit trail of changes is maintained.

User Action	VBECS
1. Select Shipments from the main menu. Select Edit Invoice Text .	<ul style="list-style-type: none">Displays options for processing blood product shipments.Displays the current text for outgoing shipment invoice.
2. Edit the text. Click OK to save the text.	<ul style="list-style-type: none">Displays entered data and changes for review.Prompts the user to save or reject the changes.
3. Review the changes. Click Yes to confirm the save and exit.	<ul style="list-style-type: none">Saves the new text and makes it available for outgoing shipment invoices.

Unit Search Screen

VBECS displays the Unit Search screen (Figure 80) when the user clicks the **find** button to select a unit through some of the Blood Units, Patients, Shipments, and Supervisor main menu items. The title of the screen may vary (e.g., Select Unit for Printing, Unit ABO/Rh Confirmation, Unit Antigen Typing, Edit Unit Information) depending on the option selected. (The user then selects criteria to continue the search.)

User Action	VBECS
<p>1. Select an item from the main menu and the desired option.</p> <p>Click additional buttons as needed.</p>	<ul style="list-style-type: none"> Displays the Unit Search screen. <p>NOTES</p> <p>Although multiple search criteria can be checked, the Unit Search screen performs independent searches not a summary search for all the selected variables.</p> <p>Users can obtain a listing of all available unit specifics by creating a custom Blood Availability Report allowing them to make the appropriate selections without looking at each unit's record.</p>
<p>2. To search for a unit, enter the full or partial unit ID and click the find button, or just click the find button (Figure 80).</p> <p>Click a check box to select a unit.</p> <p>Click OK.</p>	<ul style="list-style-type: none"> Lists units that meet the search criteria. Displays unit-associated data. <p>NOTES</p> <p>To speed the search, enter a full or partial unit ID instead of just clicking the find button (which causes VBECS to search the entire database and takes longer).</p>

VBECS checks unit IDs and product codes after the user enters the product code. When there is a discrepancy, VBECS displays the Input Error icon and allows the user to edit the data.

Figure 80: Click the Find Button

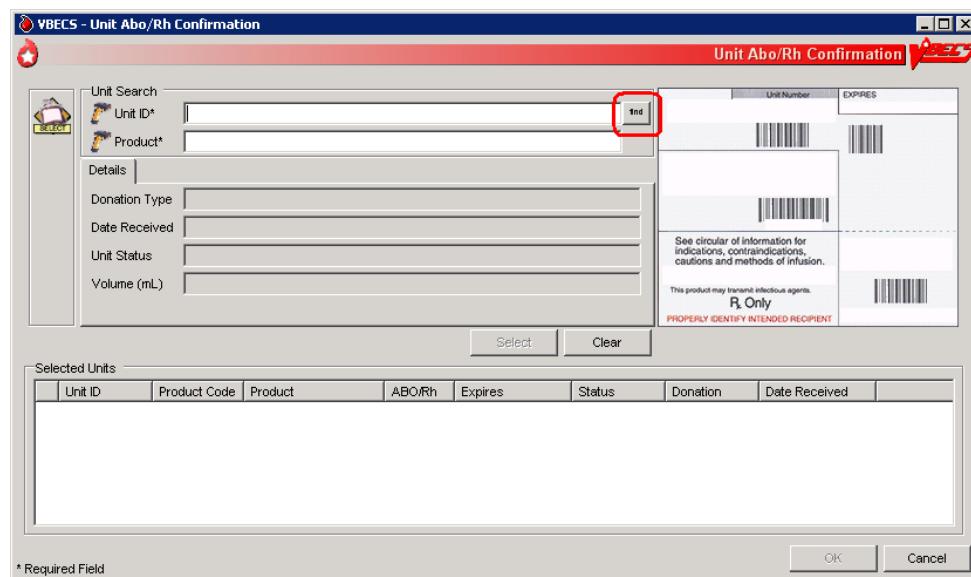


Figure 81: Display Unit Search Screen

VBECS - Select Units for Abo/Rh confirmation

Select Units for Abo/Rh confirmation

Search Criteria

Product Type

- Washed RED BLOOD CELLS
- RED BLOOD CELLS
- WHOLE BLOOD

ABO/Rh

- O Pos
- A Pos
- B Pos
- AB Pos
- O Neg
- A Neg
- B Neg

Special Testing

- CMV Negative
- Sickle Cell Negative

RBC Antigens

- D Neg
- D Pos

Received Date

- Before
- After

Expiration Date

- Before
- After

Units within Criteria : 7

Unit ID	Product	Code	ABO/Rh	Expires	Donation
<input type="checkbox"/> 1234567	RBC CPD>AS1	04210	O Neg	08/19/2011 23:59	Volunteer Allogeneic Donor
<input type="checkbox"/> 2345678	RBC CPD>AS1	04210	O Neg	08/19/2011 23:59	Volunteer Allogeneic Donor
<input type="checkbox"/> W038403000333	WB CPDA-1 PLT Red	E0065V00	B Pos	08/02/2011 12:30	Volunteer Allogeneic Donor
<input type="checkbox"/> W038412678911	RBC CPD>AS1 PA	E0338V00	A Pos	08/26/2011 12:31	Volunteer Allogeneic Donor
<input type="checkbox"/> W038412678912	RBC Wash None MNFI OPN	E0490D00	O Neg	07/28/2011 11:31	Volunteer Directed, Eligible for Cros
<input type="checkbox"/> W038412678913	RBC CPD>AS5 NMT	E0388V00	O Neg	08/26/2011 12:31	Volunteer Allogeneic Donor
<input type="checkbox"/> W040000000007	RBC CPD>AS1	E0291V00	O Neg	08/10/2011 11:43	Volunteer Allogeneic Donor

OK **Cancel**

Test Units

ABO/Rh Confirmation

Blood products that contain red blood cells (Red Blood Cells, Whole Blood, Granulocytes and Leukocytes) require ABO/Rh confirmation. The user may process ABO/Rh confirmation results by selecting an invoice or individual units of blood received in a shipment. VBECS records the ABO/Rh test results for the list of blood units tested.

Assumptions

- The division is “full service.”
- A unit’s labeled ABO/Rh was entered in the same division as the user during login of an incoming shipment.
- Only invoices processed within a user’s division may be accessed to create an ABO/Rh confirmation worklist.

Outcome

- Unit ABO/Rh observed test results and interpretations are available for selection for patient use.

Limitations and Restrictions

- This option is disabled for transfusion-only facilities.
- When building a worklist by invoice number, only invoiced units not previously tested and confirmed appear on the worklist. The user may deselect or add additional units.
- For batch processing, the user can select incoming shipment invoices from his division, processed within the previous five days.
- VBECS does not display product codes associated with unit IDs in ABO/Rh Confirmation. Therefore, do not create ABO/Rh confirmation worklists containing units with the same unit ID.
- POOLED GRANULOCYTES labeled with an Rh type of “Pooled” is considered an Rh Positive unit which does not require confirmation testing.
- POOLED GRANULOCYTES labeled with an ABO group of “Pooled” cannot be confirmed. A blood component unit received as “POOLED ABO” will have to be transfused using the Supervisory option, Document ABO Incompatible Transfusion.

Additional Information

- Within a multidivisional facility when ABO/Rh confirmation is repeated on units from a facility within the same database, both the originating facility ABO/Rh confirmation testing results and the current testing results are stored.
- The ABO/Rh confirmation worklist may be sorted by Unit ID, ABO/Rh, Product Type, Component Type, or Date/Time Received. When displaying the Date/Time Received, seconds are not displayed, but are used to sort the worklist when selecting Date/Time Received for the worklist sort criteria. VBECS uses seconds for calculations, for example in sorting items by time, even though it may not display seconds when time is shown.

User Roles with Access to This Option

All users

ABO/Rh Confirmation Testing

To confirm ABO/Rh testing, the user creates a worklist of individual units from an invoice, or creates a custom worklist. The user records observed test reactions and interpretations. VBECS uses truth tables to verify that user-entered interpretations are consistent with observed reactions. VBECS compares ABO and/or Rh interpretations to units labeled “ABO/Rh” during login and displays warnings when there are inconsistencies.



When a unit's ABO/Rh confirmation test is repeated after issue to a patient, VBECS does not compare the test results to the patient's blood type at time of issue. The user must detect discrepancies and notify clinical personnel, in accordance with local policy.

User Action	VBECS
<p>1. When processing a shipment, select Shipments from the main menu.</p> <p>When processing selected units, select Blood Units from the main menu.</p> <p>Select ABO/Rh Confirmation.</p> <p>Select an invoice or units.</p> <p>Click OK.</p>	<ul style="list-style-type: none">Displays options for processing blood product shipments. By default, displays units in the worklist in the order in which they were entered during Incoming Shipment.Lists invoices processed within the last five days in reverse chronological order.Displays fields for entering rack information.Allows units to be deleted from or added to the worklist. <p>NOTES</p> <p>The user may select only active units and those not in an inactive or final status.</p> <p>To create a customized list of units:</p> <ul style="list-style-type: none">Enter a unit ID or click the find button to search for a unit. VBECS lists units with “limited” status.Select parameters from the options on the Select Units for ABO/Rh Confirmation screen.Select Units within Criteria to work with.Click OK.Repeat this step until all units are included in the Selected Units list on the Unit ABO/Rh Confirmation screen.Click OK.
<p>2. Select a rack.</p> <p>Select a date or accept the default.</p> <p>Select a Tested By name.</p> <p>Click Using Automated Instrument, as appropriate.</p> <p>Click the Perform Anti-D Testing check box, when appropriate.</p> <p>Click OK.</p>	<ul style="list-style-type: none">Displays the racks available using the pick list provided.Displays the date and time testing was performed. The default is the current date and time. The date may be edited using the pick list provided.Displays the testing technologist’s identification. The default is the current user. Select a name from the pick list of valid division users to show a different testing technologist.Displays warnings and the option to indicate that the user is ready to begin entering results after all warnings are addressed.Generates an Exception Report entry when the user inactivates or deletes the partially completed result grids and enters a comment. <p>NOTES</p> <p>A user may disable the anti-D testing grid column in the ABO/Rh unit confirmation worksheet.</p>
<p>3. Select a unit for which to record test results.</p>	<ul style="list-style-type: none">Creates a grid for data entry of results and interpretations, displaying units in the order selected.Allows the user to select a blood test type and enter a test result for

User Action	VBECS												
<p>Select a cell and enter an observed test result from the legend at the bottom of the screen.</p> <p>Enter a valid interpretation, then press the Enter key to go to the next data grid interpretation or Comment cell.</p> <p>Select enter comments from the Unit Testing for the ABO and Rh interpretations, as needed.</p> <p>Repeat until all test results for the selected unit are entered.</p>	<p>the unit.</p> <ul style="list-style-type: none"> • Restricts data entry to valid selections for testing results. Displays testing results and interpretations entered. • An interpretation indicates that all appropriate test results were entered for the unit. <hr/> <p>NOTES</p> <p> VBECS emits an audible alert and generates a warning requiring override when the ABO and/or Rh confirmation interpretation does not match the ABO and/or Rh of the unit at initial login. When the user confirms the discrepancy, inventory units whose ABO/Rh confirmation does not match the ABO/Rh at login are marked "Quarantined." VBECS captures details for inclusion in an Exception Report (exception type: Unit ABO/Rh login vs. confirmation do not match, unit quarantined).</p> <p>ABO/Rh confirmation testing data grids accommodate uniform data entry for Forward ABO testing.</p> <p>The cursor moves with a single valid keystroke from the top left cell across the row, returning to the leftmost cell in the next row, and so on down the worksheet grid. The cursor stays in the reaction results cells by default (the technologist may use the mouse to click in the interpretation area) until the last entry is made in the reaction result portion of the grid. The cursor then goes to the first (top) interpretation cell.</p> <p>The user may enter free-text comments or comments from the Unit Testing for the ABO and Rh interpretations (Table 4).</p>												
	<p>Table 4: Valid Interpretations</p> <table border="1"> <thead> <tr> <th>ABO</th><th>Rh</th></tr> </thead> <tbody> <tr> <td>A</td><td>P, POS (positive)</td></tr> <tr> <td>B</td><td>N, NEG (negative)</td></tr> <tr> <td>AB</td><td>I, Inconclusive</td></tr> <tr> <td>O</td><td></td></tr> <tr> <td></td><td>I, Inconclusive</td></tr> </tbody> </table> <p>When a user saves a partially completed ABO/Rh testing worklist, VBECS clears unit entries with incomplete or no results, releases them from the worklist, and makes them available for inclusion in a future worklist.</p> <p>Row validation is established for the interpretation vs. reaction results pattern by system rules.</p> <p>ABO/Rh interpretations are consistent with observed test results based on system rules. Discrepancies are resolved before VBECS can verify ABO/Rh confirmation results and add them to the database. The user may reenter valid results or delete entries and start over.</p>	ABO	Rh	A	P, POS (positive)	B	N, NEG (negative)	AB	I, Inconclusive	O			I, Inconclusive
ABO	Rh												
A	P, POS (positive)												
B	N, NEG (negative)												
AB	I, Inconclusive												
O													
	I, Inconclusive												
4. Confirm that all results were reviewed and are acceptable.	<ul style="list-style-type: none"> • Repeats data validation and prompts the user to save. <hr/> <p>NOTES</p>												

User Action	VBECS
	<p>VBECS warns that it doesn't save partially completed unit ABO/Rh confirmation worksheets.</p> <p> When VBECS verifies that the ABO or Rh interpretation is "Inconclusive," it assigns the unit a "Quarantine" indicator, emits an audible alert, and captures details for inclusion in an Exception Report (exception type: ABO/Rh confirmation Inconclusive, unit quarantined).</p> <p>Upon successful completion of ABO/Rh confirmation, the unit status is updated. Successful completion means that the ABO and Rh (if required) interpretation is consistent with recorded observed results and matches the ABO/Rh at unit login and any repeat ABO/Rh tests on the unit that are not invalidated.</p>
5. Click OK to save the ABO/Rh confirmation testing.	<ul style="list-style-type: none"> Prompts the user to confirm that all results were reviewed and are acceptable.
6. Click OK .	

Unit Antigen Typing

The user performs antigen typing for one or more of the antigens in a blood unit and records the results in a worklist.

Assumptions

- The division is "full service."
- The units and the user are in the same division.
- Antiserum inventory is entered during Log In Reagents.
- Workload codes are assigned through Workload Codes.

Outcome

- Antigen typing was applied in the unit record for informational purposes in other options.
- VBECS stores the lot number, vial identifier, manufacturer, and expiration date of the positive and negative control cells with the testing record.

Limitations and Restrictions

- Antigen typing applies only to blood products that contain red blood cells. This option does not address platelet, histocompatibility locus antigens (HLAs), or Immunoglobulin A (IgA) antigenicity.
- VBECS does not check whether the lot and vial numbers selected for the positive control cells differ from those of the negative control cells.
- D antigen typing must not be combined in batch testing of Rh positive units.

Additional Information

- This option is disabled for transfusion-only facilities.
- Clinical significance and availability of typing reagents are not implied when antigens are included in the list of antigen type tests.
- A CPRS patient order is not required to perform antigen typing of units.
- A user may query the database for the presence of antigen-typed units and display or print the query results in Blood Availability.

- A user may retrieve evidence of antigen typing for any date or date range. Regulatory agencies require documentation of the reagent name, lot number, reagent expiration date, and positive and negative control.
- A user is required to complete antigen typing QC in order to save test results. If an antigen testing worklist is started but not saved and a second testing worklist is started for the same antigen, the QC grid will display and require completion because the first worklist has not been saved.

User Roles with Access to This Option

All users

Unit Antigen Typing

The user performs antigen typing for one or more of the antigens on blood units. VBECS creates a list of available antigen specificities, including weak D, that can be tested. For each antigen tested, the user must select a reagent from available inventory. VBECS will check to see if valid positive and negative controls were tested for the selected reagent on the calendar day of testing. When controls are not tested, VBECS adds tests for the positive and negative controls to the worksheet.

Red cell units prepared for patient transfusion may need to be negative for this antibody specificity based on their setting in the antibody table defined in Antibodies.

 *Deletion of a partially completed antigen typing worklist invalidates any completed tests contained on the worklist regardless of the status of the unit, including transfused units.*
It is the user's responsibility to verify that the records are being invalidated appropriately.

User Action	VBECS
<p>1. Select Blood Units from the main menu.</p> <p>Select Unit Antigen Typing.</p>	<ul style="list-style-type: none"> • Displays options for processing blood units. • Displays options for processing antisera.
<p>2. Click Create New Worklist.</p> <p>When partially completed worklists exist, click a check box to select an existing worklist and go to Step 5.</p> <p>Click Delete to delete the worklist from the list or OK to continue.</p>	<ul style="list-style-type: none"> • Accommodates selection of a partially completed testing worklist or creation of a new testing worklist. • Displays an option to select or delete a partially completed worksheet. • When selecting a partially completed worklist, the user verifies the lot numbers, manufacturer name, and expiration date of the antiserum and control cells previously selected. These may not be edited. When a discrepancy is found, the user must delete the worklist and start over. <p>NOTES —————</p> <p>VBECS displays the date and time the testing was performed (default: current date and time). These are editable to allow for retrospective entry.</p> <p>VBECS displays the testing technologist ID (default: current user). A pick list of valid division users is available for the user's selection.</p> <p>When an automated instrument was used for testing, the user enters the name of the instrument and may indicate that successful quality control testing was done for the date.</p>

User Action	VBECS
<p>3. To create a new testing worklist, select a date or accept the default.</p> <p>Select a “Tested By” name.</p> <p>Click a radio button to select Tube, Gel, or Solid Phase.</p> <p>Click Using Automated Instrument, as appropriate.</p> <p>Click OK.</p>	<ul style="list-style-type: none"> • Displays an option to select antiserum types. • Displays an option to select a lot number for each antiserum type selected. • Displays the lot number information. • Prompts the user to select the correct antiserum entry for the specificity selected. <p>NOTES</p> <p>VBECS lists available antigen type tests from which the user may select and allows the user to enter the specificity to expedite the search. The list of antigen type tests is derived from the Corresponding Antigens field from Appendix B: Table 13: Antibody and Antigen Table: Irregular Antibodies.</p> <p> VBECS displays data fields for the user to select or enter the lot number, and displays the manufacturer and expiration date of each antiserum to be used in this transaction to allow the user to choose the correct antiserum. VBECS displays in-date antisera and allows the selection of outdated antisera, if needed, emits an audible alert, and requires an override and generation of an Exception Report (exception type: expired antisera used). A comment is required.</p> <p> When a user attempts to change the phases of reactivity, VBECS emits an audible alert, warns that the user selected phases that were not used previously for this reagent, and asks whether the user wishes to continue.</p> <p>No returns the user to the screen to select the correct testing phases. Yes requires the user to enter a comment explaining why different testing phases are selected, and captures details for inclusion in an Exception Report (exception type: antigen testing phase change).</p> <p>When a user selects a reagent rack for which quality control (QC) testing was never performed in Enter Daily QC Results, VBECS displays a warning message, clears the selection, and allows the user to select a different rack.</p> <p>When the user QC'd the antigen typing reagent lot number and performed patient antigen typing on the same date under different testing phases, VBECS enables the QC grid for the antisera in the phases the user selected.</p>
<p>4. Define the antiserum used.</p> <p>Enter a lot number for each antiserum specificity.</p> <p>Select the testing phase:</p> <ul style="list-style-type: none"> • IS • IS/RT • IS/37 • AHG/CC. <p>Select the rack, if applicable.</p>	<ul style="list-style-type: none"> • Displays an option to select the testing phases for the antiserum specificity by lot number. • Displays a testing grid disabling the phases not to be entered. • Stores the antiserum information and checks whether the positive and negative controls need to be tested with this batch. • Prompts the user to enter the lot number and vial identifier for the positive and negative control used for each antiserum tested. <p>NOTES</p> <p>When VBECS displays lot numbers for selection, it checks the inventory system for the reagent name and lot numbers in a rack. The user cannot force the entry of a lot number that</p>

User Action	VBECS
Select Add . Select additional antisera, if desired. Respond to warnings and enter comments, when indicated. Select OK (Figure 82).	<p>VBECS disallowed on the selection list. The user must choose a different reagent. The user may not select a reagent with an “unsatisfactory” value in the Inspection field or with a quantity of zero.</p> <p> VBECS compares the expiration date of each reagent on the selected rack to the testing date, emits an audible alert, and warns when the reagent is expired.</p> <p>VBECS enables the rack identifier field when a user selects the weak D or AHG/CC phase and disables the field when a user selects other phases.</p> <p>VBECS automatically adds the positive and negative control tests to the selected reagent type worklist if they have not been tested on that calendar date for the selected antisera lot number. Reagent lot numbers that are selectable for positive and negative control cells:</p> <ul style="list-style-type: none"> • Reverse ABO Typing Cells (set) • A2 cell (vial) • Screening Cells (set) • Check Cells (vial) • Panel (set) • Other (set) <p>VBECS requests information about an antigen lot number only the first time it is used on a specific day.</p> <p>If there are multiple test tabs, VBECS displays a message that may not correspond to the selected tab, but will correspond to information required by one of the requested tests.</p>
5. Select one or more blood units to add to the worklist. Enter the unit ID and the product code from the inventory, or click the find button to search for a blood unit and select from the list of available units. Click OK when all units are selected.	<ul style="list-style-type: none"> • Allows the user to scan or enter units to be included in the worklist or select them from a list. • Lists the selected units.
6. Select the lot number for the reagent and the positive and negative control cells. Enter the vial identifiers, if required. Enter “0” when not using panel cells. Verify the reagent type and the positive and negative control cell lot numbers selected and vial identifiers entered, if required.	<ul style="list-style-type: none"> • Displays the reagent type and lot numbers, including the vial identifier. <p>NOTES —————</p> <p>The vial identifier, comprising one or two digits, identifies the vial contained in a kit of reagent red blood cells. A vial contains antigen positive or antigen negative blood cells to be used in quality control of the antisera.</p>
7. Respond to warnings and continue with antigen typing. Review the lot number, manufacturer name, and	<ul style="list-style-type: none"> • Displays each selected unit in a separate tab for each selected specificity with its own grids and quality control. • Displays a direct data entry (DDE) grid to enter serologic reaction results for antigen typing of units. This may already be filled with results from a previously saved partially completed worksheet.

User Action	VBECS
<p>expiration date for the antiserum and control cells. If they are not correct, start over.</p> <p>Select a control cell or unit to begin entering serologic reactions and interpretations in the data entry grid.</p>	<ul style="list-style-type: none"> Allows a user to select a unit or control cell and enter a test result. Allows only valid data entries in the grid reaction result and interpretation cells. Restricts data entry to valid selections for testing results. Keeps all testing results and entered interpretations on the screen until the user saves them. <p>NOTES</p> <p>VBECS displays the worklist for the grids by antiserum types.</p> <p>Valid observed test results entered by a user:</p> <ul style="list-style-type: none"> 1 [1+ (positive)] 2 [2+ (positive)] 3 [3+ (positive)] 4 [4+ (positive)] W [Weak (positive)] F [Mixed Field (positive)] M [Microscopic (positive)] R [rouleaux (negative)] 0 (zero) [No Agglutination (negative)] X (Not Tested) (does not indicate a disabled cell) <p>When the AHG phase is selected for testing, the AHG must be negative and the CC must be at least 1+ positive to consider the test negative for a unit or for the negative control cell.</p> <p>The column heading of the data grid specifies the phases of reactivity that may be used in a test selected for the specificity.</p> <p>The row heading of the data grid specifies the unit ID being tested. These fields cannot be edited.</p> <p>The rows for each antiserum type grid display the positive and negative control cells, then units, in the order in which a user added units to the worklist.</p> <p>A comment can be associated with any antigen test performed. A user may select canned comments required by an override from the drop-down list or select Other to enter free-text comments.</p>
8. Enter results for the control cells and for all selected units. Refer to the legend at the bottom of the screen (Figure 83).	<ul style="list-style-type: none"> Prompts the user to confirm that all results were reviewed and are acceptable. <p>NOTES</p> <p>The positive and/or negative control cell test results on the date of testing must be a valid positive or negative.</p> <p>When a weak D interpretation is positive or inconclusive and is verified for a tested blood unit that was labeled as (Rh) D negative, VBECS warns the user and adds the quarantine indicator to the unit.</p> <p>VBECS compares the antigen typing interpretation a user enters to previously assigned antigen typings for the same antigen (as assigned during Incoming Shipment, Edit Unit Information, or</p>

User Action	VBECS
	<p>previous testing). When the previous interpretation was not “inconclusive” and the current interpretation does not match the previously identified antigen type, VBECS warns the user and quarantines the unit.</p> <p>System rules establish row validation for the interpretation vs. reaction results pattern [for all user-indicated phase(s) of reactivity] for blood unit and quality control tests.</p>
9. Repeat Steps 4–8 for each antiserum type.	<ul style="list-style-type: none"> • Displays the option to save.
10. Respond to warnings and save.	<ul style="list-style-type: none"> • Saves the completed QC results with partially or fully completed unit testing data. • Reconfirms that all appropriate test results were entered, repeats all data validation, and saves the information on the worksheet displayed. <p>NOTES —————</p> <p>VBECS saves the testing of one antiserum specificity at a time.</p> <p>The Positive and Negative control cells must be complete and valid prior to the save of a unit antigen typing result. VBECS will display a tooltip message to the user if an attempt to save is made and the QC is not valid.</p> <p>VBECS stores the lot number, vial identifier, manufacturer, and expiration date of the positive and negative control cells with the testing record.</p>
11. Click OK to exit.	

Figure 82: Unit Antigen Typing

VBECS - Unit Antigen Typing

Unit Antigen Typing

Date Tested* 02/02/2006 14:00 Tested By* vhaishvbecs1

Testing Method* Tube Gel Solid Phase
Using Automated Instrument
Instrument Name*
Instrument QC'd?

Define Antiserum Specificity

Antigen Typing*
Lot Number*
Reagent Type
Manufacturer
Expires
Rack*

Testing Phases*
 IS
 IS/RT
 IS/37
 AHG/CC
 37/AHG/CC

Add Clear

Selected Specificities

	Typing	Reagent	Lot Number	Manufacturer	Expires	Phases	Rack
<input checked="" type="checkbox"/>	AGK	Anti-K	20	AAA	08/20/2006 23:59	AHG/CC	A

* Required Field

OK Cancel

Figure 83: Unit Antigen Typing Serologic Results

VBECS - Unit Antigen Typing

Unit Antigen Typing VBECS

AS TYPING

Unit Antigen Typing											
	Control	Reagent	Lot Number	Vial	IS	RT	37	AHG	CC	Interp.	Comment
✓	POS	Screening Cells	1	1				2	X	P	
✓	NEG	Screening Cells	1	2				0	2	N	

Unit ID	Product Name	IS	RT	37	AHG	CC	Interp.	Comment
✓ 1941954	RBC AS1, 04210				0	2	N	

Lot Number Manufacturer Expires

Key Valid Any Text

Valid Entries
Any Text

* Required Field OK Cancel

Modify Components

Modify Units

The user modifies a unit, which results in a different blood product (target product) with its own properties. A division's configuration determines the available modification types.

Assumptions

- Product Modifications was executed.
- The division is the active shipper for target blood products.
- A unit may be modified with or without patient assignment or restriction.
- A unit may be modified with or without a patient order.
- Units selected for modification must be in the same division as the user.
- Product labels and unit numbers for modified products are on hand.
- Supplies needed for the unit modification were entered in VBECS.

Outcome

- Data generated for an Exception Report are stored for inclusion in the next Exception Report.
- Workload-related data collected as the result of unit modification are stored for inclusion in Workload Reports.
- At the completion of a unit modification, VBECS assigns the original unit a final status of "modified."
- A permanent virtual link exists between the original unit and target units created.
- The shipper of the target unit created is the division in which the modification was performed. The cost assigned to the target unit is based on what was defined in Blood Products.

Limitations and Restrictions

- VBECS allows modification of units only in limited, assigned, crossmatched, and available statuses.
- VBECS presents the conservative expiration date for target units. VBECS does not know the specific collection details. It is the responsibility of each user to accurately enter and process any one blood unit based on its expiration date/time as labeled by the collection facility.

Additional Information

- VBECS requires the identification of the transferred-to container when using an open method or a sterile connection device (SCD).
- VBECS requires the identification of the SCD wafer data when using an SCD.
- When a user extends an expiration date beyond the default expiration date, the override message displays the original expiration date.
- Due to the variations in supplies and equipment used for other modification types, VBECS does not require supply or equipment entry, but allows entry of those data as part of the modification record.
- During the modification options (Modify Units: Pool Units; Modify Units: Split a Unit; Modify Units), a user may select an assigned blood unit for modification and transfer one patient assignment to the target units created.
- To reverse other modifications, remove the "modified" status from each unit through Remove Final Status.
- To edit the number of units in an existing pool, use Add/Remove Units from a Pool.

User Roles with Access to This Option

All users

Modify Units

The user selects the units to be modified. The user indicates the type of modification to be performed, then enters (preferably by scanning) the unit ID and product code.

VBECS displays a creation date and time (default: current date and time) for the target units, which the user may change to a date and time in the past and not before the login date and time of the original unit. VBECS calculates a new expiration date for the target unit (based on the creation date and time, modification type, method, and original unit product type) and records the date and time the original unit was modified. Some modification procedures require materials in addition to the new container. VBECS requires the entry of a new container lot number and SCD wafer, if indicated.

The modified target unit is now available for use.

User Action	VBECS
1. Select Blood Units from the main menu. Select Modify Units .	<ul style="list-style-type: none">• Displays options for processing blood units.• Displays valid modification options configured for the division and allows the user to select an option.
2. Enter the date and time of the unit modification. Enter the name of the user who modifies the unit.	<ul style="list-style-type: none">• Displays the default (current) date and time as the modification date and time.• Allows the user to change the date and time to a past date and time. <p>NOTES —————</p> <p>The user may not enter a future date and time.</p>
3. Click a radio button to select a modification type: <ul style="list-style-type: none">• Thaw• Pool (See Modify Units: Pool Units.)• Thaw/Pool Cryo (See Modify Units: Pool Units.)• Split/Divide (See Modify Units: Split a Unit.)• Irradiate• Leukoreduce• Volume Reduce• Wash• Rejuvenate• Freeze• Deglycerolize• Add/Remove Units from Pool (See Add/Remove Units from a Pool.)	<ul style="list-style-type: none">• Allows the user to enter data based on the modification type selected.• Allows the user to select a different processing technologist.• Displays only modifications enabled for the division.
4. Click a radio button to select a modification method, if required: <ul style="list-style-type: none">• Open• Closed• Sterile Connection Device (Figure 84)	<ul style="list-style-type: none">• Displays the selected method. <p>NOTES —————</p> <p>VBECS stores the date and time for the original unit when it is modified. The modification date and time are also stored as the creation time for the target unit.</p> <p>When the modification method is “Open” or “Sterile Connection</p>

User Action	VBECS
	<p>Device (SCD)," the user must select a container supply type. When the modification method is "Closed," VBECS does not allow the user to select a container supply type.</p> <p>For an SCD modification method, the user must indicate the completeness of the weld, in accordance with local policy and procedures. VBECS defines weld status as:</p> <ul style="list-style-type: none"> • Complete: the weld is complete and has no leakage or air bubbles. • Incomplete: the weld has a leak or air bubbles are present. VBECS changes the expiration date and time of the unit to those of an open system for that product type. VBECS stores this weld check with the historic record of the target unit created in the modification.
5. Click OK to accept the selected method, or edit it. 6. Scan, enter, or select the ID and product code of the unit to be modified.	<ul style="list-style-type: none"> • Allows the user to scan, enter, or select units to modify. • Displays current data for the selected unit. • Verifies the appropriateness of the selected unit based on system rules and truth tables and displays error and warning messages based on the modification type selected and the user's security level. • Prompts the user to continue with the selected unit or select a new unit based on the data in the review window. <p>NOTES</p> <p>VBECS allows the user to select a product type for the original unit based on the selected modification type.</p> <p>See Appendix B: Table 21: Allowable Product Modifications by Original Product Type.</p> <p>The user may select only active units and those not in an inactive or final status.</p> <p>When the user selects a status other than "limited," "available," "assigned," or "crossmatched" for an original unit, VBECS warns the user. There is no override.</p> <p>When the selected unit is not appropriate for the modification type selected, VBECS warns the user and instructs him to verify the selected unit and modification type. There is no override. When the unit is expired or quarantined, VBECS displays an error message.</p> <p>The user may select units containing red blood cells chosen for rejuvenation within three full days past their labeled expiration date. (For example, when a unit's expiration date is the current day, VBECS counts this day as "day 0," the next day as "day 1," etc. The last acceptable date and time for selection is day 3 at 23:59). After this three-day period, VBECS considers the unit to be expired.</p> <p>When the modification type is "volume reduce," VBECS records the discarded plasma volume on the unit record. Otherwise, the user may not edit the Discarded Plasma Volume field. The discarded plasma volume may not be "0" (zero).</p>

User Action	VBECS
7. Click OK to accept the unit for modification.	<p> See Table 5 for alerts that may occur during this option.</p> <ul style="list-style-type: none"> When the modification type is “thaw” or “irradiate,” allows the user to select additional units to process or to continue processing the selected units. Otherwise, go to Step 8. <p>NOTES</p> <p>The user may process multiple products in a batch when the units are thawed or irradiated. The creation date and time apply to all units in the batch.</p>
8. Click the Weld Complete or Weld Incomplete radio button to indicate the integrity of the SCD weld, when indicated. Click the Container, Supplies, Equipment, or SCD Wafer tab and enter another supply and/or equipment lot or ID number (Figure 85).	<ul style="list-style-type: none"> Calculates and displays a default expiration date and time of the target unit based on rules specific to the original blood product, the modification type, and processing conditions (open, closed); prompts the user to accept the date and time; and allows the user to edit them. <p>NOTES</p> <p>VBECS maintains lot and ID numbers entered, the new container lot number in the unit modification record, and the historic record of the target unit.</p> <p>VBECS checks the inventory system for the supply name, lot numbers, expiration date, and a quantity greater than “0” (zero). VBECS does not allow the user to select supplies with an “unsatisfactory” value in the Inspection field. When the user enters an unsatisfactory supply lot number, VBECS warns the user. The user must select a different supply item.</p> <p> See Table 5 for alerts that may occur during this option.</p>
9. Select a target product type from the drop-down list in the Product field.	<ul style="list-style-type: none"> Displays target component types for the modification, conditions selected, and the original unit product type, and allows the user to select a single product type. Verifies that an active target component was selected and warns the user. Allows the user to enter lot and ID numbers (with auto-completion) or select multiple items from the supply and equipment lists. Warns that the target unit is not active and allows the user to activate the target component. <p>NOTES</p> <p>The user must select a modification method (open, closed, Sterile Connection Device) for modification types, as described in Appendix B: Table 21: Allowable Product Modifications by Original Product Type.</p> <p>VBECS displays available target component types based on ISBT 128 modification rules and site parameters for the division.</p> <p>When the user selects a target component type not marked as active, VBECS warns the user, allows him to mark the target component as active, and returns him to this step in the modification process. Otherwise, VBECS returns the user to this step and allows him to select another target component type.</p> <p>Site parameters indicate which products in the product code</p>

User Action	VBECS
	<p>database are marked as active for a division. The user may define fields for a product type at a site:</p> <ul style="list-style-type: none"> Division shipper (automatically set to the user's division): <ul style="list-style-type: none"> Cost (U.S. currency) Return Credit Percentage Active? <p>These are the only fields in the blood product table that are editable.</p> <p>When the original unit type did not require ABO/Rh confirmation, the target unit, when created, does not either. The target unit is immediately available for assignment, crossmatch, issue, or transfusion on completion of the modification.</p> <p>When ABO/Rh confirmation testing was not performed on the original unit when the target unit was created, VBECS assigns "limited" status to the target unit and saves it for inclusion in an ABO/Rh confirmation worksheet.</p> <p>When the original unit is assigned or crossmatched to a patient at the time of modification, VBECS notifies the user and asks whether he wishes to automatically assign the modified units to the patient. Yes (default) applies the unit's assignment or crossmatch data to the target component. No does not apply the assignment to the target, and changes the target unit status to "available." The original unit retains patient assignment in both cases.</p> <p> See Table 5 for alerts that may occur during this option.</p>
10. Accept or edit the expiration date.	<ul style="list-style-type: none"> Verifies the validity of the expiration date and warns the user, if indicated. <p>NOTES —————</p> <p>When VBECS calculates the expiration date of a unit in days, starting from the date and time of modification, the expiration time is 23:59. When VBECS calculates the expiration date in hours, the unit expires at the specified hour and minute.</p> <p>VBECS calculates an expiration date for the target component. For all modification types except Rejuvenate and Freeze, VBECS compares the calculated expiration date with that of the original unit and displays the earlier of the two as the default. The user may edit the calculated expiration date. (Due to the nature of the process, the calculated expiration date will greatly exceed the original expiration date of the original unit for Rejuvenate and Freeze, and will supersede the unit's original expiration date.)</p>
11. Review the displayed unit data and click OK to confirm the modification. If the unit has been relabeled, click Yes to verify the label (Figure 86).	<ul style="list-style-type: none"> When the original unit is currently assigned to a patient, notifies the user and asks whether he wishes to transfer the status to the target unit. Yes updates the target unit status.

User Action	VBECS
<p>Click No to cancel and continue at Step 14.</p> <p>12. Scan the label barcodes to verify the target labels, based on local practice and labeling type of the modified unit.</p> <p>Enter a comment to complete the verification. (This step is optional.)</p>	<ul style="list-style-type: none"> • Displays the label verification window. • Indicates whether the scanned barcodes match those in the label. • For each target created, system prompts user to indicate readiness to verify the labeling of the modified product. • No new labels applied to modified product. Proceed to save the unit. • Displays a screen to accommodate data entry of the modified unit's: <ul style="list-style-type: none"> ◦ ABO/Rh ◦ Product Code ◦ Expiration Date • This step repeats until all targets have been labeled and/or are ready for continued processing. <p>NOTES</p> <p>The user may select "Verify modified unit label" to verify the modified unit label. When the user does not verify the label, he must enter a comment and then click Cancel to exit.</p> <p>When the user does not select "Verify modified unit label," he may continue at Step 14. VBECS bypasses the label verification and does not create an Exception Report.</p> <p>When a component is modified and a new label is applied, American Association of Blood Banks (AABB) standards require the labeling process to include a method to ensure the correctness of the unit number, ABO/Rh, expiration date, and component label.</p> <p>VBECS does not display the modified unit data. VBECS instructs the user to verify the labeling of the target product and notifies him that scanning the barcode is preferred. When no barcode is available, the user may enter the data.</p> <p>To verify the label, VBECS compares the entered data with the data of the target unit from the database. The user may scan or enter:</p> <ul style="list-style-type: none"> • Unit ABO/Rh Label • Unit Product Code • Unit Expiration Date
13. Scan or enter each label on the target unit.	<ul style="list-style-type: none"> • Compares the data entered in each field to the database. • Asks whether the user wishes to save the data. When the user confirms that the data are correct, notifies him that the label verification is complete. • Records the modified or restriction-released unit relabeling data as part of the unit record and stores the data for inclusion in the unit history report. • Saves only the technologist name and the date and time the label was verified. When the user indicates that a new label applied fails to complete the verification process, releases the assignment, quarantines the unit, and captures details for inclusion in an Exception Report (exception type: target product label failed checks). • When the user confirms that the comparisons are correct, notifies him and asks whether he wishes to continue. <p>NOTES</p>

User Action	VBECS
	<p> When there is no match, VBECS emits an audible alert, notifies the user, and instructs him to verify the labeling and reenter label data or select Cancel.</p> <p>When there is no previous entry for a field, VBECS notifies the user and instructs him to enter a value or select Cancel.</p> <p>VBECS requires a comment when the user selects Cancel, and captures details for inclusion in an Exception Report (exception type: target product label failed checks).</p>
14. Click OK to confirm the save.	<ul style="list-style-type: none"> • Saves the new modified unit data. • Updates the status of the original unit to “modified.” • Records the unit ID of the target unit on the record of the original unit. • Prompts the user to modify another unit. <p>NOTES —————</p> <p>The target unit automatically inherits original unit data, which the user may not edit:</p> <ul style="list-style-type: none"> • Unit ABO/Rh (from login) • ABO/Rh confirmation testing information (when performed) • Unit ID • Unit Volume (exception: when a unit is volume reduced, the target unit does not inherit the original unit's volume. The user must enter the target unit's volume.) • Donation type • Disease marker testing information, when autologous. • Restricted for patient name, ID. • Current status • Biohazard Indicator • Quarantine Indicator • Division • Additional unit testing (CMV, Sickle Cell, antigen typings)
15. Repeat Steps 3–14 to modify another unit or exit.	

Figure 84: Modify Unit(s)



Figure 85: Split/Divide Unit

VBECS - Split/Divide Unit

Split/Divide Unit VBECS

Split*

Adult Split

Pediatric Split mL

Modification

Method: Sterile Connection Device

Modification Date: 02/02/2006 14:10

Weld*

Weld Complete Weld Incomplete

Original Unit



Unit ID: 6832766
Product: 04210 - RBC AS1
ABO/Rh: O Neg Volume (mL): 350
Expiration Date: 02/28/2006 23:59
Unit Status: Available

Donation Type: For Autologous Use Only
Patient Name: VBPATIENT, EIGHTYONE
Patient ID: 000-00-0081
Disease Marker Testing: Unit fully tested and neg
Restricted For: Assigned To: Special Testing

Target Unit 1 | **Target Unit 2**

Unit ID: 6832766A
Product: 04210 - RBC AS1
ABO/Rh: O Neg Volume (mL)*: 175
Expiration Date*: 02/28/2006 23:59
Unit Status: Available

Donation Type: For Autologous Use Only
Patient Name: VBPATIENT, EIGHTYONE
Patient ID: 000-00-0081
Disease Marker Testing: Unit fully tested and neg
Restricted For: Assigned To: Special Testing

Container | **Supplies** | **Equipment** | **SCD Wafer**

Lot Number	Supply Type	Manufacturer	Expires
553466	Clean Containers	Super Supplier	03/02/2006

Original Container

*** Required Field**

Figure 86: Verify Printed Label

The screenshot shows the VBECS software interface for verifying a printed label. On the left, there's a barcode scanner icon and a comment field for skipping verification. The main area contains input fields for Unit ID (6832766A), Product Code (04210), ABO/Rh (O Neg), and Expiration Date (02/28/2006 23:59). To the right, a large window displays a printed label with the same information. The label also includes a circular logo for 'Rh NEGATIVE' and a barcode for 'LifeSource - Glenview'. A note at the bottom of the label reads 'PROPERLY IDENTIFY INTENDED RECIPIENT'.

Alerts

Table 5: Alerts That May Occur in Modify Units

Step	Alerts
6	When the unit modification date and time entered are earlier than the date and time the original unit was received in inventory, VBECS warns the user and allows him to reenter the modification date and time.
6	VBECS compares the original unit's expiration date and time to those of the modification to determine whether the selected unit was expired at the time of modification. When the user is a Blood Bank Technologist and the unit selected for modification expired, VBECS emits an audible alert, displays an error message, clears the selected unit, and allows the user to select another unit for modification. ►► When the user is an Enhanced Technologist and the unit selected for modification is expired, VBECS emits an audible alert, notifies the user, and asks him to verify to continue, or cancel to correct the entry data. Override requires a comment and VBECS captures details for inclusion in an Exception Report (exception type: expires unit modified). Cancel clears the selected unit and allows the user to select another unit for modification.
6	When ABO/Rh confirmation testing is required but not performed, VBECS notifies the user that testing must be completed before issuing a target unit. The user may continue with the modification.
6	When a unit has an indicator of "quarantined," VBECS notifies the user that it must be released before processing may continue, clears the selected unit, and allows the user to select another unit. ►► When the user is an Enhanced Technologist and a unit has an indicator of "quarantined," VBECS notifies the user that it must be released before processing may continue and asks whether the user wishes to release the unit. Cancel clears the selected unit and allows the user to select another unit. Override releases the unit through Discard or Quarantine.
6	When the selected unit has a positive biohazard Indicator, VBECS emits an audible alert, warns the user, and asks whether the user wishes to continue with the modification. Cancel clears the unit and the user may select another unit. Override requires a comment and VBECS captures details for inclusion in an Exception Report (exception type: biohazardous unit modified).

Step	Alerts
8	When a supply item is expired, VBECS emits an audible alert, warns the user, and asks whether the user still wishes to use the item. Yes requires a comment and VBECS captures details for inclusion in an Exception Report (exception type: expires supply used in unit modification). No requires the user to select another supply lot number.
10	When the user extends the system-recommended expiration date and time, VBECS warns the user, asks him to reenter the unit expiration date, and redisplays the original entry. The user may enter another expiration date and time for the target unit. ►► When an Enhanced Technologist extends the target unit expiration date and time past the system-recommended expiration date and time, VBECS emits an audible alert, warns the user, and asks whether he wishes to continue. Cancel restores the expiration date. Override allows the user to continue with this date, requires a comment, and captures details for inclusion in an Exception Report (exception type: target product outdate extended).

Modify Units: Split a Unit

The user selects the unit to be split (divided) for transfusion.

Assumptions

- See Modify Units.

Outcome

- See Modify Units.

Limitations and Restrictions

- A unit previously split in VBECS cannot be split again. A unit pooled in VBECS may not be split.
- This option appends alphabetical characters to the Codabar unit ID of split units created from original units. This prevents scanning the unit ID barcode with the split indicator.
- By ISBT 128 convention, split units created from an ISBT 128 original unit will have identical unit ID and five-digit product codes.

Additional Information

- The user may select blood units from the product types indicated in Appendix B: Table 21: Allowable Product Modifications by Original Product Type for splitting.
- To reverse other modifications, remove the “modified” status from each unit through Remove Final Status.

User Roles with Access to This Option

All users

Modify Units: Split a Unit

The user selects the unit to be split for transfusion. VBECS automatically applies directed and autologous restrictions to aliquots created. The user may or may not maintain patient assignments (usually, these units are crossmatched and available for the patient before the split is made). This section lists only rules unique to the division of a blood unit.

User Action	VBECS
1. Select Blood Units from the main menu.	<ul style="list-style-type: none"> • Displays options for processing blood units. • Displays valid modification types for selection.

User Action	VBECS
Select Modify Units . 2. Click the Split/Divide radio button. Click OK .	<ul style="list-style-type: none"> Allows the user to search for a unit.
3. Select a unit and click OK . Click Select to confirm the unit selected.	<ul style="list-style-type: none"> Displays data for the original and target units.
4. Edit the number of adult splits to be created, if required (Figure 85). Click Split* to verify the number of splits to be created and add the appropriate number of target unit tabs.	<ul style="list-style-type: none"> Displays the number of splits to be created (default: 2) for verification. Requires entry of the number of splits to be created. <p>NOTES</p> <p>When Adult Split is selected, the user must indicate the number of split units to be created. The user may split a unit into 2 to 26 subunits until the volume of the original container is consumed.</p>
5. Select or confirm the desired target component type. For a Sterile Connection Device modification method, click the Weld Complete or Weld Incomplete radio button to indicate the condition of the weld. Click the Container , Supplies , Equipment , or SCD Wafer tab and enter another supply and/or equipment lot or ID number.	<ul style="list-style-type: none"> Allows the user to enter lot and ID numbers or select multiple items from the supply and equipment lists. <p>NOTES</p> <p>When multiple target units are created, the user must associate a container lot number with each unit created. "Original container" may serve as the lot number for one of the aliquots created.</p> <p>The original blood collection facility must keep a record of the original container lot numbers so that, when a lot number is recalled, it is known which patient was exposed to it.</p> <p>The expiration date and time are the same as the original unit for each split unit created using an SCD or a closed method. This applies to all types of splits.</p>
6. Accept or edit the displayed expiration date.	<ul style="list-style-type: none"> Displays the expiration date. Calculates and displays the volume of each split unit. <p>NOTES</p> <p>The expiration date and time for each split created are 24 hours from the date and time of modification (using an open method).</p> <p>The user may edit the default volumes, which must not equal "0" (zero) or be greater than 2,000 mL. VBECS divides the volume of the original unit evenly among the splits and adds the remainder to the first split.</p>
7. Accept or edit the displayed unit volumes.	<ul style="list-style-type: none"> Displays the original and target unit data for review. <p>NOTES</p> <p>The user may not edit original unit data.</p> <p>The product code of target units created is the same as that of the original unit.</p> <p>Target units inherit the original unit data, except for volume: volume is assigned based on the conditions of modification.</p>

User Action	VBECS
8. Review the displayed data for the target units and confirm the modification type.	<ul style="list-style-type: none"> Displays the target unit's ID number. <p>NOTES</p> <p>VBECS assigns the unit ID of the target units based on the labeling of the original unit:</p> <ul style="list-style-type: none"> Codabar: the ID number comprises the original unit ID number and a suffix (A–Z) for the 2 to 26 target units, beginning with A. The five-digit product code remains the same. ISBT 128: the unit ID, five-digit product description code, and sixth-character donation code remain the same. Characters 7 and 8 of the ISBT 128 product code barcode reflect the aliquot level: when a previously undivided unit is split, character 7 of the unique eight-digit ISBT 128 product code is incremented alphabetically in uppercase letters, resulting in characters 7 and 8: A0, B0, C0, D0 ... Z0.
9. Click OK to save. Click Yes to verify the label.	<ul style="list-style-type: none"> Prompts the user to save the data, including patient assignment and restriction information, if any. Saves aliquots created from the modified unit. Records the unit ID of the target unit on the record of the original unit. Saves only the technologist name and the date and time the label was verified. When the user indicates that a new label applied fails to complete the verification process, releases the assignment, quarantines the unit, and captures details for inclusion in an Exception Report (exception type: target product label failed checks).
10. Click Yes to save. Repeat Steps 2–9 to split additional units, or return to Step 1 to use a different modification method on another unit.	<ul style="list-style-type: none"> Prompts the user to split or modify another unit.
11. Exit.	

Modify Units: Pool Units

The user selects units to be pooled.

Assumptions

- See Modify Units.

Outcome

- See Modify Units.

Limitations and Restrictions

- VBECS allows modification of units only in limited, assigned, crossmatched, and available statuses.
- When a user adds a biohazardous unit to a pool, VBECS does not designate the pool as biohazardous on the Blood Transfusion Record Form (BTRF) or Caution Tag.

Additional Information

- When autologous or directed units are pooled, the “restricted for” patient must be the same for all units. VBECS assigns the same donation type and applies the “restricted for” patient data to the pool.
- To reverse the pooling of units:
 - Inactivate the target through Edit Unit Information.
 - Remove the “modified” status from each unit in the pool through Remove Final Status.
- To edit the number of units in an existing pool, use Add/Remove Units from a Pool.

User Roles with Access to This Option

All users

Modify Units: Pool Units

When the pooling is in process or almost complete, the user records it in VBECS. Data entry may also be completed retrospectively. The user pools units and enters the ID and product code of the original unit and the IDs and product codes for each unit added to the pool.

VBECS records the date and time the units are pooled, uses the date and time to calculate a new expiration date for the pool, and assigns a unit ID and ABO/Rh, as defined by system rules. When an original unit in a pool was previously assigned to a patient, VBECS asks the user whether the entire pool is to be assigned to that patient. VBECS automatically assigns autologous and directed patient restriction data from an individual unit to the pool. VBECS determines the volume of the pool based on the sum of the volumes of the units in the pool.

The new, pooled unit needs a full set of unit properties, some of which were defined by the original units' properties. VBECS assigns new properties such as volume, date and time created, container lot number, and expiration date to the target.

VBECS uses system rules to check pertinent blood bank requirements associated with pooling and displays appropriate target unit information for review. The user reviews and edits the information, as needed, and completes the transaction. VBECS updates internal records, including the collection of workload and unit statuses, for each unit involved.

The user may now assign the units or print the BTRF and Caution Tag, and notify appropriate personnel that the ordered units are ready for issue.

User Action	VBECS
1. Select Blood Units from the main menu. Select Modify Units .	<ul style="list-style-type: none"> • Displays options for processing blood units. • Displays valid modification types for selection. <p>NOTES —————</p> <p>The >> Enhanced Technologist or above is authorized to modify expired blood products.</p>
2. Click the Pool or Thaw/Pool Cryo radio button. Click OK .	<ul style="list-style-type: none"> • Allows the user to search for a unit.
3. Scan, enter, or select the ID and product code of the most restricted (e.g., biohazardous, autologous, directed) unit to be pooled.	<ul style="list-style-type: none"> • Displays current data for the selected unit. <p>NOTES —————</p> <p>VBECS applies the restriction of the first unit selected to the pool.</p> <p> Based on the product code, VBECS determines the component class of each selected unit to verify that the component class of the selected unit matches that of the first unit selected for the pool (only units from the same component class may be pooled). When they do not match, VBECS warns the user and clears the entry. The user may enter another unit.</p> <p>VBECS automatically “thaws” frozen components when the user selects Thaw/Pool.</p>
4. Enter the ID and product code of another unit of the same donation type for inclusion in the pool.	<ul style="list-style-type: none"> • May warn the user and asks whether to release the unit from “restricted for patient.” <p>NOTES —————</p> <p>When the donation type of a unit selected for inclusion in a pool is restricted, the user may release the directed unit to the general blood supply, if indicated, and complete the pooling process.</p> <p> When a user selects an autologous unit for inclusion in a pool originally intended to be allogeneic or directed based on previously selected units, VBECS warns the user and clears the screen. The user may select another unit.</p> <p>When the first unit entered in a pool is labeled in Codabar, VBECS determines the unit ID of the pool created from Codabar units using the five-character station identifier, two-number year, and sequential number of the pool event at the division, and displays it for acceptance or edit. The user may enter an alphanumeric combination of six to 13 characters, for example, a station number of 55555, year 2002, and pool number 0520, so the VBECS-determined unit ID number is 5555502520.</p> <p>When the first unit entered for a pool is labeled in ISBT 128, VBECS determines the unit ID of the pool created from ISBT-labeled units using the country and ICCBBA-assigned facility code (division preparing the pool), two-number year, and six-number sequential number of the pool event at the division, and displays it for acceptance only.</p>

User Action	VBECS
5. Click Yes to release the unit from its “restricted for patient” status.	<ul style="list-style-type: none"> Allows the user to release the directed unit and move it to inventory. <p>NOTES</p> <p>When needed, VBECS allows the user to select Free Directed Unit For Crossover.</p>
6. Accept the previously selected unit and repeat Step 4 to select another unit to add to the pool.	<ul style="list-style-type: none"> Displays all data for all original units. <p>NOTES</p> <p>VBECS tallies the number of units pooled. The user may not edit this number. VBECS calculates the pooled unit's cost and passes that amount to the patient billing system in U.S. currency.</p> <p>VBECS does not allow a unit to be entered twice in a pool. A split unit created in VBECS may not be pooled. A unit that has a division code (7th/8th characters in the ISBT Product code) present upon receipt in incoming shipment may be pooled.</p>
7. Accept the displayed units' data for the pool.	<ul style="list-style-type: none"> Displays target component types for the modification, conditions selected, and the original unit product type, and allows the user to select a single modification type.
8. Select or confirm the desired target component type. Click the Container, Supplies, Equipment , or SCD Wafer tab and enter another supply and/or equipment lot or ID number.	<ul style="list-style-type: none"> Verifies that an active target component was selected. Requires activation, as needed. Allows the user to enter lot and ID numbers (with auto-completion) or select multiple items from the supply and equipment lists.
9. Accept the displayed expiration date or edit it.	<ul style="list-style-type: none"> Calculates and displays the expiration date of the pooled unit based on rules specific to the participant blood product. Verifies the validity of the expiration date and warns the user, when indicated. <p>NOTES</p> <p>When using the “Thaw/Pool” option to prepare a pool from frozen components, the units are simultaneously thawed and pooled..</p> <p>VBECS calculates an expiration date for the pool based on the product type, the date and time created, and the expiration dates and times of each pool participant. When one or more participants in the pool have an expiration date that is earlier than the calculated expiration date, VBECS uses the earlier date as the default expiration date for the pool.</p> <p> For each unit in the pool that was previously assigned or crossmatched to a patient, VBECS warns the user and asks whether he wishes to assign the entire pool to that patient. Yes applies the patient-assigned data to the pool. No releases the unit from the patient.</p>
10. Click OK to save.	<ul style="list-style-type: none"> Calculates the volume, total number of participants, ID, and ABO/Rh for the pool. Applies CMV negative status to the pool, when indicated. Determines whether any pool participants were previously assigned to a patient and queries the user to maintain the association, when indicated. <p>NOTES</p>

User Action	VBECS
	<p>When the ABO/Rh of each unit in a pool is identical, VBECS assigns that ABO/Rh to the pool. When the ABO and/or Rh of the units in a pool are not identical, VBECS assigns "mixed" as the ABO and/or Rh group and type.</p> <p>VBECS determines the volume of the pooled unit based on the sum of the individual volumes of the original units included in the pool and displays that as the default. The user may edit this volume (1–2,000 mL).</p> <p>VBECS labels the pool as CMV negative, Irradiated, or Leukocyte reduced only when all participant units are CMV negative, Irradiated, or Leukocyte reduced. There is no override.</p>
11. Click Yes to verify the label. Click OK to save.	<ul style="list-style-type: none"> • Saves the new pooled unit. • Updates the status of each participant unit to "modified." • Records the unit ID of the target unit on the record of the original unit. • Prompts the user to modify another unit. • Saves only the technologist name and the date and time the label was verified. When the user indicates that a new label applied fails to complete the verification process, releases the assignment, quarantines the unit, and captures details for inclusion in an Exception Report (exception type: target product label failed checks).
12. Repeat Steps 2–11 to modify another unit, or exit.	

Add and/or Remove Units from a Pool

The user selects a pooled unit to add and/or remove some of the participants.

Assumptions

- See Modify Units.

Outcome

- See Modify Units.

Limitations and Restrictions

- See Modify Units.

Additional Information

- Units must be in “available” status for data to be inactivated regardless of historic patient associations.
- To edit the number of units in an existing pool, use Add/Remove Units from a Pool.

User Roles with Access to This Option

- Enhanced Technologist (add/remove prior to patient associations)
➤➤➤➤ Traditional Supervisor (add/remove with patient associations)

Add/Remove Units from a Pool

The user adds or deletes units to or from a pool.

User Action	VBECS
1. Select Blood Units from the main menu. Select Modify Units .	<ul style="list-style-type: none">• Displays options for processing blood units.• Displays valid modification types and allows the user to select an option.
2. Click the Add/Remove Units from a Pool radio button. Click OK .	<ul style="list-style-type: none">• Allows the user to search for a pooled unit.
3. Scan, enter, or select the ID and product code for a pooled unit.	<ul style="list-style-type: none">• Displays current data for the selected unit. <p>NOTES —————</p> <p>The user must select a unit that was previously pooled within the user's division.</p> <p>The user may edit a pooled unit only once.</p> <p> VBECS does not list units (for selection) with a final status. When a user enters a unit with a final status, VBECS alerts the user and clears the screen. The user may select another unit.</p>
4. Review the existing pooled unit data. Select one or more units to add to the pool.	<ul style="list-style-type: none">• Allows the user to select one or more units to add to or remove from the pool. <p>NOTES —————</p> <p>The user may not remove all of the original units.</p>

User Action	VBECS
Clear one or more check boxes to remove units from the pool.	The user may add only units with the same product and donation types as those of the previously pooled unit.
5. Confirm to add the selected units to or remove them from the pool (Figure 87).	<ul style="list-style-type: none"> • Displays the updated pooled unit data. • Requires the user to confirm to continue. • Requires a comment. <p>NOTES</p> <p>All units in a pool must be of the same restriction type (allogeneic, autologous, directed). The specific donation type within the restriction type of the first unit selected for the pool is assigned to the target. For each individual unit removed from a pool, VBECS compares the unit's original expiration date and time with the current date and time to determine the status to be assigned to units that are no longer part of the pool.</p> <p>When the unit's original expiration date and time is in the future, VBECS assigns units removed from the pool a status of "available." When the unit's original expiration date and time is in the past, VBECS assigns units removed from the pool a status of "available."</p> <p> When the user adds a unit to a pool and the unit expires before the pool expires, VBECS emits an audible alert, warns the user, and asks the user to continue or cancel. Override requires a comment and VBECS captures details for inclusion in an Exception Report (exception type: expired unit modified). Cancel clears the selected unit data and allows the user to select another unit to modify.</p>
6. Review the entered data and save the updated pool record.	<ul style="list-style-type: none"> • Saves the modified pooled unit. Updates the status of each added unit to "modified." Updates the status of each removed unit. Records the unit ID of the target unit on the record of the original unit.
7. Repeat Steps 1–6 to modify another unit or exit.	<ul style="list-style-type: none"> • Allows the user to modify another unit.

Figure 87: Add/Remove Units from Pool

VBECS - Add/Remove Units from Pool

Add/Remove Unit(s) from Pool

Pooled Unit

Unit ID	589A6061	Expiration Date	02/02/2006 18:36	ABO/Rh	O Neg
Product	12091 - PLT Pool	Modification Date	02/02/2006 14:36	Modified By	vhaishvbecs1

Unit Search

<input type="checkbox"/> Unit ID*	<input type="text"/>	1nd
<input type="checkbox"/> Product*	<input type="text"/>	

Details

Donation Type	<input type="text"/>
Date Received	<input type="text"/>
Unit Status	<input type="text"/>
Volume (mL)	<input type="text"/>

Select **Clear**

Selected Units

	Unit ID	Product Code	Product	ABO/Rh	Expires	Status	Donation	Date Received
<input checked="" type="checkbox"/>	4654681	12000	PLT	O Neg	02/04/2006 23:59	Modified	Volunteer Allogeneic Donor	01/31/2006 10:05
<input checked="" type="checkbox"/>	34631818	12000	PLT	O Neg	02/04/2006 23:59	Modified	Volunteer Allogeneic Donor	01/31/2006 10:05

* Required Field

OK **Cancel**

Maintain Unit Records

Edit Unit Information

The user may inactivate unit tests, correct a data entry error from login, and add information about the unit that was not recorded during login.

Assumptions

- The user is logged into the division that contains the unit record.
- The connection to VistA is active.

Outcome

- The unit's new or changed information was saved with its record.

Limitations and Restrictions

- The user may correct a limited amount of unit information.
- VBECS does not limit the discarded plasma volume.
- If errors are found in the unit's product code, division or blood type at login, supplier, or donation type, the user must inactivate the unit and reenter it correctly.
- When a unit does not have a patient association, editable fields are available to an ➤➤ Enhanced Technologist.
- ABO/Rh confirmation and antigen type testing cannot be inactivated when the unit is in a final unit status.
- A user can change the restricted to patient associated with a split unit when it is not assigned to the patient. **Additional Information**
- Units must be in "available" status to have information inactivated regardless of historic patient associations.
- VBECS does not add corrected unit tests to the Pending Task List (PTL): the user simply repeats the unit test.
- Crossmatches are considered a patient test and are inactivated in Invalidate Test Results.

User Roles with Access to This Option

- Enhanced Technologist

Edit Unit Information

The user records additional unit information normally recorded during incoming shipment (login), i.e., CMV negative, sickle cell negative, and antigen typing provided by the supplier.

To inactivate a unit, remove the final status, return the unit from issue, and release it from previous assignment. No changes may be made to units in a final status.

User Action	VBECS
1. Select Blood Units from the main menu. Select Edit Unit Information .	<ul style="list-style-type: none">• Displays options for processing blood units.• Displays options for selecting and editing a unit.
2. Scan or enter a unit ID and product code. May cancel and	<ul style="list-style-type: none">• Allows the user to scan or enter a unit ID and product code or select a

User Action	VBECS
<p>select a different unit.</p> <p>When the unit entered is in a final status, click Clear to clear the Unit ID and Product fields and enter another unit.</p>	<p>unit.</p> <ul style="list-style-type: none"> • Displays the editable unit information based on the user's security level. <hr/> <p>NOTES</p> <p>VBECS does not allow selection of a unit when it:</p> <ul style="list-style-type: none"> • Is currently assigned to a patient. • Has a final unit status. • Has a previous patient assignment record and the user is a Lead or Enhanced Technologist. <p>There is no override. VBECS clears the screen and allows the user to select a new unit.</p>
<p>3. Review and enter or edit the data (Figure 88),</p> <p>or</p> <p>Go to Inactivate a Unit.</p>	<ul style="list-style-type: none"> • Displays tabs to select and edit: <ul style="list-style-type: none"> ◦ Unit information ◦ Unit ABO/Rh confirmation ◦ Antigen typing • Displays fields to add or deactivate unit information and an option to deactivate a unit based on the condition of the unit and the user's security level. <hr/> <p>NOTES</p> <p>The user may add data to a field with no previous entry.</p> <p>Users may edit unit volume, add special testing and RBC antigens, mark a unit biohazardous, and change the patient restriction.</p> <p>When a blood product code in the RED BLOOD CELLS, WHOLE BLOOD, or APHERESIS RED BLOOD CELLS ICCBBA component class is entered, VBECS enables data entry of sickle cell negative and antigen typing. The CMV negative data entry option is always enabled.</p> <p>An antigen type may be positive or negative for any antigen listed as a VBECS corresponding antigen in Appendix B: Table 13: Antibody and Antigen Table: Irregular Antibodies. An antigen may not be positive and negative for the same unit.</p> <p>Only unit antigen type results originally entered through Incoming Shipment are available for invalidation. The user must invalidate unit antigen type results originally entered through Unit Antigen Typing by selecting Inactivate Unit Testing.</p> <p>VBECS compares the antigen typing interpretation a user enters to previously assigned antigen typings for the same antigen (as assigned during Incoming Shipment, Edit Unit Information, or previous testing). When the previous interpretation was not "inconclusive" and the current interpretation does not match the previously identified antigen type, VBECS warns the user and quarantines the unit. The user may edit the discarded plasma volume only for units with previously reduced volume.</p> <p>VBECS restricts the editing of disease marker testing to autologous units.</p>

User Action	VBECS
	<p>When the user changes the disease marker testing to "POS" or "NFT (not fully tested)," VBECS enables the biohazard indicator.</p> <p>The user may not clear the biohazard indicator for a unit when:</p> <ul style="list-style-type: none"> • The disease marker testing is POS or NFT, • The ISBT 128 Donation Type codes are "X, For Autologous Use Only, Biohazardous" or "3, For Directed Donor Use Only, Biohazardous," or • The unit is quarantined. <p>When editing the "restricted for" patient name for a unit, the user may change the name of a restricted patient but may not save the transaction without a new patient name for the restriction.</p> <p>When the user restricts a unit to a patient, the patient ID is the primary identifier of the intended recipient. When VBECS searches the VistA PATIENT File (#2) and finds a match, it populates the VBECS database with the selected patient name. When the intended recipient is not in the VistA database, VBECS warns the user, clears the screen, and does not allow the unit to be further processed.</p> <p>►► The Enhanced Technologist may edit items associated with a unit and add the new data to the unit. VBECS displays only the latest entry for these fields to the Enhanced Technologist before the unit has any patient associations.</p> <p>►►► The Traditional Supervisor may edit items associated with a unit and add the new data when the unit has a history of patient associations.</p>
4. Select a test for inactivation (Figure 89).	<ul style="list-style-type: none"> • Displays the previously recorded results for the selected test. • Requires a comment to process the inactivation. <p>NOTES</p> <p>The user may view, not edit, previously entered serologic results and interpretations.</p> <p>When the user selects the unit antigen typing test, VBECS displays antigen typing reagent types and interpretations available for invalidation at the user's division. Inherited testing is not available for invalidation.</p> <p>When the user selects the ABO/Rh confirmation test, VBECS displays the unit confirmation test interpretations available for invalidation at the user's division.</p> <p>When a selected unit has current valid tests or XMs, VBECS warns the user and does not allow him to deactivate the unit. There is no override.</p>
5. Verify the test to be deactivated and enter a comment. Click OK to confirm and exit.	<ul style="list-style-type: none"> • Displays "Test Result Invalidated" in the unit record for the selected test. <p>NOTES</p> <p>When the last ABO/Rh confirmation testing is invalidated, VBECS returns the unit to a limited status and requires an ABO/Rh confirmation test according to the component class</p>

User Action	VBECS
	<p>setting in Component Classes.</p> <p>When a unit that requires ABO/Rh confirmation testing no longer has a valid test result in VBECS due to invalidation of previous unit tests, VBECS treats the unit as though no ABO/Rh confirmation testing was performed.</p> <p>VBECS flags invalidated test results as "Test Result Invalidated" and captures details for inclusion in an Exception Report.</p> <p>VBECS does not inactivate workload associated with the invalidated test.</p> <p>Replacement results are entered according to local policy.</p>
6. Return to Step 1 to edit another unit.	
7. Click OK to save.	<ul style="list-style-type: none"> Updates and saves the data.

Figure 88: Edit Unit Information

The screenshot shows the 'Edit Unit Information' dialog box. On the left, there's a 'Unit Search' section with fields for 'Unit ID*' (W0000000000005) and 'Product*' (E0291V00), and a 'Clear' button. Below this is a large barcode area with a barcode and the identifier 'W0000000000005'. To the right of the barcode is a blood type label 'O Rh NEGATIVE' with a 'VOLUNTEER DONOR' label underneath. Further down are two more barcodes, one labeled 'E0291V00' and the other '0112381233', along with the text 'RBC CPD>AS1' and '26 Aug 2011 12:33'. At the bottom of this section is a red 'X' button next to the text 'Inactivate Unit'. The main body of the dialog has three tabs at the top: 'Unit Information' (selected), 'ABO/Rh Confirmation', and 'Antigen Typing'. Under 'Unit Information', there are fields for 'Volume' (Unit Volume: 350 mL, Discarded Plasma Volume: 0 mL) and 'Special Testing' (checkboxes for CMV Negative, Sickle Cell Negative, and RBC Antigens checked). To the right of these is a scrollable list of RBC Antigens: C NEG, C POS, E NEG, E POS, c NEG, and c POS. The 'ABO/Rh Confirmation' tab shows 'O Rh NEGATIVE'. The 'Antigen Typing' tab is not visible. At the bottom right are 'OK' and 'Cancel' buttons.

* Required Field

Figure 89: Edit Unit Information Invalidate Test

The screenshot shows the VBECS software interface for editing unit information. The main title bar says "VBECS - Edit Unit Information". Below it, a red header bar says "Edit Unit Information". The left side has a "Unit Search" panel with fields for "Unit ID*" (W00000000000005) and "Product*" (E0291V00), and a "Clear" button. The center displays a barcode and unit details: "W0000000000005", "O Rh NEGATIVE", "VOLUNTEER DONOR", "E0291V00", "0112381233", "Expiration Date", and "26 Aug 2011 12:33". Below this is a "Barcode" button. At the bottom left is an "Inactivate Unit" button with a red X icon. The right side contains tabs for "Unit Information", "ABO/Rh Confirmation", and "Antigen Typing". Under "Unit Information", there are fields for "Date Tested" (07/27/2011 16:28), "ABO Interp." (O), and "Rh Interp." (N). A "Selected Test Details" section shows results for Anti-A (0), Anti-B (0), Anti-AB (blank), Anti-D (0), and D Control (blank). Below that are sections for "ABO" (O) and "Rh" (N), each with a "Comment" field. An "Invalidate" section includes a "Comment*" dropdown set to "Lab equipment failure" and a "Details" text area. At the bottom right are "Cancel" and "Invalidate Test" buttons.

Inactivate a Unit

User Action	VBECS
1. Select Blood Units from the main menu. Select Edit Unit Information .	<ul style="list-style-type: none"> Displays options for processing blood units. Displays options for selecting and editing a unit.
2. Scan or enter a unit ID and product code. May cancel and select a different unit.	<ul style="list-style-type: none"> Allows the user to scan or enter a unit ID and product code or select a unit. Displays the editable unit information based on the user's security level. <p>NOTES</p> <p>VBECS does not allow selection of a unit when it:</p> <ul style="list-style-type: none"> Is currently assigned to a patient. Has a final unit status. Has a previous patient assignment record and the user is a Lead or Enhanced Technologist. <p>There is no override. VBECS clears the screen and allows the user to select a new unit.</p>
3. Click Inactivate Unit , as desired.	<ul style="list-style-type: none"> Displays fields to add or deactivate unit information and an option to deactivate a unit, based on the condition of the unit and the user's security level. Displays unit data. Requires a comment to process the unit deactivation.

User Action	VBECS
	<p>NOTES —————</p> <p>►►► When the Traditional Supervisor attempts to deactivate a unit's record, VBECS warns the user.</p> <p>An inactivated unit is not accessible. The record is available in the Audit Trail and Unit History Reports and in division statistical and informational reports. When VBECS inactivates a unit, it also inactivates costs previously associated with that unit.</p> <p>The user may deactivate the unit with prior but no current patient associations. VBECS requires a comment to explain the deactivation.</p> <p>Unit deactivation appears on the Audit Trail Report.</p>
4. Enter a comment. Click OK to save.	<ul style="list-style-type: none"> • Marks the selected unit record "inactive." • Saves the record.

Edit Financial Data

The user enters and/or edits charges associated with a unit after it is logged into inventory.

Assumptions

- The user is logged into the division containing the unit to be edited.

Outcome

- Unit billing information was changed.

Limitations and Restrictions

- When a Technologist enters a unit ID and product code for a unit that is not in final status, VBECS does not enable the Special Test Cost field. An Enhanced Technologist or higher access can change the Special Test Cost field in the Edit Financial Data option.

Additional Information

- None

User Roles with Access to This Option

All users

Edit Financial Data

The user adds and/or edits a unit charge. The user confirms the unit selection, enters the dollar amount of the charge, and reviews and accepts the charge entered.

User Action	VBECS
1. Select Blood Units from the main menu. Select Edit Financial Data .	<ul style="list-style-type: none">• Displays options for processing blood units.• Displays fields for selecting a unit and editing its financial data.
2. Scan or enter a unit ID or product code, or click the find button to select a unit from the drop-down list.	<ul style="list-style-type: none">• Displays unit data.
3. Edit the Base Cost, Special Test Cost, and/or Return Credit fields (Figure 90).	<ul style="list-style-type: none">• Allows the user to edit the fields in the Financial Information area.• Displays the edited data. <p>NOTES —————</p> <p>►► The Enhanced Technologist may edit existing amounts.</p> <p>VBECS allows the user to enter a single special typing cost (the sum of charges from \$1.00 to \$999.00, other than the base unit cost, added to a unit by the shipper) to associate with a unit.</p> <p>When no special typing costs are entered through Incoming Shipment, a user may enter them in this option until the unit is in a final status.</p> <p>►► The Enhanced Technologist may edit a previously entered unit cost (unit base cost, return credit or special testing cost).</p> <p>VBECS displays a return credit in U.S. currency (defined by the return credit percentage multiplied by the base unit cost for the blood product type in Blood Products). A user may enter a return credit amount only when VBECS assigned the return credit percentage through Outgoing Shipment or Discard or Quarantine. Otherwise, VBECS disables the Return Credit field.</p>
4. Review the data, click OK to save, and click OK again to confirm the save.	<ul style="list-style-type: none">• Requests confirmation to update the database.• Saves and updates the unit record.
5. Repeat Steps 2–4, as needed, to update additional unit financial data.	
6. Click Cancel to exit.	

Figure 90: Edit Unit Financial Data

The screenshot shows the 'Edit Unit Financial Data' window in VBECS. On the left, there's a sidebar with a 'FINANCE' icon. The main area has a 'Unit Search' section with fields for 'Unit ID*' (W0000000000000005) and 'Product*' (E0291V00), along with a 'Clear' button. Below this is a large barcode and a unit label with the text: 'Property Intended Recipient', 'Caution: Federal Law prohibits dispensing without a prescription', 'VOLUNTEER DONOR', 'Rh NEGATIVE', 'E0291V00', 'RBC CPD>AS1', '0112381233', and '26 Aug 2011 12:33'. To the right, there are tabs for 'Details' (selected) and 'Testing'. Under 'Details', fields include 'Donation Type' (Volunteer Allogeneic Donor), 'Date Received' (07/27/2011 12:33), 'Unit Status' (Available), and 'Volume (mL)' (350). A 'Financial Information' section contains fields for 'Base Cost' (\$1.00), 'Special Test Cost' (\$0.00), and 'Return Credit' (\$0.00). At the bottom are 'OK' and 'Cancel' buttons. A note at the bottom left says '* Required Field'.

Free Directed Unit For Crossover

The user removes a patient restriction from a directed unit that is eligible for crossover, making that unit available for any patient, according to system rules.

Assumptions

- The division is “full service.”
- The unit and the user are in a blood bank in the same division.
- The unit is restricted for a patient.
- The unit is not autologous.
- The unit is not in a final status.

Outcome

- The unit has no restrictions.
- VBECS documents the restriction and its removal in the unit record.

Limitations and Restrictions

- The unit must not be in a final status to release the restriction.
- This option does not allow the removal of patient restrictions for most of the ISBT donation types that are eligible for crossover.
- This option cannot process multiple directed units simultaneously.

Additional Information

- This option is disabled for transfusion-only facilities.

User Roles with Access to This Option


►► Enhanced Technologist

Free Directed Unit For Crossover

The user releases the restriction on a directed donor unit to make the unit available in general inventory.

User Action	VBECS
1. Select Blood Units from the main menu. Select Free Directed Unit For Crossover .	<ul style="list-style-type: none">• Displays options for processing units.• Prompts the user to scan or enter the unit ID and product code.
2. Scan or enter the unit ID and product code for each unit.	<ul style="list-style-type: none">• Displays the unit information. <p>NOTES —————</p> <p>When a user adds a directed unit to an undirected pool, VBECS queries the user regarding releasing the unit's restriction.</p> <p> The unit must have a donation type of only Volunteer Directed, Eligible for Crossover. If not, VBECS warns that the unit has an invalid donation type and is in a final status."</p> <p> The unit must not be assigned to a patient to release the restriction. VBECS warns that the unit is currently assigned and that he must release the assignment before continuing.</p>
3. Enter a date and time and possibly a different user, or accept the defaults (Figure 91). Enter the required comment.	<ul style="list-style-type: none">• Displays the unit information.•  Emits an audible alert, requires a comment, and captures details for inclusion in an Exception Report (exception type: released patient restriction).
4. Click OK to process the unit and add it to available inventory.	<ul style="list-style-type: none">• Displays the return time, the editable (past date and time) field (default: current date and time), and the list of users for selection.• Saves the information to the database.• Changes the donation type of the unit to "Voluntary Allogeneic" and removes the restriction associated with the unit.
5. Click Yes to confirm the save.	

Figure 91: Free Directed Unit For Crossover

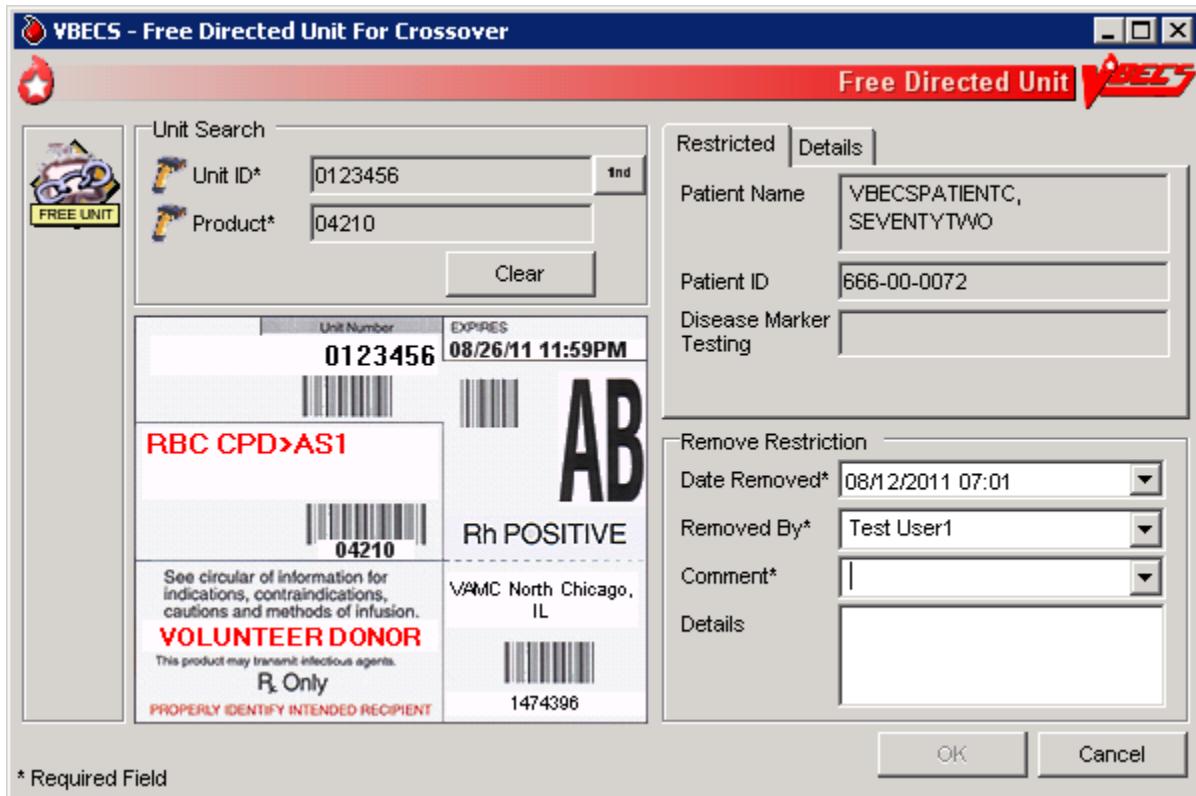

Discard or Quarantine

The user discards, quarantines, or releases from quarantine one or more units of blood.

Assumptions

- The unit of blood is in the same division as the user.

Outcome

- Discarded units are not available.
- Quarantined units are available for limited processing, not for patient assignment.
- Units released from quarantine are available.

Limitations and Restrictions

- When a user retroactively updates a unit status through this option, the Unit History Report displays the updated information, but does not display the date the change was made.

Additional Information

- None

User Roles with Access to This Option

All users

➤➤ Enhanced Technologist (release blood products from quarantine)

Discard or Quarantine

The user changes the status (discard, quarantine, release from quarantine) of one or more units of blood. The user must indicate the ID and product code of the units to be updated. VBECS also allows batch data entry to process all units from an invoice.

A quarantined unit has limited availability for continued processing. Under certain circumstances, VBECS assigns quarantine status automatically; otherwise, the user determines and enters the status. The user must specify whether a discarded unit is for credit or for waste, which determines whether or not a return credit is assigned to the selected unit.

Discard or Quarantine Individual Units

User Action	VBECS
1. Select Blood Units from the main menu. Select Discard or Quarantine .	<ul style="list-style-type: none">• Displays options for processing blood units.• Displays fields for entering unit information.• Allows the user to discard or quarantine individual units.
2. Enter an ID and product code for a blood unit.	<ul style="list-style-type: none">• When a unit is appropriate for processing, displays unit information, including the current status. <p>NOTES —————</p> <p>The user may deselect units prior to editing.</p>
3. Click Select to add the selected unit to the worklist, or click Clear to select another unit.	<ul style="list-style-type: none">• Adds the unit to a temporary list of selected units.
4. Repeat Step 3, as needed, to process multiple blood units as a single batch. Click OK to continue.	<ul style="list-style-type: none">• Displays selected units.
5. In the Activity area, click the Discard , Quarantine , or Release from Quarantine radio button (Figure 92).	<p>NOTES —————</p> <p>►► The Enhanced Technologist may release a unit from quarantine.</p>
6. Select or enter a comment in the Comment field to explain the change in status, when indicated. Enter details, as needed.	<ul style="list-style-type: none">• Displays entered data. <p>NOTES —————</p> <p>When a status of "Discard" is assigned to a unit, the user must enter a comment.</p> <p> When a blood unit is assigned an indicator of "quarantined," VBECS emits an audible alert, requires a comment, and captures details for inclusion in an Exception Report (exception type: Quarantine, Discard or Release from Quarantine, depending on the user action).</p> <p> When a unit is released from quarantine, a new status must be assigned. VBECS assigns the previous status as the default status. When the unit expiration date is in the past, VBECS emits an audible alert, warns the user, and asks whether he wishes to continue.</p> <p>Yes accepts the default status, requires a comment, and captures details for inclusion in an Exception Report (exception type: expired unit released from quarantine). No stops the</p>

User Action	VBECS
	<p>process; the unit remains marked as quarantined.</p> <p>When the user assigns a status of “Discard,” he must indicate whether the discard is waste or credit.</p> <p>When VBECS processes a unit as a discard that is waste, it does not apply a return credit to the discarded unit. (There is no refund of monies from the blood shipper.)</p> <p>When VBECS saves a unit as “discard for credit,” it applies a default return credit. A user may edit the default return credit applied to the unit before updating the database.</p> <p>When a unit has had a quarantine indicator assigned, VBECS releases any patient assignments. Patient restrictions associated with the unit remain intact.</p> <p>When a user attempts to release a unit from quarantine and the unit was quarantined because of a data conflict, VBECS checks the inconsistency before allowing the unit to be released. VBECS warns the user. There is no override. The user must invalidate the incorrect information or test to resolve the data conflict.</p>
7. Enter the date and time of the status update. Click the Waste or Credit radio button.	<ul style="list-style-type: none"> • VBECS defaults the date and time of the status update to the date and time the user opened the Discard or Quarantine option • Displays entered data. <p>NOTES —————</p> <p>The default date and time are the current date and time.</p> <p>When the user updates the status of a unit and enters a date or time in the past, VBECS compares the retrospective date and time with the last record update for the unit. VBECS allows the user to update a unit record with a date and time between the last record update and the current date and time.</p> <p> When the user attempts to update a unit record prior to the last record update, VBECS emits an audible alert, warns the Enhanced Technologist, and instructs him to check the date and time entered. There is no override.</p> <p>►►► When the user is a Lead Technologist, VBECS also requires a comment and asks whether the user wishes to continue. Yes accepts the date and time entered, requires a comment, and captures details for inclusion in an Exception Report (exception type: retrospective update prior to last record update). No allows the user to correct the entry.</p>
8. Review the unit status update information and click Save to save.	<ul style="list-style-type: none"> • Evaluates the unit's current status to determine whether it is eligible for the selected function. • Displays appropriate error messages when the unit cannot be processed. • Prompts the user to confirm the save.
9. Click Yes to confirm the save.	<ul style="list-style-type: none"> • Saves the data. <p>NOTES —————</p> <p>During batch data entry of “Discard,” “Quarantine,” and “Release</p>

User Action	VBECS
	from Quarantine," VBECS saves transaction details, including user-entered comments, to each blood unit included in the batch.

Figure 92: Discard or Quarantine Units

The screenshot shows the "VBECS - Discard or Quarantine Units" application window. The title bar includes the VBECS logo. The main area displays a table with one row of data:

Unit ID	Product Name	Code	Status	Credit	Quarantined?
1941954	RBC AS1	04210	Available	\$10.00	No

To the left of the table is a vertical toolbar with a "DISCARD" icon. Below the table is a scroll bar. At the bottom of the window are three groups of controls:

- Activity***: Radio buttons for "Discard" (selected), "Quarantine", and "Release from Quarantine".
- Discard Unit(s)**: A dropdown menu for "Comment*" set to "Incinerated", a text input field for "Details", and a dropdown menu for "Discard Date*" showing "02/02/2006 15:21".
- Reimbursement***: Radio buttons for "Waste" (selected) and "Credit".

At the bottom right are "Save" and "Cancel" buttons. A note at the bottom left indicates "* Required Field".

Discard or Quarantine Units by Invoice

User Action	VBECS
1. Select Shipments from the main menu. Select Discard or Quarantine . Select an invoice and click OK .	<ul style="list-style-type: none"> • Displays options for processing blood product shipments. • Displays fields for entering shipment information. • Allows the user to discard or quarantine units by invoice.
2. Click (clear) one or more check boxes to remove units from the invoice. Click OK .	<ul style="list-style-type: none"> • Allows the user to select individual units. • Displays unit information, including the current status included in the invoice. • Allows the user to deselect units from the batch prior to editing. <p>NOTES</p> <p>For batch processing, the user may select incoming shipment invoices from the previous five days processed in the same division as the user.</p> <p>Only incoming shipment invoices processed within the user's division may be accessed for batch processing.</p>
3. In the Activity area, click the Discard, Quarantine, or Release from Quarantine radio button.	<p>NOTES</p> <p>►► The Enhanced Technologist may release a unit from quarantine.</p>
4. Select or enter a comment in the Comment field to explain the change in status, when indicated. Enter details, as needed.	<ul style="list-style-type: none"> • Displays entered data. <p>NOTES</p> <p>When a status of "Discard" is assigned to a unit, the user must enter a comment.</p> <p>When a blood unit is assigned an indicator of "quarantined," VBECS emits an audible alert, requires a comment, and captures details for inclusion in an Exception Report (exception type: Quarantine, Discard or Release from Quarantine, depending on the user action).</p> <p> When a unit is released from quarantine, a new status must be assigned. VBECS assigns the previous status as the default status. When the unit expiration date is in the past, VBECS emits an audible alert, warns the user, and asks whether he wishes to continue.</p> <p>Yes accepts the default status, requires a comment, and captures details for inclusion in an Exception Report (exception type: expired unit released from quarantine). No stops the process; the unit remains marked as quarantined.</p> <p>When the user assigns a status of "Discard," he must indicate whether the discard is waste or credit.</p> <p>When VBECS processes a unit as a discard that is waste, it does not apply a return credit if not applied to the discarded unit. (There is no refund of monies from the blood shipper.)</p>

User Action	VBECS
	<p>When VBECS saves a unit as “discard for credit,” it applies a default return credit. A user may edit the default return credit applied to the unit before updating the database.</p> <p>When a unit has had a quarantine indicator assigned, VBECS releases any patient assignments. Patient restrictions associated with the unit remain intact.</p> <p>When a user attempts to release a unit from quarantine and the unit was quarantined because of a data conflict, VBECS checks the inconsistency before allowing the unit to be released. VBECS warns the user. There is no override. The user must invalidate the incorrect information or test to resolve the data conflict.</p>
5. Enter the date and time of the status update. Click the Waste or Credit radio button.	<ul style="list-style-type: none"> • VBECS defaults the date and time of the status update to the date and time the user opened the Discard or Quarantine option • Displays entered data. <p>NOTES —————</p> <p>The default date and time are the current date and time.</p> <p>When the user updates the status of a unit and enters a date or time in the past, VBECS compares the retrospective date and time with the last record update for the unit. VBECS allows the user to update a unit record with a date and time between the last record update and the current date and time.</p> <p> When the user attempts to update a unit record prior to the last record update, VBECS emits an audible alert, warns the Enhanced Technologist, and instructs him to check the date and time entered. There is no override.</p> <p>>>> When the user is a Lead Technologist, VBECS also requires a comment and asks whether the user wishes to continue. Yes accepts the date and time entered, requires a comment, and captures details for inclusion in an Exception Report (exception type: retrospective update prior to last record update). No allows the user to correct the entry.</p>
6. Review the unit status update information and click Save to save.	<ul style="list-style-type: none"> • Evaluates the unit's current status to determine whether it is eligible for the selected function. • Displays appropriate error messages when the unit cannot be processed. • Prompts the user to confirm the save.
7. Click Yes to confirm the save.	<ul style="list-style-type: none"> • Saves the data. <p>NOTES —————</p> <p>During batch data entry of “Discard,” “Quarantine,” and “Release from Quarantine,” VBECS saves transaction details, including user-entered comments, to each blood unit included in the batch.</p>

Select Units

A user selects a patient component order and assigns one or more products to a patient for possible transfusion.

Assumptions

- An appropriate current order was accepted through Accept Orders: Accept an Order.
- The division may or may not allow electronic crossmatch (eXM).
- When a division is “transfusion only,” eXM is not available. The user must enter a serologic crossmatch (XM) interpretation.

Outcome

- The status of the unit and the associated component order task status are updated.
- VBECS indicates that units are available for issue when units do not require XM or other additional testing.
- VBECS allows the user to print a Blood Transfusion Record Form (BTRF) and a Caution Tag for available blood products.
- When serologic XM is required, VBECS automatically adds units to the XM data entry grid when Patient Testing: Record a Crossmatch is accessed.
- When eXM is enabled for the facility, the patient is eligible for eXM in accordance with system rules (see Table 23: Rules for Electronic and Serologic Crossmatch), and the user selects a unit, VBECS assigns the unit to the patient without further testing.

Limitations and Restrictions

- When a user checks the “D Pos” or “D Neg” check box to modify a search for blood, VBECS returns no units, regardless of the availability of units.
- VBECS does not allow a user to select a pool of mixed ABO/Rh units for a patient. Instead, the user must select random units for the patient and pool the mixed ABO/Rh units to allow issue.
- VBECS displays a component order as “filled” when the number of units transfused and the number of units ordered are the same. CPRS displays a component order as “completed” when the number of units ready for issue and the number of units ordered are the same. Therefore, CPRS displays the order as “completed” before VBECS displays the order as “filled.”
- Changing a patient's antigen negative requirement or antibody identified after a unit is crossmatched does not provide an alert at Issue Unit regarding the type of crossmatch (electronic or serologic) performed on a unit available for issue.

Additional Information

- Active component orders and the user must be in the same division.
- The user may select only active units and those not in an inactive or final status.
- Multiple users may select units associated with a component order.
- When a user selects a patient order, VBECS searches for available autologous and directed units regardless of the unit’s current division within a multidivisional database, displays the unit information, and allows the user to assign the unit when it is in the user’s division.
- A user may select units for only one ordered component class at a time.
- An existing specimen and/or order remains available to select additional units until its expiration date.
- The selection of blood units does not complete an order.

- A user must accept or cancel each component task in its entirety.
- The user reserves autologous and directed units and all units processed in a “transfusion-only” facility, regardless of donation type, for a specific patient. VBECS displays this information to all divisions in a multidivisional database and accommodates immediate selection of the unit for the patient when the user is in the same division.
- When VBECS compares a patient’s current and previous ABO/Rh results, it does not evaluate previously recorded ABO/Rh results marked “Test Result Invalidated.”
- An ABO/Rh confirmation test interpretation must match the unit ABO/Rh at login and when repeat confirmation tests are not invalidated. A user may not select a unit for a patient based on ABO/Rh at login. The unit must have a satisfactory confirmation result before VBECS considers it available for selection.
- When a user accepts an order without a specimen, he may associate the specimen when he selects units for XM.
- ABO/Rh confirmation testing is not required for components in the OTHER component class that do not traditionally contain red blood cells. A user may request a XM for these components in accordance with local policy.

User Roles with Access to This Option

All users

Select Units

A user selects a patient’s component order and assigns one or more products to a patient. The user may then enter individual units or select units from a list of selectable, potentially compatible units for each component order.

User Action	VBECS
1. Select Patients from the main menu and select Patient Testing , or Select Blood Units from the main menu and select Select Units .	<ul style="list-style-type: none"> • Displays options for processing patient-related functions. • Displays the Pending Task List (PTL) and PTL search parameters in the Diagnostic Tests tab or in the Component Orders tab.
2. Enter or select PTL search parameters, if appropriate (Figure 93). Click the check boxes in the PTL to select one to four component orders. Click OK to continue (Figure 94).	<ul style="list-style-type: none"> • Displays information for each component order. • When the selected order class is RED BLOOD CELLS and the division enabled eXM, evaluates the patient for eXM unit processing. <p>NOTES —————</p> <p> See Table 6 for alerts that may occur during this option.</p>
3. Steps 3 through 8 are optional. If a specimen is already associated with the order, skip to Step 9. To associate a specimen with an emergency-issued unit requiring XM, select (highlight) an issued unit in the Order Group tree view in the left pane (Figure 95).	<ul style="list-style-type: none"> • Enables the Associate Specimen check box. <p>NOTES —————</p> <p>If eXM is enabled, complete the TAS prior to associating the specimen UID with the emergency-issued units.</p>
4. Click the Associate with Specimen check box.	<ul style="list-style-type: none"> • Opens the Select Patient Specimen window.

User Action	VBECS
Click the find button to select a specimen UID (Figure 95).	
5. Scan or enter the correct specimen UID and click Search .	<ul style="list-style-type: none"> Displays the selected patient specimen information.
6. Click OK .	<ul style="list-style-type: none"> Saves the specimen information and returns the user to the Select Units window. Displays the associated specimen UID in the Specimen UID field.
7. Repeat Steps 3–6 to associate the sample with other emergency-issued units, or click OK to save.	<ul style="list-style-type: none"> When the user clicks OK, displays a message asking the user if he wishes to complete serologic crossmatch.
8. Click Yes to complete the serologic crossmatch. Click No to return to the VBECS desktop.	<ul style="list-style-type: none"> Allows the user to enter results in the Serologic Crossmatch (XM) cells in the patient testing grid.
9. Click Add Units to add units to the component orders selected.	<ul style="list-style-type: none"> Searches the database and displays available restricted autologous and directed units that the user must select before selecting allogeneic units of the same ICCBBA component class. Displays the unit ID, product type, unit ABO/Rh, unit expiration date, specimen UID, if applicable, of units assigned to the patient in the division's blood bank on a different specimen or order. The user may release and reassign the units to the current order. When eXM was configured for the facility, displays the eXM status (eligible or not eligible) for the patient. <p>NOTES —</p> <p> See Table 6 for alerts that may occur during this option.</p> <p>Units that are selectable with an override do not appear on the list of available units and must be scanned or entered by the user.</p> <p>When no valid ABO exists for the patient, the user may select components, including autologous units:</p> <ul style="list-style-type: none"> For RBC or WB order: <ul style="list-style-type: none"> Group O, Rh positive or O Rh negative RBC WHOLE BLOOD is not selectable For FFP order only, group AB, Rh positive or Rh negative For PLT, CRYO, or OTHER order, all available units are selectable.
10. Scan, enter or select *allogeneic units. Units may be selected from the Order Group tree view or by clicking the find button for a list of compatible units. * when typing in a unit ID and Product Code, use the tab key or mouse to leave the field. Enter or select a date in the Selection Date field. Each unit may have its own selection date and time.	<ul style="list-style-type: none"> Maintains a list of selected units. The user may deselect selected units. Evaluates the selected unit: <ul style="list-style-type: none"> For the ABO/Rh compatibility of each unit and warns the user. Against patient Transfusion Requirements (TRs) and warns the user. For unit antigen negative compatibility and warns the user. Performs additional checks that include days remaining until the unit's expiration date, the current status, assignments to other patients, quarantine and biohazard indicators, when applied, and associated special testing. <p>NOTES —</p>

User Action	VBECS
Click Add Units to display the unit selection form.	<p> See Table 6 for alerts that may occur during this option.</p> <p><i>When the user selects a unit in the tree view, he may accept or edit the current date and time displayed; the date and time must be in the past.</i></p> <p>Restricted (autologous or directed) unit selection eligibility criteria:</p> <ul style="list-style-type: none"> • Restricted For Patient ID: must match • Restricted For Patient Name: must match • Unit expiration date must be in the future. • Component class (transfusion only requires restriction of all units) <p>The criteria for selecting or including an individual unit in a pick list include evaluating a unit for:</p> <ul style="list-style-type: none"> • Considered substitutable for the order • Same division as the user • Current location is in the blood bank • Has no quarantine indicator • Is not already assigned or crossmatched to the patient • Is not “restricted” for a different patient • Meets ABO/Rh compatibility requirements specific to the component class • Future expiration date <p>VBECS lists only ABO and Rh compatible units for selection.</p> <p>When ABO/Rh results were not created for the specimen associated with the component order and the order is not designated “emergency issue,” VBECS allows the user to select units using historic records, and notifies the user that ABO/Rh testing must be performed at this division on this patient before any selected units can be issued.</p> <p>VBECS displays and allows the user to complete a pending XM associated with a patient order without a final unit status.</p>
11. Repeat Step 10 until all units are selected.	<ul style="list-style-type: none"> • Requests the user to confirm the list of selected blood units. <hr/> <p>NOTES</p> <p>The user may select or deselect additional units.</p>
12. Click OK to confirm the selection of units.	<ul style="list-style-type: none"> • Determines whether the units are available for issue or whether they require additional modification and/or testing based on the component class. • Updates the database, as appropriate. • Sends selected blood units with patient details to the BCE COTS system via an HL7 message. <hr/> <p>NOTES</p> <p>When the total number of blood units eligible for issue, issued, or transfused is greater or equal to the number of ordered units, VBECS sends an order completion message to CPRS.</p> <p> See Table 6 for alerts that may occur during this option.</p>

User Action	VBECS
	<p>When eXM is enabled at the division and blood units require XM, VBECS determines whether the patient and units are eligible for eXM. All parameters listed in Appendix B: Table 23: Rules for Electronic and Serologic Crossmatch must be met. Units not eligible for eXM must have serologic XM performed.</p> <p>When eXM is not enabled, the user must enter serologic XM results for red blood cell classes or when a user orders an optional XM for OTHER component classes. VBECS does not change the unit status to "crossmatched," make units available for issue, or allow the Caution Tag or BTRF to be printed until a user enters a valid XM result, unless an emergency-issued unit is being processed.</p>
13. Select another component order to continue selecting units, or exit.	<ul style="list-style-type: none"> • Displays orders for the next patient.

Figure 93: Pending Task List Component Orders

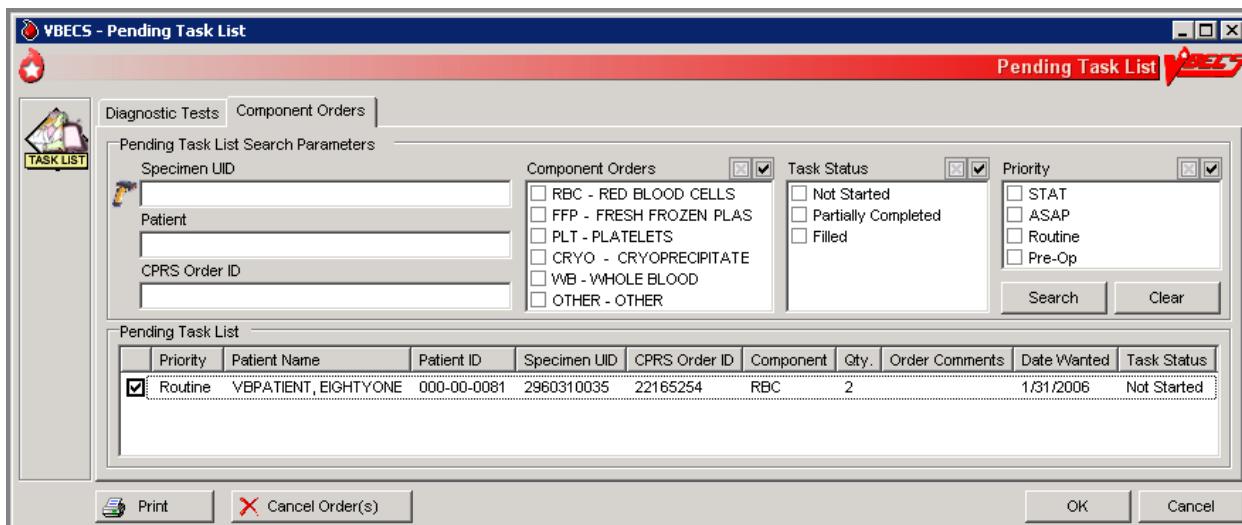


Figure 94: Select Units for a Patient


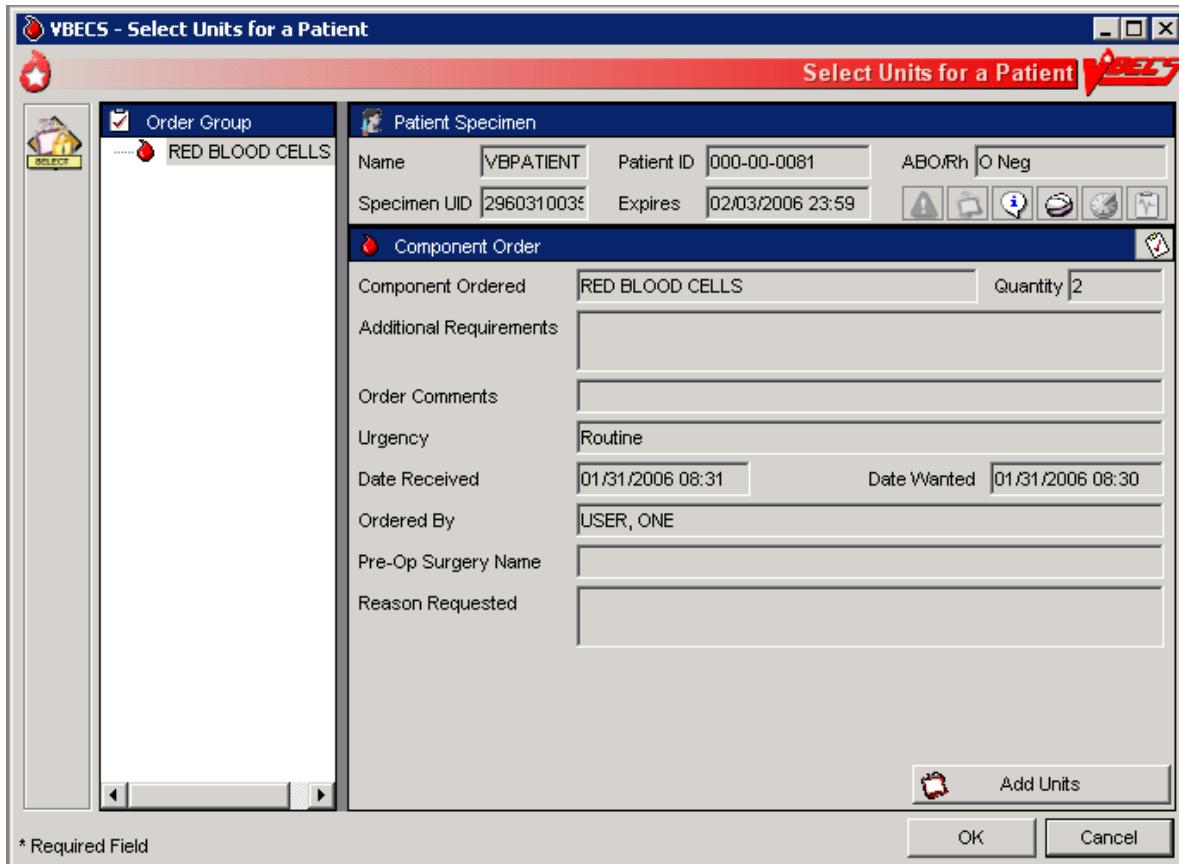


Figure 95: Associate a Specimen with an Emergency-Issued Unit

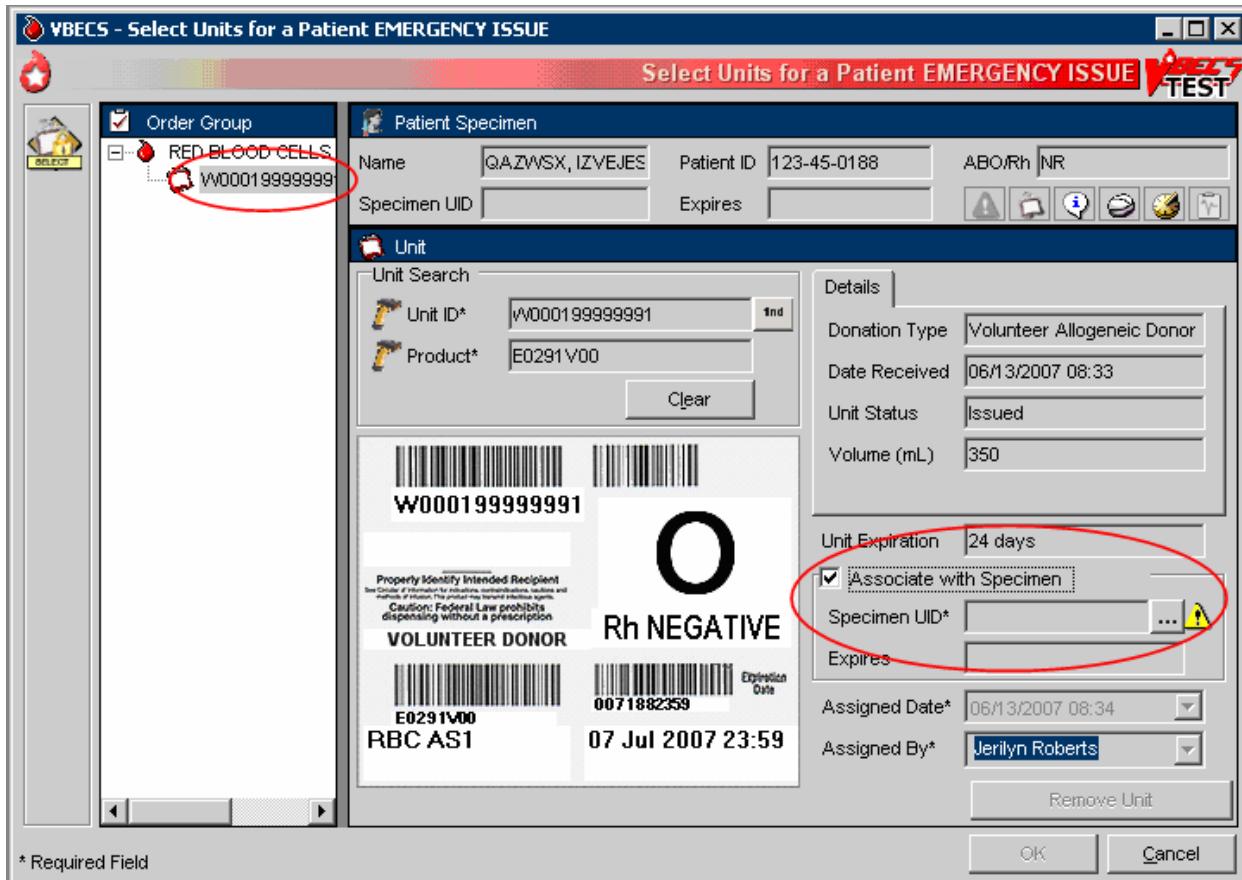

Alerts

Table 6: Alerts That May Occur in Select Units

Step	Alerts
2	VBECS displays and emits an audible alert when there are Special Instructions (SIs) or TRs for the patient.
2	When the user selects blood units that require XM, and eXM is enabled at the division but the current patient is not eligible for eXMs, VBECS displays a warning that includes the reason for the ineligibility. Selected units must have serologic XM performed. The user may not print the Caution Tag or the BTRF. VBECS does not change the unit status to "crossmatched" until the user enters a valid serologic XM result in Patient Testing: Record a Crossmatch.
2	When the user clicked the Emergency Order check box during specimen acceptance and selects blood units requiring crossmatch, VBECS does not allow the user to perform an electronic or serologic crossmatch.
9	When a current specimen is required for the selected component class based on Component Classes and one does not exist for the patient, the user may not select ABO type-specific blood products. VBECS notifies the user that he is selecting units using the emergency issue compatibility rules.
9	When the patient record is marked as having a previous ABO/Rh discrepancy that was justified, a current specimen is required. The user must select units based on the current specimen ABO/Rh. When the current specimen ABO/Rh was not performed, VBECS warns that the patient had a previous ABO/Rh discrepancy and instructs the user to perform ABO/Rh on the current specimen to continue using normal rules and policies or follow emergency issue rules and policies.

Step	Alerts
9	When there is a current specimen and ABO/Rh testing was performed, VBECS compares the ABO/Rh of the current specimen to the historical ABO/Rh. When there is a discrepancy, the ABO type is inconclusive, and there is no previous justification for the discrepancy on the patient record, VBECS notifies the user that the ABO/Rh on the current specimen does not agree with the historical ABO/Rh and instructs the user to resolve the discrepancy to continue using normal rules and policies, and recommends that the user release previously selected units or follow emergency issue rules and local policies.
9	When a repeat ABO/Rh is performed on the patient specimen and the interpretation does not match the original ABO/Rh interpretation.
9	When a user selects a blood unit that is evaluated as ABO incompatible (based on comparison of the patient's current ABO/Rh entered by the user with the ABO/Rh on record), VBECS emits an audible alert and warns the user that the unit is not ABO compatible for this patient and that the user may not select the unit. There is no override; VBECS clears the screen and allows the user to select another unit.
10	The user may override the selection of restricted units. VBECS warns that more restricted units were not selected for the patient and asks whether the user wishes to continue. The override requires a comment and captures details for inclusion in an Exception Report (exception type: less restrictive unit issued).
10	VBECS warns the user when he selects type-specific plasma products that do not require a specimen, according to the division configuration of component classes, for a patient with a history of justified ABO/Rh change.
10	Once a valid XM interpretation is verified and the unit is compatible, VBECS does not allow a user to select that unit (ID and product code) for the same specimen. The user may select the unit for a XM with the same patient's new specimen after the unit is released from the current assignment. VBECS warns that the unit has a valid current XM with this patient. The user clears the warning and may continue or enter a new unit number for selection.
10	When a selected unit is currently assigned to another patient, VBECS warns the user and asks whether the user wishes to continue to select for this patient. Yes allows the user to continue. No clears the screen and allows the user to select a new unit.
10	When a user is a Technologist and selects an expired unit, VBECS warns him. There is no override. ►► When the user is an Enhanced Technologist or above, VBECS warns that the unit is expired and asks whether the user wishes to continue to select. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: expired unit selected).
10	When a unit is selected for a patient, VBECS displays a tool tip stating that the unit may not be selected or removed.
10	When the user selects a unit and these conditions exist, VBECS warns the user. There is no override. VBECS clears the screen and the user may select another unit: <ul style="list-style-type: none"> • For XM for the same patient, VBECS warns that the unit was already selected for the patient. • That is not substitutable for the order, VBECS warns that the selected unit is not valid for order. • That is quarantined, VBECS warns that the unit is quarantined and may not be selected. • That is on a pending outgoing shipment list, VBECS warns that the unit is currently included on a pending outgoing shipment invoice and must be removed from the shipment invoice before it can be selected. • That is not in the blood bank or is issued to another patient, VBECS warns that the unit is not in the blood bank and may not be selected. • That is restricted for a different patient, VBECS warns that the unit is restricted for a different patient and may not be selected. • The confirmation is not performed or is incomplete, VBECS emits an audible alert and warns that the unit ABO/Rh confirmation was not performed or is incomplete and that the unit may not be selected.
10	When the selected unit had an incompatible crossmatch result with the current patient in the past, VBECS warns the user. The user may continue with this unit and proceed with electronic or serologic crossmatch testing.
10	When the user selects an antigen positive (or inconclusive for testing) unit and the patient has the corresponding antigen negative requirement, VBECS warns the user and requires action determined by the user's security level.

Step	Alerts
10	►► When the user is an Enhanced Technologist or below, VBECS emits an audible alert and warns that the patient has a history of an antibody, that the unit does not have a negative test for the corresponding antigen, and that the user does not have the proper security level to select the unit. There is no override.
10	►►► When the user is a Lead Technologist or above, VBECS emits an audible alert, warns that the patient has a history of an antibody and that the unit is positive or inconclusive for the corresponding antigen, and asks whether the user wishes to continue to select the unit. No allows the user to continue to select units for the patient. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: antigen positive units selected).
10	When a unit was not antigen typed for the corresponding antigen, VBECS emits an audible alert, warns that the unit must be typed and found negative to be available for this patient, and asks whether the user wishes to continue to select this unit. The user may continue to select units for the patient. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: antigen positive units selected). See Appendix B: Table 13: Antibody and Antigen Table: Irregular Antibodies.
10	When an Rh negative patient has an active anti-D antibody or a D antigen negative requirement set in Special Instructions and Transfusion Requirements, VBECS verifies the presence or absence of the corresponding antigen based on the ABO/Rh of the selected blood unit. When a user selects an Rh positive blood unit for a patient with anti-D, VBECS warns the user and requires action based on the user's security level. The antigen negative requirement for anti-D is D, not AGwD. AGwD is not a selectable antigen negative requirement in Special Instructions & Transfusion Requirements: Antigen Negative RBCs.
10	For an Rh negative, unknown, or inconclusive patient: when the selection of Rh positive or Rh pooled units from the WHOLE BLOOD, RED BLOOD CELLS, PLATELETS, and OTHER VBECS component classes are individually entered and confirmed as required, VBECS emits an audible alert, warns that the patient is Rh negative and the selected unit is Rh positive, instructs the user to check the patient sex and hospital policy before continuing, and asks whether the user wishes to continue to select the unit. No clears the unit data and allows the user to select a different unit. Yes requires additional comment (unit selection canned comment context) and captures details for inclusion in an Exception Report (exception type: Rh positive units selected).
10	When a user selects blood units to fill an OTHER component order, VBECS allows him to order an optional XM. When the user requests an optional XM order, he must specify a specimen UID to associate with the XM. When the selection process is complete and eXM is enabled at the site or added to the serologic XM data grid, VBECS evaluates the units for eXM. See Appendix B: Table 20: Enable Crossmatch Option for OTHER Product Types.
10	When a patient has TRs other than CMV negative or Sickle Cell negative, VBECS verifies whether the selected unit satisfies the requirement based on system rules TT_3.04 Blood Unit Substitution and the Fulfillment of Transfusion Requirements. When the TR is not satisfied, VBECS emits an audible alert, warns that the unit does not satisfy all of the patient's TRs, and asks whether the user wishes to continue to select the unit. No clears the screen. Yes allows the user to continue without comment.
10	When a patient has a TR of CMV negative or Sickle Cell negative, VBECS verifies the selected unit for compliance. When the selected unit does not fit the criteria, VBECS emits an audible alert and warns the user: <ul style="list-style-type: none">• ►► Enhanced Technologist or below: VBECS warns that the patient requires CMV negative or Sickle Cell negative blood, that the selected unit is not CMV negative or Sickle Cell negative, and that the user does not have the proper security level to select the unit. There is no override. VBECS clears the screen, and the user must select another unit.• ►►► Lead Technologist or above: VBECS warns that the patient requires CMV negative or Sickle Cell negative blood, that the selected unit is not CMV negative or Sickle Cell negative, and asks whether the user wishes to continue to select. No clears the screen and the user must select a new blood unit. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: Transfusion Requirement incompatible unit selected).

Step	Alerts
10	<p>When a user selects a unit with a biohazardous indicator for a patient, VBECS emits an audible alert, warns the user, and asks whether the user still wishes to select the unit.</p> <p>No clears the entry and allows the user to select a new unit. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: biohazardous unit selected).</p>
10	<p><input checked="" type="checkbox"/> When the user selects blood components from the frozen product types for a patient, VBECS emits an audible alert and warns that the selected blood unit must be thawed and may require additional processing prior to issue. The unit is not eligible for printing the BTRF or Caution Tag. VBECS updates the blood unit status to "assigned" and the unit is no longer available for issue.</p>
10	<p>When a user selects a unit requiring XM, eXM is enabled, and the patient is eligible, the selected unit ABO/Rh must be compatible with the patient ABO/Rh. VBECS displays a message based on Appendix B: Table 23: Rules for Electronic and Serologic Crossmatch. If the user continues to select that unit in which VBECS requires a serologic XM, VBECS does not allow the Caution Tag or BTRF to be printed, and does not update the unit status to "crossmatched" under normal system rules and policies. When the selected unit is available for emergency issue to the patient, VBECS allows the Caution Tag and BTRF to be printed.</p>
12	<p>When eXM is enabled and the patient and selected units are eligible, VBECS notifies the user and asks whether he wishes to continue with eXM.</p> <p>Yes assigns a XM result of "Compatible: Electronically Crossmatched." VBECS updates the database with the XM result, date and time of XM, XM technologist, division, specimen UID, specimen phlebotomist, and specimen expiration date and time. VBECS allows a Caution Tag and BTRF to be printed and the unit is ready for issue and transfusion. VBECS changes the unit status to "crossmatched."</p> <p>No keeps the unit assigned to the patient and available for data entry of serologic XM. VBECS does not allow the Caution Tag or BTRF to be printed and does not change the unit status to "crossmatched." When a division is "transfusion only," the unit remains assigned to the patient and is available for the user to enter the serologic XM interpretation, VBECS does not allow the Caution Tag or BTRF to be printed, and the unit is not crossmatched.</p>
12	<p>When the user selects a unit requiring XM, the division is defined as "full service," eXM is enabled, and the patient is not eligible for eXM, VBECS warns the user that the selected units require serologic XM and asks whether he wishes to continue.</p> <p>Yes adds the units to the XM data entry grid and transfers the user to Patient Testing: Record Patient Test Results to begin the XM process. No records that a serologic XM test is needed. The user may not print a Caution Tag or Blood Transfusion Record Form; the unit status does not change to "crossmatched."</p>

Issue Blood Components

The user issues one or more patient-assigned blood units to a transporter for storage in a remote location or for possible transfusion.

Assumptions

- The facility has a procedure for the emergency issue of blood products.
- The unit and the user are in the same division.
- The unit is assigned and/or crossmatched for the patient currently selected for issue.
- VistA is available for the user to choose a physician and valid hospital locations.
- All issue-to locations must be valid entries in the VistA HOSPITAL LOCATION (#44), including home transfusion.

Outcome

- The unit status was updated to “issued.”
- The unit’s current location was updated to the “issue-to” location.

Limitations and Restrictions

- VBECS does not check for restricted units (not assigned to the patient) received between selection and issue.
- All issue-to locations must be valid entries in the VistA hospital location file, including home transfusion.
- If VistA link is not available, the user may proceed with the issue process but the current issue-to hospital location is not recorded. All of the associated report records display “VistALink Down” in the Issue To location field.
- ISBT 128 tables include blood product codes labeled as not for manufacturing or transfusion. VBECS does not restrict such products from being entered, issued, or transferred (local policy controls their use) and does not accommodate blood modification for these products.
- When a site has configured subservient VistA institutions to the VBECS division, users can only issue units to locations in the primary VistA. Sites may add locations to VistA to account for the relocations to the subservient VistA institutions.
- BCE information applies ONLY when a BCE COTS product is available for communication and is configured.
- Changing a patient's transfusion requirements (component requirements, antigen negative requirement or antibody identified) after a unit is crossmatched does not provide an alert at Issue Unit regarding the type of crossmatch (electronic or serologic) performed on a unit available for issue.

Additional Information

- When no VistA order is available, the user must process the request manually until VistA provides them to VBECS. The user must print blank Caution Tags and Blood Transfusion Record Forms (BTRFs) to accommodate this rare occurrence. When VistA does not provide a pseudonym or associated pseudo-ID for the patient that allows placement of a VistA order, VBECS cannot create a patient order for emergency issue of blood.
- “Issue with Medical Director Approval” components must be approved prior to this option and do not require special handling by VBECS.
- Emergency issue processing accommodates the issue of blood to patients without completed specimen testing, as required.

- Sites may enter the lot number of blood administration sets and leukoreducing filters when they enter the name of the individual picking up the blood product. VBECS retains this information in the unit history record.
- User must date all signature lines on the BTRF dating of the signatures [i.e., Inspected and issued by, Issued to, First Identifier, Second Identifier, Transfusion Data section completed by, Transfusion Reaction Data section completed by (as needed)].

User Roles with Access to This Option

All users

To correct issue time and other errors, return the unit from issue and reissue it with the correct information.

Issue Blood Components (Routine)

The user issues assigned or crossmatched blood to a transporter for patient transfusion at valid inpatient and outpatient hospital locations, and records issue details, including the transporter's name or identifier.

User Action	VBECS
1. Select Patients from the main menu. Select Issue Blood Components .	<ul style="list-style-type: none"> • Displays options for processing patient-related functions. • Displays an option to scan the Caution Tag or enter a patient name or ID. <p>NOTES _____</p> <p>The user must verify the correct patient selection.</p> <p>The >> Enhanced Technologist or above is authorized to issue expired blood products.</p> <p>The >>> Lead Technologist or above is authorized to issue antigen positive and untyped red blood cells.</p>
2. Select a patient.	<ul style="list-style-type: none"> • Displays patient data, including entries that meet the entered patient criteria, and allows the user to select the correct patient. <p>NOTES _____</p> <p> See Table 7 for alerts that may occur during this option.</p> <p>Not all possible error messages may display when multiple messages may apply to the unit/patient. Correction of any one of the patient/unit problems and subsequent re-attempts to issue the unit to the patient will display another of the possible messages until all errors are corrected.</p>
3. Click OK to verify that the correct patient was selected.	<ul style="list-style-type: none"> • Displays fields for the user to enter data for the issue transaction.
4. Enter the required data in the transporter name, issue-to location, and date and time fields (Figure 96). Click OK .	<p>NOTES _____</p> <p>The user may edit the transfusion-ordering clinician, which does not change the blood-ordering clinician.</p> <p>When the user updates the transfusion-ordering-clinician field,</p>

User Action	VBECS
	<p>VistALink queries VistA for valid providers only. The user may narrow the search by entering the provider's name.</p> <p>The user may edit the issue date and time (default: current date and time) to allow retrospective data entry; no future issue date or time is allowed. When the issue time is incorrect, the user must undo the issue and reissue the unit with the correct date and time.</p> <p>A remote storage indicator is optional (default: No). The transporter name field allows free text and does not display a selection list. The issue-to location lists current hospital locations provided by VistALink, including the inpatient or outpatient location indicator (VistA refers to this as a "patient class"). This is related to billing reimbursement and is not an enabled field in the issue option.</p> <p>A user could return the unit from issue and reissue it with a hospital location when VistALink is up again.</p> <p>When a remote storage location indicator is checked, the user may select from a list of previously entered locations or enter a new, free-text location. Each location entered and saved is available for future selection as a remote location.</p>
5. Scan the unit ID and product code, or click one or more check boxes to select units for issue to the selected patient (Figure 97).	<ul style="list-style-type: none"> • Displays the option to select units from the Assigned Units or Emergency Issue Units tabs. • Lists the selected units, allows the user to select additional units or deselect units, and prompts the user to indicate that the selections are complete. • Allows the user to select the unsatisfactory unit and edit its issue selections. <p>NOTES —————</p> <p> See Table 7 for alerts that may occur during this option.</p> <p>When the user attempts to issue a unit with an inconclusive crossmatch result, VBECS warns the user that the repeat crossmatch is required before the unit can be issued.</p> <p>When the user attempts to issue a unit that is not assigned or crossmatched to the patient, VBECS warns the user that the unit is not assigned to the patient and allows the user to enter another unit number.</p> <p>When required by the component, the unit must have a satisfactory ABO/Rh confirmation result before it is considered available for issue (emergency or routine) and the Caution Tag and BTRF may be generated.</p> <p>When a current specimen is required, VBECS compares the current and previous ABO/Rh results. It does not evaluate results marked "Test Result Invalidated."</p> <p>To issue blood components that require a crossmatch, the TAS must be completed and the unit status must be crossmatched to the patient.</p> <p>Whole blood is not allowed for the emergency issue process. Completed Type and Screen and compatible crossmatch tests are required for whole blood issuance.</p>

User Action	VBECS
6. Click OK to save and complete the unit issue.	<ul style="list-style-type: none"> • Displays entries for review and editing. <p>NOTES</p> <p>When the user selects multiple units for a patient and creates a batch for issue, all units in the batch are issued simultaneously to the same location. VBECS assigns the issue information to all units in the batch. The user may not enter different information, including any override comment, for an individual unit in the batch.</p>
7. Click Yes when the unit is satisfactory to save, or No when the unit is unsatisfactory.	<ul style="list-style-type: none"> • Saves the data. • Displays the unit inspection indicator (satisfactory or unsatisfactory) listing each selected unit. • Sends updated patient and blood unit information to the BCE COTS system via an HL7 message. • Sends updated patient and blood unit information to the BCE COTS system via an HL7 message. The message contains the date time when VBECS creates the message and not the date time the user selects the unit for issue. <p>NOTES</p> <p> See Table 7 for alerts that may occur during this option.</p> <p>Sites using BCE COTS must issue blood components in real-time for BCE COTS to be able to perform transfusion verification.</p>
8. Exit.	

Figure 96: Issue Blood Components (1)


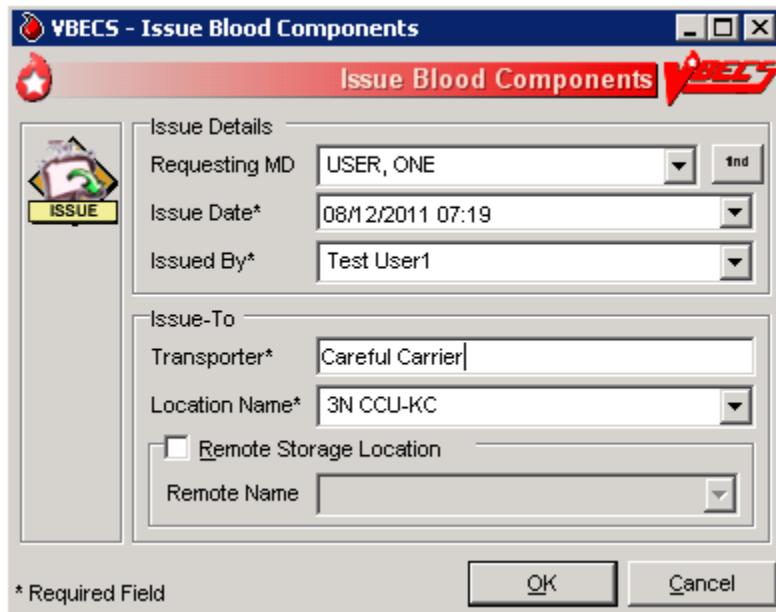


Figure 97: Issue Blood Components (2)

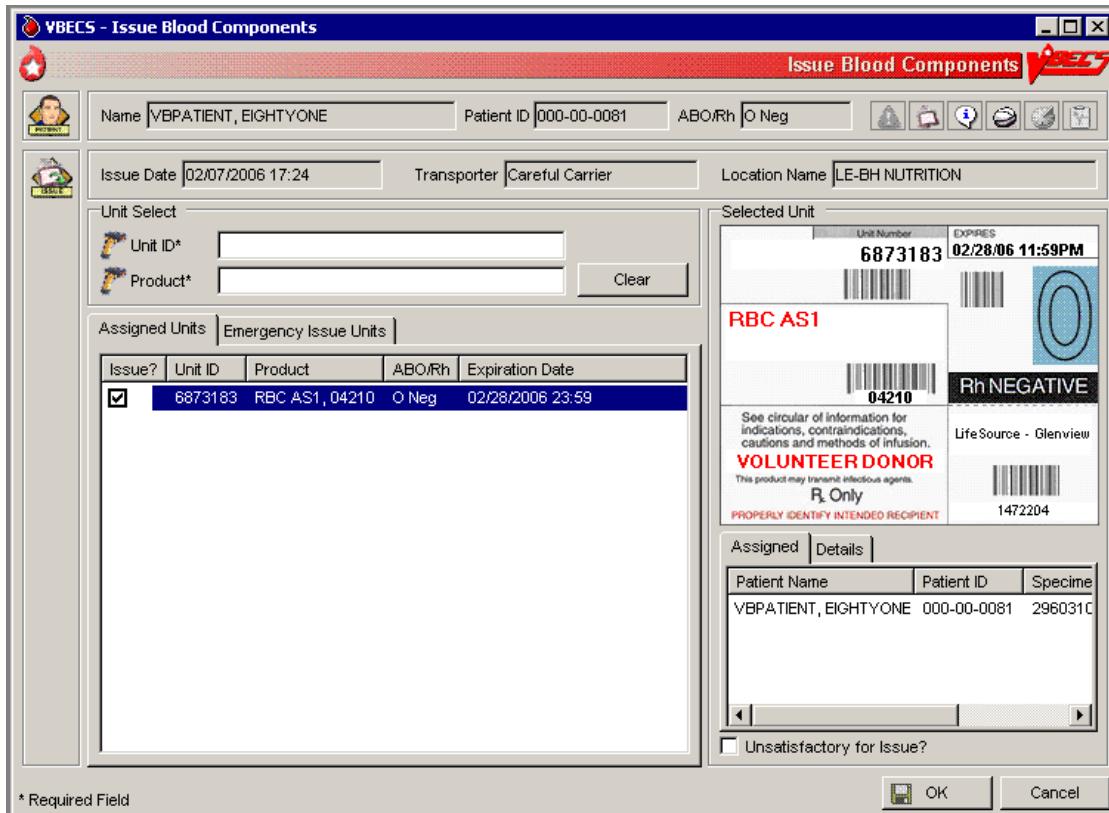


Issue Blood Components (Emergency)

The user issues assigned or crossmatched blood to a transporter for patient transfusion that does not meet routine issue conditions at valid inpatient and outpatient hospital locations, and records issue details, including the transporter's name or identifier.

User Action	VBECS
1. Select Patients from the main menu. Select Issue Blood Components .	<ul style="list-style-type: none"> Displays options for processing patient-related functions. Displays an option to scan the Caution Tag or enter a patient name or ID.
2. Select a patient.	<ul style="list-style-type: none"> Displays patient data, including entries that meet the entered patient criteria, and allows the user to select the correct patient. <p>NOTES</p> <p> See Table 8 for alerts that may occur during this option.</p>
3. Click OK to verify that the correct patient was selected.	<ul style="list-style-type: none"> Displays fields for the user to enter the data for the issue transaction.
4. Enter the required data.	<ul style="list-style-type: none"> Displays the option to select units from the Assigned Units or Emergency Issue Units tabs.
5. Click the Emergency Issue Units tab.	<ul style="list-style-type: none"> Lists the selected units, allows the user to select additional units or deselect units, and prompts the user to indicate that the selections are complete.

User Action	VBECS
Scan or click one or more check boxes to select units for issue to the selected patient (Figure 98).	<ul style="list-style-type: none"> Allows the user to select the unsatisfactory unit and edit its issue selections. <p>NOTES —</p> <p> See Table 8 for alerts that may occur during this option.</p>
Click OK .	
6. Click OK to accept the entries and complete the unit issue.	<ul style="list-style-type: none"> Displays entries for review and editing.
7. Click Yes when the unit is satisfactory or No when the unit is unsatisfactory.	<ul style="list-style-type: none"> Saves the data. Displays the unit inspection indicator (satisfactory or unsatisfactory) listing each selected unit. Sends updated patient and blood unit information to the BCE COTS system via an HL7 message. <p>NOTES —</p> <p> See Table 8 for alerts that may occur during this option.</p>
8. Exit.	

Figure 98: Issue Blood Components Emergency Issue Units Tab

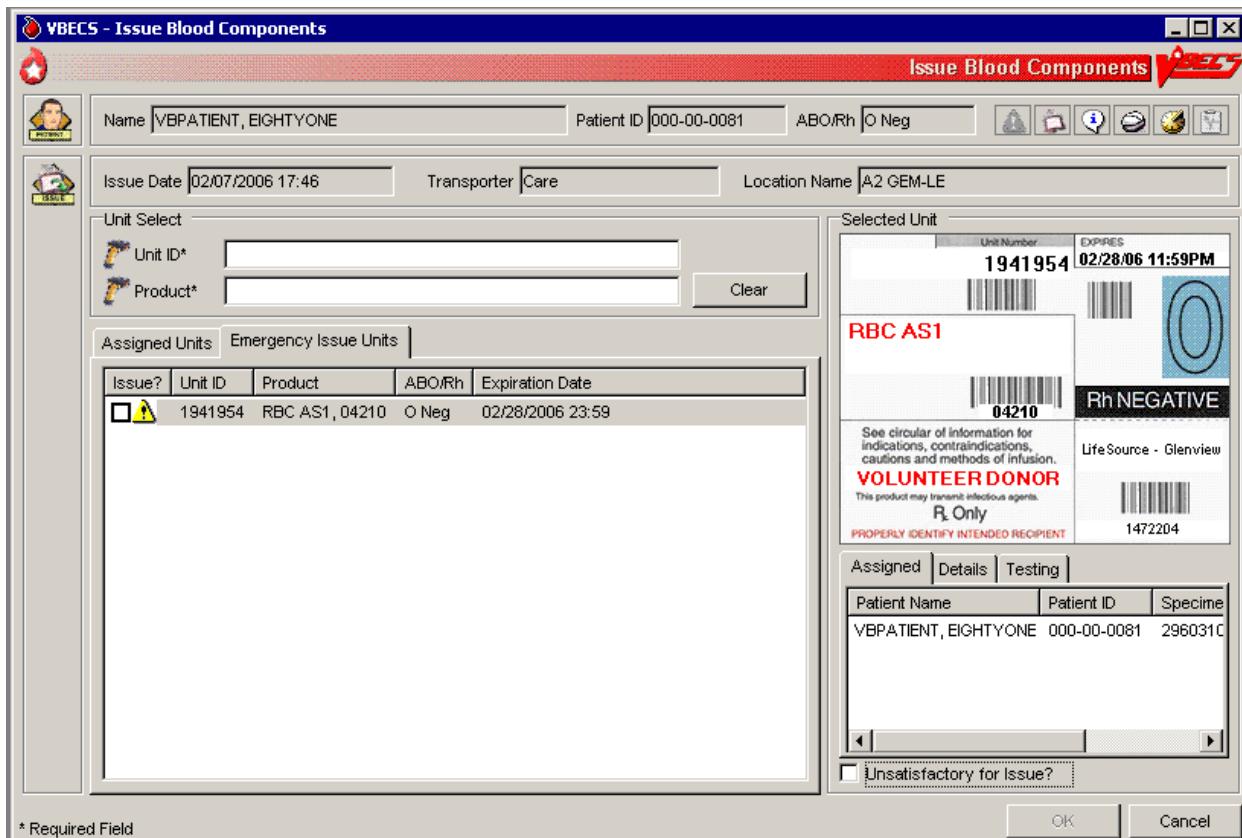

The screenshot shows the VBECS software interface for issuing blood components. The main title bar reads "VBECS - Issue Blood Components". The window is titled "Issue Blood Components".

On the left, there's a vertical toolbar with icons for patient selection, issue date, transporter, location, and other functions. The main area has several input fields and dropdown menus:


- Name: VBPATIENT, EIGHTYONE
- Patient ID: 000-00-0081
- ABO/Rh: O Neg
- Issue Date: 02/07/2006 17:46
- Transporter: Care
- Location Name: A2 GEM-LE
- Unit Select section: Unit ID* (1941954), Product* (RBC AS1)
- Assigned Units: Emergency Issue Units
- Table of Assigned Units:

Issue?	Unit ID	Product	ABO/Rh	Expiration Date
<input checked="" type="checkbox"/>	1941954	RBC AS1, 04210	O Neg	02/28/2006 23:59
- Selected Unit card (for unit 1941954):

Unit Number: 1941954	EXPIRES: 02/28/06 11:59PM
RBC AS1	
04210	Rh NEGATIVE
See circular of information for indications, contraindications, cautions and methods of infusion.	
VOLUNTEER DONOR	
This product may transmit infectious agents.	
R Only	
PROPERLY IDENTIFY INTENDED RECIPIENT	
- Details tab (Patient Name: VBPATIENT, EIGHTYONE, Patient ID: 000-00-0081, Specimen ID: 296031C)
- Testing tab (checkbox: Unatisfactory for Issue?)
- Buttons at the bottom: OK and Cancel

* Required Field

Figure 99: Emergency Issue

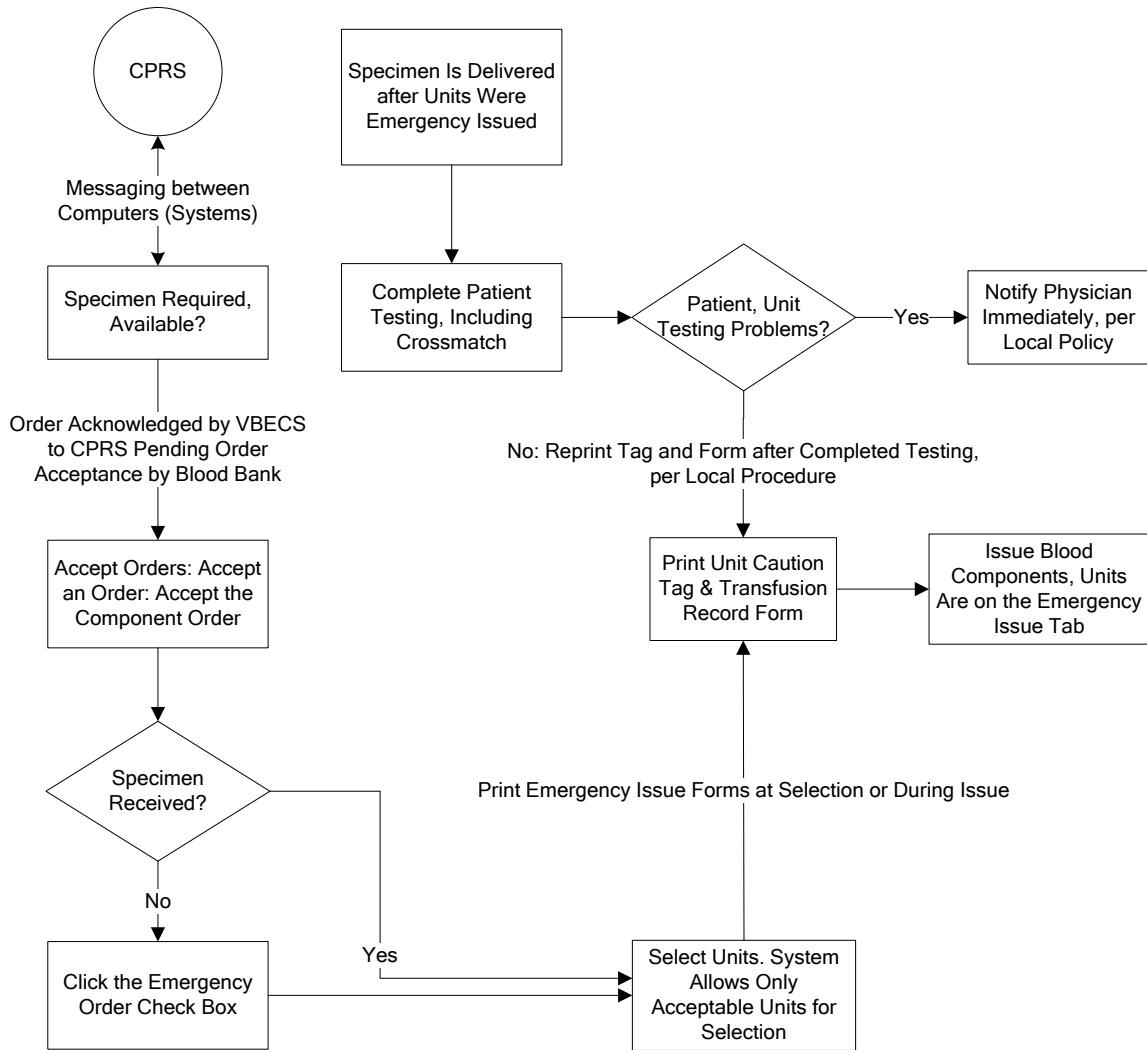

Alerts

Table 7: Alerts That May Occur in Issue Blood Components (Routine)

User Action Step	Alerts
2	VBECS displays  and emits an audible alert when there are Special Instructions (SIs) or Transfusion Requirements (TRs) for the patient.
2	When no patient is found, VBECS warns the user and does not allow the user to continue.
2	When the patient has no components available for issue, VBECS warns the user and does not allow the user to continue.
2	When required testing is incomplete at issue, VBECS emits an audible alert, warns the user, and asks whether to emergency issue the unit in accordance with hospital policy. No stops the user from issuing the unit and requires selection of another unit. Yes creates an Exception Report entry for the issue, indicates on the unit record that the unit was issued emergently, and asks whether the user wishes to print the Caution Tags and Blood Transfusion Record Forms. No allows the user to continue with the emergency issue. Yes prints a Caution Tag and/or Blood Transfusion Record Form and proceeds with the issue.
5	When the user attempts to issue a unit located in a different division, VBECS warns him that the unit cannot be issued and does not allow him to continue.
5	When autologous or directed, or designated or dedicated blood units of the same component class are available in the same division as the user, VBECS requires blood units to be issued in this order: 1. Autologous units 2. Restricted units 3. Voluntary allogeneic units When the user attempts to issue a unit out of order, VBECS warns him that more restrictive units are available for issue and asks whether he wishes to continue. No allows the user to enter another unit number. Yes allows the user to continue and captures transaction details for inclusion in an Exception Report (exception type: less restrictive unit issued).
5	When a user indicates that the inspected unit is unsatisfactory, VBECS does not allow the unit to be issued, releases it from assignment, and warns the user that the unit is unacceptable for transfusion. VBECS allows the user to enter another unit number.
5	When the user attempts to select frozen blood components, VBECS warns the user. Units selected in their frozen state must be thawed before issue.
5	When the patient has Component Requirements (CRs) and Antigen Negative Requirements, the user attempts to issue (including emergency issue) red blood cells that do not satisfy the requirements, and the user is: ►► An Enhanced Technologist or below, VBECS emits an audible alert and warns the user that the unit is not selectable and he is not authorized to issue the unit. The unsatisfied requirements are not detailed. VBECS does not allow an override.
5	When the patient has Component Requirements, the user attempts to issue (including emergency issue) red blood cells that do not satisfy the requirement, and the user is: • ►► An Enhanced Technologist or below, VBECS emits an audible alert and warns the user that the patient's component requirement is not satisfied and he is not authorized to issue the unit. VBECS does not allow an override. • ►►► A Lead Technologist or above, VBECS emits an audible alert and warns the user that the patient's component requirement is not satisfied and allows for an override. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: unit issued with unsatisfied TR).
5	When a repeat ABO/Rh performed on the patient specimen does not match the original ABO/Rh interpretation and the unit doesn't qualify for emergency issue, VBECS warns the user and instructs him to resolve the discrepancy before issuing the unit. There is no override. Units that qualify for emergency issue may continue and VBECS will display a message that the ABO/Rh on the current specimen does not agree with the historic ABO/Rh.

User Action Step	Alerts
5	When a unit is issued for a patient with a specimen expiration date more than three days in the future, VBECS changes the expiration date to three days from the date of issue of the unit.
5	When a unit has a biohazard indicator, VBECS emits an audible alert, warns the user, and asks whether he still wishes to issue the unit. No clears the entry and allows the user to select a new unit. Yes requires a comment and captures details of the transaction for inclusion in an Exception Report (exception type: biohazardous unit issued).
5	For an Rh negative, unknown, or inconclusive patient: the user may enter and confirm Rh positive or Rh pooled units from the WHOLE BLOOD, RED BLOOD CELLS, PLATELETS, and OTHER VBECS component classes. When VBECS component classes are individually entered and confirmed by the user, VBECS emits an audible alert, warns the user, and asks whether he wishes to issue the unit. No clears the unit data and allows the user to select a different unit. Yes requires additional comment and captures details for inclusion in an Exception Report (exception type: antigen positive/untested units issued).
5	When a user attempts to issue units assigned and/or crossmatched or that require an ABO/Rh associated with an expired specimen when compared to the issue date and time, VBECS warns the user. There is no override. ►► When the user is an Enhanced Technologist or above, VBECS emits an audible alert, warns the user, and asks whether he wishes to continue to issue. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: unit issued on expired specimen).
5	When the user attempts to issue an expired unit, VBECS warns him. There is no override. ►► When the user is an Enhanced Technologist or above, VBECS emits an audible alert, warns him, and asks whether he wishes to continue to issue. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: issued expired unit).
5	When the user selects units from the WHOLE BLOOD, RED BLOOD CELLS, or APHERESIS RED BLOOD CELLS ICCBBA component classes, VBECS determines when the patient has a calculated or a persistent antigen negative requirement. This does not apply when selecting autologous (ISBT 128 units with an autologous donation type). When the selected unit was antigen typed and is positive for the corresponding antigen or is untested for the antigen, VBECS warns the user and requires action based on his security level. ►► When the user is an Enhanced Technologist or below, VBECS warns that he does not have the proper security level to issue the unit. There is no override. ►►► When the user is a Lead Technologist or above, VBECS emits an audible alert, warns him, and asks whether he wishes to continue to issue. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report entry (exception type: antigen-positive/untested units issued).
5	VBECS compares the ABO/Rh of the current specimen to the historic ABO/Rh. When there is a discrepancy, there is no previous justification for the discrepancy on the patient record, as described in Justify ABO/Rh Change, and the unit does not qualify for emergency issue, VBECS warns the user and instructs him to resolve the discrepancy before issuing the unit. There is no override. VBECS exits the option, clears the screen, and returns the user to the point of entry. Units that qualify for emergency issue may continue and VBECS will display a message that the ABO/Rh on the current specimen does not agree with the historic ABO/Rh.
5	When no valid ABO exists for the patient, the user may select components, including autologous units: <ul style="list-style-type: none">• For RBC or WB order: group O, Rh positive or O Rh negative RBC. WHOLE BLOOD is not selectable.• For FFP order only: group AB, Rh positive or Rh negative.• For PLT, CRYO, or OTHER order, all available units are selectable.
5	When the division is “transfusion only” and a blood unit is ABO incompatible at issue based on the interpretation of the current patient specimen’s ABO test, VBECS warns the user that the unit is not ABO compatible for the patient’s current specimen and may not be issued. The user must resolve pending ABO discrepancies before issuing the unit.

User Action Step	Alerts
7	When the division is “full service” and a blood unit is ABO incompatible at issue based on the system interpretation of the current patient specimen’s ABO test result, VBECS emits an audible alert, warns that the unit is not ABO compatible with the patient’s current specimen results, and asks whether he wishes to continue. No clears the screen and allows the user to select another unit. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: ABO discrepant unit issue).

Table 8: Alerts That May Occur in Issue Blood Components (Emergency)

User Action Step	Alerts
2	VBECS displays  and emits an audible alert when there are Special Instructions (SIs) or Transfusion Requirements (TRs) for the patient.
2	When no patient is found, VBECS warns the user and does not allow the user to continue.
2	When the patient has no components available for issue, VBECS warns the user and does not allow the user to continue.
2	When required testing is incomplete at issue, VBECS emits an audible alert, warns the user, and asks whether to emergency issue the unit in accordance with hospital policy. No stops the user from issuing the unit and requires selection of another unit. Yes creates an Exception Report entry for the issue, indicates on the unit record that the unit was issued emergently, and asks whether the user wishes to print the Caution Tags and Blood Transfusion Record Forms. No allows the user to continue with the emergency issue. Yes prints a Caution Tag and/or Blood Transfusion Record Form and proceeds with the issue.
5	When the user attempts to issue a unit located in a different division, VBECS warns him that the unit cannot be issued and does not allow him to continue.
5	When a user indicates that the inspected unit is unsatisfactory, VBECS does not allow the unit to be issued, releases it from assignment, and warns the user that the unit is unacceptable for transfusion. VBECS allows the user to enter another unit number.
5	When the user attempts to select frozen blood components, VBECS emits an audible alert and warns the user. Units selected in their frozen state must be thawed before issue.
5	When the patient has Component Requirements and Antigen Negative Requirements, the user attempts to issue (including emergency issue) red blood cells that do not satisfy the requirements, and the user is: ►► An Enhanced Technologist or below, VBECS emits an audible alert and warns the user that the unit is not selectable and he is not authorized to issue the unit. The unsatisfied requirements are not detailed. VBECS does not allow an override.
5	When the patient has Component Requirements, the user attempts to issue (including emergency issue) red blood cells that do not satisfy the TR, and the user is: • ►► An Enhanced Technologist or below, VBECS warns the user that he is not authorized to issue untyped red blood cells (for the clinically significant antibodies in the table). VBECS does not allow an override. • ►►► A Lead Technologist or above, VBECS emits an audible alert and warns the user. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: unit issued with unsatisfied TR).
5	When a repeat ABO/Rh performed on the patient specimen does not match the original ABO/Rh interpretation and the unit does not qualify for emergency issue, VBECS warns the user and instructs him to resolve the discrepancy before issuing the unit. There is no override. Emergency issue may continue for qualifying units.
5	When a unit is issued for a patient with a specimen expiration date more than three days in the future, VBECS changes the expiration date to three days from the date of issue of the unit.
5	When a unit has a biohazard indicator, VBECS emits an audible alert, warns the user, and asks whether he still wishes to issue the unit. No clears the entry and allows the user to select a new unit. Yes requires a comment and captures details of the transaction for inclusion in an Exception Report (exception type: biohazardous unit issued).
5	For an Rh negative, unknown, or inconclusive patient, the user may enter and confirm Rh positive and

User Action Step	Alerts						
	Rh pooled units from the WHOLE BLOOD, RED BLOOD CELLS, PLATELETS, and OTHER VBECS component classes. VBECS emits an audible alert, warns that the patient is Rh negative and the selected unit is Rh positive, instructs the user to check the patient sex and hospital policy before continuing, and asks whether the user wishes to continue to select the unit. No clears the unit data and allows the user to select a different unit. Yes requires additional comment and captures details for inclusion in an Exception Report (exception type: antigen positive/untested units issued).						
5	<p>When the user attempts to issue an expired unit, VBECS emits an audible alert and warns him. There is no override.</p> <p>►► When the user is an Enhanced Technologist or above, VBECS emits an audible alert, warns him, and asks whether he wishes to continue to issue. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: issued expired unit).</p>						
5	<p>When the user selects units from the WHOLE BLOOD, RED BLOOD CELLS, or APHERESIS RED BLOOD CELLS ICCBBA component classes, VBECS determines when the patient has a calculated or a persistent antigen negative requirement. This does not apply when selecting autologous (ISBT 128 units with an autologous donation type). When the selected unit was antigen typed and is positive for the corresponding antigen or is untested for the antigen, VBECS warns the user and requires action based on his security level.</p> <p>►► When the user is an Enhanced Technologist or below, VBECS emits an audible alert and warns that he does not have the proper security level to issue the unit. There is no override.</p> <p>►►► When the user is a Lead Technologist or above, VBECS emits an audible alert, warns him, and asks whether he wishes to continue to issue. No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report entry (exception type: antigen-positive/untested units issued).</p>						
5	VBECS compares the ABO/Rh of the current specimen to the historic ABO/Rh. When the unit was electronically crossmatched, the current and historic ABO/Rh results must have been recorded in the login user's division. When there is a discrepancy, there is no previous justification for the discrepancy on the patient record, as described in Justify ABO/Rh Change, and the unit does not qualify for emergency issue, VBECS warns the user and instructs him to resolve the discrepancy before issuing the unit. There is no override. VBECS exits the option, clears the screen, and returns the user to the point of entry. Emergency issue may continue for qualifying units.						
5	<p>When no valid ABO exists for the patient, the user may select components, including autologous units:</p> <ul style="list-style-type: none"> • For RBC or WB order: group O, Rh positive or O Rh negative RBC. WHOLE BLOOD is not selectable. • For FFP order only: group AB, Rh positive or Rh negative. • For PLT, CRYO, or OTHER order, all available units are selectable. 						
5	<input checked="" type="checkbox"/> When the user selects blood components from the frozen product types for a patient, VBECS emits an audible alert and warns that the selected blood unit must be thawed and may require additional processing prior to issue. The unit is not eligible for printing the BTRF or Caution Tag. VBECS updates the blood unit status to "assigned" and the unit is no longer available for issue.						
5	<p>When a Caution Tag or BTRF is printed for an emergency-issued unit, "Emergency Issue" is included on the form.</p> <p>The recipient ABO/Rh prints as "Unk" or "Unknown" during emergency issue when no current specimen results are on file; VBECS may not use a historic record blood type (from any source) for the typing printed on the tag/form.</p> <table border="1"> <thead> <tr> <th>Compatibility Interpretation</th><th>Text on BTRF</th></tr> </thead> <tbody> <tr> <td>Crossmatch</td><td>TESTING INCOMPLETE</td></tr> <tr> <td>Crossmatch Not Required</td><td>Not Required</td></tr> </tbody> </table> <p>VBECS prints "Emergency Issue" in the TR section of the form.</p>	Compatibility Interpretation	Text on BTRF	Crossmatch	TESTING INCOMPLETE	Crossmatch Not Required	Not Required
Compatibility Interpretation	Text on BTRF						
Crossmatch	TESTING INCOMPLETE						
Crossmatch Not Required	Not Required						
5	When blood products are emergency issued (units has an emergency issue" order indicator, or the						

User Action Step	Alerts
	<p>testing is incomplete or inconclusive on the current specimen), the remarks section of the Blood Transfusion Record Form includes “This unit is issued for transfusion BEFORE the completion of required compatibility testing. The blood bank staff will continue to test the compatibility of the unit using the current patient specimen. If you are alerted to a unit’s incompatibility by the blood bank, stop the transfusion immediately.”</p> <p>When “****” prints for the patient ABO/Rh and crossmatch, the form includes “Specimen testing incomplete at time of issue.”</p>
5	<p>When a unit was emergency issued and a conflicting test result found that indicates a potential unit incompatibility with the patient, such as a patient ABO/Rh discrepancy, the unit XM is incompatible with the patient, the patient ABS is positive, the patient DAT is positive, or the unit is antigen positive where antigen negative is required on the tested specimen, VBECS warns the user and instructs him to notify the physician of this potential problem immediately, according to hospital policy. VBECS warns the user during the execution of the various tests related to the compatibility of a blood unit and captures details for inclusion in an Exception Report (exception type: unit ER issued, testing problem).</p>
7	<p>When the division is “full service” and a blood unit is ABO incompatible at issue based on the system interpretation of the current patient specimen’s ABO test result, VBECS emits an audible alert and warns that the unit is not ABO compatible with the patient’s current specimen results and asks whether he wishes to continue. No clears the screen and allows the user to select another unit. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: ABO discrepant unit issue).</p>

Return Issued Units To Blood Bank

An issued unit was returned to the blood bank in satisfactory condition after issue to a valid hospital location.

Assumptions

- The user is logged into the division that contains the unit.
- The blood unit is in the “issued” status.
- The connection to VistA is active.

Outcome

- The unit is returned to a hospital location of a blood bank.
- The unit retains its patient assignment at the time of issue when returned in satisfactory condition.
- The user may mark the unit “unsatisfactory” and quarantine the unit to generate an Exception Report, if needed per hospital policy.

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option

All users

Return Issued Units To Blood Bank

An issued unit was returned to the blood bank in satisfactory condition after issue to a valid hospital location. The user updates VBECS to reflect the change in location.

User Action	VBECS
1. Select Blood Units or Patients from the main menu. Select Return Issued Units To Blood Bank .	<ul style="list-style-type: none"> Displays options for processing blood units or for processing patient-related functions. Displays the return date and time and the name of the user who returned the issued unit. <p>NOTES</p> <p>The user may edit the return date and time (default: current date and time) to a past date and time.</p> <p>VBECS assigns the same return date and time to all units processed in a batch.</p>
2. Enter or edit the date and time and/or select a user, or click OK to accept the default data.	
3. Enter the patient data.	<ul style="list-style-type: none"> Prompts the user to enter the unit ID and product code or enter patient data to retrieve the records for the units issued. Displays an option to search for or enter a patient or unit based on the original menu option selected.
4. Select a patient or unit.	<ul style="list-style-type: none"> Lists the units in "issued" status that are assigned to the patient.
5. Select one or more units to return from issue.	<ul style="list-style-type: none"> Displays the units selected for return for review and allows the user to select or deselect additional units.
6. Scan the unit ID and product code, or click one or more check boxes to select units to return from issue.	<ul style="list-style-type: none"> Displays the data entered and allows the user to enter additional units and indicate that he completed the entry for this transaction. Displays data for each unit entered. <p>NOTES</p> <p>When a selected unit is not in "issued" status, VBECS warns the user, clears the unit field, and allows the user to enter a different unit number.</p> <p>VBECS allows the user to select multiple units to create a worklist for simultaneous processing.</p> <p>VBECS requires the user to enter a "returned-from" location (default: "issue-to" location) and allows the user to change this location to a different valid hospital location. To do this, click the find button and select a hospital location.</p>
7. Respond to the warnings and messages (Figure 100).	<ul style="list-style-type: none"> Allows the user to indicate that each unit or batch of units was visually inspected and marked satisfactory or unsatisfactory. <p>NOTES</p> <p> When the user indicates that the unit is unsatisfactory, VBECS emits an audible alert and quarantines the unit, relocates it to the hospital location of the blood bank, and captures details for inclusion in an Exception Report (exception type: unit unsatisfactory upon return from issue). Upon quarantining the unit, VBECS releases it from assignment and requires a comment regarding the return.</p>

User Action	VBECS
	<p>When a unit has had a quarantine indicator assigned, VBECS releases any patient assignments. Patient restrictions associated with the unit remain intact.</p> <p>VBECS places units returned from issue to the blood bank in satisfactory condition in the status they had prior to issue.</p> <p>When the unit returned and/or the specimen used for assignment or crossmatch is expired, VBECS warns the user.</p>
8. Click OK to accept the entries, and complete the unit return. Exit.	<ul style="list-style-type: none"> • Saves the data. • Sends update to the BCE COTS system via an HL7 message to remove the unit from the system. <p>NOTES</p> <p>VBECS returns units to the blood bank location within the division.</p>

Figure 100: Return Issued Units

The screenshot shows the 'VBECS - Return Issued Units' window. At the top, there are fields for Name (VBPATIENT, EIGHTYONE), Patient ID (000-00-0081), and ABO/Rh (O Neg). Below these are buttons for Home, Logout, and Help. The main area is divided into sections:

- Select Units to Return:** A table with columns: Return?, Unit ID, Product, ABO/Rh, and Expires. One row is selected, showing: 6873183 RBC AS1, 04210 O Neg 2/28/2006 11:59:00 PM.
- Selected Unit:** Details for the selected unit. Fields include: Donation Type (Volunteer Allogeneic Donor), Date Received (01/31/2006 10:05), Unit Status (Issued), and Volume (mL) (350).
- Unit Issue Details:** Fields include: Issue Date (2/7/2006 5:24:56 PM), Transporter (Careful Carrier), Issued To (LE-BH NUTRITION), and Remote Storage Location (Identity).
- Unit Return Details:** Fields include: Returned From* (LE-BH NUTRITION) and a checkbox for Unit Unsatisfactory for Return.

At the bottom left is a note: * Required Field. At the bottom right are OK and Cancel buttons.

Release Units From Patient Assignment

The user releases crossmatched and assigned blood units not issued after a period determined by regulatory standards and by internal policy.

Assumptions

- The blood unit is assigned or crossmatched to at least one patient.

Outcome

- Patient restrictions still apply to a unit even after it is released from assignment and made available.

Limitations and Restrictions

- This option does not apply to the release of patient restrictions (autologous and directed).
- When a patient, specimen ,or blood unit expires, VBECS does not automatically release an assigned blood unit: a user must release the unit.

Additional Information

- A user may release units from assignment as a group or individually from a patient's list of units.

User Roles with Access to This Option

All users


Release Units From Patient Assignment

A user may select one or more blood units for release from patient assignment and return to “available” status. A user may release one of multiple patient assignments without affecting other patient assignments on a unit.

User Action	VBECS
<p>1. To release units by patient, select Patients from the main menu.</p> <p>To release units by unit, select Blood Units from the main menu.</p> <p>Select Release Units From Patient Assignment.</p>	<ul style="list-style-type: none">• Displays options for processing patient-related functions or for processing units.
<p>2. Scan or enter a unit ID or select a patient to release units from patient assignment and/or crossmatch (XM).</p>	<ul style="list-style-type: none">• Allows the user to edit the date and time and the user who released the units.• Lists currently assigned units.• Allows the user to list assigned units by patient, by unit ID, or by unit expiration date. <p>NOTES —</p> <p>VBECS indicates which specimens and units are expired.</p> <p>VBECS displays a unit's donation type only when it is restricted.</p> <p>VBECS displays all currently assigned/crossmatched units located in the blood bank. When the user enters a unit ID that is not currently in the blood bank location, VBECS warns the user, clears the unit ID entry, and allows the user to enter another unit</p>

User Action	VBECS
3. Indicate which units to release. Click OK . Indicate a reason for release of units (Figure 101).	<p>or patient.</p> <ul style="list-style-type: none"> Allows the user to enter a reason for releasing for an individual unit. Allows the user to review the selections, select additional units, and deselect units before saving. <p>NOTES</p> <p>VBECS allows the user to select no units or one or more units for each patient on the list.</p>
4. Review the entered information and click OK to save. Click OK to exit.	<ul style="list-style-type: none"> Updates the patient record, the unit record, and the status. Sends update to the BCE COTS system via an HL7 message to remove the unit from the system. <p>NOTES</p> <p>When a unit has a future expiration date and holds no other assignments, VBECS updates the unit to "available" status. When a unit is assigned to a different patient, it remains assigned to that patient.</p> <p>VBECS alerts the user when the unit is expired and reminds him to remove additional patient assignments and the unit from the available inventory according to local policies. VBECS instructs the user to release the unit from all assignments and remove it from the available inventory per hospital policy.</p>

Figure 101: Confirm Release

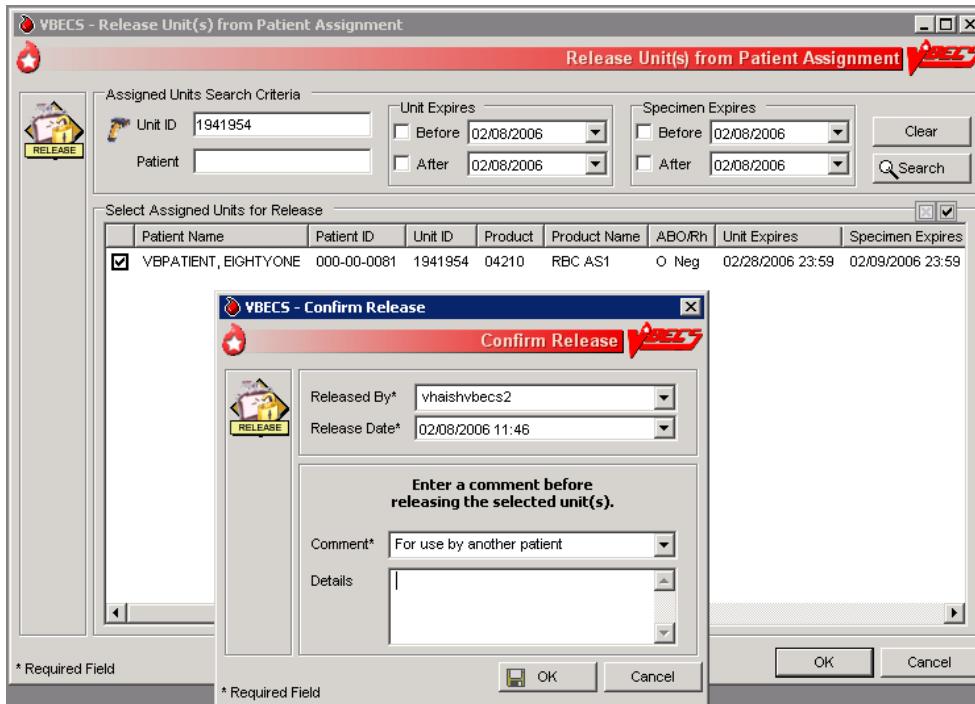

Print Unit Caution Tag & Transfusion Record Form

The user may print a Caution Tag and/or a Blood Transfusion Record Form (BTRF) when crossmatch compatibility, assignment, or emergency issue are completed or as an independent action.

Assumptions

- The blood unit, patient order, and user are in the same division.
- There is a current component class or emergency issue order for the patient, and the blood unit printed on the Caution Tag fulfills the patient order requirements.
- When the blood unit belongs to a component class that does not require crossmatch, it has a status of “assigned” (non-emergency issue only).
- For non-emergency issue, the blood unit meets all patient Transfusion Requirements (TRs).

Outcome

- The barcode generated on the Caution Tag may be used in bedside verification procedures, in accordance with local policy and procedures.
- A user may reprint a BTRF and Caution Tag.

Limitations and Restrictions

- VBECS cannot reprint a Caution Tag on an expired unit or expired patient order.
- VBECS does not print a unit’s biohazardous status on the BTRF or Caution Tag.

Additional Information

- The informational barcode printed on the Caution Tag includes only the patient and unit IDs and the unit product code.
- A user may print multiple Caution Tags and BTRFs for a blood unit when it is assigned to multiple patients.
- Attach the Caution Tag labels to blood bags in accordance with local policy and procedures.

User Roles with Access to This Option

All users


Print Backup Caution Tags and Blood Transfusion Record Forms

It is recommended that blank Caution Tags and BTRFs be kept on hand in case the system or printer is not available to generate them.

User Action	VBECS
<p>1. Select Tools from the main menu.</p> <p>Select Print Backup Forms.</p> <p>Enter or select the number of blank tags or forms to print (Figure 102).</p> <p>Click OK to print and exit.</p>	<ul style="list-style-type: none">• Displays options for processing administrative functions.• Prompts the user to enter the number of Caution Tags to be printed.• Prints the blank Caution Tags. <p>NOTES —————</p> <p>VBECS displays an option to print a Caution Tag or a BTRF when the user selects a unit or a patient.</p> <p>The user may print 1 to 99 blank Caution Tags and BTRFs to use in case of downtime.</p>

User Action	VBECS
	VBECS interrupts or cancels a print job to accommodate the printing of an emergency issue tag. VBECS notifies the user of the number of tags printed and instructs the user to enter a new request to print more tags.

Figure 102: Print Backup Forms

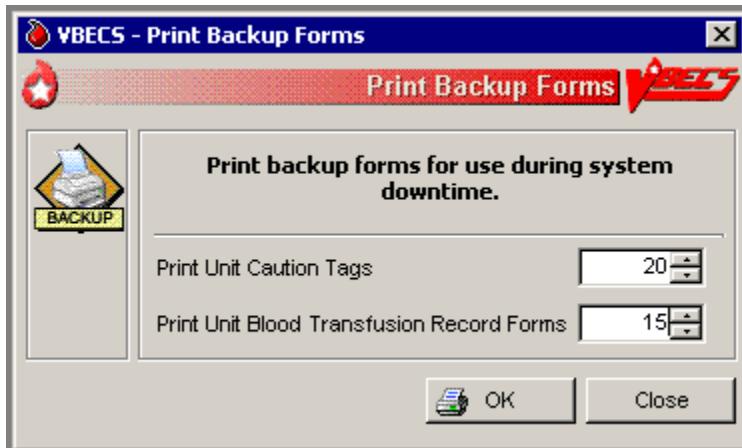

Print Unit Caution Tag & Transfusion Record Form

The user may need to generate Caution Tags when they were not printed immediately after unit selection and crossmatch data entry, emergency issue, or when there is a problem with the printed Caution Tag. According to local policy, the user selects a current patient order and one or more units to print from the VBECS-generated list of units eligible for immediate transfusion.


VBECS does not allow a user to print a Caution Tag unless the unit is satisfactory for immediate issue.

User Action	VBECS
1. Select Patients or Blood Units from the main menu. Select Print Unit Caution Tag & Transfusion Record Form .	<ul style="list-style-type: none"> Displays options for processing patient-related functions or for processing units. Displays the option to select a patient, if one is not already selected.
2. Enter the patient name or ID and click Search to select a patient, or Click Search by Unit (Figure 103). Enter the unit ID and/or product.	<ul style="list-style-type: none"> Requests that the user select a Caution Tag and/or a BTRF to print by unit ID or by patient. Displays an informational message, when appropriate.
3. Select a patient, or Select or scan a blood unit.	<ul style="list-style-type: none"> Compiles and displays the information associated with the user's selection. Displays the patient's name, ID, date of birth, the unit ID, and the long product name of units to appear on the Caution Tags. Prompts the user to select the units when selecting by patient. Allows the user to select and deselect blood units to print. See Table 9 for alerts that may occur during this option. <p>NOTES —</p> <p>VBECS retrieves the most recent crossmatch result (original or</p>

User Action	VBECS						
	<p>repeat) to determine unit availability and compatibility status to allow printing the Caution Tag and/or BTRF and to allow the unit to be issued.</p> <p>When no valid ABO exists for the patient, the user may select components, including autologous units:</p> <ul style="list-style-type: none"> • For RBC order only: Group O, Rh positive or O Rh negative. • For FFP order only: Group AB, Rh positive or Rh negative, or not specified. • For PLT, or Cryo order: All available units are issuable. • For OTHER order: All available units are issuable when the optional crossmatch has not been selected. When the optional crossmatch has been selected, OTHER are handled as RBC. • WHOLE BLOOD IS NOT SELECTABLE. <p>When a user selects a Caution Tag and/or BTRF from the Blood Unit Menu and the units are assigned or have been assigned and/or crossmatched to multiple patients, the selected form may be printed for only one patient during that session. The message will display that the form cannot be printed again for another assigned patient. VBECS must be closed and re-opened to print a tag for one of the other patients. If you need to print tags on the same unit for multiple assigned patients, access the option from the Patient Menu.</p> <p>When VBECS searches for units available for a patient name, it makes available for printing only units considered available for issue (standard processing or emergency issue).</p> <p>See Appendix B: Table 23: Rules for Electronic and Serologic Crossmatch. The user may select only active units and those not in an inactive or final status. See Appendix B: Figure 136: Unit Status Flowchart.</p>						
<p>4. Acknowledge the selected patient and blood units for printing.</p> <p>Click Print Unit Caution Tag, Print Transfusion Record Form, or Print Both.</p> <p>Click OK to confirm and exit.</p>	<ul style="list-style-type: none"> • Displays an option to print the Caution Tag and BTRF for units with a patient assignment. • Generates a separate Caution Tag and BTRF for each blood unit. <hr/> <p>NOTES</p> <p>When a Caution Tag or BTRF is printed for an emergency-issued unit, "Emergency Issue" is included on the form.</p> <p>The recipient ABO/Rh prints as "Unk" or "Unknown" during emergency issue when no current specimen results are on file; VBECS may not use a historic record blood type (from any source) for the typing printed on the tag/form.</p> <table border="1" style="margin-left: auto; margin-right: auto;"> <thead> <tr> <th style="text-align: center;">Compatibility Interpretation</th> <th style="text-align: center;">Text on BTRF</th> </tr> </thead> <tbody> <tr> <td style="text-align: center;">Crossmatch</td> <td style="text-align: center;">TESTING INCOMPLETE</td> </tr> <tr> <td style="text-align: center;">Crossmatch Not Required</td> <td style="text-align: center;">Not Required</td> </tr> </tbody> </table> <p>VBECS prints "Emergency Issue" in the TR section of the form.</p> <p>The text that appears in the compatibility section is based on the crossmatch interpretation applied.</p>	Compatibility Interpretation	Text on BTRF	Crossmatch	TESTING INCOMPLETE	Crossmatch Not Required	Not Required
Compatibility Interpretation	Text on BTRF						
Crossmatch	TESTING INCOMPLETE						
Crossmatch Not Required	Not Required						

User Action	VBECS										
	<table border="1"> <thead> <tr> <th>Compatibility Interpretation</th> <th>Text on BTRF</th> </tr> </thead> <tbody> <tr> <td>Crossmatch Compatible</td> <td>Compatible</td> </tr> <tr> <td>Compatible: Electronically Crossmatched</td> <td>Compatible</td> </tr> <tr> <td>Crossmatch Not Required</td> <td>Not Applicable</td> </tr> <tr> <td>Crossmatch Incompatible: Give only with Medical Director Approval</td> <td>Incompatible: Given with Medical Director Approval</td> </tr> </tbody> </table>	Compatibility Interpretation	Text on BTRF	Crossmatch Compatible	Compatible	Compatible: Electronically Crossmatched	Compatible	Crossmatch Not Required	Not Applicable	Crossmatch Incompatible: Give only with Medical Director Approval	Incompatible: Given with Medical Director Approval
Compatibility Interpretation	Text on BTRF										
Crossmatch Compatible	Compatible										
Compatible: Electronically Crossmatched	Compatible										
Crossmatch Not Required	Not Applicable										
Crossmatch Incompatible: Give only with Medical Director Approval	Incompatible: Given with Medical Director Approval										
5. Print and exit.											

Figure 103: Print Unit Caution Tag and Transfusion Record Form

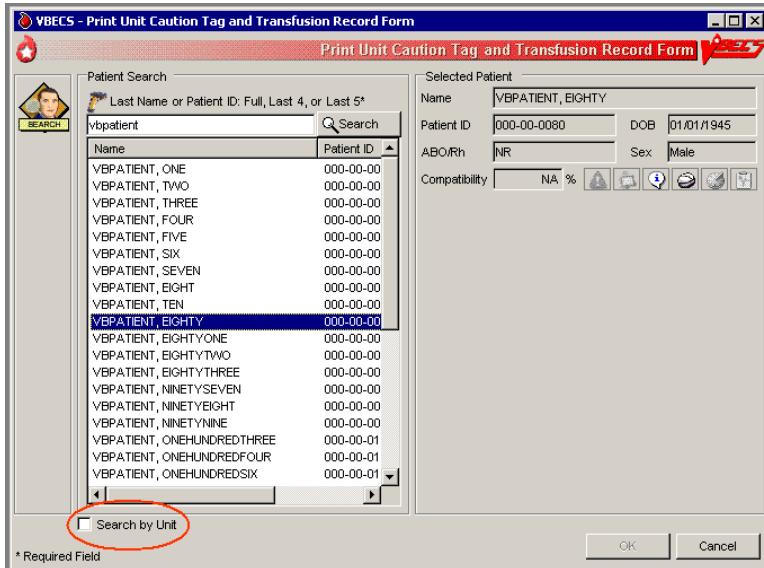

 **Alerts**
Table 9: Alerts That May Occur in Print Unit Caution Tag & Transfusion Record Form (Step 3)

User Action Step 3: When This Condition Exists  VBECS Emits an Audible Alert and Warns the User:	User Role
VBECS determines whether the patient has a calculated or persistent antigen negative requirement (except when selecting autologous units). The selected unit was antigen typed and is positive for the corresponding antigen or is untested for the antigen.	The patient has a history of a clinically significant antibody, the unit is positive or untested for the corresponding antigen, the user does not have the proper security level. There is no override.	Enhanced Technologist or below
	The patient has a history of an antibody(s) and the unit ID is positive or untested for the corresponding antigen(s). Continue to print Caution Tag or BTRF? No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: Caution Tag or BTRF printed for antigen positive or untested unit).	Lead Technologist
When the patient has TRs, VBECS evaluates the selected unit to see whether a TR was satisfied based on system rules, including emergency issue.	The patient has at least one TR that is not satisfied, the user does not have the proper security level to print the Caution Tag or BTRF for this unit. Further processing is needed and must be documented in VBECS. There is no override.	Enhanced Technologist or below
	The patient has at least one TR that is not satisfied. Continue to print Caution Tag or BTRF? No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: Caution Tag printed with unsatisfied TR).	Lead Technologist
The patient specimen is expired.	The patient specimen is expired. There is no override.	Technologist
	The patient specimen is expired. Continue to print the Caution Tag and/or the BTRF? No clears the screen. Yes requires a comment and captures details for inclusion in an Exception Report (exception type: Caution Tag printed on expired specimen).	Enhanced Technologist
<p>The user attempts to print Caution Tags for emergency issue and:</p> <ul style="list-style-type: none"> • The current specimen has an "Emergency Issue" order. • The patient has no record of ABO/Rh typing. • The required ABO/Rh specimen testing is incomplete. <p>The user may select:</p> <ul style="list-style-type: none"> • For WB or RBC order only: group O, Rh positive or O Rh negative (WHOLE BLOOD is not available). • For FFP order only: group AB, Rh positive or Rh negative. • For PLT order: all available units are issuable. 	The unit does not meet these requirements. Testing is incomplete and this unit is not eligible for emergency issue. The Caution Tag cannot be printed. There is no override.	All users
VBECS determines whether a blood unit that		
All of the conditions for the unit are met but		All users

User Action Step 3: When This Condition Exists  VBECS Emits an Audible Alert and Warns the User:	User Role
<p>does require crossmatch (electronic or serologic) is ready for issue based on the following rules:</p> <ul style="list-style-type: none"> • Unit is not in a frozen state. • The unit is crossmatched to the patient (eXM, or serologic XM). • TAS Specimen testing is complete and valid (ABO/Rh and ABS results are entered). <p>One of these crossmatch interpretations was entered:</p> <ul style="list-style-type: none"> • Crossmatch Compatible • Compatible: Electronically crossmatched • Crossmatch Incompatible: Give only with Medical Director approval 	<p>an emergency issue order indicator is in place or the testing is incomplete or inconclusive on the current specimen and units were selected. Do you want to print emergency issue tags/forms for the unit/patient? Yes prints the emergency issue tags and forms. No does not allow the user to continue.</p>	
For blood units that do not require crossmatch when the unit is assigned to the patient and is not in a frozen state, VBECS determines when a blood unit is ready for issue.	All of the conditions for the unit are not met. The Caution Tag and/or the BTRF cannot be printed at this time and the reason they cannot be printed. Additional patient testing or unit processing must be performed.	All users

Figure 104: Caution Tag

Product Code	2D Barcode
Division Name City, State	
Patient Name:	
Patient ID:	Recipient Blood Type:
Unit/Pool ID:	Donor Blood Type:
Compatibility Interpretation:	
CAUTION: Before transfusing this blood component, perform bedside verification of the patient's identity by matching the patient identification on this tag with the patient wristband and the Blood Transfusion Record Form.	
Assigned Date/Time:	Assigning Tech/ Crossmatch Tech:

Figure 105: Blood Transfusion Record Form

Blood Transfusion Record Form		
VA HEARTLAND - WEST, VISN 15, Kansas City, MO		Transfusion Requirements
Patient Name: VBECSPATIENTC, TWOHUNTWEIGHT		
Patient ID: 666-00-0228	Patient Blood Type: AB Neg	Blood Component Information
Unit/Pool ID: W3416 00 000020 Unit Blood Type: AB Neg		Component: Washed RED BLOOD CELLS None/XX/refg Open Irradiated Plasma added
Donation: Volunteer Allogeneic Donor		Product code: E5654V00
Compatibility Interpretation: TESTING INCOMPLETE		Component Expiration Date/Time: 04/29/2014 16:56
Assigned Date/Time: 04/29/2014 09:47		Number of units in pool: 0
		Additional Patient Information
		Specimen UID, Expiration: 2941190224, 05/02/2014 16:21
		Location: VBEC KANSAS CITY INPT
Remarks		
<small>*This unit is issued for transfusion BEFORE the completion of required compatibility testing. The blood bank staff will continue to test the compatibility of the unit using the current patient specimen. If you are alerted to a unit's incompatibility by the blood bank, stop the transfusion immediately.</small>		
EMERGENCY ISSUE		
Pretransfusion Data		
Inspected and Issued by:	PRINT	Sign:
Issued to:	PRINT	Date and time issued: Issued to location:
<small>I have verified and compared, AT THE BEDSIDE, the transfusion recipient's identity (i.e., wristband), the unit ID tag, the blood component container label, and this form. I verify that all information matches and is consistent ITEM for ITEM. I verify that the intended recipient is the same person named on this form and on the unit ID tag. Furthermore, I verify that there is a current, valid INFORMED CONSENT, and a PHYSICIAN'S ORDER for this transfusion.</small>		
1st identifier (Sign and Date):	2nd identifier (Sign and Date):	

Transfusion Data					
	Date/Time	Temperature	Pulse	Blood Pressure	Respiration
Start					
15 minutes					
Mid (optional)					
Stop					
Amount given (mL):	<input type="text"/>	<input type="checkbox"/> Completed <input type="checkbox"/> Interrupted <input type="checkbox"/> No reaction <input type="checkbox"/> Reaction suspected (Complete Transfusion Reaction Data section)			
Other difficulties (equipment, clots, etc.):	<input type="checkbox"/> No <input type="checkbox"/> Yes (specify):	Transfusion Data section completed by (Sign and Date):			
Transfusion Reaction Data					
Patient identity match with unit compatibility tag and this form reverified by (Sign and Date):					
Description of reaction: <input type="checkbox"/> Urticaria <input type="checkbox"/> Chill <input type="checkbox"/> Pain <input type="checkbox"/> Fever <input type="checkbox"/> Dyspnea <input type="checkbox"/> Other (specify):					
Transfusion Reaction Data section completed by (Sign and Date):					
Suspected Transfusion Reaction: Immediate Response					
If a reaction is suspected, IMMEDIATELY: 1. Interrupt the transfusion. Give emergency treatment. Keep the intravenous line open. 2. Reverify patient and blood product identification. 3. Notify physician and Transfusion Service. 4. Follow transfusion reaction procedures. 5. Do NOT discard the component. Return the blood component bag, filter set, and attached solutions to Transfusion Service. Collect patient specimens as required. 6. Order a Transfusion Reaction Workup.					

Patient Name: VBECS PATIENTC, TWOHUNTWEIGHT Patient ID: 666-00-0228 Patient DOB: 2/3/1980 BTRF, Version. 3.0

Invalidate Test Results

The user invalidates incorrectly entered patient test results and reenters correct patient test results and interpretations.

Assumptions

- A completed patient order exists in VBECS.

Outcome

- Test results are invalidated.
- VBECS may assign “correction” status to a task to allow the user to repeat a test and add it to the Pending Task List (PTL).

Limitations and Restrictions

- A user may not deactivate a component task in this option.
- This option does not release units from patient assignment. The user must release “assigned” units before invalidating the test results.
- To maintain a clear record and system checks, the user must invalidate both original and repeat ABO/Rh tests.
- Do not invalidate part of the configured test. Doing this will put the specimen into a state where it cannot be processed. Invalidate all tests configured as part of the TAS regardless of which sub test (ABO/Rh, ABS or DAT) is the problem.

Additional Information

- VBECS invalidates the entire test group during this option. The user must reenter results to replace the invalidated results.
- Workload and billing are inactivated or rescinded from the inactivated testing.
- VBECS invalidates transfusion reaction workups (TRWs) in total and crossmatches (XMs) by unit.
- Invalidation of an antibody identification (ABID) task invalidates all antibody specificities entered in that task..
- VBECS displays results and interpretations of the tests but not any QC performed with the tests.
- VBECS displays “Test Result Invalidated” to indicate invalidated test results.

User Roles with Access to This Option

➤➤ Enhanced Technologist

Invalidate Test Results

The user may invalidate one or more previously entered test results and replace them with corrected results.

User Action	VBECS
1. Select Patients from the main menu. Select Invalidate Test Results .	<ul style="list-style-type: none">• Displays options for processing patient-related functions. <p>NOTES —————</p> <p>Click the Component Orders tab to view its PTL.</p>
2. Enter the patient name and/or specimen UID.	<ul style="list-style-type: none">• Displays an option to enter a specimen UID or a patient for a specific date or for a range of dates.

User Action	VBECS
<p>Edit the "expires after" date, if required.</p> <p>Click Search, if needed.</p>	<ul style="list-style-type: none"> • Displays the complete diagnostic tests associated with the entry.
<p>3. Select a test result from the Specimen Test Results list.</p>	<ul style="list-style-type: none"> • Displays the original results for the user to review. <p>NOTES</p> <p>When the user chooses to invalidate an ABID, VBECS warns that the invalidation of the antibody specificity also invalidates the corresponding calculated antigen negative requirements. No user action is required.</p> <p>►►► When the user is a Traditional Supervisor and the selected specimen is the most recent specimen tested for the selected patient, the user may invalidate a previously justified ABO/Rh test result. When the specimen ABO/Rh test associated with a justification is invalidated, VBECS removes the justification information from the patient record.</p> <p>When the user attempts to invalidate a non-XM test for a patient specimen and blood units are currently associated with the specimen (assigned, electronically or serologically crossmatched, issued or transfused), VBECS warns that blood components are associated to this specimen and that the associated units must be removed before testing can be invalidated. The user may not continue to invalidate the test.</p> <p>When the user attempts to invalidate a XM test for a patient specimen and the associated blood units are currently issued or transfused, VBECS warns the user that the selected test cannot be invalidated because the unit was issued based on the XM result. The user may not continue to invalidate the XM test.</p>
<p>4. Verify the accuracy of the selected patient specimen and patient information (Figure 106).</p> <p>Click OK to accept the results.</p>	<ul style="list-style-type: none"> • Marks these test results as "invalid" and asks the user whether he needs to replace the invalidated results with new test results. <ul style="list-style-type: none"> o Yes (test result replacement is needed) changes the task status to "correction" to return the task to the PTL for completion. o No (test replacement is not needed) marks the erroneous test results as "Test Result Invalidated"; the task status does not change. • When the user selects all the diagnostic tests associated with a specimen, displays an option to mark the specimen unacceptable and invalidate the tests.
<p>5. Select the test to invalidate from the Patient Test Results list.</p> <p>Click Invalidate Results (Figure 107).</p> <p>Enter a comment.</p> <p>Click the Add test to the Pending Task List to allow for replacement results check box, if appropriate.</p> <p>Click OK.</p>	<ul style="list-style-type: none"> • Asks the user whether he wishes to change the specimen status to "unacceptable." <p>NOTES</p> <p>To maintain a clear record and system checks, invalidate both original and repeat ABO/Rh tests.</p> <p>When the user invalidated all of the diagnostic tests associated with a given specimen, VBECS asks whether the user wishes to change the specimen to "unacceptable." No allows the user to invalidate the test results without making the specimen unacceptable. Yes changes the specimen to "unacceptable" and requires a comment.</p>

User Action	VBECS
6. Click OK to exit.	

Figure 106: Invalidate Patient Test Results

VBECS - Invalidate Patient Test Results

Specimen Search Criteria*

Patient: vbpa

Specimen UID:

Expires After: 02/08/2006

Specimen Search Results

Specimen UID	Expires	Patient Name	Patient ID
2960310037	02/07/2006 23:59	VBPATIENT, EIGHTYTWO	000-00-0082
2960310035	02/09/2006 23:59	VBPATIENT, EIGHTYONE	000-00-0081

Include Unacceptable Specimens

Selected Patient Specimen

Specimen UID: 2960310037

Collection Date: 02/07/2006 08:42

Received Date: 02/07/2006 10:58

Expiration Date: 02/10/2006 23:59

Phlebotomist: Kevin

Visual Inspection: Acceptable

Comment:

Selected Patient

Patient Name: VBPATIENT, EIGHTYTWO

Patient ID: 000-00-0082

Patient Sex: Male


Date of Birth: 01/01/1945

ABO/Rh: O Neg

Buttons: OK, Cancel

* Required Field

Figure 107: Confirm Invalidation

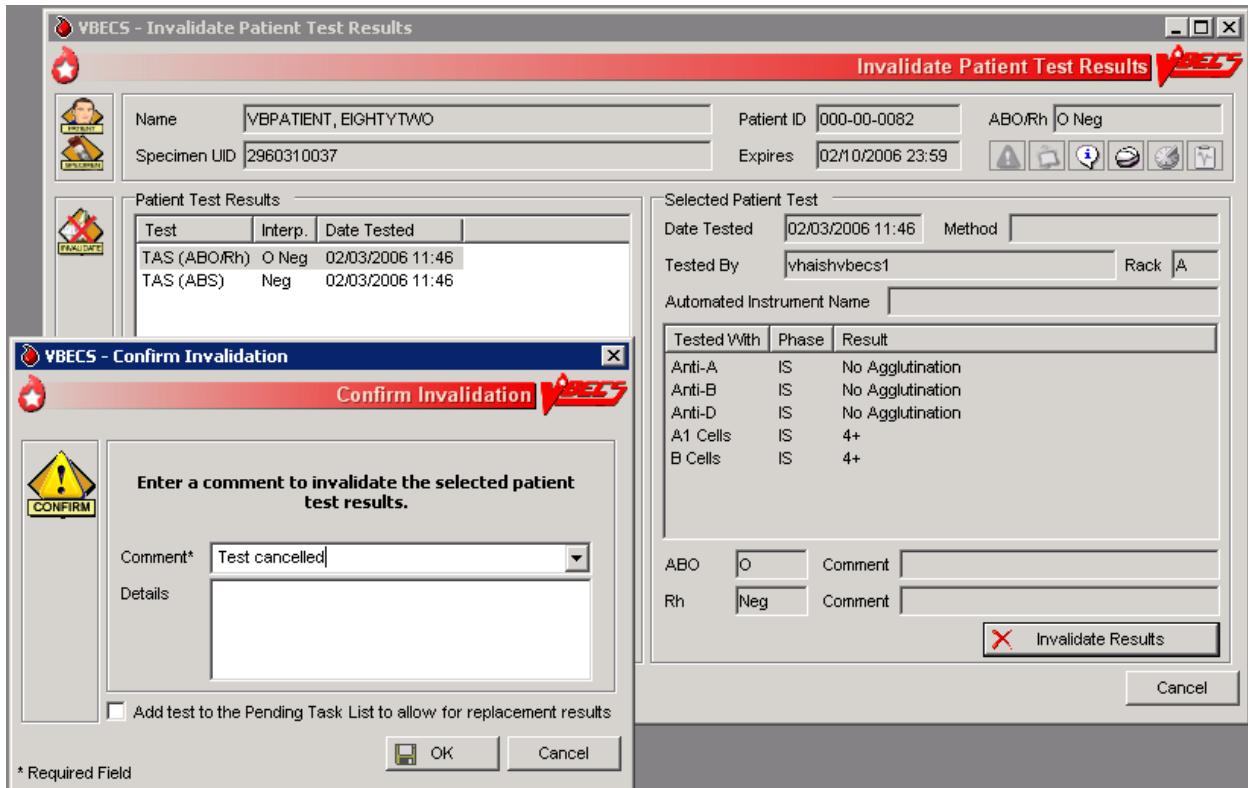

Sub-Window (Bottom):

- Title Bar:** VBECS - Confirm Invalidation
- Message:** Enter a comment to invalidate the selected patient test results.
- Comment:** Test cancelled
- Details:** (Empty text area)
- Checkboxes:**
 - Add test to the Pending Task List to allow for replacement results
- Buttons:** OK, Cancel

Remove Final Status

The user removes a final status from a unit. All other statuses are manipulated through normal system rules and practices.

Assumptions

- The user is logged into the division containing the unit record.
- The unit is in a final (transfused, transferred, modified, discarded) status.
- The user must invalidate the modification status of the split or pooled units, and then the original unit or pool.

Outcome

- The new information was saved with the unit's historic record.
- The unit is in the status immediately prior to the final status.

Limitations and Restrictions

- This option applies only to the release of units in final statuses.

Additional Information

- This option allows a user to restore a unit to inventory from a final status. It may be used to begin undoing a transfusion episode (step by step in reverse order) or correct a unit modification record.
- “Issued” is not a final unit status.
- This option is not used to add units to or remove them from a pool.

User Roles with Access to This Option

►►► Traditional Supervisor

Remove Final Status


The user corrects a clerical error: a unit was processed and moved to a final status in the computer incorrectly.

 *Removing final status of a unit reported as transfused by the Transfusion Verification application will create a discrepancy between the systems. Changes or updates to a transfusion record must be made in the Transfusion Verification application first.*

User Action	VBECS
1. Select Supervisor from the main menu. Select Remove Final Status .	<ul style="list-style-type: none">• Displays options available to supervisors.• Displays fields for selecting a unit and removing its final status.

<p>2. Scan or enter a unit ID or product code, or click the find button to select a unit from the drop-down list.</p>	<ul style="list-style-type: none"> • Displays unit data, including the previous and current unit statuses. • Warns the user to proceed with caution. <p>NOTES</p> <p>VBECS warns the user when the unit does not have a final status or was not modified (other than pool, thaw/pool, and split). The user may select another unit.</p> <p>When the final status is “modified” (pooled, thaw/pooled, or split), VBECS verifies that the target unit is invalidated. When the target unit was not invalidated, the user may not continue.</p>
<p>3. Verify that the unit selected is the desired unit.</p> <p>Click OK.</p> <p>Enter a comment (Figure 108).</p>	<ul style="list-style-type: none"> • Asks the user if he would like to remove the final status. • Requires the user to enter an explanatory comment. <p>NOTES</p> <p>VBECS updates the unit status to the previous status and maintains restrictions on the original unit.</p>
<p>4. Review the data, click OK to save, and click OK again to confirm the save.</p>	<ul style="list-style-type: none"> • Displays the unit data and asks the user to verify the return to the previous status. • Allows the user to review his entries and choices and save the correction to the unit record or exit without saving. • Updates the unit record. <p>NOTES</p> <p>When VBECS removes a final status of “discard” or “transfer,” it also removes return credits previously assigned to the unit.</p>

Figure 108: Remove Final Unit Status

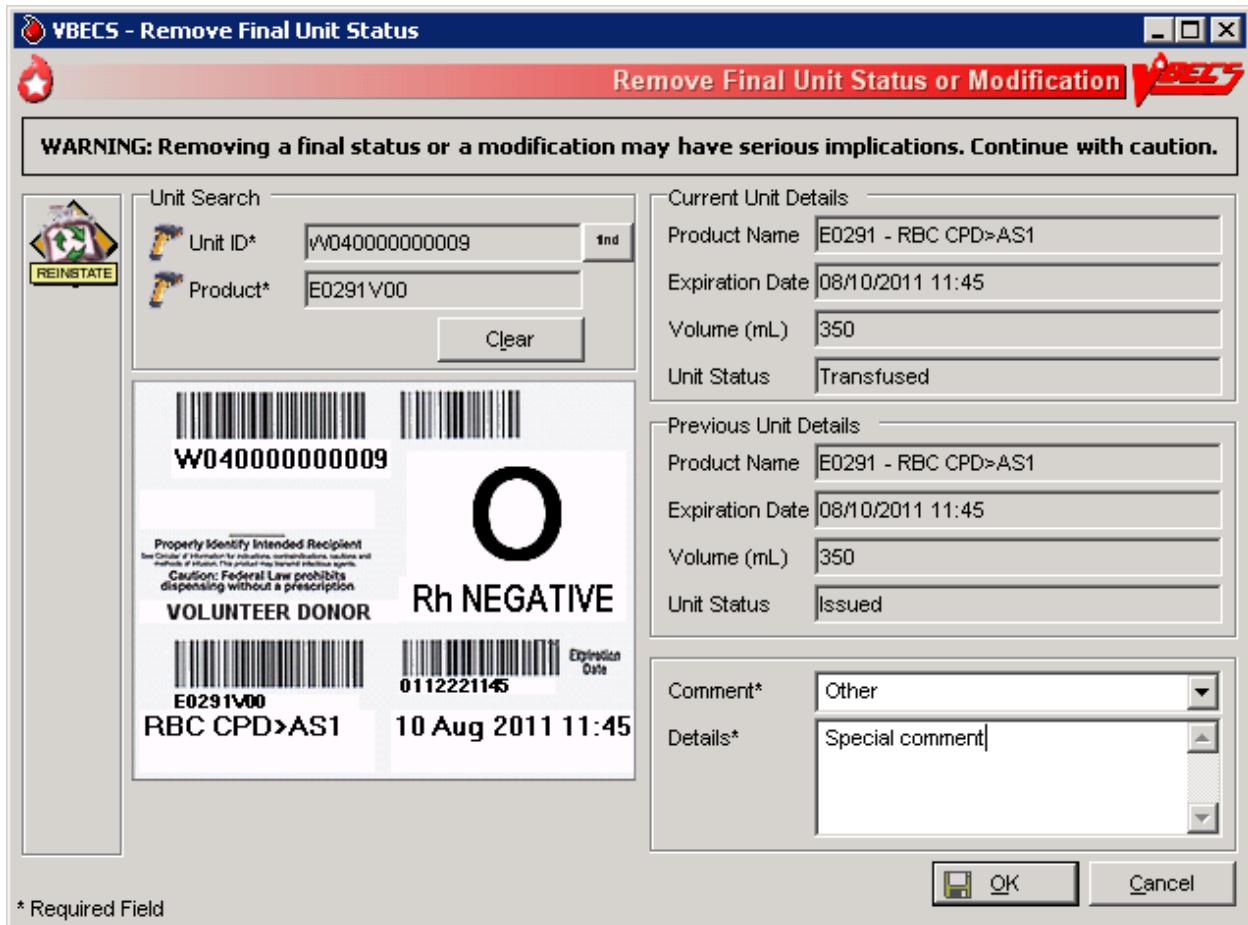

* Required Field

This page intentionally left blank.

Processing Orders

Accept Orders: Accept an Order

The user accepts the receipt of an order.

Assumptions

- Orders were placed in CPRS for a patient and a VBECS order number was generated.
- When a specimen is required, the user located it.
- Inappropriate Transfusion Request parameters are defined for orderable component classes.
- The connection to VistA is active.

Outcome

- Ordered tests may be processed and are available to various worklists.

Limitations and Restrictions

- An order for TAS must be accepted before the component order if the user wants to include the component request to the Inappropriate Transfusion Request Report. Accepting the TAS after the component order is marked inappropriate and accepted clears the user's indication that this is an inappropriate request.
- When no CPRS order exists, the user processes the specimen offline.
- The user must accept the Type & Screen (TAS) before he accepts other diagnostic tests or component orders.
- The user cannot click  (Restricted or Assigned Units) to access information on autologous or directed units that may be restricted to a patient.
- Orders placed for patients with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in VistA are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.
- A VBECS order on a Collection List has to be received with a date and time in VistA to make the specimen available in VBECS.

Additional Information

- VBECS automatically cancels previously accepted orders the day after they expire.
- To maintain a specimen with a TAS and diagnostic test orders, the TAS must be processed first, before any diagnostic test order (see Maintain Specimen).
- The user may investigate patient clinical information using VistA and CPRS.
- A facility may capture the phlebotomist's signature and the date and time of specimen collection on a printed order from CPRS. [Local IRM personnel may configure CPRS to print the blood bank order: refer to *Setup Guide* (August 2000) in VistA Document Library, Computerized Patient Record System (CPRS)].

User Roles with Access to This Option

All users

Accept Orders: Accept an Order

The user reviews the list of pending orders and the specimen, determines that local policy requirements are satisfied, and accepts all or part of a routine order. The user accepts the VBECS order and associates the Laboratory generated specimen UID number with the order.

The user uses the Pending Order List (POL) to view additional details, print outstanding orders, and to select, accept, and add orders to the Pending Task List (PTL).

 Changes made in the VBECS Administrator option mapping orders to another VBECS division do not affect delivered orders. Orders delivered to a VBECS division must be completed, rejected, or canceled in that division. Resubmit orders after mapping is completed to send an order to another VBECS division.

User Action	VBECS
1. Select Orders from the main menu. Select Accept Orders .	<ul style="list-style-type: none">• Displays options for processing orders.• Displays the POL and POL search criteria.
2. Select an order to accept from the POL (Figure 109), or In the Pending Order List Search Criteria area, enter the known data, including the date wanted. Click Search to find an order.	<ul style="list-style-type: none">• Displays options for selecting, reviewing, and processing patient data and pending orders.• Allows the user to associate a specimen with an order.• Displays the Patient Information Toolbar that allows the user to review. <p>NOTES</p> <p>Select an order for TAS before selecting a corresponding order for diagnostic tests or components.</p> <p> VBECS emits an audible alert when there are patient Special Instructions or Transfusion Requirements.</p> <p>VBECS alerts the user when the component class ordered requires a specimen. The status may be REQ (specimen required) or NRQ (no specimen required). The user may not continue until specimen information is added.  When modifiers are included in the physician's component order, VBECS provides visual and audible alerts to remind the user to evaluate the request. VBECS does not monitor subsequent action.</p> <p>VBECS alerts the user when units (autologous, directed, or assigned/crossmatched) with a future expiration date are associated with the patient in any division in the database.</p>
3. Scan or enter the specimen UID in the Specimen UID field. (See Maintain Specimen.) (Figure 110) When the component requires a current specimen and a valid	<ul style="list-style-type: none">• Allows the user to enter new specimen information or correct problems with an existing specimen with the current VBECS order.• Allows the user to process an order without an associated specimen for emergency issue and retroactively associate a specimen with the order.

User Action	VBECS
<p>sample is not available, click the Emergency Order check box.</p> <p>Click Accept to accept the order without an associated specimen.</p>	<p>NOTES</p> <p>Alphabetical characters in specimen UIDs must be capitalized: VBECS recognizes only uppercase letters.</p> <p>Specimens collected for diagnostic tests only are to be used for performing the specifically ordered tests and are not available for processing blood component orders.</p> <p>For all orders requiring a specimen, VBECS queries VistALink to verify that the lab order numbers are consistent with the specimen. If they are not consistent, VBECS notifies the user that the specimen's lab number does not match the number associated with the order. VBECS does not accept the order until the user identifies a specimen with a matching lab order number.</p> <p>When the order status is "In Blood Bank," a specimen was associated with the order; when the status is "Specimen Required," the user must identify a specimen by its UID and enter specimen details according to Maintain Specimen.</p> <p><input type="checkbox"/> When an order for a component class requires a specimen (WHOLE BLOOD and RED BLOOD CELLS component classes, and, optionally, based on Setting Component Classes) the user must select this function and replace the existing specimen with a new one. However, when VBECS displays the "emergency issue" order indicator, the user may continue without identifying an associated specimen. VBECS captures details for inclusion in an Exception Report [exception type: Emergency order accepted (no specimen)].</p> <p>Within an order group containing orders for blood components and the Type & Screen (TAS) test, the specimen associated with the TAS test applies to the TAS test and to all blood component orders contained in the order group.</p>
<p>4. Click  to review order details.</p> <p>5. When the urgency is "Pre-Op," click  in the Patient Information Toolbar to review the MSBOS information for the division.</p> <p>View and/or search entries in the MSBOS for possible matches with the patient surgery.</p>	<ul style="list-style-type: none"> Displays order details that may not be edited. Allows the user to review current MSBOS settings. Lists the current MSBOS site parameters, based on generic surgery types. <p>NOTES</p> <p>When the pre-op indicator is "Yes," the user may access the MSBOS to determine whether to include details of the order in the Inappropriate Request Report.</p>
<p>6. Based on the review of the MSBOS parameter associated with the surgery type and the details of the current order, click the Include In Inappropriate Transfusion Request Report check box, according to local policy.</p>	<ul style="list-style-type: none"> Adds the details of the order to the Inappropriate Transfusion Request Report when the user completes Accept Orders: Accept an Order. <p>NOTES</p> <p>To display a list of available autologous and/or directed units in the system for the selected patient, see Blood Availability.</p> <p>A user must accept the corresponding TAS order before including a component order in Inappropriate Transfusion Request Report.</p>
<p>7. Click  in the Patient Information</p>	<ul style="list-style-type: none"> Displays read-only TR and SI patient VBECS data.

User Action	VBECS
Toolbar to review patient data, Transfusion Requirements and Special Instructions.	
8. Click Accept or Reject for each order.	<ul style="list-style-type: none"> • For each order accepted, updates the: <ul style="list-style-type: none"> ◦ Status to “accepted” and adds it to the PTL for further processing. ◦ Existing lab accession number to collect data associated with every order. <p>NOTES</p> <p>When a user rejects an order, VBECS acknowledges the order and cancels it through Accept Orders: Cancel a Pending Order.</p> <p>When no specimen is required, VBECS assigns an order expiration date to three days in the future at 23:59 (the date of acceptance of the component order acceptance is day “0”).</p> <p>When a specimen is required, the order expiration date is the same as the associated specimen. VBECS automatically expires the associated component order when the specimen is expired. (See Appendix B: Figure 135: Order Status Flowchart.)</p> <p>The user evaluates each diagnostic test and component order separately.</p> <p>When the pre-op indicator is “No,” VBECS compares the Laboratory test values included in the order message from CPRS to the threshold values defined for the ordered component class requirements in Configure ICCBBA Component Classes Inappropriate Request Parameters setting (see Component Classes). When the threshold values were exceeded for one or more ordered component classes, VBECS notifies the user that the Laboratory test value exceeds site parameters for the component class orders and asks whether to add the order to the Inappropriate Transfusion Request report. Yes adds the details of the order to the Inappropriate Request Report. No requires a comment to be associated with the order.</p> <p>The user may cancel one or more components of an order through Accept Orders: Cancel a Pending Order to cancel while working in Accept Orders: Accept an Order. When this occurs, VBECS sends a message to CPRS and notifies the VistA Laboratory application that the order was canceled.</p>
9. Repeat Steps 2 through 8 to accept another order, or click Close to exit.	

Figure 109: Accept Orders

The screenshot shows the 'VBECS - Accept Orders' window. At the top, there is a search criteria section with fields for Specimen UID, Patient, Date Wanted (02/03/2006), To (02/03/2006), Lab Order No., and search/clear buttons. Below this is a 'Pending Order List' table with columns: Patient Name, Patient ID, Lab Order No., Urgency, Specimen, Request, Ordered By, and Patient Location. Two rows are listed:

Patient Name	Patient ID	Lab Order No.	Urgency	Specimen	Request	Ordered By	Patient Location
VBPATIENT, EIGHTYTWO	000-00-0082	65028	Routine	Required	TAS	USER, ONE	VBEC LEAVENW
VBPATIENT, EIGHTYTHREE	000-00-0083	65060	Routine	Required	TAS	USER, ONE	VBEC LEAVENW

Below the table is a toolbar with icons for Print, Accept, and Reject. On the left, there is a tree view under 'Order Group' showing 'VBPATIENT, EIGHTYTWO' expanded to show 'TAS', 'RED BLOOD CELLS', and 'FRESH FROZEN PL'. On the right, there are input fields for Patient Name, ID, Ordered By, Request, Comments, Specimen UID, ABO/Rh, Patient Location, Lab Order No., Order Status, Requirements, and Quantity. Buttons for 'Accept' and 'Reject' are at the bottom right. A note at the bottom left says '* Required Field'.

Figure 110: Maintain Specimen

The screenshot shows the 'VBECS - Maintain Specimen' window. At the top, it displays 'Name: VBECS PATIENTC, THIRTYSE', 'Patient ID: 666-00-0037', and 'ABO/Rh: A Pos'. Below this is a toolbar with icons for various actions. The main interface is divided into two sections: 'Patient Specimen' on the left and 'Associated Orders' on the right.

Patient Specimen:

- Specimen UID*: 2912220073
- Collection Date*: 08/10/2011 13:03
- Received Date*: 08/10/2011 13:06
- Expiration Date*: 08/12/2011 23:59
- Phlebotomist ID*: [empty]
- Processing Tech*: Test User1
- Specimen Status: Radio buttons for 'Acceptable Specimen' (selected) and 'Unacceptable Specimen'.
- Comment*: [text input field]
- Details: [text input field]

Associated Orders:

Order	Status	Date Wanted
[empty]	[empty]	[empty]

At the bottom, there are 'OK' and 'Cancel' buttons. A note at the bottom left says '* Required Field'.

Accept Orders: Pending Order List

The Pending Order List (POL) lists pending blood bank orders.

Assumptions

- A CPRS order is acknowledged by VBECS.
- An order is printed or an alert is issued per site parameters.
- The connection to VistA is active.

Outcome

- A list of pending orders displays and may be sorted and searched.
- A list of pending orders may be printed in the default sort presentation.

Limitations and Restrictions

- Orders placed for patients with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in VistA are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator “Interface Failure Alert Recipient” field.

Additional Information

- This option is used for diagnostic and blood component test orders.

User Roles with Access to This Option

All users

Accept Orders: Pending Order List

VBECS displays the POL with search criteria, which allows the user to limit the number of viewable items. The user may review details of an order without taking action on it.

User Action	VBECS
1. Select Orders from the main menu. Select Accept Orders , or Click the Pending Orders alert icon , if available, in the status bar at the bottom of the screen (Figure 109).	<ul style="list-style-type: none">• Displays options for processing orders.• Displays the Pending Order List and Pending Order List search criteria.• Allows the user to search for specific orders.• Provides the criteria to limit the search and subsequent presentation and processing.• Allows acceptance of an order. <p>NOTES _____</p> <p>VBECS displays orders directed to the division by a user during system configuration.</p> <p>See Appendix B: Figure 135: Order Status Flowchart.</p> <p>When searching for pending orders based on patient name, the system displays all pending and current accepted orders for the patient. A current order is a previously accepted order with a future order/specimen expiration date.</p>
2. Change the selection criteria, as needed, to prepare the list or to accept the default presentation.	<ul style="list-style-type: none">• Displays the POL requested by user's selection(s). <p>NOTES _____</p>

User Action	VBECS
	<p>The POL can be sorted by Patient, Procedure, Urgency, Date/time received or date/time wanted.</p> <p>The user may select an entry on the POL to view additional details for the order, including:</p> <ul style="list-style-type: none"> • Patient Name • Patient ID • Additional component class or diagnostic tests included in the same order group as the selected order. (These are displayed whether they are "in progress" or still on the POL. There should be some indication to distinguish those that are already in progress.) • Order comments (additional requirements from CPRS) • Quantity (component class orders only) • Patient hospital location • Lab order # • Patient ABO/Rh (if known) • The CPRS Order ID, Group Order ID for all diagnostic tests and component orders placed in one CPRS order is the CPRS Order ID used by VBECS on the Order History Report and the Pending Task Lists.
3. To continue processing the POL, review the displayed list, print the list, or go to Accept Orders: Accept an Order.	<ul style="list-style-type: none"> • Allows the user to return to Step 1 to change the presentation of the list, review details of pending orders, print the list, or select an order or proceed to Accept Orders: Accept an Order. <hr/> <p>NOTES</p> <p>When a user requests a POL report and VBECS does not find any information, VBECS notifies the user and asks whether he wishes to continue to print. The user may cancel or print the report, which will include "No pending orders found."</p>
4. Exit.	

Display Order Alerts

VBECS notifies the user of new and updated patient orders.

Assumptions

- VBECS receives a new order or an order update message from VistA for the user's division.

Outcome

- VBECS displays notification icons only to users logged into VBECS.

Limitations and Restrictions

- VBECS issues no notifications for order cancellations.

Additional Information

- The user may set VBECS to automatically print an order alert with limited information on the designated printer when VBECS acknowledges the order.
- The alert displays the individual test's (child) CPRS ID. The individual test's (child) CPRS ID is not displayed on the Pending Order List, elsewhere in VBECS or on reports and is not usable in search functions.

User Roles with Access to This Option

All users

Display Order Alerts

VBECS notifies the user that order updates occurred and allows him to accept the new orders.

User Action	VBECS
<p>1. Hold the mouse over the icon to view a tool tip, or</p> <p>Click  to view a panel showing new or updated orders and hide the icon.</p> <p>Click  again to hide the panel.</p>	<ul style="list-style-type: none">• Acknowledges a new or updated patient order.• Displays an order notification icon in the status bar when new orders are received, according to local settings.• Displays information regarding new orders and order updates and allows the user to accept them. <p>NOTES</p> <p>When the division is configured to alert the user of orders received by VBECS, VBECS notifies the user of ordered components or ordered tests that were:</p> <ul style="list-style-type: none">• Directed to the user's division.• Added or updated in VistA. <p>When the division is configured to alert the user of orders needed in the next specified number of hours, VBECS also notifies the user of ordered components or ordered tests with a required time less than a specified number of hours in the future.</p>

Display Patient Update Alerts

VBECS notifies the user of patient update, death, and merge events.

Assumptions

- A patient demographic update is received from VistA for the user's division for a VBECS patient with an active order.

Outcome

- The patient demographics are updated in VBECS.

Limitations and Restrictions

- VBECS does not notify users of additional VistA patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed VistA patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor updates the patient records and clears the icon.
- Patient updates involving patient's with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in VistA are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.
- Patient update alerts only display when active component or diagnostic test orders are pending or active³ at the time the update occurs.

Additional Information

- None

User Roles with Access to This Option

All users

Display Patient Update Alerts

VBECS notifies the user that patient updates occurred and allows him to view additional details.

User Action	VBECS
<p>1. Hold the mouse over the Patient Alert icon to view a tool tip, or</p> <p>Click  to view a menu and select Updated Patients, Deceased Patients, or Patient Merge to view the updates.</p> <p>Click  again to hide the menu.</p>	<ul style="list-style-type: none">• Displays a Patient Alert icon in the status bar when it receives a patient update or death event message and the patient has active or pending orders in the user's VBECS division.• Displays a Patient Alert icon in the status bar when it receives a merge event message from VistA. <p>NOTES —————</p> <p>The user may also select Patients from the main menu, Patient Updates, and Updated Patients, Deceased Patients, or Patient Merge to view the updates.</p>

³ Active: VBECS order that is not rejected, canceled or completed. When the patient update is received and the order is 'active', the update is displayed regardless of order cancelation until it is viewed by a VBECS user.

Display Patient Merge Alerts

VBECS notifies the user of VistA patient merges.

Assumptions

- A patient demographic merge alert is received from VistA for the user's division for a VBECS patient.

Outcome

- VBECS notified the user that patient records were merged in VistA and that the user must take action to match the records in VBECS.

Limitations and Restrictions

- VBECS does not notify users of additional VistA patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed VistA patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor updates the patient records and clears the icon.

Additional Information

- VBECS notifies users of VistA patient merges until a Traditional Supervisor clears the Patient Alert icon. See Patient Merge for instructions.
- VBECS displays an alert to users for a VistA patient merge event if either the merge-from or merge-to patient exists in the VBECS system when the merge message is received.

User Roles with Access to This Option

All users

Display Patient Merge Alerts

VBECS notifies the user that patient merges occurred in VistA and allows him to view additional update details.

User Action	VBECS
<p>1. Hold the mouse over the Patient Alert icon to view a tool tip, or</p> <p>Click  to view a menu and select Updated Patients, Deceased Patients, or Patient Merge to view the updates.</p> <p>Click  again to hide the menu.</p>	<ul style="list-style-type: none">• When it receives a patient update or death event message and the patient has active or pending orders in the user's VBECS division, displays a Patient Alert icon in the status bar.• When it receives a merge event message from VistA, displays a Patient Alert icon in the status bar. <p>NOTES</p> <p>The user may also select Patients from the main menu, Patient Updates, and Patient Merge to view the updates.</p> <p>VBECS displays a patient update or merge icon alert to all users until all merge events are cleared and all updates, including death notices are viewed.</p>

Maintain Specimen

A specimen for testing is received in the blood bank and is processed as part of initial order acceptance, as described in Accept Orders: Accept an Order.

Assumptions

- Pending CPRS orders and active orders require an acceptable specimen.
- A specimen was collected and previously accessioned in the VistA Laboratory system, which generated the specimen's universal ID (UID).

Outcome

- The user accepted or rejected a specimen.

Limitations and Restrictions

- When a search by UID causes VBECS to time out, repeat the search by UID and patient name or by patient name only.
- Specimen recollection requests must be made offline.
- Designating a specimen as “unacceptable” after accepting a CPRS order closes that order and requires a new CPRS order and sample.

Additional Information

- A specimen may be marked unsuitable for use after its acceptance.
- A specimen expiration date may be edited after acceptance.
- A specimen may be recollected without canceling or reordering the VistA order. This practice is dictated by local policy and procedures.
- Specimens marked as originally associated with diagnostic tests only or unacceptable are not considered available specimens when CPRS is queried for specimen availability.

User Roles with Access to This Option

All users

Maintain Specimen

A specimen for testing is received in the blood bank and is processed as part of initial order acceptance, as described in Accept Orders: Accept an Order.

Information about a previously accepted specimen may be updated. The user may verify whether an active specimen is currently in the blood bank or update the information concerning a previously accepted specimen. After the user checks the physician's order against the specimen requirement, the user verifies the specimen's condition, labeling, and accompanying paperwork, according to local policy. The user verifies the specimen information and acceptability and prepares the specimen for testing.

The user may edit specimen information, such as extending the expiration time according to local policy, and marking a specimen as unsatisfactory or expired so that future orders require that a new specimen be collected.

When a specimen is required, the user must find that specimen before acknowledging that order, as described in Accept Orders: Accept an Order.

User Action	VBECS
1. Select Orders from the main menu. Select Maintain Specimen .	<ul style="list-style-type: none"> • Displays options for processing orders. • Displays fields for entering patient specimen information.
2. Scan or enter the specimen's UID.	<ul style="list-style-type: none"> • Displays a field to enter the specimen UID and allows the user to search for existing specimens. • Displays previously entered information for the specimen and accommodates data entry of new specimen information. <p>NOTES</p> <p>The user may:</p> <ul style="list-style-type: none"> • Create and/or maintain a specimen associated with an order from Accept Orders: Accept an Order. • Update specimen information from a previously received specimen. <p>Specimens associated with blood bank diagnostic tests are not available in Maintain Specimen. Specimens collected for diagnostic tests only are to be used for performing the specifically ordered tests and are not available for processing blood component orders.</p> <p>The orders and associated specimens are maintained for three days before the order or specimen expires: Type & Screen (TAS), component orders (RBC, FFP, CRYO, PLT, WB, OTHER).</p> <p>The list of orderable and reflex tests below may be added to a specimen maintained for a component order, but may not have a component order added to them: they do not have maintained specimens.</p> <p>These tests have a specimen drawn and used for their ordered test. The specimens are not maintained for component preparation:</p> <ul style="list-style-type: none"> • DAT • ABS • ABO/Rh • Transfusion Reaction Workup (TRW) <p>VBECS allows the user to search all existing patient specimens associated with a specific patient. The user may enter the patient name or a specimen UID to access the list.</p>
3. Enter new specimen information or edit previously entered information (Figure 110).	<ul style="list-style-type: none"> • Displays uneditable patient demographic information. • Displays an option to edit. <p>NOTES</p> <p>When the user extends the specimen expiration date and time, VBECS calculates the number of days by which the specimen was extended (based on the originally calculated expiration date) and the newly entered expiration date. VBECS compares the difference to the site parameter for the maximum number of extension days allowed for a specimen. When the entered expiration date exceeds the parameter, VBECS warns the user. There is no override.</p>

User Action	VBECS
	<p>The default expiration date is calculated as three days in the future with the day of the original specimen collection as "0." For example, a specimen collected on September 16, 2007 (regardless of the time), is day "0" and the calculated expiration date and time are September 19, 2007 @ 23:59 (midnight).</p> <p>When a new or existing specimen is being associated with a pending component order during Accept Orders: Accept an Order, the specimen must be acceptable.</p> <p>When there is a specimen requirement for orders, each order must be evaluated for a valid specimen before it can be accepted.</p> <p>VBECS displays specimen UIDs for selection by CPRS-associated maintainable specimen and then by acceptable and unexpired maintainable specimens, in reverse chronological order.</p> <p>VBECS allows the user to include expired and unacceptable specimens.</p> <p> When the user attempts to save an expiration date past the default of three days, VBECS checks the patient's record for evidence of transfusion within the previous three months and the history of antibody identification, regardless of the division. If either is found, VBECS emits an audible alert, warns that the patient has evidence of transfusion within previous three months and/or antibody formation, and instructs the user to enter a comment to extend the specimen expiration date beyond three days. The user must enter a canned comment (with the specimen extension) and/or free-text comment to continue. VBECS captures details for inclusion in an Exception Report (exception type: expiration date extended).</p> <p>When the specimen is deemed unacceptable, the status of the orders associated with the specimen remains "pending" or "active."</p> <p> When a user attempts to reactivate a previously expired or unacceptable specimen, VBECS emits an audible alert, warns that this is a previously expired or unacceptable specimen, and instructs the user to enter a comment to make it available. A comment is required to continue. VBECS captures details for inclusion in an Exception Report (exception type: unacceptable/expired specimen used).</p> <p>VBECS requires the date and time a specimen was collected and the phlebotomist ID.</p> <p>When a specimen is processed, VBECS displays the details of associated orders (read only).</p> <p>The user must enter a comment when marking a specimen as unacceptable.</p>

User Action	VBECS
4. Enter or verify specimen-related information from Step 3.	<ul style="list-style-type: none"> Allows the user to review and edit previously entered information.
5. Review the entered information and click OK to save.	<ul style="list-style-type: none"> Saves the data for use in further specimen order processing and subsequent retrieval. <p>NOTES —————</p> <p>When saving or updating an expired specimen, VBECS warns the user and asks whether he wishes to continue. The user may continue without entering a comment.</p> <p>The UID initially used to identify the specimen becomes the VBECS specimen number, which VBECS retains and uses to process subsequent orders during the life of the specimen, including CPRS, VBECS, and Laboratory orders, and UIDs associated with the specimen. A historical record lists changes made to the specimen and is available for review.</p> <p>A permanent link remains between an exhausted or unacceptable specimen UID and its order, which allows the user to document the recollection of a new specimen, according to local policy.</p>
6. Click OK to exit.	

Patient Testing: Cancel an Active Order

The user cancels an accepted order on the Pending Task List (PTL) and adds a comment to justify the cancellation.

Assumptions

- The order and the user are in the same division.
- A user accepted an order that originated in CPRS through Accept Orders: Accept an Order, or created a reflex test order through Order Reflex Tests.
- The connection to VistA is active.

Outcome

- The order was removed from the PTL.
- An order update message was sent to CPRS.

Limitations and Restrictions

- Cancellation of one order does not affect other orders placed with it.
- Before a user may cancel a partially completed task, he must invalidate the results for that task, thereby changing its status to “not started.”

Additional Information

- None

User Roles with Access to This Option

All users

Patient Testing: Cancel an Active Order

The user selects one or more tasks to cancel from the PTL, cancels the accepted order from within VBECS, and adds a comment regarding the cancellation.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the Pending Task List in the Diagnostic Tests tab. <p>NOTES _____</p> <p>Click the Component Orders tab to view its PTL.</p>
2. Click one to four check boxes in the PTL in the Diagnostic Tests tab or the Component Orders tab to select tasks (Figure 111). Click OK to continue, or click Cancel Order(s) to cancel the selected tasks. Select or enter a comment. Click OK to confirm the cancellation.	<ul style="list-style-type: none">• Deletes the selected tasks from the PTL and sends a message to CPRS with an order status update. <p>NOTES _____</p> <p>Select a canned comment and/or free-text comment VBECS allows a user to cancel only tasks that are not started. See Appendix B: Figure 135: Order Status Flowchart.</p>

Figure 111: Pending Task List Diagnostic Tests

The screenshot shows the 'Pending Task List' window for VBECS. At the top, there are tabs for 'Diagnostic Tests' and 'Component Orders'. Below the tabs are 'Pending Task List Search Parameters' which include fields for 'Specimen UID', 'Patient', and 'CPRS Order ID', along with dropdown menus for 'Diagnostic Tests', 'Task Status', and 'Priority'. A 'Search' button and a 'Clear' button are located at the bottom of these parameters. The main area displays a 'Pending Task List' grid with columns: Priority, Patient Name, Patient ID, Specimen UID, CPRS Order ID, Diagnostic Test, Order Comments, Date Wanted, and Task Status. Two rows of data are visible:

Priority	Patient Name	Patient ID	Specimen UID	CPRS Order ID	Diagnostic Test	Order Comments	Date Wanted	Task Status
Routine	VBPATIENT, EIGHTY	000-00-0080	2960310033	22165251	TAS		1/31/2006	Not Started
Routine	VBPATIENT, EIGHTYTWO	000-00-0082	2960310037	22165259	TAS		1/31/2006	Not Started

At the bottom of the window are buttons for 'Print', 'Cancel Order(s)', 'OK', and 'Cancel'.

Patient Testing: Pending Task List

The user accesses a Pending Task List (PTL) to select tasks (accepted orders), organize a batch of tests to be performed in a session, process patient component requests, and view and/or print a list of outstanding blood bank tasks.

Assumptions

- CPRS orders were accepted or a reflex order was placed within the VBECS division.

Outcome

- The task is available for further processing.

Limitations and Restrictions

- This option is for diagnostic, blood component, and reflex test orders placed in or directed to the user's division.
- The PTL does not include daily quality control, equipment maintenance schedules, or pending ABO/Rh confirmation and auxiliary tests on red blood cell units.
- Users may perform testing for no more than four patients in one VBECS session.

Additional Information

- This is the main entry point for other testing options.
- The PTL does not display completed diagnostic tests.
- When a specimen is not associated with a component order, partially completed and filled component orders remain on the Pending Task List for 10 days from date of order acceptance.
- When a specimen is associated with a diagnostic test or component order, partially completed diagnostic tests and filled component orders remain on the Pending Task List until the specimen expires.
- VBECS cancels accepted orders (tasks) the day after they expire.
- Multiple users are prohibited from performing the same task on the same unit of blood for the same patient.
- The task may be selected and assigned to a technologist, preventing duplicate selection by another technologist.

User Roles with Access to This Option

All users

Patient Testing: Pending Task List

The user accesses a list of pending tasks. VBECS displays the PTL with search criteria to allow the user to limit the number of items viewed.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing (Figure 111).	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the PTL and PTL search parameters in the Diagnostic Tests tab or in the Component Orders tab.
2. Click the check boxes in the PTL to select one to four tasks for processing and click OK , or	<ul style="list-style-type: none">• Allows the user to limit the search and display.• Prompts the user to select a task for continued processing, print a task list, or return to Step 1 to change the display of the list.

User Action	VBECS
Click Print to print the list.	<p>NOTES</p> <p>Additional information is available in Appendix H: Frequently Asked Questions.</p> <p>The user may click Cancel Order(s), select a comment, and click OK to exit.</p> <p>The user may select one to four tasks for testing at a time. The tasks may be component class order tasks, diagnostic test tasks, or a combination of both.</p> <p>VBECS displays “NRQ” (not required) when the task does not require a specimen for processing.</p> <p>For each order, VBECS displays the specific diagnostic test or component order. The default searches for the current date's tasks, then yesterday's tasks, and so on until the task status is “completed.” The display includes the task's status and a user name, when applicable.</p> <p>VBECS alerts the user when he attempts to select multiple tasks of the same type for a patient, and does not allow the selection of two or more tasks of the same type ordered for a patient during a testing session. Tasks are considered to be of the same type when they need to accommodate the use of the same testing grid in Patient Testing: Record Patient Test Results. For example, a user may not select a TAS task and one of its subtests that is ordered in a testing session. VBECS provides for one testing grid display per patient per test, not per specimen UID. If the user selected a TAS and an ABS, VBECS displays only the first one selected for result entry.</p> <p>Task and order statuses and functions that are displayed on the PTL:</p> <ul style="list-style-type: none"> • Not Started: the order was not selected. • Partially Completed <user ID>: allows selection by another user for a diagnostic test to complete, invalidate, or delete the work on this specimen. <i>Component orders remain in this status until the specimen and order expire.</i> • Filled: This order status is assigned when the total number of units ordered is updated to a unit status of “issued” or “transfused.” • Correction: previous results were invalidated; a corrected entry is required. <p>Task and order statuses that are not displayed for selection on the PTL:</p> <ul style="list-style-type: none"> • Completed: applies to diagnostic tests • Expired: applies to component orders at expiration of the associated specimen. <p>See Appendix B: Figure 135: Order Status Flowchart for the transmission of orders between CPRS and VBECS.</p> <p>When a user prints a PTL and no pending tasks are found, VBECS asks whether the user wishes to continue to print. The user may cancel or print the report, which states that no pending</p>

User Action	VBECS
	orders were found.
3. Click OK or Cancel to exit.	

Accept Orders: Cancel a Pending Order

The user cancels an unaccepted order on the Pending Order List (POL).

Assumptions

- An order was not accepted by a user.
- The connection to VistA is active.

Outcome

- An order update message was sent to CPRS.
- The order was removed from the POL.

Limitations and Restrictions

- Cancellation of one order does not affect other orders placed with it.

Additional Information

- Each pending order is managed individually.
- VBECS does not automatically clear or cancel pending or unaccepted orders from the POL.

User Roles with Access to This Option


All users

Accept Orders: Cancel a Pending Order

The user selects an unaccepted order for cancellation.

User Action	VBECS
1. Select Orders from the main menu. Select Accept Orders .	<ul style="list-style-type: none"> • Displays options for processing orders. • Displays the Pending Order List and Pending Order List search criteria.
2. Select an order to cancel from the Pending Order List (Figure 112). Click Reject .	<ul style="list-style-type: none"> • Requires the user to enter a comment. • Requires the user to confirm the cancellation of the pending order.
3. Enter or select a comment in the Comment field. Enter details in the Details field, if desired. Click OK and Close to confirm the order cancellation.	<ul style="list-style-type: none"> • Removes the order from the Pending Order List and sends a message to CPRS with an order status update. <p>NOTES —————</p> <p>VBECS does not cancel a lab order if the lab order has already been accessioned. See CPRS VBECS Order Details and VBECS Order Cancellation for more information.</p>

Figure 112: Accept Orders: Reject an Order

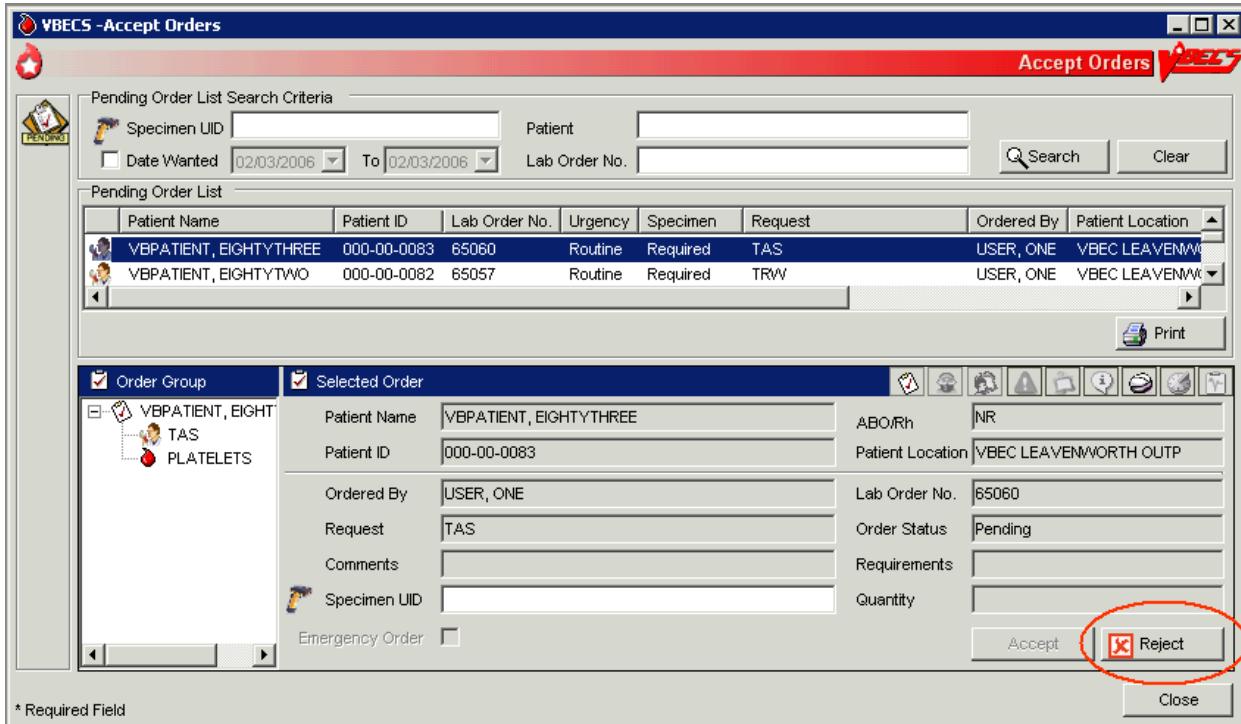

Order Reflex Tests

The user may order additional blood bank tests for a complete or incomplete accepted patient specimen to add to the Pending Task List (PTL). The ordered test or test group limits the type of reflex and repeat test order that may be associated with a specimen.

Assumptions

- A patient's test result requires further investigation or repetition of a diagnostic test.

Outcome

- One or more reflex orders are in the PTL with a "routine" urgency.

Limitations and Restrictions

- A crossmatch (XM) and a repeat XM cannot be ordered unless they are associated with a specimen with a future expiration date and a phlebotomist is associated with that specimen.
- When a user selects weak D for the antigen typing test and the selected patient is Rh positive, VBECS warns that weak D testing cannot be performed. There is no override and VBECS does not add the test.

Additional Information

- Unlimited repeat tests may be ordered and completed, but are not available to CPRS or to billing.

User Roles with Access to This Option

All users

Order Reflex Tests

The user orders additional blood bank tests for a patient specimen that was already accepted and that may or may not be completed. The ordered reflex test(s) are added to the PTL.

User Action	VBECS
1. Select Orders from the main menu. Select Order Reflex Tests .	<ul style="list-style-type: none">• Displays options for processing orders.• Displays specimen search criteria and search results. NOTES ————— The user may also access Order Reflex Tests by clicking Order Reflex Tests when executing Patient Testing or by right-clicking the testing grid that generates the reflexive order.
2. Enter the patient name and/or the specimen UID. Click Search , if necessary.	NOTES ————— The user may identify specific specimen for reflex test ordering by scanning or by manual data entry of the UID on the specimen.
3. Select an item from Specimen Search Results. Review the information in the Selected Patient Specimen area and in the Patient Details and Order Details tabs (Figure 113). Click OK to accept the order.	<ul style="list-style-type: none">• Displays an option to search for specimens currently available (in-date) for a patient by standard patient search conventions.• Accommodates direct selection of a specimen by scanning or by entering a specimen UID.• Displays available (in-date) patient specimens or orders in reverse chronological order within a date range defined by the site parameter. Only one patient specimen may have reflex tests added in one session. The expired and unacceptable specimens will follow.• Displays a message directing the user to complete Patient Testing

User Action	VBECS
Respond to messages, as required.	<p>when the selected test is not completed. User returns to Patient Testing.</p> <ul style="list-style-type: none"> • Displays a message when a specimen is expired or unacceptable. User may click Yes to continue with the order or No to select or request another specimen, as appropriate. • Displays the date/time/division of the specimen collection when selecting a specimen for ordering reflex tests.
	NOTES ————— <p> VBECS emits an audible alert, warns that the specimen is expired, and instructs the user to enter a comment to continue. When a user enters or scans an expired specimen. VBECS requires a comment and captures details for inclusion in an Exception Report (exception type: reflex test ordered on an expired specimen).</p>
4. When a specimen is already selected, go to Step 5. Select or enter a patient specimen that requires the additional test.	<ul style="list-style-type: none"> • Displays all orders associated with a specimen that can be used to request reflex tests for an existing diagnostic test order within an order group. • Automatically associates the reflex test with the order when there is only one order associated with the specimen.
5. Click one or more check boxes to select reflex tests from the Order Reflex Test(s) list.	<ul style="list-style-type: none"> • Displays the selected tests for review. • Allows the user to select or to deselect reflex tests.
6. Click OK to save and exit. Enter a comment, if required. cancellation.	<ul style="list-style-type: none"> • Adds the ordered patient reflex tests to the PTL.

Figure 113: Order Reflex Test(s)

VBECS - Order Reflex Test(s)

Order Reflex Test(s)

Name	VBPATIENT, EIGHTYONE	Patient ID	000-00-0081	ABO/Rh	O Neg
Specimen UID	2960310035	Expires	02/03/2006 23:59		

Order Reflex Test(s)

Order	Test Code	Reflex Test Name
<input type="checkbox"/>	Repeat ABO/Rh	Repeat ABO/Rh
<input type="checkbox"/>	Repeat ABS	Repeat ABS
<input checked="" type="checkbox"/>	DAT	DAT
<input type="checkbox"/>	ABID	ABID
<input type="checkbox"/>	AGwD	Weak D Typing
<input type="checkbox"/>	AGC	C Typing
<input type="checkbox"/>	AGE	E Typing
<input type="checkbox"/>	AGc	c Typing
<input type="checkbox"/>	AGe	e Typing
<input type="checkbox"/>	AGK	K Typing
<input type="checkbox"/>	AGk	k Typing
<input type="checkbox"/>	AGKp(a)	Kp(a) Typing
<input type="checkbox"/>	AGKp(b)	Kp(b) Typing
<input type="checkbox"/>	AGJk(a)	Jk(a) Typing
<input type="checkbox"/>	AGJk(b)	Jk(b) Typing

* Required Field

OK Cancel

This page intentionally left blank.

Patient Services

Patient Search Screen

VBECS displays the Patient Search screen to search for a specific patient in VBECS.

Assumptions

- The patient record exists in VBECS with the exception of specific searches that return patient information from VistA. For the Edit Unit Information, Incoming Shipment, and Document ABO Incompatible Transfusion options, VBECS searches VistA.

Outcome

- VBECS presents the user with a list of patients for selection based on their search input.

Limitations and Restrictions

- None.

Additional Information

- The Patient Search Screen is used by options throughout VBECS to select patients.
- The Patient Location displayed in Figure 114 is either a representation of the location provided by the subscription to Admit, Discharge and Transfer events or in the event that a patient has been discharged, the location displayed could be the location where a CPRS order was placed for the patient. Once a patient is admitted, the location will be updated with the location at time of admission. If a patient is discharged and no orders are placed for the patient post discharge, no location will be provided.

User Roles with Access to This Option


All users

Patient Search Screen

User Action	VBECS
<p>1. To access the Patient Search screen, open one of the options listed in Notes. To search for a patient, enter the patient's full name (last name first) or ID (Figure 114), or Enter the last four digits of the ID, or the initial of the last name and the last four digits of the ID, and click Search.</p>	<p>• Displays the Patient Search screen. • Lists patients that meet the search criteria.</p> <p>NOTES —————</p> <p>Options that allow the user to search for a patient:</p> <ul style="list-style-type: none">• Blood Units from the main menu. Select Edit Unit Information (for restricted units only).• Patients from the main menu. Select:<ul style="list-style-type: none">○ Recent Orders○ Blood Availability○ Recent Transfusions/Issued Units○ Patient Testing Worklist Report○ Print Unit Caution Tag & Transfusion Record Form○ Issue Blood Components○ Post-Transfusion Information○ Return Issued Units To Blood Bank and click OK.○ Special Instructions & Transfusion Requirements○ Transfusion Reaction History

User Action	VBECS
	<ul style="list-style-type: none"> ○ Patient History Report ● Shipments from the main menu. Select Incoming Shipment (for restricted units only). ● Orders from the main menu. Select Recent Patient Orders. ● Reports from the main menu. Select: <ul style="list-style-type: none"> ○ Blood Availability (for individual patient only) ○ Medication Profile ○ Patient History Report ○ Transfusion Effectiveness Report (for individual patient report only) ● Supervisor from the main menu. Select: <ul style="list-style-type: none"> ○ Document ABO Incompatible Transfusion ○ Justify ABO/Rh Change
2. Select a patient from the list.	<ul style="list-style-type: none"> ● Displays patient-associated data. ● Displays the Patient Information Toolbar. ● Displays Special Instructions (SIs) and Transfusion Requirements (TRs), when appropriate. <p>NOTES —————</p> <p>Displays a message when the user selects a patient with a blank first or last name, or the name is too long. The user is instructed to contact the VBECS Administrator for assistance.</p>

Figure 114: Patient Search

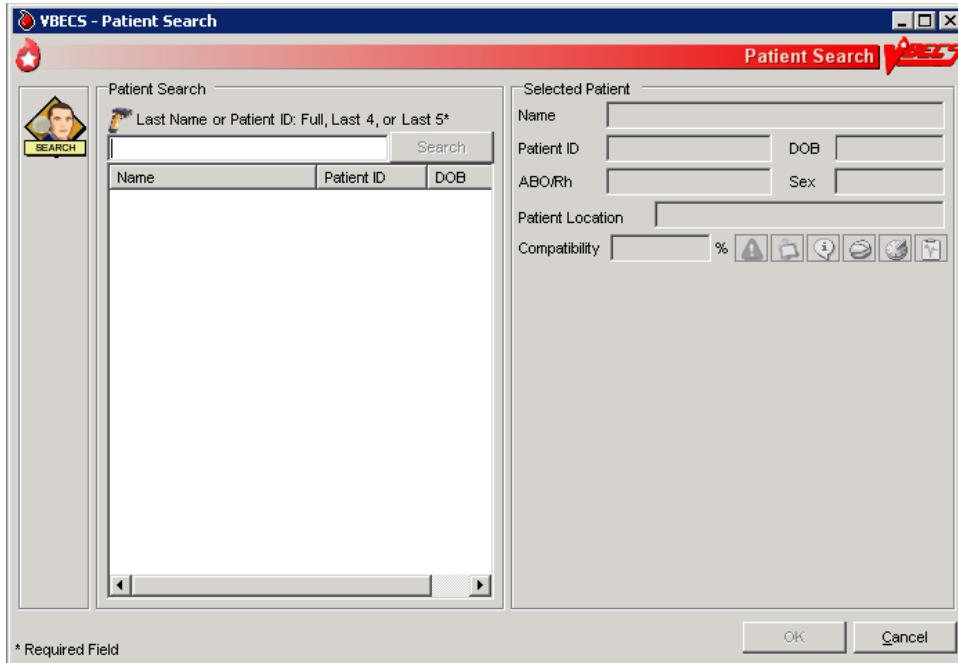

Patient Merge

VBECS notifies the user of VistA patient merges and displays details for his review.

Assumptions

- The connection to VistA is active.
- VBECS received the message about a VistA patient merge.
- VBECS accepts merge messages when either or both the merge from or merge to patients exist in VBECS.

Outcome

- VBECS notified the user that patient records were merged in VistA.

Limitations and Restrictions

- The Traditional Supervisor or above must update patient records manually.
- VBECS does not notify users of additional VistA patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor completes the update and clears the icon.
- Patient merges involving patient's with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in VistA are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.

Additional Information

- If more than two records are merged in VistA, VBECS will receive a corresponding number of merge messages.

User Roles with Access to This Option

All users

Patient Merge

The VBECS user examines records that were merged in VistA and decides whether to update them.

 VBECS does not create links for merged records. Users must update patient information and Special Instructions to show the merged-patient relationship.


 Patient merges that do not change the patient name or social security number (e.g. patient record moved to a new DFN) are displayed to the VBECS user. However, the records displayed in the Merge window will not appear to have any differences. Users are advised to contact customer support if a record merge of this type occurs.

User Action	VBECS
1. Select Patients from the main menu.	<ul style="list-style-type: none">• Displays  and notifies the user of the VistA patient merge.• Lists the patient records to be updated.• Displays limited merge-from and merge-to patient information for

User Action	VBECS
Select Patient Updates. Select Patient Merge, or Click  in the Patient Information Toolbar.	<p>evaluation, including:</p> <ul style="list-style-type: none"> o Patient name o Patient ID, if available o Date and time the merge message was received from VistA <p>NOTES</p> <p>VBECS displays  at the first instance of a VistA patient update or merge. The icon remains on the screen until the Traditional Supervisor clears the patient merge alert.</p> <p>The list of patient records is sorted in chronological order: the oldest received updated record appears at the top of the list.</p> <p>The user may update both patients' VBECS files in accordance with local policy and procedure.</p>
2. Select a VistA patient merge event from the list and review the data (Figure 115).	<ul style="list-style-type: none"> • Displays data about merged VistA patient records: <ul style="list-style-type: none"> o Patient name o Patient ID, if available o Patient sex, if available o Date and time the merge message was received from VistA o Patient ABO/Rh, if available o Patient Transfusion Requirements (TRs) and Special Instructions (SIs), if available o Recent orders and recent transfusion records for the patient (available only by clicking Patient Information Toolbar icons), if available <p>NOTES</p> <p>"Patient has no VBECS records at the receipt of the merge message" displays when the merge-from or merge-to patient is not found in VBECS.</p>
Merge-to Patients 3. Update the merge-to patient record to include SIs and TRs from the merge-from patient record. Insert a comment in the Special Instruction field to note that a merge occurred.	<p>NOTES</p> <p>A user cannot update SIs or TRs if the merge to patient does not have a VBECS record.</p> <p>Local policy must define how to handle records when the merge to patient does not exist in VBECS.</p> <p>The merge-to patient record is the active record in VBECS.</p> <p>Use Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement, and Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions to enter, update, or delete relevant component requirements, antibodies, or antigen requirements.</p> <p>Enter an SI comment in the merge-to patient record, for example, "This patient record was merged with <insert merge-from patient name> in VistA and is used to process orders after <insert merge date and time>. Refer to the records of <insert merge-from patient name> for testing and transfusion information prior to the merge event."</p>
4. Resolve ABO/Rh discrepancies	

User Action	VBECS
among patients. Determine which blood type is correct for use.	NOTES Use Justify ABO/Rh Change to finalize blood type changes for the merge-to patient.
5. Complete pending orders for the merge-to patient.	
Merge-from Patients 6. Update the merge-from patient record to include SIs and TRs from the merge-to patient record Insert a comment in the Special Instruction field to note that a merge occurred.	NOTES A user cannot update SIs or TRs if the merge from patient does not have a VBECS record. Local policy must define how to handle records when the merge from patient does not exist in VBECS. The merge-from patient record is the inactive record in VBECS. Use Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement, and Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions to enter, update, or delete relevant component requirements, antibodies, or antigen requirements. Enter an SI comment in the merge-from patient record, for example, "This patient record was merged with <insert merge-to patient name> in VistA and must not be used to process orders after <insert merge date and time>. Refer to the records of <insert merge-to patient name> for testing and transfusion information after the merge event."
7. Cancel pending orders on the merge-from patient. Request replacement orders for the merge-to patient, as appropriate.	NOTES Local policy must define how to handle orders when the merge from patient or merge to patient does not exist in VBECS.
8. To confirm that the merge-from and merge-to patient SIs and TRs match on the Patient Merge display.	NOTES Local policy must define how to handle records when the merge from patient or merge to patient does not exist in VBECS.
9. Click Clear Patient Merge Alert .	<ul style="list-style-type: none"> When there are no pending patient updates, clears  (the Patient Alert icon) from the main status bar. Displays an information message stating that the user viewed the record. NOTES >>> The Clear Patient Merge Alert function is available only to the Traditional Supervisor and above.
10. Click OK to acknowledge the message.	<ul style="list-style-type: none"> Clears the icon when the user acknowledges the message.
11. Click Close to exit, or return to Step 2.	

Figure 115: Patient Merge

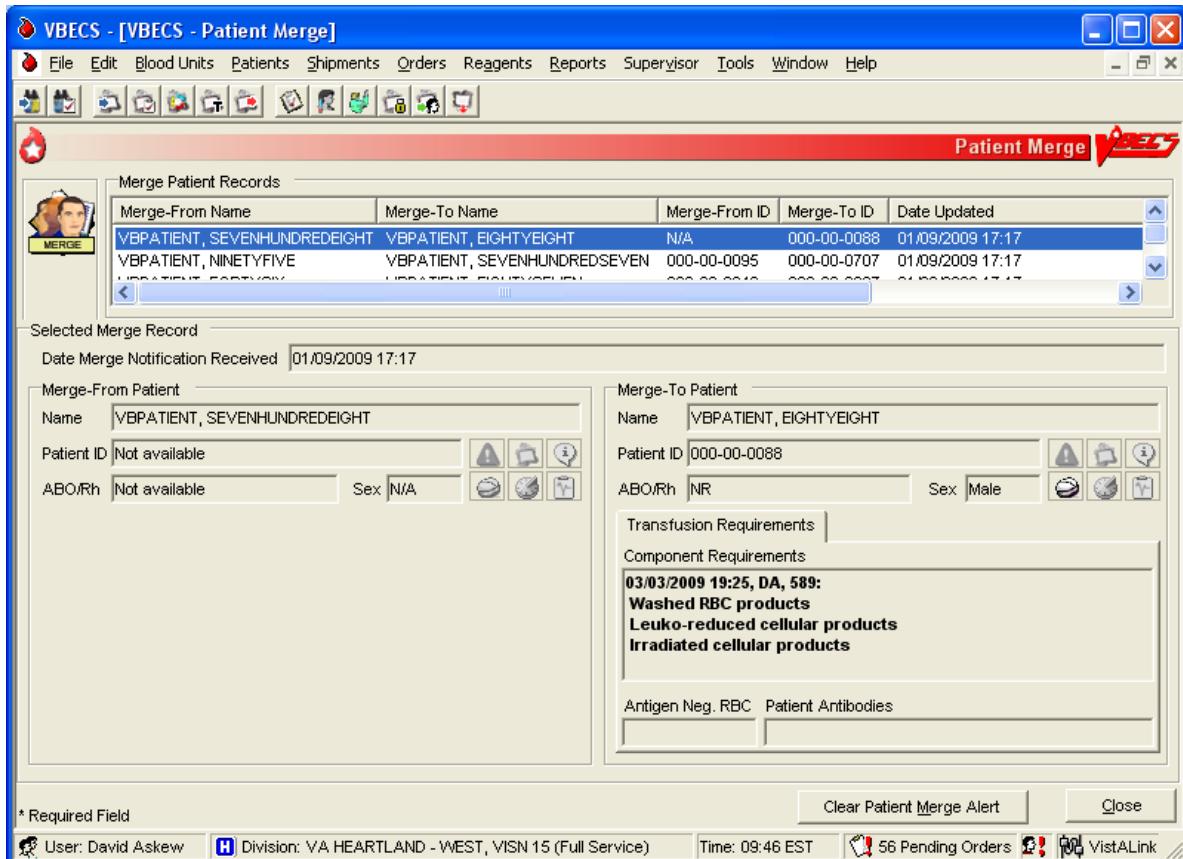

Updated Patients and Deceased Patients

The user views patient updates and death events for patients with pending or active orders in VBECS.

Assumptions

- The connection to VistA is active.
- The patient record was updated in VistA.
- VBECS received a message from VistA about a patient record update or death.
- The affected patient's records are in the VBECS database.
- The affected patient has a pending or active order in VBECS for updates and/or death event.

Outcome

- The user reviews updated patient records.

Limitations and Restrictions

- Demographic and death event updates do not display information for patients without pending or active orders in VBECS. Updates will occur without user input.
- VBECS does not notify users of additional VistA patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor completes the update and clears the icon.
- When VBECS notifies a user of an update to a patient's prefix or suffix, VBECS does not display the updated information: the user must access VistA to view the updated information.

Additional Information

- None

User Roles with Access to This Option

All users

Updated Patients and Deceased Patients

The user views updated patient records.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Updates . Select Updated Patient or Deceased Patient , or Click  in the Patient Information Toolbar, when it is active.	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Lists updated or deceased patients (also lists these data after the update if they changed), including:<ul style="list-style-type: none">◦ Patient name◦ Patient ID◦ Patient date of birth◦ Date and time of the update, or date of death
2. Click Close to exit, or return to Step 1 to view another patient's data.	

Patient Testing: General Instructions

The user records a patient's serologic test results.

Assumptions

- The patient specimen and order were processed and found satisfactory for testing.
- The user selected one to four diagnostic and/or component tasks for testing from the Pending Task List (PTL).
- The patient may or may not have a previous blood bank record. Previously saved VBECS data related to this patient are available for comparison with the user's test results and interpretations.

Outcome

- The test record was filed, is unique, and is retrievable for the life of VBECS.
- Partially completed testing is available for completion or invalidation but is not available to the patient record or reports.
- Completed, current specimen test results are available for inclusion in various reports.
- VBECS changes the task status to "completed" and triggers the generation of a message to CPRS that the order is complete.
- The test record includes the details of the verified, completed testing.

Limitations and Restrictions

- This test record is not editable after the user completes a test and saves it to the database. When a correction needs to be filed, a corrected report must be generated, as described in Invalidating Patient Test Results.
- User selection of partially completed testing grids with different phases of testing (All Phases or AHG Only) as previously selected creates a testing discrepancy in VBECS.
- Limited batch testing to partially completed orders with the same testing grid selection, either All Phases or AHG Only.
- VBECS cannot compare test results for a patient when two users are entering test results for the same patient at the same time.

Additional Information

- Partially completed tests may be invalidated within the testing option.

User Roles with Access to This Option

All users

Patient Testing: General Instructions

A user records a patient's test results: Type & Screen (TAS), ABO/Rh, Antibody Screen Test (ABS), Direct Antiglobulin Test (DAT), and crossmatches (XMs) on units selected for a patient. This option describes the configuration of reflex or repeat testing and patient antigen typing.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the Pending Task List (PTL) in the Diagnostic Tests tab. NOTES _____

User Action	VBECS
	<p>A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click one to four check boxes in the Pending Task List to select tasks for testing (Figure 111). Click OK .	<ul style="list-style-type: none"> Displays patient names and associated task information.
3. Click a check box to select a reagent rack. In the ABS or XM Testing Phases area, click the All Phases or AHG Only radio button (Figure 116). Click OK to continue with testing.	<ul style="list-style-type: none"> Allows the user to select the All Phases or AHG Only radio buttons only when applicable. Displays reagent rack names and associated information. Requires the user to select the reagent rack name. The testing grid configuration for serologic testing defaults to the setting selected in the Configure Division option. The user is allowed to indicate whether an automated instrument was used. Displays the date and time testing was performed. Displays the testing technologist name. <p>NOTES</p> <p>When the testing method selected is AHG only, VBECS enables specific phases for the ABS and XM data grid boxes for:</p> <ul style="list-style-type: none"> Antibody Screen (ABS): AHG phase Serologic crossmatch (XM): IS and AHG phases <p>Rack configuration of AHG applies only to the antibody screen and crossmatch testing grids.</p> <p>The user may select another user's name from the list, as appropriate. The logged on user's name appears on the Testing Worklist Report as the Processing Tech in addition to the selected user name.</p> <p>The user may indicate that testing was performed at a past date and time.</p> <p> See Table 10 for alerts that may occur during this option.</p>
4. Respond to warnings, which may include entering a comment and details, and click OK to continue.	<ul style="list-style-type: none">  Displays  and emits an audible alert when there are SIs or TRs for the patient. Prepares a direct data entry (DDE) grid for the user to enter the serologic reaction results for testing associated with the selected tasks. Displays testing results and entered interpretations until the user saves them.
5. Select a patient to begin entering his serologic reactions (Figure 117). Repeat this step for each selected patient until all test results and interpretations are entered in the test grid.	<ul style="list-style-type: none"> Allows the user to select one patient at a time and enter a test result for that patient. Allows only valid data entries in the grid reaction result and interpretation cells. Allows a user to save partially completed tests. When the user enters interpretation data, allows the save with or without warnings and overrides based on the test results entered for the patient and the user's security level. <p>NOTES</p> <p>The user may order a reflex test and may update the patient's</p>

User Action	VBECS
	<p>TRs and SIs.</p> <p>The ABO and Rh tests are considered one test and must be resulted together before a save is permitted.</p>
<p>6. Respond to warnings for the selected patient and confirm that appropriate test results were entered.</p> <p>Repeat Steps 5–6 for each patient until all patient work is saved.</p>	<ul style="list-style-type: none"> When a user enters a valid interpretation and saves it for a test, updates the order status to “completed.” Displays each patient’s data for review before allowing the user to save them. Saves patient testing results for retrieval from completed tests. Saves only the patient testing displayed. <p>NOTES</p> <p>The user may save part or all of one patient’s testing results at a time. When all data for a patient are saved, the patient tab is no longer visible.</p> <p>When a user saves a test with a “correction” task status, VBECS saves the results and interpretations and adds a comment to the Testing Worklist Report.</p> <p>VBECS sends an order completion message to CPRS when the user enters and successfully saves all test results and interpretations for a diagnostic test.</p> <p> May display multiple warning and override messages for the user’s response.</p>
7. When the selected patient testing is complete, click OK and Yes to confirm changes and continue to the next patient, or exit.	

Figure 116: Testing Details


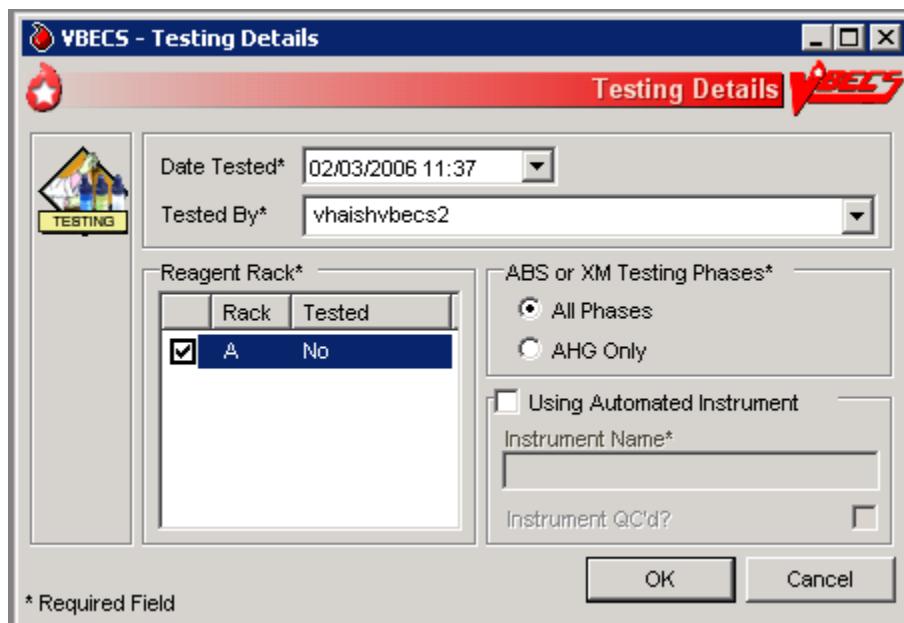


Figure 117: Patient Testing

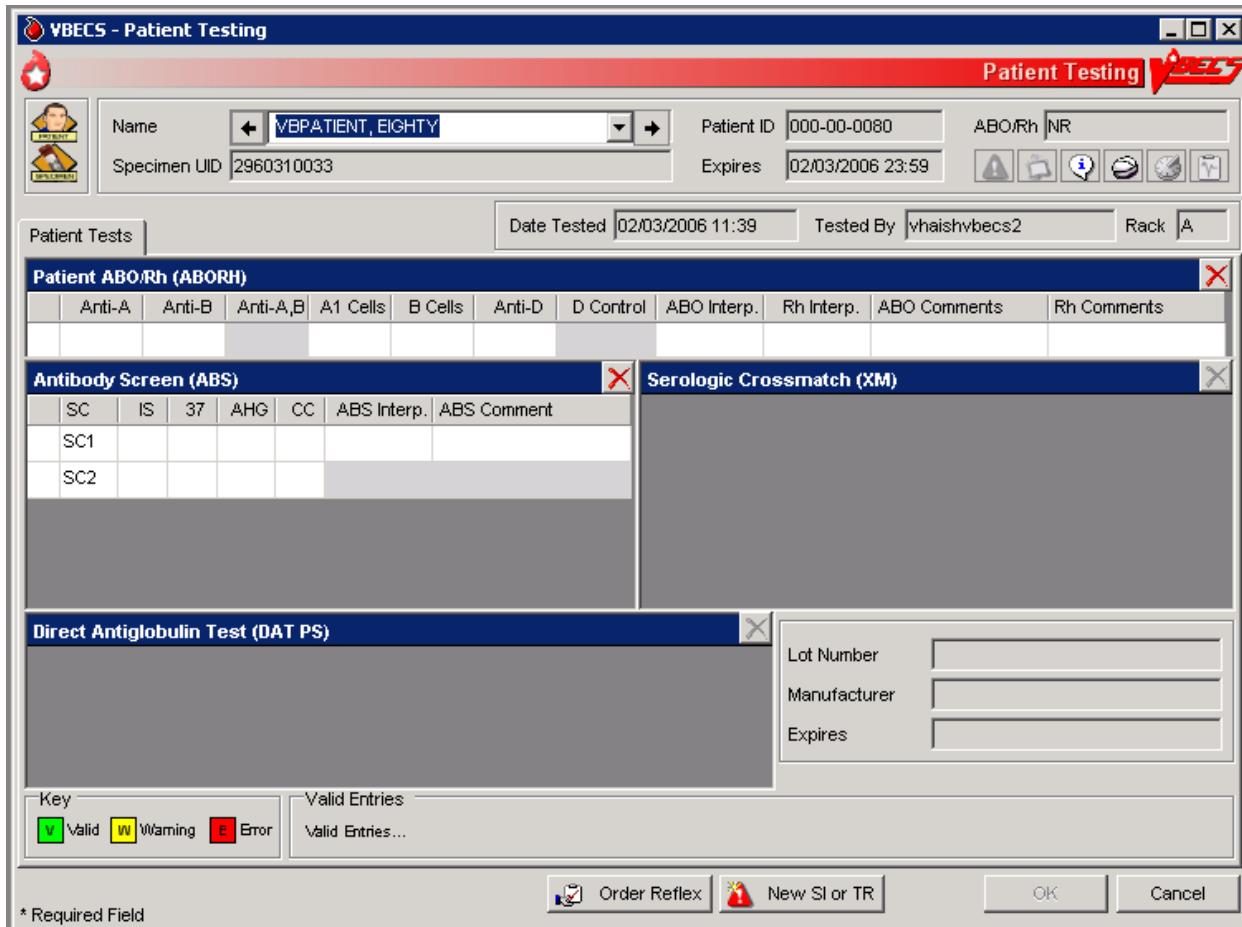

Alerts

Table 10: Alerts That May Occur in Patient Testing: General Instructions

User Action Step	Alerts
4	When the division is configured as “full service” in Configure Division, VBECS displays data fields for the user to select or enter the lot number, and displays the manufacturer and expiration date of each antiserum to be used in this transaction to allow the user to choose the correct antiserum. VBECS displays in-date antisera and allows the user to select outdated antisera, if needed, emits an audible alert, and requires an override. VBECS captures details for inclusion in an Exception Report (exception type: expired antisera used).
5	The user may cancel or delete a partially completed grid. VBECS: <ul style="list-style-type: none"> Emits an audible alert and requires a comment. Clears previously entered observed results from the data grid and allows the user to enter new data. Captures details for inclusion in an Exception Report (exception type: previously recorded results invalidated).
6	Warns that the user is saving the first instance of a historic ABO/Rh.
6	When the ABO/Rh interpretation matches observed results, is not inconclusive, and does not match the patient's historic ABO/Rh (if any, including from database conversion), VBECS emits an audible alert, warns that the current ABO/Rh does not match previous results, and instructs the user to enter a

User Action Step	Alerts
	comment to continue, or click Cancel to clear the ABO/Rh testing results and interpretation from the screen. VBECS captures details for inclusion in an Exception Report (exception type: ABO/Rh discrepancy).
6	When the interpretation of a repeat ABO/Rh does not match the original interpretation, VBECS: <ul style="list-style-type: none"> Verifies the retest results. Emits an audible alert and warns the user. Instructs the user to resolve the discrepancy before the specimen may be used. Asks whether the user wishes to continue to save. No clears the ABO/Rh testing results and interpretation from the screen. Yes requires a comment and captures details for inclusion in an Exception Report (Exception type: ABO/Rh discrepancy).
6	When a user documents an ABO/Rh discrepancy, VBECS releases units with "assigned" and "crossmatched" statuses associated with this specimen from the patient.
6	When a unit was emergency issued and a conflicting test result found that indicates a potential unit incompatibility with the patient, such as a patient ABO/Rh discrepancy, the unit XM is incompatible with the patient, the patient ABS is positive, the patient DAT is positive, or the unit is antigen positive where antigen negative is required on the tested specimen, VBECS warns the user and instructs him to notify the physician of this potential problem immediately, according to hospital policy. VBECS warns the user during the execution of the various tests related to the compatibility of a blood unit and captures details for inclusion in an Exception Report (exception type: unit ER issued, testing problem).

Patient Testing: Record a Patient ABO/Rh

The user records observed results when performing ABO/Rh blood typing.

Assumptions

- A CPRS order exists that includes an ABO/Rh typing [ABO/Rh Test or Type & Screen (TAS)].
- The user accesses patient ABO/Rh grids through Patient Testing: Record Patient Test Results.

Outcome

- The ABO/Rh result is used to determine overall blood component compatibility.
- When ABO/Rh is ordered as a stand-alone test, VBECS changes the task status to “completed” and triggers the generation of a message to CPRS that the order is completed.

Limitations and Restrictions

- None

Additional Information

- When an ABO/Rh discrepancy exists, VBECS releases assigned and crossmatched units. These units are then available for emergency issue based on system rules.
- VBECS does not enforce local policies that require running a D control for AB positive patients.
- The most recent instance of ABO/Rh testing in the database serves as the historic ABO/Rh record.

User Roles with Access to This Option

All users

Patient Testing: Record a Patient ABO/Rh

This option describes a test to determine the patient’s ABO/Rh typing from a blood specimen. The user enters the serologic reaction results in the data entry grid using the instructions and rules in this option. The grid is used for the Confirm Patient ABO/Rh reflex test as part of a TAS and crossmatch (XM) test battery. VBECS monitors the user’s entries and interpretations and maintains the records for patient reports.

These steps describe the entry of data for a patient ABO/Rh test when the test is ordered individually or as part of a TAS.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the Pending Task List (PTL) in the Diagnostic Tests tab. <p>NOTES —————</p> <p>A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click one to four check boxes in the Pending Task List to select tasks for testing. Click OK .	<ul style="list-style-type: none">• Displays patient names and associated task information.
3. Enter patient ABO/Rh testing results.	<ul style="list-style-type: none">• Displays a direct data entry (DDE) grid for entering serologic reaction results for ABO/Rh patient testing.

User Action	VBECS
	<p>NOTES —————</p> <p>The user may reenter valid results or delete entries and start over before saving the results. The user may invalidate, but not edit, data in a grid from a previous session by clicking the red "X" in the upper right corner of the testing grid, as described in Invalidate Test Results.</p> <p>VBECS determines the validity of the user-entered patient ABO/Rh interpretation by comparing the reaction results pattern of the row with the truth table and with the previous ABO/Rh entries for the patient.</p>
4. Select a patient to begin entering his serologic reactions.	<ul style="list-style-type: none"> • Allows the user to select one patient at a time and enter a test result for that patient. • Allows only valid data entries in the grid reaction result and interpretation cells. <p>NOTES —————</p> <p>When A,B or D control is not part of the daily QC configuration, VBECS disables it in the ABO/Rh test grid.</p> <p>When an AB Positive patient requires the performance of a control cell, the user must order a repeat ABO/Rh test.</p> <p>When A,B is enabled, the user enters a result for Anti-A,B as part of the patient test, it is included in determining the interpretation for that patient test. "X" is an acceptable entry.</p> <p>A,B and D control are always enabled for the repeat ABO/Rh test grid, regardless of the QC configuration.</p> <p>Reflex testing grid fields are available for data entry regardless of QC rack configuration.</p> <p>See Table 4: Valid Interpretations. Row validation is established for the interpretation vs. reaction results pattern defined by system rules.</p>
5. Repeat Steps 2–4 until some or all portions of the test grid for the selected patient are entered. Enter the interpretation for the patient ABO/Rh typing, if appropriate.	<ul style="list-style-type: none"> • Prompts the user to confirm that all results were reviewed and are acceptable. <p>NOTES —————</p> <p>When the ABO/Rh interpretation matches observed results, is not inconclusive, and does not match the patient's historic ABO/Rh (if any, including from database conversion), VBECS emits an audible alert, warns that the current ABO/Rh does not match previous results, and instructs the user to enter a comment to continue, or click Cancel to clear the ABO/Rh testing results and interpretation from the screen. VBECS captures details for inclusion in an Exception Report (exception type: ABO/Rh discrepancy).</p> <p> >> VBECS allows the Enhanced Technologist to enter an ABO/Rh interpretation that matches the patient's previously recorded (historic) ABO/Rh interpretation when the system</p>

User Action	VBECS
	interpretation is inconclusive. This allows unit selection to continue based on the patient's ABO/Rh interpretation entered by the user. VBECS requires the user to document the findings of any offline investigation, alerting the reviewer to appropriately follow up, according to local policy. VBECS emits an audible alert and warns the user. VBECS allows an override with comment and captures details for inclusion in an Exception Report (exception type: discrepant ABO/Rh override) that includes all interpretations and historic patient information.
6. When all required fields are complete, click OK to save and continue to the next patient, or exit.	<ul style="list-style-type: none"> • Saves the data.

Patient Testing: Record a Patient Antibody Screen

The user records the serologic results when performing the Antibody Screen Test (ABS).

Assumptions

- The user accesses patient ABO/Rh grids through Patient Testing: Record Patient Test Results.
- A CPRS order exists that includes an ABO/Rh typing [ABO/Rh Test or Type & Screen (TAS)]. An ABS may also be ordered as a reflex patient test.
- Configure Daily QC and Configure Testing are complete.

Outcome

- When an ABS is ordered as a stand-alone test, VBECS changes the task status to “completed” and triggers the generation of a message to CPRS that the order is complete.

Limitations and Restrictions

- This option is for patient testing only.

Additional Information

- The user performs lab tests and enters data simultaneously.

User Roles with Access to This Option

All users

Patient Testing: Record a Patient Antibody Screen

The user tests the patient’s specimen for the presence of antibodies-to-irregular red cell antigens. The user may perform this test alone or as part of a battery of tests.

The user accesses patient antibody screen grids through Patient Testing: Record Patient Test Results and enters observations of the serologic reactions in the data entry grid during testing. VBECS makes the interpretation available to CPRS and patient reports.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the Pending Task List (PTL) in the Diagnostic Tests tab. <p>NOTES —————</p> <p>A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click one to four check boxes in the Pending Task List to select tasks for testing. Click OK .	<ul style="list-style-type: none">• Displays patient names and associated task information.
3. Click a check box to select a reagent rack. In the ABS or XM Testing Phases area, click the All Phases or AHG Only radio button. Click OK to continue with testing.	<p>NOTES —————</p> <p>When the testing method selected is AHG only, VBECS disables the data grid boxes for phases IS, 37, and CC.</p>

User Action	VBECS
4. When processing more than one patient, select a name from the drop-down list in the Name field (or click the right or left arrow) to record testing for that patient.	<ul style="list-style-type: none"> Allows the user to select a patient and enter a test result for the patient. Allows only valid data entries in the grid reaction result and interpretation cells. <p>NOTES —————</p> <p>The logic described in Appendix B: Table 15: Antibody Screen Test Interpretation is an integral part of the configuration of the reaction result grid developed for Patient Testing: Record a Patient Antibody Screen.</p> <p>There may be two to five rows of screening cells to be validated as a group (the ABS). The user enters only one interpretation for the combination of cells tested.</p> <p>Valid interpretations of an ABS Test are positive or negative.</p>
5. Repeat Step 4 until results are entered in some or all portions of the test grid for the selected patient. Enter the interpretation for the patient ABS Test, if desired.	<ul style="list-style-type: none"> Verifies that appropriate tests are entered. Prompts the user to confirm that results were reviewed and are acceptable. <p>NOTES —————</p> <p>When the ABS Test interpretation must be consistent with observed test results based on Table 15: Antibody Screen Test Interpretation, discrepancies must be resolved before VBECS verifies the results and adds them to the database. The user may reenter valid results or delete the entries and start over before saving.</p>
6. Enter test results in preparation for updating the database.	<ul style="list-style-type: none"> Saves the data.
7. Click OK to save, and OK to confirm and exit.	

Patient Testing: Record a Direct Antiglobulin Test

The user records observed results when performing a Direct Antiglobulin Test (DAT).

Assumptions

- A CPRS order exists that includes a DAT, which the user may also order as a reflex patient test.
- The user is testing a patient specimen in Patient Testing: Record Patient Test Results.
- VBECS is configured to accommodate the DAT as part of a Type & Screen (TAS).

Outcome

- When a DAT is ordered as a stand-alone test, VBECS changes the task status to “completed” and triggers the generation of a message to CPRS that the order is complete.

Limitations and Restrictions

- VBECS does not check whether the lot and vial numbers selected for the positive control cells are different from those of the negative control cells.

Additional Information

- The site settings determine which reagent types and lot numbers are available.
- VBECS checks the rack for the DAT reagents in use. When the rack testing does not include Poly AHG, Anti-IgG, or anti-complement, VBECS displays additional positive and negative control rows with the test grid for the DAT sera.
- The user may select reagent lot numbers for positive and negative control cells: Reverse ABO Typing Cells (Set), A2 cell (vial), Screening Cells (set), Check Cells (vial), Panel (set), and Other (set).

User Roles with Access to This Option

All users

Patient Testing: Record a Direct Antiglobulin Test

The user performs a DAT on a patient sample.

The DAT may generate additional reflex test orders such as elution and antibody identification. VBECS saves the reactions and interpretation in the database as part of the patient’s record.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the Pending Task List (PTL) in the Diagnostic Tests tab. <p>NOTES —————</p> <p>A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click one to four check boxes in the Pending Task List to select tasks for testing. Click OK .	<ul style="list-style-type: none">• Displays patient names and associated task information.
3. When QC was not performed, enter a lot number for the antiserum.	<ul style="list-style-type: none">• Displays reaction result grids for DAT tests for each patient.• Lists the AHG antisera related to the test.• Displays information about the lot number, including the manufacturer

User Action	VBECS
	<p>and expiration date.</p> <ul style="list-style-type: none"> Prompts the user to select the correct antisera and possibly the control cell entry. <p>NOTES</p> <p>Reagent lot numbers that are selectable for positive and negative control cells:</p> <ul style="list-style-type: none"> Reverse ABO Typing Cells (set) A2 cell (vial) Screening Cells (set) Check Cells (vial) Panel (set) Other (set) <p>VBECS lists the AHG antisera related to the test from previously entered inventory. When VBECS displays multiple DAT tests or the full battery, the user selects each DAT serum separately by repeating this step until all are selected.</p> <p>Each positive and negative control cell requires the user to select a lot number.</p> <p>VBECS displays options to select the lot number of the antisera and check cells used for QC of the antiglobulin sera. The user must select complement check cells for Anti-Complement AHG tests.</p> <p> VBECS compares the control cell results entered to those of the previous day for the same rack number and alerts the user when there is a decrease in reactivity of 2 or more within the acceptable limits for each reagent type by lot number. VBECS emits an audible alert, warns that there is a decrease in reagent reactivity of two or more, and asks whether the user wishes to continue to use this vial.</p> <p>Yes allows the user to continue, with comment. VBECS captures details for inclusion in an Exception Report (exception type: QC decrease reagent reactivity ≥ 2).</p> <p>No allows the user to reenter the lot number for this reagent type and allows the approval of this same lot number (different vial) or the entry of a new lot number, returning to the testing worksheet where the user left off. (The user may replace an individual vial that is losing its reactivity with one of the same lot number.)</p> <p>When the user-entered test results indicate a satisfactory interpretation, VBECS indicates that testing is satisfactory and allows the user to save the results. VBECS considers the affected QC incomplete until it is retested with satisfactory results. VBECS generates warnings through other options when QC testing is incomplete.</p> <p> VBECS compares the expiration date of the reagent on the selected reagent to the testing date. When the reagent is expired, VBECS emits an audible alert, warns the user and asks whether he still wishes to use this reagent. Yes requires a</p>

User Action	VBECS
	<p>comment and captures details for inclusion in an Exception Report (exception type: expired reagent QC'd).</p> <p>No removes the unsatisfactory reagent and requires the user to enter a new one. Unsatisfactory lot numbers are never displayed. VBECS captures details for inclusion in an Exception Report (exception type: Expired reagent QC'd).</p> <p>When the antisera QC is required, the user verifies the lot numbers, manufacturer name, and expiration date and saves them as part of the QC record for the associated test.</p>
4. Review and accept the selected reagent information.	<ul style="list-style-type: none"> Displays the previously entered information for user review.
5. Respond to warnings and continue with DAT testing.	<ul style="list-style-type: none"> Displays a direct data entry (DDE) grid for the user to enter serologic reaction results for DAT testing.
6. Enter valid test serologic reactions in the data entry grid.	<ul style="list-style-type: none"> Allows the user to enter valid data and test results in the interpretation cells. <p>NOTES —————</p> <p>Row validation is established for the interpretation vs. reaction results pattern as defined by system rules.</p>
7. Repeat Step 4 until some or all portions of the test grid for the selected patient are entered. When appropriate, enter the interpretation for the patient DAT testing.	<ul style="list-style-type: none"> Prompts the user to confirm that all results were reviewed and are acceptable.
8. Enter test results in preparation for updating the database.	
9. Click OK to exit.	

Patient Testing: Record a Crossmatch

The user records observed results and interpretations when performing a serologic crossmatch (XM) between a patient and one or more units of blood.

Assumptions

- A CPRS order exists for a blood component that requires a serologic XM. The patient order and the selected units are in the same division as the user.
- A blood unit that contains red blood cells, selected in Select Units, is acceptable for testing, according to the patient's current or historic records and ABO/Rh.
- A current specimen was accepted and is in process for ABO/Rh and Antibody Screen Test (ABS), or was previously completed and is not expired.

Outcome

- When patient Transfusion Requirements (TRs), if any, are fulfilled and a user enters a compatible XM interpretation, VBECS makes units available for generating the Blood Transfusion Record Form (BTRF), for printing the Caution Tag, and for issue to the patient.
- When a user enters an unacceptable XM interpretation, VBECS releases the selected units from the patient.

Limitations and Restrictions

- The specimen associated with testing must have an expiration date and time that are later than the testing date and time entered in Patient Testing: Record Patient Test Results. When this is not the case, VBECS warns that the patient specimen expired. There is no override.
- This option does not include electronic crossmatch (eXM).

Additional Information

- When the user saves a partially completed XM testing worklist, VBECS does not clear patient entries without results or release them from the worklist, which VBECS saves for completion later.
- When VBECS saves a partially completed XM worksheet, it maintains a status of "assigned" for selected units until the serologic XM is complete. When VBECS saves the XM test, including interpretation, to the database, it assigns "crossmatched" status to the unit. The user may correct verified patient data only through Invalidate Test Results.
- Units that require antigen typing and/or modification to fulfill a component requirement and/or Sickle Cell testing may be crossmatched pending the completion of further processing of the unit.
- VBECS does not allow users in a transfusion-only facility to save partially completed worksheets. The user must enter interpretations for all selected units before VBECS enables the OK button.
- XM compatibility is one aspect of blood unit compatibility. VBECS takes into account the known requirements for compatibility is taken into account when the user prints Caution Tags and BTRFs, and includes such issues as antigen negative requirements and component requirements.

User Roles with Access to This Option

All users

Patient Testing: Record a Crossmatch

During Select Units, the user selects a blood unit to fill a patient order. If the unit selected contains red blood cells and eXM is not enabled, or the patient or unit does not qualify, VBECS directs the user to enter serologic XM results. (When eXM is enabled, the user may select serologic XM.)

The user performs the serologic XM test and records results in VBECS. VBECS monitors the user's entries and interpretations to determine (with the user) compatibility for each crossmatched unit.

User Action	VBECS
1. Select Patients from the main menu and select Patient Testing , or Select Blood Units from the main menu and select Select Units .	<ul style="list-style-type: none"> Displays options for processing patient-related functions. Displays the Pending Task List (PTL) and PTL search parameters in the Diagnostic Tests tab or in the Component Orders tab.
2. Enter or select PTL search parameters, if appropriate. Click the check boxes in the PTL to select one to four component orders. Click OK to continue.	<ul style="list-style-type: none"> Displays information for each component order. When the selected order class is RED BLOOD CELLS and the division enabled eXM, evaluates the patient for eXM unit processing. <p>NOTES —————</p> <p> See Table 11 for alerts that may occur during this option.</p>
3. Click Add Units to add units to the component orders selected (Figure 94).	<ul style="list-style-type: none"> Searches the database and displays available restricted autologous and directed units that the user must select before selecting allogeneic units of the same ICCBBA component class. Displays the unit ID, product type, unit ABO/Rh, unit expiration date, specimen UID, if applicable, of units assigned to the patient in the division's blood bank on a different specimen or order. The user may release and reassign the units to the current order. When eXM was configured for the facility, displays the eXM status (eligible or not eligible) for the patient. <p>NOTES —————</p> <p> See Table 11 for alerts that may occur during this option.</p> <p>When no valid ABO exists for the patient, the user may select components, including autologous units:</p> <ul style="list-style-type: none"> For RBC or WB order: <ul style="list-style-type: none"> Group O, Rh positive or O Rh negative RBC WHOLE BLOOD is not selectable For FFP order only, group AB, Rh positive or Rh negative, or not specified. For PLT, CRYO, or OTHER order, all available units are selectable.
4. Scan or enter allogeneic units, select compatible units from the Order Group list, or click the find button to select units. Enter or select a date in the Selection Date field. Each unit may have its own selection date and time. Click Add Units to display the	<ul style="list-style-type: none"> Maintains a list of selected units. The user may deselect selected units. Evaluates the selected unit: <ul style="list-style-type: none"> For the ABO/Rh compatibility of each unit and warns the user. Against patient Transfusion Requirements (TRs) and warns the user. For unit antigen negative compatibility and warns the user. Performs additional checks that include days remaining until the unit's expiration date, the current status, assignments to other patients, quarantine and biohazard indicators, when applied, and associated special testing.

User Action	VBECS
unit selection form. To change the Assigned Date and Assigned By fields, click a unit in the Order Group list.	<p>NOTES _____</p> <p> See Table 11 for alerts that may occur during this option.</p> <p>The user may accept or edit the current date and time displayed; the date and time must be in the past.</p> <p>Restricted (autologous or directed) unit selection eligibility criteria:</p> <ul style="list-style-type: none"> • Restricted For Patient ID: must match • Restricted For Patient Name: must match • Unit expiration date must be in the future. • Component class (transfusion only requires restriction of all units) <p>The criteria for selecting or including an individual unit in a pick list include evaluating a unit for:</p> <ul style="list-style-type: none"> • Considered substitutable for the order • Same division as the user • Current location is in the blood bank • Has no quarantine indicator • Is not already assigned or crossmatched to the patient • Is not “restricted” for a different patient • Meets ABO/Rh compatibility requirements specific to the component class • Future expiration date <p>VBECS lists only ABO/Rh compatible units for selection.</p> <p>When ABO/Rh results were not created for the specimen associated with the component order and the order is not designated “emergency issue,” VBECS allows the user to select units using historic records, and notifies the user that ABO/Rh testing must be performed at this division on this patient before any selected units can be issued.</p> <p>VBECS displays and allows the user to complete a pending XM associated with a patient order without a final unit status.</p>
5. Repeat Step 4 until all units are selected.	<ul style="list-style-type: none"> • Requests the user to confirm the list of selected blood units. <p>NOTES _____</p> <p>The user may select or deselect additional units.</p>
6. Click OK to confirm the selection of units.	<ul style="list-style-type: none"> • Determines whether the units are available for issue or whether they require additional modification and/or testing based on the component class. • Updates the database, as appropriate. <p>NOTES _____</p> <p> See Table 11 for alerts that may occur during this option.</p> <p>When eXM is enabled at the division and blood units require XM, VBECS determines whether the patient and units are eligible for</p>

User Action	VBECS
	<p>eXM. All parameters listed in Appendix B: Table 23: Rules for Electronic and Serologic Crossmatch must be met. Units not eligible for eXM must have serologic XM performed.</p> <p>When eXM is not enabled, the user must enter serologic XM results for red blood cell classes or when a user orders an optional XM for OTHER component classes. VBECS does not change the unit status to "crossmatched," make units available for issue, or allow the Caution Tag or BTRF to be printed until a user enters a valid XM result, unless an emergency-issued unit is being processed.</p>
7. Select another component order to continue selecting units, or exit.	<ul style="list-style-type: none"> Displays orders for the next patient.
8. Select a unit to begin entering serologic reactions. Enter valid text, then press the Enter key to move to the next data grid comment cell.	<ul style="list-style-type: none"> Allows only valid data entries in the grid reaction result and interpretation cells. Restricts data entry to valid selections for testing results. All testing results and entered interpretations remain on the screen. <p>NOTES</p> <p>During XM result data entry within a given row, the IS phase result must be entered first. Once the IS phase has been entered, data entry may progress and an observed result or NT entered in each phase before the XM interpretation is validated by VBECS.</p> <p>XM interpretations must be consistent with observed test results based on system rules for the testing method. Discrepancies must be resolved before VBECS can verify results and add them to the database. The user may reenter valid results or delete entries and start over.</p> <p>If the patient has a history of a clinically significant antibody or a persistent antigen negative requirement, VBECS warns that the patient is not eligible for IS XM and instructs the user to perform the antiglobulin transaction XM on the patient. There is no override. The user must complete all phases of the XM test.</p>
9. Repeat Step 3 until the test grid for the selected patient is complete. Enter the interpretation for the patient XM test, if desired. Respond to warnings and enter comments, when indicated.	<ul style="list-style-type: none"> Verifies that all appropriate test results have been entered for the patient and the selected unit(s). Displays various warnings based on rules. Prompts the user to confirm that results were reviewed and are acceptable. <p>NOTES</p> <p> When VBECS emits an audible alert and requires a comment, it captures details for inclusion in an Exception Report (exception type: inconclusive crossmatch).</p> <p>A user may enter a "Compatible" XM interpretation only after entering an observed test result for the IS phase (other phases are NT) when the antibody screen on the current specimen is complete and negative <i>and</i> when the patient has no calculated or persistent antigen negative requirement.</p> <p>When observed test result patterns are invalid, VBECS warns that test results are invalid and that the user must repeat tests.</p>

User Action	VBECS
	<p>VBECS stores these results as documentation and clears the grid for the user to begin again.</p> <p> VBECS emits an audible alert and requires a comment documenting medical director approval, according to site policy. VBECS captures details for inclusion in an Exception Report (exception type: Crossmatch Incompatible: Give Only with Medical Director Approval).</p>
10. Enter test results in preparation for saving.	<ul style="list-style-type: none"> Saves data. Prompts the user to generate a BTRF and a Caution Tag during Patient Testing: Record Patient Test Results when the crossmatched unit is acceptable for issue.
11. Click OK to exit.	

Crossmatch Interpretations

In all cases, VBECS saves the crossmatch information to the database.

Table 11: Crossmatch Interpretations for Print Unit Caution Tag & Transfusion Record Form

User Entry	Crossmatch Interpretation	VBECS Updates the Unit Status to:	Is Information Available to CPRS?	May User Print the BTRF and Caution Tag?	Other
C	Crossmatch Compatible	Crossmatched	Yes (unit is available)	Yes	
D	Crossmatch Compatible: Don't Transfuse	Available	No	No	
G	Crossmatch Incompatible: Give Only with Medical Director Approval	Crossmatched	Yes (unit is available)	Yes	
I	Crossmatch Incompatible	Available	No	No	
V	Inconclusive	Available	No	No	

Patient Testing: Record a Patient Antigen Typing

The user records observed results when performing a patient antigen typing test other than the ABO/Rh test.

Assumptions

- Patient antigen typing is ordered through Order Reflex Tests.
- The user accesses patient antigen typing grids through Patient Testing: Record Patient Test Results.
- Reagent inventory was entered in Maintain Minimum Levels.

Outcome

- The testing record includes the patient's antigen typing and the control cell information.

Limitations and Restrictions

- Antigen typing applies only to red blood cell antigens. (This option does not address platelet, HLA, or IgA antigens.)
- VBECS does not check:
 - The patient antigen typing with the specificity of any patient's antibodies entered.
 - Whether the lot and vial numbers selected for the positive control cells differ from those of the negative control cells.
- The patient antigen typing result is not available to CPRS.

Additional Information

- Clinical significance and availability of typing reagents are not implied when antigens are included in the list of antigen type tests.
- Unit antigen tests may be added through Incoming Shipment and Edit Unit Information.
- The patient antigen typing test is associated with a CPRS patient order.
- The patient's antigen typing information does not appear as part of the Transfusion Requirements or Special Instructions in the Patient Information Toolbar.
- When VBECS displays lot numbers for selection, it checks the inventory for the reagent name and lot numbers in use in a rack. VBECS does not offer reagents with an "unsatisfactory" value in the Inspection field or with a quantity of zero as a selectable reagent. The user may not enter a lot number that VBECS disallowed on the selection list. The user must choose a different reagent.
- When a user selects a reagent rack that was never QC'd, VBECS warns the user, clears the selection, and allows the user to select a different rack.
- VBECS automatically adds the positive and negative control tests to the selected reagent type worklist if they have not been tested on that calendar date for the selected antisera lot number.
- VBECS stores the lot number, vial identifier, manufacturer, and expiration date of the positive and negative control cells with the testing record.
- The user may not enter "H" (hemolysis) in antigen typing testing.

User Roles with Access to This Option

All users

Patient Testing: Record a Patient Antigen Typing

This option describes a test to determine the patient antigen typing from a blood specimen.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none"> • Displays options for processing patient-related functions. • Displays the Pending Task List (PTL) in the Diagnostic Tests tab. <p>NOTES ——————</p> <p>A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click one to four check boxes in the Pending Task List to select tasks for testing. Click OK .	<ul style="list-style-type: none"> • Displays patient names and associated task information. • Displays a direct data entry (DDE) grid for entering serologic reaction results for antigen typing testing. This grid may be partially populated from a previous session. • Allows the user to select a unit and enter a test result for the unit. • Displays options for finding and selecting tests to process. <p>NOTES ——————</p> <p>When the user selects a partially completed worksheet, he must verify the lot numbers, manufacturer name, and expiration date of the antiserum and control cells previously selected. These fields may not be edited. (A user who finds a discrepancy must delete the worksheet and start over.)</p>
3. Confirm the selected patient. When control cells were tested for the antiserum, go to Step 7.	<ul style="list-style-type: none"> • Displays an option to select a lot number for each of the antiserum types selected. <p>NOTES ——————</p> <p>VBECS displays grids and quality control cells for each antiserum type.</p>
4. Select a lot number for each antiserum specificity.	<ul style="list-style-type: none"> • Allows the user to select or enter the lot number and displays the manufacturer and expiration date of each antiserum to be used in this transaction so the user can select the correct antiserum. VBECS displays in-date antisera. When the user selects outdated antisera, VBECS emits an audible alert and requires a comment and an override. VBECS captures details for inclusion in an Exception Report (exception type: expired antisera used). • Prompts the user to select the correct antiserum entry for the specificity selected. • Displays the lot number information. <p>NOTES ——————</p> <p>VBECS requests information about an antigen lot number only the first time it is used on a specific day.</p> <p>If there are multiple test tabs, VBECS displays a message that may not correspond to the selected tab, but will correspond to information required by one of the requested tests.</p>
5. Verify the antiserum lot number selected (Figure 118).	<ul style="list-style-type: none"> • Allows the user to select the testing phases for the antiserum specificity by lot number. • Displays a testing grid and disables the phases that are not to be completed. • Stores the antiserum information and checks whether the positive and negative controls are needed. <p>NOTES ——————</p>

User Action	VBECS
	<p>The user must indicate the phases of reactivity for each antiserum specificity:</p> <ul style="list-style-type: none"> • IS • IS/RT • IS/37 • AHG/CC <p>VBECS enables the Rack Identifier field when the user selects the weak D or AHG/CC phase. VBECS disables the rack identifier field when the user selects other phases.</p> <p>When the user QC'd the antigen typing reagent lot number and performed patient antigen typing on the same date under different testing phases, VBECS enables the QC grid for the antisera in the phases the user selected.</p> <p> When the user attempts to change the phases of reactivity previously recorded for the lot number, VBECS emits an audible alert, warns that the user selected phases that were not used previously for this reagent, and asks whether the user wishes to continue. No requires the user to select the correct testing phases. Yes requires the user to enter a comment explaining why he is selecting different testing phases. VBECS captures details for inclusion in an Exception Report (exception type: antigen testing phase change). VBECS requires a comment and captures details for inclusion in an Exception Report. Antigen testing phase change</p>
6. Select the lot number for the reagent and the positive and negative control cells. Enter the vial identifiers, if required. Enter "0" when not using panel cells. Verify the reagent type and the positive and negative control cell lot numbers selected and vial identifiers entered, if required (Figure 119).	<ul style="list-style-type: none"> • Displays the lot number information, including the vial identifier. <p>NOTES —————</p> <p>Selectable reagent lot numbers for positive and negative control cells are:</p> <ul style="list-style-type: none"> • Reverse ABO Typing Cells • A2 cell • Screening Cells • Check Cells • Panel • Other (set) <p>The vial identifier, comprising one or two digits, identifies the vial contained in a kit of reagent red blood cells. A vial contains antigen positive or antigen negative blood cells to be used in quality control of the antisera.</p>
7. Respond to warnings and continue with antigen typing.	<ul style="list-style-type: none"> • Displays a DDE grid for entering serologic reaction results for antigen typing testing. This grid may be partially populated with results from a previous session. <p>NOTES —————</p> <p>VBECS organizes the worklist for the grids by antiserum types.</p> <p>The row title of the data grid specifies the positive and negative control cells; the column title specifies the phases of reactivity that may be used in the test. VBECS fills in this information; the user may not edit it.</p>

User Action	VBECS																																	
8. Select one control cell or patient to begin entering serologic reactions and interpretations.	<ul style="list-style-type: none"> Allows the user to select a patient or control cell and enter test results. <p>NOTES</p> <p>Valid observed test results include:</p> <table border="1"> <thead> <tr> <th>Entry</th><th>VBECS Translation</th><th></th></tr> </thead> <tbody> <tr> <td>1</td><td>1+</td><td>Positive</td></tr> <tr> <td>2</td><td>2+</td><td>Positive</td></tr> <tr> <td>3</td><td>3+</td><td>Positive</td></tr> <tr> <td>4</td><td>4+</td><td>Positive</td></tr> <tr> <td>W</td><td>Weak</td><td>Positive</td></tr> <tr> <td>F</td><td>Mixed Field</td><td>Positive</td></tr> <tr> <td>M</td><td>Microscopic</td><td>Positive</td></tr> <tr> <td>X</td><td>Not Tested X (Not Tested) (does not indicate a disabled cell)</td><td>Not Tested</td></tr> <tr> <td>0</td><td>No Agglutination</td><td>Negative</td></tr> <tr> <td>R</td><td>Rouleaux, No Agglutination</td><td>Negative</td></tr> </tbody> </table> <p>System rules establish row validation for the interpretation vs. reaction results pattern [for all user-indicated phase(s) of reactivity] for blood unit and quality control tests.</p> <p>When completing previously saved partially completed testing, the user may not change the antiserum lot number, manufacturer name, or expiration date. The user must invalidate the antiserum testing and start over.</p>	Entry	VBECS Translation		1	1+	Positive	2	2+	Positive	3	3+	Positive	4	4+	Positive	W	Weak	Positive	F	Mixed Field	Positive	M	Microscopic	Positive	X	Not Tested X (Not Tested) (does not indicate a disabled cell)	Not Tested	0	No Agglutination	Negative	R	Rouleaux, No Agglutination	Negative
Entry	VBECS Translation																																	
1	1+	Positive																																
2	2+	Positive																																
3	3+	Positive																																
4	4+	Positive																																
W	Weak	Positive																																
F	Mixed Field	Positive																																
M	Microscopic	Positive																																
X	Not Tested X (Not Tested) (does not indicate a disabled cell)	Not Tested																																
0	No Agglutination	Negative																																
R	Rouleaux, No Agglutination	Negative																																
9. Enter results for selected patients and control cells.	<ul style="list-style-type: none"> Prompts the user to confirm that results were reviewed and are acceptable. <p>NOTES</p> <p> VBECS searches the current order and the patient's historical record for previously entered red cell antigen results for each specificity within the user's division and compares the results to prevent a discrepancy with previous test entries. When a discrepancy is found, VBECS emits an audible alert, notifies the user that the current specimen and a previously entered antigen typing on this patient do not match, and asks whether the user wishes to save the current result.</p> <p>No closes the form and returns to the grid. Yes allows the user to update the conflicting results in the database and/or enter a comment explaining why he is accepting this result. VBECS captures details for inclusion in an Exception Report (exception type: discrepant patient antigen typing).</p>																																	
10. Repeat Steps 5–9 for each antiserum type.	<ul style="list-style-type: none"> Allows the user to save. <p>NOTES</p> <p>The positive and/or negative control cell test results on the date of testing must be a valid positive or negative.</p>																																	
11. Respond to warnings and click OK to save the data.	<ul style="list-style-type: none"> Reconfirms that all appropriate test results are entered, repeats all data validation, and saves the data. 																																	
12. Click OK to exit.	<ul style="list-style-type: none"> Saves the data. 																																	

Figure 118: Patient Testing Antigen Typing (1)

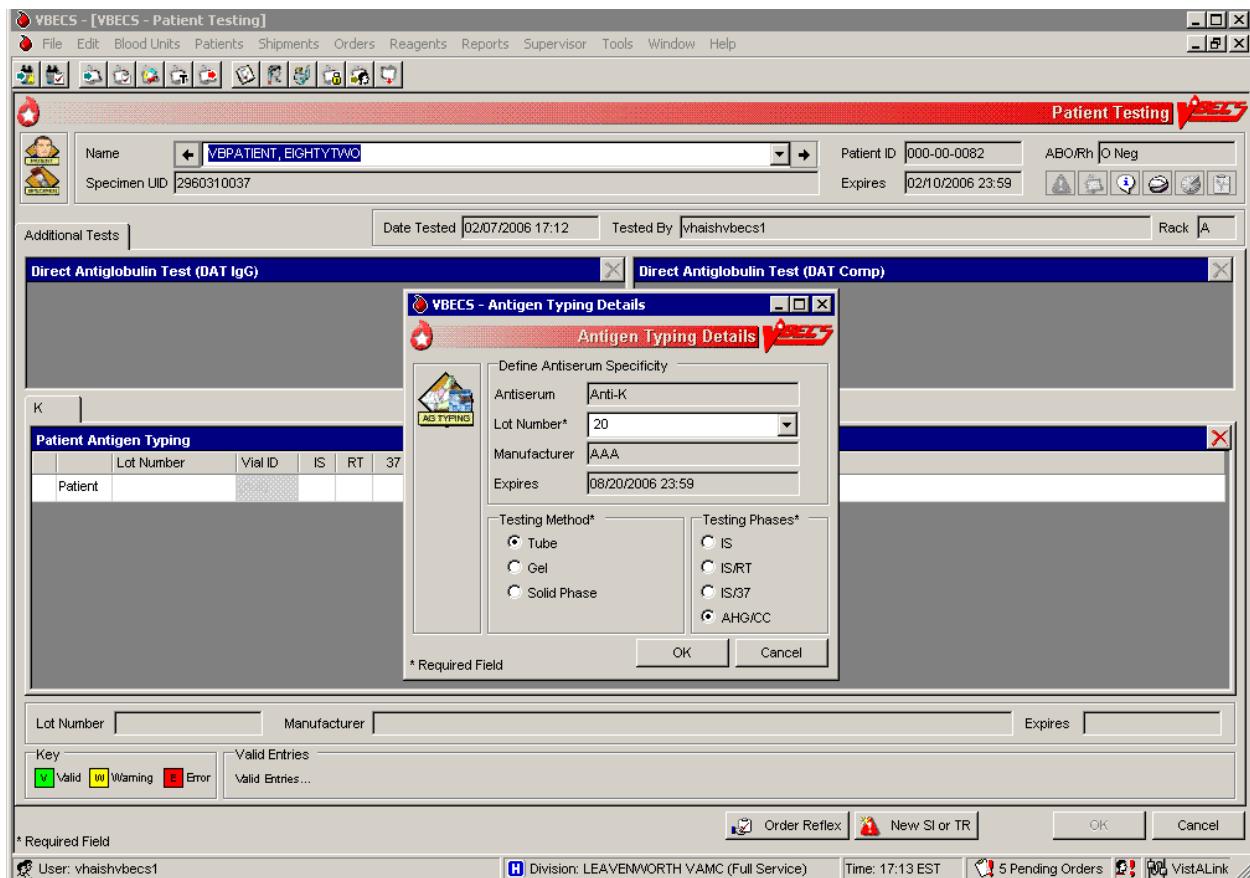


Figure 119: Patient Testing Antigen Typing (2)

The screenshot shows the VBECS - Patient Testing software interface. At the top, there is a header bar with the title "VBECS - Patient Testing" and a logo. Below the header, the patient information is displayed: Name (VBPATIENT, EIGHTYTWO), Patient ID (000-00-0082), ABO/Rh (O Neg), Specimen UID (2960310037), and Expires (02/03/2006 23:59). There are also buttons for additional tests, date tested (02/03/2006 11:48), tested by (vhaishvbecs2), and rack (A).

Two tabs are open: "Direct Antiglobulin Test (DAT IgG)" and "Direct Antiglobulin Test (DAT Comp)".

A modal window titled "Patient Antigen Typing" is displayed. It contains a table with columns: Lot Number, Vial ID, IS, RT, 37, AHG, CC, Interp., and Comment. The rows show:

	Lot Number	Vial ID	IS	RT	37	AHG	CC	Interp.	Comment
POS									
NEG									
Patient	20								

Below the table, there are fields for Lot Number, Manufacturer, and Expires. A key at the bottom left indicates: V Valid, W Warning, E Error. Buttons at the bottom right include Order Reflex, New SI or TR, OK, and Cancel. A note at the bottom left says "* Required Field".

Patient Testing: Enter Antibody Identification Results

The user records results of reflex tests needed to complete a patient blood bank antibody workup.

Assumptions

- An order for an antibody identification (ABID) was placed and is opened in the Pending Task List (PTL).

Outcome

- VBECS displays the ABID in the patient's record, which the user may review by clicking an icon in the Patient Information Toolbar.
- Entry of antibody specificity through this option may automatically set an antigen negative requirement [Transfusion Requirement (TR)] based on the division-specific parameters.

Limitations and Restrictions

- When an antigen negative requirement is not configured for an antibody, the user may enter the requirement for a patient through Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement.

Additional Information

- Correction to verified patient data may be accessed only through Invalidate Test Results.
- Only one technologist at a time can enter test results for a selected PTL task.
- Patient antigen testing is not included in this option.

User Roles with Access to This Option

All users

Patient Testing: Enter Antibody Identification Results

The user may select workload codes for work performed on-site during the ABID workup. When the ABID is performed off-site, the user enters the associated charges.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the Pending Task List (PTL) in the Diagnostic Tests tab. <p>NOTES —————— A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click one to four check boxes in the Pending Task List to select tasks for testing. Click OK .	<ul style="list-style-type: none">• Displays patient names and associated task information.
3. Click a check box to select a reagent rack. Click OK to continue with testing.	<p>NOTES —————— No serologic testing is recorded for ABID.</p>
4. Select antibodies. Enter associated testing details.	<ul style="list-style-type: none">• Displays the list of antibodies and allows the user to select from the identified specificities and to add antibodies or to indicate that no antibody was identified.

User Action	VBECS
Respond to comments and continue. Verify Patient information. Click OK .	NOTES _____ Multiple workload codes may apply when work is performed on-site.
5. Review the selections displayed. Return to Step 1 to change selections when necessary. Accept the reviewed entries.	<ul style="list-style-type: none"> • Displays applicable data fields for the ordered ABID. • Allows the user to indicate the testing site and associated billing from outside sources.
6. Enter the location where testing was performed and cost, when off-site.	<ul style="list-style-type: none"> • Requires the user to indicate whether the ABID was performed at the division and/or off-site and saves that information with the ABID. • Allows the user to enter a charge, but not the workload, when the ABID work was performed off-site. • Allows the user to indicate multiples of any processes when a test was performed on-site. • Allows the user to enter comments.
7. Review the antibody tests entered and comments selected and, when necessary, change them before accepting (Figure 120).	<ul style="list-style-type: none"> • Displays the requirements for antigen negative blood components, if any, based on the antibody entered.
8. Click OK to save the information to the database.	<ul style="list-style-type: none"> • Saves the data.
9. Click OK to exit.	

Figure 120: Patient Testing Antibody Identification

The screenshot shows the VBECS - Patient Testing interface. At the top, the menu bar includes File, Edit, Blood Units, Patients, Shipments, Orders, Reagents, Reports, Supervisor, Tools, Window, and Help. Below the menu is a toolbar with various icons. The main window displays patient information: Name (VBPATIENT_EIGHTYTWO), Specimen UID (2960310037), Patient ID (000-00-0082), ABO/Rh (O Neg), and Expires (02/10/2006 23:59). The 'Antibody ID' section lists identified antibodies, with 'Anti-K' checked and highlighted in blue. Other listed antibodies include Anti-C, Anti-E, Anti-c, Anti-e, Anti-f, Anti-G, Anti-N, Anti-S, Anti-s, Anti-U, Anti-P1, Anti-Le(a), Anti-Le(b), Anti-Js(a), Anti-Js(b), Anti-Lu(a), Anti-Lu(b), Anti-Lu3, Anti-C(w), Anti-Ce, Anti-C(x), Anti-E(w), Anti-V, Anti-Xg(a), Anti-M(g), Anti-Mi(a), Anti-P, Anti-A1, Anti-H, Anti-I, Anti-i, Warm auto-antibody, Cold auto-antibody, Anti-Le(ab), Anti-Di(a), Anti-Di(b), Anti-Do(a), and Anti-Do(b). To the right, there are sections for 'Test Performed Off-Site' (unchecked) and 'Test Performed On-Site' (checked), which includes a dropdown for 'Workload Process' (Antibody Identification panel, standa) and a multiplier input field (set to 1). At the bottom, there are buttons for Order Reflex, New SI or TR, OK, and Cancel. A note at the bottom left indicates a required field. The status bar at the bottom shows the user (vhaishvbechs1), division (LEAVENWORTH VAMC (Full Service)), time (17:16 EST), and pending orders (5 Pending Orders).

Patient Testing: Record a Transfusion Reaction Workup

The user documents a transfusion reaction workup (TRW).

Assumptions

- A TRW is ordered and is “not started” or “in progress” on the Pending Task List (PTL).

Outcome

- The user views or prints a preliminary TRW report of findings.
- Preliminary report findings are available for display in CPRS.

Limitations and Restrictions

- Adverse reactions to transfusion detected at a later time (disease transmission) are not recorded through this option. Refer to Transfusion Complications Report.

Additional Information

- All implicated units are added to data grids.
- Testing is based on local policies.
- VBECS allows the user to invalidate previously entered and saved information. When the user invalidates specimen data entered, VBECS invalidates the TRW.
- A delayed hemolytic transfusion reaction is typically identified about 7 to 10 days post-transfusion. It may be addressed in this option or as an antibody identification. (See Patient Testing: Enter Antibody Identification Results.)
- This test record is not editable after the user saves it to the database. When a correction needs to be filed, a corrected report must be generated, as described in Invalidate Test Results.
- See the printed chart copy for the medical director’s signature and TRW interpretation, as described in Finalize/Print TRW.
- A user may order a TRW prior to the completion of a unit’s transfusion record.
- The user must complete the unit’s transfusion record prior to implicating a unit in a TRW.
- These data remain encapsulated as part of the TRW and are not saved, compared to the patient record, or associated with other pre-transfusion grid testing.
- TRW ABO/Rh typing results are not saved as historical ABO/Rh typing results. These test records are maintained as part of the TRW only.

User Roles with Access to This Option

All users

Patient Testing: Record a Transfusion Reaction Workup

The user accepts the TRW order, inspects the sample for hemolysis, and performs implicated clerical checks. The TRW may require entry of patient data, including vital signs and serologic test results. A signed, printed TRW report must be filed with the patient's chart.

Transfusion reactions may be worked up in stages depending on patient symptoms and serologic findings. The first stage includes the initial clerical and specimen checks. The second stage includes serologic test interpretations, including pre- and post-patient testing, ABO/Rh, ABS, XM, and DAT.

User Action	VBECS
1. Select Patients from the main menu. Select Patient Testing .	<ul style="list-style-type: none"> Displays options for processing patient-related functions. Displays the Pending Task List (PTL) in the Diagnostic Tests tab. <p>NOTES</p> <p>A user may search for specific tasks using search criteria detailed in Pending Task List.</p>
2. Click a check box in the PTL to select a task. Click OK .	<ul style="list-style-type: none"> Displays TRW patient order information and allows the user to indicate which portion of the TRW to perform (initial clerical and specimen checks). <p>NOTES</p> <p>The TRW patient name and patient ID are visible at all times.</p> <p>Initial clerical and specimen checks must be performed before TRW serologic testing is available.</p>
3. Enter the date the transfusion reaction was noted and the date it was investigated.	<ul style="list-style-type: none"> Allows the user to document observations about the pre- and post-transfusion specimens.
4. Click a radio button in the Specimen Checks tab to indicate whether pre- and post-transfusion specimens are available (Figure 121). Scan or enter specimen UIDs in the Specimen UID fields and enter comments, as required. Select an item from the drop-down menu in the Hemolysis Inspection field.	<ul style="list-style-type: none"> Allows the user to enter the date and time the transfusion reaction was noted (provided by the clinician) and the date and time of the investigation (default: current date and time). <p>NOTES</p> <p>When the pre- and/or post-transfusion samples are available to perform the clerical check, VBECS requires the user to indicate whether the clerical check is okay.</p> <p>VBECS allows the user to scan or enter the pre- and post-transfusion specimen UIDs. VBECS compares the patient name and ID associated with each specimen UID to those associated with the TRW workup.</p> <p> When the patient name does not match VBECS data for the TRW patient, VBECS emits an audible alert, warns the user, and instructs him to notify the physician immediately. VBECS allows the user to reenter the specimen UID or verify that the data are correct. VBECS requires a comment when there is no match and captures details for inclusion in an Exception Report (exception type: Transfusion Reaction Workup Specimen).</p> <p>The user may enter "icteric," "absent," "1+," "2+," "3+," or "4+" in the Hemolysis Inspection field to identify the presence or absence of hemolysis in the pre- and post-transfusion patient</p>

User Action	VBECS
5. Click the Clerical Check OK or Clerical Check Failed radio button to perform initial clerical and identification checks (Figure 122).	specimens or blood container, when available.
6. Click OK to save, and Yes to confirm the save.	<p>NOTES —</p> <p>When the user exits and saves the data, VBECS changes the task status of the TRW to “partially completed.”</p> <p>A partially completed TRW may be corrected only by invalidating and reentering the TRW.</p>
7. Enter the date and time of the transfusion reaction noted by the clinician and the date and time of the reaction investigation.	<p>NOTES —</p> <p>VBECS allows the user to enter the date and time of the reaction noted by the clinician.</p> <p>When the implicated unit is unknown, VBECS allows the user to continue without identifying a unit or entering implicated unit information.</p> <p>VBECS requires the user to enter the date and time of the TRW investigation (default: current date and time), including a date and time in the past for retrospective data entry. The user may not enter a future date and time.</p> <p>If the user indicates that implicated units are unknown, he continues at Step 10.</p>
8. In the Implicated Units tab, click the Implicated Unit(s) Identified radio button to select an implicated unit. Scan or enter the unit ID and product code. Click the Blood Bag Returned or Blood Bank Not Returned radio button. Select an item from the drop-down menu in the Hemolysis Inspection field. Click the Clerical Check OK or Clerical Check Failed radio button to perform initial clerical and identification checks. Enter a comment and details, as required (Figure 123). Click Add to add the unit. or	<ul style="list-style-type: none"> • Allows the user to: <ul style="list-style-type: none"> ○ Enter implicated unit information. ○ Indicate whether the bag and infusion set were returned. ○ Review and/or edit the entered information. <p>NOTES —</p> <p>When the user indicates that the blood bag was not returned, VBECS disables the hemolysis and clerical checks.</p> <p>All implicated units entered are available for TRW ABO/Rh and TRW crossmatching. Actual data entry of test results for any entries is optional, based on a facility's local policies and procedures.</p> <p>The user may enter a free-text comment to document site-specific information such as other solutions or additives, urinary hemoglobin checks, and whether the product was warmed during transfusion.</p>

User Action	VBECS
<p>Click the No Implicated Units radio button.</p> <p>Click OK to save, and Yes to confirm the save.</p>	
9. Repeat Step 8 for each implicated unit.	<ul style="list-style-type: none"> Allows the user to enter the patient's symptoms. <p>NOTES</p> <p>VBECS accepts data for multiple implicated units and stores and reports them as part of the clerical check. When the user enters implicated unit data, VBECS adds the unit to TRW serologic testing.</p>
<p>10. In the Symptoms tab, click the Patient Symptoms Identified radio button.</p> <p>Click the check boxes to select one or more patient symptoms (Figure 124),</p> <p>or</p> <p>Click the No Patient Symptoms Identified radio button.</p> <p>Click OK to save, and Yes to confirm the save.</p>	<ul style="list-style-type: none"> Allows the user to review and/or edit the entered data before it saves the record. Allows the user to indicate that the TRW is complete.
<p>11. In the Serologic Tests tab, click a check box to select which tests to perform.</p> <p>Select a rack from the drop-down list in the Rack field.</p> <p>Click the Using Automated Instrument check box, as required.</p> <p>Enter the name of the instrument in the Name field.</p> <p>Click the Instrument QC'd? check box, as appropriate.</p> <p>Enter interpretations for TRW testing (Figure 125).</p>	<ul style="list-style-type: none"> Allows the user to enter one or more additional TRW test interpretations: <ul style="list-style-type: none"> TRW ABO/Rh testing TRW ABS testing TRW DAT testing TRW XM testing Requires the user to enter the rack name and test method to be applied to the specimens selected in this testing group for serologic testing. The user must also indicate whether an automated instrument was used in testing. Displays the date and time the testing was performed (default: current date and time), which the user may edit retrospectively. Displays the testing technologist's identification (default: current user). Lists valid division users for the user's selection. <p>NOTES</p> <p>When no units were implicated or no pre- or post-transfusion specimens are available, VBECS disables those portions of the grid.</p> <p>VBECS displays the standard TRW tests for the user to select to perform some, none or all.</p> <p>When a facility is "full service" and an automated instrument was used for testing, a user must enter the name of the instrument and may indicate that successful quality control testing was done for the testing date.</p>

User Action	VBECS
	<p>For implicated units, the ABO/Rh TRW, ABS TRW, DAT TRW, and Crossmatch TRW testing data grids are pre-populated with pre- and post-specimen data. Data entry of test results is optional, based on a facility's local policies and procedures.</p>
12. Review the entries. Click OK to save, and Yes to confirm the save.	<ul style="list-style-type: none"> • Stores data for retrieval. • Makes the findings viewable through a CPRS query. <p>NOTES —————</p> <p>VBECS allows the user to indicate that the TRW is complete, then updates the TRW task status to "complete."</p> <p>Once the TRW is completed, the summary of results and observations includes the comment, "Final report to follow in hard copy with medical director evaluation of the transfusion reaction workup."</p> <p>VBECS indicates on the implicated unit's patient transfusion record that it was implicated in a transfusion reaction and notes the patient's name and ID and the date and time.</p> <p> VBECS saves the record and one or more clerical checks are marked "failed," VBECS emits an audible alert, requires a comment, and captures details for inclusion in an Exception Report (exception type: Transfusion Reaction Workup).</p>
13. Click Finalize/Print TRW and continue at Finalize/Print TRW to finalize the TRW (Figure 126), or Click Cancel to complete the TRW and exit.	<p>NOTES —————</p> <p>A completed TRW indicates that serologic testing was completed. A finalized TRW indicates that the medical director reviewed the TRW, classified the reaction, and may make recommendations for future transfusions.</p>

Figure 121: Transfusion Reaction Workup: Specimen Check

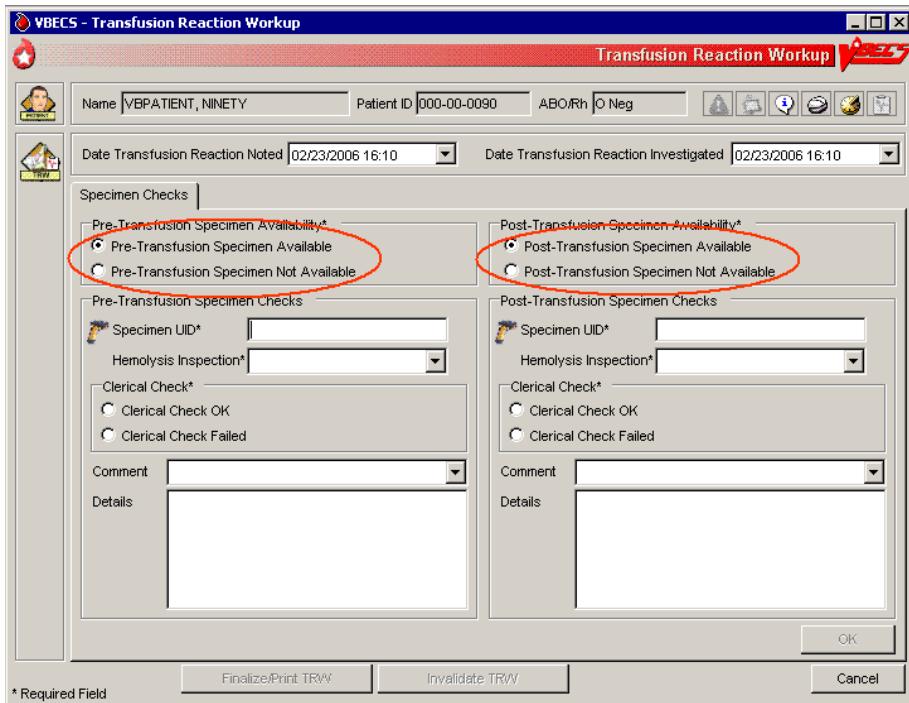


Figure 122: Transfusion Reaction Workup: Clerical Check

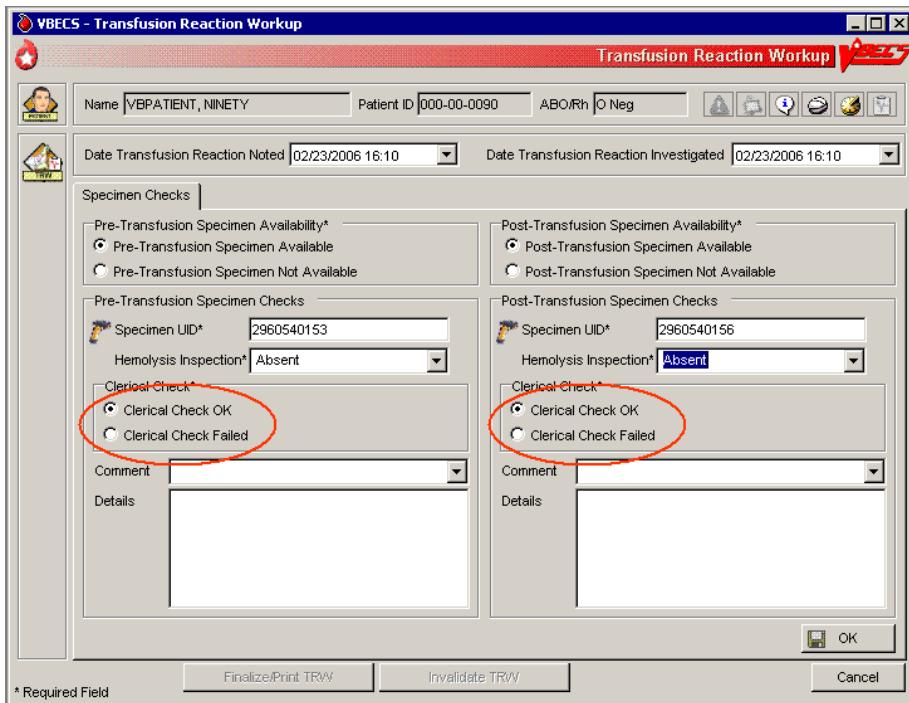


Figure 123: Transfusion Reaction Workup: Implicated Units

This screenshot shows the 'Implicated Units' tab of the Transfusion Reaction Workup interface. At the top, there are fields for Name (VBPATIENT, NINETY), Patient ID (000-00-0090), and ABO/Rh (O Neg). Below these are two date pickers: 'Date Transfusion Reaction Noted' (02/23/2006 16:10) and 'Date Transfusion Reaction Investigated' (02/23/2006 16:10). The 'Specimen Checks' tab is selected. Under 'Implicated Units?'*, the radio button 'Implicated Unit(s) Identified' is selected. In the 'Implicated Unit Checks' section, there are fields for Unit ID (6167616) and Product (04210). The 'Blood Bag Checks' section includes 'Hemolysis Inspection*' (Absent) and 'Clerical Check*' (Clerical Check OK). A 'Comment' field and a 'Details' scrollable area are also present. At the bottom, there are buttons for 'OK', 'Finalize/Print TRW', 'Invalidate TRW', and 'Cancel'. A note at the bottom left indicates '* Required Field'.

Figure 124: Transfusion Reaction Workup: Symptoms

This screenshot shows the 'Symptoms' tab of the Transfusion Reaction Workup interface. At the top, there are fields for Name (VBPATIENT, NINETY), Patient ID (000-00-0090), and ABO/Rh (O Neg). Below these are two date pickers: 'Date Transfusion Reaction Noted' (02/23/2006 16:10) and 'Date Transfusion Reaction Investigated' (02/23/2006 16:10). The 'Specimen Checks' tab is selected. Under 'Symptoms Identified?*', the radio button 'Patient Symptoms Identified' is selected. In the 'Patient Symptoms*' section, there is a list of symptoms with checkboxes. Some checked symptoms include Wheezing, shortness of breath, Flushing, and Pain at infusion site. Other symptoms listed include Urticaria (Rash), Pruritis (Itching), Dyspnea (difficulty breathing), Tachypnea (rapid respiration), Hypoxemia, Chills (rigors), Fever (>1C or 2F increase), Pain in Flanks, Pain in Chest, Pain in Abdomen, and Localized edema (angioedema). There is also a 'Other Details' field and a scrollable 'Other' section. At the bottom, there are buttons for 'OK', 'Finalize/Print TRW', 'Invalidate TRW', and 'Cancel'. A note at the bottom left indicates '* Required Field'.

Figure 125: Transfusion Reaction Workup: Serologic Tests

The screenshot shows the 'Transfusion Reaction Workup' window with the following details:

- Name:** VBPATIENT, NINETY
- Patient ID:** 000-00-0090
- ABO/Rh:** O Neg
- Date Transfusion Reaction Noted:** 02/23/2006 16:10
- Date Transfusion Reaction Investigated:** 02/23/2006 16:10
- Specimen Checks:** Implicated Units, Symptoms, Serologic Tests
- Tested By:** Test User1
- Rack:** A
- Using Automated Instrument:** Checked
- Name:** SeroSerious X422
- Instrument QC'd?**:
- Comment:**
- ABO/Rh:** ABO Interp. Rh Interp.
Pre (2960540153) O Neg
Post (2960540156) O Neg
6167616 O Neg
- XM:** Pre Interp. Post Interp.
6167616 Compatible - Sa Compatible - Sa
- ABS:** Interp.
Pre (2960540153) Neg
Post (2960540156) Neg
- DAT:** Interp.
Pre (2960540153)
Post (2960540156)

* Required Field Finalize/Print TRW Invalidate TRW OK Cancel

Figure 126: Transfusion Reaction Workup Report

The screenshot shows the 'Transfusion Reaction Workup Report' window with the following details:

- Patient Name:** VBPATIENT, NINETY
- Patient ID:** 000-00-0090
- Reaction Type:** Undetermined
- Specimen Checks:** Specimen USD, Hemolysis, Checks OK?, Comments
- Reported Symptom:** Wheezing, shortness of breath, Flushing, Anxiety
- Serologic Test:** Specimen USD, Unit ID, Product, Interp., Comments

Serologic Test Name	Specimen USD	Unit ID	Product	Interp.	Comments
Pre ABO Interp	2960540153			O	
Pre Rh Interp	2960540153			Neg	
Pre Antibody Screen Interp	2960540153			Neg	
Post ABO Interp	2960540156			O	
Post Rh Interp	2960540156			Neg	

Current Page No: 1 Total Page No: 1 Zoom Factor: Page Width
 Finalize Transfusion Reaction Workup Report Schedule Print Print Now Close

* Required Field

Finalize/Print TRW

The user finalizes and views and/or prints the Transfusion Reaction Workup Report.

Assumptions

- Patient Testing: Record a Transfusion Reaction Workup was completed.

Outcome

- The Transfusion Reaction Workup Report is available for printing and signing.
- Final report findings are available for display in CPRS.

Limitations and Restrictions

- The user cannot invalidate or change a finalized TRW.
- Transfusion reaction reports are not available to clinicians through CPRS: they must be printed and filed with patient charts.
- VBECS displays these default values for transfusion reactions converted from VistA to VBECS. They do not reflect actual interpretations or results:

Converted Transfusion Reaction	Default
Pre-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Post-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Symptoms	No symptoms identified.
Implicated units	No units implicated in reaction.
Serologic tests	No serologic tests performed.

Additional Information

- The medical director's signature is on the patient's printed chart.
- If a report needs to be replaced, the user may retrieve this report after it is complete. The final report is a printed document with the medical director's signature.

User Roles with Access to This Option

All users

Finalize Transfusion Reaction Workup Report

►►► The Lead Technologist reviews the TRW details, identifies the type of reaction, and adds comments, according to the medical director's instructions. The user prints a copy for the medical director's signature and the Lead Technologist finalizes the TRW, generating a preliminary report. The final printed report requires the signature of the medical director.

User Action	VBECS
1. Select Reports from the main menu.	<ul style="list-style-type: none">Lists report names.Lists the completed TRWs for selection.
Select Finalize/Print TRW .	NOTES _____
Click the Include Non-Finalized TRW Reports check box (Figure 127).	When the user enters and saves the portion recorded in Patient Testing: Record a Transfusion Reaction Workup, the TRW status is "completed." When the Transfusion Reaction Workup Report is complete, the TRW status is "finalized."
Select a completed TRW to finalize.	

User Action	VBECS
Click OK .	
2. Click the Finalize Transfusion Reaction Workup Report check box.	<ul style="list-style-type: none"> Prepares the completed report for review. Displays the name of the blood bank medical director to sign the final report. VBECS recreates the final report for viewing or printing in the future.
3. Select the Transfusion Reaction Type and complete the Transfusion Reaction Comments required fields (Figure 128).	<ul style="list-style-type: none"> Displays the selected Transfusion Reaction Workup Report data. Requires the user to select a transfusion reaction type listed in Table 12: Transfusion Reaction Types. Allows the user to enter a TRW comment (maximum 1,000 characters). Allows the user to save the data entered. <p>NOTES</p> <p>The user must indicate one transfusion reaction type. Entry of a comment is optional for all types except O (other), which requires details to be entered.</p>
4. Click OK and Yes to save.	<ul style="list-style-type: none"> Updates the transfusion reaction to "finalized" in the database and prints a report for the medical director's signature. <p>NOTES</p> <p>VBECS notifies the user when DSS updates cannot be sent because VistALink is not active.</p> <p>VBECS will attempt to resend the DSS information when VistALink is available after any user logs into the VBECS application's VistALink connection.</p>
5. Select Schedule Print or Print Now and follow the instructions to print.	<ul style="list-style-type: none"> Prints the report on the selected printer.
6. Click Close to exit.	

Figure 127: Select Transfusion Reaction Workup


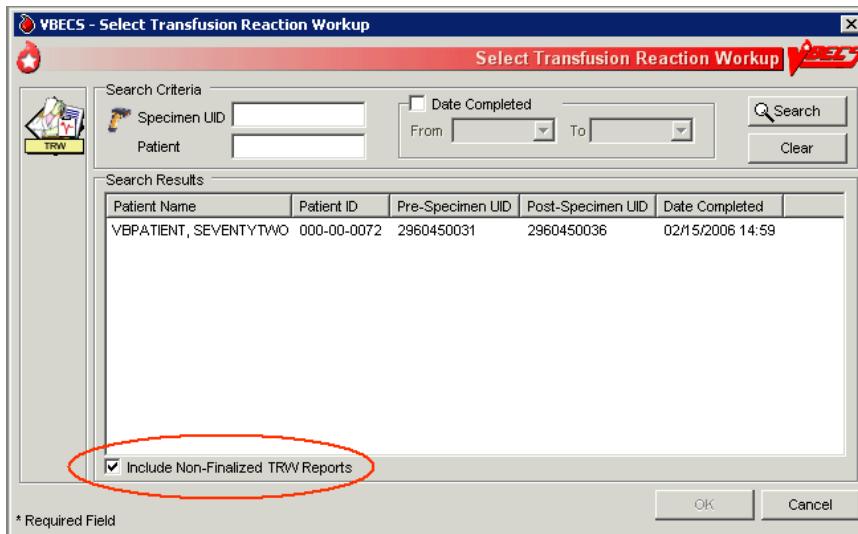


Figure 128: Finalize TRW Report

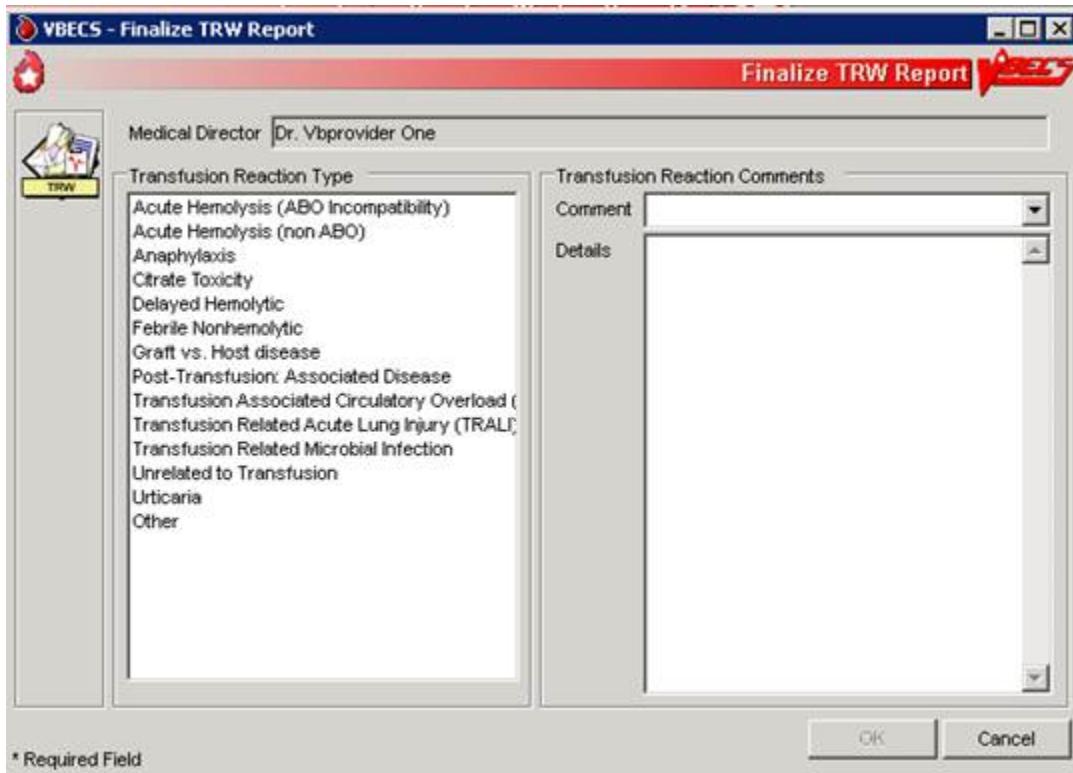

View/Print a Finalized Transfusion Reaction Workup Report

The user views and/or prints the finalized Transfusion Reaction Workup Report.

User Action	VBECS
<p>1. Select Reports from the main menu.</p> <p>Select Finalize/Print TRW.</p> <p>Select a completed TRW (Figure 126).</p>	<ul style="list-style-type: none"> Lists report names. Displays the selected Transfusion Reaction Workup Report data. Displays an option to view and/or print a finalized TRW by selecting a patient. <p>NOTES —————</p> <p>The user finalizing the report may not edit previously entered TRW data.</p>
2. Select a patient.	<ul style="list-style-type: none"> Compiles the report and displays an option to print the report or exit. <p>NOTES —————</p> <p>Details included in the finalized TRW report:</p> <ul style="list-style-type: none"> Patient name Patient ID Specimen UID Reported symptoms Specimen clerical checks (pre-, post-, and unit) Serologic testing (if performed, interpretations only) Transfusion reaction type Comment Implicated units (include from Patient Testing: Record a

User Action	VBECS
	<p>Transfusion Reaction Workup)</p> <ul style="list-style-type: none"> • TRW report completed by: <ul style="list-style-type: none"> ◦ Date and time the report was completed ◦ Signature line for the medical director ◦ Line for the date and time the medical director completed the review ◦ Finalized by ◦ Finalized date and time
3. Select Schedule Print or Print Now and follow the instructions to print.	<ul style="list-style-type: none"> • Prints the report on the selected printer.
4. Click Close to exit.	

Table 12: Transfusion Reaction Types

Transfusion Reaction Type Text
Anaphylaxis
Citrate Toxicity
Delayed Hemolytic
Febrile Nonhemolytic
Graft vs. Host Disease
Acute Hemolysis (non ABO)
Acute Hemolysis (ABO Incompatibility)
Transfusion Associated Circulatory Overload (TACO)
Transfusion Related Acute Lung Injury (TRALI)
Transfusion Related Microbial Infection
Post-Transfusion: Associated Disease
Unrelated to Transfusion
Urticaria
Other: free-text data required on new entries

Maintain Patient Records

Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement

The user enters or inactivates the Transfusion Requirements (TRs) in the patient's record by division.

Assumptions

- The patient exists in the VBECS database.

Outcome

- VBECS displays TRs to all divisions within a multidivisional database. The division in which they were entered maintains and controls them.
- VBECS maintains a historic record of changes to the patient's TRs.

Limitations and Restrictions

- Inactivation is defined and restricted by security level.
- VistA Special Instructions (SIs) from database conversion do not become TRs: they remain SIs. SIs are information only and are not enforced unless entered in the Component Requirements tab in Special Instructions & Transfusion Requirements.

Additional Information

- No current patient specimen or order is required to access this option.
- When TRs are not met, VBECS prohibits the printing of the Caution Tag and Blood Transfusion Record Form (BTRF) and prevents issue of the blood unit.
- User-entered antigen negative requirements are persistent. VBECS enforces them even when they are in conflict with the division's antibody table settings.
- Patient antibody identifications from VistA are included in the database conversion.
- By default, inactivation of an antibody removes requirements for antigen negative blood derived from the antibody entry.

User Roles with Access to This Option

All users

Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement

The user assigns patient blood component requirements or indicates the need for antigen negative blood to be applied during blood component selection and issue.

VBECS displays TR data and applies requirements during blood unit selection.

All items in these TR categories indicate the need for pre-transfusion processing of a unit before the unit may be issued to the patient. All users may activate entries. ►►► Only a Traditional Supervisor may inactivate a component requirement for a patient. A patient may have one or more TRs:

- Irradiate cellular products*
- Leukoreduce cellular products
- Washed RBC products

- Washed PLT products
- Sickle Cell Negative RBC products
- CMV negative cellular products*
- Antibodies Identified
- Antigen Negative RBC: A user or VBECS may create an antigen negative requirement.

*All product types in OTHER are evaluated for the requirement including those that are non-cellular.

User Action	VBECS
1. Select Patients from the main menu. Select Special Instructions & Transfusion Requirements , or Click  in the Patient Information Toolbar, when it is active.	<ul style="list-style-type: none"> • Displays options for processing patient-related functions. • Displays the option to select a patient, if one is not already selected. <p>NOTES</p> <p>The user may also click New SI or TR in Patient Testing to access this option.</p>
2. Select the Antigen Negative RBCs, Component Requirements, or Antibodies Identified tab (Figure 129).	<ul style="list-style-type: none"> • Displays the patient's TRs. • Allows the user to add or inactivate information based on his security level. <p>NOTES</p> <p>When a user accesses antigen negative requirements, VBECS displays and enforces the antigen negative requirement for the patient, regardless of the division setting for the antibody.</p> <p>When a user accesses the Antibodies Identified tab, VBECS compares the patient's antibodies with the division's antibody table setting. When the antibody is set to require antigen negative blood at issue, VBECS displays and enforces the antigen negative requirement.</p>
3. Select a TR from the pick list, or enter an antigen or antibody, as the tab allows. Click Add to move the selection to the patient's record. Click OK when all TRs are selected to save them to the patient's record. Click Yes to confirm changes and exit.	<ul style="list-style-type: none"> • Allows multiple selections from the Component Requirements, Antibodies Identified, and Antigen Negative tabs. <p>NOTES</p> <p>VBECS applies the TR to component classes, as indicated.</p> <p>The user may not add duplicate SIs or TRs in the same session, but may duplicate them in another session.</p> <p>The division code and date of entry are associated with the entered component or antigen negative requirement, or antibody specificity so that other divisions will know whom to contact for further information, when needed. The original entries may have been made through Patient Testing (associated with an antibody identification) or TRs from any division within a multidivisional database.</p>
4. Select the TR to be inactivated (Figure 130). Select and/or enter an explanatory comment and click OK .	<ul style="list-style-type: none"> • Allows the user to enter a comment. <p>NOTES</p> <p>When a Transfusion Requirement is inactivated in any division in a multidivisional configuration, a comment (free text or from a canned comment with the context: Transfusion Requirement Inactivation) must be added that will be saved and printed on the audit trail of originating division alerting the originating division of</p>

User Action	VBECS
	<p>the change to the patient's TR setting if they were not consulted regarding that change.</p> <p>Note: It is best practice that inactivating a TR in the non-originating division involves a consultation with the originating division and allowing the originating to make the change with any locally required documentation.</p> <p>For inactivation of an item:</p> <ul style="list-style-type: none"> • >>> A Traditional Supervisor from any division may inactivate TRs and persistent antigen negative requirements in a multidivisional database regardless of the source of the TR. • When an antibody specificity is inactivated, VBECS inactivates the corresponding calculated (VBECS-generated) antigen negative requirement. • For antibodies identified, VBECS displays a summary of antibodies currently associated with the patient and identifies their source. Additional antibodies may be added and previously identified antibodies may be inactivated.
5. Return to Step 2 to continue editing the patient's TRs.	
6. Accept entries and exit, or Return to Step 1 to make entries for a different patient.	

Figure 129: Enter Patient Special Instructions and Transfusion Requirements


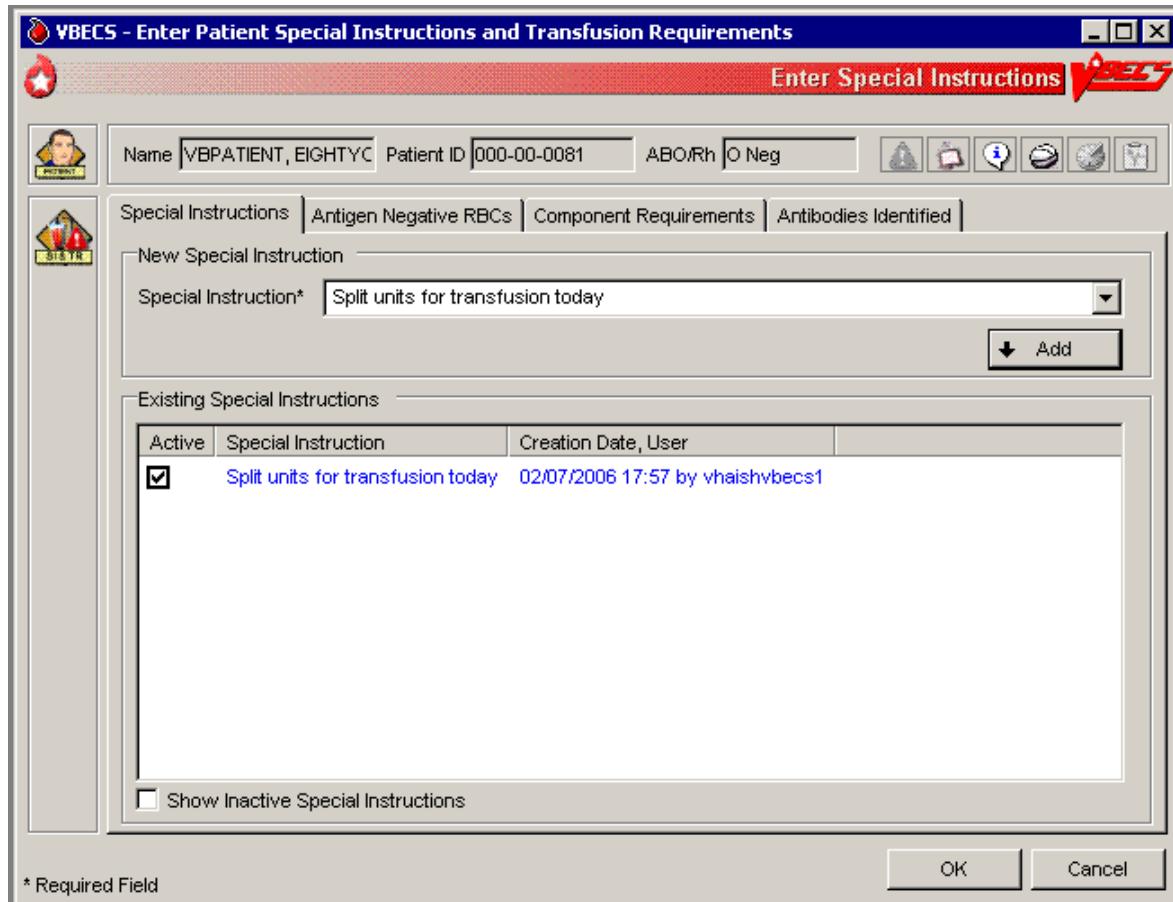


Figure 130: Inactivate Component Requirements

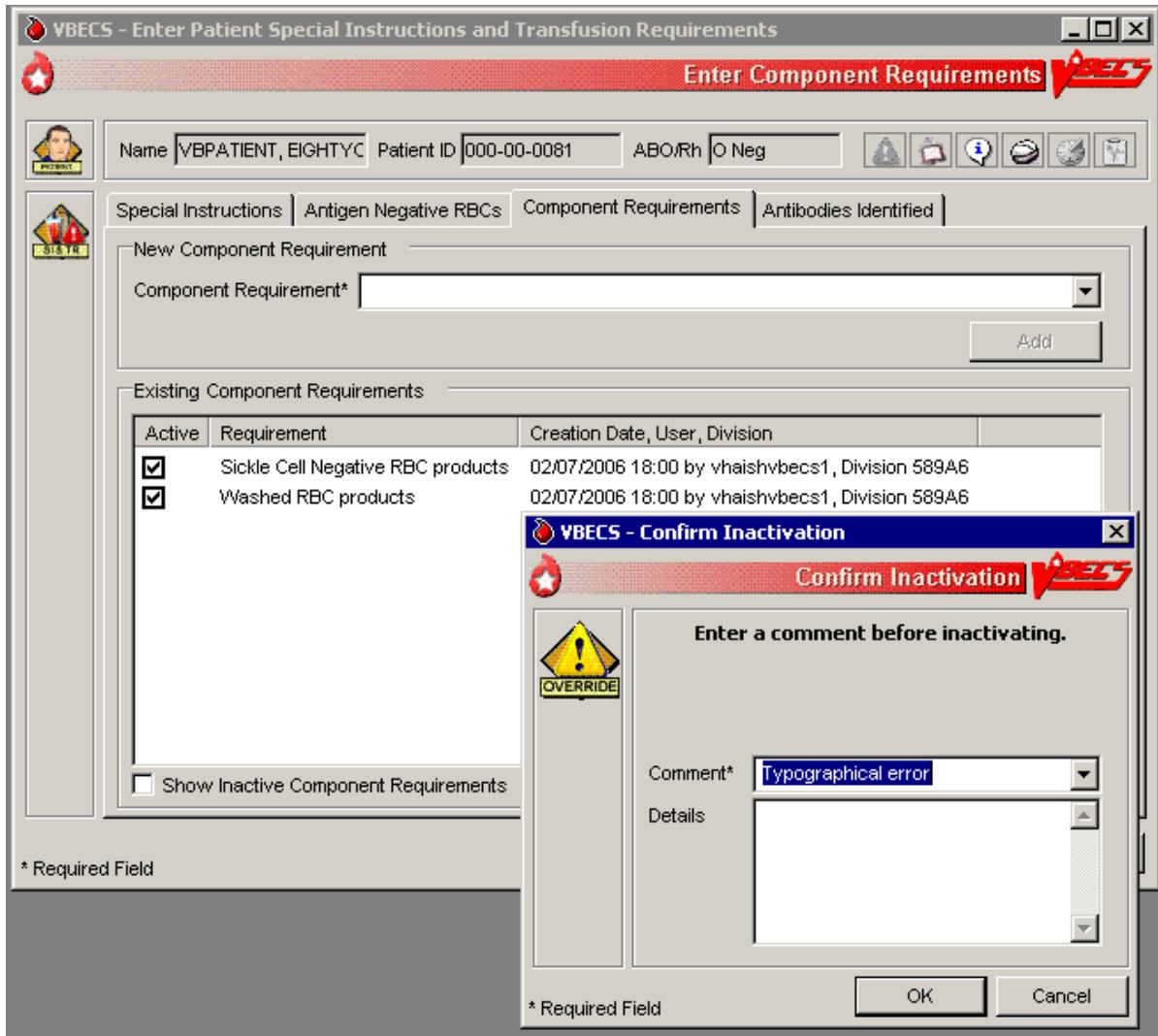

Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions

The user enters comments in a patient's file that are not enforced by VBECS.

Assumptions

- Fields containing Special Instructions (SIs) may be populated during database conversion.

Outcome

- When the user searches for a patient or when the Patient Information Toolbar is displayed, VBECS displays  to alert the user to view the patient's SI entries.
- When a user retrieves a patient record and is in a multidivisional setting, he may view SIs created in all divisions within the database.

Limitations and Restrictions

- VBECS does not enforce Transfusion Requirements (TRs) or antigen negative requirements entered through this option.
- SIs moved to VBECS during the database conversion may contain information that must be entered as TRs in VBECS. Inactivate the old SI and create a TR and new updated SI: the database conversion software cannot automatically translate the SI from database conversion into a rule-based TR.
- VBECS will truncate SI text from the VistA database conversion.

Additional Information

- Inactivation of patient SI entries is allowed only in the division of entry.
- Only currently active SIs at a division will automatically display via the database to all divisions.
- Inactive SIs are not deleted from the VBECS database. They may not be reactivated but may be viewed in the maintenance fields when requested.
- Canned comments in the Special Instructions context must be entered by the division before they can be selected for use.
- Previously saved Special Instructions that require editing due to context or spelling errors must be inactivated and reentered in total.
- Database conversion entries are labeled as such.
- SIs from database conversion may include component or antigen negative requirements that must be entered in TRs so that VBECS can enforce them.

User Roles with Access to This Option

All users


Special Instructions & Transfusion Requirements: Enter a Special Instruction

A user may add to and deactivate SI comments in a patient's blood bank record. Users in all divisions in a multidivisional database may view these comments.

User Action	VBECS
1. Select Patients from the main menu. Select Special Instructions & Transfusion Requirements , or	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays the option to select a patient, if one is not already selected. NOTES _____

<p>Click  in the Patient Information Toolbar, when displayed.</p>	<p>The user may also click New SI or TR in Patient Testing to access this option.</p>
<p>2. Select a patient. Click OK.</p>	<ul style="list-style-type: none"> • Displays options to select a patient. • Displays the option to enter a new comment or remove an existing comment in the division. • Displays the active SI entries in reverse chronological order with all active SIs first, along with date, division code, and text of the entry. <p>NOTES —————</p> <p>The need for SIs is usually temporary and may not be required for all transfusions.</p> <p>The user may choose to display inactivated SI comments.</p>
<p>3. Select the Special Instructions tab (Figure 131). Click the Show inactive special instructions check box to display inactivated SI entries.</p>	<ul style="list-style-type: none"> • Allows the user to enter or deactivate comments. • Displays inactive Special Instructions but does not allow the user to reactivate them. <p>NOTES —————</p> <p>A user may add or deactivate a patient-specific comment at any time in the division in which it was entered.</p> <p>►►►► A Traditional Supervisor may deactivate SIs from database conversion.</p>
<p>4. Confirm the activation of the entry and enter a comment. Click OK. Return to Step 2 to add more SIs for this patient. Select another patient and return to Step 1, if desired. Exit.</p>	<ul style="list-style-type: none"> • Saves the entry in the database.

Figure 131: Enter Patient Special Instructions and Transfusion Requirements

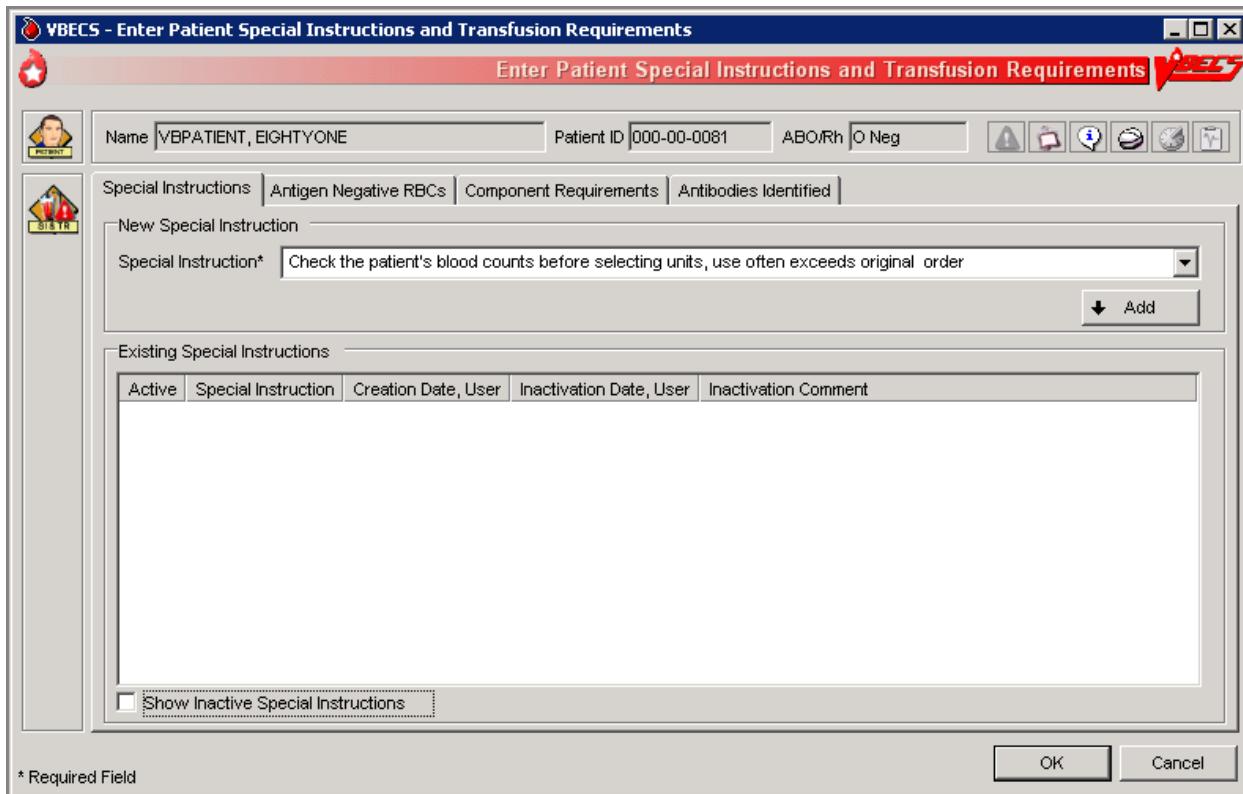

Post-Transfusion Information

The user enters post-transfusion details from the Blood Transfusion Record Form (BTRF) or VBECS saves post-transfusion details received from the Transfusion Verification System.

Assumptions

- The BTRF is returned and completed for each unit.
- The unit status is “issued” or “presumed transfused” in VBECS.
- Interface connections to the BCE COTS system are enabled to receive transfusion update messages.

Outcome

- The unit status is updated to “transfused.”
- The VBECS patient record includes the blood transfusion information.

Limitations and Restrictions

- VBECS cannot update a unit record to presumed transfused when a user has that unit record open (locked) in Enter Post Transfusion Information. VBECS will attempt to update the unit the next time the presumed transfused background routing runs if the unit remains in an issued status.
- If post-transfusion information is not transmitted to VBECS by the BCE COTS interface, the user will need to enter that information to update the record from presumed transfused to transfused.
- The transfusion comments and volume transfused entered in the Transfusion Verification application as an update to a completed transfusion are not sent to VBECS.
- Transfusion reaction symptoms are not provided to VBECS by the Transfusion Verification application unless a transfusion reaction has also been indicated for that event.

Additional Information

- A user with the appropriate security level may enter or edit bedside verification information unless it is received from BCE COTS interface.
- The user may indicate that a transfusion reaction and symptoms were noted on the BTRF as free text. Noting a transfusion reaction does not create an order for a transfusion reaction workup (TRW).
- Remove a “transfused” status through Remove Final Status unless that unit was transfused by the Transfusion Verification application. See limitations in Remove Final Status.

User Roles with Access to This Option

All users

Post-Transfusion Information

The user verifies the date and time transfused with a returned BTRF, adds other transfusion details, initiates a system check of the bedside verification information, enters that information, and indicates whether a reaction occurred during or immediately after the transfusion, as documented on the BTRF. VBECS does not require entry of post-transfusion information: once a day, without user input, VBECS moves units 48 hours after issue to “presumed transfused” status.

User Action	VBECS
1. Select Patients from the main menu. Select Post-Transfusion	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays an option to scan the Caution Tag or enter a patient name or ID.

User Action	VBECS
Information. 2. Select a patient and click OK .	<ul style="list-style-type: none"> Lists the units issued to the selected patient. Displays the issue date and time. Displays the units one at a time for update in the order in which the user selected them. <p>NOTES</p> <p>►►► VBECS displays an option to obtain a list of transfused (those not received by BCE COTS) units for the patient when the user is a Traditional Supervisor or above.</p> <p>When post-transfusion information is not entered into VBECS within 48 hours from the date/time the unit was recorded as "issued", VBECS will change the unit status to "transfused" with a presumed indicator after the execution of a scheduled background job. The unit will then be considered transfused by the system regarding report generation. If post-transfusion information is entered later, the presumed indicator is deactivated.</p>
3. Click one or more check boxes to select units (Figure 132). Click OK .	<ul style="list-style-type: none"> Allows the user to select units for which to enter post-transfusion information. <p>NOTES</p> <p> VBECS does not list units (for selection) with a final status. When a user enters a unit with a final status, VBECS alerts the user and clears the screen. The user may select another unit.</p>
4. Enter transfusion details for a blood unit. Verify the unit selection.	<ul style="list-style-type: none"> Displays the unit information entered by the bedside verification check and allows the user to enter or edit the information. Displays an option to enter or edit unit data. <p>NOTES</p> <p>►►► A user with a security setting of Traditional Supervisor may edit user entered bedside verification, transfusion reaction, or post-transfusion detail information.</p> <p>VBECS displays an option to enter information for each unit. This information is available for patient record retrieval.</p>
5. Select one or more units. Enter or edit the bedside verification check information with the proper security level.	<ul style="list-style-type: none"> Displays the entries for review, including the comment, until all details are entered. Displays an option to save the data and exit. <p>NOTES</p> <p>Upon data entry and saving of the post-transfusion information, VBECS assigns a status of "transfused" for the unit of blood assigned to the indicated patient.</p>
6. Enter post-transfusion details, including: <ul style="list-style-type: none"> Transfusion start and end dates and times Transfusionist ID Location of transfusion Whether the transfusion was 	<ul style="list-style-type: none"> Displays an option to enter unit data. Warns the user or allows the user to continue to the next selected unit. <p>NOTES</p> <p>When a user indicates that a transfusion reaction occurred, VBECS does not automatically order a TRW.</p>

User Action	VBECS
<ul style="list-style-type: none"> completed or interrupted If interrupted, the amount transfused Whether there was a transfusion reaction If there was a reaction, symptoms (Figure 133) <p>Click OK to save.</p>	<p>VBECS notifies the user when the transfusion end date and time is earlier than the start or issue time. The user may correct the entry.</p> <p>When the unit was not issued to a remote storage location, VBECS compares the date and time the unit was issued to the start transfusion time entered.</p> <p>When the difference is greater than 30 minutes, VBECS notifies the user that this is a delayed start time and allows the user to edit the transfusion start date and time. No allows the user to continue. Yes allows the user to edit the date and time the transfusion ended. When the user edits the date and time transfused, VBECS repeats this process until the user stops editing and chooses to continue.</p> <p>When the unit was not issued to a remote storage location, VBECS compares the date and time recorded as the transfusion end time to the start date, calculates the transfusion time in minutes, and compares this period to the maximum transfusion time for the component class. When the calculated transfusion time for this unit exceeds the defined maximum transfusion time for the component class, VBECS notifies the user that this is a prolonged transfusion and allows the user to edit the transfusion end date and time. No allows the user to continue. Yes allows the user to edit the date and time the transfusion ended. When the user edits the date and time transfused, VBECS repeats this process until the user stops editing and chooses to continue.</p>
7. Respond to warnings or continue to enter data for the next unit.	<ul style="list-style-type: none"> Allows the user to repeat for subsequent selected units until all details are entered. Displays an option to save the data and exit.
	<p>NOTES —————</p> <p>VBECS saves all transfusion information for future retrieval and reporting as part of the patient and unit historic records.</p> <p>VBECS notifies the user when DSS updates cannot be sent because VistALink is not active.</p> <p>VBECS will attempt to resend the DSS information when VistALink is available after any user logs into the VBECS application's VistALink connection</p>
8. Exit.	

Figure 132: Post-Transfusion Information

The screenshot shows the VBECS - Post-Transfusion Information window. On the left, a list of transfused units is displayed:

Unit ID	Product	ABO/Rh	Expires	Date Issued	Transfusion End Date
12030802	RBC CPD>AS1, 04210	A Pos	12/17/2008 11:59:00 PM	12/3/2008 3:46:01 PM	12/3/2008 3:46:41 PM
12030801	PLT None, 12000	A Pos	12/4/2008 11:59:00 PM	12/3/2008 3:46:01 PM	
W039911111111	RBC CPD>AS1, E0291D00	A Pos	11/21/2008 2:23:00 PM	10/27/2008 3:16:38 PM	10/27/2008 3:28:35 PM
W08888888888881	RBC CPD>AS1, E0291V00	A Pos	9/30/2008 11:59:00 PM	9/15/2008 2:41:47 PM	9/15/2008 2:42:11 PM
W036477777777	WB ACD-A, E0001V00	O Pos	9/11/2008 11:59:00 PM	9/4/2008 1:34:54 PM	9/4/2008 1:35:33 PM
0822002	PLT None, 12000	A Pos	8/23/2008 11:59:00 PM	8/22/2008 5:50:50 PM	11/13/2008 12:26:43 PM

A checkbox at the bottom left allows selecting "Show Transfused Units (those not received by BCE COTS)". On the right, a detailed view of a selected unit (Unit Number 12030801, Expiry 12/04/08 11:59PM) is shown in a modal window. The details include:

- PLT None**
- A Rh POSITIVE**
- Davis Blood Center**
- VOLUNTEER DONOR**
- R. Only**
- PROPERLY IDENTIFY INTENDED RECIPIENT**
- See circular of information for indications, contraindications, cautions and methods of infusion.**
- Disease Marker Testing**: Unit fully tested and negative for all disease markers.

Figure 133: Post-Transfusion Information Details

The screenshot shows the VBECS - Post-Transfusion Information Details window. The top header includes fields for Name (VBPATIENT, EIGHTYONE), Patient ID (000-00-0081), and ABO/Rh (O Neg). The main form contains the following sections:

- Issue Details**: Date Issued (02/09/2006 10:45), Issued By (vhaishvbecs1).
- Post-Transfusion Information**: Tranfusion End Date (02/13/2006 16:22), Patient Location at Transfusion (A2MED IM-LE).
- Tranfusion Completed?**: Radio button selected for "Tranfusion Completed".
- Amount Transfused?**: Radio buttons for 1/4, 1/2, 3/4, and Volume (0 mL).
- Tranfusion Reaction Noted?**: Radio button selected for "Reaction Noted". Symptom listed: goofy smile.
- Comments**: A text input field.

* Required Field

Automated Entry of Post-Transfusion Details from BCE COTS Application

VBECS receives and saves the details of a unit transfusion episode as recorded in the BCE COTS system.

User Action	VBECS
1. BCE COTS system user completes a transfusion for a patient which triggers a HL7 message containing transfusion information to VBECS for a blood unit.	<ul style="list-style-type: none">VBECS rejects a BCE COTS HL7 message when data elements are not acceptable per the data entry rules for manual entry. <p>NOTES</p> <p>Units that are not updated by a BCE COTS message will be marked presumed transfused 48 hours after issue by a VBECS background process.</p> <p>Transfusion information is available in the Unit History Report.</p>
2. VBECS updates the patient transfusion record with the details received in the successfully acknowledged HL7 message.	<ul style="list-style-type: none">VBECS generates an exception report entry: Electronic Transfusion Record Update when a BCE COTS update to a previously saved manually or electronically entered transfusion record is received.When an update to a previously saved transfusion record is made by BCE COTS or a user, VBECS will display the original entry and the updated details in the Patient History and Unit History Reports. <p>NOTES</p> <p>The BCE COTS system delivers these patient and unit specific transfusion data elements to VBECS for update</p> <ul style="list-style-type: none">• Bedside verification date and time• Transfusion verification manual Unit ID entry reason• Transfusion Start date and time• Transfusionist ID #1• Transfusionist ID #2 (may be blank if so configured in BCE COTS application)• Transfusion End Date and time• Patient location at Transfusion• Transfusion interrupted indicator• Transfused volume, (ml not in fractions by the BCE COTS)• Transfusion Reaction noted• Transfusion Reactions symptoms entered, if any• Transfusion comments, if any

Document ABO Incompatible Transfusions

The user enters post-transfusion details from the Blood Transfusion Record Form (BTRF) to document ABO incompatible and other inadvertent transfusions.

Assumptions

- A unit was transfused to a patient other than the intended recipient.
- The connection to VistA is active.

Outcome

- The unit status is updated to “transfused.”
- The patient record includes the blood transfusion information.

Limitations and Restrictions

- VBECS does not enforce compatibility checks between the patient and the unit.
- VBECS does not update CPRS component order status when the Document ABO Incompatible Transfusion option is used.

Additional Information

- FDA guidances govern when to file blood product deviation reports.
- Follow local policy and procedures for additional reporting requirements.
- VBECS requires a specimen and an order placed through CPRS to process transfusions. When a patient was transfused and no specimen was collected or order placed for blood components, the user may record the transfusion using this option.

User Roles with Access to This Option

►►► Traditional Supervisor

Document ABO Incompatible Transfusions

When a unit is transfused to the incorrect patient, the user enters transfusion details for that unit.

 *VBECS does not enforce compatibility checks to allow consistent and accurate online availability of a unit's final disposition. Use of this option is restricted to supervisory blood bank personnel who have thoroughly investigated the transfusion and documented the inability to use normal VBECS business processes to enter the transfusion documentation.*

User Action	VBECS
1. Select Supervisor from the main menu. Select Document ABO Incompatible Transfusion .	<ul style="list-style-type: none">• Displays options for processing patient-related functions.• Displays an option to scan the Caution Tag or enter a patient name or ID.
2. Select a patient and click OK .	<p>Displays an option to select a patient that may not be available from the VBECS files.</p> <p>NOTES _____</p> <p>Displays a message when the user selects a patient with a blank</p>

User Action	VBECS
3. Select a unit to document a transfusion that does not follow system rules.	<p>first or last name, or the name is too long. The user is instructed to contact the VBECS Administrator for assistance.</p> <ul style="list-style-type: none"> Displays an option to enter a unit ID and product code or search for one. <p>NOTES _____</p> <p>VBECS presents all units in the system for selection when the search button is used. This search will require time to complete as the number of units to display increases. VBECS does not enforce compatibility checks in this option.</p> <p> When a user enters a unit with a final status, VBECS alerts the user and clears the screen. The user may select another unit.</p>
4. Enter post-transfusion details and a comment.	<ul style="list-style-type: none"> Displays an option to enter or edit unit data. See Post-Transfusion Information.  VBECS emits an audible alert, requires a comment, and captures the transaction details for inclusion in an Exception Report (exception type: ABO incompatible unit transfused). Displays the entries for review, including the comment, until all details are entered.
5. Click OK and Yes to save.	<ul style="list-style-type: none"> Saves the data. Allows the user to enter comments and details.
6. Enter or select a comment and details, in accordance with local policy and procedures. Click OK .	<p>NOTES _____</p> <p>VBECS saves transfusion information for future retrieval and reporting as part of the patient and unit historic records.</p> <p>Upon data entry and saving of the post-transfusion information, VBECS assigns a status of "transfused" to the unit of blood for the indicated patient.</p> <p>VBECS captures details of this override in an Exception Report.</p>
7. Exit.	

Justify ABO/Rh Change

The user changes a patient's historical blood type only when there is a clinical change in blood type, not to correct an error. (See Invalidate Test Results.)

Assumptions

- The patient has a previous ABO/Rh on record from the user's current division, from a different division within the database, or from the database conversion.
- The patient has a current specimen test result with an ABO/Rh interpretation that does not match the historical record. Both previous and current specimen test results are valid.

Outcome

- VBECS permanently marks the patient record and the specimen test to indicate that the patient has a justified ABO/Rh change, regardless of the user's division.
- The justified ABO/Rh on the patient record forces the user's selection to be used as the new historical record, replacing the current test result as the historic record. VBECS does not check for ABO/Rh discrepancies with tests entered for this patient prior to the justified test result.
- VBECS permanently disallows electronic crossmatch (eXM) for this patient.

Limitations and Restrictions

- VBECS does not allow the selection or issue of certain blood components that may appear ABO/Rh compatible when discrepant ABO/Rh results are not justified.

Additional Information

- When the discrepant ABO/Rh tests are on the same specimen, the user may not justify the patient's ABO/Rh.
- Do not use this option for ABO/Rh discrepancies that occur because the patient had a recent transfusion of a component of a type other than his own.
- ABO/Rh discrepancies that occur because the previously tested or current specimen was tested in error are handled by inactivating the previously recorded test results through Invalidate Test Results.
- Subsequent specimens may need to be justified in the presence of continued discrepant ABO results.
- This option is used in rare clinical situations such as following bone marrow transplant.

User Roles with Access to This Option

►►► Traditional Supervisor


Justify ABO/Rh Change

The user determines that the patient's current ABO/Rh testing results are valid, although they differ from previously recorded results, and that the patient requires a justification indicator on his record.

This option is used only when both the previous specimen testing and the current specimen testing are considered valid. There is a valid explanation for why the ABO/Rh interpretations do not match and invalidation of the test result is not appropriate. Once the current patient ABO/Rh result change is justified, VBECS marks the patient's record so that special handling for component orders for the patient apply from this point forward and remain in effect until the patient record is inactivated by the recording of his death in VBECS or by removing the justification.

User Action	VBECS
1. Select Supervisor from the main menu. Select Justify ABO/Rh Change .	<ul style="list-style-type: none"> • Displays options available to supervisors.
2. Enter a patient name in the Patient field and click Search to select a patient. Select a patient. Click OK .	<ul style="list-style-type: none"> • Displays a screen for patient selection. • Disallows the justification or allows the user to continue. <p>NOTES</p> <p>VBECS compares the results of the current ABO/Rh with the previous ABO/Rh results. When there is no discrepancy, VBECS warns that no justification is required. There is no override. VBECS clears the screen and the user may select another patient.</p> <p>VBECS warns the user to use this option only when both the current and previous specimen tests are valid, otherwise, use Invalidate Patient Test Results.</p>
3. Select an ABO/Rh interpretation as the patient's justified historical record (Figure 134).	<ul style="list-style-type: none"> • Displays details of the historic and current ABO/Rh testing for comparison and a list of available ABO and Rh types for user selection. • Requires a comment to continue. <p>NOTES</p> <p>The user may select one ABO group and one Rh type for the patient to serve as the justified ABO/Rh typing in the historic record. VBECS displays the current ABO and Rh interpretations as the defaults.</p> <p>VBECS requires a comment to continue.</p>
4. Confirm the update of the database with the justified ABO/Rh. Exit.	<ul style="list-style-type: none"> • Updates the patient record.

Figure 134: Justify ABO/Rh Change

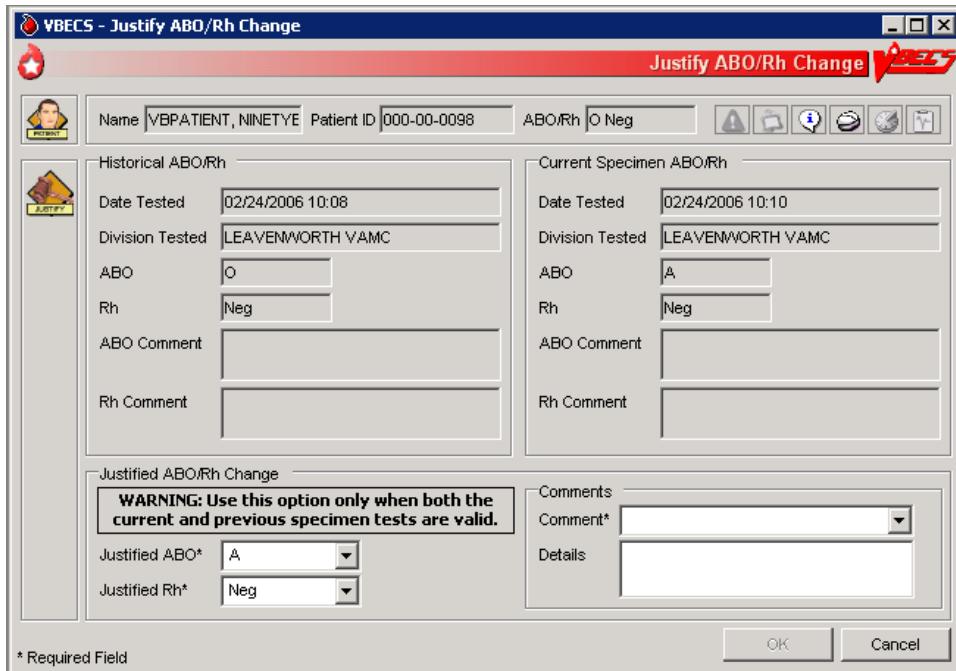

Reports

The report information that is presented on-screen, not the print preview, can be sorted by clicking on each column header. The sort changes made do not apply to the printed reports. Reports print in a default sort order.

Note: Switching to Normal layout view from Print layout view so the user can search the report may cause the error message “**Error. Thread was being aborted.**” to appear for a large report. The report proceeds to load in the different layout view chosen by the user. This is a normal behavior of a report engine and no harm is done.

Administrative Data Report

The user views and/or prints the Administrative Data Report.

Assumptions

- Data are available for inclusion in the report.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- None

Additional Information

- The report is available as a single-division or multidivisional (all divisions) report.

User Roles with Access to This Option

All users

Administrative Data Report

The user views and/or prints the Administrative Data Report, which includes blood product transactions for a given date range.

User Action	VBECS
1. Select Reports from the main menu. Select Administrative Data Report .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. Click the Current Divisions or All Divisions radio button.	<ul style="list-style-type: none">• Allows the user to select one or all divisions. <p>NOTES —————</p> <p>For each component class and its subdivisions, VBECS tallies:</p> <ul style="list-style-type: none">• Products received• Products modified/prepared• Units transfused• Units shipped• Outdated units

User Action	VBECS
	<ul style="list-style-type: none"> Discarded credit units Discarded waste units
3. Enter or select start and end dates.	
4. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none"> Compiles and displays the report.
5. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> Prints the report on the selected printer.

Audit Trail

The user views and/or prints an Audit Trail Report.

Assumptions

- Data are available for inclusion in the report.
- The user must be logged into the division in which the report is to be generated.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- VBECS does not accommodate online review (verification, signature) of reports.

Additional Information

- Audit trail entries reflect the actual date and time the change was made.

User Roles with Access to This Option

All users

Audit Trail

The user views and/or prints an Audit Trail Report, which lists changes to specific patient, unit, and division configuration settings. It is a regulatory requirement to maintain documentation proving that these changes were reviewed. Modify local record retention policy to instruct users to print, review, document, and save these reports. The Audit Trail Report contains only data changes generated in the user's login division.

User Action	VBECS
1. Select Reports from the main menu. Select Audit Trail .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. Enter or select start and end dates.	<p>NOTES —————</p> <p>The parameter table and Transfusion Requirements Report sections include changes from no previous entry to an entry and from an entry to no entry.</p> <p>Each entry includes an area for supervisor acknowledgement (signature or initials and date of review) and comment.</p>
3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none">• Compiles and displays the report.
4. Select Print to print the report and	<ul style="list-style-type: none">• Prints the report on the selected printer.

User Action	VBECS
exit, or click Close to exit without printing.	<ul style="list-style-type: none">• See Appendix B: Table 18: Details in Audit Trail Report.

Blood Availability

The user views and/or prints the Blood Availability Report.

Assumptions

- A blood unit is logged into a division in the database.

Outcome

- VBECS users may view available blood unit data in several formats.

Limitations and Restrictions

- None.

Additional Information

- VBECS displays the defined reports and allows the user to view an individual Patient Blood Availability Report or to create a custom report.
- The report is available as a single-division or multidivisional (all divisions) report.

User Roles with Access to This Option

All users

Blood Availability

The user may check the blood availability for a specific patient and view and/or print a variety of Blood Inventory Reports from the user's and other divisions, as available. Blood Availability Reports include the product's unit status, indicators, patient information for units in the "assigned" and "issued" statuses, restrictions, and special testing or handling of a unit.

User Action	VBECS
<p>1. For a non-individual report, select Reports from the main menu.</p> <p>Select Blood Availability, or</p> <p>For a patient-specific report, select Patients from the main menu.</p> <p>Select Blood Availability.</p>	<p>• Displays options for generating, viewing, and printing blood reports.</p> <p>• Displays non-individual report type and compilation criteria.</p> <p>• Displays options for generating a patient-specific report.</p> <p>NOTES</p> <p>The "issued to" location for an issued unit, special testing associated with a unit, disease marker testing, and comments are printed in the additional information area of the report.</p>
<p>2. Click a radio button to select a report format:</p> <ul style="list-style-type: none">• Daily Inventory (All Units)• Daily Inventory (Available Units)• Antigen Typed Inventory (Available Units)• Units with No Disposition (Available Expired)• Units with Final Disposition (Not Transfused)• Custom Report• Individual Patient Blood	<p>• Displays a print preview of the report.</p> <p>NOTES</p> <p>The user may click the Custom Report radio button and OK, and then click the appropriate check boxes in the Component Class, Testing, Donation Type, Unit Status, Date Range, ABO/Rh, Product Attributes, and Division tabs to customize the report.</p> <p>The user may view individual patient Blood Availability Reports, including only blood units currently assigned, crossmatched, or issued for a patient, online only.</p>

User Action	VBECS
Availability	<p>When a user selects Units with Final Disposition, the report includes units with a status of "discarded," "modified," or "transferred."</p> <p>When a user selects Units with No Disposition, the report includes expired units with a status of "limited" or "available."</p>
<p>3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or</p> <p>Click the Schedule Print Job radio button and OK, then enter or select a print date and printer name. Click OK to view and print the report.</p>	<ul style="list-style-type: none"> • Compiles and displays the report.
4. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Cost Accounting Report

The user views and/or prints the Cost Accounting Report.

Assumptions

- Data for the report were saved in VBECS.

Outcome

- The report is viewable or printed.

Limitations and Restrictions

- None

Additional Information

- The report is for a single division.

User Roles with Access to This Option

All users

Cost Accounting Report

The user views and/or prints the Cost Accounting Report for a division. The Cost Accounting Report itemizes incoming and outgoing shipments, modifications of blood products, and patient testing costs for a given period. The report includes base cost, special testing, modifications, and return credits for units, as applicable.

User Action	VBECS
1. Select Reports from the main menu. Select Cost Accounting Report .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. Enter or select start and end dates.	<ul style="list-style-type: none">• Compiles and displays the report. <p>NOTES —————</p> <p>Final calculations appear at the end of the report:</p> <ul style="list-style-type: none">• Subtotal of base cost per invoice (Incoming Shipment section only)• Subtotal of special testing cost (Incoming Shipment section only)• Subtotal of return credit for an invoice (Outgoing Shipment section only)• Total base cost of all units• Total return credits of all units• Total special testing cost of all units• Total off-site antibody identification testing cost• Total number of units received• Total number of units transferred• Total number of units modified• Total cost of modifications• Total number of units discarded
3. In the Report Compilation area,	<ul style="list-style-type: none">• Compiles and displays the report.

User Action	VBECS
<p>click the Create Report Now radio button. Click OK to view the report, or</p> <p>Click the Schedule Print Job radio button and OK, then enter or select a print date and printer name. Click OK to view and print the report.</p>	
4. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

C:T Ratio Report

The user views and/or prints the C:T Ratio Report.

Assumptions

- Data for the report are available.

Outcome

- The report is viewed or printed.

Limitations and Restrictions

- None

Additional Information

- This calculation does not include component classes that do not require XMs (FRESH FROZEN PLASMA, PLATELETS, CRYOPRECIPITATE) and component types in the OTHER component class (when no crossmatch is performed).
- The report is available as a single-division or multidivisional (all divisions) report.

User Roles with Access to This Option

All users

C:T Ratio Report

This option compiles and calculates the crossmatch-to-transfusion ratio. The user may select one of five report formats.

User Action	VBECS
1. Select Reports from the main menu. Select C:T Ratio Report .	<ul style="list-style-type: none">• Lists report names.• Displays report type and compilation criteria.• Displays a crossmatched date range for selection.
2. Enter a crossmatched date range.	<ul style="list-style-type: none">• Displays an option to compile the report for one or all divisions in the database.
3. Click the Current Division or All Divisions radio button.	<ul style="list-style-type: none">• Allows the user to select a report format. <p>NOTES —————</p> <p>The user may request the Detailed Report, Physician Summary Report, Treating Specialty Summary Report, Detailed Report with Physician Summary, or Detailed Report with Treating Specialty Summary.</p>
4. Click a radio button to select a report format: <ul style="list-style-type: none">• Detailed Report• Physician Summary Report• Treating Specialty Summary Report• Detailed Report with Physician Summary• Detailed Report with Treating Specialty Summary	<ul style="list-style-type: none">• Compiles the report and displays an option to print the report or exit. <p>NOTES —————</p> <p>The physician included in this report is the ordering physician. VBECS compiles the C:T ratio report as a summary or detailed report. All versions include a single C:T ratio for the division.</p> <p>VBECS determines the C:T ratio for a division by dividing the number of units crossmatched (electronically or serologic) during the period selected by the number of units transfused at the</p>

User Action	VBECS
	<p>division.</p> <p>This calculation does not include units with a XM result of "Compatible - Don't Transfuse" or "Incompatible - Unsafe to Transfuse."</p> <p>The detailed C:T ratio report includes a list of units crossmatched (electronically or serologic) for a patient during the selected date range. The report includes the unit ID, physician name, treating specialty, XM interpretation, date and time the XM was completed, product short name, and an indication as to whether the unit was transfused or presumed transfused.</p> <p>VBECS determines the C:T ratio for a physician by dividing the number of units that were crossmatched (electronically or serologic) during the period selected by the number of transfused units.</p> <p>VBECS determines the C:T ratio for a treating specialty by dividing the number of units crossmatched (electronically or serologic) for all patients listed in that treating specialty by the number of units transfused to all patients listed for that specialty.</p>
<p>5. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or</p> <p>Click the Schedule Print Job radio button and OK, then enter or select a print date and printer name. Click OK to view and print the report.</p>	<ul style="list-style-type: none"> • Compiles and displays the report.
6. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Division Transfusion Report

The user views and/or prints the Division Transfusion Report.

Assumptions

- Transfused unit data for the report were saved in VBECS.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- None

Additional Information

- VBECS compiles data for the report from the user's login division.
- Appendix B: Table 22: National Treating Specialty Table defines the treating specialty.

User Roles with Access to This Option

All users

Division Transfusion Report

The user views and/or prints a Division Transfusion Report, which lists details of units transfused during a selected period.

User Action	VBECS
1. Select Reports from the main menu. Select Division Transfusion Report .	<ul style="list-style-type: none">• Lists report names.• Displays fields for entering a report date range.
2. Enter or select start and end dates. In the Report Type area, select the Issuing Physician Report , Treating Specialty Report , or Issued-To Location Report radio button.	<ul style="list-style-type: none">• Allows the user to compile the report by issuing physician, by treating specialty, or by issued-to location.• Allows the user to include some or all physicians, treating specialties, or issued-to locations. <p>NOTES ——————</p>
3. In the Report Compilation area, click the Create Report Now radio button. Click OK to review the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none">• Compiles and displays the report. <p>NOTES ——————</p> <p>The report displays tallies of units transfused by component class within the date range, by:</p> <ul style="list-style-type: none">• Physician name• Treating specialty• Issued-to location• Number of units transfused• Number of units presumed transfused• Total number of patients transfused• Total number of units transfused and presumed transfused

User Action	VBECS
	<p>Reports include:</p> <ul style="list-style-type: none"> • Patient name • Patient ID • Unit ID • Product short name • If pooled, number of participants in pool • Date and time Issued • Issued-to location • Physician • Treating specialty • Transfused or presumed transfused indicator, if present • Transfusion comments, if any
4. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Division Workload Report

The user views and/or prints a workload report.

Assumptions

- VistA workload codes were assigned to VBECS processes through Workload Codes.
- Data for the report were saved in VBECS.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- None

Additional Information

- The report is available as a single-division or multidivisional (all divisions) report.

User Roles with Access to This Option

All users

Division Workload Report

The user views and/or prints a workload report, which includes totals of workload types, for one or all divisions.

User Action	VBECS
1. Select Reports from the main menu. Select Division Workload Report .	<ul style="list-style-type: none">• Lists report names.• Displays report division and compilation criteria.• Displays a date range for selection.
2. Click the Current Divisions or All Divisions radio button.	<ul style="list-style-type: none">• Allows the user to select one or all divisions.• Displays a date range for selection.
3. Enter or select start and end dates.	<ul style="list-style-type: none">• Compiles and displays the report.• Allows the user to print the report.
	<p>NOTES —————</p> <p>A multiple-division workload report includes totals for each type of workload for each division requested in the report.</p> <p>VBECS searches its workload file for entries within the selected date range.</p>
4. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none">• Compiles and displays the report.

User Action	VBECS
5. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none">Prints the report on the selected printer.

Exception Report

The user views and/or prints an Exception Report.

Assumptions

- Data are available for inclusion in the report.
- The user must be logged into the division in which the report is to be generated.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- VBECS does not accommodate online review (verification, signature) of reports.
- Only the first twenty patient antigens will print on the report when an exception is generated for a unit issued with unsatisfied antigen negative requirements.

Additional Information

- A supervisor may use a printed report to document the review.
- Exception Report entries reflect the actual date and time an exception occurred.
- Issue Blood Component exceptions that include the Testing Status field displays numeric values 0, 1, 2 or 3. Those numbers represent the unit's disease marker testing.
 - (0) – N/A for non-autologous units,
 - (1) – NEG – Unit fully tested and negative for all disease markers,
 - (2) – POS – Unit tested and positive for one or more disease markers, Biohazard,
 - (3) – NFT – Unit not fully tested for one or more disease markers, Biohazard.

User Roles with Access to This Option

All users

Exception Report

The user views and/or prints an Exception Report. Exception Reports list details of transactions that involve overrides to warnings and selected transactions of interest to the supervisor (for example, quarantined units).

User Action	VBECS
1. Select Reports from the main menu. Select Exception Report .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. Enter or select start and end dates.	<p>NOTES —————</p> <p>The report is sorted by exception type. Exception types are in chronological order.</p>
3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job	<ul style="list-style-type: none">• Compiles and displays the report.

User Action	VBECS
radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	
4. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer. • See Appendix B: Table 19: Details in Exception Report.

Finalize/Print TRW (See Patient Services, Finalize/Print TRW)

Inappropriate Transfusion Request Report

The user views and/or prints the Inappropriate Transfusion Request Report.

Assumptions

- The order is flagged as an inappropriate transfusion request during Accepting an Order.
- Setting Component Requirements was processed and there are inappropriate test indicators associated with the component classes in the login division.
- Data for the report were saved in VBECS.

Outcome

- The user view or prints the report.

Limitations and Restrictions

- None

Additional Information

- None

User Roles with Access to This Option

All users

Inappropriate Transfusion Request Report

The user views and/or prints the Inappropriate Transfusion Request Report, which includes inappropriate transfusion orders determined by the division's configuration of each component class.

User Action	VBECS
1. Select Reports from the main menu. Select Inappropriate Transfusion Request Report .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. Enter or select start and end dates.	<ul style="list-style-type: none">• Compiles and displays the report. <p>NOTES —————</p> <p>See Appendix B: Table 22: National Treating Specialty Table.</p>
3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none">• Compiles and displays the report.
4. Select Print to print the report and exit, or click Close to exit without	<ul style="list-style-type: none">• Prints the report on the selected printer.

User Action	VBECS
printing.	

Issued/Returned Unit Report

The user views and/or prints the Issued/Returned Unit Report.

Assumptions

- Data for the report were saved in the VBECS database.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- None

Additional Information

- The report is compiled for one division.
- Pooled units include the number of participants in the pool.

User Roles with Access to This Option

All users

Issued/Returned Unit Report

The user views and/or prints a report of units issued and/or returned to patients within a specified date range for a division.

User Action	VBECS
1. Select Reports from the main menu. Select Issued/Returned Report .	<ul style="list-style-type: none">• Lists report names.• Displays fields for entering a report date range.
2. Enter or select start and end dates.	<ul style="list-style-type: none">• Allows the user to sort the report by issued date range or by patient name.
3. In the Report Sort Order area, select the Issue/Return Date or Patient Name radio button.	<ul style="list-style-type: none">• Compiles the report and, for each unit issued and when units are returned, includes return information. <p>NOTES —</p> <p>The report includes a unit that is issued or returned within the entered date range. When a unit is identified as issued or returned during the entered date range, the report also includes unit return and issuance data, respectively, even when they fall outside the date range of the report.</p>
4. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none">• Compiles and displays the report.
5. Select Print to print the report and	<ul style="list-style-type: none">• Prints the report on the selected printer.

User Action	VBECS
exit or click Close to exit without printing.	

Medication Profile (See Patient Information Toolbar)

Order History Report

The user views and/or prints a Single-Order History Report or an Order Summary Report.

Assumptions

- Data for the report are available.

Outcome

- The user may view and/or print a report.

Limitations and Restrictions

- VBECS does not display the identity of the CPRS user who discontinues a signed order in CPRS. The CPRS order details option contains all pertinent order information.

Additional Information

- Test counts and workload information for other VBECS processes are in the Workload Report.

User Roles with Access to This Option

All users

Order History Report

The user views and/or prints an Order History Report. The user may retrieve a report for a single order or a summary report for all VBECS orders within a date range. The user may sort the Order Summary Report by treating specialty or by patient location.

User Action	VBECS
1. Select Reports from the main menu. Select Order History Report .	<ul style="list-style-type: none">• Lists report names.• Displays options for generating, viewing, and printing reports.• Displays report type and compilation criteria.
2. Click the Order Summary Report or Single-Order History Report radio button to select a report type.	
3. For an Order Summary Report, create a report date range: Select or edit the dates and times in the Start Date and End Date fields. In the Report Sort Order area, click the Treating Specialty or Ordering Location radio button and continue at Step 5.	<ul style="list-style-type: none">• Displays the option to enter a date range for the report based on the date and time VBECS acknowledged the order.• Displays an option to sort the report by treating specialty or by ordering location. The report tallies orders in a given date range by ordering location and treating specialty.
4. For a Single-Order History Report: Enter the search criteria and/or	<ul style="list-style-type: none">• Displays the Select Patient Order window.

NOTES —————

User Action	VBECS
click the find button to select an order from the list.	<p>The user enters a date or a specimen UID, or selects a patient to locate an order. VBECS displays the following information based on the selected search option:</p> <ul style="list-style-type: none"> • When a user enters a date, VBECS displays orders accepted or canceled by VBECS on that date. • When the specimen UID is entered, VBECS displays only orders associated with that specimen. • When a patient name is selected, VBECS displays the patient's orders listed in reverse chronological order for the last 90 days with an option to display all of the patient's accepted and canceled orders.
5. In the Report Compilation area, click the Create Report Now radio button and OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name, then click OK to view and print the report.	<ul style="list-style-type: none"> • Compiles and displays the report. <p>NOTES —————</p> <p>The Order Summary Report includes a section for each treating specialty or ordering location, including accepted and canceled VBECS orders and the Diagnostic Test or Component Orders for the selected date range. The report lists individual orders by order type, including the Patient Name, ID, CPRS and Lab Order Numbers, and Date.</p> <p>The Single-Order History Report contains details of an individual order.</p>
6. Click Close to exit.	

Patient History Report

The user views and/or prints the Patient History Report.

Assumptions

- VBECS displays the report header even when there is no information in the report.

Outcome

- The user views and prints a patient's VBECS record.

Limitations and Restrictions

- When units are restricted for a patient during Incoming Shipment, VBECS does not display them on the Patient History Report..
- VBECS displays these default values for transfusion reactions converted from VistA to VBECS. They do not reflect actual interpretations or results:

Converted Transfusion Reaction	Default
Pre-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Post-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Symptoms	No symptoms identified.
Implicated units	No units implicated in reaction.
Serologic tests	No serologic tests performed.

Additional Information

- The patient record report is divided into sections for report retrieval.
- This report includes:
 - o Post-transfusion details.
 - o Data from all divisions within the database.

User Roles with Access to This Option

All users

Patient History Report

The user views and/or prints the patient's VBECS record.

User Action	VBECS
1. Select Patients or Reports from the main menu. Select Patient History Report .	<ul style="list-style-type: none">• Displays options for processing patient-related functions, or for generating, viewing, and printing reports.• Displays report sections and compilation criteria.
2. Enter the patient name and/or ID, or click the find button to search for patient data. Click OK .	<ul style="list-style-type: none">• Allows the entry or selection of a patient.• Displays the patient report section types for selection:<ul style="list-style-type: none">o Patient Demographic Historyo Transfusion Requirementso Testing Detailso Unit Assignment Historyo Transfusion Historyo Specimen Detailso Special Instructionso Transfusion Reaction Detailso Exception Report Entries

User Action	VBECS
	<ul style="list-style-type: none"> o Assignments/Releases o Audit Trail Report
3. Select one or more sections to include in the report.	
4. Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none"> • Compiles and displays the report. • Allows the user to print the report. <p>NOTES —————</p> <p>When the patient has an ABO/Rh as the result of a conversion, VBECS displays a message in the Transfusion History section instructing the user to check the legacy records for additional information.</p>
5. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Prolonged Transfusion Time Report

The user views and/or prints the Prolonged Transfusion Time Report.

Assumptions

- Data for the report were saved in VBECS.

Outcome

- The user views and prints the report.

Limitations and Restrictions

- This report does not include presumed transfused units or units issued to a remote storage location.

Additional Information

- This option uses the Appendix B: Table 22: National Treating Specialty Table, not a local list.

User Roles with Access to This Option

All users

Prolonged Transfusion Time Report

The user views and/or prints the Prolonged Transfusion Time Report, which shows transfusions started more than 30 minutes after issue and units transfused for more than the maximum transfusion time set by each division.

User Action	VBECS
1. Select Reports from the main menu. Select Prolonged Transfusion Time Report .	<ul style="list-style-type: none">• Displays options for generating, viewing, and printing reports.• Displays a date range for selection.• Displays report sort order and compilation criteria.
2. Enter a date range. Click the Transfusion Date (Reverse Chronological) or Issue-To Location radio button to select a sort order.	
3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none">• Compiles and displays the report. <p>NOTES —————</p> <p>The report includes:</p> <ul style="list-style-type: none">• Issue-to location• Unit ID• Product short name• Issue Date/time• Start Date/time• End Date/time• Number of minutes elapsed between unit issue and start time, when exceeds 30 minutes• Number of minutes elapsed between the unit start and end dates and times, when they exceed the maximum transfusion setting for the component class.

User Action	VBECS
	<ul style="list-style-type: none"> • Patient Name • Patient ID • Treating Specialty <p>The report includes instances of a prolonged start time and/or the prolonged transfusion time.</p> <p>A prolonged transfusion time is the number of minutes elapsed between the transfusion start and end times and the maximum transfusion time setting for the component class. When the number of minutes exceeds the maximum transfusion time setting for the component class, the unit and its transfusion data are included in the report.</p> <p>A delayed start of transfusion is the number of minutes elapsed between the issue time and the transfusion start time. When the start time of the unit's transfusion is greater than 30 minutes from the unit's issue time, the unit and its transfusion data are included in the report. When no start time is entered, VBECS records the issue time as the start time; it does not calculate a delayed start time.</p> <p>The total transfusion time is the number of minutes elapsed between the issue and end times.</p> <p>The report includes:</p> <ul style="list-style-type: none"> • Issue-To Location • Unit ID • Product Short Name • Issue Date/Time • Start Date/Time • End Date/Time • Delayed start interval (the number of minutes elapsed between unit issue and start times, when it exceeds 30 minutes) • Prolonged transfusion interval (the number of minutes elapsed between unit start and end dates and times, when they exceed the maximum transfusion setting for the component class) • Patient Name • Patient ID • Treating Specialty
4. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Patient Testing Worklist and Testing Worklist Reports

The user views and/or prints the Patient Testing Worklist and Testing Worklist Reports.

Assumptions

- The user must be logged into the division in which the report is to be generated.

Outcome

- The user views and/or prints a report, or schedules a report for future compilation and printing.

Limitations and Restrictions

- Reports do not include Transfusion Reaction Workup (TRW) test results.
- VBECS does not accommodate online review (verification, signature) of reports.

Additional Information

- In a transfusion-only facility, VBECS does not populate the test cell or phase results sections.
- VBECS collates complete, incomplete, and invalidated unit test results entered in the user's division when the date tested falls within the supplied date range for the unit testing, reagent QC, and patient testing sections of the report.
- VBECS includes test results on all units, regardless of status, for the unit testing section. CMV and Sickle Cell test results are not included in the report.
- VBECS includes the results of the daily rack QC and patient and unit control testing for the reagent QC section.
- For the patient testing section:
 - VBECS collates as above for the selected patient specimens within a specified date range.
 - The report does not include TRW test results.
 - BR_48.08 VBECS includes XMs only in the patient testing section.
 - VBECS displays individual tests within the TAS and a statement in the Additional Information Column to this effect.
 - When Antibody Identification (ABID) testing was performed at an outside facility, VBECS prints the name of the outside facility in the Additional Information column.
- If a rack has only been tested offline (never tested within VBECS) and the QC has expired, the message "Rack was never QC'd and cannot be used." is displayed.
- The Testing Worklist Report does not display any information for Rack QC when QC Documentation Offline has been checked for that date.

User Roles with Access to This Option

All users

Patient Testing Worklist and Testing Worklist Reports

The user views and/or prints a Patient Testing Worklist Report for a specific patient or a Testing Worklist Report for a division. The reports contain both partially completed and completed testing data entries on patients, units, and reagents created in the same division as the user generating the report.

User Action	VBECS
<p>1. For a Patient Testing Worklist Report (for a single patient, select Patients from the main menu. Select Patient Testing Worklist Report. For a Testing Worklist Report for a division, go to Step 5.</p>	<ul style="list-style-type: none"> Displays options for processing patient-related functions. Displays patient search criteria.
<p>2. Enter the patient name or identification number. Click Search, as required. Select a patient and click OK.</p>	<ul style="list-style-type: none"> Requires the user to select a patient. Displays specimen search criteria.
<p>3. Select or edit the dates and times in the Specimen Received Start Date and End Date fields. Click Search, as required.</p>	<ul style="list-style-type: none"> Allows the user to indicate a date range and patient specimens to include in the report. Displays specimens with associated testing for the selected patient. <p>NOTES</p> <p>The individual Patient Testing Worklist Report includes VBECS specimen UIDs from the last 90 days, by default. The user may remove one or more specimen UIDs and/or edit the default date range.</p> <p>When the patient record is not found, VBECS notifies the user that no matching patients were found.</p> <p>When the selected patient has no specimen on file within the division, VBECS warns the user and asks the user to select a different patient.</p> <p>Requires the user to indicate the format of the report:</p> <ul style="list-style-type: none"> Patient Testing Unit Testing Daily QC Miscellaneous Reagent QC
<p>4. Click one or more check boxes to select specimens. Click OK. Go to Step 8.</p>	<ul style="list-style-type: none"> Displays the report. Allows the user to print the report.
<p>5. For a Testing Worklist Report for a division, select Reports from the main menu. Select Testing Worklist Report.</p>	<ul style="list-style-type: none"> Lists report names. Displays the report date range, sections (Patient Testing, Unit Testing, Rack QC, and Miscellaneous Reagent QC), and compilation criteria.
<p>6. Select or edit the dates and times in the Start Date and End Date</p>	<ul style="list-style-type: none"> Displays the option to enter a date range for the report based on the date and time VBECS acknowledged the order.

User Action	VBECS
fields.	
7. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report.	<ul style="list-style-type: none"> Compiles and displays the report.
8. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> Prints the report on the selected printer.

Transfusion Complications Report

The user views and/or prints the Transfusion Complications Report.

Assumptions

- VistA tests are defined in Transfusion Complications.
- The connection to VistA is active.
- Transfusion Complications was processed and post-transfusion complication tests are defined.
- Units have a status of presumed transfused or transfused for the patient.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- The report does not include units in an issued status.

Additional Information

- The report searches for evidence of a transfusion in VBECS only. If the user enters a start date that requires a search of VistA records to check back a full 12 months, the report includes a message indicating that only transfusions documented in VBECS are included in the report and that the user needs to run a supplementary report in VistA.
- The report searches for VistA tests defined in Transfusion Complications when the report is requested.
- If threshold values are words, they are case insensitive.
- A patient's test results and transfusions are not directly linked by date, episode, or admission event.
- Lab test data retrieved for the report are not saved in VBECS.
- The report is not saved in VBECS.

User Roles with Access to This Option

All users

Transfusion Complications Report

The user views and/or prints the Transfusion Complications Report, which includes patients transfused in his division with a test result that exceeds the threshold set in Transfusion Complications for that test. Tests associated with possible post-transfusion complications are gathered by searching for results in a user-defined date range for the report. The report also includes the documented patient transfusions from 12 months prior to the start date of the report through the end date of the report.

User Action	VBECS
1. Select Reports from the main menu. Select Transfusion Complications Report .	<ul style="list-style-type: none">• Lists report names.• Displays fields for entering a report date range.
2. Enter or select start and end dates.	<ul style="list-style-type: none">• Compiles and displays the report and allows the user to print it. <p>NOTES —————</p> <p>VBECS searches the Laboratory files (VistALink) and gathers instances of tests that exceed the threshold result of the tests set</p>

User Action	VBECS
	<p>in Transfusion Complications within the specified date range. Retrieved test results are then filtered so that the report includes results only for patients who have transfusions documented in VBECS within the last 12 months prior to the start date of the report through the end date.</p> <p>VBECS includes the evaluated tests and their threshold values in the report criteria section.</p> <p>Details of the report include:</p> <ul style="list-style-type: none"> • Patient name • Patient ID • Current hospital location • Date of death (in lieu of location if patient died) • Unit ID • Unit short name • End date and time of the unit's transfusion, if available <p>From VistA for each patient's tests:</p> <ul style="list-style-type: none"> • Test name • Associated test result • Date and time recorded
<p>3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or</p> <p>Click the Schedule Print Job radio button and OK, then enter or select a print date and printer name. Click OK to view and print the report.</p> <p>4. Select Print to print the report and exit, or click Close to exit without printing.</p>	<ul style="list-style-type: none"> • Compiles and displays the report. <p>NOTES —————</p> <ul style="list-style-type: none"> • This report requires a logged on VistA user to retrieve the lab results. It cannot be scheduled to print. <ul style="list-style-type: none"> • Prints the report on the selected printer.

Transfusion Effectiveness Report

The user views and/or prints a Transfusion Effectiveness Report.

Assumptions

- VistA tests are defined in Transfusion Effectiveness.
- The connection to VistA is active.
- Units have a status of presumed transfused or transfused for the patient.
- Data for the report were saved in VBECS and in VistA.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- The report does not include VBECS diagnostic tests.

Additional Information

- Only verified Laboratory test results are available from VistA.
- Laboratory test data retrieved for the report from VistA are not saved in VBECS.
- A new page is started for each patient in the report.
- Test results in the report do not include units of measure or normal ranges.

User Roles with Access to This Option

All users

Transfusion Effectiveness Report

The user views and/or prints the Transfusion Effectiveness Report, which contains transfusion episodes and Laboratory results for patients for a selected period.

User Action	VBECS
1. Select Reports from the main menu. Select Transfusion Effectiveness Report .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. Enter or select start and end dates.	
3. In the Report Type area, select the All Patients radio button, or Click the Individual Patient Report radio button and the find button to search for a patient.	<ul style="list-style-type: none">• Compiles and displays the report.
4. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer	<ul style="list-style-type: none">• Compiles and displays the report. <p>NOTES —</p> <ul style="list-style-type: none">• This report requires a logged on VistA user to retrieve the lab results. It cannot be scheduled to print.

User Action	VBECS
name. Click OK to view and print the report.	
5. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Transfusion Reaction Count Report

The user views and/or prints a Transfusion Reaction Count Report.

Assumptions

- Data for the report were saved in VBECS.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- None

Additional Information

- The report is available as a single-division or multidivisional (all divisions) report.

User Roles with Access to This Option

All users

Transfusion Reaction Count Report

The user views and/or prints the Transfusion Reaction Count Report for a single division or for all divisions in the database.

User Action	VBECS
1. Select Reports from the main menu. Select Transfusion Reaction Count Report .	<ul style="list-style-type: none">• Lists report names.• Displays report division, format, and compilation criteria. <p>NOTES —————</p> <p>The report includes transfusion reactions that are ordered, pending, completed, or finalized, and those indicated during Post-Transfusion Information.</p>
2. Enter or select start and end dates.	<ul style="list-style-type: none">• Compiles and displays the report.• Allows the user to print the report.
3. Click the Current Divisions or All Divisions radio button.	<ul style="list-style-type: none">• Allows the user to select one or all divisions.• Displays a date range for selection. <p>NOTES —————</p> <p>The Transfusion Reaction Count Report by type includes specific patient information indicated by the user. The default report setting excludes patient names.</p>
4. Click the Detailed Report or Summary Report radio button.	<ul style="list-style-type: none">• Displays an option to compile the report as detailed or summary.• Compiles and displays the report.• Allows the user to print the report. <p>NOTES —————</p> <p>The Transfusion Reaction Count Report prints the implicated unit numbers or "no implicated units," as appropriate.</p> <p>VBECS reports the total number of patient transfusion reactions</p>

User Action	VBECS
	in the selected date range.
5. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or Click the Schedule Print Job radio button and OK , then enter or select a print date and printer name. Click OK to view and print the report, or click Close to exit without printing.	<ul style="list-style-type: none"> Compiles and displays the report.
6. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> Prints the report on the selected printer.

Transfusion Requirements Report

The user views and/or prints a Transfusion Requirements Report.

Assumptions

- The database conversion of VistA patient files is complete.
- Data for the report were saved in VBECS.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- Deceased patient records do not appear on the report when they are generated after the patient has been marked as deceased in the VBECS database via a VistA Death Notification Update.

Additional Information

- A Transfusion Requirements report is compiled from all VBECS divisions defined in the VBECS database.
- VBECS does not allow a report to be printed when no data are available.
- When the date range includes a future date and time, VBECS displays “Partial” in the report header on each page of the report.
- Each page containing patient information includes the patient name.
- Divisions in a multidivisional database may view each other’s Special Instructions (SI) comments.
- Important transfusion related information for a patient may appear on more than one page. Always check the next page for additional patient information.

User Roles with Access to This Option

All users

Transfusion Requirements Report

The user views and/or prints the Transfusion Requirements Report for a division. This option describes the continuous compilation of the complete file of patients in the VBECS database for historic record checking during downtime.

User Action	VBECS
1. Select Reports from the main menu. Select Transfusion Requirements Report .	<ul style="list-style-type: none">• Lists report names.• Displays report compilation criteria.• Displays a date range for selection.
2. In the Report Type area, select the Cumulative Report radio button, or Click the Partial Report radio button and enter or select start and end dates.	<ul style="list-style-type: none">• Compiles and displays the report. <p>NOTES —————</p> <p>Patients who do not have a blood type, SI, TR, or Transfusion Reaction in their VBECS record do not display on the Transfusion Requirements Report in either the Cumulative or Partial format.</p> <p>The Cumulative Report displays patients that have a blood type,</p>

User Action	VBECS
	<p>SI, TR, or Transfusion Reaction.</p> <p>The Partial Report displays patients that have a change to their VBECS record within the date range of the requested report. Changes are defined as:</p> <ul style="list-style-type: none"> • New order • Record of testing • Patient Update event • New SI entered • New TR entered • New Transfusion Reaction entered
<p>3. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or</p> <p>Click the Schedule Print Job radio button and OK, then enter or select a print date and printer name. Click OK to view and print the report.</p> <p>4. Select Print to print the report and exit, or click Close to exit without printing.</p>	<ul style="list-style-type: none"> • Compiles and displays the report.
	<ul style="list-style-type: none"> • Prints the report on the selected printer.

Unit History Report

The user views and/or prints a Unit History Report.

Assumptions

- The unit exists in the user's division.

Outcome

- The user views and/or prints the report.

Limitations and Restrictions

- The report is not multidivisional.
- When a user retroactively updates a unit status through Discard or Quarantine, the Unit History Report displays the updated information, but does not display the date the change was made.

Additional Information

- There may be multiple instances of the same unit number and product code in the same division. Inactivated unit records may have a unit ID and product code that duplicate an active unit ID and product code.
- VBECS always displays the current unit status when the user selects a blood unit, so there is no need for a separate report option, as is possible in VistA.
- Unit testing details are included in Testing Worklist Report.

User Roles with Access to This Option

All users

Unit History Report

The user views and/or prints the complete record for a blood unit, from incoming shipment through final status, including unit activities saved in the VBECS database. VBECS compiles data for that unit for one division.

User Action	VBECS
1. Select Reports from the main menu. Select Unit History Report .	<ul style="list-style-type: none">• Lists report names.• Allows the user to scan or enter a unit ID and product code.• Displays compilation criteria.
2. Scan or enter the unit ID and product code, or click the find button to search for and select a unit.	<p>NOTES —</p> <p>The user may select any unit (active or inactive) in a division. VBECS marks the units so that the active ones are distinct from the inactive ones.</p>
3. In the Include Report Sections area, click one or more check boxes to select report sections.	<ul style="list-style-type: none">• Displays the unit report section types for selection:<ul style="list-style-type: none">o Incoming Shipmento Outgoing Shipmento Unit Quarantine Informationo Unit Financialso Patient Associations Historyo Issue Informationo Testing

User Action	VBECS
	<ul style="list-style-type: none"> o Crossmatch o Modification o Transfused To o Exception Report
<p>4. In the Report Compilation area, click the Create Report Now radio button. Click OK to view the report, or</p> <p>Click the Schedule Print Job radio button and OK, then enter or select a print date and printer name. Click OK to view and print the report.</p>	<ul style="list-style-type: none"> • Compiles and displays the report for the user's login division. <p>NOTES —————</p> <p>The page header includes the standard VBECS header, the unit ID, product, donation type, expiration date, restricted for patient, status, ABO/Rh, division, product long name, and user.</p> <p>VBECS displays entries from only the related reports (Audit Trail and Exception Reports) specific to the selected unit or patient.</p> <p>Audit Trail and Exception Reports that record only unit override information are linked to the unit's historic record. Those that record only patient override information are linked to the patient's historic record. Those that record both are accessible in unit and patient historic records.</p> <p>See Appendix B: Table 18: Details in Audit Trail Report and Table 19: Details in Exception Report.</p> <p>The report displays the product long name.</p> <p>The report displays "CMV neg" and/or "SC neg" in the "Special Testing" column when a unit qualifies.</p> <p>The report displays the unit's quarantine indicator in the "Quarantined" column: blank when never quarantined, "Yes" when quarantined, and "No" when released from quarantine.</p> <p>The report displays the unit's biohazardous indicator in the "Biohazard" column: blank when never marked as biohazardous, "Yes" when biohazardous, and "No" when the biohazard indicator was removed.</p>
5. Select Print to print the report and exit, or click Close to exit without printing.	<ul style="list-style-type: none"> • Prints the report on the selected printer.

This page intentionally left blank.

References

- Association for the Advancement of Medical Instrumentation, Arlington, Virginia. *The Quality System Compendium: GMP Requirements & Industry Practice*, 1998.
- Association for the Advancement of Medical Instrumentation. *Medical device software—Software life cycle processes (ANSI/AAMI SW68:2001)*, 2001.
- Brandt, Gwyneth. ClinicalTrialJobs.com. “Preparing for FDA Inspections.” July 18, 2002.
- *Clinical Laboratory Technical Procedure Manuals; Approved Guideline—Fourth Edition*, NCCLS document GP2-A4 (Wayne, PA: National Committee for Clinical Laboratory Standards, 2002).
- CMMI Product Development Team, Carnegie Mellon Software Engineering Institute. *CMMISM-SE/SW, V1.0, Capability Maturity Model®—Integrated for Systems Engineering/Software Engineering, Version 1.0* (Pittsburgh, PA: Carnegie Mellon Software Engineering Institute, August 2000).
- Commission on Laboratory Accreditation, Laboratory Accreditation Program. *Laboratory General Checklist* (Northfield, IL: College of American Pathologists, 2005).
- Commission on Laboratory Accreditation, Laboratory Accreditation Program. *Transfusion Medicine Checklist* (Northfield, IL: College of American Pathologists, 2005).
- *Comprehensive Accreditation Manual for Hospitals: The Official Handbook (CAMH)*. (Oakbrook Terrace, Illinois: Joint Commission on Accreditation of Healthcare Organizations, 2004).
- Joint Commission on Accreditation of Healthcare Organizations. *Comprehensive Accreditation Manual for Hospitals: The Official Handbook (CAMH)*. (Oakbrook Terrace, Illinois: Joint Commission on Accreditation of Healthcare Organizations, 2002).
- Cotliar, William and Larry Bauer. *Complying with FDA Good Manufacturing Practice Requirements: How to Develop Your Quality System Manual*. 5th edition. (Arlington, VA: Association for the Advancement of Medical Instrumentation, 1996).
- Department of Veterans Affairs, Software Engineering Process Group (SEPG) and Software Quality Assurance (SQA), *National Documentation Standards* (October 2005).
- Department of Veterans Affairs, Veterans Health Administration. *VHA Directive 99-053, VistA Blood Bank Software*, November 10, 1999.
- Department of Veterans Affairs, Veterans Health Administration. *VHA Directive 2005-029, Transfusion Verification and Identification Requirements for All Sites*, July 1, 2005.
- Department of Veterans Affairs: Office of Acquisition and Materiel Management. *Veterans Affairs Acquisition Regulation* (Washington, DC: 1997).
- Galloway, Dianne. *Mapping Work Processes*. (Milwaukee, Wisconsin: ASQ Quality Press, 1994).
- *General Principles of Software Validation; Final Guidance for Industry and FDA Staff* (U.S. Department Of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health, Center for Biologics Evaluation and Research, January 11, 2002).
- *Guidance for Industry and FDA Premarket and Design Control Reviewers, Medical Device Use-Safety: Incorporating Human Factors Engineering into Risk Management* (U.S. Department of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health, Division of Device User Programs and Systems Analysis, Office of Health and Industry Programs, July 18, 2000).
- *Guidance for Industry, FDA Reviewers and Compliance on Off-The-Shelf Software Use in Medical Devices* (U.S. Department Of Health And Human Services, Food and Drug Administration, Center for Devices and Radiological Health, Office of Device Evaluation, September 9, 1999).
- *Guidelines for Implementing the Electronic Crossmatch* (Bethesda, MD: American Association of Blood Banks, 2003).

- *Guidelines for the Labeling of Specimens for Compatibility Testing*. Ed. James P. AuBuchon, MD. (Bethesda, MD: American Association of Blood Banks, 2002).
- *Guidelines for the Laboratory Evaluation of Transfusion Reactions* (Bethesda, MD: American Association of Blood Banks, 2003).
- *IEEE Standard for Software Configuration Management Plans* (IEEE Std 828-1998) (New York: Institute of Electrical and Electronic Engineers, Inc., 1998).
- *ISBT 128 Bar Code Symbology and Application Specification for Labeling of Whole Blood and Blood Components*, Version Number 1.4.0 [York, PA: International Council for Commonality in Blood Banking Automation (ICCBBA), June 2001].
- Jacobson, Ivar, Grady Booch, and James Rumbaugh. *The Unified Software Development Process* (Boston: Addison-Wesley, 1999).
- *Medical Device Reporting—Baseline Report*, Form FDA 3417.
- *Medical devices—Application of risk management to medical devices* (ANSI/AAMI/ISO 14971:2000), 2001. Association for the Advancement of Medical Instrumentation.
- Nevalainen, D. E. and M. F. Callery. *Applying Self-Assessment and Process Analysis to Create a Quality System* (Basel, Switzerland: S. Karger Publishing, 1995).
- Nozick, Robin F. “Demystifying the Information Systems Validation Process,” *AABB News* 7:6 (2005): 50–51.
- *Object-Oriented Analysis and Design Using the UML: Student Manual, Book 1, Version 4.2* (Cupertino, CA: Rational Corporation, 1999).
- *Rational ClearQuest® Introduction. Version: 2002.05.00. Part Number: 800-025123-000*. Windows/UNIX. ©1997–2002 Rational Software Corporation.
- *Reviewer Guidance for a Premarket Notification Submission for Blood Establishment Computer Software* (U.S. Department Of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health, Center for Biologics Evaluation and Research, Office of Blood research and Review, Division of Blood Applications (January 13, 1997).
- *Software Documentation Style Standards*, SOP 11.
- *Standards for Blood Banks and Transfusion Services*. Twenty-third edition (Bethesda, MD: American Association of Blood Banks, 2005).
- *Technical Manual. Fifteenth edition* (Bethesda, MD: American Association of Blood Banks, 2005).
- *The Chicago Manual of Style*, 15th ed., (Chicago: The University of Chicago Press, 2003).
- *Title 21 Code of Federal Regulations*, Part 11, Electronic Records; Electronic Signatures.
- *Title 21 Code of Federal Regulations*, Part 803, Medical Device Reporting.
- *Title 21 Code of Federal Regulations*, Part 806, Medical Devices; Reports of Corrections and Removals.
- *Title 21 Code of Federal Regulations*, Part 807, Establishment Registration and Device Listing for Manufacturers and Initial Importers of Devices.
- *Title 21 Code of Federal Regulations*, Part 810, Medical Device Recall Authority.
- *Title 21 Code of Federal Regulations*, Part 820, Quality System Regulation.
- *Title 21 Code of Federal Regulations*, Parts 200 to 299, Food and Drugs.
- *Title 21 Code of Federal Regulations*, Parts 600 to 799, Food and Drugs.
- *Title 21 Code of Federal Regulations*, Parts 800 to 1299, Food and Drugs.
- Trautman, Kimberly A. *The FDA and Worldwide Quality System Requirements Guidebook for Medical Devices* (Milwaukee, Wisconsin: ASQ Quality Press, 1997).
- U.S. Department Of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health. *Design Control Guidance for Medical Device Manufacturers*, March 11, 1997.

- U.S. Department Of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health. *Inspection of Medical Device Manufacturers, Final Guidance for Industry and FDA*, February 7, 2001.
- U.S. Department Of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health. *Medical Device Quality Systems Manual: A Small Entity Compliance Guide*. First edition, December 1996.
- U.S. Department Of Health And Human Services, Food and Drug Administration, Center for Devices and Radiological Health, Office of Device Evaluation. *Guidance for the Content of Premarket Submissions for Software Contained in Medical Devices*, May 29, 1998.
- U.S. Department Of Health And Human Services, Food and Drug Administration, Center for Devices and Radiological Health, Office of Regulatory Affairs. *Glossary of Computerized System and Software Development Terminology*, August 1995.
- U.S. Department Of Health and Human Services, Food and Drug Administration, Center for Devices and Radiological Health. *Premarket Approval Manual*. January 1998.
- U.S. Food and Drug Administration, Office of Regulatory Affairs and the Center for Devices and Radiological Health. *Guide to Inspections of Medical Device Manufacturers*, December 1997.
- U.S. Food and Drug Administration, ORA Inspectional References, Quality Inspection System Technique. *Guide to Inspections of Quality Systems*, August 1999.
- *United States Industry Consensus Standard for the Uniform Labeling of Blood and Blood Components Using ISBT 128, Version 1.2.0*, prepared by the International Council for Commonality in Blood Banking Automation, Inc. (U.S. Department of Health and Human Services, Food and Drug Administration, Center for Biologics Evaluation and Research, November 1999).
- VHA Office of Information, System Design & Development, Software Engineering Process Group, *RequisitePro Requirements Plan*, October 2000.
- *Veterans Health Administration Software Documentation Handbook*.
- VistA Data Systems & Integration Service (VDSI), VHA Office of Information—System Design & Development. *VistA Data Systems & Integration Review Program*. VistA Data Systems & Integration Service (VDSI), VHA Office of Information—System Design & Development, April 2002].
- Wilson, Paul F., Larry D. Dell, and Gaylord F. Anderson. *Root Cause Analysis: A Tool for Total Quality Management* (Milwaukee, Wisconsin: ASQ Quality Press, 1993).
- *ZPL II Programming Guide*, Rev. 1, Customer order # 46530L, Manufacturer part # 46530LB (Zebra Technologies Corporation, 1998).

This page intentionally left blank.

Glossary

Acronym, Term	Definition
AABB	American Association of Blood Banks.
ABID	Antibody identification.
ABO	A group for classifying human blood, based on the presence or absence of specific antigens in the blood, which contains four blood types: A, B, AB, and O. The ABO group is the most critical of the human blood systems. It is used to determine general compatibility of donor units to a recipient.
ABS	Antibody screen, Antibody Screen Test.
Access Code	A field in the VistA New Person file used to uniquely identify a user on the VistA system.
Accessioned	Describes an order when a specimen has been collected or has been placed on a current collection list by the VistA Laboratory software. <i>A mechanism for accessioning orders for blood components that do not require a specimen needs to be designed.</i>
Active order	A patient order for blood product(s) that is not completed or expired. Depending on the component class, an accompanying specimen may be required.
Active user	A user with a setting in VistA that allows him to access VBECS. Inactive users may no longer access VBECS.
Activities	Business processes that the product supports.
Activity	An identified verification and validation task.
ADPAC	Automated Data Processing Application Coordinator.
AHG	Anti-Human Globulin.
Aliquot	A sample/portion divided off or separated as a part of a whole.
Allogeneic	See "Homologous."
Antibody	A protein in the serum of some people that will react to a specific antigen on the blood cells. <i>In the case of a red blood cell transfusion, a patient with a clinically significant (leading to a transfusion reaction and possible patient harm) identified antibody must receive only red cells that are typed and found negative for the associated antigen. Once a clinically significant antibody is identified, antigen-negative cells must always be transfused even if a current specimen no longer shows the presence of the antibody. The antibody table in VBECS addresses the antibodies associated with red cell antigens. Though many antibodies are considered clinically significant, some are more of a nuisance to testing than a problem for transfusion. Many antibodies have documentation that defines their clinical significance; that information is used to define the VBECS table system warnings and overrides.</i>
Antigen	A substance on the surface of a red cell that stimulates an immune response (formation of an antibody).
Antigen negative requirement	A prerequisite that specified blood product categories must not contain stated characteristics.

Acronym, Term	Definition
Antigen typing	Testing for the presence or absence of an antigen on a red blood cell.
Antiglobulin	An antibody that combines with a specific protein on a different type of antibody.
Anti-Human Globulin	A commercial reagent that binds to immunoglobulin molecules attached to the surface of a red blood cell.
Antisera	Plural of “antiserum.”
Antiserum	Blood serum (human, animal, monoclonal) containing antibodies that are specific for one or more antigens.
Assigned	Describes a unit of blood in a state that can be issued for transfusion and assigned to a patient. See “Selected.”
Attribute	ISBT 128 term that describes a property of a collected blood product. <i>Examples include “Irradiated” and “Residual Leukocyte Content” (Leukopoor).</i>
Audit Trail Report	The report of changes made to previously verified data. <i>Regulations require that a mechanism be in place to prevent “electronic white-out.” An audit report contains details of the change, before and after data, date, and who made the change.</i>
Autocontrol	A test using a patient’s serum (plasma) and cells under the same conditions as the accompanying test.
Autodiscriminate	Describes the ability of a barcode scanner to identify and interpret more than one barcode symbology without having to adjust any software or hardware parameters.
Autologous	Describes blood donated in advance for use by the donor in the future, often for elective surgery. <i>Regulations allow for processing patients with diseases that otherwise would preclude them from donating blood to donate for themselves autologously. It is imperative that software guarantee that such blood is never available for use by anyone other than the intended recipient.</i>
Automated problem-reporting system (also referred to as the problem-reporting system)	See “Remedy.”
Barcode scanner	An electronic device that uses a finely focused beam of light to sense and transmit the data in a barcode image.
BCE	Bar Code Enhancement.
BCE-PPI	Bar Code Enhancement – Positive Patient Identification.
Billable	A characteristic of an LMIP workload code that determines if statistics gathered as a result of that workload code should be included in the monthly report to the Austin Automation Center.
Blood product	See “Component Type.”
Blood Transfusion Record Form (BTRF)	A form printed in the blood bank used to record the actual events of a transfusion. <i>Examples: start/stop time, transfusionists, vitals, reaction information. Formerly “SF 518 equivalent” (VistA nomenclature).</i>
Blood type (blood group)	The ABO and Rh type of blood determined by testing.

Acronym, Term	Definition
Canned comment	Defined comments, entered in abbreviated form and displayed as expanded text, used to enter supporting information at a variety of data entry points.
Category type	When used in relation to the file of canned comments, the category of the comment that determines in which functionality the comment may be selected and associated with other data.
Caution Tag	A blood bank-generated form attached to a unit of blood when the unit is assigned or crossmatched.
CC	Check Cell.
Check Cell (CC)	A type of reagent red cell used in the blood bank to verify the reactivity of antiglobulin antisera when the results of an antiglobulin test are negative. <i>Also referred to as "Coombs Control cells."</i>
Check character	A character used to check the integrity of a barcode symbol.
Clinically significant antibodies	If there is a chance of an adverse reaction to a transfusion, the antibodies associated with the antigens in that unit are considered clinically significant. <i>Only clinically significant antibodies require that the patient receive red cells that have been phenotyped as negative for the associated antigen.</i>
CMS	Centers for Medicare and Medicaid Services. The federal agency responsible for administering the Medicare and Medicaid programs as well as certain other healthcare entitlement programs. "CMS" is the new name for HCFA (Health Care Financing Administration) under the Bush administration.
CMV	Cytomegalovirus. Some patients should receive only blood units that have been tested and found negative for CMV.
COTS	Commercial-Off-The-Shelf.
Codabar	Standard symbology used to define the format of barcode labels used to label all blood units. <i>ISBT 128 will replace Codabar.</i>
Compatibility testing	Tests performed in the blood bank to determine the suitability of a unit of blood for a patient. <i>Depending on the component class, compatibility testing can be as limited as comparing donor and patient ABO/Rh or as complex as requiring crossmatch.</i>
Completed TRW	A completed TRW indicates that serologic testing was completed.
Component	Any raw material, substance, piece, part, software, firmware, labeling, or assembly that is intended to be included as part of the finished, packaged, and labeled device.
Component class	The general form of a specific product produced by a blood supplier. <i>Example: Red Blood Cells, Fresh Frozen Plasma, and Platelets. (Product Type in VBECS model)</i>
Coombs test	See "Anti-Human Globulin."
CPRS	Computerized Patient Record System. A VistA software application that provides an integrated patient record system for use by clinicians, managers, quality assurance staff, and researchers.

Acronym, Term	Definition
CPT	Current Procedural Terminology. A proprietary terminology system owned by the AMA that enumerates medical and surgical services typically provided by physicians. <i>Although CPT is proprietary, CMS (formerly HCFA) includes it in HCPCS and requires its use for outpatient claims submitted for Medicare and Medicaid beneficiaries.</i>
Crossmatch	A procedure that demonstrates ABO compatibility and the absence of clinically significant antibody reactions between a patient and a unit of blood. <i>Traditionally, a crossmatch consists of a serological test between the patient's blood and a sample of the blood unit to be transfused. New regulations allow for what is called "computer/electronic crossmatch" in defined circumstances.</i>
Crossover	A process by which an autologous or directed donor unit is released from the intended recipient and made available for general transfusion. <i>VBECS supports only directed unit crossover.</i>
D	The major antigen of the Rh system. Also "Rh _o ." <i>The presence of this antigen on red blood cells determines if a patient or donor is Rh Pos or Rh Neg. A weak expression of the D antigen, referred to as "Du" or "weak D," may require additional testing to determine its presence. "D negative" is another term for "Rh negative"; "D control" is another term for "Rh control." The convention for VBECS GUIs is D control.</i>
Data entry grid	A grid that contains observations and interpretations of observations.
Database	A collection of data arranged for ease and speed of retrieval that may or may not contain data from multiple VBECS divisions.
Day of use	The calendar date times (0000–2359) a reagent is used.
DDE	Direct data entry.
Decision support system (DSS)	A financial/cost tracking system for VA hospitals.
Delete	Remove a VBECS record from user access and keep it in the database.
Direct Antiglobulin Test (DAT)	A laboratory test performed in a blood bank in conjunction with pre-transfusion testing or as a diagnostic test. <i>A positive DAT indicates that an antibody-antigen reaction may have occurred in vivo, which may be related to medications or previous transfusions. "Direct Coombs" and "DAT" are used interchangeably.</i>
Direct data entry (DDE)	The process of entering test results in a computer rather than recording them on paper.
Directed	Blood that is donated in advance to be set aside for a specific patient. <i>Directed donation units must meet all of the requirements of any regular unit of donated blood (all disease marker testing must be negative). If not needed by the intended recipient, a directed unit may be released for general use.</i>
Disposition	The recorded final state of a blood unit. <i>A blood unit with a disposition is no longer available for routine processing.</i>
Edit	To enter or change previously saved information. <i>VistA overwrites the data and does not create an audit trail. VBECS overlays the data and creates an audit trail.</i>

Acronym, Term	Definition
Elution	A procedure performed by a blood bank technologist to investigate a positive DAT. <i>A positive eluate is further tested to determine the specificity of the antibody, which has been eluted from the patient's red cells.</i>
Emergency issue	The release of blood prior to pre-transfusion testing and crossmatch.
Exception Report	A report that details a user deviation from a standard process. <i>Example: a user overrides a warning and is required to enter a justification for continuing.</i>
eXM	Electronic crossmatch.
Expiration date	The last date during which a blood component can be considered safe for transfusion. <i>If the expiration date does not include a time, the time is assumed to be midnight.</i>
Eye-readable prefix	Two non-scannable numbers appended to a Codabar Unit ID label by some collection facilities. <i>Most Red Cross blood suppliers use eye-readable prefixes; not all private blood suppliers do. Eye-readable prefixes correlate with collection facilities.</i>
Facility (see "Site")	A blood bank laboratory where a blood unit is located.
FDA	Food and Drug Administration.
File	Source code, data, and documentation.
Finalized TRW	A finalized TRW indicates that the medical director reviewed the TRW, classified the reaction, and may make recommendations for future transfusions.
FOIA	Freedom of Information Act.
HCPCS	Healthcare Common Procedure Coding System. Codes used to enumerate medical and surgical services and supplies in the outpatient setting. <i>CMS (formerly HCFA) uses HCPCS for Medicare and Medicaid reimbursement. HCPCS are the required code set for medical and surgical services and supplies for the outpatient setting under the final code set standards for HIPAA. There are three levels of HCPCS codes: level I are CPT codes from AMA; level II are additional codes developed by HCFA for services and supplies not covered by CPT; level III are local codes assigned by individual state Medicare carriers. Level III codes are scheduled to be phased out under HIPAA.</i>
HIPAA	Health Insurance Portability and Accountability Act of 1996.
Historical ABO/Rh	The most recent valid entry entered in the same division as the current order/treatment record, not from the database conversion.
Historical record	Cumulative, redundant data, referenced continually throughout VBECS.
HL7	Health Level 7. An international messaging standard used predominantly in the healthcare industry.
HLA (histocompatibility locus antigens)	Antigens associated with cellular immunity.
Homologous (also "allogeneic")	The standard type of blood donation: a volunteer donor reports to a blood collection facility and donates a unit of blood. The blood is tested and released for general use.

Acronym, Term	Definition
ICCBBA	International Council for Commonality in Blood Banking Automation, Inc. The governing body that developed and maintained the database that governs the ISBT 128 blood labeling system.
IEN	Internal Entry Number.
IgA (Immunoglobulin A)	Immune protein usually associated with tissue response (hypersensitivity).
In-date	Describes the state of a unit of blood when the expiration date/time of the unit is in the future.
Interpretation	The evaluation of serologic test reactions.
IRM	Information Resource Management.
Irradiate	To expose blood products to radiation to damage white cells and avoid graft-versus-host disease.
ISBT	International Society of Blood Transfusion.
ISBT 128	New standardized method for labeling blood units. <i>It includes much more information within the barcodes than is possible with Codabar. Its rules are much more structured than Codabar. It is not backwards compatible with Codabar.</i>
Leukopoor	See "leukoreduce."
Leukoreduce	To remove white cells from a cellular product through filtration. The resulting product is considered "leukoreduced" or "leukopoor."
ListServ	An email distribution list to manage discussion groups.
LMIP	Laboratory Management Index Program. Coding system used by the College of American Pathologists (CAP) to describe various procedures performed in a clinical laboratory with the purpose of assisting in calculating productivity.
Medical device	An instrument, apparatus, implement, machine, contrivance, implant, in vitro reagent, or other similar or related article, including a component part, or accessory that is intended for use in the diagnosis of disease or other conditions, or in the cure, mitigation, treatment, or prevention of disease.
Modifier	An ISBT 128 term that describes a characteristic of a donated unit of blood that distinguishes it from other members of the component class. <i>The VBECS application is concerned with the modifiers "Washed," "Thawed," "Liquid," "Frozen," "Rejuvenated," and "Deglycerolized."</i>
Modifier (order)	Component requirements associated with a specific component order. <i>Example: a component order where red blood cells must be leukoreduced. When a component order is accompanied by special requirements, units issued for transfusion must meet specific criteria that translate into "Modifiers" or "Attributes" in ISBT 128.</i>
MSBOS	Maximum Surgical Blood Ordering Schedule. The appropriate number of specific blood components to be ordered for selected surgeries at a facility. <i>Preoperative blood component orders are compared with this list and details of orders falling outside of these established criteria are recorded and reportable.</i>
Multidivisional database	More than one separate transfusion services share a database.

Acronym, Term	Definition
Negative control	The negative control cell is known to lack the antigen the antiserum type is to detect. This cell should demonstrate negative serologic reaction, and, therefore, have a valid negative interpretation.
	The positive and negative control cells must demonstrate a valid positive or negative set of serologic reactions and interpretations for the test of the unit or patient to be considered valid, regardless of the row validation of the unit test itself. (See “Positive control.”)
NFT	Not fully tested.
NLT	National Laboratory Test code. Coding schema used by the VistA Laboratory software to assign numeric identifier code to specific tests. Also known as Workload Code (WKLD CODE).
OI NT&EO	Office of Information National Training and Education Office.
Order	Request for a service. <i>The request may specify when the order is to be completed, e.g., immediately (STAT). VBECS acts on the order when it has been “accessioned.”</i>
Original unit	A unit before modification.
Outdated	Describes the state of a unit of blood that has passed the expiration date/time.
Patch	Software developed and released to correct a defect, or change or enhance existing functionality; usually a modified subset or small addition to existing software.
Patient care encounter (PCE)	A contact with a patient documented in the patient's record.
PCE	Patient care encounter.
Pending Order List (POL)	A list of orders transmitted from CPRS waiting to be processed.
Pending Task List (PTL)	A list of tasks to be processed in the blood bank.
Phenotype	Describes the presence or absence of specific antigens.
Phenotyping	Testing for the presence or absence of specific antigens. See “Antigen typing.”
Pick list	A system-generated list that allows a user to select one or more entries.
POL	Pending Order List.
Pool	A blood product from many donors stored for later use.
Positive control	The positive control cell is known to have the antigen that the reagent antiserum type is to detect. This cell should demonstrate positive serologic reaction and, therefore, have a valid positive interpretation.
	The positive and negative control cells must demonstrate a valid positive or negative set of serologic reactions and interpretations for the test of the unit or patient to be considered valid, regardless of the row validation of the unit test itself. (See “Negative control.”)
Potential hazard	An event that compromises the safety of the patient or blood supply.
Pre-transfusion testing	A group of tests, including ABO/Rh and an antibody screen, that must be performed on a sample of blood prior to crossmatching red cells for transfusion.
Product	Components, manufacturing materials, in-process devices, finished devices, and returned devices.

Acronym, Term	Definition
Product type	See "Component class."
Presumed transfused	Unit status that VBECS assigns 48 hours after issue of a blood product.
QC	Quality control. Testing routinely performed to ensure the proper functioning of materials, equipment, and methods.
Quality	The totality of features and characteristics that affect the ability of a device to satisfy fitness-for-use, including safety and performance.
Reagent rack (also "rack")	A group of reagents used in routine blood bank testing.
Reagent red blood cells	Non-patient testing materials, usually from a commercial source.
Reflex order	An order for additional testing to be performed resulting from observations of previous testing.
Reflex test	A laboratory test not ordered by a physician but performed in response to results of other laboratory tests.
Remedy	A computer application for reporting, logging, and tracking problems associated with the daily operation of computer systems within VistA.
Rh (also "Rh factor")	Rhesus factor. The major antigen of the Rh system. Also "D." <i>The presence of this antigen on red blood cells determines if a patient or donor is Rh Pos or Rh Neg. A weak expression of the D antigen, referred to as "Du" or "weak D," may require additional testing to determine its presence.</i>
SC	Screening cells.
SCD	Sterile connection device.
Screening cell (SC)	Antigenically defined cells used to detect alloantibodies.
Selected	Status of a unit of red cells when the blood bank technician selects it from available inventory to be crossmatched with a patient specimen. <i>A unit in "selected" state is not available to be issued for transfusion. By entering a crossmatch result of "compatible," the state changes from "selected" to "assigned."</i>
Serologic crossmatch	Performance of a laboratory test to predict post-transfusion survival of transfused blood.
Serum	The liquid portion of blood used to detect antibodies.
Server	A computer system in a network shared by multiple users.
SF 518 equivalent	See "Blood Transfusion Reaction Form."
Site (see also "Facility")	A blood bank laboratory where a blood unit is located.
SME	Subject matter expert.
Software	Computer programs and associated documentation pertaining to the operation of a system.
Special Instructions (SIs)	Free-text or canned comments entered by the Blood Bank Technologist to communicate problems, restrictions, and considerations not otherwise documented in a patient's record.
Status	The disposition of a unit of blood.
Sterile connection device	An instrument that connects two plastic tubes while maintaining sterility.
SSRS	SQL Server Reporting Services.

Acronym, Term	Definition
Target unit	A unit formed as a result of a pool or split.
TAS	Type & Screen.
TCP/IP	Transmission Control Protocol/Internet Protocol.
Test (in general and laboratory contexts)	A means of examination to determine a property. In the context of the VBECS database, the contents of a single grid entry (test cell).
Test cell	In the context of the VBECS database: A grid space. In a laboratory: Red cells, in a reagent vial, used in a test.
Testing phase	1. The portion of a laboratory test in which an interpretation is recorded. 2. The phase of the development process in which the functionality is available for testing.
TIU	Text Integration Utility.
Transfusion Requirements (TRs)	Mandatory blood product attributes applied to every transfusion.
Truth table	A table, specific to a test or action, listing all possible combinations of user inputs and the corresponding acceptable results and/or interpretations. It allows the system to verify user input based on established parameters.
TRW	Transfusion reaction workup.
UBS	United Blood Services.
UID	Universal ID. An eye-readable and barcode-scannable coded representation of the Accession Number assigned by the laboratory package. <i>UID is unique to the generating information system where the specimen was received.</i>
Unit	A single example of a blood component type. Example: a single unit of red cells, a unit of plasma, a unit of platelets.
Unit ID	The label identifier, consisting of a string of characters, of the unit of blood, barcoded on the blood bag label and scanned when a blood unit is received in the blood bank in a shipment.

Acronym, Term	Definition
Unit status	<p>The state of a unit of blood at any given time.</p> <ul style="list-style-type: none"> • Available: Required processing is complete. • Crossmatched/Assigned: The unit is immediately available for transfusion. All required compatibility testing, including crossmatch if indicated, is performed. • Discarded: The unit has been destroyed and is no longer available. • Issued: The unit is not currently located in the blood bank. • Limited: Required processing is complete, except for ABO/Rh confirmation. • Quarantine: The unit cannot be selected. • Selected: The unit is reserved for one or more patients. Additional compatibility testing is required, such as antigen phenotype, crossmatch, and unit ABO/Rh confirmation. • Transferred: The unit has been returned to the original supplier or transferred to a different facility. • Transfused: The unit has been infused into a patient and is no longer available.
Validation	Establishment by objective evidence that a specific process or device specification conforms to user needs and intended use.
VASD	VA Service Desk (formerly National Help Desk).
VBECS specimen number	A blood bank-specific identifier for specimens.
VDL	VistA Documentation Library.
Verification	Confirmation by examination and provision of objective evidence that requirements were fulfilled.
Verify Code	A field in the VistA New Person file used to verify the identity of a user associated with an Access Code.
Vial identifier	A one- or two-character number that identifies a specific vial contained within a kit of reagents. <i>A vial contains antigen positive or antigen negative blood cells to be used in quality control of antisera. Also referred to as "cell identifier."</i>
VISN	Veterans Integrated Services Network.
VistA	Veterans Health Information Systems and Technology Architecture.
VistALink	A real-time communication link that provides connectivity and utilities for data exchange between VistA and VBECS.
WAN	Wide Area Network.
Worklist	A list detailing work to be performed.
Workload	Statistical data used to calculate productivity.
XM	Crossmatch.
XML	Extensible Markup Language.

Appendices

Appendix A: Downtime Forms and Instructions

Each site establishes its own policies, procedures, and processes to address unexpected VBECS downtime issues such as:

- The interface is unavailable. It is recommended that blood bank staff notify clinicians when messaging between CPRS and VBECS fails during order completion: clinicians may need to know that blood units are ready for pickup and transfusion.
- Test order handling and completion when the serologic testing is performed off-site or in a non-VA facility. The CPRS test order must be processed and completed even when testing is done off-site.
- The physical handling of a blood unit when it cannot be logged into inventory because VistALink is down and the patient is not in the VBECS database. The unit does not exist in VBECS until it is logged in, but is physically available.
- Manual blood labeling, verification of labels, and tracking documentation of the downtime labeling.
- Some Special Instructions (SIs) that were moved to VBECS during the database conversion must be entered as Transfusion Requirements (TRs) in VBECS. Inactivate the old SI and create a TR and new updated SI: the database conversion software cannot automatically translate the SI from database conversion into a rule-based TR.
- Creating the Audit Trail, Exception, and Testing Worklist Reports.
- Mitigation of the risk of patient identification conflicts due to patient merge issues (changing patient name or ID). It is recommended that staff submit the patient name and ID in writing to the blood bank when picking up blood products.
- VBECS is unavailable. To ensure accurate entry of information when VBECS is running, users must record the unit number and ID on their downtime form.

It is recommended that each site have available:

- Antibody Workup Forms
- Antigen Typing Forms (for patient)
- Antigen Typing Forms (for unit)
- Blood Transfusion Record Forms
- Caution Tags
- Daily QC Forms
- Emergency Issue Forms: developed by each site to reflect site-specific emergency issue policies and procedures
- Reagent Receipt Log
- Supply Receipt Log
- Transfusion Reaction Forms
- Transfusion Requirements Report: current printed version must be available at all times for users to check patient historic records and Transfusion Requirements. (See Appendix J: Recommended Report Usage: Table 39: Recommended Report Usage.)

Each site may use and modify (included in this appendix):

- Patient Testing Forms
- Unit ABO/Rh Confirmation Forms
- Unit Issue and Inspection Logs

- Unit Modification Forms

The information on these forms matches information to be input in VBECS.

Patient Testing Form

Facility Name:						
Address:						
City, State, Zip Code:						
Patient Name	Patient Identification Number	Location	CPRS Order Number	Lab Order Number	Urgency	Appropriate? (Yes, No)
Specimen UID	Date and Time Collected	Phlebotomist Name	Date and Time Received	Processed By:	Acceptable? (Yes, No)	
Ordered Tests:						
Comments:				Entered By:	Date and Time:	

Reagent Rack		Date and Time	Tech Name					
Historic ABO/Rh Check	Forward ABO Grouping			Rh Typing		Reverse ABO Grouping		ABO/Rh Interpretation
	Anti-A	Anti-B	Anti-A,B	Anti-D	Rh Control	A1 Cell	B Cell	

Antibody Screening																	
Screen Cell 1				Screen Cell 2				Screen Cell 3				Screen Cell 4				Screen Interpretation	
RT/IS ⁴	37°	AHG ⁵	CC ⁶	RT/IS	37°	AHG	CC	RT/IS	37°	AHG	CC	RT/IS	37°	AHG	CC		
Comments:																	

Patient Name		Patient Identification Number				Location		Specimen UID			
Compatibility Testing											
Unit Identification Number		Product Code	ABO/Rh	RT/IS	37°	AHG	CC	Interpretation		Tech Name	

⁴ Room Temperature/Immediate Spin

⁵ Anti-Human Globulin

⁶ Coated cell

Direct Antiglobulin Test					
Reagent Antiserum	Lot Number	AHG	CC	Interpretation	Tech Name
Polyspecific					
Anti-IgG					
Anti-C3d					

Direct Antiglobulin QC									
Reagent Antiserum	Lot Number	Expiration Date	Reagent Red Cells	Lot Number	Expiration Date	AHG	CC	Interpretation	Tech Name
Polyspecific									
Anti-IgG									
Anti-C3d									

Test Results: Positive, Graded = M, F, W, 1+, 2+, 3+, 4+; Negative = 0; NT = Not Tested

Unit ABO/Rh Confirmation Form

Facility Name:				
Address:				
City, State, Zip Code:				
Invoice Number	Source (Shipper) Name	Reagent Rack	Date and Time	Tech Name

Unit Identification Number	Product Code	Unit ABO/Rh	Anti-A	Anti-B	Anti-A,B	Anti-D	Rh Control	ABO/Rh Interpretation

Unit Issue and Inspection Log

Unit Modification Form

Facility Name:				
Address:				
City, State, Zip Code:				
Modification Date and Time	Modification Type	Method (Open, Closed, SCD)	Modified By:	
Unit Assigned? (Yes, No)	Patient Name	Patient Identification Number	Specimen UID	

Unit(s) to Be Modified				
Unit Identification Number	Product Code	ABO/Rh	Container Lot Number	Expiration Date and Time

Unit(s) to Be Modified				
Unit Identification Number	Product Code	ABO/Rh	Container Lot Number	Expiration Date and Time

Unit Assigned? (Yes, No)	Patient Name	Patient Identification Number	Specimen UID

Modified Unit(s)				
Unit Identification Number	Product Code	ABO/Rh	Container Lot Number	Expiration Date and Time

This page intentionally left blank.

Appendix B: Database Table Information

Antibody and Antigen Tables

Antibody and Antigen Table: Irregular Antibodies

VBECS creates a division-specific copy in each division in the database. Setting Antibody Parameters allows the division to reset the calculated antigen negative requirement of antibodies to its medical director's policy. Each orderable test name has a corresponding repeat test.

Clinical significance of antibody specificity is set using the 14th edition of the AABB Technical Manual. A specific instance of an antibody for certain patients or different treatment conditions may behave quite differently, requiring antigen negative or crossmatch-compatible blood for the patient. The user may edit these settings for an individual patient.

The antigen negative required settings apply to all instances of the antibody. Individual hospital policy may indicate that a setting of "OPT" should be made "REQ." A policy of a global OPT may remain in the event of a specific clinically significant instance of a specificity in a particular patient. The individual patient's antigen negative requirement may be set for the specificity in Entering and Removing Transfusion Requirements. In the Antigen Negative Requirement column:

- "REQ" sets an antigen negative requirement that is not editable because these are generally considered to require the provision of antigen negative blood for patient transfusion.
- "NO" does not set an antigen negative requirement and is not editable.
- "OPT" in the Antibody sets an antigen negative requirement is editable and may be reset repeatedly on the division's copy of the table.

Table 13: Antibody and Antigen Table: Irregular Antibodies

Antibody Name		Clinically Significant?§	Does Antibody Set an Antigen Negative Requirement?	Associated Antigen		Antigen Negative Compatibility Percentage	VBECS Antigen Typing Orderable Test
VBECS*	VistA†			VBECS	VistA†		
Anti-D	Anti-D	Yes	Yes	D	D	13	AGD
N/A	N/A	N/A	N/A	N/A	N/A	N/A	AGwD
Anti-C	Anti-C	Yes	Yes	C	C	53	AGC
Anti-E	Anti-E	Yes	Yes	E	E	75	AGE
Anti-c	Anti-c	Yes	Yes	c	c	12	AGc
Anti-e	Anti-e	Yes	Yes	e	e	2	AGe
Anti-f	Anti-f	Yes	Yes	c and e	f	N/A	AGc AGe
Anti-G	Anti-G	Yes	Yes	C and D	G	N/A	AGD AGC
Anti-K	Anti-K	Yes	Yes	K	K	91	AGK
Anti-k	Anti-k	Yes	Yes	k	k	1	AGk
Anti-Kp(a)	Anti-Kp(a)	Yes	Optional	Kp(a)	Kp(a)	98	AGKp(a)
Anti-Kp(b)	Anti-Kp(b)	Yes	Optional	Kp(b)	Kp(b)	2	AGKp(b)
Anti-Jk(a)	Anti-Jk(a)	Yes	Yes	Jk(a)	Jk(a)	15	AGJk(a)
Anti-Jk(b)	Anti-Jk(b)	Yes	Yes	Jk(b)	Jk(b)	42	AGJk(b)
Anti-Jk3	Anti-Jk(a)Jk(b)	Yes	Yes	Jk(a) and Jk(b)	None	N/A	AGJk(a) AGJk(b)
Anti-Fy(a)	Anti-Fy(a)	Yes	Yes	Fy(a)	Fy(a)	63	AGFy(a)
Anti-Fy(b)	Anti-Fy(b)	Yes	Yes	Fy(b)	Fy(b)	47	AGFy(b)
Anti-M	Anti-M	No	Optional	M	M	26	AGM
Anti-N	Anti-N	No	Optional	N	N	27	AGN
Anti-S	Anti-S	Yes	Yes	S	S	58	AGS
Anti-s	Anti-s	Yes	Yes	s	s	7	AGs
Anti-U	Anti-U	Yes	Optional	U	U	0	AGU
Anti-P1	Anti-P1	No	Optional	P1	P1	21	AGP1
Anti-Le(a)	Anti-Le(a)	No	Optional	Le(a)	Le(a)	78	AGLe(a)
Anti-Le(b)	Anti-Le(b)	No	Optional	Le(b)	Le(b)	37	AGLe(b)
Anti-Js(a)	Anti-Js(a)	Yes	Optional	Js(a)	Js(a)	90	AGJs(a)
Anti-Js(b)	Anti-Js(b)	Yes	Optional	Js(b)	Js(b)	10	AGJs(b)
Anti-Lu(a)	Anti-Lu(a)	No	Optional	Lu(a)	Lu(a)	92	AGLu(a)
Anti-Lu(b)	Anti-Lu(b)	No	Optional	Lu(b)	Lu(b)	2	AGLu(b)
Anti-Lu3	Anti-Lu(a)Lu(b)	No	Optional	Lu (a) and Lu(b)	None	N/A	AGLu(a) AGLu(b)
Anti-C(w)	Anti-Cw	Yes	Optional	C(w)	Cw	99	AGC(w)
Anti-Ce	Anti-Ce	Yes	Yes	C and e	Ce	N/A	AGC AGe
Anti-C(x)	Anti-Cx	Yes	Optional	C(x)	Cx	99	AGC(x)
Anti-E(w)	Anti-E(w)	Yes	Optional	E(w)	E(w)	99	AGE(w)
Anti-V	Anti-V	Yes	Optional	V	V	85	AGV
Anti-Xg(a)	Anti-Xg(a)	No	Optional	Xg(a)	Xg(a)	25	AGXg(a)
Anti-M(g)	Anti-M(g)	No	Optional	M(g)	M(g)	99	AGM(g)
Anti-Mi(a)	Anti-Mi(a)	No	Optional	Mi(a)	Mi(a)	99	AGMi(a)
Anti-P	Anti-P	No	Optional	P	P (not otherwise specified)	25	AGP

VBECS*	VistA†	Clinically Significant?§	Does Antibody Set an Antigen Negative Requirement?	Associated Antigen		Antigen Negative Compatibility Percentage	VBECS Antigen Typing Orderable Test
				VBECS	VistA†		
Anti-A1	Anti-A-1	No	Optional	A1	A-1	10	AGA1
Anti-H	Anti-H	No	Optional	H	H	0	AGH
Anti-I	Anti-I	No	No	None	None	0	None
Anti-i	Anti-i	No	No	None	None	99	
Warm autoantibody	Warm autoantibody	Unknown§	No	None	None	0	
Cold autoantibody	Cold autoantibody	No	No	None	None	0	
Anti-Le(ab)	Anti-Le(x)	No	Optional	Le(a) and Le(b)	None	N/A	AGLe(a) AGLe(b)
Anti-Di(a)	Anti-Di(a)	No	Optional	Di(a)	Di(a)	0	AGDi(a)
Anti-Di(b)	Anti-Di(b)	No	Optional	Di(b)	Di(b)	99	AGDi(b)
Anti-Do(a)	Anti-Do(a)	No	Optional	Do(a)	Do(a)	33	AGDo(a)
Anti-Do(b)	None	No	Optional	Do(b)	None	17	AGDo(b)
Anti-Co(a)	None	No	Optional	Co(a)	None	33	AGCo(a)
Anti-Co(b)	None	No	Optional	Co(b)	None	17	AGCo(b)
Anti-LW(a)	None	No	Optional	LW(a)	None	0	AG LW(a)
Anti-LW(b)	None	No	Optional	LW(b)	None	99	AG LW(b)
Anti-Sc1	None	No	Optional	Sc1	None	0	AGSc1
Anti-Sc2	None	No	Optional	Sc2	None	99	AGSc2
Anti-Sd(a)	None	No	Optional	Sd(a)	None	1	AGSd(a)
Anti-Au(a)	Anti-Au(a)	No	Optional	Au(a)	Au(a)	20	AGAu(a)
Anti-Au(b)	None	No	Optional	Au(b)	None	50	AGAu(b)
Anti-Ch	None	No	Optional	Ch	None	2	AGCh
Anti-Rg	None	No	Optional	Rg	None	3	AGRg
Anti-Cs(a)	None	No	Optional	Cs(a)	None	4	AGCs(a)
Anti-Yk(a)	None	No	Optional	Yk(a)	None	4	AGYk(a)
Anti-Kn(a)	None	No	Optional	Kn(a)	None	1	AGKn(a)
Anti-McC(a)	None	No	Optional	McC(a)	None	1	AGMcC(a)
Anti-JMH	None	No	Optional	JMH	None	1	AGJMH
Anti-Yt(a)	None	No	Optional	Yt(a)	None	1	AGYt(a)
Anti-Yt(b)	None	No	Optional	Yt(b)	None	92	AGYt(b)
HTLA (probable)	None	No	No	None	None	0	None
Antibody, No Specificity Identified	Antibody (not otherwise specified)	Unknown§	No	None	None	0	
Antibody to Low-Incidence Antigen	Antibody, private (not otherwise specified)	Unknown§	No	None	None	99	

Antibody Name	VistA†	Clinically Significant?§	Does Antibody Set an Antigen Negative Requirement?	Associated Antigen		Antigen Negative Compatibility Percentage	VBECS Antigen Typing Orderable Test
				VBECS	VistA†		
VBECS*	VistA†	Unknown§	No	None	None	0	
Antibody to High-Incidence Antigen	Antibody, public (not otherwise specified)						
Serologic Crossmatch Test Required	N/A	No	No	None	None	N/A	None

*Specificities written as superscript or subscript letters in blood bank literature are shown as letters enclosed in parentheses in VBECS (the GUI cannot display superscript or subscript letters). Specificities written as subscript characters in blood bank literature are shown as numbers or letters not enclosed in parentheses.

†Examples of common VistA antibody and antigen terminology. Individual sites may vary. Refer to division-specific mapping documents.

‡The user will order the tests for each antigen separately; VBECS does not display a duplicate antigen typing for this combination antibody specificity.

§When the clinical significance is unknown, it is displayed as "No" due to the limitations of the system.

Antibody and Antigen Table: Database Conversion

Users cannot set antigen negative requirements for these antibodies in VBECS. Clinical significance of antibody specificity is set using the 14th edition of the AABB Technical Manual.

Table 14: Antibody and Antigen Table: Database Conversion

Antibody Name		Associated Antigen	
VBECS*	VistA†	VBECS*	VistA†
Anti-P1+P+P(k)	Anti-PP1	P1+P+P(k)	PP1
Anti-Pk	Anti-p(k)	P(k)	p(k)
Anti -M, other	Anti-M, other	None	None
Anti-LW	Anti-LW	None	LW
Anti-Lu	Anti-Lu	None	Lu
Anti-Le, other	Anti-Le, other	None	None
Anti-N, other	Anti-N, other	None	None
Anti-I(int)	Anti-I(int)	None	I(int)
Anti-rhesus (not otherwise specified)	Anti-rhesus (not otherwise specified)	None	None

*Specificities written as superscript or subscript letters in blood bank literature are shown as letters enclosed in parentheses in VBECS (the GUI cannot display superscript or subscript letters). Specificities written as subscript characters in blood bank literature are shown as numbers or letters not enclosed in parentheses.

†Examples of common VistA antibody and antigen terminology. Individual sites may vary. Refer to division-specific mapping documents.

Antibody Screen Test Interpretation

For an Antibody Screen Test (ABS) to be correctly interpreted, the system must take the row validation (Table 15) for each screening cell and test method and determine whether the screening in total is positive or negative. This table is the input from the Blood Bank Technologist. Interpreting the ABS: The Blood Bank Technologist does not enter an interpretation for each row/cell of the antibody screen. The technologist enters an interpretation for the combination of screening cells including the autocontrol when tested. The system must correlate the information from the rows in the test, as described in system rules , to determine the proper interpretation for the ABS. Valid interpretations are positive or negative. A technologist may enter an antibody screen comment with the interpretation.

- SR 1: Valid interpretation. Update database.
- SR 2: Warning message. No override. Do not update database at this time.
- SR 3: Valid interpretation. Update database. Collect comment and exception.
- EM 1: "Interpretation does not match your results."
- EM 2: "You are attempting to file an invalid antibody screen. This is not allowed."

Table 15: Antibody Screen Test Interpretation

Antibody Screen Test Interpretation					
Antibody Screen Interpretation: includes all cells tested; might be SC 1 and SC 2 only, or SC 1, 2, 3, or SC 1, 2, 3, 4, as indicated by site parameter in Configure Daily QC. Method of Tube, Gel, or Solid Phase	SC 1, SC 2, SC 3, SC 4, Autocontrol Interpretation	Tech Entry	System Interpretation of Screening Cell Test	System Response	Error Message
System interpretation if all individual row validations are negative.	Negative	Negative (N)	SR 1	N/A	
	Positive	Negative (N)	SR 2	EM 1	
System interpretation if all individual row validations are negative and indicates that a comment is required.	Negative	Negative (N)	SR 3	N/A	
	Positive	Positive (P)	SR 1	N/A	
System interpretation if one, several, or all rows are positive (P).	Negative	Positive (P)	SR 2	EM 1	
	Positive	Positive (P)	SR 3	N/A	
System interpretation if one, several, or all rows are invalid (I).	Positive or Negative	Invalid (I)	SR 2	EM2	

Canned Comment Category Types and Text

Canned comments are defined comments that can be inserted in a record during execution of a variety of functions. The canned comments available in a function are restricted to the category type associated with that function. Table 16 lists the category types available when defining canned comments.

Table 16: Canned Comment Category Types and Text

Text	Category Types																																					
	ABO Incompatible Issue	Crossmatch	Incoming Shipment	Justify ABO/Rh Change	Justify ABO/Rh Discrepancy	Order Cancellation	Order Processing	Parameter Setting Change	Patient Testing	Quarantine/Discard	Reagent and Supply	Release from Quarantine	Reason for Release	Special Instruction Inactivation	Special Instructions	Specimen Extension	Specimen General	Specimen Processing	Specimen Unacceptable	Transfusion Processing	Transfusion Reaction Workup	Specimen Required	Pre-Transfusion Specimen Checks	Post-Transfusion Specimen Checks	Implicated Units	Symptom Details	Transfusion Blood Bag Not Returned	Transfusion Clerical Checks Failed	Transfusion Requirement Inactivation	Unit Inactivation	Unit Data Validation	Unit Issue	Unit Modification	Unit Return	Unit Restriction	Unit Selection	Unit Status Removal	Unit Testing
Remote storage cooler A																																						
Resetting due to policy change						X																																
Rh mismatch selection necessary																																						
Ruptured during thaw							X																													X		
Specimen and paperwork do not match																								X	X	X												
Specimen expired								X																														
Specimen hemolyzed																		X																				
Specimen misplaced																	X	X					X	X														
Specimen not received																	X	X																				
Specimen recollected																																						

CPRS Orderable Blood Components (Component Classes) Mapped to ICCBBA Component Classes

Table 17 is used to map the CPRS/VistA orderables available to clinicians to the ICCBBA component classes.

Table 17: CPRS Orderable Blood Components (Component Classes) Mapped to ICCBBA Component Classes

VBECS Component Class	Code	Short Name	Product Type Name (ICCBBA Component Class)
RED BLOOD CELLS	E002	RBC	RED BLOOD CELLS
	E003	RBC_Wash	Washed RED BLOOD CELLS
	E004	RBC_Froz	Frozen RED BLOOD CELLS
	E005	RBC_Froz_Rej	Frozen Rejuvenated RED BLOOD CELLS
	E006	RBC_Deg	Deglycerolized RED BLOOD CELLS
	E007	RBC_Deg_Rej	Deglycerolized Rejuvenated RED BLOOD CELLS
	E008	RBC_Rej	Rejuvenated RED BLOOD CELLS
	E009	RBC_AFR	Apheresis RED BLOOD CELLS
	E049	RBC_AFR_Wash	Washed Apheresis RED BLOOD CELLS
	E050	RBC_AFR_Froz	Frozen Apheresis RED BLOOD CELLS
	E051	RBC_AFR_Deg	Deglycerolized Apheresis RED BLOOD CELLS
	E052	RBC_AFR_Rej	Rejuvenated Apheresis RED BLOOD CELLS
	E053	RBC_AFR_Froz_Rej	Frozen Rejuvenated Apheresis RED BLOOD CELLS
	E054	RBC_AFR_Deg_Rej	Deglycerolized Rejuvenated Apheresis RED BLOOD CELLS
WHOLE BLOOD	E001	WB	WHOLE BLOOD
PLATELETS	E020	PLT	PLATELETS
	E021	PLT_Wash	Washed PLATELETS
	E022	PLT_Pool	POOLED PLATELETS
	E023	PLT_Pool_Wash	Washed POOLED PLATELETS
	E024	PLT_AFR	Apheresis PLATELETS
	E025	PLT_AFR_Froz	Frozen Apheresis PLATELETS
	E026	PLT_AFR_Thaw	Thawed Apheresis PLATELETS
	E027	PLT_AFR_Wash	Washed Apheresis PLATELETS
	E060	PLT_Pool_Froz	Frozen POOLED PLATELETS
	E061	Plt_Pool_Thaw	Thawed POOLED PLATELETS
	E064	PLT_Pool_Thaw_Wash	Washed Thawed POOLED PLATELETS
	E065	PLT_AFR_Thaw_Wash	Washed Thawed Apheresis PLATELETS
	E068	PLT_AFR_Pool_Wash	Washed Apheresis POOLED PLATELETS
	E069	PLT_Froz	Frozen PLATELETS
FRESH FROZEN PLASMA	E010	FFP	FRESH FROZEN PLASMA
	E011	FFP_Thaw	Thawed FRESH FROZEN PLASMA
	E012	FFP_AFR	Apheresis FRESH FROZEN PLASMA
	E013	FFP_AFR_Thaw	Thawed Apheresis FRESH FROZEN PLASMA
	E014	PLASMA_AFR	Apheresis PLASMA
	E015	PLASMA_AFR_Thaw	Thawed Apheresis PLASMA
	E016	PLASMA_Liq	Liquid PLASMA
	E017	PLASMA	PLASMA
	E018	PLASMA_Thaw	Thawed PLASMA
	E041	PLASMA_AFR_Liq	Liquid Apheresis PLASMA
	E055	PLASMA_Pool_Thaw	Thawed Pooled Plasma
	E056	FFP_Pool_Thaw	Thawed POOLED FRESH FROZEN PLASMA

VBECS Component Class	Code	Short Name	Product Type Name (ICCBBA Component Class)
	E058	PLASMA_Liq_Pool	Liquid POOLED PLASMA
	E059	FFP_Pool	POOLED FRESH FROZEN PLASMA
	E062	FFP_Afr_Pool_Thaw	Thawed Apheresis POOLED FRESH FROZEN PLASMA
	E063	Plasma_Afr_Pool_Thaw	Thawed Apheresis POOLED PLASMA
	E066	FFP_Lyfo	Lyophilized FRESH FROZEN PLASMA
	E067	FFP_Recon	Reconstituted FRESH FROZEN PLASMA
CRYOPRECIPITATE	E028	CRYO	CRYOPRECIPITATE
	E029	CRYO_Thaw	Thawed CRYOPRECIPITATE
	E030	CRYO_Pool	POOLED CRYOPRECIPITATE
	E031	CRYO_Pool_Thaw	Thawed POOLED CRYOPRECIPITATE
	E032	CRYO_AFR	Apheresis CRYOPRECIPITATE
	E033	CRYO_AFR_Thaw	Thawed Apheresis CRYOPRECIPITATE
OTHER	E019	PLASMA_Pltrch	PLATELET RICH PLASMA
	E034	GRAN	GRANULOCYTES
	E035	GRAN_AFR	Apheresis GRANULOCYTES
	E036	GRAN_Pool	POOLED GRANULOCYTES
	E037	GRAN-PLT_AFR	Apheresis GRANULOCYTES-PLATELETS
	E038	LEUK	LEUKOCYTES
	E039	LEUK_AFR	Apheresis LEUKOCYTES
	E040	PLASMA_Pool	POOLED PLASMA
	E042	BUFFY_Pltrch	PLATELET RICH BUFFY COAT
	E043	BUFFY_Pltrch_Pool	POOLED PLATELET RICH BUFFY COAT
	E044	LYMPH_AFR	Apheresis LYMPHOCYTES
	E045*	MONO_AFR	Apheresis MONOCYTES
	E046	SERUM	SERUM
	E047*	SERUM_Pool	POOLED SERUM
	E048	SERUM_Pool_Froz	Frozen POOLED SERUM
	E057	GRAN_Wash	Washed GRANULOCYTES

*This product type currently has no ICCBBA assigned product codes.

Details in Audit Trail Report

Changes include the date and time of the change, ID of the user who made the change, and required comments.

Table 18: Details in Audit Trail Report

Report Section	Generating Options	GUI Fields Monitored	Accessible in:	
			Unit History Report?	Patient History Report?
Blood Unit Changes	Modify Units: Pool Units (when a previously pooled unit is edited)	For each pooled unit edited (unit ID, long name, product code): <ul style="list-style-type: none"> • Added or Removed Unit IDs and Product Codes • Comment 	Yes	No
	Editing A Unit	For each unit edited (unit ID, long name, product code): <ul style="list-style-type: none"> • Unit Antigen Information • Special Testing Information • Discarded Plasma Volume • Unit Volume 	Yes	No

Report Section	Generating Options	GUI Fields Monitored	Accessible in:	
			Unit History Report?	Patient History Report?
		<ul style="list-style-type: none"> • Biohazardous Indicator • Unit ABO/Rh Confirmation Test, Invalidation Indicator • Unit Antigen Typing Test, Invalidation Indicator • Unit Record Inactivation Indicator • Restricted For Patient Name and Patient ID • Comment 		
Blood Unit Financial Changes	Edit Financial Data	For each unit edited (unit ID, long name, product code): <ul style="list-style-type: none"> • Unit Base Cost • Special Test Cost • Return Credit 	Yes	No
Patient Testing Changes	Invalidate Test Results	For each patient (name, ID): <ul style="list-style-type: none"> • Ordered Test Original Results, Invalidation Indicator • Comment 	No	Yes
	Justify ABO/Rh Change	For each patient (name, ID): <ul style="list-style-type: none"> • Patient ABO/Rh • Comment 	No	Yes
Patient Transfusion Requirements	Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement Patient Testing: Record Patient Test Results (new TR created from ABID data entry)	For each Transfusion Requirement category or Special Instruction: <ul style="list-style-type: none"> • Status Indicator • Comment 	No	Yes
Maintenance: Local Facilities	Local Facilities	For each facility name and FDA Registration Number: <ul style="list-style-type: none"> • Eye-readable prefix • ICCBBA Registration Number • Alpha Characters Indicator • Active Facility Indicator • Collection Facility Indicator • Testing Facility Indicator • Address Line 1 • Address Line2 • Address Line 3 • City • State • Zip • Phone • Fax 	No	No
Maintenance: Blood Product	Blood Products	For each blood product (name, code): <ul style="list-style-type: none"> • HCPCS Code and Text For each shipper of this product code (name, ID):	No	No

Report Section	Generating Options	GUI Fields Monitored	Accessible in:	
			Unit History Report?	Patient History Report?
		<ul style="list-style-type: none"> • Cost • Return Credit Percentage • Status Indicator 		
Maintenance: Division Configuration	Configure Division	<p>For each Division (Name and Code):</p> <ul style="list-style-type: none"> • Full Service, Transfusion Only Indicator • Facility Name • Active Status • Full-Face-Label Printer Status • Full-Face-Label Printer Port Number • Full-Face-Label Printer IP address • Full-Face-Label Printer COM • Report Printer Name • Accession Area Name • Lock Inactivity Timeout Setting • Time Zone Setting • Daylight savings time Setting • Daylight savings time start date • Daylight savings time end date <p>For each Associated Institution (Division Code)</p> <ul style="list-style-type: none"> • Status indicator 	No	No
Maintenance: Division Configuration	Configure Division	<p>For each division (name, code):</p> <ul style="list-style-type: none"> • Electronic Crossmatch Enabled Indicator • ICCBBA Registration Number • Maximum Specimen Expiration Days • Blood Bank MD 	No	No
Maintenance: Division Configuration	Edit Invoice Text	<p>For each division (name, code):</p> <ul style="list-style-type: none"> • Invoice Text 	No	No
Maintenance: Division Configuration	Configure Testing	<p>For each division (name, code):</p> <ul style="list-style-type: none"> • Antibody Screen Indicator • ABS and XM Testing Phases 	No	No
Maintenance: Division Configuration	Order Alerts	<p>For each division (name, code):</p> <ul style="list-style-type: none"> • Login Message 	No	No
Maintenance: Modifications Configuration	Product Modifications	<p>For each modification type:</p> <ul style="list-style-type: none"> • Enabled or Disabled Indicator • Cost 	No	No
Maintenance: User Role Changes	Update User Roles	<p>For each user:</p> <ul style="list-style-type: none"> • User Role Security Setting 	No	No
Maintenance: User Role Changes	Configure Users	<p>For each NT user (name, ID):</p> <ul style="list-style-type: none"> • VistA DUZ • VistA User Name • Email Address • User Initials • Active VBECS User Indicator 	No	No

Report Section	Generating Options	GUI Fields Monitored	Accessible in:	
			Unit History Report?	Patient History Report?
		<ul style="list-style-type: none"> • Active Division Indicator 		
Maintenance: Antibody Configuration	Antibodies	For each antibody specificity: <ul style="list-style-type: none"> • Compatibility Percentage • Higher Level Override Required Indicator 	No	No
Maintenance: Workload Configuration	Workload Codes	For each workload process and for each LMIP procedure and LMIP/NLT code: <ul style="list-style-type: none"> • Start Date • Stop Date • Status Indicator 	No	No
Maintenance: MSBOS Configuration	MSBOS	For each surgery name: <ul style="list-style-type: none"> • Surgery Status Indicator For each MSBOS recommendation: <ul style="list-style-type: none"> • TAS Indicator • No Blood Required Indicator For each selected component class: <ul style="list-style-type: none"> • Number of Units • Recommendation Status Indicator 	No	No
Maintenance: Transfusion Effectiveness Changes	Transfusion Effectiveness	For each Laboratory test name, for each VistA Laboratory test Internal Entry Number (IEN), and for each specimen type: <ul style="list-style-type: none"> • Status Indicator 	No	No
Maintenance: Transfusion Complications Changes	Transfusion Complications	For each Laboratory test name, for each VistA Laboratory test IEN, and for each specimen type: <ul style="list-style-type: none"> • Threshold Result • Threshold Setting Status • Status indicator 	No	No
Maintenance: Component Class	Component Classes	For each component class configured: <ul style="list-style-type: none"> • Current ABO/Rh Specimen Results Required Indicator • Maximum Transfusion Time • For each Laboratory Test Name, for each VistA Laboratory Test Internal Entry Number (IEN), and for each Specimen Type: Lab Test Status • Threshold Result 	No	No
Maintenance: Daily QC Configuration	Configure Daily QC	For each daily QC setup changed: <ul style="list-style-type: none"> • Commercial or Non-Commercial QC Indicator • Rack Name Indicator (A through Z or 1 through 26) • Number of Racks to Test • Daily Alert Time • Primary Enhancement Media • Secondary Enhancement Media • Comment For each reagent type and test with combination: <ul style="list-style-type: none"> • Status Indicator 	No	No

Report Section	Generating Options	GUI Fields Monitored	Accessible in:	
			Unit History Report?	Patient History Report?
		<ul style="list-style-type: none"> • Minimum Reaction Entry 		
Maintenance: Minimum Reagent Inventory Levels	Maintain Minimum Levels	For each reagent type: <ul style="list-style-type: none"> • Minimum Inventory Level 	No	No
Maintenance: User Alert Configuration	Order Alerts	For each division (name, code): <ul style="list-style-type: none"> • Printer Alerts On • Printer Name • (Printer) Order Alert Type • (Printer) Orders Needed In • Icon Alerts On • Refresh Rate • (Icon) Order Alert Type • (Icon) Orders Needed In 	No	No
Maintenance: Interface Control	Configure Interfaces	<ul style="list-style-type: none"> • VistA TCP/IP Address • VistA TCP Port Number • VistA Domain 	No	No
Maintenance: Interface Control	Configure Interfaces	For each interface defined: <ul style="list-style-type: none"> • Interface IP Address • Interface Domain • Interface Port Number • Interface Facility ID • VBECS IP Address • VBECS Port Number • VBECS Facility ID • Completed Messages • Messages in Error • Log Events and HL7 Messages to Event Log Indicator • Ack Timeout • Retransmit Attempts • Interface Administrator Email • Interface Status 	No	No

Details in Exception Report

Table 19 lists exception types and the reports in which they appear.

Table 19: Details in Exception Report

Exception Type	Available to Role	Available in Option	Accessible in:	
			Unit History Report?	Patient History Report?
ABO discrepant unit issue	All	Issue Unit	No	Yes
ABO incompatible unit selected	None	Select Unit	No	Yes
ABO Incompatible Unit Transfused	\geq Traditional Supervisor	Document ABO Incompatible Transfusion	No	Yes
Antigen positive unit crossmatched	\geq Lead Tech	Select Unit (when crossmatched)	No	Yes
Antigen positive unit selected	> Lead Tech	Select Unit	No	Yes
Antigen testing phase change	All	Unit or Patient Antigen typing	Yes (Antigen Typing)	Yes (Patient Testing: Record a Patient Antigen Typing)
Antigen-positive/untested units issued	\geq Lead Tech	Issue Blood	No	Yes
Biohazardous unit issued	All	Issue Blood	No	Yes
Biohazardous unit modified	All	Modification	Yes	No
Biohazardous unit selected	All	Select Unit	No	Yes
Caution Tag or Blood Transfusion Record Form printed for antigen positive or untested unit	\geq Lead Tech		No	Yes
Caution Tag or Blood Transfusion Record Form printed on expired specimen	\geq Enhanced Tech		No	Yes
Caution Tag or Blood Transfusion Record Form printed on expired unit	\geq Enhanced Tech		No	Yes
Caution Tag or Blood Transfusion Record Form printed with unsatisfied Transfusion Requirement	\geq Lead Tech		No	Yes
Crossmatch incompatible: Give only with MD approval	All		No	Yes
Deletion of partially completed QC	All	Enter Daily QC	No	No
Discrepant ABO/Rh override	\geq Enhanced Tech	Patient Testing	No	Yes
Discrepant patient antigen typing	All	Patient Testing	No	Yes
Electronic Transfusion Record Update	BCE COTS system	Post -Transfusion Information	No	No
Emergency issue request (no specimen)	All	Accept Order	No	Yes
Emergency issue, require testing incomplete	All	Patient Testing	No	Yes

Exception Type	Available to Role	Available in Option	Accessible in:	
			Unit History Report?	Patient History Report?
Expired antisera used	All	Enter Rack QC, Patient and Unit Antigen typing	Yes (Antigen Typing)	Yes (Patient Testing: Record a Patient Antigen Typing)
Expired reagent QC'd	All	Enter Daily QC	No	No
Expired supply used in modification	All	Modify Units	No	No
Expired task processed	All	Pending Task List	No	Yes
Expired unit modified	> Enhanced Tech	Modify Units	Yes	No
Expired unit received	All	Incoming Shipment	Yes	No
Expired unit released from quarantine	> Enhanced Tech	Release from Quarantine	Yes	No
Expired unit selected	> Enhanced Tech	Select Unit	No	Yes
Final unit status or modification removed	> Traditional Supervisor	Remove Final Status	Yes	No
Inconclusive crossmatch	All	Patient Testing, serologic crossmatch grid	No	Yes
Issued expired unit	> Enhanced Tech	Issue Blood	No	Yes
Less restrictive unit issued	All	Issue Blood	No	Yes
Less restrictive unit selected	All	Select Unit	No	Yes
Non-standard test procedure	All	Patient Testing, Antibody screen, serologic crossmatch	No	Yes
(Patient) ABO/Rh discrepancy	All	Patient testing, ABO/Rh test	No	Yes
Previously recorded results inactivated	All	All testing grids where the red x is used to clear a grid.	No	Yes
Previously recorded results invalidated	All	Invalidate Patient Testing	Yes	No
QC decrease reagent reactivity ≥ 2	All	All testing where Rack is selected	No	No
QC not performed on rack used for testing	All	All testing where Rack is selected	No	Yes
Reflex test ordered on expired specimen	All	Order Reflex Test	No	Yes
Released patient restriction	> Enhanced Tech	Free Directed Donor	No	Yes
Retrospective update prior to last record update	> Lead Tech	Discard/Quarantine	Yes	No
Rh positive units selected	All	Select Units	No	Yes
Specimen expiration date extended	All	Maintain Specimen	No	Yes
Target product label failed checks	All	All modification processes	Yes	No
Target product outdate extended	> Enhanced Tech	All modification processes	Yes	No
Transfusion Reaction Workup	All	Transfusion Reaction Workup	No	Yes
Transfusion Requirement incompatible unit selected	> Lead Tech	Select Unit	No	Yes

Exception Type	Available to Role	Available in Option	Accessible in:	
			Unit History Report?	Patient History Report?
Unacceptable/expired specimen used	All	Patient Testing	No	Yes
Unit ABO/Rh confirmation inconclusive, unit quarantined	All	ABO/Rh Confirmation	Yes	No
Unit ABO/Rh log-in vs. confirmation do not match, unit quarantined	All	ABO/Rh Confirmation	Yes	No
Unit discarded	All	Discard	Yes	No
Unit ER issued, testing problem	All	Patient testing, unit antigen typing	No	Yes
Unit issued on expired specimen	> Enhanced Tech	Issue Unit	No	Yes
Unit issued with unsatisfied Transfusion Requirement	> Lead Tech	Issue Unit	No	Yes
Unit quarantined	All	Quarantine	Yes	No
Unit released from quarantine	> Enhanced Tech	Release Unit from Quarantine	Yes	No
Unit unsatisfactory upon return from issue	All	Return Issued Units to Blood Bank	Yes (quarantined)	No
			Yes (quarantined with patient assignment)	Yes

Enable Crossmatch Option for OTHER Product Types

During the execution of Select Units, the user has the option to request a crossmatch test for the OTHER component class. Table 20 indicates which of the OTHER product types must have the crossmatch test enabled.

Table 20: Enable Crossmatch Option for OTHER Product Types

Component Class	ICCBBA Component Class Allowed	Allow XM Order
OTHER	PLASMA	No
OTHER	APHERESIS PLASMA	No
OTHER	PLATELET-RICH PLASMA	No
OTHER	GRANULOCYTES	Yes
OTHER	APHERESIS GRANULOCYTES-PLATELETS	Yes
OTHER	LEUKOCYTES	Yes
OTHER	POOLED PLASMA	No
OTHER	PLATELET-RICH BUFFY COAT	Yes
OTHER	APHERESIS LYMPHOCYTES	Yes
OTHER	APHERESIS MONOCYTES	Yes

Allowable Product Modifications by Original Product Type

ICCBBA provides a comprehensive database with every known blood product type that includes a five-digit product code, Product Type classification, and the Product Attribute for each blood product. The VBECS Blood Product Table combines this ICCBBA ISBT 128 blood product database and the known Codabar blood products. Each Codabar blood product was assigned a Product Type and the Product Attribute Form column populated with a core condition and some basic attributes.

Calculation of the Expiration Date

VBECS calculates the expiration date of the new unit based on system rules specific to the blood product and modification types, and on whether the system is open or closed. The user may accept or edit this expiration date.

Table 21 shows which kind of modifications are allowed for a specific product type. The first column contains product types; the second set of columns contains allowable modifications. Blank cells indicate that a modification type cannot be performed on an original product type. Modifications are not allowed for product types not listed.

Table 21: Allowable Product Modifications by Original Product Type

Allowable Product Modifications by Original Product Type												
This table defines the allowable Modification using the Original/Parent Unit's ICCBBA Product Type*												
Product Type Name	ICCBBA Product Type*	Split/Divide	Thaw Pool	Thaw/Pool	Irradiate	Leuko-reduce	Volume Reduce	Wash	Rejuve-nate	Freeze	Deglyce-rolize	
WHOLE BLOOD	E001	X				X	X	X	X	X	X	
RED BLOOD CELLS	E002	X				X	X	X	X	X	X	
Washed RED BLOOD CELLS	E003	X				X	X			X	X	
Frozen RED BLOOD CELLS	E004											X
Frozen Rejuvenated RED BLOOD CELLS	E005											X
Deglycerolized RED BLOOD CELLS	E006	X				X	X					
Deglycerolized Rejuvenated RED BLOOD CELLS	E007	X				X	X					
Rejuvenated RED BLOOD CELLS	E008	X				X	X				X	
APHERESIS RED BLOOD CELLS	E009	X				X	X		X	X	X	

Allowable Product Modifications by Original Product Type												
This table defines the allowable Modification using the Original/Parent Unit's ICCBBA Product Type*												
Product Type Name	ICCBBA Product Type*	Split/Divide	Thaw	Pool	Thaw/Pool	Irradiate	Leuko-reduce	Volume Reduce	Wash	Rejuvenate	Freeze	Deglycerolize
FRESH FROZEN PLASMA	E010		X									
Thawed FRESH FROZEN PLASMA	E011	X				X	X					
APHERESIS FRESH FROZEN PLASMA	E012		X		X							
Thawed APHERESIS FRESH FROZEN PLASMA	E013	X		X		X	X					
APHERESIS PLASMA	E014		X		X							
Thawed APHERESIS PLASMA	E015	X		X		X	X					
Liquid PLASMA	E016	X				X	X					
PLASMA	E017		X									
Thawed PLASMA	E018	X				X	X					
PLATELET-RICH PLASMA	E019	X										
PLATELETS	E020	X		X		X	X		X			
Washed PLATELETS	E021	X		X		X	X					
POOLED PLATELETS	E022			X		X	X		X			
Washed POOLED PLATELETS	E023			X		X	X					
APHERESIS PLATELETS	E024	X				X	X		X			
Frozen APHERESIS PLATELETS	E025		X									
Thawed APHERESIS PLATELETS	E026	X				X	X					

Allowable Product Modifications by Original Product Type												
This table defines the allowable Modification using the Original/Parent Unit's ICCBBA Product Type*												
Product Type Name	ICCBBA Product Type*	Split/Divide	Thaw	Pool	Thaw/Pool	Irradiate	Leuko-reduce	Volume Reduce	Wash	Rejuvenate	Freeze	Deglycerolize
Washed APHERESIS PLATELETS	E027	X					X	X				
CRYOPRECIPITATE	E028		X		X							
Thawed CRYOPRECIPITATE	E029	X		X			X					
POOLED CRYOPRECIPITATE	E030		X									
Thawed POOLED CRYOPRECIPITATE	E031			X			X					
APHERESIS CRYOPRECIPITATE	E032		X									
Thawed APHERESIS CRYOPRECIPITATE	E033	X					X					
GRANULOCYTES	E034	X		X			X					
APHERESIS GRANULOCYTES	E035	X					X					
POOLED GRANULOCYTES	E036			X			X					
APHERESIS GRANULOCYTES - PLATELETS	E037	X					X					
LEUKOCYTES	E038	X					X					
APHERESIS LEUKOCYTES	E039	X					X					
POOLED PLASMA	E040		X				X					
Liquid APHERESIS PLASMA	E041	X					X					

Allowable Product Modifications by Original Product Type												
This table defines the allowable Modification using the Original/Parent Unit's ICCBBA Product Type*												
Product Type Name	ICCBBA Product Type*	Split/Divide	Thaw	Pool	Thaw/Pool	Irradiate	Leuko-reduce	Volume Reduce	Wash	Rejuvinate	Freeze	Deglycerolize
PLATELET - RICH BUFFY COAT	E042	X				X						
POOLED PLATELET - RICH BUFFY COAT	E043					X						
APHERESIS LYMPHOCYTES	E044	X				X						
APHERESIS MONOCYTES* **	E045	X				X						
SERUM	E046											
POOLED SERUM**	E047											
Frozen POOLED SERUM	E048											
Washed APHERESIS RED BLOOD CELLS	E049	X				X	X			X		
Frozen APHERESIS RED BLOOD CELLS	E050											X
Deglycerolized APHERESIS RED BLOOD CELLS	E051	X				X	X			X		
Rejuvenated APHERESIS RED BLOOD CELLS	E052	X				X	X				X	
Frozen Rejuvenated APHERESIS RED BLOOD CELLS	E053											X
Deglycerolized Rejuvenated APHERESIS RED BLOOD CELLS	E054	X				X	X					

Allowable Product Modifications by Original Product Type												
This table defines the allowable Modification using the Original/Parent Unit's ICCBBA Product Type*												
Product Type Name	ICCBBA Product Type*	Split/Divide	Thaw	Pool	Thaw/Pool	Irradiate	Leuko-reduce	Volume Reduce	Wash	Rejuvenate	Freeze	Deglycerolize
Thawed POOLED PLASMA	E055					X	X					
Thawed POOLED FRESH FROZEN PLASMA	E056					X	X					
Washed GRANULOCYTES	E057	X				X						
Liquid POOLED PLASMA	E058					X	X					
POOLED FRESH FROZEN PLASMA	E059		X									
Frozen POOLED PLATELETS	E060		X									
Thawed POOLED PLATELETS	E061					X	X					
Thawed Apheresis POOLED FRESH FROZEN PLASMA	E062					X	X					
Thawed Apheresis POOLED PLASMA	E063					X	X					
Washed Thawed POOLED PLATELETS	E064											
Washed Thawed Apheresis PLATELETS	E065											
Lyophilized FRESH FROZEN PLASMA	E066		X									

Allowable Product Modifications by Original Product Type												
This table defines the allowable Modification using the Original/Parent Unit's ICCBBA Product Type*												
Product Type Name	ICCBBA Product Type*	Split/Divide	Thaw	Pool	Thaw/Pool	Irradiate	Leuko-reduce	Volume Reduce	Wash	Rejuvenate	Freeze	Deglycerolize
Reconstituted FRESH FROZEN PLASMA	E067											
Washed Apheresis POOLED PLATELETS	E068											
Frozen PLATELETS	E069		X									

**This product type has no ICCBBA assigned product codes.

National Treating Specialty Table

Table 22: National Treating Specialty Table


PT_105.01 National Treating Specialty Table	
Treating Specialty Code	Treating Specialty Name
VA1	ALLERGY
VA2	CARDIOLOGY
VA3	PULMONARY, TUBERCULOSIS
VA4	PULMONARY, NON-TB
VA5	GERONTOLOGY
VA6	DERMATOLOGY
VA7	ENDOCRINOLOGY
VA8	GASTROENTEROLOGY
VA9	HEMATOLOGY/ONCOLOGY
VA10	NEUROLOGY
VA11	EPILEPSY CENTER
VA12	MEDICAL ICU/CCU
VA14	METABOLIC
VA15	GENERAL (ACUTE MEDICINE)
VA16	CARDIAC-STEP DOWN UNIT
VA17	TELEMTRY
VA18	NEUROLOGY OBSERVATION
VA19	STROKE UNIT
VA20	REHABILITATION MEDICINE
VA21	BLIND REHAB
VA22	SPINAL CORD INJURY
VA23	SPINAL CORD INJURY OBSERVATION
VA24	MEDICAL OBSERVATION
VA25	PSYCH RESID REHAB TRMT PROG
VA26	PTSD RESIDENTIAL REHAB PROG
VA27	SUBSTANCE ABUSE RES TRMT PROG
VA28	HOMELESS CWT/TRANS RESID
VA29	SUBST ABUSE CWT/TRANS RESID
VA31	GEM ACUTE MEDICINE
VA32	GEM INTERMEDIATE CARE

PT_105.01 National Treating Specialty Table	
Treating Specialty Code	Treating Specialty Name
VA33	GEM PSYCHIATRIC BEDS
VA34	GEM NEUROLOGY
VA35	GEM REHABILITATION MEDICINE
VA36	BLIND REHAB OBSERVATION
VA37	DOMICILIARY CHV
VA38	PTSD CWT/TR
VA39	GENERAL CWT/TR
VA40	INTERMEDIATE MEDICINE
VA41	REHAB MEDICINE OBSERVATION
VA42	NH LONG STAY DEMENTIA CARE
VA43	NH LONG STAY SKILLED NURSING
VA44	NH LONG STAY MAINTENANCE CARE
VA45	NH LONG STAY PSYCHIATRIC CARE
VA46	NH LONG STAY SPINAL CORD INJ
VA47	NH RESPITE CARE (NHCU)
VA50	GENERAL SURGERY
VA51	GYNECOLOGY
VA52	NEUROSURGERY
VA53	OPHTHALMOLOGY
VA54	ORTHOPEDIC
VA55	OTORHINOLARYNGOLOGY
VA56	PLASTIC SURG, INC HEAD/NECK
VA57	PROCTOLOGY
VA58	THORACIC SURGERY, INC CARDIAC
VA59	UROLOGY
VA60	ORAL SURGERY
VA61	PODIATRY
VA62	PERIPHERAL VASCULAR
VA63	SURGICAL ICU
VA64	NH SHORT STAY REHABILITATION
VA65	SURGICAL OBSERVATION
VA66	NH SHORT STAY RESTORATIVE
VA67	NH SHORT STAY MAINTENANCE
VA68	NH SHORT STAY PSYCHIATRIC CARE
VA69	NH SHORT STAY DEMENTIA CARE
VA70	ACUTE PSYCHIATRY (<45 DAYS)
VA71	LONG-TERM PSYCHIATRY (>45 DAYS)
VA72	ALCOHOL DEPENDENCE TRMT UNIT
VA73	DRUG DEPENDENCE TRMT UNIT
VA74	SUBSTANCE ABUSE TRMT UNIT
VA75	HALFWAY HOUSE
VA76	PSYCHIATRIC MENTALLY INFIRM
VA77	PRRTP
VA79	SIPU (SPEC INPT PTSD UNIT)
VA80	NHCU
VA81	NH GEM NURSING HOME CARE
VA83	RESPITE CARE (MEDICINE)
VA84	SUBSTANCE ABUSE INTERMED CARE
VA85	DOMICILIARY
VA86	DOMICILIARY SUBSTANCE ABUSE
VA87	GEM DOMICILIARY
VA88	DOMICILIARY PTSD
VA89	STAR I, II & III
VA90	SUBST ABUSE STAR I, II & III

PT_105.01 National Treating Specialty Table	
Treating Specialty Code	Treating Specialty Name
VA91	EVAL/BRF TRMT PTSD UNIT(EBTPU)
VA92	GEN INTERMEDIATE PSYCH
VA93	HIGH INTENSITY GEN PSYCH INPAT
VA94	PSYCHIATRIC OBSERVATION
VA95	NH SHORT STAY SKILLED NURSING
VA96	NH HOSPICE
VA98	NON-DOD BEDS IN VA FACILITY
VA99	DOD BEDS IN VA FACILITY

Order Status Flowchart

Figure 135: Order Status Flowchart

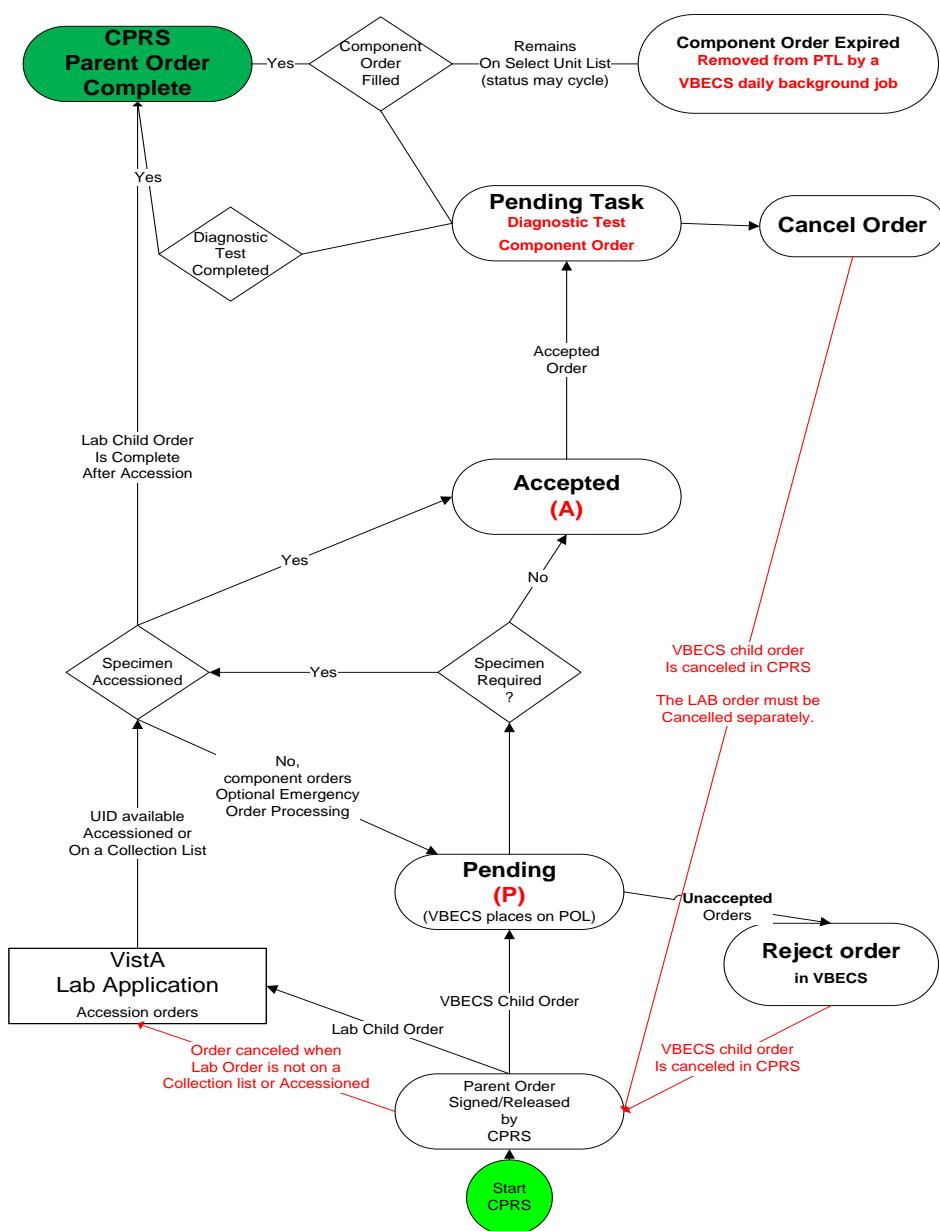

Rules for Electronic and Serologic Crossmatch

The WHOLE BLOOD, RED BLOOD CELLS, and APHERESIS RED BLOOD CELLS ICCBBA component classes require crossmatch evaluation, serologic or electronic (when enabled). For the OTHER ICCBBA component class, a reflex crossmatch test was ordered based on the technologist's evaluation of the unit itself. Each collection is different and the technologist must evaluate the red blood cell content of the unit at the time of selection.

When electronic crossmatch is enabled at a division and units that require crossmatch (WHOLE BLOOD, RED BLOOD CELLS, and APHERESIS RED BLOOD CELLS ICCBBA component classes or a reflex crossmatch test was ordered) are selected, the system performs a comprehensive set of data validation, which does not include inactivated patient or unit tests. This patient and unit data validation is to ensure that the combination of patient and blood unit testing allows the unit to be eligible for electronic crossmatch. When the combination does not allow for electronic crossmatch, VBECS warns the user that the electronic crossmatch will not be applied.

Table 23 provides the mapping for the actual parameters verified and associated error messages. All parameters must comply in order for the unit to be eligible for electronic crossmatch: there are no exceptions or overrides. The parameters are checked in the order listed and once a parameter fails, VBECS displays the appropriate error message, no further checking occurs, and the selected unit must have serologic crossmatch performed.

-  To use the electronic crossmatch function, sites must perform additional validation of electronic crossmatch to meet regulatory requirements. See *Guidelines for Implementing the Electronic Crossmatch* (Bethesda, MD: American Association of Blood Banks, 2003).
-  Changing a patient's antigen negative requirement or antibody identified after a unit is crossmatched does not provide an alert at Issue Unit regarding the type of crossmatch (electronic or serologic) performed on a unit available for issue.

Table 23: Rules for Electronic and Serologic Crossmatch

Parameter Checked	Error Message When Parameter Fails
Current specimen ABO/Rh and Antibody Screen testing must be complete.	Patient not eligible for eXM. Current specimen not fully tested.
Current specimen Antibody Screen interpretation must be negative. VBECS checks only the last Antibody Screen entered.	Patient not eligible for eXM due to positive Antibody Screen.
Patient must have no history of previously identified antibodies, regardless of clinical significance based on system rules or division (or from conversion) of the antibody ID.	Patient not eligible for eXM due to previous antibody history.
Patient cannot have a persistent antigen negative requirement regardless of the division of the requirement entry.	Patient not eligible for eXM due to transfusion antigen negative requirement.
The patient current ABO/Rh must match the historical record when there is a historical record, and there must be no previous justifications. An ABO/Rh marked as "entered in error" is not considered part of the patient's record.	Patient not eligible for eXM due to ABO/Rh discrepancy from previous record.
The patient cannot have a documented instance of current or previous Inconclusive ABO/Rh typing.	Patient not eligible for eXM due to ABO/Rh typing difficulty.
There must be at least two valid instances of ABO/Rh typing performed in the same division for the patient. One instance must be a valid ABO/Rh on the current specimen and not "unknown" or "inconclusive"; the other may be from a previous typing episode or a repeat ABO/Rh performed on the current specimen. The patient ABO/Rh interpretations must match in all instances. The user must	Patient not eligible for eXM: he must have two instances of ABO/Rh typing at the division for eXM, including the current specimen.

Parameter Checked	Error Message When Parameter Fails
enter valid forward and reverse typing for the system to use as the second test performed in the division.	
Selected unit must be ABO compatible based on system rules.	Unit not eligible for eXM. Unit is ABO incompatible with the patient.
Selected Rh negative units must have ABO and Rh confirmation results entered. The ABO and Rh confirmation interpretation must match the unit ABO/Rh at login.	Unit not eligible for eXM. ABO/Rh confirmation not performed.
Selected Rh positive units must have ABO confirmation results entered. The ABO confirmation interpretation must match the unit ABO at login.	Unit not eligible for eXM. ABO confirmation not performed.
Selected unit was previously serologically crossmatch incompatible with this patient.	Unit not eligible for eXM. Unit was previously serologically crossmatch incompatible with the patient.

Unit Status Flowchart


Types of Unit Statuses

- Temporary: awaiting ABO/Rh confirmation (lime green on the flowchart)
- Adjustable: may move in and out of the status as part of normal processing
- Issued: pending an update to transfusion or return (bright green in the flowchart)
- Final: the unit is no longer available for most processes in the division; the unit record remains active (pink and stop signs in the flowchart)

Indicators Associated with Units

- Quarantined: attaches to the unit record and status when the unit is quarantined. Removed by releasing the unit from quarantined.
- Presumed: attaches to the unit record only when the unit is in an issued unit status for more than 48 hours. Remains on the unit record if no transfusion information is entered through Post-Transfusion Information and Document ABO Incompatible Transfusion.

Figure 136: Unit Status Flowchart

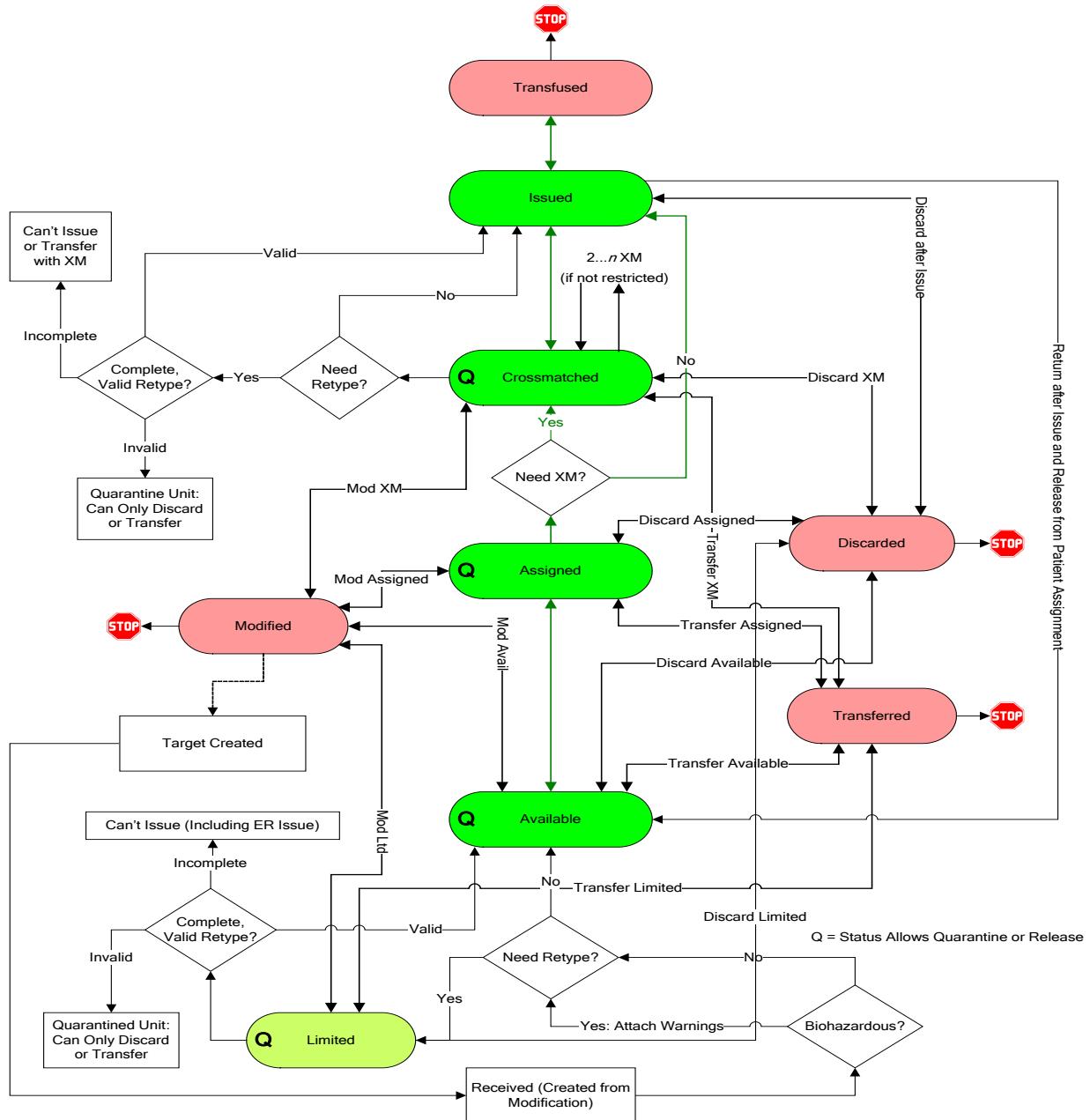

VBECS Maximum Surgical Blood Order Schedule

Table 24: VBECS Maximum Surgical Blood Order Schedule (MSBOS)

VBECS Maximum Surgical Blood Order Schedule (MSBOS)	
Abdominal-Perineal resection	Keratectomy
Adrenalectomy	Hernia repair
Amputation (A/K, B/K)	Hip Replacement, total
Aneurysm, aortic	Hip Replacement, revision
Aneurysm, abdominal	Knee Replacement, total
Aneurysm, thoracic	Laminectomy, cervical, thoracic, or lumbar
Aneurysm, cranial	Laminectomy, disk repair
Aorto-Femoral Bypass	Laparotomy, exploratory
Appendectomy	Laryngectomy
Coronary Bypass Graft (CABG)	Lobectomy, pulmonary
Coronary Bypass Graft (CABG), redo	Mandibulectomy
Carotid Endarterectomy	Mastectomy, radical
Cervical Discectomy	Maxillectomy
Cervical Fusion	Nephrectomy
Cholecystectomy	Pancreatectomy
Cholostomy	Parathyroidectomy
Colectomy	Portocaval Shunt
Colon Resection	Prostatectomy
Craniotomy	Prostatectomy, total
Cystectomy	Radical Neck Dissection
Cystolithotomy	Renal Transplant
Cystoscopy	Skin Flap
Embolectomy	Skin Graft
Esophageal Resection	Splenectomy
Esophogectomy	Thyroidectomy
Femoral-Popliteal Bypass	Tracheostomy
Fracture, open reduction	Transurethral Resection (TURP)
Fusion, cervical or lumber	Vagotomy
Gastrectomy	Valve replacement, aortic or mitral
Gastric Bypass	Other*
Glossectomy	
Hematoma Evacuation, sub- or epidural	

*Other surgeries entered in the Surgery Name field during order entry. (Local facilities may add surgeries to this list.) Entries made during order entry will not have MSBOS component order recommendations.

VBECS does not display “Other” in MSBOS.

BCE COTS Message Updates

Table 25: BCE COTS Message Updates

This table identifies the specific subsequent processing of an assigned unit for one patient unit association that may or may not result in a change to the unit status or other quality and how that may or may not trigger an status change that requires update to the BCE COTS product.

Unit status code legend:

A	Available	S	Assigned
C	Crossmatched	T	Transfused
D	Discarded	X	Transferred
I	Issued	DS	Dispense Unit (Issued)
L	Limited	RE	Release Unit
M	Modified	RS	Reserve Unit

BCE COTS Message Updates					
Row #	Initial Unit Status	Processing Option	Additional Information About the Processing Option	Unit Status Exiting the Option	Unit Status Reported to BCE COTS
1	A	Select Unit for Patient	Normal flow	S	RS
2	A	Select Unit for Patient	Emergency issue flow	S	RS
3	A	Select Unit for Patient	Electronic Crossmatch processed	C	RS
4	C	Process Outgoing Shipment	Maintain Assignment	X	RE upon confirmation
5	C	Process Outgoing Shipment	Release Assignment upon confirmation	A (then X when confirmed)	RE sent when release processed
6	S	Process Outgoing Shipment	Maintain Assignment	X	RE upon confirmation
7	S	Process Outgoing Shipment	Release Assignment upon confirmation	A (then X when confirmed)	RE sent when release processed
8	D	Remove Final Unit Status	Unit was previously Crossmatched	C	RS
9	D	Remove Final Unit Status	Unit was previously Assigned	S	RS
10	S	Remove Final Unit Status	Remove Modification	A	No message resent
11	X	Remove Final Unit Status	Unit was previously Crossmatched	C	RS
12	X	Remove Final Unit Status	Unit was previously Assigned	S	RS
13	T	Remove Final Unit Status	Normal flow	I	DS
14	S	Pool Units	Maintain Assignment	M	RE-original
15	S	Pool Units	Maintain Assignment	S	RS-target

BCE COTS Message Updates					
Row #	Initial Unit Status	Processing Option	Additional Information About the Processing Option	Unit Status Exiting the Option	Unit Status Reported to BCE COTS
16	S	Pool Units	Release Assignment	M	RE-original
17	S	Pool Units	Release Assignment	A	Target, no message sent
18	C**	Pool Units	Maintain Assignment	M	RE-original(s), may be multiple messages.
19	C**	Pool Units	Maintain Assignment	C	RS-target
20	C**	Pool Units	Release Assignment	M	RE-original
21	C**	Pool Units	Release Assignment	A	Target, no message sent
22	S	Split Units	Maintain Assignment	M	RE-original
23	S	Split Units	Maintain Assignment	S	RS-targets, one message for each target created.
24	S	Split Units	Release Assignment	M	RE-original
25	S	Split Units	Release Assignment	A	Target, no message sent
26	C	Split Units	Maintain Assignment	M	RE-original
27	C	Split Units	Maintain Assignment	C	RS-target
28	C	Split Units	Release Assignment	M	RE-original
29	C	Split Units	Release Assignment	A	Target, no message sent
30	S	Modify Units	Maintain Assignment	M	RE-original(when processing in batch, multiple messages are sent)
31	S	Modify Units	Maintain Assignment	S	RS-target (when processing in batch, multiple messages are sent)
32	S	Modify Units	Release Assignment	M	RE-original(when processing in batch, multiple messages are sent)
33	S	Modify Units	Release Assignment	A	Target, no message sent
34	C	Modify Units	Maintain Assignment	M	RE-original(when processing in batch, multiple messages are sent)
35	C	Modify Units	Maintain Assignment	C	RS-target (when processing in batch, multiple messages are sent)
36	C	Modify Units	Release Assignment	M	RE-original(when processing in batch, multiple messages are sent)
37	C	Modify Units	Release Assignment	A	Target, no message sent

BCE COTS Message Updates					
Row #	Initial Unit Status	Processing Option	Additional Information About the Processing Option	Unit Status Exiting the Option	Unit Status Reported to BCE COTS
38	C	Discard or Quarantine Unit	Discarded or Wasted	D	RE
39	S	Discard or Quarantine Unit	Discarded or Wasted	D	RE
40	S	Discard or Quarantine Unit	Quarantine	A	RE
41	C	Discard or Quarantine Unit	Quarantine	A	RE
42	S	Record Testing when saving Crossmatch Results	Compatible (C)	C	No message is sent
43	S	Record Testing when saving Crossmatch Results	Incompatible (I)	A	RE
44	S	Record Testing when saving Crossmatch Results	Crossmatch Compatible: Don't Transfuse (D)	A	RE
45	S	Record Testing when saving Crossmatch Results	Compatible-Give with MD approval (G)	C	No message is sent
46	S	Record Phenotype Results	Unit is antigen positive for a patient's antigen negative requirement	S	No message is sent
47	C	Record Phenotype Results	Unit is antigen positive for a patient's antigen negative requirement	C	No message is sent
48	S	Record Phenotype Results	Antigen typing result creates a discrepant antigen typing.	A	RE
49	C	Record Phenotype Results	Antigen typing result creates a discrepant antigen typing.	A	RE
50	S	Release Unit from Patient Assignment	Normal flow, or any overrides processed	A	RE
51	C	Release Unit from Patient Assignment	Normal flow, or any overrides processed	A	RE
52	I	Return Issued Units to Blood Bank	Unit without a crossmatch	S	RS
53	I	Return Issued Units to Blood Bank	Crossmatched unit	C	RS
54	I	Return Issued Units to Blood Bank	Unit without a crossmatch, unit is unsatisfactory	A	RE

BCE COTS Message Updates					
Row #	Initial Unit Status	Processing Option	Additional Information About the Processing Option	Unit Status Exiting the Option	Unit Status Reported to BCE COTS
55	I	Return Issued Units to Blood Bank	Unit crossmatched, unit is unsatisfactory	A	RE
56	S	Issue Blood Components	Normal flow with or without overrides	I	DS
57	C	Issue Blood Components	Normal flow with or without overrides	I	DS
58	S	Issue Blood Components	Unit is marked as unsatisfactory for issue	A	RE
59	C	Issue Blood Components	Unit is marked as unsatisfactory for issue	A	RE
60	C	Invalidate Patient Test Results	Invalidate Crossmatch results for the unit-patient association	S	No message is sent as the unit has not been released from assignment from the patient.

**Pooling of products that also require crossmatch is not currently possible.

Workload Process to Use Case Mapping Table

Table 26: Workload Process to Use Case Mapping Table

*Transaction type: Patient (P), Unit (U), Misc (M), no donor workload types in VBECS at this time.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
ABO Forward and reverse typing (patient)	Transfusion Reaction Workup	Positive	P	An ABO/Rh test for pre or post is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for Pre and Post entries.
	Patient Testing: ABO/Rh	Positive	P	Accrue workload when an CPRS ordered ABO/Rh test performed from the ordered test ABO/Rh only (Record Patient ABO-Rh is completed and saved by Record Testing and was not performed as part of the Type & Screen) (when completed and saved via Record Testing)
	Invalidate Patient Testing: Patient ABO/Rh	Negative	P	Accrue workload when a completed ABO/Rh test is invalidated.
ABO Forward and reverse typing (patient) Repeat Test	Patient Testing: ABO/Rh	Positive	M	Accrue workload when a reflex or repeat ABO/Rh test is performed, completed and saved via Record Testing. Note: ABO/Rh ordered in Order Reflex Test from main menu or via use of Order Reflex Test from Record Testing
	Invalidate Patient Testing: Patient ABO/Rh	Negative	M	Accrue workload when a reflex or repeat ABO/Rh test is invalidated.
ABO forward typing (unit)	Unit ABO/RH confirmation test	Positive	U	An ABO confirmation test is performed. When multiple units are selected in a batch, each unit in the batch accrues a workload event. Note: Workload generated when Anti-D testing is not included in the unit's confirmation test. Workload is not accrued when an ABO or Rh discrepancy override is processed and VBECS releases all patient assignments. Workload is not accrued when VBECS quarantines the unit due to a discrepancy. No special handling for workload collection for additional confirmation tests on a unit.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
	Invalidate Unit ABO/RH Confirmation test	Negative	U	Accrue workload when an ABO confirmation test is invalidated
ABO/Rh forward typing (unit)	Unit ABO/RH confirmation test	Positive	U	An ABO/Rh confirmation test is performed. When multiple units are selected in a batch, each unit in the batch accrues a workload event. Note: Workload generated when Ant-D is part of the unit's confirmation test. Workload is not accrued when an ABO or Rh discrepancy override is processed and VBECS releases all patient assignments. Workload is not accrued when VBECS quarantines the unit due to a discrepancy. Any unit successfully confirmed accrues workload. For modification Split: Workload is not inherited by split units. A split unit that requires confirmation accrues confirmation workload at the time of testing. No special handling for workload collection for additional confirmation tests on a unit.
	Invalidate Unit ABO/RH Confirmation test	Negative	U	Accrue workload when an ABO/Rh confirmation test is invalidated
Accept Order	Accept Order	Positive	M	Accrue workload when an order is accepted. When a multiple orders are selected, each order accrues workload.
Antibody identification Work-Up	ABID Test	Positive	P	User enters additional workload associated with the individual reflex ordered ABID. The selected VBECS multiplier will multiply against the VistA multiplier and display the (multiplication) product on the report in Division Workload Report.
	Invalidate ABID Test	Negative	P	Accrue workload when the ABID is invalidated. All workload processes associated in Enter Reflex Test Results, including the selected VBECS multiplier, are invalidated
Antibody Screen (patient)	Transfusion reaction workup	Positive	P	An ABS test for pre or post is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for Pre and Post entries.
	Patient Testing ABS test	Positive	P	Accrue workload when an ordered ABS test performed from the CPRS ordered test ABS only (Record Patient Antibody Screen is completed and saved by Record Testing. The workload accrues when the ABS is completed and saved via Record Testing no matter how many cells (rows) are required by the division configuration.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
	Invalidate Patient Testing ABS test	Negative	P	Accrue workload when a completed ABS test is invalidated.
Antibody Screen (patient) Repeat Test	Patient Testing ABS test	Positive	M	Accrue workload when a reflex or repeat ABS test is performed, completed and saved via Record Testing. Note: ABS ordered in Order Reflex Test from main menu or via use of Order Reflex Test from Record Testing.
	Invalidate Patient Testing ABS test	Negative	M	Accrue workload when a reflex or repeat ABS test is invalidated.
Antigen phenotyping, Single Test phase (QC)	QC for unit or patient antigen typing	Positive	M	Accrue workload when Antiserum QC is included in Record Phenotype Results or Record Patient Antigen Types includes the testing of both the positive and negative control cells, per specificity by lot number, when only one phase of reactivity is chosen for the test grid, IS or AHG. One workload event is collected per completed tab for regular or repeat antigen tests.
Antigen phenotyping, Multiple Test phases (QC)	QC for unit or patient antigen typing	Positive	M	Accrue workload when Antiserum QC is included in Record Phenotype Results or Record Patient Antigen Types includes the testing of both the positive and negative control cells, per specificity by lot number, when only multiple phases of reactivity is chosen for the test grid, IS/RT, RT/37 or weak D. One workload event is collected per completed tab for regular or repeat antigen tests. When weak D is the selected test, QC may not be accrued for the rack selection. QC is accrued when positive and negative cells must be tested for the lot number.
Cancel Order	Accept Order, Reject Order	Positive	M	Accrue workload when an order on the pending order list is canceled. When a multiple orders are canceled, each order accrues workload.
Cancel Order	Patient Testing List, Cancel order	Positive	M	Accrue workload when an order on the pending task list is canceled. When a multiple orders are canceled, each order accrues workload.
Crossmatch unit, electronic	Select Unit	Positive	P	This process is invoked when an individual unit is selected for patient assignment and the unit is electronically crossmatched. When a multiple units are selected, each unit accrues workload.

Workload Process to Use Case Mapping Table				
VBECs Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
Crossmatch unit, serologic immediate spin	Patient Testing, serologic crossmatch	Positive	P	Accrue workload when an individual unit (row) crossmatch, selected to include only the IS phase, is completed and is saved by Record Testing. When a multiple units are selected, each unit (row) accrues workload.
	Invalidate Patient Testing, serologic crossmatch	Negative	P	Accrue workload when a completed crossmatch test is invalidated. The workload originally saved with the test, serologic immediate spin.
Crossmatch unit, serological Coombs	Transfusion Reaction Workup, serologic crossmatch	Positive	P	A crossmatch test for pre or post is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for Pre and Post entries. When a multiple units are selected, each unit accrues workload.
	Patient Testing, serologic crossmatch	Positive	P	Accrue workload when an individual unit (row) crossmatch, selected to include all phases or only AHG phase, is completed and is saved by Record Testing. When a multiple units are selected, each unit (row) accrues workload.
	Invalidate Patient Testing, serologic crossmatch	Negative	P	Accrue workload when a completed crossmatch test is invalidated. The workload originally saved with the test, serological Coombs.
Crossmatch, Repeat Test	Patient Testing, serologic crossmatch	Positive	M	Accrue workload when an individual unit (row) crossmatch, selected to include all phases or IS or only AHG phase, is completed and is saved by Record Testing. When a multiple units are selected, each unit (row) accrues workload. Note: ordered via Order Reflex Test from main menu or the use of Order Reflex Test by Record Testing.
	Invalidate Patient Testing, serologic crossmatch	Negative	M	Accrue workload when an individual unit crossmatch is invalidated.
Daily Rack Quality Control (QC)	Enter Daily Reagent Rack QC	Positive	M	Accrue workload when Daily QC rack completed for one individual rack includes all rows in configured QC. When multiple racks are tested, each completed and saved tab accrues a workload event.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
DAT (QC)	Enter Daily Reagent Rack QC or DAT QC	Positive	M	Accrue workload when Reagent QC completed or in Record Patient Direct Antiglobulin Test includes the testing of both the positive and negative control cells, per specificity per lot number, when only one phase of reactivity is chosen for the test grid, IS or AHG. One workload event is collected per completed tab for regular or repeat antiglobulin tests (PS, IgG, Comp) tests.
Direct Antiglobulin Test (DAT)	Transfusion Reaction Workup, DAT test	Positive	P	A DAT test for pre or post is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for Pre and Post entries.
	Patient Testing, DAT	Positive	P	Accrue workload when a DAT completed and saved (Record Patient Direct Antiglobulin Test is saved by Record Testing as part of a TAS configuration or individual test ordered. This count is used for all antiglobulin tests (PS, IgG, Comp) when ordered from CPRS or Reflex testing.
	Invalidate Patient Testing, DAT	Negative	P	Accrue workload when a completed DAT, PS, IgG, or Comp is invalidated.
Direct Antiglobulin Test (DAT) Repeat test	Patient Testing, DAT	Positive	M	Accrue workload when a reflex or repeat DAT test is performed, completed and saved via Record Testing. This applies to all repeat antiglobulin tests (PS, IgG, Comp) Note: Any DAT ordered in Order Reflex Test from main menu or via use of Order Reflex Test from Record Testing.
	Invalidate Patient Testing, DAT	Negative	M	Accrue workload when a completed Repeat DAT, PS, IgG, or Comp is invalidated.
Deglycerolize unit	Modify Unit, Deglycerolize	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Deglycerolize. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Deglycerolize	Negative	U	An individual blood unit's status is invalidated when the original modification process was Deglycerolize.
Discard unit	Discard Unit	Positive	U	Accrue workload when an individual blood unit's status is invalidated. When a batch of units is selected, each unit accrues workload.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
	Invalidate Unit Status, Discard	Negative	U	Accrue workload when a unit is discarded for waste or credit. When a batch of units is selected, each unit accrues workload.
Freeze unit	Modify Unit, Freeze	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Freeze. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Freeze	Negative	U	An individual blood unit's status is invalidated when the original modification process was Freeze.
Irradiate unit	Modify Unit, Irradiate	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Irradiate. When a batch of units is irradiated, each unit accrues workload. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Irradiate	Negative	U	An individual blood unit's status is invalidated when the original modification process was Irradiate.
Leukoreduce unit	Modify Unit, Leukoreduce	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Leukoreduce. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Leukoreduce	Negative	U	An individual blood unit's status is invalidated when the original modification process was Leukoreduce.
Split unit	Modify Unit, Split	Positive	U	Accrue workload when a unit modification of Split and a single workload event is recorded regardless of the number of units created by the modification. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit. Workload accrued by splitting is invalidated in Remove Final Unit Status with the inactivation of the modification</i>

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
	Invalidate Unit Status, Split	Negative	U	A Split Unit has its unit status invalidated. A single workload event is recorded regardless of the number of units originally created by the modification.
Rejuvenate unit	Modify Unit, Rejuvenate	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Rejuvenate. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Rejuvenate	Negative	U	An individual blood unit's status is invalidated when the original modification process was Rejuvenate.
Thaw	Modify Unit, Thaw	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Thaw. When a batch of units is thawed, each unit accrues workload. This applies to Thaw FFP and Thaw Cryo. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Thaw	Negative	U	An individual blood unit's status is invalidated when the original modification process was Thaw. This modification type is applicable to Thaw FFP and Thaw Cryo.
Wash unit	Modify Unit, Wash	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Wash. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Wash	Negative	U	An individual blood unit's status is invalidated when the original modification process was Wash.
Volume Reduce	Modify Unit, Volume Reduce	Positive	U	Accrue workload when an individual blood unit is processed by the modification type: Volume Reduce. <i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i>
	Invalidate Unit Status, Volume Reduce	Negative	U	An individual blood unit's status is invalidated when the original modification process was Volume Reduce.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
Issue unit	Issue Blood Unit	Positive	P	Accrue workload when a unit (row) is issued or is returned from issue to a patient. When a batch of units is processed, each unit invokes one workload process.
Justification	Justify ABO/Rh Change	Positive	M	Workload is accrued when a patient's ABO or Rh typing is justified. One workload event per patient justification.
Login Equipment	Equipment, Login	Positive	M	Accrue workload when a lot number of any type of equipment is logged into the system.
Login Reagent	Reagent, Login	Positive	M	Accrue workload when a lot number of any type of reagent is logged into the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
Login Supply	Supply, Login	Positive	M	Accrue workload when a lot number of any type of supply is logged into the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
Maintain Specimen	Maintain Specimen	Positive	M	Accrue workload when a specimen is maintained during order acceptance when a specimen is requirement for the order. Note: This event is accrued when the specimen is maintained in Maintain Specimen. This is collected in addition to the accept order workload accrued by Accept Order. Marking a specimen unacceptable does not create a negative workload event.
Patient antigen phenotype (multiple phases)	Patient Antigen typing	Positive	P	Accrue workload when a patient antigen phenotype test [specificity (tab)] with IS/RT or IS/37 phases is completed and saved (when completed and saved via Record Testing). One workload event is collected per completed tab for repeat or regular antigen tests.
	Invalidate, Patient Antigen typing	Negative	P	Accrue workload when a patient antigen phenotype test as defined by the antiserum specificity tested with any phases is invalidated.
Patient antigen phenotype (single phase)	Patient Antigen typing	Positive	P	Accrue workload when a patient antigen phenotype test [specificity (tab)] with AHG or IS phase is completed and saved (when completed and saved via Record Testing). One workload event is collected per completed tab for repeat or regular antigen tests.
	Invalidate, Patient Antigen typing	Negative	P	Accrue workload when a patient antigen phenotype test as defined by the antiserum specificity tested with a single phases is invalidated.

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
Pool unit	Modify, Pool	Positive	U	<p>Accrue workload when a pooled unit is created and a single workload event is recorded regardless of the number of units included in the pooled unit. This applies to modification type: Pool. Add/Remove unit from a pool does not accrue any workload.</p> <p><i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i></p> <p><i>Workload accrued by Pooling is invalidated in Edit Unit Processing Information with the inactivation of the pooled unit.</i></p>
	Edit Unit, Inactive Unit (that was in a Pool)	Negative	U	Accrue workload when a unit is inactivated (unit record inactivated) when the pooled unit was created in VBECS.
	Invalidate unit status, Pool	None	N/A	No effect on workload accrual when a unit is removed from a modified status that was included in a pool.
Quarantine unit	Discard/Quarantine, Quarantine unit	Positive	U	Accrue workload when a unit is marked for quarantine. When a batch of units is selected, each unit accrues workload.
Release directed unit	Free Directed Unit	Positive	U	Accrue workload when an individual blood unit with the restriction type of directed is released for use as an allogeneic unit.
Release unit from patient back to inventory	Release Unit from Patient Assignment	Positive	U	Accrue workload when an individual unit is released from patient assignment. When a multiple units are selected, each unit accrues workload.
Release unit from Quarantine	Discard/Quarantine, Release Unit from Quarantine	Positive	U	Accrue workload when a unit is released from quarantine. When a batch of units is selected, each unit accrues workload.
Return Issued unit	Return Issued unit to Inventory	Positive	U	Accrue workload when a unit is returned from issue status.
See Serologic Crossmatch, Patient ABO/Rh, Patient ABS, Patient DAT, Patient Antigen Typing	Patient Testing	Positive	Various	<i>Patient Testing encapsulates Serologic Crossmatch, Patient ABO/Rh, Patient ABS, Patient DAT, Patient Antigen Typing. Patient Testing does not in and of itself accrue workload, but accrues workload for the various testing grids for which it controls the data save. See the various options in this list for more information.</i>

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
Thaw/pool Cryo	Modify Unit, Thaw/Pool	Positive	U	<p>Accrue workload when an individual unit modification of Thaw/Pool Cryo. A single workload event is recorded regardless of the number of units included in the pooled unit.</p> <p><i>Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.</i></p>
	Edit, Unit, Inactivate Unit (that was in a Thaw/Pool)	Negative	U	Accrue workload when a unit is inactivated (unit record inactivated) when the pooled unit was created in VBECS.
	Invalidate Unit status, Thaw/Pool	None	N/A	No effect on workload accrual when a unit is removed from a modified status that was included in a Thaw/pool Cryo pool.
Transfuse Unit	Enter Post Transfusion Data	Positive	U	Accrue workload when an individual blood unit's status is updated to transfused by Enter Post-Transfusion Details basic flow or Alt Flow A.
	Invalidate Unit Status, Transfused	Negative	U	An individual blood unit's status is invalidated when the unit was in a status of transfused.
Transfusion Reaction Investigation	Patient Testing, TRW	Positive	P	<p>Accrue workload when a transfusion reaction investigation is saved by Record Transfusion Reaction Workup. <i>This does not include workload accrued by the optional TRW serologic testing.</i></p> <p><i>Note: The Division Workload Report workload report reports the sum of Workload Events recorded in Record Transfusion Reaction Workup for an individual TRW.</i></p>
	Invalidate Patient Testing, TRW	Negative	P	<p>Accrue workload when a transfusion reaction investigation previously saved by Record Transfusion Reaction Workup is invalidated. The workload for the transfusion reaction, as well as all of the workload accrued by the optional TRW serologic testing (pre- and post-ABO/Rh, pre- and post-ABS, pre- and post-DAT) are invalidated, as included in selected TRW (Record Transfusion Reaction Workup).</p> <p><i>Note: Any pre- or post-crossmatch tests performed as part of the TRW must be invalidated individually.</i></p>
Unit Antigen phenotyping, Multiple Test phases	Unit Antigen typing	Positive	U	<p>Accrue workload when a unit antigen phenotype test [specificity (tab)] with IS/RT or IS/37 phases is selected and completed for an individual blood unit (row).</p> <p>No special handling for workload collection for additional repeat antigen typing tests on a unit.</p>

Workload Process to Use Case Mapping Table				
VBECS Process	Record Saved in This Use Case	Positive or Negative Workload Accrued	Transaction Type*	Explanation
	Edit Unit, Invalidate Unit Antigen Typing	Negative	U	Accrue workload when a unit antigen phenotype test (specificity) with Multiple Test phases is invalidated for an individual blood unit.
Unit Antigen phenotyping, Single Test phase	Unit Antigen typing	Positive	U	A unit (row) antigen phenotype test [specificity (tab)] with AHG or IS phase is selected and completed for an individual blood unit (row). No special handling for workload collection for additional repeat antigen typing tests on a unit.
	Edit Unit, Invalidate Unit Antigen Typing	Negative	U	Accrue workload when a unit antigen phenotype test (specificity) with Single Test phase is invalidated for an individual blood unit.
Unit login	Incoming Shipment	Positive	U	An individual unit record is activated as saved to an incoming shipment invoice. When a multiple units are entered, each unit added to the database accrues workload.
	Edit Unit, Invalidate Unit Antigen Typing	Negative	U	Accrue workload when a unit is inactivated (unit record inactivated) when logged in thru incoming shipment or is a pooled unit created in VBECS. When the unit was created by split modification no workload is invalidated in this option.
Unit logout	Outgoing Shipment	Positive	U	An individual unit's status is updated to <i>transferred</i> on a confirmed outgoing shipment invoice. When a multiple units are selected, each unit accrues workload. Accrue workload upon the confirmation of the invoice, not the addition of a unit to a temporary outgoing shipment invoice as an invoice may be confirmed only once.
	Invalidate Unit Status, Transferred	Negative	U	An individual unit status is invalidated when the unit had a previous unit status of transferred.
Update Equipment Record	Equipment, Update	Positive	M	Accrue workload when a lot number of any type of equipment is updated in the system.
Update Reagent Inventory	Reagents, Update	Positive	M	Accrue workload when a lot number of any type of reagent is updated in the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
Update Supply Inventory	Supplies, Update	Positive	M	Accrue workload when a lot number of any type of supply is updated in the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.

Appendix C: VBECS Work Flow

Each of the six user roles is associated with a security level (see Table 1). Privileges accumulate as the security level increases. For example, a Lead Technologist's privileges include those of a Blood Bank Technologist and Enhanced Technologist.

Figure 137: VBECS Work Flow (left)


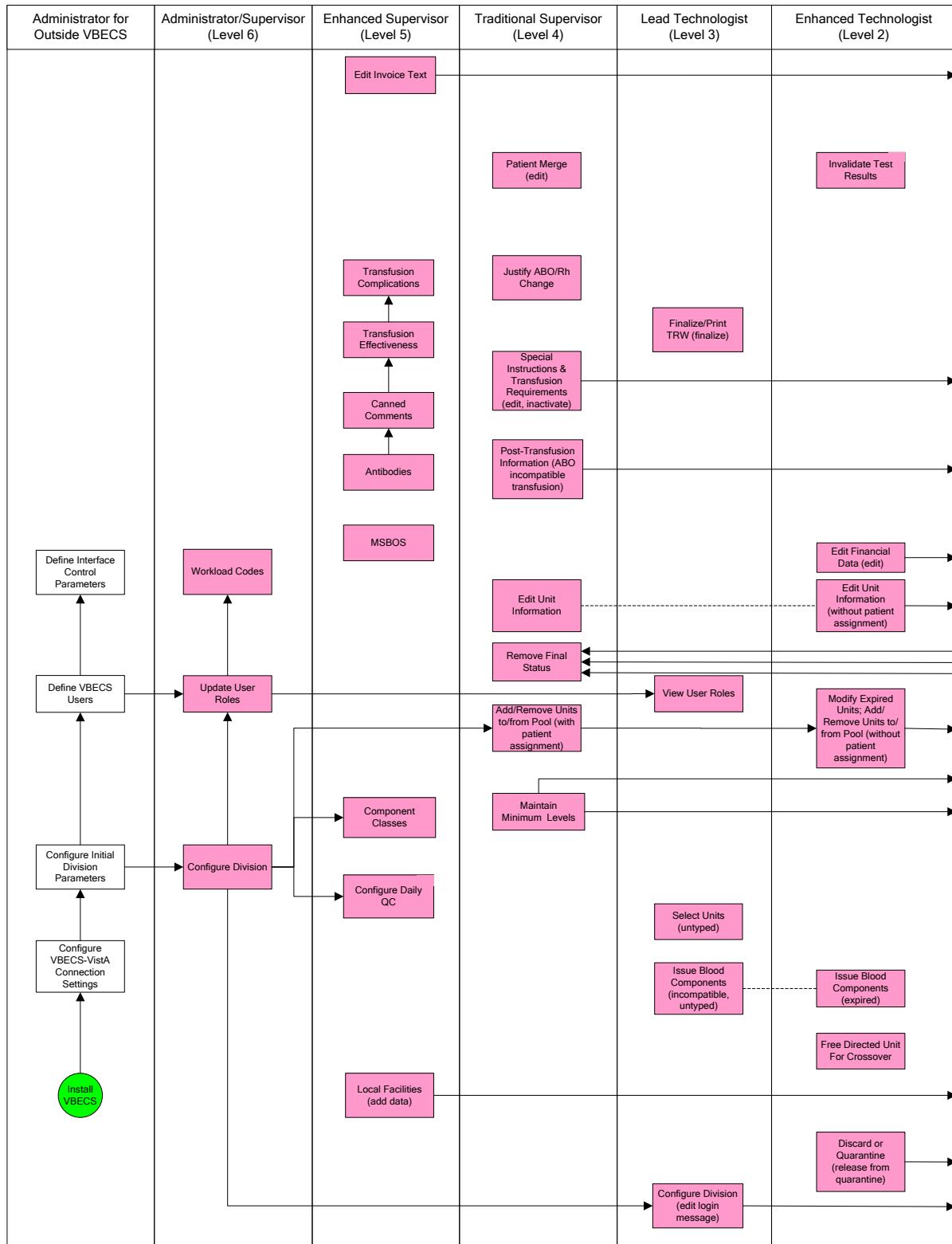


Figure 138: VBECS Work Flow (right)

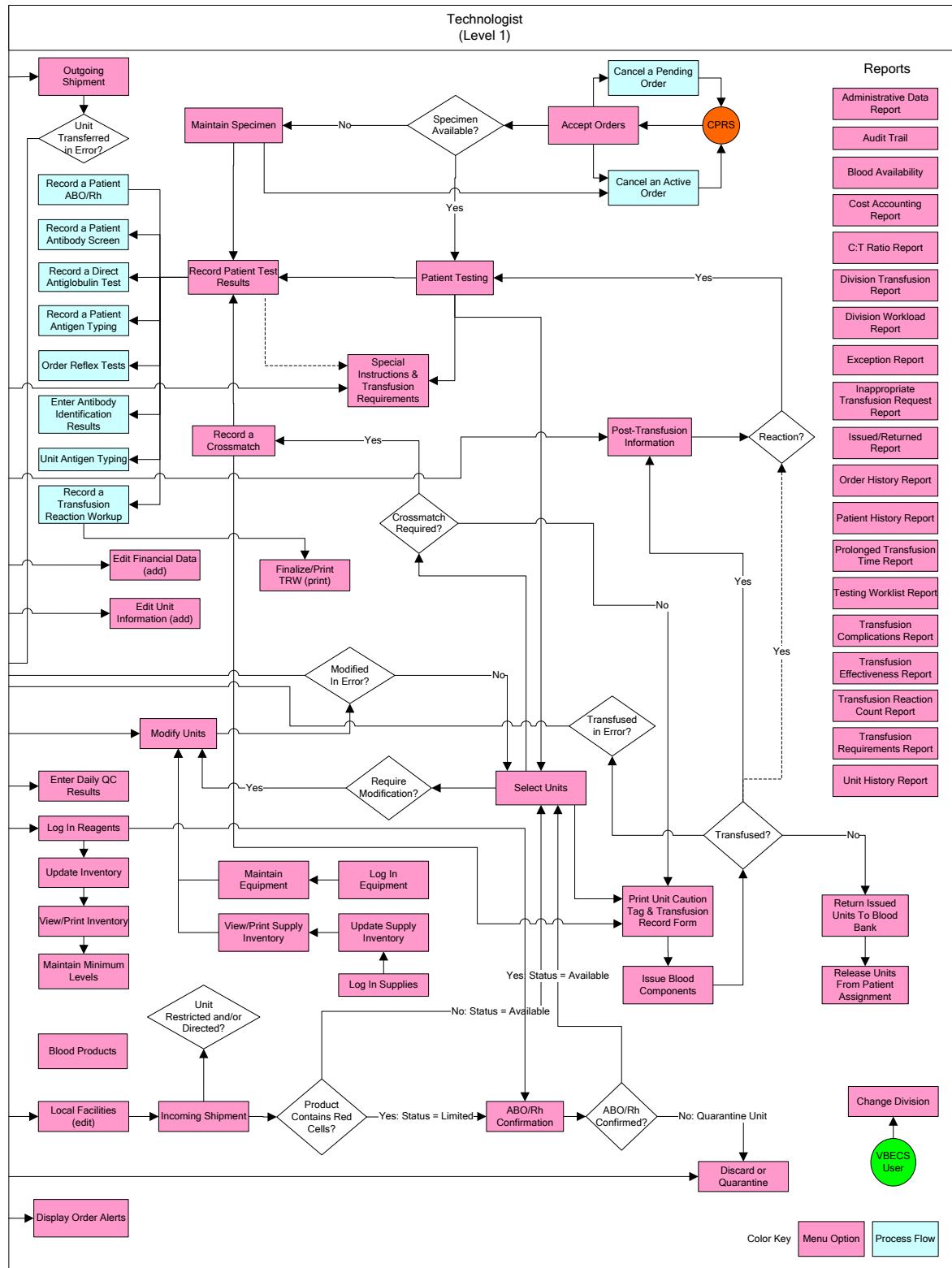

This page intentionally left blank.

Appendix D: Limitations and Restrictions

Limitations and restrictions are listed by the options in which they occur.

Table 27: Limitations and Restrictions

Functions Not Supported
Blood donor collection activities (donor module)
Direct data entry of antibody identification
Direct data entry of transfusion reaction workups
Patient testing, transfusions, and blood components recorded in Legacy VistA are not available in VBECS
Pediatric blood product preparation
Printing of blood product labels
User defined customized queries or reports
User editing blood product, antigen, antibody, or transfusion reaction tables
Validation record storage

Throughout VBECS
After entering a canned or free-text comment in the testing comment field of the testing grid, the user must press Enter to store the comment.
Although the copy-and-paste function is enabled, do not use it to enter patient identification information. It is highly recommended that the user scan this information when possible; otherwise, enter it .
Data retrieval from VBECS is available only by approved database integration agreement.
Information may not be transferred from one window to another without refreshing the screen. When a window is partially displayed or freezes; or when you move from one window to another, enter information, return to the original window, and the entered information is not available, minimize then maximize the window to refresh it.
Initial blood orders must be placed in CPRS. Tests and products may be added to existing orders in VBECS.
Maximize option windows to display safety-critical information in full.
The date and time displayed in VBECS depend on the date and time set on the server. Sites are not notified when an administrator changes the date or time on the server. Users must contact their System Administrator when the date or time is incorrect.
VBECS depends on a connection to VistA and CPRS to retrieve patient information and orders, respectively.
When a message window requires a user response and VBECS displays it behind the active window, the user must minimize the active window or click and drag dialog boxes to expose the message window.
When the default screen resolution is changed, division or VistALink information may not be entirely visible.
When the user does not enter a full or partial unit ID and clicks the find button in the Unit Search Screen, response time will vary based on the size of the database.

ABO/Rh Confirmation
When building a worklist by invoice number, only invoiced units not previously tested and confirmed appear on the worklist. The user may deselect or add additional units.

ABO/Rh Confirmation

For batch processing, the user can select incoming shipment invoices from his division, processed within the previous five days.

VBECS does not display product codes associated with unit IDs in ABO/Rh Confirmation. Therefore, do not create ABO/Rh confirmation worklists containing units with the same unit ID.

Accept Orders: Accept an Order

An order for TAS must be accepted before the component order if the user wants to include the component request to the Inappropriate Transfusion Request Report. Accepting the TAS after the component order is marked inappropriate and accepted clears the user's indication that this is an inappropriate request.

When no CPRS order exists, the user must process the specimen offline.

The user must accept the Type & Screen (TAS) before he accepts other diagnostic tests or component orders.

The user cannot click  (Restricted or Assigned Units) to access information on autologous or directed units that may be restricted to a patient.

Orders placed for patients with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in VistA are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.

A VBECS order on a Collection List has to be received with a date and time in VistA to make the specimen available in VBECS.

Accept Orders: Cancel a Pending Order

Cancellation of one order does not affect other orders placed with it.

Accept Orders: Pending Order List

Orders placed for patient's with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in VistA are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.

Add/Remove Units from a Pool

See Modify Units.

Audit Trail

VBECS does not accommodate online review (verification, signature) of reports.

Blood Products

Blood product records distributed through VBECS are national standards.

The division must have activated the valid target blood products for modification.

Canned Comments

The user may not edit comments provided with VBECS, but may create additional comments.

Site personnel must enter comments in accordance with local policy and procedures.

Component Classes
The user may not change the specimen requirement for a component class when there are incomplete orders for that component class.
There is no retrospective data entry or user selection for this option.
Configure Daily QC
This option pertains only to routine reagents.
When there are partially or currently QC'd racks in the system, the user may not change the QC template.
Configure Division
Changes made to the division setting are not applied to a user until they re-connect to VBECS after the change is made. It is highly recommended that all users log off VBECS prior to any configuration setting change.
Discard or Quarantine
When a user retroactively updates a unit status through this option, the Unit History Report displays the updated information, but does not display the date the change was made.
Display Order Alerts
VBECS issues no notifications for order cancellations.
Display Patient Merge Alerts
VBECS does not notify users of additional Vista patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed Vista patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor updates the patient records and clears the icon.
Display Patient Update Alerts
VBECS does not notify users of additional Vista patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed Vista patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor updates the patient records and clears the icon.
Patient updates involving patient's with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in Vista are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.
Document ABO Incompatible Transfusions
VBECS does not enforce compatibility checks between the patient and the unit.
Edit Financial Data
When a Technologist enters a unit ID and product code for a unit that is not in final status, VBECS does not enable the Special Test Cost field. An Enhanced Technologist or higher access can change the Special Test Cost field in the Edit Financial Data option.

Edit Unit Information
The user may correct a limited amount of unit information.
VBECS does not limit the discarded plasma volume.
If errors are found in the unit's product code, division or blood type at login, supplier, or donation type, the user must inactivate the unit and reenter it correctly.
When a unit does not have a patient association, editable fields are available to an ►► Enhanced Technologist.
ABO/Rh confirmation and antigen type testing cannot be inactivated when the unit is in a final unit status.

Enter Daily QC Results
Additional user-defined reagent types not available in Configure Daily QC or Enter Daily QC Results are not part of a reagent rack's daily QC testing.

Equipment
VBECS does not check for duplicate entries of equipment names and identifiers.
There is no retrospective data entry or user selection for this option.

Exception Report
VBECS does not accommodate online review (verification, signature) of reports.

Finalize/Print TRW
The user cannot invalidate or change a finalized TRW.
Transfusion reaction reports are not available to clinicians through CPRS: they must be printed and filed with patient charts.
VBECS displays these default values for transfusion reactions converted from VistA to VBECS. They do not reflect actual interpretations or results:

Converted Transfusion Reaction	Default
Pre-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Post-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Symptoms	No symptoms identified.
Implicated units	No units implicated in reaction.
Serologic tests	No serologic tests performed.

Free Directed Unit For Crossover
The unit must not be in a final status to release the restriction.
This option does not allow the removal of patient restrictions for most of the ISBT donation types that are eligible for crossover.
This option cannot process multiple directed units simultaneously.

Functions Not Supported
Blood donor collection activities (donor module)

Functions Not Supported
Direct data entry of antibody identification
Direct data entry of transfusion reaction workups
Patient testing, transfusions, and blood components recorded in Legacy VistA are not available in VBECS
Pediatric blood product preparation
Printing of blood product labels
User defined customized queries or reports
User editing blood product, antigen, antibody, or transfusion reaction tables
Validation record storage

Incoming Shipment
VBECS does not compare the ABO/Rh of an autologous unit or display a message when the ABO/Rh is incompatible with the restricted for patient.
When a product type, code, or local code is not in VBECS or a blood supplier creates a product code, contact the Help Desk. See Customer Support for contact information.
The unit expiration [Maximum Storage Time (MST)] is set to the longest possible storage time based on product type, core condition, attributes, and available references as a guide. The MST is set at a higher level than an individual collection product code. VBECS does not know the specific collection details. It is the responsibility of each user to accurately enter and process any one blood unit based on its expiration date/time as labeled by the collection facility.
The user must enter the donation type for Codabar-labeled units.
When a user inadvertently enters a unit ID that is longer than seven characters for a Codabar unit, or 13 characters for an ISBT 128 unit, VBECS accepts the entry but will not associate it with the scanned unit ID.
The display of Special Testing for a unit is limited to 253 characters in Incoming Shipment. All of the special testing is saved in VBECS and is available on reports.
The display label that is built when logging a blood product into VBECS does not change from volunteer donor to autologous when the donation type is changed.
When blood products intended for quarantine are scanned, VBECS does not automatically quarantine them: the user must do so.
VBECS does not accept scanned information for some modified ISBT 128 blood units. VBECS accepts entered information in these instances.
The user may not scan labels for para-Bombay or Bombay ISBT units.
VBECS does not check whether blood products with the same unit ID and different product codes are the same blood type.
A site may receive duplicate unit IDs. VBECS cannot ensure that these IDs represent products from the same donor.
VBECS reads Codabar prefixes and translates them into digits. When this results in identical Codabar unit IDs, keep the original unit and return the other unit to the supplier.
When a user scans a product code or expiration date barcode in the Unit ID field, VBECS does not display the entered information or an error icon.

Incoming Shipment
When a user scans an ISBT product ID with "0" as the donation type (sixth character), VBECS indicates that the donation type is not specified. When the user selects a donation type from the drop-down menu, VBECS replaces "0" with a recognizable character in the product code text box and stores it as such. When the user rescans the product ID, he must manually replace "0" or search using the unit ID alone, or the unit ID and the first five characters of the product ID.
Blood product code entries in the blood product table include those from ICCBBA for ISBT 128 and from the AABB Codabar table. Product codes created after the release of VBECS cannot be entered until the updated tables are incorporated through patches to VBECS. Refuse acceptance of the unit, if possible. Users may not edit the blood product table.
ISBT 128 tables include blood product codes labeled as not for manufacturing or transfusion. VBECS does not restrict such products from being entered, issued, or transferred (local policy controls their use) and does not accommodate blood modification for these products.

Invalidate Test Results
A user may not inactivate a component task in this option.
This option does not release units from patient assignment. The user must release "assigned" units before invalidating the test results.
To maintain a clear record and system checks, the user must invalidate both original and repeat ABO/Rh tests.
Do not invalidate part of the configured test. Doing this will put the specimen into a state where it cannot be processed; invalidate all tests configured as part of the TAS regardless of which sub test (ABO/Rh, ABS or DAT) is the problem.

Issue Blood Components
VBECS does not check for restricted units (not assigned to the patient) received between selection and issue.
All issue-to locations must be valid entries in the VistA hospital location file, including home transfusion.
If VistA link is not available, the user may proceed with the issue process but the current issue-to hospital location is not recorded. All of the associated report records display "VistALink Down" in the Issue To location field.
ISBT 128 tables include blood product codes labeled as not for manufacturing or transfusion. VBECS does not restrict such products from being entered, issued, or transferred (local policy controls their use), but does not accommodate blood modification for these products.
When a site has configured subservient VistA institutions to the VBECS division, users can only issue units to locations in the primary VistA. Sites may add locations to VistA to account for the relocations to the subservient VistA institutions.
BCE information applies ONLY when a BCE COTS product is available for communication and is configured.
Changing a patient's transfusion requirements (component requirements, antigen negative requirement or antibody identified) after a unit is crossmatched does not provide an alert at Issue Unit regarding the type of crossmatch (electronic or serologic) performed on a unit available for issue.

Justify ABO/Rh Change
VBECS does not allow the selection or issue of certain blood components that may appear ABO/Rh compatible when discrepant ABO/Rh results are not justified.

Local Facilities
A user may:
<ul style="list-style-type: none"> • Activate and edit collection facility records distributed with VBECS. • Add and locally edit active collection facilities. • Activate and/or deactivate local collection facilities. • Edit a division name, which may result in a mismatch between VistA and VBECS division names. • The Local Facility Address maximum length is 45 characters including spaces to allow the address to display on VBECS reports. The third line of the facility address is not displayed on VBECS reports.

Log In Reagents, Update Inventory, View/Print Inventory, Maintain Minimum Levels
VBECS does not display a reagent report by invoice number.
The unsatisfactory portion of a shipment must be recorded individually.
There is no retrospective data entry or user selection for these options.
Each site must set a standard for recording quantity: it may use the number of cases or the number of units within the cases.
When a user enters a lot number in the Lot # field and clicks the magnifying glass , VBECS displays only the first instance of a duplicate lot number.

Log Into VBECS and VistA
VBECS automatically logs a user with access to only one division into that division.

Maintain Specimen
When a search by UID causes VBECS to time out, repeat the search by UID and patient name or by patient name only.
Specimen recollection requests must be made offline.
Designating a specimen as "unacceptable" after accepting a CPRS order closes that order and requires a new CPRS order and sample.

Modify Units
VBECS allows modification of units only in limited, assigned, crossmatched, and available statuses.
VBECS presents the conservative expiration date for target units. VBECS does not know the specific collection details. It is the responsibility of each user to accurately enter and process any one blood unit based on its expiration date/time as labeled by the collection facility.

Modify Units: Pool Units
VBECS allows modification of units only in limited, assigned, crossmatched, and available statuses.
When a user adds a biohazardous unit to a pool, VBECS does not designate the pool as biohazardous on the Blood Transfusion Record Form (BTRF) or Caution Tag.

Modify Units: Split a Unit
A unit previously split in VBECS cannot be split again. A unit pooled in VBECS may not be split.
This option appends alphabetical characters to the Codabar unit ID of split units created from original units. This prevents scanning the unit ID barcode with the split indicator.

Modify Units: Split a Unit

By ISBT 128 convention, split units created from an ISBT 128 original unit will have identical unit ID and five-digit product codes.

Order Alerts

This option does not create a mechanism to accept or reject orders or specimens or define reports for viewing orders.

The user may select the option to display printer alerts only when the designated printer is configured and available.

The printer must be configured by the System Administrator.

Order Reflex Tests

A crossmatch (XM) and a repeat XM cannot be ordered unless they are associated with a specimen with a future expiration date and a phlebotomist is associated with that specimen.

When a user selects weak D for the antigen typing test and the selected patient is Rh positive, VBECS warns that weak D testing cannot be performed. There is no override and VBECS does not add the test.

Outgoing Shipment

A user may edit only unconfirmed invoices.

The number of tests displayed on an invoice is limited to 8,000 characters per unit.

When the user creates an invoice and cancels before saving it, VBECS excludes the invoice number from use.

Patient History Report

When units are restricted for a patient during Incoming Shipment, VBECS does not display them on the Patient History Report.

VBECS displays these default values for transfusion reactions converted from VistA to VBECS. They do not reflect actual interpretations or results:

Converted Transfusion Reaction	Default
Pre-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Post-Transfusion Specimen Checks, Hemolysis Checks OK?	Yes
Symptoms	No symptoms identified.
Implicated units	No units implicated in reaction.
Serologic tests	No serologic tests performed.

Order History Report

VBECS does not display the identity of the CPRS user who discontinues a signed order in CPRS.

Patient Merge

The Traditional Supervisor or above must update patient records manually.

VBECS does not notify users of additional VistA patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor completes the update and clears the icon.

Patient Merge

Patient merges involving patient's with names that exceed the database standard for File #2 (30 characters) or do not have a first or last name in Vista are rejected and do not appear in VBECS. The VBECS System Administrator is notified of each rejection via an email to the address entered in VBECS Administrator "Interface Failure Alert Recipient" field.

Patient Testing: Cancel an Active Order

Cancellation of one order does not affect other orders placed with it.

Before a user may cancel a partially completed task, he must invalidate the results for that task, thereby changing its status to "not started."

Patient Testing: Enter Antibody Identification Results

When an antigen negative requirement is not configured for an antibody, the user may enter the requirement for a patient through Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement.

Patient Testing: General Instructions

This test record is not editable after the user completes a test and saves it to the database. When a correction needs to be filed, a corrected report must be generated, as described in Invalidating Patient Test Results.

User selection of partially completed testing grids with different phases of testing creates a testing discrepancy in VBECS. Batch testing is limited to testing at the same phases.

Patient Testing: Pending Task List

This option is for diagnostic, blood component, and reflex test orders placed in or directed to the user's division.

The PTL does not include daily quality control, equipment maintenance schedules, or pending ABO/Rh confirmation and auxiliary tests on red blood cell units.

Users may perform testing for no more than four patients in one VBECS session.

Patient Testing: Record a Crossmatch

The specimen associated with testing must have an expiration date and time that are later than the testing date and time entered in Patient Testing: Record Patient Test Results. When this is not the case, VBECS warns that the patient specimen expired. There is no override.

This option does not include electronic crossmatch (eXM).

Patient Testing: Record a Direct Antiglobulin Test

VBECS does not check whether the lot and vial numbers selected for the positive control cells are different from those of the negative control cells.

Patient Testing: Record a Patient Antibody Screen

This option is for patient testing only.

Patient Testing: Record a Patient Antigen Typing

Antigen typing applies only to red blood cell antigens. [This option does not address platelet, histocompatibility locus antigen (HLA), or Immunoglobulin A (IgA) antigens.]

VBECS does not check:

- The patient antigen typing with the specificity of any patient's antibodies entered.
- Whether the lot and vial numbers selected for the positive control cells differ from those of the negative control cells.

The patient antigen typing result is not available to CPRS.

Patient Testing: Record a Transfusion Reaction Workup

Adverse reactions to transfusion detected at a later time (disease transmission) are not recorded through this option. Refer to Transfusion Complications Report.

Patient Testing Worklist and Testing Worklist Reports

Reports do not include Transfusion Reaction Workup (TRW) test results.

VBECS does not accommodate online review (verification, signature) of reports.

Post-Transfusion Information

VBECS cannot update a unit record to presumed transfused when a user has that unit record open (locked) in Enter Post Transfusion Information. VBECS will attempt to update the unit the next time the presumed transfused background routing runs if the unit remains in an issued status.

If post-transfusion information is not transmitted to VBECS, the user will need to enter that information to update the record from presumed transfused to transfused.

The transfusion comments and volume transfused entered in the Transfusion Verification application as an update to a completed transfusion are not sent to VBECS.

Transfusion reaction symptoms are not sent to VBECS by the Transfusion Verification application unless a transfusion reaction has also been indicated for that event.

Print Unit Caution Tag & Transfusion Record Form

VBECS cannot reprint a Caution Tag on an expired unit or expired patient order.

VBECS does not print a unit's biohazardous status on the BTRF or Caution Tag.

Prolonged Transfusion Time Report

This report does not include presumed transfused units or units issued to a remote storage location.

Reagents

VBECS does not display a reagent report by invoice number.

The unsatisfactory portion of a shipment must be recorded individually.

There is no retrospective data entry or user selection for these options.

Each site must set a standard for recording quantity: it may use the number of cases or the number of units within the cases.

Reagents
When a user enters a lot number in the Lot # field and clicks the magnifying glass, VBECS displays only the first instance of a duplicate lot number.
Free text details can be added only when "Other" is selected as the canned comment.
Release Units From Patient Assignment
This option does not apply to the release of patient restrictions (autologous and directed).
When a patient, specimen, or blood unit expires, VBECS does not automatically release an assigned blood unit: a user must release the unit.
Remove Final Status
This option applies only to the release of units in final statuses.
Screen Settings
VBECS was designed to operate with the settings described in this section to display complete screens. Failure to maintain these settings may result in partial display of safety-critical information.
Select Units
When a user checks the "D Pos" or "D Neg" check box to modify a search for blood, VBECS returns no units, regardless of the availability of units.
VBECS does not allow a user to select a pool of mixed ABO/Rh units for a patient. Instead, the user must select random units for the patient and pool the mixed ABO/Rh units to allow issue.
VBECS displays a component order as "filled" when the number of units transfused and the number of units ordered are the same. CPRS displays a component order as "completed" when the number of units ready for issue and the number of units ordered are the same. Therefore, CPRS displays the order as "completed" before VBECS displays the order as "filled."
Special Instructions & Transfusion Requirements: Enter a Transfusion Requirement
Inactivation is defined and restricted by security level.
VistA Special Instructions (SIs) from database conversion do not become Transfusion Requirements (TRs): they remain SIs. SIs are information only and are not enforced unless entered in the Component Requirements tab in Special Instructions & Transfusion Requirements.
Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions
VBECS does not enforce Transfusion Requirements (TRs) or antigen negative requirements entered through this option.
SIs moved to VBECS during the database conversion may contain information that must be entered as TRs in VBECS. Inactivate the old SI and create a TR and new updated SI: the database conversion software cannot automatically translate the SI from database conversion into a rule-based TR.
VBECS will truncate SI text longer than 255 characters from the VistA database conversion.

Supplies
The division's modification processes may allow or require the user to select previously entered supplies for association with a blood unit's modification. An inventory of containers (transfer bags) and sterile connection device (SCD) wafers is <i>required</i> for various modification processes.
There is no retrospective data entry or user selection for this option.
The unsatisfactory portion of a shipment must be recorded individually.
Each site must set a standard for recording quantity: it may use the number of cases or the number of units within the cases.

Transfusion Complications Report
The report does not include units in an issued status.

Transfusion Effectiveness Report
The report does not include VBECS diagnostic tests.

Transfusion Requirements Report
Deceased patient records do not appear on the report when they are generated after the patient has been marked as deceased in the VBECS database via a VistA Death Notification Update.

Unit Antigen Typing
Antigen typing applies only to blood products that contain red blood cells. This option does not address platelet, HLAs, or IgA antigenicity.
VBECS does not check whether the lot and vial numbers selected for the positive control cells differ from those of the negative control cells.
D antigen typing must not be combined in batch testing of Rh positive units.

Unit History Report
The report is not multidivisional.
When a user retroactively updates a unit status through Discard or Quarantine, the Unit History Report displays the updated information, but does not display the date the change was made.

Update User Roles
One role at a time may be assigned to a user at a division. The role may be reset.
VBECS does not allow users to edit or delete the sole Administrator/Supervisor.

Updated Patients and Deceased Patients
Demographic and death event updates do not display information for patients without pending or active orders in VBECS. Updates will occur without user input.
VBECS does not notify users of additional VistA patient updates or merges while the icon is displayed. When a user views an update and there is an unprocessed patient merge, the Patient Alert icon remains on the screen until the Traditional Supervisor completes the update and clears the icon.

Updated Patients and Deceased Patients

When VBECS notifies a user of an update to a patient's prefix or suffix, VBECS does not display the updated information: the user must access VistA to view the updated information.

Workload Codes

Transfusion-only sites must not map workload codes to processes such as ABO/Rh testing not performed on site.

This page intentionally left blank.

Appendix E: Known Defects and Anomalies

Copies of *Known Defects and Anomalies* may be obtained from the VA Software Document Library Web site.

This page intentionally left blank.

Appendix F: System Validation

Introduction

Each facility is responsible for validating the VBECS installation, operation, and performance before using VBECS in a live (production) environment. This includes evaluating local intended uses of the software and changing settings to conform to the facility's procedures, processes, and policies. Use internal operating procedures and policies to validate that the software performs in accordance with your site's parameters. Determine the level of validation required for changes to hardware and software.

This section describes only basic system validation. Sample test scripts are available from the National Training & Education Office (NT&EO).

Validation Plan

Write a validation plan with clearly stated goals and strategies to demonstrate that the system hardware, application, new features, and software patches function as expected in your facility's environment. Include descriptions of documents, hardware and software Verification & Validation (V&V) tools, techniques, methods, and operating and test environments to be used. Determine and document when an option is not applicable to your facility and does not require testing.

You may separate this plan into three major sections:

Installation Qualification

The Installation Qualification (IQ) section includes system identification, definition of physical properties, interfaces and peripheral devices, and general functionality. Determine what types of testing must be done. Create test cases in accordance with your internal procedures and intended use.

Apply review, analysis, and testing techniques to determine whether VBECS hardware and peripheral components comply with requirements, including functionality and attributes such as consistency, safety, and security.

Operational Qualification

The Operational Qualification (OQ) section includes testing the software (the primary validation activity). Use your OQ to challenge VBECS and the connection to VistA. Base the number and types of tests on the criticality of the software function. The focus of this validation is to test the limits of what the system allows you to do, not how the system will be used. Each facility will have a different level of validation based on local practices.

Analyze risks to ensure that the scope and degree of validation are appropriate. (A system error that may result in direct and severe patient harm requires more detailed and rigorous testing than one that may have minor consequences.)

Create a testing environment as much like the production environment as possible. Design test cases to:

- Address critical software functions.
- Detect failures and flaws.
- Verify and challenge system functions and limitations, and user documentation.

Verify data from VistA displayed in VBECS and data from VBECS displayed in VistA. (Results will vary based on your VBECS settings.) Evaluate the output of the testing and determine whether the test passed or failed.

Maintain evidence of software testing, as part of the validation record, in accordance with standards established by accrediting agencies associated with the blood bank and laboratory. Include, at minimum, your inputs, outputs, and evidence that tests performed were satisfactory or unsatisfactory. (You may record input manually, make copies of unit labels or patient information, capture screen shots⁷ from VistA and VBECS, print various reports, etc.)

Performance Qualification

The Performance Qualification (PQ) section focuses on the integration of the application and processes into your facility. This begins with the revision of Standard Operation Procedures (SOPs) and policies to include the proper use of the new system. Once the SOPs are in place and validated, train the users. Ensure that training materials and records are complete before allowing users to operate the system. Document the validation, including:

- Personnel participating in validation activities
- Hardware and software change controls
- Testing results
- SOP reviews, training plans, installation and implementation
- Evaluation and acceptance
- Review and approval authority and dates

Prerequisites

Perform these tasks for test and production environments in VBECS before you test the application:

- Install required patches in your VistA test and production environments.
- Install VBECS.
- Complete the database conversion.
- Record system settings.

⁷ To make a copy of the active window, press **Alt+Print Screen**. To copy the entire screen as it appears on your monitor, press **Print Screen**. To paste the image into a file in an Office program, open the file you want to paste the image into, click **Edit**, and then click **Paste**.

Record System Settings

These sections describe how to capture system settings in preparation for validation.

Configure Division

A user defines site parameters that affect how VBECS behaves at a specific division. VBECS may be configured to function differently at each division within a consolidated database. Site parameters may be redefined.

User Roles with Access to This Option

►►►►► Administrator/Supervisor

Preparation for Testing

- 1) Select **Configure Division**.
- 2) Record:
 - ICCBBA registration number: _____
 - Blood Bank Medical Director's name as you want it to appear on reports: _____
 - Latest expiration date (default: three days) for specimens, in accordance with local policies and procedures: _____
- 3) Is your division full service (FS) or transfusion only (TO)? _____
- 4) If your site does not perform on-site testing, continue at Order Alerts.
- 5) If your division is full service, do you want to enable electronic crossmatch? _____
- 6) Print, sign, date, and save an Audit Trail Report for your validation record.

Configure Testing: Full Service

The user sets testing preferences in accordance with a division's medical policies.

User Roles with Access to This Option

►►►►► Administrator/Supervisor

Preparation for Testing

- 1) Set testing preferences.
- 2) Capture and save a screen shot of each piece of information.
- 3) Print, sign, date, and save an Audit Trail Report for your validation record.

Product Modifications: Full Service

The user enables the types of modifications made at a division.

User Roles with Access to This Option

►►►►►► Administrator/Supervisor

Preparation for Testing

Modification processes are disabled by default. You must activate modification types performed at your facility. You may select none, some, or all. Enter a processing cost for each modification type, as applicable. You may update these settings later.

Table 28: Modification Types and Costs

Modification Type	Modification Processes		Processing Cost
	Enabled	Disabled	
Thaw			
Pool			
Thaw/Pool Cryo			
Split/Divide			
Irradiate			
Leukoreduce			
Volume Reduce			
Wash			
Rejuvenate			
Freeze			
Deglycerolize			

- 1) Capture and save a screen shot of each piece of information.
- 2) Print, sign, date, and save an Audit Trail Report for your validation record.

Order Alerts

The user configures diagnostic test or component order alerts that appear at login, and turns them on or off. The alerts notify users when new orders are entered, updated, or canceled in VistA.

User Roles with Access to This Option

►►►►► Administrator/Supervisor

Preparation for Testing

- 1) Capture and save a screen shot of each piece of information.
- 2) Print, sign, date, and save an Audit Trail Report for your validation record.

Configure Daily QC: Full Service

The user sets up the routine reagent racks and QC template for the division before performing daily reagent QC.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

- 1) If you perform daily reagent rack quality control using a commercial kit, select the commercial kit template. If not, select the non-commercial kit template. Template selected: _____
- 2) Daily alert time at which users will be prompted to perform QC: _____
- 3) Do you want to use letters or numbers as rack names? _____
- 4) Maximum number of reagent racks for daily use: _____
- 5) Primary enhancement medium: _____
- 6) Secondary enhancement medium: _____

The template contains the minimum number of tests required to test the reagents. Select additional tests that are part of your daily reagent rack QC, in accordance with local policies and procedures.

- 7) Capture and save a screen shot of each piece of information.
- 8) Print, sign, date, and save an Audit Trail Report for your validation record.

Maintain Minimum Levels: Full Service

The user sets minimum stock levels for reagent types used in the division.

User Roles with Access to This Option

►►►► Traditional Supervisor

Preparation for Testing

You may enter the minimum number of vials for each reagent type that your blood bank stocks to alert staff when it reaches a critical level.

- 1) Capture and save a screen shot of each piece of information.
- 2) Print, sign, date, and save an Audit Trail Report for your validation record.

Log In Reagents: Full Service

A user records the receipt of specific reagents by type. The reagent record includes the date the shipment was received, reagent name, lot number, expiration date, manufacturer, and number of vials received in a shipment.

User Roles with Access to This Option

All users

Preparation for Testing

- 1) Enter reagent types and quantities used. (For QC validation, enter a lot number that was received and found unsatisfactory, and a lot number that expires at 23:59 on the day that you test the reagent rack QC.)
- 2) Capture and save a screen shot of all entries.
- 3) Print, sign, date, and save the View/Print Inventory (Reagent) Report for your validation record.

Blood Products

The user activates, deactivates, and edits a limited number of variables for an existing blood product site parameter. VBECS provides a Codabar and an ISBT 128 Blood Product reference table to facilitate these and other actions.

User Roles with Access to This Option

All users

Preparation for Testing

Activate blood products here or during incoming shipment. It may be helpful to set up your most common products before entering blood inventory in the system.

Healthcare Common Procedure Coding System (HCPCS) codes are provided by your hospital's VistA file. Select the code that is associated with the blood product.

- 1) Use the sample worksheet to record information.
- 2) Capture and save a screen shot of each piece of information before continuing.
- 3) Print, sign, date, and save an Audit Trail Report for your validation record.

Table 29: Sample Blood Products Worksheet

Blood Product		HCPCS	Product Shippers			
Product Code	Short Product Name	HCPCS Code	Shipper Name	FDA Reg. #	Cost (\$)	Return Credit %
04020	RBC ACD-A	P2028	Abington Memorial Hospital	2577690	\$150.00	0

Table 30: Blood Products Worksheet

Blood Product		HCPCS	Product Shippers			
Product Code	Short Product Name	HCPCS Code	Shipper Name	FDA Reg. #	Cost (\$)	Return Credit %

Local Facilities

The user customizes collection facility information.

User Roles with Access to This Option

All users (activate a facility)

►►►► Enhanced Supervisor (add a facility)

Preparation for Testing

- 1) Use the sample worksheet to record information.
 - 2) Capture and save a screen shot of each piece of information.
 - 3) Print, sign, date, and save an Audit Trail Report for your validation record.

Table 31: Sample Local Facilities Worksheet

FDA Reg. No.	Facility Name	ICCBBA Reg. No.	Facility Type		Active Facility?
			Collection Facility?	Testing Facility?	
2577690	Abington Memorial Hospital	W0002	✓		✓

Table 32: Local Facilities Worksheet

Component Classes

The user configures the component class parameters for the selected division.

User Roles with Access to This Option

>>>> Enhanced Supervisor

Preparation for Testing

- 1) Use the sample worksheet to record information.
- 2) Capture and save a screen shot of each piece of information.
- 3) Print, sign, date, and save an Audit Trail Report for your validation record.

Table 33: Component Class Requirements Worksheet

Component Class	Require Current Specimen ABO/Rh Results?	Maximum Transfusion Time ⁸	Laboratory Test Name	Specimen Type	Threshold Result
RED BLOOD CELLS	✓ (uneditable)				
FRESH FROZEN PLASMA					
PLATELETS					
CRYOPRECIPITATE					
WHOLE BLOOD	✓ (uneditable)				
OTHER					

⁸ Default is 270 minutes (4 hours, 30 minutes).

Transfusion Complications

The user defines Laboratory tests and threshold values used to generate a report listing patients transfused within a specified period who may have developed a transfusion-associated disease.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

- 1) Use the sample worksheet to record information.
- 2) Capture and save a screen shot of all entries.
- 3) Print, sign, date, and save an Audit Trail Report for your validation record.

Table 34: Transfusion Complications Worksheet

Laboratory Test Name	Specimen Type	Threshold Result

Transfusion Effectiveness

The user defines Laboratory tests to monitor the effectiveness of transfusions.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

- 1) Use the sample worksheet to record information.
- 2) Capture and save a screen shot of all entries.
- 3) Print, sign, date, and save an Audit Trail Report for your validation record.

Table 35: Transfusion Effectiveness Worksheet

Laboratory Test Name	Specimen Type	Threshold Result

Workload Codes

The user assigns workload Laboratory Management Index Program/National Laboratory Test (LMIP/NLT) tests to processes performed through VBECS. VBECS test names are associated with VistA workload codes, used when compiling and storing data to provide workload-related reports to VistA.

User Roles with Access to This Option

►►►►► Administrator/Supervisor

Preparation for Testing

Print, sign, date, and save an Audit Trail Report for your validation record.

Canned Comments

The user defines canned comments, by category, for a division.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

Print, sign, date, and save an Audit Trail Report for your validation record.

MSBOS

The user defines parameters for comparing blood component orders placed by a clinician with the Maximum Surgical Blood Ordering Schedule (MSBOS).

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

Print, sign, date, and save an Audit Trail Report for your validation record.

Edit Invoice Text

The user may edit the explanatory text above the signature on the last page of a shipping invoice.

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

Print, sign, date, and save an Audit Trail Report for your validation record.

Antibodies

The user may edit certain fields in the supplied antibody table (Table 13). Use Table 36 (extracted from Table 13) to document local decisions on antigen-negative settings.

Table 36: User-Defined Antigen Negative Requirements

Antibody Name		VBECS Clinically Significant?	VBECS Antigen Negative Setting Default	Local comments and decisions on antigen- negative settings.
VBECS	VistA			
Anti-Kp(a)	Anti-Kp(a)	Yes	No	
Anti-Kp(b)	Anti-Kp(b)	Yes	No	
Anti-M	Anti-M	No	No	
Anti-N	Anti-N	No	No	
Anti-U	Anti-U	Yes	No	
Anti-P1	Anti-P1	No	No	
Anti-Le(a)	Anti-Le(a)	No	No	
Anti-Le(b)	Anti-Le(b)	No	No	
Anti-Js(a)	Anti-Js(a)	Yes	No	
Anti-Js(b)	Anti-Js(b)	Yes	No	
Anti-Lu(a)	Anti-Lu(a)	No	No	
Anti-Lu(b)	Anti-Lu(b)	No	No	
Anti-Lu3	Anti-Lu(a)Lu(b)	No	No	
Anti-C(w)	Anti-Cw	Yes	No	
Anti-C(x)	Anti-Cx	Yes	No	
Anti-E(w)	Anti-E(w)	Yes	No	
Anti-V	Anti-V	Yes	No	
Anti-Xg(a)	Anti-Xg(a)	No	No	
Anti-M(g)	Anti-M(g)	No	No	
Anti-Mi(a)	Anti-Mi(a)	No	No	
Anti-P	Anti-P	No	No	
Anti-A1	Anti-A-1	No	No	
Anti-H	Anti-H	No	No	
Anti-Le(ab)	Anti-Le(x)	No	No	
Anti-Di(a)	Anti-Di(a)	No	No	
Anti-Di(b)	Anti-Di(b)	No	No	
Anti-Do(a)	Anti-Do(a)	No	No	
Anti-Do(b)	None	No	No	
Anti-Co(a)	None	No	No	
Anti-Co(b)	None	No	No	
Anti-LW(a)	None	No	No	
Anti-LW(b)	None	No	No	
Anti-Sc1	None	No	No	
Anti-Sc2	None	No	No	
Anti-Sd(a)	None	No	No	
Anti-Au(a)	Anti-Au(a)	No	No	
Anti-Au(b)	None	No	No	
Anti-Ch	None	No	No	
Anti-Rg	None	No	No	
Anti-Cs(a)	None	No	No	
Anti-Yk(a)	None	No	No	
Anti-Kn(a)	None	No	No	
Anti-McC(a)	None	No	No	
Anti-JMH	None	No	No	
Anti-Yt(a)	None	No	No	
Anti-Yt(b)	None	No	No	

User Roles with Access to This Option

►►►► Enhanced Supervisor

Preparation for Testing

- 1) Change the antigen negative requirement setting on those currently not invoked. Change the percent compatible for your donor population.
- 2) Print, sign, date, and save an Audit Trail Report for your validation record.

Update User Roles

The user reviews and changes the user roles and associated security levels that determine which options and functions are available to a user in the division.

User Roles with Access to This Option

►►►►►► Administrator/Supervisor (change user roles)
►►► Lead Technologist (view user roles)

Preparation for Testing

- 1) It is recommended to assign one Lead Technologist on each shift to respond to system overrides.
- 2) Print, sign, date, and save an Audit Trail Report for your validation record.

Testing

Sample test scripts will be provided during the implementation of VBECS. Each site must evaluate the intended uses of the software and make appropriate changes (additions or deletions) to the facility's procedures, processes, policies, and test scripts.

You are required to validate your process and how you use the software. You must also validate the software itself. Training users is an integral part of the validation process.

VBECS Forms and Reports


Print and review these reports and forms at the end of validation testing.

- Administrative Data
- Audit Trail
- Blood Availability
- Cost Accounting
- C:T Ratio
- Division Transfusion
- Division Workload
- Exception
- Finalize/Print TRW
- Inappropriate Transfusion Request
- Issued/Returned
- Maintain Equipment (Print)
- Order History
- Patient History
- Print Unit Caution Tag & Transfusion Record Form
- Prolonged Transfusion Time
- Testing Worklist
- Transfusion Complications
- Transfusion Effectiveness
- Transfusion Reaction Count
- Transfusion Requirements
- Unit History
- View/Print Equipment Inventory
- View/Print Inventory (Reagent)
- View/Print Supply Inventory


This page intentionally left blank.

Appendix G: Work Process Flowcharts


Configure Division Settings: Full-Service Divisions

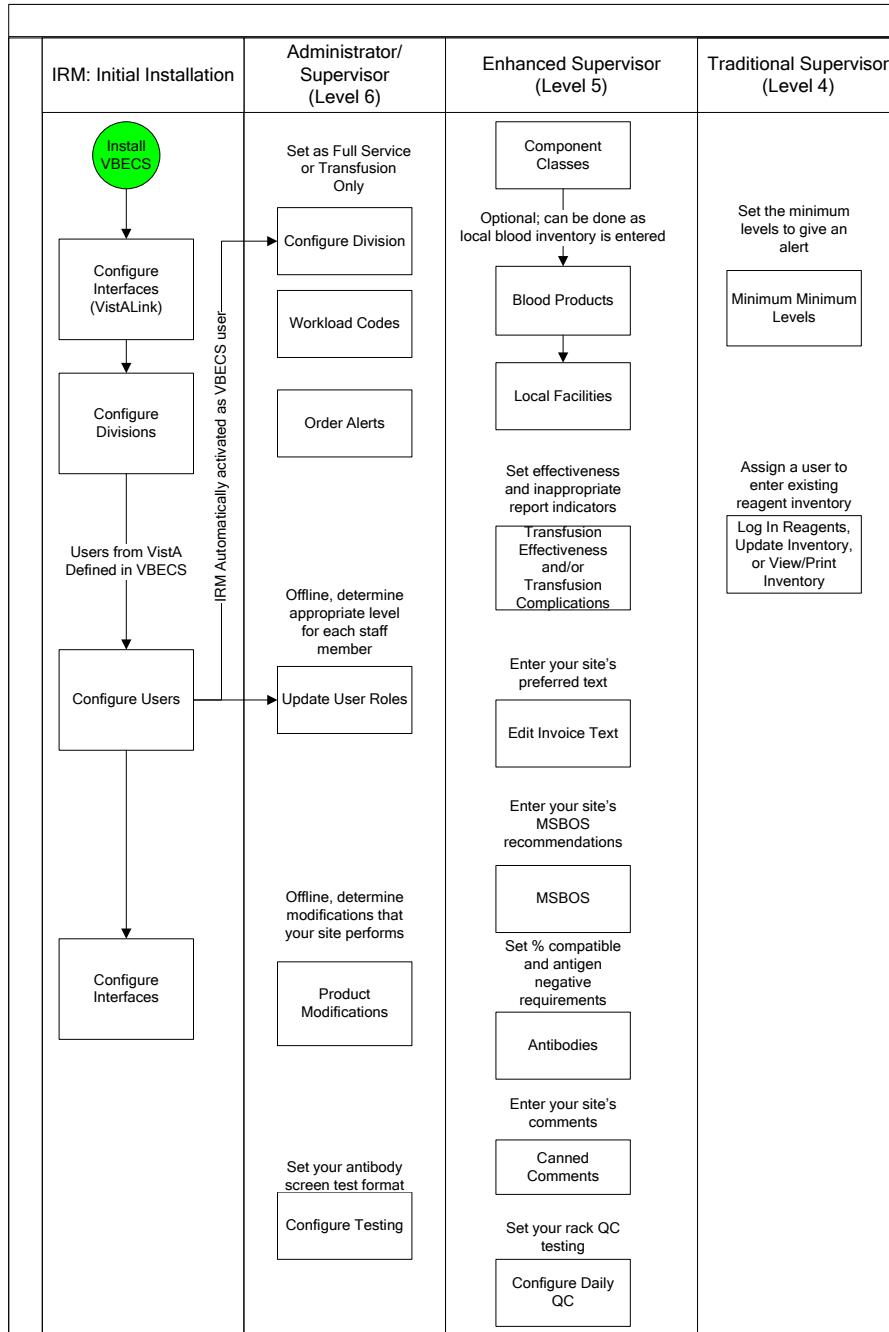


Configure Division Settings: Transfusion-Only Divisions

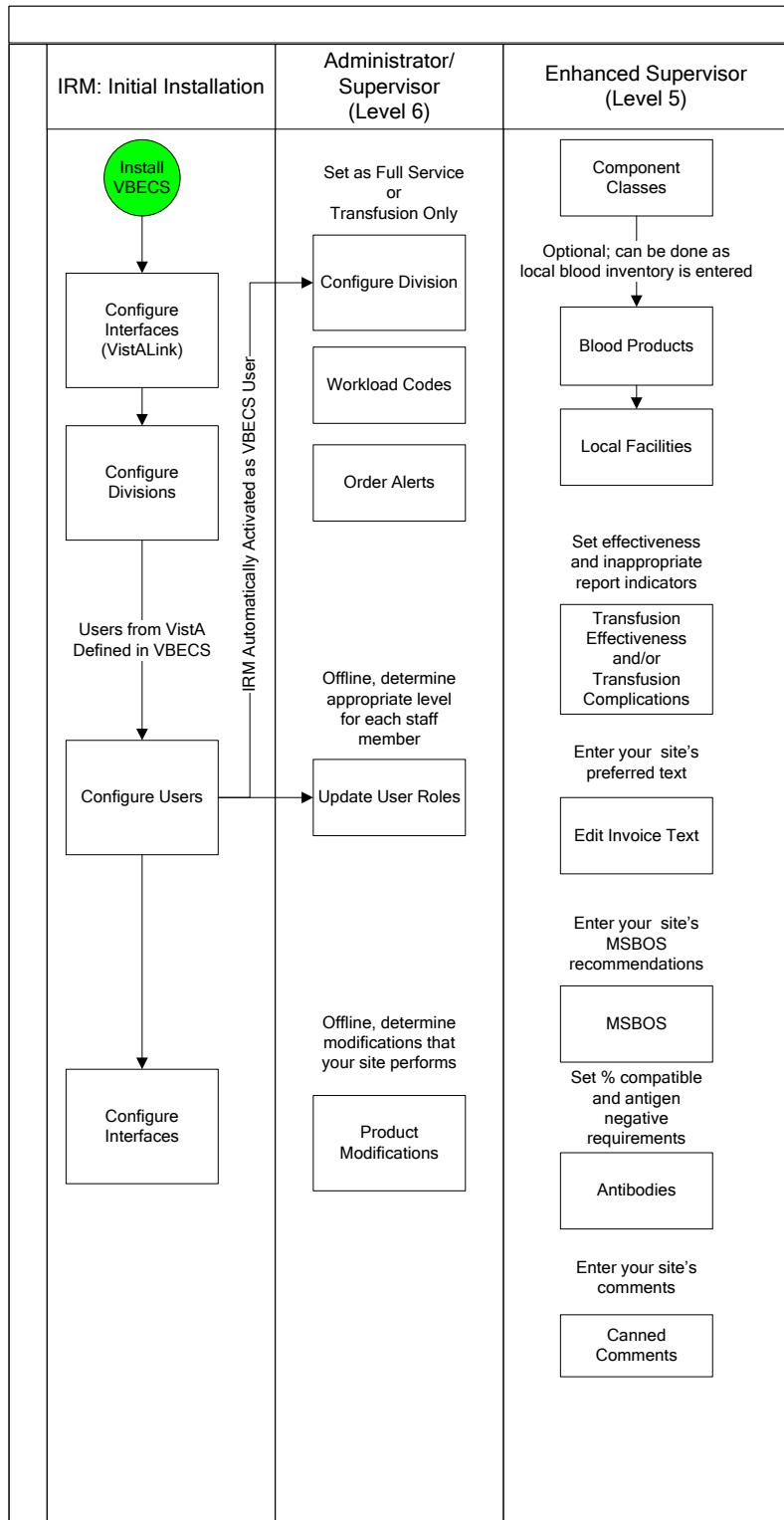


Maintain Reagents and Daily QC: Full-Service Divisions

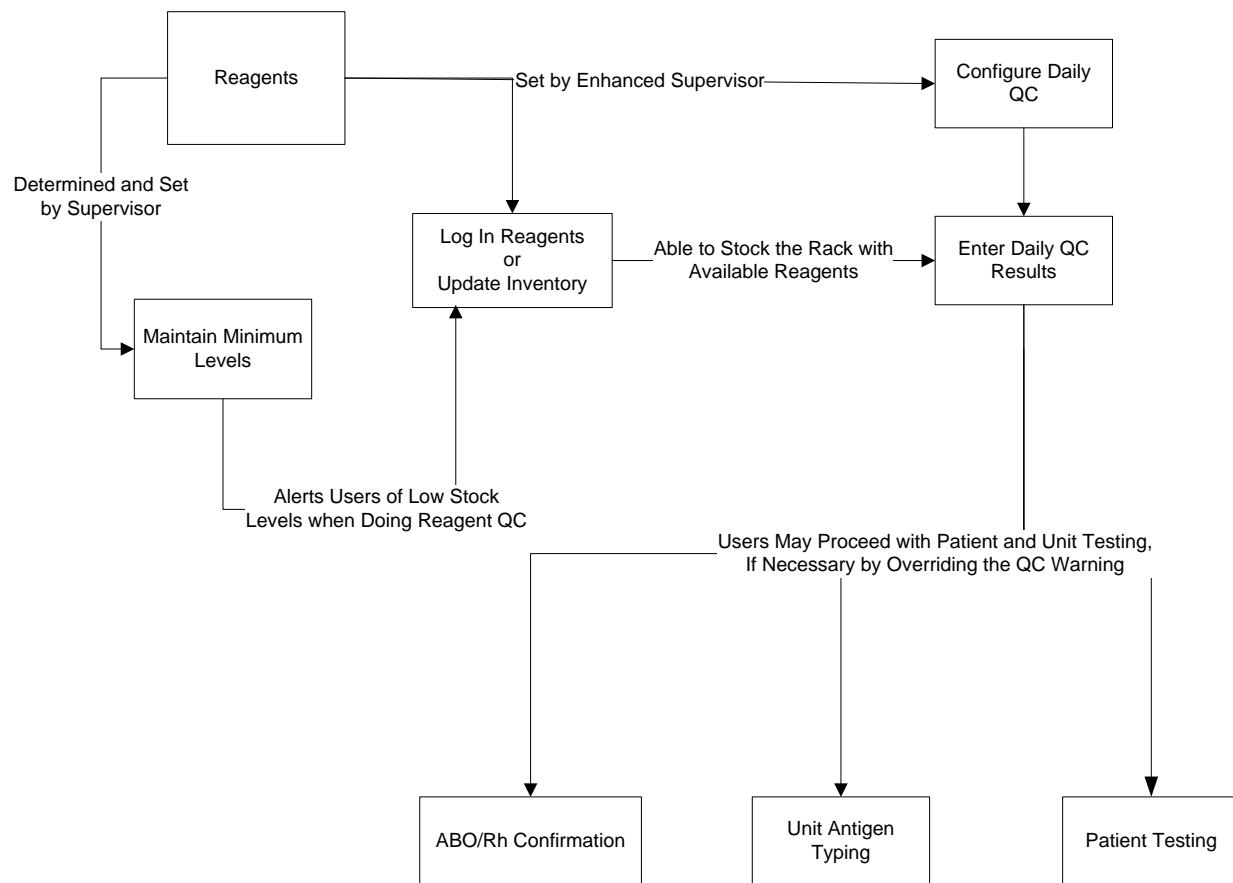


Process Incoming Shipment: Full-Service Divisions

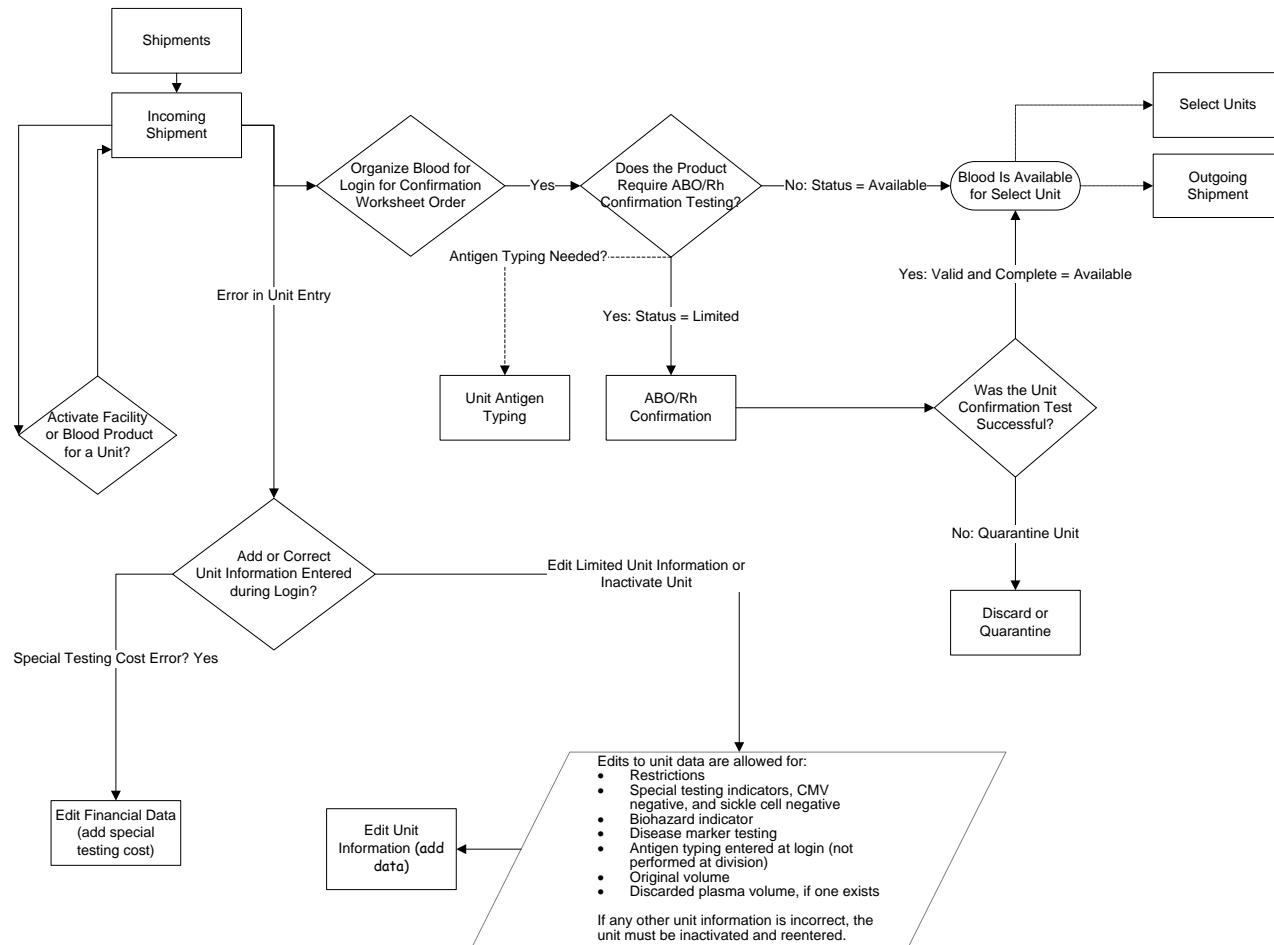


Process Incoming Shipment: Transfusion-Only Divisions

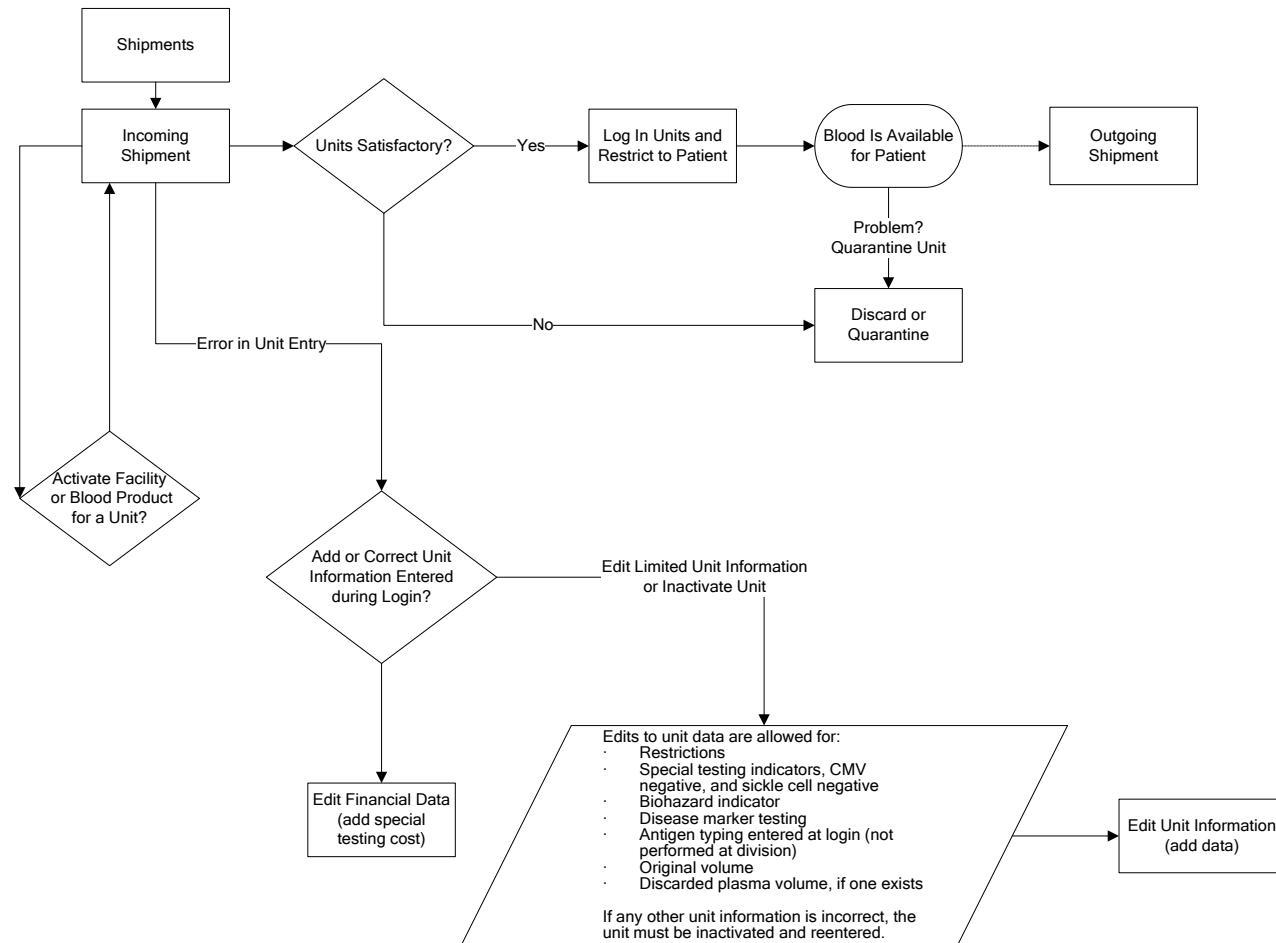


Modify Units

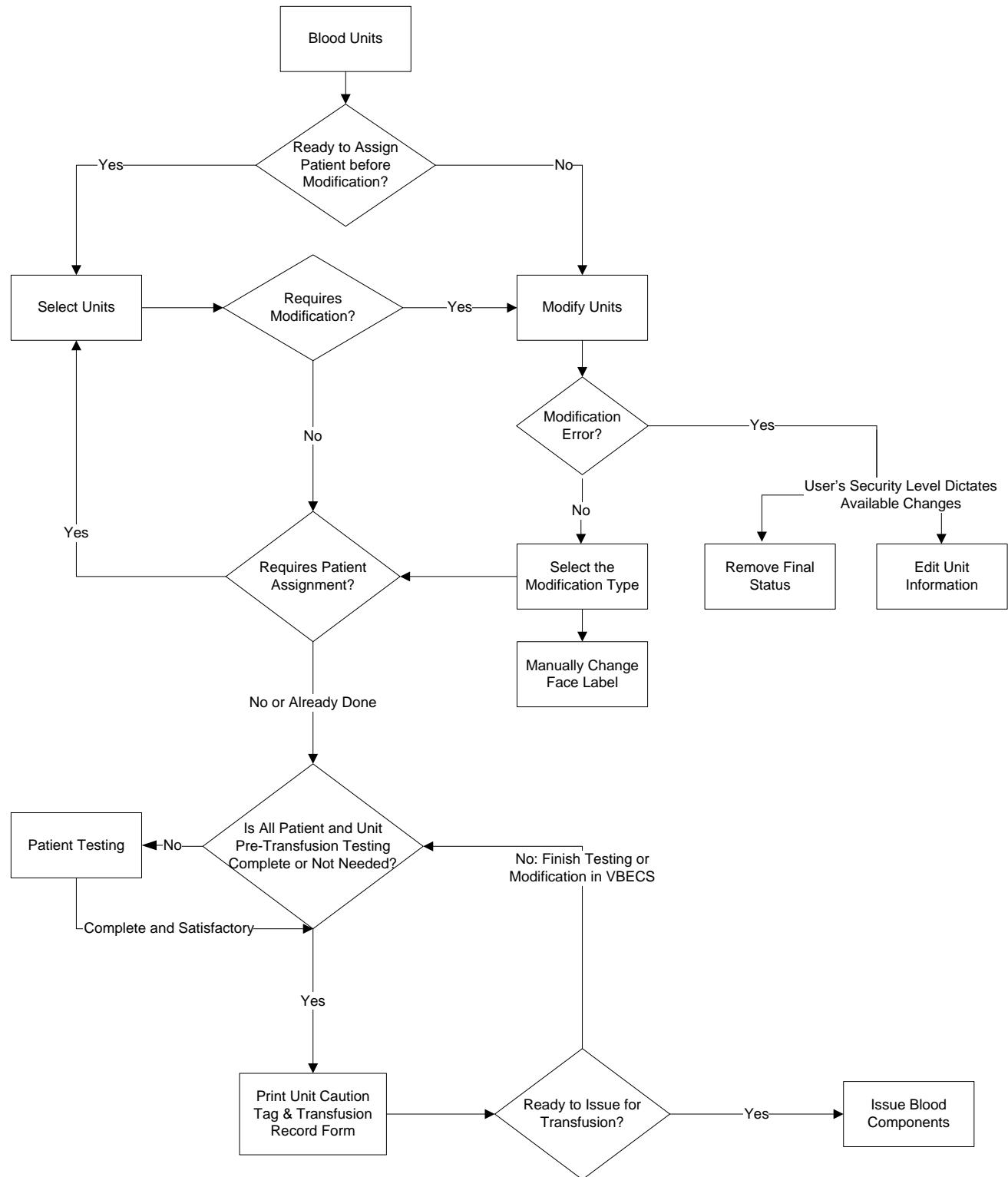


Issue, Return, or Transfuse Units: Full-Service Divisions

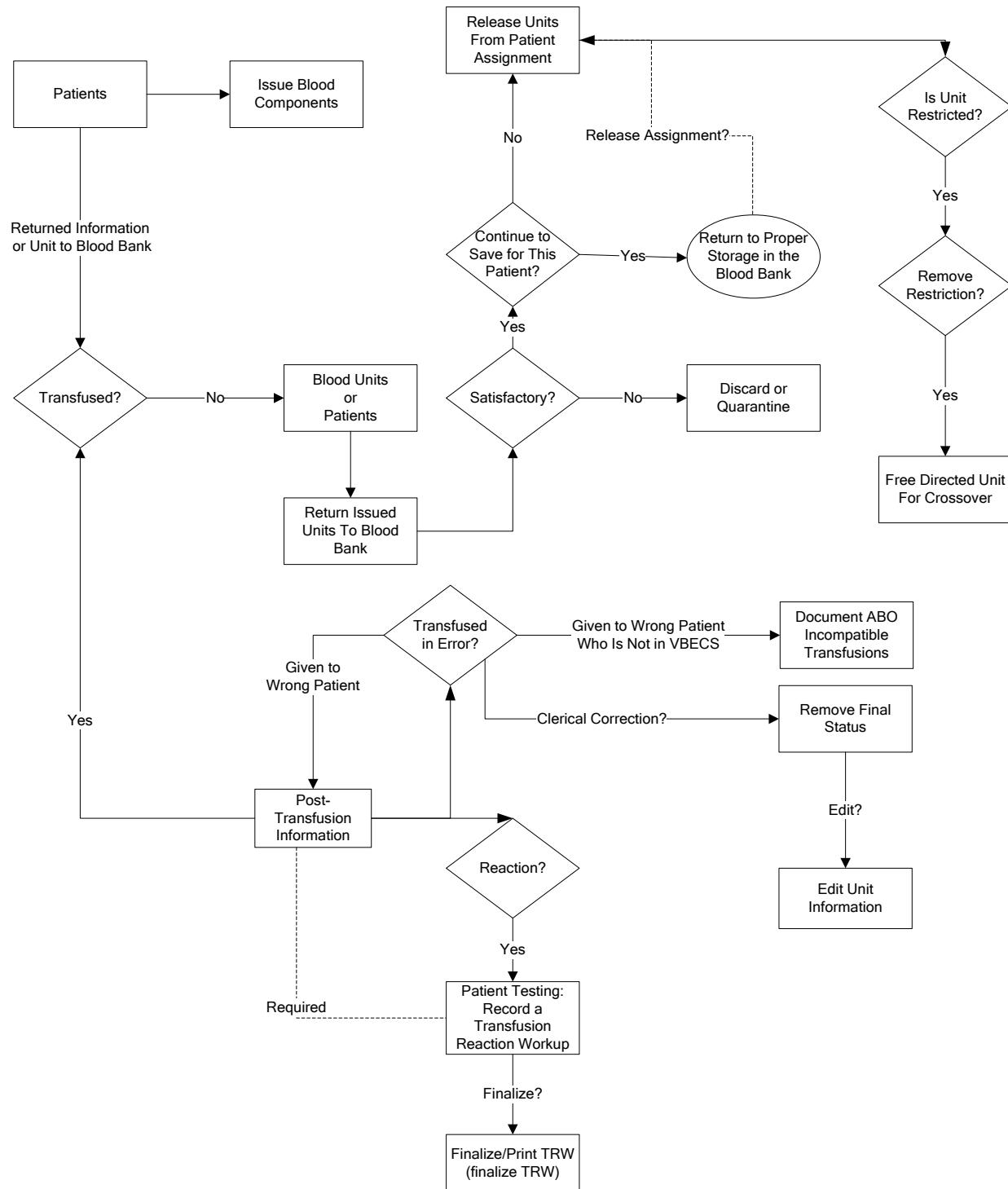


Issue, Return, or Transfuse: Transfusion-Only Divisions

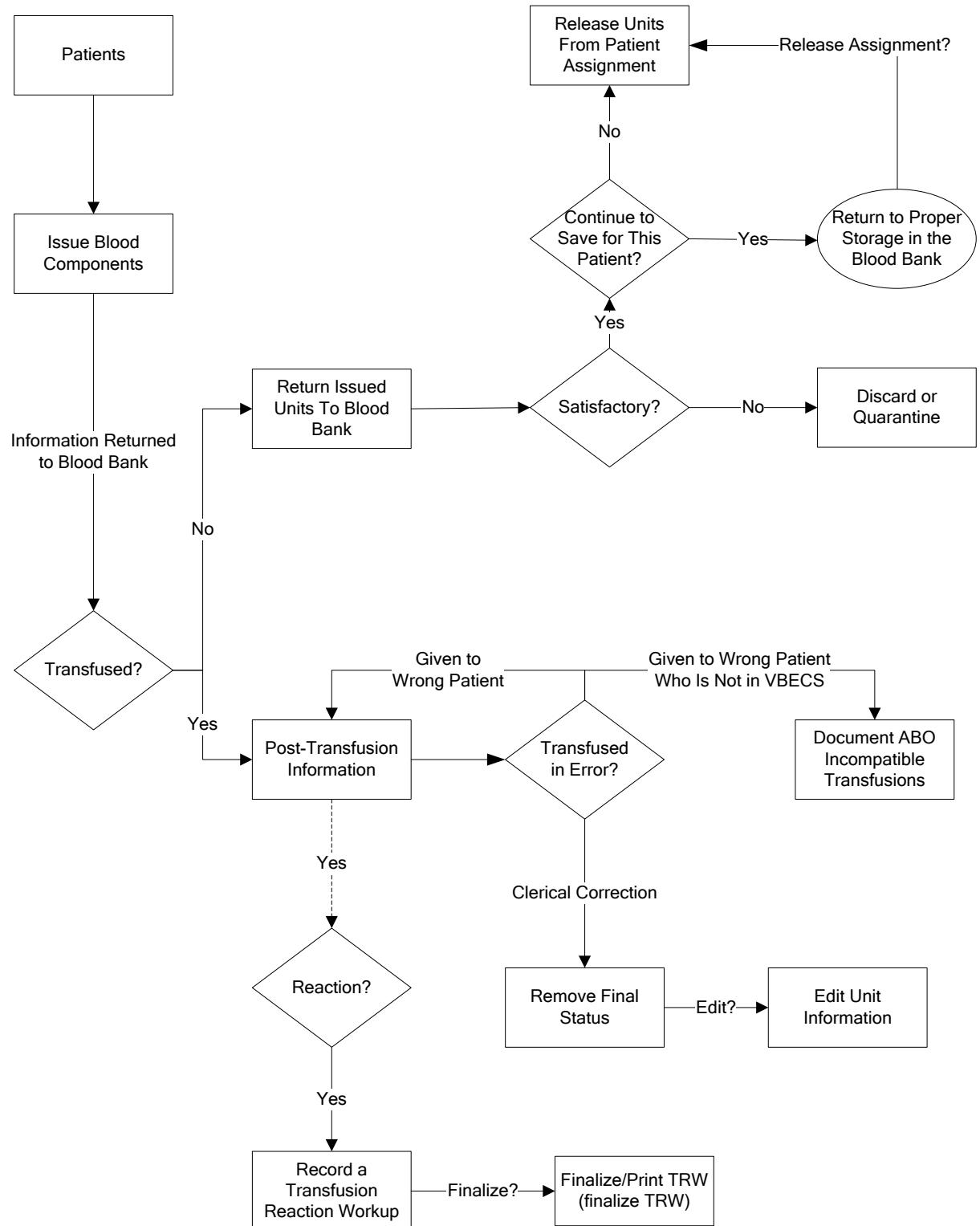

Process a Patient Order

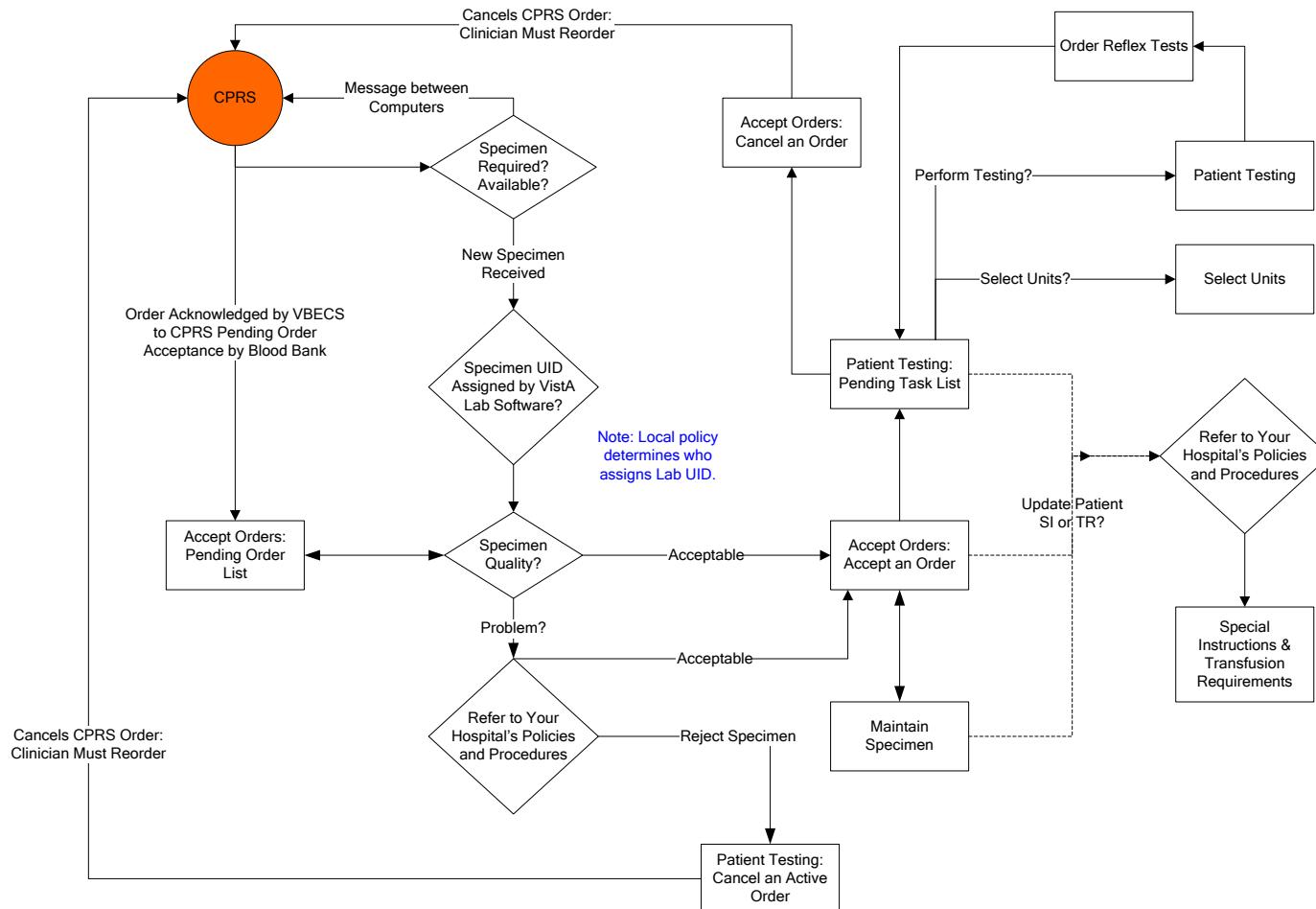

Patient Testing

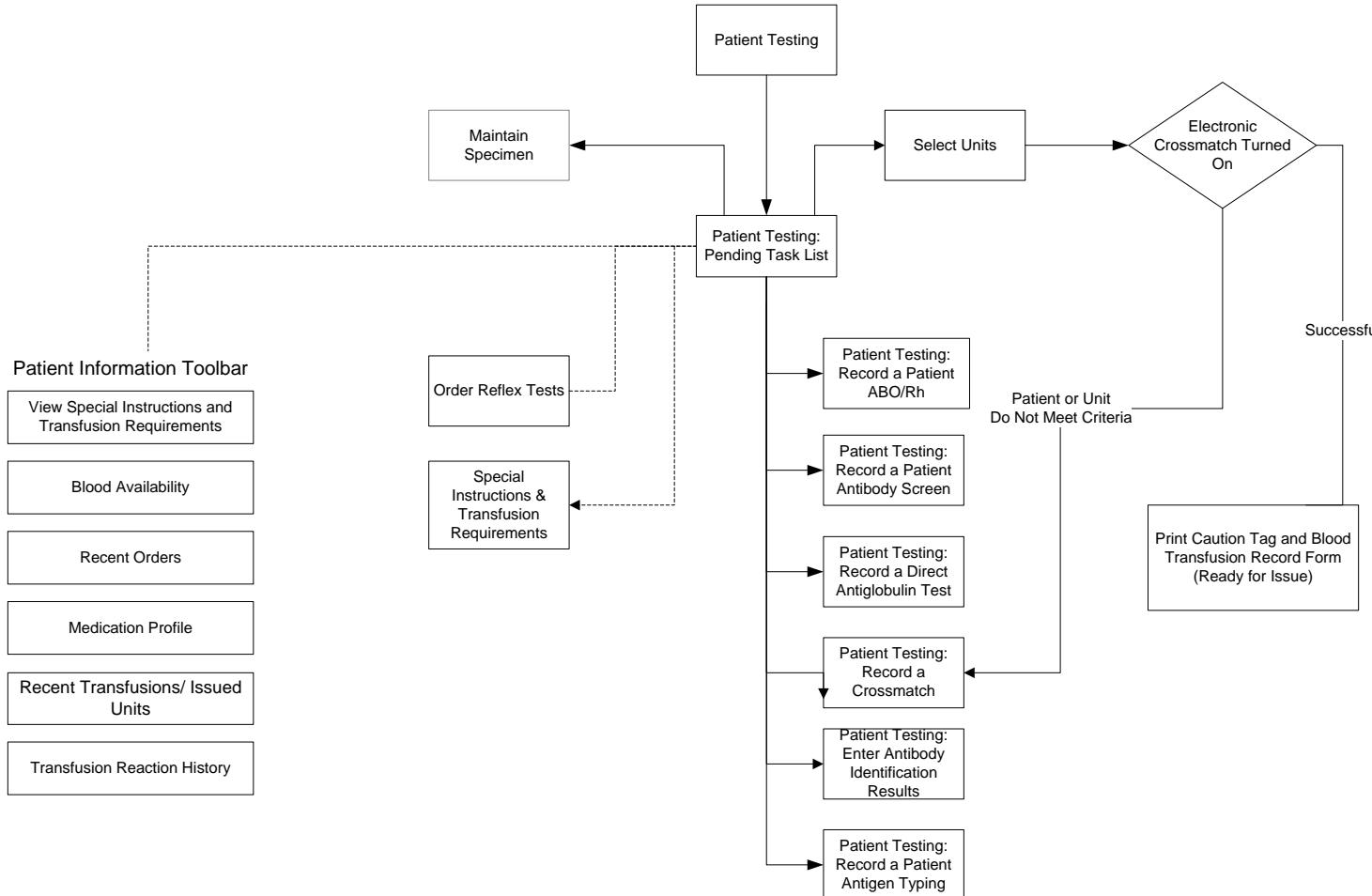

Appendix H: Frequently Asked Questions

Frequently Asked Questions (FAQs) are released as standalone documents. Topics are explained in a question and answer format. FAQs are available at the following SharePoint site:

<http://vaww.oed.portal.va.gov/projects/vbecls/VBEC%20FAQs/Forms/AllItems.aspx>

This page intentionally left blank.

Appendix I: Common Error Corrections

Correct Patient Test Results

Define the problem:

Specifically what data needs correction?

Is it a patient test?

Are any units assigned to the specimen UID associated with this patient order?

The path to transfusion is linear and follows the path from assignment to issue to transfusion. VBECS requires the unit-patient connection to be severed in order to change patient testing used to evaluate compatibility.

All the units must be released from assignment, no matter when the error occurred and needs correction. VBECS requires undoing all of the unit associations so the computer can re-evaluate the corrected information entered and present any warnings or overrides based on that new information.

Information to gather before you start trying to correct any test results:

- Actual dates/times of processing for retrospective data entry
 - Testing Worklist Report: Test result entry information details from the testing grids, reaction results, tech who did the work, any comments they entered
 - Unit History Report: Unit modification, Issue and Transfusion details
 - Patient History Report: Other patient information
 - Order History report: See the search option only for the CPRS Order Number if you have to retrieve an expired order. You do not need to actually open the report.
 - Any offline results that may need to be entered.
 - VBECS CR, HD (Remedy Ticket) and/or CPRS CQ numbers from the Known Anomalies and Defects document related to this problem. Have the user enter this number in the comment details field for all overrides processed.
 - Any local incident number that they may need to enter as a comment.

Invalidate Test Results

Important nuances about invalidating tests:

- When a Type and Screen (TAS) test needs to be invalidated, all the tests in the TAS test configuration must be invalidated to be able to select a test on the Pending Task List (PTL) (It is a loop of a partially completed test).
- When an Antibody Identification (ABID) test is invalidated:
 - The antibody specificities will be removed from the Special Instructions/Transfusion requirements Antibodies Identified at the same time. This also clears any antigen negative requirements that are created by the antibody identified.
 - If separate antigen negative requirements for each of these were entered in SI/TR, they will require update separately.

- If you need to enter corrected results for the ABID, be sure you check the box to add the order for the PTL. If not, you will have to reorder the reflex test.

A finalized Transfusion Reaction Workup (TRW) cannot be invalidated. Invalidation requires a Remedy ticket to have the development team intervene and make that correction. A TRW that has not been finalized can be invalidated and re-entered as any other VBECS test.

Where to start...

1. Are the units transfused?

Go to Supervisor, Remove Final Status to get the unit back to an issued status.

2. Are the units issued?

Go to Return Unit from Issue by either the patient or blood unit menu.

Return the unit from issue at the current date and time.

3. Are the units selected for a patient?

Go to Release Units from Assignment by either the patient or blood unit menu.

When the units are no longer assigned to the patient, go into Invalidate Patient Testing by the Patient Menu and invalidate the patient test.

STOP and make sure the user checks the box in the lower left corner if they want the test to appear on the PTL for correction.

UN-Pooling a Pooled Unit (in the computer only)

This may happen when a pool is created in the computer and the physical pool was not created after all. If there is a patient assignment, issue or transfusion, that must be removed first as described above in the **Where to start** section.

1. Go to Blood Unit, Edit Unit Information and INACTIVATE the Pooled Unit.
2. Go to Supervisor, Remove Final Status, select each original unit removing the MODIFIED status and making each available.
3. Repeat Step 2 for each original unit included in the pool.

UN-Modifying a Unit

Go to Supervisor, Remove Final Status.

If there is a patient assignment, issue or transfusion, that must be removed first as described above in the **Where to start** section.

Appendix J: Recommended Report Usage

Table 39 lists reports, how often to run them, and how to use them. These are recommendations for meeting regulatory requirements. Reports, including the non-individual Blood Availability Report, are available through the Reports menu item. The patient-specific Blood Availability Report is available through the Patients menu item.

Table 37: Recommended Report Usage

Run This Report:	Report Option	Purpose	To Meet Regulatory Requirements:
Daily	Audit Trail	Quality assurance (QA) report details data changes in system.	Print, review, sign, date, and save hard copy.
	Blood Availability	Inventory report	This report may not be needed by transfusion-only or other sites.
	Exception Report	QA report details overrides of system rules processed by users.	Print, review, sign, date, and save hard copy.
	Testing Worklist Report	QA report displays all testing within a division.	Print Patient Testing, Unit testing, Rack QC and miscellaneous reagent QC. Print, review, sign, date, and save hard copy.
Daily or Weekly	Issued/Returned Report	QA report displays all relocation transactions.	Review for anomalies and take corrective action.
Daily, Weekly, and Monthly	Transfusion Requirements Report	A compendium of patient blood types, SIs, and TRs are printed for use during downtime.	<ul style="list-style-type: none"> • Days 1–6: print report daily. • Day 7: print report for previous week; discard daily reports. • Start of month: print report for previous month; discard weekly reports. • Start of year: print report for previous year; discard monthly reports. <p>Generate cumulative reports starting with the date of first data entry through the current date to include all patient records. Store reports so that they are easily accessible.</p>
Weekly	Prolonged Transfusion Time Report	Details units that were not infused within the required time frame.	Generate these reports for blood utilization committee and other quality management activities.
	Transfusion Effectiveness Report	Compares transfusions with user-configured lab data indicators to monitor appropriate blood usage.	
Weekly or Monthly	C:T Ratio Report	Details crossmatch-to-transfusion ratio by physician or treating specialty.	
	Transfusion Complications Report	Compares transfusions with user configured lab data indicators to monitor adverse affects from transfusion.	

Run This Report:	Report Option	Purpose	To Meet Regulatory Requirements:
	Inappropriate Transfusion Request Report	Details requests for blood products that did not meet the facility's transfusion criteria.	
Monthly	Administrative Data Report	Management report details products received, transfused, and returned.	Generate these reports for administrative purposes.
	Cost Accounting Report	Management report of expenditure by invoice number.	
	Division Transfusion Report	Details of units transfused during a selected period.	
	Division Workload Report	Totals workload types for one or all divisions.	Generate this report for blood utilization committee and other quality management activities.
	Transfusion Reaction Count Report	Details transfusion reactions that are ordered, pending, completed, or finalized for one or all divisions.	
As needed	Finalize/Print TRW	Patient Transfusion Reaction Report details a transfusion reaction episode.	Print, review, sign, date, and save hard copy.

Appendix K: Barcode Scanner Configuration and Troubleshooting

Make sure the cord is plugged into the scanner and the computer. Point the scanner at a piece of white paper and squeeze the trigger to verify that the laser works.

For correct barcode scanning results, the Caps Lock key must not be engaged.



Never look directly at the laser. Never direct the laser beam toward yourself or another person. Avoid mirror-like surfaces such as glass, metal, watches, jewelry, and other highly reflective materials when directing the laser beam.

To configure a barcode scanner, scan the barcode in Figure 140. Repeat this process if scanned data are incorrectly displayed on the terminal. Consult the product reference guide or call the VA Service Desk for further assistance.

If your site uses barcode symbology for specimens not covered in these instructions, use the configuration barcodes in the manual included with the scanner.

Figure 139: Barcode Scanner Configuration



This page intentionally left blank.

Appendix L: United Blood Services (UBS) Facility Barcodes

The FDA Registration Number barcode on blood unit bags is specific to Scottsdale. When a facility receives blood units from a UBS blood center other than Scottsdale, scan the corresponding barcode in the FDA Registration Number field to record the appropriate information.

Albuquerque (1671462)



Las Vegas (2971457)



Billings (1771461)



Bismarck (1776389)



El Paso (1671451)



Cheyenne (1771470)



Ft. Smith (2371453)



Lafayette (2371456)



Lubbock (1671458)



McAllen (1671454)



Meridian (1071460)



Rapid City (1771468)



Reno (2971465)



Tupelo (1044534)



Ventura (2077574)



Appendix M: System Responses to Active Transfusion Requirements in Select Unit and Issue Blood Component

The following series of tables, Figures 141 through 152, detail VBECS response to an active transfusion requirement, component class of the unit selected and the user's VBECS role at unit selection and blood component issue.

- Figure 140: Irradiate Cellular Products Only and Select Blood Unit
- Figure 141: Irradiate Cellular Products Only and Issue Blood Component
- Figure 142: Leuko-reduce Cellular Products and Select Blood Unit
- Figure 143: Leuko-reduce Cellular Products and Issue Blood Component
- Figure 144: Washed Red Blood Cell (RBC) Products and Select Blood Unit
- Figure 145: Washed Red Blood Cell (RBC) Products and Issue Blood Component
- Figure 146: Washed PLATELET (PLT) Products and Select Blood Unit
- Figure 147: Washed PLATELET (PLT) Products and Issue Blood Component
- Figure 148: Sickle Cell Negative RBC Products and Select Blood Unit
- Figure 149: Sickle Cell Negative RBC Products and Issue Blood Component
- Figure 150: CMV Negative Cellular Products* and Select Blood Unit
- Figure 151: CMV Negative Cellular Products* and Issue Blood Component

Figure 140: Irradiate Cellular Products Only and Select Blood Unit

Select Unit: IRRADIATED				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
Irradiate Cellular Products Only	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Irradiated	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Not Irradiated	ALL VBECS Users	<p>The Transfusion Requirement is not satisfied; the system emits an audible alert and displays a warning message: <i>This unit does not currently satisfy all of the patient's Transfusion Requirements. Further processing is needed. Continue to select this unit?</i></p> <p>If the user clicks "No," the screen is cleared.</p> <p>If the user clicks "Yes," selection can continue to select the unit.</p>

Figure 141: Irradiate Cellular Products Only and Issue Blood Component

Issue Blood Component: IRRADIATED				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
Irradiate Cellular Products Only	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Irradiated	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Not Irradiated	Technologist Enhanced Technologist	The system displays the warning message: <i>Patient's Transfusion Requirement is not satisfied. You do not have the proper security to issue this unit. Further processing is needed and must be documented in VBECS.</i> The tech may not issue this unit in the computer.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Not Irradiated	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's Irradiate Cellular Products Only Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .

Figure 142: Leuko-reduce Cellular Products and Select Blood Unit

Select Unit: Leuko-reduce Cellular Products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
Leuko-reduce Cellular Products	RED BLOOD CELLS WHOLE BLOOD PLATELETS	Leuko-reduced	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS	Not Leuko-reduced	ALL VBECS Users	The Transfusion Requirement is not satisfied; the system emits an audible alert and displays a warning message: <i>This unit does not currently satisfy all of the patient's Transfusion Requirements. Further processing is needed. Continue to select this unit?</i> If the user clicks "No," the screen is cleared. If the user clicks "Yes," selection can continue to select the unit.

Figure 143: Leuko-reduce Cellular Products and Issue Blood Component

Issue Blood Component: Leuko-reduce Cellular Products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	Selected User Role	Expected System Response
Leuko-reduce Cellular Products	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Leuko-reduced	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS	Not Leuko-reduced	Technologist Enhanced Technologist	The system displays the warning message: <i>Patient's Transfusion Requirement is not satisfied. You do not have the proper security to issue this unit. Further processing is needed and must be documented in VBECS.</i> The tech may not issue this unit in the computer.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS	Not Leuko-reduced	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's Leuko-reduce Cellular Products Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .

Figure 144: Washed Red Blood Cell (RBC) Products and Select Blood Unit

Select Unit: Washed Red Blood Cell (RBC) products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
Washed RBC products	RED BLOOD CELLS WHOLE BLOOD	Washed	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD	Not Washed	ALL VBECS Users	The Transfusion Requirement is not satisfied; the system emits an audible alert and displays a warning message: <i>This unit does not currently satisfy all of the patient's Transfusion Requirements. Further processing is needed. Continue to select this unit?</i> If the user clicks "No," the screen is cleared. If the user clicks "Yes," selection can continue to select the unit.

Figure 145: Washed Red Blood Cell (RBC) Products and Issue Blood Component

Issue Blood Component: Washed RBC products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	Selected User Role	Expected System Response
Washed RBC products	RED BLOOD CELLS WHOLE BLOOD	Washed	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD	Not Washed	Technologist Enhanced Technologist	The system displays the warning message: <i>Patient's Transfusion Requirement is not satisfied. You do not have the proper security to issue this unit. Further processing is needed and must be documented in VBECS.</i> The tech may not issue this unit in the computer.
	RED BLOOD CELLS WHOLE BLOOD	Not Washed	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's Washed RBC Products Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .

Figure 146: Washed PLATELET (PLT) Products and Select Blood Unit

Select Unit: Washed PLT products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
Washed PLT products	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Washed	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	Not Washed	ALL VBECS Users	The Transfusion Requirement is not satisfied; the system emits an audible alert and displays a warning message: <i>This unit does not currently satisfy all of the patient's Transfusion Requirements. Further processing is needed. Continue to select this unit?</i> If the user clicks "No," the screen is cleared. If the user clicks "Yes," selection can continue to select the unit.

Figure 147: Washed PLATELET (PLT) Products and Issue Blood Component

Issue Blood Component: Washed PLT products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Platelet Unit is	Selected User Role	Expected System Response
Washed PLT products	PLATELETS	Washed	ALL VBECS Users	No warning message will display.
	PLATELETS	Not Washed	Technologist Enhanced Technologist	The system displays the warning message: <i>Patient's Transfusion Requirement is not satisfied. You do not have the proper security to issue this unit. Further processing is needed and must be documented in VBECS.</i> The tech may not issue this unit in the computer.
	PLATELETS	Not Washed	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's Washed PLT products Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .

Figure 148: Sickle Cell Negative RBC Products and Select Blood Unit

Select Unit: Sickle Cell Negative RBC products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
Sickle Cell Negative RBC products	RED BLOOD CELLS WHOLE BLOOD	Sickle Cell Negative	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD	Not Sickle Cell Negative	Technologist Enhanced Technologist	The Transfusion Requirement is not satisfied; the system emits an audible alert and displays a warning message: <i>This unit does not currently satisfy all of the patient's Transfusion Requirements. Further processing is needed. Continue to select this unit?</i> If the user clicks "No," the screen is cleared. If the user clicks "Yes," selection can continue to select the unit.
	RED BLOOD CELLS WHOLE BLOOD	Not Sickle Cell Negative	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The following warning(s) is informational: <i>This unit does not currently satisfy the patient's Sickle Cell Negative Transfusion Requirements. Further processing is needed. Continue to select this unit?</i> <i>Do you wish to proceed with the selection of unit <insert unit ID number>?:?</i> If "No," is selected, the screen is cleared and a new blood unit must be selected. If "Yes," is selected, a comment is required and the details of the transaction are captured for later inclusion in an exception report entry: <i>Transfusion Requirement incompatible unit selected.</i>

Figure 149: Sickle Cell Negative RBC Products and Issue Blood Component

Issue Blood Component: Sickle Cell Negative RBC Products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	Selected User Role	Expected System Response
Sickle Cell Negative RBC products	RED BLOOD CELLS WHOLE BLOOD	Sickle Cell negative	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD	NOT Sickle Cell Negative	Technologist Enhanced Technologist	The system displays the warning message: <i>Patient's Transfusion Requirement is not satisfied. You do not have the proper security to issue this unit. Further processing is needed and must be documented in VBECS.</i> The tech may not issue this unit in the computer.
	RED BLOOD CELLS WHOLE BLOOD	NOT Sickle Cell Negative	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's Sickle Cell Negative RBC products Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .

Figure 150: CMV Negative Cellular Products* and Select Blood Unit

Select Unit: CMV Negative Cellular Products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	User Role	VBECS Response
CMV negative cellular products*	RED BLOOD CELLS	CMV Negative*	ALL VBECS Users	No warning message will display.
	WHOLE BLOOD			
	PLATELETS			
	OTHER			
	RED BLOOD CELLS	Not CMV negative	Technologist Enhanced Technologist	The Transfusion Requirement is not satisfied; the system emits an audible alert and displays a warning message: <i>This unit does not currently satisfy all of the patient's Transfusion Requirements. Further processing is needed. Continue to select this unit?</i> If the user clicks "No," the screen is cleared. If the user clicks "Yes," selection can continue to select the unit.
	WHOLE BLOOD			
	PLATELETS			
	OTHER*		Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's CMV Negative cellular products Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .
	RED BLOOD CELLS			
	WHOLE BLOOD			
	PLATELETS			
	OTHER*			

*VistA test results are not evaluated by VBECS. CMV Negative status is recorded during receipt at Incoming Shipment or added to the unit record in Edit Unit Information. All product types in OTHER are evaluated for the requirement including those that are non-cellular.

Figure 151: CMV Negative Cellular Products* and Issue Blood Component

Issue Blood Component: CMV Negative Cellular Products				
Patient Active Transfusion Requirement	Ordered Component Class*	Selected Unit is	Selected User Role	Expected System Response
CMV negative cellular products*	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER	CMV Negative*	ALL VBECS Users	No warning message will display.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER*	Not CMV negative	Technologist Enhanced Technologist	The system displays the warning message: <i>Patient's Transfusion Requirement is not satisfied. You do not have the proper security to issue this unit. Further processing is needed and must be documented in VBECS.</i> The tech may not issue this unit in the computer.
	RED BLOOD CELLS WHOLE BLOOD PLATELETS OTHER*	Not CMV negative	Lead Technologist Traditional Supervisor Enhanced Supervisor Administrator/Supervisor	The system displays the warning message: <i>Patient's CMV Negative Cellular Products Transfusion Requirement is not satisfied</i> and prompts the user whether to proceed. If the user responds "No", the system will clear the selected unit. If the user responds "Yes", the system will require a comment (free text or canned, "unit issue" context) and capture the transaction details for inclusion in an Exception Report entry: <i>Unit issued with unsatisfied Transfusion Requirement</i> .

*VistA test results are not evaluated by VBECS. CMV Negative status is recorded during receipt at Incoming Shipment or added to the unit record in Edit Unit Information. All product types in OTHER are evaluated for the requirement including those that are non-cellular.

This page intentionally left blank.

Appendix N: Workload Process Mapping to Application Option Table

Table 41 associates record saves with workload processes. The data fields identified for transmission at the completion of a Workload Event are based on current VistA workload-related files and fields. VBECS will transmit information to a new flat file. There are no donor workload types in VBECS.

Table 38: Workload Process Mapping to Application Option

Record Save Option	VBECS Process	Transaction Type [P (Patient), U (Unit), M (Miscellaneous)]	Explanation
Record a Transfusion Reaction Workup	ABO Forward and reverse typing (patient)	P	An ABO/Rh test for "pre" or "post" is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for "Pre" and "Post" entries.
Record Patient ABO/Rh		P	Accrue workload when a CPRS-ordered ABO/Rh test is performed.
Invalidate Patient Test Results*		P	Accrue workload when a completed ABO/Rh test is invalidated.
Record Patient ABO/Rh	ABO Forward and reverse typing (patient) Repeat Test	M	Accrue workload when a reflex or repeat ABO/Rh test is performed, completed, and saved.
Invalidate Patient Test Results*		M	Accrue workload when a reflex or repeat ABO/Rh test is invalidated.
ABO/Rh Confirmation	ABO forward typing (unit)	U	An ABO confirmation test is performed. When multiple units are selected in a batch, each unit in the batch accrues a workload event. Note: Workload generated during Anti-D testing is not included in the unit's confirmation test. Workload is not accrued when an ABO or Rh discrepancy override is processed and VBECS releases all patient assignments. Workload is not accrued when VBECS quarantines the unit due to a discrepancy. There is no special handling for workload collection for additional confirmation tests on a unit.
Edit Unit Information*		U	Accrue workload when an ABO confirmation test is invalidated.
ABO/Rh Confirmation	ABO/Rh forward typing (unit)	U	An ABO/Rh confirmation test is performed. When multiple units are selected in a batch, each unit in the batch accrues a workload event. Note: Workload generated during Anti-D is part of the unit's confirmation test. Workload is not accrued when an ABO or Rh discrepancy override is processed and VBECS releases all patient assignments. Workload is not accrued when VBECS quarantines the unit due to a discrepancy. Any unit successfully confirmed accrues workload. For split modifications: workload is not inherited by split units. A split unit that requires confirmation accrues confirmation workload at the time of testing. There is no special handling for workload collection for additional confirmation tests on a unit.
Edit Unit Information*		U	Accrue workload when an ABO/Rh confirmation test is invalidated.

Accept Order	Accept Order	M	Accrue workload when an order is accepted. When a multiple orders are selected, each order accrues workload.
Enter Antibody Identification Results	Antibody identification Work-Up	P	User enters additional workload associated with the individual reflex-ordered ABID. The selected VBECS multiplier will multiply against the Vista multiplier and display the (multiplication) product on the Division Workload Report.
Invalidate Patient Test Results*		P	Accrue workload when the ABID is invalidated.
Record a Transfusion Reaction Workup	Antibody Screen (patient)	P	An ABS test for "pre" or "post" is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for "Pre" and "Post" entries.
Record Patient Antibody Screen		P	Accrue workload when an ordered ABS test is performed.
Invalidate Patient Test Results*		P	Accrue workload when a completed ABS test is invalidated.
Record Patient Antibody Screen	Antibody Screen (patient) Repeat Test	M	Accrue workload when a reflex or repeat ABS test is performed, completed, and saved.
Invalidate Patient Test Results*		M	Accrue workload when a reflex or repeat ABS test is invalidated.
Unit Antigen Typing / Patient Antigen Typing	Antigen phenotyping, Single Test phase (QC)	M	Accrue workload when Antiserum QC in Unit or Patient Antigen Typing includes the testing of both the positive and negative control cells, per specificity by lot number, when only one phase of reactivity is chosen for the test grid (IS or AHG). One workload event is collected per completed tab for regular or repeat antigen tests.
Unit Antigen Typing / Patient Antigen Typing	Antigen phenotyping, Multiple Test phases (QC)	M	Accrue workload when Antiserum QC in Unit or Patient Antigen Typing includes the testing of both the positive and negative control cells, per specificity by lot number, when only multiple phases of reactivity are chosen for the test grid, IS/RT, RT/37, or weak D. One workload event is collected per completed tab for regular or repeat antigen tests. When weak D is the selected test, QC may not be accrued for the rack selection. QC is accrued when positive and negative cells must be tested for the lot number.
Cancel Pending Order	Cancel Order	M	Accrue workload when an order on the pending order list is canceled. When multiple orders are canceled, each order accrues workload.
Cancel Active Order	Cancel Order	M	Accrue workload when an order on the pending task list is canceled. When multiple orders are canceled, each order accrues workload.
Select Units for Crossmatch	Crossmatch unit, electronic	P	This process is invoked when an individual unit is selected for patient assignment and the unit is electronically crossmatched. When multiple units are selected, each unit accrues workload.
Enter Crossmatch Results	Crossmatch unit, serologic immediate spin	P	Accrue workload when an individual unit crossmatch is selected to include only the IS phase, is completed, and is saved. When multiple units are selected, each unit accrues workload.

Invalidate Patient Test Results*		P	Accrue workload when a completed crossmatch test is invalidated. This applies to the workload originally saved with the serologic immediate spin test.
Record a Transfusion Reaction Workup	Crossmatch unit, serological Coombs	P	A crossmatch test for "pre" or "post" is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for "Pre" and "Post" entries. When multiple units are selected, each unit accrues workload.
Enter Crossmatch Results		P	Accrue workload when an individual unit crossmatch is selected to include all phases or only the AHG phase, is completed, and is saved. When multiple units are selected, each unit accrues workload.
Invalidate Patient Test Results*		P	Accrue workload when a completed crossmatch test is invalidated. This applies to the workload originally saved with the test, serological Coombs.
Enter Crossmatch Results	Crossmatch, Repeat Test	M	Accrue workload when an individual unit crossmatch is selected to include all phases or IS or only the AHG phase, is completed, and is saved. When multiple units are selected, each unit accrues workload.
Invalidate Patient Test Results*		M	Accrue workload when an individual unit crossmatch is invalidated.
Enter Daily QC Results	Daily Rack Quality Control (QC)	M	Accrue workload when Daily QC rack completed for one individual rack includes all rows in configured QC. When multiple racks are tested, each completed and saved tab accrues a workload event.
Record Patient Direct Antiglobulin Test	DAT (QC)	M	Accrue workload when Reagent QC completed in Patient DAT testing includes the testing of both the positive and negative control cells, per specificity per lot number, when only one phase of reactivity is chosen for the test grid (IS or AHG). One workload event is collected per completed tab for regular or repeat antiglobulin tests (PS, IgG, Comp).
Record a Transfusion Reaction Workup	Direct Antiglobulin Test (DAT)	P	A DAT test for "pre" or "post" is enabled and a valid interpretation other than Not Tested is selected. A workload event is accrued separately for "Pre" and "Post" entries.
Record Patient Direct Antiglobulin Test		P	Accrue workload when a DAT is completed and saved. This count is used for all antiglobulin tests (PS, IgG, Comp) when ordered from CPRS or Reflex testing.
Invalidate Patient Test Results*		P	Accrue workload when a completed DAT, PS, IgG, or Comp is invalidated.
Record Patient Direct Antiglobulin Test	Direct Antiglobulin Test (DAT) Repeat test	M	Accrue workload when a reflex or repeat DAT test is performed, completed, and saved. This applies to all repeat antiglobulin tests (PS, IgG, Comp).
Invalidate Patient Test Results*		M	Accrue workload when a completed Repeat DAT, PS, IgG, or Comp is invalidated.
Modify Units	Deglycerolize unit	U	Accrue workload when an individual blood unit is processed by the Deglycerolize modification type. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.

Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Deglycerolize."
Discard or Quarantine Unit	Discard unit	U	Accrue workload when an individual blood unit's status is invalidated. When a batch of units is selected, each unit accrues workload.
Remove Final Status*		U	Accrue workload when a unit is discarded for waste or credit. When a batch of units is selected, each unit accrues workload.
Modify Units	Freeze unit	U	Accrue workload when an individual blood unit is processed by the Freeze modification type. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Freeze."
Modify Units	Irradiate unit	U	Accrue workload when an individual blood unit is processed by the Irradiate modification type. When a batch of units is irradiated, each unit accrues workload. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Irradiate."
Modify Units	Leukoreduce unit	U	Accrue workload when an individual blood unit is processed by the Leukoreduce modification type. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Leukoreduce."
Split a Unit	Split unit	U	Accrue workload when a unit modification of Split and a single workload event is recorded regardless of the number of units created by the modification. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	A Split Unit has its unit status invalidated. A single workload event is recorded regardless of the number of units originally created by the modification.
Modify Units	Rejuvenate unit	U	Accrue workload when an individual blood unit is processed by the Rejuvenate modification type. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Rejuvenate."

Modify Units	Thaw	U	Accrue workload when an individual blood unit is processed by the Thaw modification type. When a batch of units is thawed, each unit accrues workload. This applies to Thaw FFP and Thaw Cryo. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Thaw." This modification type is applicable to Thaw FFP and Thaw Cryo.
Modify Units	Wash unit	U	Accrue workload when an individual blood unit is processed by the Wash modification type. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was "Wash."
Modify Units	Volume Reduce	U	Accrue workload when an individual blood unit is processed by the Volume Reduce modification type. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Remove Final Status*		U	An individual blood unit's status is invalidated when the original modification process was Volume Reduce.
Issue Blood Components	Issue unit	P	Accrue workload when a unit is issued to a patient. When a batch of units is processed, each unit invokes one workload process.
Justify Patient ABO/Rh Change	Justification	M	Workload is accrued when a patient's ABO or Rh typing is justified. One workload event is accrued per patient justification.
Login Equipment	Login Equipment	M	Accrue workload when a lot number of any type of equipment is logged into the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
Login Reagent	Login Reagent	M	Accrue workload when a lot number of any type of reagent is logged into the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
Login Supply	Login Supply	M	Accrue workload when a lot number of any type of supply is logged into the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
Maintain Specimen	Maintain Specimen	M	Accrue workload when a specimen is maintained during order acceptance and is required for acceptance of the order. Note: This is collected in addition to the accept order workload accrued by accepting an order. Marking a specimen unacceptable does not create a negative workload event.

Patient antigen phenotype	Patient antigen phenotype (multiple phases)	P	Accrue workload when a patient antigen phenotype test with IS/RT or IS/37 phases is completed and saved. One workload event is collected per completed tab for repeat or regular antigen tests.
Invalidate Patient Test Results*		P	Accrue workload when a patient antigen phenotype test as defined by the antiserum specificity tested with any phases is invalidated.
Patient antigen phenotype	Patient antigen phenotype (single phase)	P	Accrue workload when a patient antigen phenotype test with AHG or IS phase is completed and saved. One workload event is collected per completed tab for repeat or regular antigen tests.
Invalidate Patient Test Results*		P	Accrue workload when a patient antigen phenotype test as defined by the antiserum specificity tested with a single phases is invalidated.
Pool Units	Pool unit	U	Accrue workload when a pooled unit is created and a single workload event is recorded regardless of the number of units included in the pooled unit. This applies to the Pool modification type. Add/Remove unit from a pool does not accrue any workload. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Edit Unit Information*		U	Accrue workload when a unit is inactivated if the pooled unit was created in VBECS.
Remove Final Status		N/A	No effect on workload accrual when a unit is removed from a modified status that was included in a pool.
Discard or Quarantine Unit	Quarantine unit	U	Accrue workload when a unit is marked for quarantine. When a batch of units is selected, each unit accrues workload.
Free Directed Unit For Crossover	Release directed unit	U	Accrue workload when an individual blood unit with the restriction type of "directed" is released for use as an allogeneic unit.
Release Unit from Patient Assignment	Release unit from patient back to inventory	U	Accrue workload when an individual unit is released from patient assignment. When multiple units are selected, each unit accrues workload.
Discard or Quarantine Unit	Release unit from Quarantine	U	Accrue workload when a unit is released from quarantine. When a batch of units is selected, each unit accrues workload.
Return Issued Unit	Return Issued unit	U	Accrue workload when a unit is returned from issue status.
Modify Units	Thaw/pool Cryo	U	Accrue workload when an individual unit has a modification of Thaw/Pool Cryo. A single workload event is recorded regardless of the number of units included in the pooled unit. Note: Workload is not accrued when a patient assignment is processed and VBECS releases all other patient assignments. Workload is not accrued when VBECS is required to quarantine the unit.
Edit Unit Information*		U	Accrue workload when a unit is inactivated (unit record inactivated) when the pooled unit was created in VBECS.

Remove Final Status		N/A	There is no effect on workload accrual when a unit is removed from a modified status that was included in a Thaw/pool Cryo pool.
Enter Post-Transfusion Details	Transfuse Unit	U	Accrue workload when an individual blood unit's status is updated to "transfused."
Remove Final Status*		U	An individual blood unit's status is invalidated when the unit was in a status of "transfused."
Record a Transfusion Reaction Workup	Transfusion Reaction Investigation	P	Accrue workload when a transfusion reaction investigation is saved. This does not include workload accrued by the optional TRW serologic testing.
Invalidate Patient Test Results*		P	Accrue workload when a transfusion reaction investigation previously saved is invalidated.
Unit Antigen Typing	Unit Antigen phenotyping, Multiple Test phases	U	Accrue workload when a unit antigen phenotype test with IS/RT or IS/37 phases is selected and completed for an individual blood unit. There is no special handling for workload collection for additional repeat antigen typing tests on a unit.
Edit Unit Information*		U	Accrue workload when a unit antigen phenotype test with Multiple Test phases is invalidated for an individual blood unit.
Unit Antigen Typing	Unit Antigen phenotyping, Single Test phase	U	A unit antigen phenotype test with AHG or IS phase is selected and completed for an individual blood unit. There is no special handling for workload collection for additional repeat antigen typing tests on a unit.
Edit Unit Information*		U	Accrue workload when a unit antigen phenotype test with Single Test phase is invalidated for an individual blood unit.
Incoming Shipment	Unit login	U	An individual unit record is activated as "saved" to an incoming shipment invoice. When multiple units are entered, each unit added to the database accrues workload.
Edit Unit Information*		U	Accrue workload when a unit is inactivated and logged in through incoming shipment or is a pooled unit created in VBECS. When the unit was created by split modification, no workload is invalidated in this option.
Outgoing Shipment	Unit logout	U	An individual unit's status is updated to "transferred" on a confirmed outgoing shipment invoice. When multiple units are selected, each unit accrues workload. Accrue workload on confirmation of the invoice, not the addition of a unit to a temporary outgoing shipment invoice: an invoice may be confirmed only once.
Remove Final Status*		U	An individual unit status is invalidated when the unit had a previous unit status of "transferred."
Update Equipment Record	Update Equipment Record	M	Accrue workload when a lot number of any type of equipment is updated in the system.
Update Reagent Inventory	Update Reagent Inventory	M	Accrue workload when a lot number of any type of reagent is updated in the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.

Update Supply Inventory	Update Supply Inventory	M	Accrue workload when a lot number of any type of supply is updated in the system. When multiple lot numbers are processed in a batch, each lot number's workload is counted.
-------------------------	-------------------------	---	--

*Accumulates negative workload when it is associated with inactivation of a unit or removal of a final status.

Index

A

Accept Orders: Accept an Order.....	225
Accessing the System	35
Administrative Data Report	313
Appendices	4
Application Architecture.....	6
Assumptions	147
Audit Trail	315
Automated Entry of Post-Transfusion Details from BCE COTS Application	307

B

Background and Basic Knowledge	3
Barcode Scanner Configuration and Troubleshooting	481
Blood Availability	317

C

C:T Ratio Report.....	321
check character	116
Collecting Workload.....	55
Common Error Corrections.....	477
Commonly Used System Rules	12
Component Processing	115
Configure Division	65
Configuring Site Parameters	57
Contents Tab.....	28
Cost Accounting Report.....	319
Creating and Viewing Reports.....	16
Crossmatch Interpretations	273
Customer Support	10

D

Database Table Information.....	377
Date and Time	12
Discard or Quarantine Individual Units	177
Discard or Quarantine Units by Invoice.....	180
Display Order Alerts	232
Display Patient Merge Alerts.....	234
Display Patient Update Alerts.....	233
Division Transfusion Report	323
Division Workload Report	325
Document ABO Incompatible Transfusions	308
Downtime Forms and Instructions	367
Drop-Down Menus	47

E

Enter Daily QC Results.....	92
Equipment.....	109

Equipment: Log In Equipment.....	110
Equipment: Maintain Equipment.....	111
Exception Report	327

F

Figures

Barcode Scanner Configuration.....	481
Blood Transfusion Record Form.....	215
Caution Tag	214
Order Status Flowchart.....	407
System Schematic.....	9
Unit Status Flowchart	409
VBECS Work Flow.....	427
Figures and Tables	3
Finalize Transfusion Reaction Workup Report.....	291
Finalize/Print TRW.....	291
Forms	
Patient Testing Form	369
Unit ABO/Rh Confirmation Form.....	372
Unit Issue and Inspection Log	373
Unit Modification Form.....	374
Freeware Disclaimer	1
Frequently Asked Questions.....	475

G

Glossary	357
Glossary Tab.....	31

H

Hardware and Infrastructure Architecture.....	6
How This User Guide Is Organized	3

I

Icons and Buttons	48
Inactivate a Unit.....	171
Inappropriate Transfusion Request Report	329
Index Tab.....	29
Introduction	1
Issue Blood Components	192
Issue Blood Components (Emergency).....	196
Issue Blood Components (Routine)	193
Issued/Returned Unit Report.....	331

J

Justify ABO/Rh Change	310
-----------------------------	-----

K

Keyboard Shortcuts for Windows.....	46
Known Defects and Anomalies.....	445

L

Limitations and Restrictions	431
Local Facilities.....	75
Local Machine Screen Resolution	39
Locking and Time-Outs	12
Log in Reagents	99
Log in Supplies.....	105
Login Message.....	73

M

Main Status Bar	55
Main Toolbar	48
Maintain Minimum Levels	103
Maintain Patient Records	295
Maintain Specimen	235
Maintain Unit Records.....	167
Medication Profile	50
Modify Component.....	147
MSBOS.....	88

O

Online Help F1	42
Options	3
Order Alerts	71
Order History Report	333
Order Reflex Tests	245
Orientation	3
Other VBECS Functions.....	55
Outgoing Shipment.....	123

P

Patient Blood Availability.....	49
Patient History Report	335
Patient Information Toolbar.....	49
Patient Merge.....	251
Patient Search Screen.....	249
Patient Services.....	249
Patient Testing Form.....	369
Patient Testing Worklist and Testing Worklist Reports.....	339
Patient Testing: Cancel an Active Order.....	239
Patient Testing: Enter Antibody Identification Results.....	280
Patient Testing: General Instructions	256
Patient Testing: Pending Task List	241
Patient Testing: Record a Crossmatch	269
Patient Testing: Record a Direct Antiglobulin Test	266
Patient Testing: Record a Patient ABO/Rh	261
Patient Testing: Record a Patient Antibody Screen	264
Patient Testing: Record a Patient Antigen Typing	274
Patient Testing: Record a Transfusion Reaction Workup	283
Post-Transfusion Information	303
Print Backup Caution Tags and Blood Transfusion Record Forms	209
Problems?	10

Processing Orders	225
Product Modifications.....	69
Prolonged Transfusion Time Report.....	337

Q

Quick K Calculator	116
--------------------------	-----

R

Reagents	99
Reagents and Supplies	92
Recommended Report Usage.....	479
References	353
Related Manuals and Materials.....	2
Remote Desktop Connections Screen Size	40
Restore VistA Connection	38

S

Screen Display.....	42
Screen Settings	39
Search Tab	30
Searching the Database.....	13
Security.....	4
Setting Transfusion Parameters	85
Shipments	115
Software Basics	10
Sound.....	43
Special Instructions & Transfusion Requirements: Enter and Remove Special Instructions	300
Supplemental Information	4
Supplies	105
System Responses to Active Transfusion Requirements in Select Unit and Issue Blood Component.....	485
System Validation	
Antibodies	461
Blood Products	455
Canned Comments.....	460
Component Classes	457
Configure Daily QC: Full Service	452
Configure Division	449
Configure Testing: Full Service.....	449
Edit Invoice Text	460
Installation Qualification	447
Introduction	447
Local Facilities	456
Log In Reagents: Full Service	454
Maintain Minimum Levels: Full Service.....	453
MSBOS	460
Operational Qualification	447
Order Alerts	451
Performance Qualification.....	448
Prerequisites	448
Product Modifications: Full Service	450
Record System Settings.....	449
Testing.....	462
Transfusion Complications.....	458

Transfusion Effectiveness.....	459
Update User Roles.....	462
Validation Plan.....	447
VBECS Forms and Reports.....	463
Workload Codes.....	460
System Validation.....	447

T

Tables

Alerts That May Occur in Issue Blood Components (Emergency).....	201
Alerts That May Occur in Issue Blood Components (Routine).....	199
Alerts That May Occur in Modify Units.....	156
Alerts That May Occur in Patient Testing	
General Instructions	259
Alerts That May Occur in Print Unit Caution Tag & Transfusion Record Form.....	213
Alerts That May Occur in Select Units.....	188
Allowable Product Modifications by Original Product Type	399
Antibody and Antigen Table, Database Conversion.....	380
Antibody and Antigen Table, Irregular Antibodies	377
Antibody Screen Test Interpretation.....	381
BCE COTS Message Updates	412
Blood Products Worksheet	455
Canned Comment Category Types and Text	382
Component Class Requirements Worksheet.....	457
CPRS Orderable Blood Components (Component Classes) Mapped to ICCBBA Component Classes.....	390
Crossmatch Interpretations for Print Unit Caution Tag & Transfusion Record Form	273
Details in Audit Trail Report	391
Details in Exception Report.....	396
Enable Crossmatch Option for OTHER Product Types.....	398
Local Facilities Worksheet	456
Modification Types and Costs	450
National Treating Specialty Table	404
Rules for Electronic and Serologic Crossmatch	408
Sample Blood Products Worksheet	455
Sample Local Facilities Worksheet	456
Transfusion Complications Worksheet.....	458
Transfusion Effectiveness Worksheet	459
Transfusion Reaction Types	294
User Roles and Functions	4
Valid Interpretations.....	137
Valid Observed Test Results Entries	14
VBECS Maximum Surgical Blood Order Schedule (MSBOS)	411
Terms	3
Test Units.....	136
Testing	14
Throughout VBECS.....	55
Tool Tips	47
Toolbars	48
Transfusion Complications Report	342
Transfusion Effectiveness	86
Transfusion Effectiveness Report	344
Transfusion Reaction Count Report.....	346
Transfusion Requirements Report	348

U

Unit ABO/Rh Confirmation Form	372
Unit History Report	350
Unit Issue and Inspection Log	373
Unit Modification Form.....	374
Unit Search Screen	134
United Blood Services (UBS) Facility Barcodes	483
Update Reagent Inventory	100
Update Supply Inventory	107
Updated Patients and Deceased Patients.....	255
Using Online Help	26
Using the Software	10
Using VBECS.....	25

V

VA Service Desk.....	10
VBECS User Report	58
VBECS Work Flow	427
View/Print a Finalized Transfusion Reaction Workup Report.....	293
View/Print Reagent Inventory	102
View/Print Supply Inventory	108
Viewing the VBECS Version Number	32
VistA Records in VBECS.....	34
Volume	45

W

Work Process Flowcharts	465
Configure Division Settings: Full-Service Divisions.....	465
Configure Division Settings: Transfusion-Only Divisions	466
Issue, Return, or Transfuse Units: Full-Service Divisions	471
Issue, Return, or Transfuse: Transfusion-Only Divisions.....	472
Maintain Reagents and Daily QC: Full-Service Divisions	467
Modify Units	470
Patient Testing	474
Process a Patient Order	473
Process Incoming Shipment: Full-Service Divisions.....	468
Process Incoming Shipment: Transfusion-Only Divisions	469
Working with Data	12
Workload Codes	90
Workload Process Mapping to Application Option Table	495
Workload Process to Use Case Mapping Table	416

This is the last page of *VistA Blood Establishment Computer Software (VBECS) 2.1.0 User Guide*.